
Allen SDK Documentation

Release dev

Allen Institute for Brain Science

May 20, 2021

Contents

1	Install Guide	1
1.1	Quick Start Using Anaconda	1
1.2	Quick Start Using Pip	1
1.3	Other Distribution Formats	2
1.4	Installation with Docker (Optional)	2
2	Data Resources	3
2.1	Brain Observatory	3
2.2	Cell Types	6
2.3	Mouse Connectivity	9
2.4	Reference Space	11
2.5	API Access	11
2.6	Visual Coding – Neuropixels	14
2.7	Visual Behavior - Optical Physiology	22
3	Models	45
3.1	Generalized LIF Models	45
3.2	Biophysical Models	51
4	Examples	61
5	Authors	63
6	allensdk package	65
6.1	Subpackages	65
6.2	Submodules	371
6.3	Module contents	371
7	Allen Brain Observatory	373
8	Allen Cell Types Database	375
9	Allen Mouse Brain Connectivity Atlas	377
10	What’s New - 2.11.0	379
11	What’s New - 2.11.0	381

12 What's New - 2.10.3	383
13 What's New - 2.10.2	385
14 What's New - 2.10.1	387
15 What's New - 2.9.0	389
16 What's New - 2.8.0	391
17 What's New - 2.7.0	393
18 What's New - 2.6.0	395
19 What's New - 2.5.0 (January 29, 2021)	397
20 What's New - 2.4.0 (December 21, 2020)	399
21 What's New - 2.3.2 (October 19, 2020)	401
22 What's New - 2.3.1 (October 13, 2020)	403
23 What's New - 2.3.0 (October 9, 2020)	405
24 What's New - 2.2.0 (September 3, 2020)	407
25 What's New - 2.1.0 (July 16, 2020)	409
26 What's New - 2.0.0 (June 11, 2020)	411
27 Previous Release Notes	413
Bibliography	415
Python Module Index	417
Index	423

CHAPTER 1

Install Guide

This guide is a resource for using the Allen SDK package. It is maintained by the [Allen Institute for Brain Science](#).

Attention: As of October 2019, we have dropped Python 2 support. The Allen SDK is developed and tested with Python 3.6 and 3.7. We do not guarantee consistent behavior with other Python versions.

1.1 Quick Start Using Anaconda

1. From the [Anaconda downloads page](#), download the Python 3.7 version for your operating system and run the installer.
2. After the installation is complete, open up a terminal (in Windows open Anaconda3 Command Prompt).
3. Install the AllenSDK using PIP:

```
pip install allensdk
```

4. Download one of our many [Jupyter Notebook examples](#) to a new folder.
5. In your terminal, navigate to the directory where you downloaded the Jupyter Notebook example and run the following command:

```
jupyter notebook
```

6. Your browser should open and you should see the Jupyter Notebook example. Enjoy using the Allen SDK!

1.2 Quick Start Using Pip

First ensure you have [pip](#) installed. It is included with the Anaconda distribution.

```
pip install allensdk
```

To uninstall the SDK:

```
pip uninstall allensdk
```

1.3 Other Distribution Formats

The Allen SDK is also available from the Github source repository.

1.4 Installation with Docker (Optional)

[Docker](#) is an open-source technology for building and deploying applications with a consistent environment including required dependencies. The AllenSDK is not distributed as a Docker image, but example Dockerfiles are available.

1. Ensure you have Docker installed.
2. Use Docker to build the image:

```
docker pull alleninstitute/allensdk
```

Other docker configurations are also available under docker directory in the source repository.

3. Run the docker image:

```
docker run -i -t -p 8888:8888 alleninstitute/allensdk /bin/bash
```

4. Run the SDK tests:

```
cd allensdk
make test
```

5. Start a Jupyter Notebook:

```
cd allensdk/doc_template/examples_root/examples/nb
jupyter notebook --ip=* --no-browser --allow-root
```


Using the browser on your host machine, navigate to the path provided by the output from the jupyter notebook command.

The Allen SDK features Python code to support data and model access for the Allen Cell Types Database. Resources for other Allen Brain Atlas data resources will come in future updates.

2.1 Brain Observatory

The [Allen Brain Observatory](#) is a database of the visually-evoked functional responses of neurons in mouse visual cortex based on 2-photon fluorescence imaging. Characterized responses include orientation tuning, spatial and temporal frequency tuning, temporal dynamics, and spatial receptive field structure.

The data is organized into experiments and experiment containers. An experiment container represents a group of experiments with the same targeted imaging area, imaging depth, and Cre line. The individual experiments within an experiment container have different stimulus protocols, but cover the same imaging field of view.

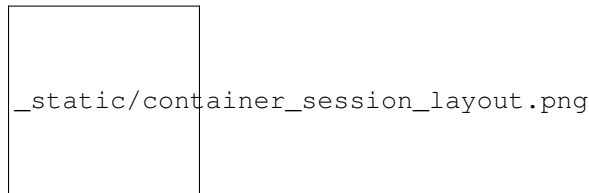

Note: Version 1.3 of scipy fixed an error in its 2 sample Kolmogorov-Smirnoff test implementation. The new version produces more accurate p values for small and medium-sized samples. This change impacts speed tuning analysis p values (as returned by *StimulusAnalysis.get_speed_tuning*). If you access precalculated analysis results via *BrainObservatoryCache.get_ophys_experiment_analysis*, you will see values calculated using an older version of scipy's *ks_2samp*. To access values calculated from the new version, install *scipy* $\geq 1.3.0$ in your environment and construct a *StimulusAnalysis* object from a *BrainObservatoryNwbDataSet* (as returned by *BrainObservatoryCache.get_ophys_experiment_data*).

Note: Data collected after September 2016 uses a new session C stimulus designed to better-characterize spatial receptive fields in higher visual areas. The original locally sparse noise stimulus used 4.65 visual degree pixels. Session C2 broke that stimulus into two separate stimulus blocks: one with 4.65 degree pixels and one with 9.3 degree pixels.

Note that the `stimulus_info` module refers to these as `locally_sparse_noise_4deg` and `locally_sparse_noise_8deg`, respectively.

For more information on experimental design and a data overview, please visit the [Allen Brain Observatory data portal](#).

2.1.1 Data Processing

For all data in Allen Brain Observatory, we perform the following processing:

1. Segment cell masks from each experiment's 2-photon fluorescence video
2. Associate cells from experiments belonging to the same experiment container and assign unique IDs
3. Extract each cell's mean fluorescence trace
4. Extract mean fluorescence traces from each cell's surrounding neuropil
5. Demix traces from overlapping ROIs
6. Estimate neuropil-corrected fluorescence traces
7. Compute dF/F
8. Compute stimulus-specific tuning metrics

All traces and masks for segmented cells in an experiment are stored in a Neurodata Without Borders (NWB) file. Stored traces include the raw fluorescence trace, neuropil trace, demixed trace, and dF/F trace. Code for extracting neuropil-corrected fluorescence traces, computing dF/F, and computing tuning metrics is available in the SDK.

New in June 2017: Trace demixing is a new addition as of June 2017. All past data was reprocessed using the new demixing algorithm. We have also developed a new module to better characterize a cell's receptive field. Take a look at the [receptive field analysis example notebook](#)

For more information about data processing, please [read the technical whitepapers](#).

2.1.2 Getting Started

The Brain Observatory [Jupyter notebook](#) has many code samples to help get started with the available data:

- [Download experimental metadata by visual area, imaging depth, and Cre line](#)
- [Find cells with specific response properties, like direction tuning](#)
- [Download data for an experiment](#)
- [Plot raw fluorescence traces, neuropil-corrected traces, and dF/F](#)
- [Find the ROI mask for a given cell](#)
- [Run neuropil correction](#)
- [Get pupil location and size](#)

The code used to analyze and visualize data in the [Allen Brain Observatory data portal](#) is available as part of the SDK. Take a look at this [Jupyter notebook](#) to find out how to:

- [Plot cell's response to its preferred stimulus condition](#)
- [Compute a cell's on/off receptive field based on the locally sparse noise stimulus](#)

More detailed documentation is available demonstrating how to:

- [Read and visualize the stimulus presentation tables in the NWB files](#)
- [Understand the layout of Brain Observatory NWB files](#)

- Map previous cell specimen IDs to current cell specimen IDs

2.1.3 Precomputed Cell Metrics

A large table of precomputed metrics are available for download to support population analysis and filtering. The table below describes all of the metrics in the table. The `get_cell_specimens()` method will download this table as a list of dictionaries which can be converted to a pandas DataFrame as shown in this [Jupyter notebook](#).

Stimulus	Metric	Field Name
drifting gratings	orientation selectivity	osi_dg
	direction selectivity	dsi_dg
	preferred direction	pref_dir_dg
	preferred temporal frequency	pref_tf_dg
	response p value	p_dg
	global ori. selectivity	g_osi_dg
	global dir. selectivity	g_dsi_dg
	response reliability	reliability_dg
	running modulation	run_mod_dg
	running modulation p value	p_run_mod_dg
	pref. condition mean df/f	peak_dff_dg
	TF discrimination index	tfdi_dg
static gratings	orientation selectivity	osi_sg
	preferred orientation	pref_ori_sg
	preferred spatial frequency	pref_sf_sg
	preferred phase	pref_phase_sg
	mean time to peak response	time_to_peak_sg
	response p value	p_sg
	global ori. selectivity	g_osi_sg
	reponse reliability	reliability_sg
	running modulation	run_mod_sg
	running modulation p value	p_run_mod_sg
	pref. condition mean df/f	peak_dff_ns
	SF discrimination index	sfdi_sg
natural scenes	mean time to peak response	time_to_peak_ns
	preferred scene index	pref_scene_ns
	response p value	p_ns
	image selectivity	image_sel_ns
	running modulation	run_mod_ns
	running modulation p value	p_run_mod_ns
	pref. condition mean df/f	peak_dff_ns
natural movie 1	response reliability (session A)	reliability_nm1_a
	response reliability (session B)	reliability_nm1_b
	response reliability (session C)	reliability_nm1_c
natural movie 2	response reliability	reliability_nm2
natural movie 3	response reliability	reliability_nm3
locally sparse noise	RF area (on subunit)	rf_area_on_lsn
	RF area (off subunit)	rf_area_off_lsn
	RF center (on subunit)	rf_center_on_x, rf_center_on_y
	RF center (off subunit)	rf_center_off_x, rf_center_off_y
	RF chi^2	rf_chi2_lsn
	RF on-off subunit distance	rf_distance_lsn
	RF on-off subunit overlap index	rf_overlap_lsn

2.2 Cell Types

The Allen Cell Types data set is a database of mouse and human neuronal cell types based on multimodal characterization of single cells to enable data-driven approaches to classification and is fully integrated with other Allen Brain Atlas resources. The database currently includes:

- **electrophysiology**: whole cell current clamp recordings made from Cre-positive neurons
- **morphology**: 3D bright-field images of the complete structure of neurons from the visual cortex

This page describes how the SDK can be used to access data in the Cell Types Database. For more information, please visit the Cell Types Database [home page](#) and the [API documentation](#).

2.2.1 Examples

The Cell Types [Jupyter notebook](#) has many code samples to help get started with analysis:

- Download and plot stimuli and responses from an NWB file for a cell
- Download and plot a cell's morphological reconstruction
- Download and plot precomputed electrophysiology features
- Download precomputed morphology features to a table
- Compute electrophysiology features for a single sweep

2.2.2 Cell Types Cache

The `CellTypesCache` class provides a Python interface for downloading data in the Allen Cell Types Database into well known locations so that you don't have to think about file names and directories. The following example demonstrates how to download meta data for all cells with 3D reconstructions, then download the reconstruction and electrophysiology recordings for one of those cells:

```
from allensdk.core.cell_types_cache import CellTypesCache

ctc = CellTypesCache(manifest_file='cell_types/manifest.json')

# a list of cell metadata for cells with reconstructions, download if necessary
cells = ctc.get_cells(require_reconstruction=True)

# open the electrophysiology data of one cell, download if necessary
data_set = ctc.get_ephys_data(cells[0]['id'])

# read the reconstruction, download if necessary
reconstruction = ctc.get_reconstruction(cells[0]['id'])
```

`CellTypesCache` takes care of knowing if you've already downloaded some files and reads them from disk instead of downloading them again. All data is stored in the same directory as the `manifest_file` argument to the constructor.

2.2.3 Feature Extraction

The `EphysFeatureExtractor` class calculates electrophysiology features from cell recordings. `extract_cell_features()` can be used to extract the precise feature values available in the Cell Types Database:

```

from allensdk.core.cell_types_cache import CellTypesCache
from allensdk.ephys.extract_cell_features import extract_cell_features
from collections import defaultdict

# initialize the cache
ctc = CellTypesCache(manifest_file='cell_types/manifest.json')

# pick a cell to analyze
specimen_id = 324257146

# download the ephys data and sweep metadata
data_set = ctc.get_ephys_data(specimen_id)
sweeps = ctc.get_ephys_sweeps(specimen_id)

# group the sweeps by stimulus
sweep_numbers = defaultdict(list)
for sweep in sweeps:
    sweep_numbers[sweep['stimulus_name']].append(sweep['sweep_number'])

# calculate features
cell_features = extract_cell_features(data_set,
                                     sweep_numbers['Ramp'],
                                     sweep_numbers['Short Square'],
                                     sweep_numbers['Long Square'])

```

2.2.4 File Formats

This section provides a short description of the file formats used for Allen Cell Types data.

Morphology SWC Files

Morphological neuron reconstructions are available for download as SWC files. The SWC file format is a white-space delimited text file with a standard set of headers. The file lists a set of 3D neuronal compartments, each of which has:

Column	Data Type	Description
id	string	compartment ID
type	integer	compartment type
x	float	3D compartment position (x)
y	float	3D compartment position (y)
z	float	3D compartment position (z)
radius	float	compartment radius
parent	string	parent compartment ID

Comment lines begin with a '#'. Reconstructions in the Allen Cell Types Database can contain the following compartment types:

Type	Description
0	unknown
1	soma
2	axon
3	basal dendrite
4	apical dendrite

The Allen SDK comes with a `swc` Python module that provides helper functions and classes for manipulating SWC files. Consider the following example:

```
import allensdk.core.swc as swc

# if you ran the examples above, you will have a reconstruction here
file_name = 'cell_types/specimen_485909730/reconstruction.swc'
morphology = swc.read_swc(file_name)

# subsample the morphology 3x. root, soma, junctions, and the first child of the root
# are preserved.
sparse_morphology = morphology.sparsify(3)

# compartments in the order that they were specified in the file
compartment_list = sparse_morphology.compartment_list

# a dictionary of compartments indexed by compartment id
compartments_by_id = sparse_morphology.compartment_index

# the root soma compartment
soma = morphology.soma

# all compartments are dictionaries of compartment properties
# compartments also keep track of ids of their children
for child in morphology.children_of(soma):
    print(child['x'], child['y'], child['z'], child['radius'])
```

Neurodata Without Borders

The electrophysiology data collected in the Allen Cell Types Database is stored in the [Neurodata Without Borders](#) (NWB) file format. This format, created as part of the [NWB initiative](#), is designed to store a variety of neurophysiology data, including data from intra- and extracellular electrophysiology experiments, optophysiology experiments, as well as tracking and stimulus data. It has a defined schema and metadata labeling system designed so software tools can easily access contained data.

The Allen SDK provides a basic Python class for extracting data from Allen Cell Types Database NWB files. These files store data from intracellular patch-clamp recordings. A stimulus current is presented to the cell and the cell's voltage response is recorded. The file stores both stimulus and response for several experimental trials, here called "sweeps." The following code snippet demonstrates how to extract a sweep's stimulus, response, sampling rate, and estimated spike times:

```
from allensdk.core.nwb_data_set import NwbDataSet

# if you ran the examples above, you will have a NWB file here
file_name = 'cell_types/specimen_485909730/ephys.nwb'
data_set = NwbDataSet(file_name)

sweep_numbers = data_set.get_sweep_numbers()
sweep_number = sweep_numbers[0]
sweep_data = data_set.get_sweep(sweep_number)

# spike times are in seconds relative to the start of the sweep
spike_times = data_set.get_spike_times(sweep_number)

# stimulus is a numpy array in amps
stimulus = sweep_data['stimulus']
```

(continues on next page)

(continued from previous page)

```
# response is a numpy array in volts
reponse = sweep_data['response']

# sampling rate is in Hz
sampling_rate = sweep_data['sampling_rate']

# start/stop indices that exclude the experimental test pulse (if applicable)
index_range = sweep_data['index_range']
```

HDF5 Overview

NWB is implemented in [HDF5](#). HDF5 files provide a hierarchical data storage that mirrors the organization of a file system. Just as a file system has directories and files, and HDF5 file has groups and datasets. The best way to understand an HDF5 (and NWB) file is to open a data file in an HDF5 browser. [HDFView](#) is the recommended browser from the makers of HDF5.

There are HDF5 manipulation libraries for many languages and platforms. MATLAB and Python in particular have strong HDF5 support.

2.3 Mouse Connectivity

The Allen Mouse Brain Connectivity Atlas consists of high-resolution images of axonal projections targeting different anatomic regions or various cell types using Cre-dependent specimens. Each data set is processed through an informatics data analysis pipeline to obtain spatially mapped quantified projection information.

This page describes how to use the SDK to access experimental projection data and metadata. For more information, please visit the Connectivity Atlas [home page](#) and the [API documentation](#)

2.3.1 Structure-Level Projection Data

All AAV projection signal in the Allen Mouse Connectivity Atlas has been registered to the expert-annotated Common Coordinate Framework (CCF) and summarized to structures in the adult mouse structure ontology. Most commonly used for analysis are measures of the density of projection signal in all brain areas for every experiment. This data is available for download and is described in more detail on the [structure unionizes page](#).

2.3.2 Voxel-Level Projection Data

The CCF-registered AAV projection signal is also available for download as a set of 3D volumes for each experiment. The following data volumes are available for download:

- **projection density:** sum of detected projection pixels / sum of all pixels in voxel
- **injection_fraction:** fraction of pixels belonging to manually annotated injection site
- **injection_density:** density of detected projection pixels within the manually annotated injection site
- **data_mask:** binary mask indicating if a voxel contains valid data. Only valid voxels should be used for analysis.

2.3.3 Code Examples

The Mouse Connectivity Jupyter notebook has many code samples to help get started with analysis:

- Download experimental metadata by injection structure and transgenic line
- Download projection signal statistics at a structure level
- Build a structure-to-structure matrix of projection signal values
- Download and visualize gridded projection signal volumes

2.3.4 Mouse Connectivity Cache

The `MouseConnectivityCache` class saves all of the data you can download via the `MouseConnectivityApi` in well known locations so that you don't have to think about file names and directories. It also takes care of knowing if you've already downloaded some files and reads them from disk instead of downloading them again. The following example demonstrates how to download meta data for all experiments with injections in the isocortex and download the projection density volume for one of them:

```
from allensdk.core.mouse_connectivity_cache import MouseConnectivityCache

# tell the cache class what resolution (in microns) of data you want to download
mcc = MouseConnectivityCache(resolution=25)

# use the structure tree class to get information about the isocortex structure
structure_tree = mcc.get_structure_tree()
isocortex_id = structure_tree.get_structures_by_name(['Isocortex'])[0]['id']

# a list of dictionaries containing metadata for non-Cre experiments
experiments = mcc.get_experiments(file_name='non_cre.json',
                                  injection_structure_ids=[isocortex_id])

# download the projection density volume for one of the experiments
pd = mcc.get_projection_density(experiments[0]['id'])
```

2.3.5 File Formats

This section provides a short description of the file formats used for data in the Allen Mouse Connectivity Atlas.

NRRD Files

All of the volumetric data in the connectivity atlas are stored as **NRRD (Nearly Raw Raster Data)** files. A NRRD file consists of a short ASCII header followed by a binary array of data values.

To read these in Python, we recommend the `pynrrd` package. Usage is straightforward:

```
import nrrd

file_name = 'mouse_connectivity/experiment_644250774/projection_density_25.nrrd'
data_array, metadata = nrrd.read(file_name)
```

2.4 Reference Space

Allen Institute atlases and data are registered, when possible, to one of several common reference spaces. Working in such a space allows you to easily compare data across subjects and experimental modalities.

This page documents how to use the Allen SDK to interact with a reference space. For more information and a list of reference spaces, see the [atlas drawings and ontologies API documentation](#) and the [3D reference models API documentation](#). For details about the construction of the Common Coordinate Framework space, see the [CCFv3 whitepaper](#).

2.4.1 Structure Tree

Brain structures in our reference spaces are arranged in trees. The leaf nodes of the tree describe the very fine anatomical divisions of the space, while nodes closer to the root correspond to gross divisions. The `StructureTree` class provides an interface for interacting with a structure tree.

To download a structure tree, use the `allensdk.core.reference_space_cache.ReferenceSpaceCache` class as seen in [this example](#)

2.4.2 Annotation Volumes

An annotation volume is a 3d raster image that segments the reference space into structures. Each voxel in the annotation volume is assigned an integer value that describes the finest structure to which that point in space definitely belongs.

To download a nrrd formatted annotation volume at a specified isometric resolution, use the `allensdk.core.reference_space_cache.ReferenceSpaceCache` class. There is [an example](#) in the notebook.

2.4.3 ReferenceSpaceCache Class

The `allensdk.core.reference_space_cache.ReferenceSpaceCache` class provides a Python interface for downloading structure trees and annotation volumes. It takes care of knowing if you've already downloaded the files and reads them from disk instead of downloading them again.

The class contains methods for working with our reference spaces. Some use cases might include:

- Building an indicator mask for one or more structures
- Viewing the annotation
- Querying the structure graph

Please see the [example notebook](#) for more code samples.

2.5 API Access

The `allensdk.api` package is designed to help retrieve data from the [Allen Brain Atlas API](#). `api` contains methods to help formulate API queries and parse the returned results. There are several pre-made subclasses available that provide pre-made queries specific to certain data sets. Currently there are several subclasses in Allen SDK:

- `CellTypesApi`: data related to the Allen Cell Types Database
- `BiophysicalApi`: data related to biophysical models
- `GlifApi`: data related to GLIF models

- *AnnotatedSectionDataSetsApi*: search for experiments by intensity, density, pattern, and age
- *GridDataApi*: used to download 3-D expression grid data
- *ImageDownloadApi*: download whole or partial two-dimensional images
- *MouseConnectivityApi*: common operations for accessing the Allen Mouse Brain Connectivity Atlas
- *OntologiesApi*: data about neuroanatomical regions of interest
- *ConnectedServices*: schema of Allen Institute Informatics Pipeline services available through the RmaApi
- *RmaApi*: general-purpose HTTP interface to the Allen Institute API data model and services
- *SvgApi*: annotations associated with images as scalable vector graphics (SVG)
- *SynchronizationApi*: data about image alignment
- *TreeSearchApi*: list ancestors or descendants of structure and specimen trees

2.5.1 RMA Database and Service API

One API subclass is the *RmaApi* class. It is intended to simplify constructing an RMA query.

The *RmaApi* is a base class for much of the `allensdk.api.queries` package, but it may be used directly to customize queries or to build queries from scratch.

Often a query will simply request a table of data of one type:

```
from allensdk.api.queries.rma_api import RmaApi

rma = RmaApi()

data = rma.model_query('Atlas',
                       criteria="[name$il'*Mouse*']")
```

This will construct the RMA query url, make the query and parse the resulting JSON into an array of Python dicts with the names, ids and other information about the atlases that can be accessed via the API.

Using the criteria, include and other parameter, specific data can be requested.

```
associations = ''.join(['[id$eq1]',
                       'structure_graph(ontology)', ', ',
                       'graphic_group_labels'])

atlas_data = rma.model_query('Atlas',
                             include=associations,
                             criteria=associations,
                             only=['atlases.id',
                                   'atlases.name',
                                   'atlases.image_type',
                                   'ontologies.id',
                                   'ontologies.name',
                                   'structure_graphs.id',
                                   'structure_graphs.name',
                                   'graphic_group_labels.id',
                                   'graphic_group_labels.name'])
```

Note that a ‘class’ name is used for the first parameter. ‘Association’ names are used to construct the include and criteria parameters nested using parentheses and commas. In the only clause, the ‘table’ form is used, which is

generally a plural lower-case version of the class name. The only clause selects specific ‘fields’ to be returned. The schema that includes the classes, fields, associations and tables can be accessed in JSON form using:

```
# http://api.brain-map.org/api/v2/data.json
schema = rma.get_schema()
for entry in schema:
    data_description = entry['DataDescription']
    clz = list(data_description.keys())[0]
    info = list(data_description.values())[0]
    fields = info['fields']
    associations = info['associations']
    table = info['table']
    print("class: %s" % (clz))
    print("fields: %s" % (','.join(f['name'] for f in fields)))
    print("associations: %s" % (','.join(a['name'] for a in associations)))
    print("table: %s\n" % (table))
```

2.5.2 Using Pandas to Process Query Results

When it is difficult to get data in exactly the required form using only an RMA query, it may be helpful to perform additional operations on the client side. The pandas library can be useful for this.

Data from the API can be read directly into a pandas [Dataframe](#) object.

```
import pandas as pd

structures = pd.DataFrame(
    rma.model_query('Structure',
                    criteria='[graph_id$eq1]',
                    num_rows='all'))
```

Indexing subsets of the data (certain columns, certain rows) is one use of pandas: specifically `.loc`:

```
names_and_acronyms = structures.loc[:, ['name', 'acronym']]
```

and Boolean indexing

```
mea = structures[structures.acronym == 'MEA']
mea_id = mea.iloc[0,:].id
mea_children = structures[structures.parent_structure_id == mea_id]
print(mea_children['name'])
```

`Concatenate`, `merge` and `join` are used to add columns or rows:

When an RMA call contains an include clause, the associated data will be represented as a python dict in a single column. The column may be converted to a proper Dataframe and optionally dropped.

```
criteria_string = "structure_sets[name$eq'Mouse Connectivity - Summary']"
include_string = "ontology"
summary_structures = \
    pd.DataFrame(
        rma.model_query('Structure',
                        criteria=criteria_string,
                        include=include_string,
                        num_rows='all'))

ontologies = \
    pd.DataFrame(
```

(continues on next page)

(continued from previous page)

```
list(summary_structures.ontology)).drop_duplicates()
flat_structures_dataframe = summary_structures.drop(['ontology'], axis=1)
```

Alternatively, it can be accessed using normal python dict and list operations.

```
print(summary_structures.ontology[0]['name'])
```

Pandas Dataframes can be written to a CSV file using `to_csv` and read using `load_csv`.

```
summary_structures[['id',
                    'parent_structure_id',
                    'acronym']].to_csv('summary_structures.csv',
                                      index_label='structure_id')
reread = pd.read_csv('summary_structures.csv')
```

Iteration over a Dataframe of API data can be done in several ways. The `.itertuples` method is one way to do it.

```
for id, name, parent_structure_id in summary_structures[['name',
                                                         'parent_structure_id']].
    .itertuples():
    print("%d %s %d" % (id, name, parent_structure_id))
```

2.5.3 Caching Queries on Disk

`wrap()` has several parameters for querying the API, saving the results as CSV or JSON and reading the results as a pandas dataframe.

```
from allensdk.api.warehouse_cache.cache import Cache

cache_writer = Cache()
do_cache=True
structures_from_api = \
    cache_writer.wrap(rma.model_query,
                     path='summary.csv',
                     cache=do_cache,
                     model='Structure',
                     criteria='[graph_id$eq1]',
                     num_rows='all')
```

If you change `to_cache` to `False` and run the code again it will read the data from disk rather than executing the query.

2.6 Visual Coding – Neuropixels

The Visual Coding – Neuropixels project uses high-density extracellular electrophysiology (**Ecephys**) probes to record spikes from a wide variety of regions in the mouse brain. Our experiments are designed to study the activity of the visual cortex and thalamus in the context of passive visual stimulation, but these data can be used to address a wide variety of topics.

Spike-sorted data and metadata are available via the AllenSDK as [Neurodata Without Borders](#) files. However, if you're using the AllenSDK to interact with the data, no knowledge of the NWB data format is required.

2.6.1 Getting Started

To jump right in, check out the [quick start guide \(download .ipynb\)](#), which will show you how to download the data, align spikes to a visual stimulus, and decode natural images from neural activity patterns. For a quick summary of experimental design and data access, see the [cheat sheet](#).

If you would like more example code, the [full example notebook \(download .ipynb\)](#) covers all of the ways to access data for each experiment.


Additional tutorials are available on the following topics:

1. [Data access \(download .ipynb\)](#)
2. [Unit quality metrics \(download .ipynb\)](#)
3. [LFP data analysis \(download .ipynb\)](#)
4. [Receptive field mapping \(download .ipynb\)](#)
5. [Optotagging \(download .ipynb\)](#)

For detailed information about the experimental design, data acquisition, and informatics methods, please refer to our [technical whitepaper](#). AllenSDK API documentation is [available here](#).

A note on terminology: Throughout the SDK, we refer to neurons as “units,” because we cannot guarantee that all the spikes assigned to one unit actually originate from a single cell. Unlike in two-photon imaging, where you can visualize each neuron throughout the entire experiment, with electrophysiology we can only “see” a neuron when it fires a spike. If a neuron moves relative to the probe, or if it’s far away from the probe, some of its spikes may get mixed together with those from other neurons. Because of this inherent ambiguity, we provide a variety of quality metrics to allow you to find the right units for your analysis. Even highly contaminated units contain potentially valuable information about brain states, so we didn’t want to leave them out of the dataset. But certain types of analysis require more stringent quality thresholds, to ensure that all of the included units are well isolated from their neighbors.

2.6.2 Data Processing



`_static/neuropixels_data_processing.png`

Neuropixels probes contain 374 or 383 channels that continuously detect voltage fluctuations in the surrounding neural tissue. Each channel is split into two separate data streams, or “bands,” on the probes. The “spike band” is digitized at 30 kHz, and contains information about action potentials fired by neurons directly adjacent to the probe. The “LFP band” is digitized at 2.5 kHz, and records the low-frequency (<1000 Hz) fluctuations that result from synchronized neural activity over a wider area.

To go from the raw spike-band data to NWB files, we perform the following processing steps:

1. Median-subtraction to remove common-mode noise from the continuous traces
2. High-pass filtering (>150 Hz) and whitening across blocks of 32 channels
3. Spike sorting with [Kilosort2](#), to detect spikes and assign them to individual units
4. Computing the mean waveform for each unit
5. Removing units with artifactual waveforms
6. Computing quality metrics for every unit
7. Computing stimulus-specific tuning metrics

For the LFP band, we:


1. Downsample the signals in space and time (every 4th channel and every 2nd sample)
2. High-pass filter at 0.1 Hz to remove the DC offset from each channel
3. Re-reference to channels outside of the brain to remove common-mode noise

The packaged NWB files contain:

1. Spike times, spike amplitudes, mean waveforms, and quality metrics for every unit
2. Information about the visual stimulus
3. Time series of the mouse's running speed, pupil diameter, and pupil position
4. LFP traces for channels in the brain
5. Experiment metadata

All code for data processing and packaging is available in the [ecephys_spike_sorting](#) and the `ecephys` section of the AllenSDK.

2.6.3 Visual Stimulus Sets



`_static/neuropixels_stimulus_sets.png`

A central aim of the Visual Coding – Neuropixels project is to measure the impact of visual stimuli on neurons throughout the mouse visual system. To that end, all mice viewed one of two possible stimulus sets, known as “Brain Observatory 1.1” or “Functional Connectivity”. Both stimulus sets began with a Gabor stimulus flashed at 81 different locations on the screen, used to map receptive fields of visually responsive units. Next, the mice were shown brief flashes of light or dark, to measure the temporal dynamics of the visual response.

The remainder of the visual stimulus set either consisted of the same stimuli shown in the two-photon experiments


(“Brain Observatory 1.1”), or a subset of those stimuli shown with a higher number of repeats. We also added a dot motion stimulus, to allow us to measure the speed tuning of units across the mouse visual system.

2.6.4 AllenSDK 2.0 and Data Compatability

AllenSDK version 2.0 marks a major update to released Visual Coding Neuropixels datasets. Due to newer versions of pynwb/hdmf having issues reading previously released Visual Coding Neuropixels NWB files and due to the significant reorganization of updated NWB file contents, this release contains breaking changes that necessitate a major version revision. NWB files released prior to 6/11/2020 are not guaranteed to work with the 2.0.0 version of AllenSDK. If you cannot or choose not to re-download the updated NWB files, you can continue using a prior version of AllenSDK (< 2.0.0) to access them. However, no further features or bugfixes for AllenSDK (< 2.0.0) are planned. Data released for other projects (Cell Types, Mouse Connectivity, etc.) are *NOT* affected and will *NOT* need to be re-downloaded.

When using the Visual Coding **EcephysProjectCache** from this updated AllenSDK version, if a **ManifestError** is encountered, this indicates that previously downloaded cached data files need to be removed and re-downloaded. The location these files as well as manifest are user defined and are set when instantiating an **EcephysProjectCache**.

2.6.5 Quality Metrics



`_static/neuropixels_quality_metrics.png`

Every NWB file includes a table of quality metrics, which can be used to assess the completeness, contamination, and stability of units in the recording. By default, we won't show you units below a pre-determined quality threshold; we hide any units that are not present for the whole session (`presence_ratio < 0.95`), that include many contaminating spikes (`isi_violations > 0.5`), or are likely missing a large fraction of spikes (`amplitude_cutoff > 0.1`). However, even contaminated or incomplete units contain information about brain states, and may be of interest to analyze. Therefore, the complete units table can be accessed via special flags in the AllenSDK.

In general, we do not make a distinction between ‘single-unit’ and ‘multi-unit’ activity. There is no obvious place to draw a boundary in the overall distributions of quality metrics, and setting a strict cutoff (e.g. `isi_violations = 0`) will remove a lot of potentially valuable data. We prefer to leave it up to the end user to decide what level of contamination is tolerable. But that means you need to be aware that different units will have different levels of cleanliness.

It should also be noted that all of these metrics assume that the spike waveform is stable throughout the experiment. Given that the probe drifts, on average, about 40 microns over the course of the ~3 hour recordings, this assumption is almost never valid. The resulting changes in waveform shape can cause a unit’s quality to fluctuate. If you’re unsure about a unit’s quality, it can be helpful to plot its spike amplitudes over time. This can make it obvious if it’s drifting below threshold, or if it contains spikes from multiple neurons.

Documentation on the various quality metrics can be found in the [ecephys_spike_sorting](#) repository.

For a detailed discussion of the appropriate way to apply each of these metrics, please check out [this tutorial \(download .ipynb\)](#)

2.6.6 Precomputed Stimulus Metrics

Tables of precomputed metrics are available for download to support population analysis and filtering. The table below describes all of the available metrics. The `get_unit_analysis_metrics()` method will load this table as a [pandas DataFrame](#).

Stimulus	Metric	Field Name
drifting gratings	preferred orientation	pref_ori_dg
	preferred temporal frequency	pref_tf_dg
	global ori. selectivity	g_osi_dg
	global dir. selectivity	g_dsi_dg
	running modulation	run_mod_dg
	running modulation p-value	p_run_mod_dg
	firing rate	firing_rate_dg
	fano factor	fano_dg
	modulation index	mod_idx_dg
	f1/f0	f1_f0_dg
	lifetime sparseness	lifetime_sparseness_dg
	c50 (contrast tuning stimulus)	c50_dg
static gratings	preferred orientation	pref_ori_sg
	preferred spatial frequency	pref_sf_sg
	preferred phase	pref_phase_sg
	global ori. selectivity	g_osi_sg
	running modulation	run_mod_sg
	running modulation p-value	p_run_mod_sg
	firing rate	firing_rate_sg
	fano factor	fano_sg
	lifetime sparseness	lifetime_sparseness_sg
natural scenes	preferred image index	pref_image_ns
	image selectivity	image_selectivity_ns
	running modulation	run_mod_ns
	running modulation p-value	p_run_mod_ns
	firing rate	firing_rate_ns
	fano factor	fano_factor_ns
	lifetime sparseness	lifetime_sparseness_ns
dot motion	preferred speed	pref_speed_dm
	preferred direction	pref_dir_dm

Continued on next page

Table 2 – continued from previous page

Stimulus	Metric	Field Name
	running modulation	run_mod_dm
	running modulation p-value	p_run_mod_dm
	firing rate	firing_rate_dm
	fano factor	fano_factor_dm
	lifetime sparseness	lifetime_sparseness_dm
full-field flashes	on/off ratio	on_off_ratio_fl
	running modulation	run_mod_fl
	running modulation p-value	p_run_mod_fl
	firing rate	firing_rate_fl
	fano factor	fano_factor_fl
	lifetime sparseness	lifetime_sparseness_fl
gabors	RF area	area_rf
	RF elevation	elevation_rf
	RF azimuth	azimuth_rf
	RF p-value	p_value_rf
	running modulation	run_mod_rf
	running modulation p-value	p_run_mod_rf
	firing rate	firing_rate_rf
	fano factor	fano_factor_rf
	lifetime sparseness	lifetime_sparseness_rf

2.7 Visual Behavior - Optical Physiology

The Visual Behavior 2P project used in vivo 2-photon calcium imaging (also called optical physiology, or “ophys”) to measure the activity of populations of genetically identified neurons in the visual cortex of mice performing a visually guided behavioral task (image change detection, described below). We used single- and multi-plane imaging approaches to record the activity of populations of excitatory neurons and two inhibitory classes, Somatostatin (Sst) and Vasoactive Intestinal Peptide (Vip) expressing interneurons, across multiple cortical depths and two visual areas (VSIp and VISl). Each population of neurons was imaged repeatedly over multiple days under different sensory and behavioral contexts, including with familiar and novel stimuli, as well as active behavior and passive viewing conditions. This dataset can be used to evaluate the influence of experience, expectation, and task engagement on neural coding and dynamics.

`_static/visual_behavior_2p/datasets.png`

While 2-photon imaging data was acquired in well-trained mice, the full behavioral training history of all imaged mice is also provided, allowing investigation into task learning, behavioral strategy, and inter-animal variability.

Overall, the dataset includes neural and behavioral measurements from 82 mice, including 3021 behavior training sessions and 551 in vivo 2-photon imaging sessions, resulting in longitudinal recordings from 34,619 cortical cells.

The table below describes the numbers of mice, sessions, and unique recorded neurons for each transgenic line and imaging platform in the dataset:

`_static/visual_behavior_2p/variants_table.png`

2.7.1 GETTING STARTED

First, [install or update the AllenSDK](#), our Python based toolkit for accessing and working with Allen Institute datasets.

Data is provided in [NWB](#) format and can be downloaded using the AllenSDK, or accessed directly via an S3 bucket (instructions provided in notebook #1 below). Regardless of which method of file download you choose, we recommend that you load and interact with the data using the tools provided in the AllenSDK, which have been designed to simplify data access and subsequent analysis. No knowledge of the NWB file format is required.

Specific information about how Visual Behavior Optical Physiology data is stored in NWB files and how AllenSDK accesses NWB files can be found [here](#).

To get started, check out these jupyter notebooks to learn how to:

- 1) [Download data using the AllenSDK or directly from our Amazon S3 bucket \(download .ipynb\)](#)
- 2) [Identify experiments of interest using the dataset manifest \(download .ipynb\)](#)
- 3) [Load and visualize data from a 2-photon imaging experiment \(download .ipynb\)](#)
- 4) [Examine the full training history of one mouse \(download .ipynb\)](#)
- 5) [Compare behavior and neural activity across different trial types in the task \(download .ipynb\)](#)

For a description of available AllenSDK methods and attributes for data access, see this [further documentation](#).

For detailed information about the experimental design, data acquisition, and informatics methods, please refer to our [technical whitepaper](#).


If you have questions about the dataset that aren't addressed by the whitepaper or any of our tutorials, please reach out by posting at <https://community.brain-map.org/>

2.7.2 CHANGE DETECTION TASK

`_static/visual_behavior_2p/task.png`

We trained mice to perform a go/no-go visual change detection task in which they learned to lick a spout in response to changes in stimulus identity to earn a water reward. Visual stimuli are continuously presented over a 1-hour session, with no explicit cue to indicate the start of a trial. Mice are free to run on a circular disk during the session.


We used a standardized procedure to progress mice through a series of training stages, with transitions between stages determined by specific advancement criteria. First, mice learned to detect changes in the orientation of full field static grating stimuli. Next, a 500ms inter stimulus interval period with mean luminance gray screen was added between the 250ms stimulus presentations, incorporating a short-term memory component to the task. Once mice successfully and consistently performed the orientation change detection with flashed gratings, they moved to the image change detection version of the task. During image change detection, 8 natural scene images were presented during each behavioral session, for a total of 64 possible image transitions. When behavioral performance again reached criterion, mice were transitioned to the 2-photon imaging stage in which they performed the task under a microscope to allow simultaneous measurement of neural activity and behavior.



`_static/visual_behavior_2p/automated_training.png`

Behavioral training data for mice progressing through these stages of task learning is accessible via the **BehaviorSession** class of the AllenSDK or the `get_behavior_session()` method of the **VisualBehaviorOphysProjectCache**. Each **BehaviorSession** contains the following data streams, event times, and metadata:

- Running speed
- Lick times
- Reward times
- Stimulus presentations
- Behavioral trial information
- Mouse metadata (age, sex, genotype, etc)



`_static/visual_behavior_2p/behavior_timeseries.png`

2.7.3 2-PHOTON IMAGING DATASET

Once mice are well-trained on the image change detection task, they transition to performing the behavior under a 2-photon microscope. During the imaging phase, mice undergo multiple **session types**, allowing measurement of neural activity across different sensory and behavioral contexts.

`_static/visual_behavior_2p/experiment_design.png`

Mice initially perform the task under the microscope with the same set of images they observed during training, which have become highly familiar (each image is viewed thousands of times during training). Mice also undergo several sessions with a novel image set that they had not seen prior to the 2-photon imaging portion of the experiment. Interleaved between active behavior sessions, are passive viewing sessions where the mice are given their daily water before the session (and are thus satiated) and view the task stimuli with the lick spout retracted so they are unable to earn water rewards. This allows investigation of the impact of motivation and attention on patterns of neural activity. Finally, stimuli were randomly omitted with a 5% probability, resulting in an extended gray screen period between two presentations of the same stimulus, and disrupting the expected cadence of stimulus presentations. Stimuli were only omitted during the 2-photon imaging sessions (not during training), and change stimuli were never omitted.

We used both single- and multi-plane 2-photon imaging to record the activity of GCaMP6 expressing cells in populations of excitatory (Slc17a7-IRES2-Cre;Camk2a-tTA;Ai93(TITL-GCaMP6)) and inhibitory (Vip-IRES-Cre;Ai148(TIT2L-GC6f-ICL-tTA2) & Sst-IRES-Cre;Ai148(TIT2L-GC6f-ICL-tTA2)) neurons. Imaging took place between 75-400um below the cortical surface.

`_static/visual_behavior_2p/cre_lines.png`

The data collected in a single continuous recording is defined as a **session**. For single-plane imaging experiments, there is only one imaging plane (referred to as an **experiment**) per session. For multi-plane imaging experiments, there can be up to 8 imaging planes (aka 8 experiments) per session. Due to our strict QC process, described below, not all multi-plane imaging sessions have exactly 8 experiments, as some imaging planes did not meet our data quality criteria.

We aimed to track the activity of single neurons across the session types described above by targeting the same population of neurons over multiple recording sessions, with only one session recorded per day for a given mouse. The collection of imaging sessions for a given population of cells, belonging to a single imaging plane measured across days, is called a **container**. A container can include between 3 and 11 separate sessions for a given imaging plane. Mice imaged with the multi-plane 2-photon microscope can have multiple containers, one for each imaging plane recorded across multiple sessions. The session types available for a given container can vary, due to our selection criteria to ensure data quality (described below).

Thus, each mouse can have one or more **containers**, each representing a unique imaging plane (**experiment**) that has been targeted across multiple recording **sessions**, under different behavioral and sensory conditions (**session types**).

`_static/visual_behavior_2p/data_structure.png`

The **BehaviorOphysExperiment** class in the AllenSDK (or the `get_behavior_ophys_experiment()` method of the **VisualBehaviorOphysProjectCache**) provides all data for a single imaging plane, recorded in a single session, and contains the following data streams in addition to the behavioral data described above:

- Max intensity image
- Average intensity image
- Segmentation masks
- dF/F traces (baseline corrected, normalized fluorescence traces)
- Corrected fluorescence traces (neuropil subtracted and demixed, but not normalized)
- Events (detected with an L0 event detection algorithm)
- Pupil position
- Pupil area

`_static/visual_behavior_2p/behavior_and_ophys_timeseries.png`

2.7.4 DATA PROCESSING

Each 2-photon movie is processed through a series of steps to obtain single cell traces of baseline-corrected fluorescence (dF/F) and extracted events, that are packaged into NWB files along with stimulus and behavioral information, as well as other metadata.

`_static/visual_behavior_2p/data_processing.png`

Detailed descriptions of data processing steps can be found in the technical white paper, as well as our [data processing repository](#).

2.7.5 QUALITY CONTROL

Every 2-photon imaging session was carefully evaluated for a variety of quality control criteria to ensure that the final dataset is of the highest quality possible. Sessions or imaging planes that do not meet our criteria are excluded from the dataset in this release. These are a few of the key aspects of the data that are evaluated:

- intensity drift
- image saturation or bleaching
- z-drift over the course of a session
- accuracy of session-to-session field of view matching
- excessive or uncorrectable motion in the image
- uncorrectable crosstalk between simultaneously recorded multiscope planes
- errors affecting temporal alignment of data streams
- hardware or software failures
- brain health
- animal stress

2.7.6 SUMMARY OF AVAILABLE DATA

Behavior	Physiology	Metadata
Running speed	Max intensity projection image	Mouse genotype, age, sex
Licks	Average projection image	Date of acquisition
Rewards	Segmentation mask image	Imaging parameters
Pupil area	Cell specimen table	Task parameters
Pupil position	Cell ROI masks	Session type
Stimulus presentations table	Corrected fluorescence traces	Stimulus images
Trials table	dF/F activity traces	Performance metrics
Stimulus timestamps	Detected events	
	Ophys timestamps	

2.7.7 DATA FILE CHANGELOG

v0.3.0

13 sessions were labeled with the wrong session_type in v0.2.0. We have corrected that error. The offending sessions were

behavior_session_id	ophys_session_id	session_type_v0.2.0	session_type_v0.3.0	
875020233		OPHYS_3_images_A	OPHYS_2_images_A_passive	
902810506		TRAINING_4_images_B_training	TRAINING_3_images_B_10uL_reward	
914219174		OPHYS_0_images_B_habituation	TRAINING_5_images_B_handoff	ready
863571063		TRAINING_5_images_A_handoff	TRAINING_1_gratings	ready
974330793		OPHYS_0_images_B_habituation	TRAINING_5_images_B_handoff	ready
863571072		OPHYS_5_images_B_passive	TRAINING_4_images_A_training	
1010972317		OPHYS_4_images_A	OPHYS_3_images_B	
1011659817		OPHYS_5_images_A_passive	OPHYS_4_images_A	
1003302686	1003277121	OPHYS_6_images_A	OPHYS_5_images_A_passive	
863571054		OPHYS_7_receptive_field_mapping	TRAINING_5_images_A_epilogue	
974282914	974167263	OPHYS_6_images_B	OPHYS_5_images_B_passive	
885418521		OPHYS_1_images_A	TRAINING_5_images_A_handoff	lapsed
915739774		OPHYS_1_images_A	OPHYS_0_images_A_habituation	

The Allen SDK currently focuses on models generated from electrophysiology data in the Allen Cell Types Database. There are two classes of models available for download: biophysical models and generalize leaky integrate-and-fire models.

3.1 Generalized LIF Models

The Allen Cell Types Database contains Generalized Leaky Integrate and Fire (GLIF) models that simulate the firing behavior of neurons at five levels of complexity. Review the GLIF technical [white paper](#) for details on these models and how their parameters were optimized.

The Allen SDK GLIF simulation module is an explicit time-stepping simulator that evolves a neuron's simulated voltage over the course of an input current stimulus. The module also tracks the neuron's simulated spike threshold and registers action potentials whenever voltage surpasses threshold. Action potentials initiate reset rules that update voltage, threshold, and (optionally) trigger afterspike currents.

The GLIF simulator in this package has a modular architecture that enables users to choose from a number of dynamics and reset rules that update the simulation's voltage, spike threshold, and afterspike currents during the simulation. The GLIF package contains a built-in set of rules, however developers can plug in custom rule implementations provided they follow a simple argument specification scheme.

The Allen SDK GLIF simulator was developed and tested with Python 2.7.9, installed as part of [Anaconda Python](#) distribution version 2.1.0.

The rest of this page provides examples demonstrating how to download models, examples of simulating these models, and general GLIF model documentation.

Note: the GLIF simulator module is still under heavy development and may change significantly in the future.

3.1.1 Downloading GLIF Models

There are two ways to download files necessary to run a GLIF model. The first way is to visit <http://celltypes.brain-map.org> and find cells that have GLIF models available for download. The electrophysiology details page for a cell has a neuronal model download link. Specifically:

1. Click ‘More Options +’ and filter for GLIF models.
2. Click the electrophysiology thumbnail for a cell on the right hand panel.
3. Choose a GLIF model from the ‘Show model responses’ dropdown.
4. Scroll down to the model response click ‘Download model’.

One such link (for a simple LIF neuronal model, ID 566302806), would look like this:

```
http://api.brain-map.org/neuronal_model/download/566302806
```

This link returns .zip archive containing the neuron configuration file and sweep metadata required to simulate the model with stimuli applied to the cell. Specifically, the .zip archive will contain:

- **472423251_neuron_config.json**: JSON config file for the GLIF model
- **ephys_sweeps.json**: JSON with metadata for sweeps presented to the cell
- **neuronal_model.json**: JSON with general metadata for the cell

If you would like to reproduce the model traces seen in the Cell Types Database, you can download an NWB file containing both the stimulus and cell response traces via a ‘Download data’ link on the cell’s electrophysiology page. See the [NWB](#) description section for more details on the NWB file format.

You can also download all of these files, including the cell’s NWB file, using the *GlifApi* class:

```
from allensdk.api.queries.glif_api import GlifApi
from allensdk.core.cell_types_cache import CellTypesCache
import allensdk.core.json_utilities as json_utilities

neuronal_model_id = 566302806

# download model metadata
glif_api = GlifApi()
nm = glif_api.get_neuronal_models_by_id([neuronal_model_id])[0]

# download the model configuration file
nc = glif_api.get_neuron_configs([neuronal_model_id])[neuronal_model_id]
neuron_config = glif_api.get_neuron_configs([neuronal_model_id])
json_utilities.write('neuron_config.json', neuron_config)

# download information about the cell
ctc = CellTypesCache()
ctc.get_ephys_data(nm['specimen_id'], file_name='stimulus.nwb')
ctc.get_ephys_sweeps(nm['specimen_id'], file_name='ephys_sweeps.json')
```

3.1.2 Running a GLIF Simulation

To run a GLIF simulation, the most important file you need is the `neuron_config` JSON file. You can use this file to instantiate a simulator and feed in your own stimulus:

```

import allensdk.core.json_utilities as json_utilities
from allensdk.model.glif.glif_neuron import GlifNeuron

# initialize the neuron
neuron_config = json_utilities.read('neuron_config.json')['566302806']
neuron = GlifNeuron.from_dict(neuron_config)

# make a short square pulse. stimulus units should be in Amps.
stimulus = [ 0.0 ] * 100 + [ 10e-9 ] * 100 + [ 0.0 ] * 100

# important! set the neuron's dt value for your stimulus in seconds
neuron.dt = 5e-6

# simulate the neuron
output = neuron.run(stimulus)

voltage = output['voltage']
threshold = output['threshold']
spike_times = output['interpolated_spike_times']

```

Note: The GLIF simulator does not simulate during action potentials. Instead it inserts NaN values for a fixed number of time steps when voltage surpasses threshold. The simulator skips `neuron.spike_cut_length` time steps after voltage surpasses threshold.

To reproduce the model's traces displayed on the Allen Cell Types Database web page, the Allen SDK provides the `allensdk.core.model.glif.simulate_neuron` module for simulating all sweeps presented to a cell and storing them in the NWB format:

```

import allensdk.core.json_utilities as json_utilities

from allensdk.model.glif.glif_neuron import GlifNeuron
from allensdk.model.glif.simulate_neuron import simulate_neuron

neuron_config = json_utilities.read('neuron_config.json')['566302806']
ephys_sweeps = json_utilities.read('ephys_sweeps.json')
ephys_file_name = 'stimulus.nwb'

neuron = GlifNeuron.from_dict(neuron_config)

sweep_numbers = [ s['sweep_number'] for s in ephys_sweeps if s['stimulus_units'] ==
↳ 'Amps' ]
sweep_numbers = sweep_numbers[:1] # for the sake of a speedy example, just run the_
↳ first one
simulate_neuron(neuron, sweep_numbers, ephys_file_name, ephys_file_name, 0.05)

```

Warning: These stimuli are sampled at a very high resolution (200kHz), and a given cell can have many sweeps. This process can take over an hour.

The `simulate_neuron` function call simulates all sweeps in the NWB file. Because the same NWB file is being used for both input and output, the cell's response traces will be overwritten as stimuli are simulated. `simulate_neuron` optionally accepts a value which will be used to overwrite these NaN values generated during action potentials (in this case 0.05 Volts).

If you would like to run a single sweep instead of all sweeps, try the following:

```
import allensdk.core.json_utilities as json_utilities
from allensdk.model.glif.glif_neuron import GlifNeuron
from allensdk.core.nwb_data_set import NwbDataSet

neuron_config = json_utilities.read('neuron_config.json')['566302806']
ephys_sweeps = json_utilities.read('ephys_sweeps.json')
ephys_file_name = 'stimulus.nwb'

# pull out the stimulus for the current-clamp first sweep
ephys_sweep = next( s for s in ephys_sweeps
                    if s['stimulus_units'] == 'Amps' )
ds = NwbDataSet(ephys_file_name)
data = ds.get_sweep(ephys_sweep['sweep_number'])
stimulus = data['stimulus']

# initialize the neuron
# important! update the neuron's dt for your stimulus
neuron = GlifNeuron.from_dict(neuron_config)
neuron.dt = 1.0 / data['sampling_rate']

# simulate the neuron
output = neuron.run(stimulus)

voltage = output['voltage']
threshold = output['threshold']
spike_times = output['interpolated_spike_times']
```

Note: The dt value provided in the downloadable GLIF neuron configuration files does not correspond to the sampling rate of the original stimulus. Stimuli were subsampled and filtered for parameter optimization. Be sure to overwrite the neuron's dt with the correct sampling rate.

If you would like to plot the outputs of this simulation using numpy and matplotlib, try:

```
import numpy as np
import matplotlib.pyplot as plt

voltage = output['voltage']
threshold = output['threshold']
interpolated_spike_times = output['interpolated_spike_times']
spike_times = output['interpolated_spike_times']
interpolated_spike_voltages = output['interpolated_spike_voltage']
interpolated_spike_thresholds = output['interpolated_spike_threshold']
grid_spike_indices = output['spike_time_steps']
grid_spike_times = output['grid_spike_times']
after_spike_currents = output['AScurrents']

# create a time array for plotting
time = np.arange(len(stimulus))*neuron.dt

plt.figure(figsize=(10, 10))

# plot stimulus
plt.subplot(3,1,1)
plt.plot(time, stimulus)
plt.xlabel('time (s)')
```

(continues on next page)

(continued from previous page)

```

plt.ylabel('current (A)')
plt.title('Stimulus')

# plot model output
plt.subplot(3,1,2)
plt.plot(time, voltage, label='voltage')
plt.plot(time, threshold, label='threshold')

if grid_spike_indices is not None:
    plt.plot(interpolated_spike_times, interpolated_spike_voltages, 'x',
             label='interpolated spike')

    plt.plot((grid_spike_indices-1)*neuron.dt, voltage[grid_spike_indices-1], '.',
             label='last step before spike')

plt.xlabel('time (s)')
plt.ylabel('voltage (V)')
plt.legend(loc=3)
plt.title('Model Response')

# plot after spike currents
plt.subplot(3,1,3)
for ii in range(np.shape(after_spike_currents)[1]):
    plt.plot(time, after_spike_currents[:,ii])
plt.xlabel('time (s)')
plt.ylabel('current (A)')
plt.title('After Spike Currents')

plt.tight_layout()
plt.show()

```

Note: There both interpolated spike times and grid spike times. The grid spike is the first time step where the voltage is higher than the threshold. Note that if you try to plot the voltage at the grid spike indices the output will be NaN. The interpolated spike is the calculated intersection of the threshold and voltage between the time steps.

3.1.3 GLIF Configuration

Instances of the *GlifNeuron* class require many parameters for initialization. Fixed neuron parameters are stored directly as properties on the class instance:

Parameter	Description	Units	Type
EI	resting potential	Volts	float
dt	time duration of each simulation step	seconds	float
R_input	input resistance	Ohms	float
C	capacitance	Farads	float
asc_vector	afterspike current coefficients	Amps	np.array
spike_cut_length	spike duration	time steps	int
th_inf	instantaneous threshold	Volts	float
th_adapt	adapted threshold	Volts	float

Some of these fixed parameters were optimized to fit Allen Cell Types Database electrophysiology data. Optimized coefficients for these parameters are stored by name in the `neuron.coeffs` dictionary. For more details on which

parameters were optimized, please see the technical [white paper](#).

The `GlifNeuron` class has six methods that can be customized: three rules for updating voltage, threshold, and afterspike currents during the simulation; and three rules for updating those values when a spike is detected (voltage surpasses threshold).

Method Type	Description
<code>voltage_dynamics_method</code>	Update simulation voltage for the next time step.
<code>threshold_dynamics_method</code>	Update simulation threshold for the next time step.
<code>AScurrent_dynamics_method</code>	Update afterspike current coefficients for the next time step.
<code>voltage_reset_method</code>	Reset simulation voltage after a spike occurs.
<code>threshold_reset_method</code>	Reset simulation threshold after a spike occurs.
<code>AScurrent_reset_method</code>	Reset afterspike current coefficients after a spike occurs.

The GLIF neuron configuration files available from the Allen Brain Atlas API use built-in methods, however you can supply your own custom method if you like:

```
# define your own custom voltage reset rule
# this one linearly scales the input voltage
def custom_voltage_reset_rule(neuron, voltage_t0, custom_param_a, custom_param_b):
    return custom_param_a * voltage_t0 + custom_param_b

# initialize a neuron from a neuron config file
neuron_config = json_utilities.read('neuron_config.json')['566302806']
neuron = GlifNeuron.from_dict(neuron_config)

# configure a new method and overwrite the neuron's old method
method = neuron.configure_method('custom', custom_voltage_reset_rule,
                                { 'custom_param_a': 0.1, 'custom_param_b': 0.0 })
neuron.voltage_reset_method = method

output = neuron.run(stimulus)
```

Notice that the function is allowed to take custom parameters (here `custom_param_a` and `custom_param_b`), which are configured on method initialization from a dictionary. For more details, see the documentation for the `GlifNeuron` and `GlifNeuronMethod` classes.

3.1.4 Built-in Dynamics Rules

The job of a dynamics rule is to describe how the simulator should update the voltage, spike threshold, and afterspike currents of the simulator at a given simulation time step.

Voltage Dynamics Rules

These methods update the output voltage of the simulation. They all expect a voltage, afterspike current vector, and current injection value to be passed in by the `GlifNeuron`. All other function parameters must be fixed using the `GlifNeuronMethod` class. They all return an updated voltage value.

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_forward_euler()
```

Threshold Dynamics Rules

These methods update the spike threshold of the simulation. They all expect the current threshold and voltage values of the simulation to be passed in by the `GlifNeuron`. All other function parameters must be fixed using the `GlifNeuronMethod` class. They all return an updated threshold value.

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_three_components_exact()
```

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_spike_component()
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_inf()
```

Afterspike Current Dynamics Rules

These methods expect current afterspike current coefficients, current time step, and time steps of all previous spikes to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated afterspike current array.

```
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_exp()
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_none()
```

3.1.5 Built-in Reset Rules

The job of a reset rule is to describe how the simulator should update the voltage, spike threshold, and afterspike currents of the simulator after the simulator has detected that the simulated voltage has surpassed threshold.

Voltage Reset Rules

These methods update the output voltage of the simulation after voltage has surpassed threshold. They all expect a voltage to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated voltage value.

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_zero()
allensdk.model.glif.glif_neuron_methods.reset_voltage_v_before()
```

Threshold Reset Rules

These methods update the spike threshold of the simulation after a spike has been detected. They all expect the current threshold and the reset voltage value of the simulation to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated threshold value.

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_inf()
allensdk.model.glif.glif_neuron_methods.reset_threshold_three_components()
```

Afterspike Reset Rules

These methods expect current afterspike current coefficients to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated afterspike current array.

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_none()
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_sum()
```

3.2 Biophysical Models

The Allen Cell Types Database contains biophysical models that characterize the firing behavior of neurons measured in slices through current injection by a somatic whole-cell patch clamp electrode. These models contain a set of 10 active conductances placed at the soma and use the reconstructed 3D morphologies of the modeled neurons. The biophysical modeling [technical white paper](#) contains details on the specific construction of these models and the optimization of the model parameters to match the experimentally-recorded firing behaviors.

The biophysical models are run with the [NEURON](#) simulation environment. The Allen SDK package contains libraries that assist in downloading and setting up the models available on the Allen Institute web site for users to run using NEURON. The examples and scripts provided run on Linux using the bash shell.

3.2.1 Prerequisites

You must have NEURON with the Python interpreter enabled and the Allen SDK installed.

The Allen Institute perisomatic biophysical models were generated using NEURON [version v7.4.rel-1370](#). Instructions for compiling NEURON with the Python interpreter are available from the NEURON team under the heading [Installation with Python as an alternative interpreter](#). The Allen SDK is compatible with Python version 2.7.9, included in the Anaconda 2.1.0 distribution.

Instructions for optional [Docker installation](#) are also available.

Note: Building and installing NEURON with the Python wrapper enabled is not always easy. This page targets users that have a background in NEURON usage and installation.

3.2.2 Downloading Biophysical Models

There are two ways to download files necessary to run a biophysical model. The first way is to visit <http://celltypes.brain-map.org> and find cells that have biophysical models available for download. The electrophysiology details page for a cell has a neuronal model download link. Specifically:

1. Click ‘More Options+’
2. Check ‘Models -> Biophysical - perisomatic’ or ‘Biophysical - all active’
3. Use the Filters, Cell Location and Cell Feature Filters to narrow your results.
4. Click on a Cell Summary to view the Mouse Experiment Electrophysiology.
5. Click the “download data” link to download the NWB stimulus and response file.
6. Click “show model response” and select ‘Biophysical - perisomatic’ or ‘Biophysical - all active’.
7. Scroll down and click the ‘Biophysical - perisomatic’ or ‘Biophysical - all active’ “download model” link.

This may also be done programmatically. The neuronal model id can be found to the left of the corresponding ‘Biophysical - perisomatic’ or ‘Biophysical - all active’ “download model” link.

```
from allensdk.api.queries.biophysical_api import BiophysicalApi

bp = BiophysicalApi()
bp.cache_stimulus = True # change to False to not download the large stimulus NWB file
neuronal_model_id = 472451419 # get this from the web site as above
bp.cache_data(neuronal_model_id, working_directory='neuronal_model')
```

More help can be found in the [online help](#) for the Allen Cell Types Database web application.

3.2.3 Directory Structure

The structure of the directory created looks like this. It includes stimulus files, model parameters, morphology, cellular mechanisms and application configuration.

```
neuronal_model
|-- manifest.json
|-- 472451419_fit.json
```

(continues on next page)

(continued from previous page)

```

|-- Nr5a1-Cre_Ai14_IVSCC_-169248.04.02.01.nwb
|-- Nr5a1-Cre_Ai14_IVSCC_-169248.04.02.01_403165543_m.swc
|-- modfiles
|   |--CaDynamics.mod
|   |--Ca_HVA.mod
|   |--Ca_LVA.mod
|   |--Ih.mod
|   `--...etc.
|
|--x86_64
`---work

```

3.2.4 Running the Simulation (Linux shell prompt)

All of the sweeps available from the web site are included in manifest.json and will be run by default. This can take some time.

```

cd neuronal_model
nrnivmodl ./modfiles # compile the model (only needs to be done once)
python -m allensdk.model.biophysical.runner manifest.json # perisomatic models
python -m allensdk.model.biophysical.runner manifest.json # legacy all-active models
# new all-active models (axon replaced by a 60 micron long 1 micron diameter stub)
python -m allensdk.model.biophysical.runner manifest.json --axon_type stub

```

3.2.5 Selecting a Specific Sweep

The sweeps are listed in manifest.json. You can remove all of the sweep numbers that you do not want run.

3.2.6 Simulation Main Loop

The top level script is in the `run()` method of the `allensdk.model.biophysical.runner` module. The implementation of the method is discussed here step-by-step:

First configure NEURON based on the configuration file, which was read in from the command line at the very bottom of the script.

`run()`:

```

# configure NEURON -- this will infer model type (perisomatic vs. all-active)
utils = Utils.create_utils(description)
h = utils.h

```

The next step is to get the path of the morphology file and pass it to NEURON.

```

# configure model
manifest = description.manifest
morphology_path = description.manifest.get_path('MORPHOLOGY')
utils.generate_morphology(morphology_path.encode('ascii', 'ignore'))
utils.load_cell_parameters()

```

Then read the stimulus and recording configuration and configure NEURON

```
# configure stimulus and recording
stimulus_path = description.manifest.get_path('stimulus_path')
nwb_out_path = manifest.get_path("output")
output = NwbDataSet(nwb_out_path)
run_params = description.data['runs'][0]
sweeps = run_params['sweeps']
junction_potential = description.data['fitting'][0]['junction_potential']
mV = 1.0e-3
```

Loop through the stimulus sweeps and write the output.

```
# run sweeps
for sweep in sweeps:
    utils.setup_iclamp(stimulus_path, sweep=sweep)
    vec = utils.record_values()

    h.finitialize()
    h.run()

    # write to an NWB File
    output_data = (numpy.array(vec['v']) - junction_potential) * mV
    output.set_sweep(sweep, None, output_data)
```

3.2.7 Customization

Much of the code in the perisomatic simulation is not core Allen SDK code. The `runner.py` script largely reads the configuration file and calls into methods in the `Utils` class. `Utils` is a subclass of the `HocUtils` class, which provides access to objects in the NEURON package. The various methods called by the `runner.py` script are implemented here, including: `generate_morphology()`, `load_cell_parameters()`, `setup_iclamp()`, `read_stimulus()` and `record_values()`.

`Utils`:

```
from allensdk.model.biophys_sim.neuron.hoc_utils import HocUtils

.....

class Utils(HocUtils):
    .....

    def __init__(self, description):
        super(Utils, self).__init__(description)
    ....
```

To create a biophysical model using your own software or data, simply model your directory structure on one of the downloaded simulations or one of the examples below. Add your own `runner.py` and `utils.py` module to the simulation directory.

Compile the `.mod` files using NEURON's `nrnivmodl` command (Linux shell):

```
nrnivmodl modfiles
```

Then call your runner script directly, passing in the manifest file to your script:

```
python runner.py manifest.json
```

The output from your simulation and any intermediate files will go in the work directory.

3.2.8 Examples

A minimal example (`simple_example.tgz`) and a multicell example (`multicell_example.tgz`) are available to download as a starting point for your own projects.

Each example provides its own `utils.py` file along with a main script (Linux shell) and supporting configuration files.

`simple_example.tgz`:

```
tar xvzf simple_example.tgz
cd simple
nrnivmodl modfiles
python simple.py
```

`multicell_example.tgz`:

```
tar xvzf multicell_example.tgz
cd multicell
nrnivmodl modfiles
python multi.py
python multicell_diff.py
```

3.2.9 Exporting Output to Text Format or Image

This is an example of using the AllenSDK to save a response voltage to other formats.

```
from allensdk.core.dat_utilities import \
    DatUtilities
from allensdk.core.nwb_data_set import \
    NwbDataSet
import numpy as np
import matplotlib
matplotlib.use("Agg")
import matplotlib.pyplot as plt

nwb_file = '313862020.nwb'
sweep_number = 52
dat_file = '313862020_%d.dat' % (sweep_number)

nwb = NwbDataSet(nwb_file)
sweep = nwb.get_sweep(sweep_number)

# read v and t as numpy arrays
v = sweep['response']
dt = 1.0e3 / sweep['sampling_rate']
num_samples = len(v)
t = np.arange(num_samples) * dt

# save as text file
data = np.transpose(np.vstack((t, v)))
with open (dat_file, "w") as f:
    np.savetxt(f, data)

# save image using matplotlib
fig, ax = plt.subplots(nrows=1, ncols=1)
ax.plot(t, v)
```

(continues on next page)

(continued from previous page)

```
ax.set_title("Sweep %s" % (sweep_number))
fig.savefig('out.png')
```

3.2.10 Model Description Files

Basic Structure

A model description file is simply a JSON object with several sections at the top level and an array of JSON objects within each section.

```
{
  "cell_section": [
    {
      "name": "cell 1",
      "shape": "pyramidal",
      "position": [ 0.1, 0.2, 0.3 ]
    },
    {
      "name": "cell 2",
      "shape": "glial",
      "position": [ 0.1, 0.2, 0.3 ]
    }
  ],
  "extra": [
    { "what": "wood",
      "who": "woodchuck"
    }
  ]
}
```

Even if a section contains no objects or only one object the array brackets must be present.

Objects Within Sections

While no restrictions are enforced on what kinds of objects are stored in a section, some rules of thumb make the file easier to work with.

1. All objects within a section are the same structure. Common operations on a section are to display it as a table, iterate over it, load from or write to a spreadsheet or csv file. These operations are all easier if the section is fairly homogeneous.
2. Objects are not deeply nested. While some shallow nesting is often useful, deep nesting such as a tree structure is not recommended. It makes interoperability with other tools and data formats more difficult.
3. Arrays are allowed, though they should not be deeply nested either.
4. Object member values should be literals. Do not use pickled classes, for example.

Comment Lines

The JSON specification does not allow comments. However, the Allen SDK library applies a preprocessing stage to remove C++-style comments, so they can be used in description files.

Multi-line comments should be surrounded by `/* */` and single-line comments start with `//`. Commented description files will not be recognized by strict json parsers unless the comments are stripped.

commented.json:

```
{
  /*
   * multi-line comment
   */
  "section1": [
    {
      "name": "simon" // single line comment
    }
  ]
}
```

Split Description Files by Section

A model description can be split into multiple files by putting some sections in one file and other sections into another file. This can be useful if you want to put a topology of cells and connections in one file and experimental conditions and stimulus in another file. The resulting structure in memory will behave the same way as if the files were not split. This allows a small experiment to be described in a single file and large experiments to be more modular.

cells.json:

```
{
  "cell_section": [
    {
      "name": "cell 1",
      "shape": "pyramidal",
      "position": [ 0.1, 0.2, 0.3 ]
    },
    {
      "name": "cell 2",
      "shape": "glial",
      "position": [ 0.1, 0.2, 0.3 ]
    }
  ]
}
```

extras.json:

```
{
  "extra": [
    {
      "what": "wood",
      "who": "woodchuck"
    }
  ]
}
```

Split Sections Between Description Files

If two description files containing the same sections are combined, the resulting description will contain objects from both files. This feature allows sub-networks to be described in separate files. The sub-networks can then be composed into a larger network with an additional description of the interconnections.

network1.json:

```
/* A self-contained sub-network */
{
  "cells": [
    { "name": "cell1" },
    { "name": "cell2" }
  ],
  /* intra-network connections */
  "connections": [
    { "source": "cell1", "target" : "cell2" }
  ]
}
```

network2.json:

```
/* Another self-contained sub-network */
{
  "cells": [
    { "name": "cell3" },
    { "name": "cell4" }
  ],
  "connections": [
    { "source": "cell3", "target" : "cell4" }
  ]
}
```

interconnect.json:

```
{
  // the additional connections needed to
  // connect the network1 and network2
  // into a ring topology.
  "connections": [
    { "source": "cell2", "target": "cell3" },
    { "source": "cell4", "target": "cell1" }
  ]
}
```

3.2.11 Resource Manifest

JSON has many advantages. It is widely supported, readable and easy to parse and edit. As data sets get larger or specialized those advantages diminish. Large or complex models and experiments generally need more than a single model description file to completely describe an experiment. A manifest file is a way to describe all of the resources needed within the Allen SDK description format itself.

The manifest section is named “manifest” by default, though it is configurable. The objects in the manifest section each specify a directory, file, or file pattern. Files and directories may be organized in a parent-child relationship.

A Simple Manifest

This is a simple manifest file that specifies the BASEDIR directory using “.”, meaning the current directory:

```
{
  "manifest": [
```

(continues on next page)

(continued from previous page)

```

    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    }
  ] }
}
```

Parent Child Relationships

Adding the optional “parent_key” member to a manifest object creates a parent-child relation. In this case WORKDIR will be found in “./work”:

```

{
  "manifest": [
    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    },
    {
      "key": "WORKDIR",
      "type": "dir",
      "spec": "/work",
      "parent_key": "BASEDIR"
    }
  ] }
}
```

File Spec Patterns

Files can be specified using the type “file” instead of “dir”. If a sequence of many files is needed, the spec may contain patterns to indicate where the sequence number (%d) or string (%s) will be interpolated:

```

{
  "manifest": [
    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    },
    {
      "key": "voltage_out_cell_path",
      "type": "file",
      "spec": "v_out-cell-%d.dat",
      "parent_key": "BASEDIR"
    }
  ] }
}
```

Split Manifest Files

Manifest files can be split like any description file. This allows the specification of a general directory structure in a shared file and specific files in a separate configuration (i.e. stimulus and working directory)

Extensions

To date, manifest description files have not been used to reference URLs that provide model data, but it is a planned future use case.

3.2.12 Further Reading

- [NEURON](#)
- [Python](#)
- [JSON](#)

CHAPTER 4

Examples

Take a look at the table below for a list of SDK example notebooks and scripts.

Description	Link
Introduction to the Mouse Connectivity Atlas	Jupyter notebook (download .ipynb)
Introduction to the Cell Types Database	Jupyter notebook (download .ipynb)
Introduction to the Brain Observatory	Jupyter notebook (download .ipynb)
Brain Observatory Stimulus Manipulation	Jupyter notebook (download .ipynb)
Brain Observatory Tuning Analysis	Jupyter notebook (download .ipynb)
Brain Observatory Receptive Field Analysis	Jupyter notebook (download .ipynb)
Brain Observatory Cell Specimen ID Mapping	Jupyter notebook (download .ipynb)
Brain Observatory Monitor	Jupyter notebook (download .ipynb)
Dynamic Brain Workshop 2015 experiment detail	Jupyter notebook (download .ipynb)
Stimulating a biophysical model with a square pulse	Jupyter notebook (download .ipynb)
Using a Reference Space	Jupyter notebook (download .ipynb)
Downloading Images	Jupyter notebook (download .ipynb)
Visual Coding Neuropixels Quick Start	Jupyter notebook (download .ipynb)
Visual Coding Neuropixels Reference	Jupyter notebook (download .ipynb)

CHAPTER 5

Authors

Michael Buice @mabuice
Nicholas Cain @nicain
Tom Chartrand @tmchartrand
Kael Dai @kaeldai
Saskia de Vries @saskiad
David Feng @dyf
Tim Fliss @tfliss
Marina Garrett @matchings
Richard Gerkin @rgerkin
Nathan Gouwens @gouwens
Nile Graddis @nilegraddis
Sergey Gratiy @sgratiy
@jennan @jennan
Xiaoxuan Jia @jiaxx
Wesley Jones @wesley-jones
Justin Kiggins @neuromusic
Joe Knox @jknox13
Peter Ledochowitsch @ledochowitsch
Nicholas Mei @njmei
Chris Mochizuki @mochic
Ani Nandi @anirban6908
Gabe Ocker @gocker

Michael Oliver @the-moliver

Doug Ollerenshaw @dougollerenshaw

Jed Perkins @jfperkins

Alex Piet @alexpier

Nick Ponvert @nickponvert

Kat Schelonka @kschelonka

Shu Shi @shus2018

Josh Siegle @jsiegle

Corinne Teeter @corinneteeter

Shreejoy Tripathy @stripathy

Werner Van Geit @wvangeit

Michael Wang @aimichaelwang

Isaak Willett @isaak-willett

Derric Williams @derricw

6.1 Subpackages

6.1.1 allensdk.api package

Subpackages

allensdk.api.cloud_cache package

Submodules

allensdk.api.cloud_cache.cloud_cache module

allensdk.api.cloud_cache.file_attributes module

```
class allensdk.api.cloud_cache.file_attributes.CacheFileAttributes (url: str,  
ver-  
sion_id:  
str,  
file_hash:  
str, lo-  
cal_path:  
str)
```

Bases: `object`

This class will contain the attributes of a remotely stored file so that they can easily and consistently be passed around between the methods making up the remote file cache and manifest classes

Parameters

url: str The full URL of the remote file

version_id: str A string specifying the version of the file (probably calculated by S3)

file_hash: str The (hexadecimal) file hash of the file

local_path: pathlib.Path The path to the location where the file's local copy should be stored (probably computed by the Manifest class)

file_hash
local_path
url
version_id

allensdk.api.cloud_cache.manifest module

class allensdk.api.cloud_cache.manifest.**Manifest** (*cache_dir: Union[str, pathlib.Path], json_input*)

Bases: object

A class for loading and manipulating the online manifest.json associated with a dataset release

Each Manifest instance should represent the data for 1 and only 1 manifest.json file.

Parameters

cache_dir: str or pathlib.Path The path to the directory where local copies of files will be stored

json_input: A ".read()" -supporting file-like object containing a JSON document to be deserialized (i.e. same as the first argument to json.load)

data_file_attributes (*self, file_id*) → allensdk.api.cloud_cache.file_attributes.CacheFileAttributes
Return the CacheFileAttributes associated with a data file

Parameters

file_id: The identifier of the data file whose attributes are to be returned. Must be a key in self._data['data_files']

Raises

RuntimeError If you try to run this method when self._data is None (meaning you haven't yet loaded a manifest.json file)

ValueError If the file_id is not a valid option

file_id_column

The column in the metadata files used to uniquely identify data files

file_id_values

List of valid file_id values

metadata_file_attributes (*self, metadata_file_name: str*) → allensdk.api.cloud_cache.file_attributes.CacheFileAttributes
Return the CacheFileAttributes associated with a metadata file

Parameters

metadata_file_name: str Name of the metadata file. Must be in self.metadata_file_names

Raises

RuntimeError If you try to run this method when self._data is None (meaning you haven't yet loaded a manifest.json)

ValueError If the `metadata_file_name` is not a valid option

metadata_file_names

List of metadata file names associated with this dataset

project_name

The name of the project whose data and metadata files this manifest tracks.

version

The version of the dataset currently loaded

allensdk.api.cloud_cache.utils module

`allensdk.api.cloud_cache.utils.bucket_name_from_url(url: str) → Union[str, None-Type]`

Read in a URL and return the name of the AWS S3 bucket it points towards.

Parameters

URL: str A generic URL, suitable for retrieving an S3 object via an HTTP GET request.

Returns

str An AWS S3 bucket name. Note: if 's3.amazonaws.com' does not occur in the URL, this method will return None and emit a warning.

`allensdk.api.cloud_cache.utils.file_hash_from_path(file_path: str) → str`

Return the hexadecimal file hash for a file

Parameters

file_path: str path to a file

Returns

str: The file hash (Blake2b; hexadecimal) of the file

`allensdk.api.cloud_cache.utils.relative_path_from_url(url: str) → str`

Read in a url and return the relative path of the object

Parameters

url: str The url of the object whose path you want

Returns

str: Relative path of the object

Notes

This method returns a str rather than a `pathlib.Path` because it is used to get the S3 object Key from a URL. If using `Pathlib.path` on a Windows system, the '/' will get transformed into '\\', confusing S3.

Module contents

allensdk.api.queries package

Submodules

allensdk.api.queries.annotated_section_data_sets_api module

class allensdk.api.queries.annotated_section_data_sets_api.**AnnotatedSectionDataSetsApi** (base...
Bases: *allensdk.api.queries.rma_api.RmaApi*

See: [Searching Annotated SectionDataSets](#)

get_annotated_section_data_sets (*self*, *structures*, *intensity_values=None*, *density_values=None*, *pattern_values=None*, *age_names=None*)

For a list of target structures, find the SectionDataSet that matches the parameters for intensity_values, density_values, pattern_values, and Age.

Parameters

structure_graph_id [dict of integers] what to retrieve
intensity_values [array of strings, optional] 'High', 'Low', 'Medium' (default)
density_values [array of strings, optional] 'High', 'Low'
pattern_values [array of strings, optional] 'Full'
age_names [array of strings, options] for example 'E11.5', '13.5'

Returns

data [dict] The parsed JSON response message.

Notes

This method uses the non-RMA Annotated SectionDataSet endpoint.

get_annotated_section_data_sets_via_rma (*self*, *structures*, *intensity_values=None*, *density_values=None*, *pattern_values=None*, *age_names=None*)

For a list of target structures, find the SectionDataSet that matches the parameters for intensity_values, density_values, pattern_values, and Age.

Parameters

structure_graph_id [dict of integers] what to retrieve
intensity_values [array of strings, optional] intensity values, 'High', 'Low', 'Medium' (default)
density_values [array of strings, optional] density values, 'High', 'Low'
pattern_values [array of strings, optional] pattern values, 'Full'
age_names [array of strings, options] for example 'E11.5', '13.5'

Returns

data [dict] The parsed JSON response message.

Notes

This method uses the RMA endpoint to search annotated SectionDataSet data.

get_compound_annotated_section_data_sets (*self*, *queries*, *fmt='json'*)

Find the SectionDataSet that matches several annotated_section_data_sets queries linked together with a Boolean 'and' or 'or'.

Parameters**queries** [array of dicts] dicts with args like build_query**fmt** [string, optional] 'json' or 'xml'**Returns****data** [dict] The parsed JSON response message.**allensdk.api.queries.biophysical_api module****class** allensdk.api.queries.biophysical_api.**BiophysicalApi** (*base_uri=None*)Bases: *allensdk.api.queries.rma_template.RmaTemplate***BIOPHYSICAL_MODEL_TYPE_IDS** = (491455321, 329230710)**build_rma** (*self, neuronal_model_id, fmt='json'*)

Construct a query to find all files related to a neuronal model.

Parameters**neuronal_model_id** [integer or string representation] key of experiment to retrieve.**fmt** [string, optional] json (default) or xml**Returns****string** RMA query url.**cache_data** (*self, neuronal_model_id, working_directory=None*)

Take a an experiment id, query the Api RMA to get well-known-files download the files, and store them in the working directory.

Parameters**neuronal_model_id** [int or string representation] found in the neuronal_model table in the api**working_directory** [string] Absolute path name where the downloaded well-known files will be stored.**create_manifest** (*self, fit_path="", model_type="", stimulus_filename="", swc_morphology_path="", marker_path="", sweeps=[]*)

Generate a json configuration file with parameters for a a biophysical experiment.

Parameters**fit_path** [string] filename of a json configuration file with cell parameters.**stimulus_filename** [string] path to an NWB file with input currents.**swc_morphology_path** [string] file in SWC format.**sweeps** [array of integers] which sweeps in the stimulus file are to be used.**get_neuronal_models** (*self, specimen_ids, num_rows='all', count=False, model_type_ids=None, **kwargs*)

Fetch all of the biophysically detailed model records associated with a particular specimen_id

Parameters**specimen_ids** [list] One or more integer ids identifying specimen records.**num_rows** [int, optional] how many records to retrieve. Default is 'all'.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

model_type_ids [list, optional] One or more integer ids identifying categories of neuronal model. Defaults to all-active and perisomatic biophysical_models.

Returns

List of dict Each element is a biophysical model record, containing a unique integer id, the id of the associated specimen, and the id of the model type to which this model belongs.

get_well_known_file_ids (*self*, *neuronal_model_id*)

Query the current RMA endpoint with a neuronal_model id to get the corresponding well known file ids.

Returns

list A list of well known file id strings.

is_well_known_file_type (*self*, *wkf*, *name*)

Check if a structure has the expected name.

Parameters

wkf [dict] A well-known-file structure with nested type information.

name [string] The expected type name

See also:

[*read_json*](#) where this helper function is used.

read_json (*self*, *json_parsed_data*)

Get the list of well_known_file ids from a response body containing nested sample,microarray_slides,well_known_files.

Parameters

json_parsed_data [dict] Response from the Allen Institute Api RMA.

Returns

list of strings Well known file ids.

```
rma_templates = {'model_queries': [{'name': 'models_by_specimen', 'description': 's
```

allensdk.api.queries.brain_observatory_api module

```
class allensdk.api.queries.brain_observatory_api.BrainObservatoryApi (base_uri=None,  
                                                                    dat-  
                                                                    acube_uri=None)
```

Bases: *allensdk.api.queries.rma_template.RmaTemplate*

CELL_MAPPING_ID = 590985414

NWB_FILE_TYPE = 'NWBophys'

OPHYS_ANALYSIS_FILE_TYPE = 'OphysExperimentCellRoiMetricsFile'

OPHYS_EVENTS_FILE_TYPE = 'ObservatoryEventsFile'

dataframe_query (*self*, *data*, *filters*, *primary_key*)

Given a list of dictionary records and a list of filter dictionaries, filter the records using Pandas and return the filtered set of records.

Parameters

data: list of dicts List of dictionaries

filters: list of dicts Each dictionary describes a filtering operation on a field in the dictionary. The general form is { 'field': <field>, 'op': <operation>, 'value': <filter_value(s)> }. For example, you can apply a threshold on the "osi_dg" column with something like this: { 'field': 'osi_dg', 'op': '>', 'value': 1.0 }. See `_QUERY_TEMPLATES` for a full list of operators.

dataframe_query_string (*self, filters*)

Convert a list of cell metric filter dictionaries into a Pandas query string.

filter_cell_specimens (*self, cell_specimens, ids=None, experiment_container_ids=None, include_failed=False, filters=None*)

Filter a list of cell specimen records returned from the `get_cell_metrics` method according some of their properties.

Parameters

cell_specimens: list of dicts List of records returned by the `get_cell_metrics` method.

ids: list of integers Return only records for cells with cell specimen ids in this list

experiment_container_ids: list of integers Return only records for cells that belong to experiment container ids in this list

include_failed: bool Whether to include cells from failed experiment containers

filters: list of dicts Custom query used to reproduce filter sets created in the Allen Brain Observatory web application. The general form is a list of dictionaries each of which describes a filtering operation based on a metric. For more information, see `dataframe_query`.

filter_experiment_containers (*self, containers, ids=None, targeted_structures=None, imaging_depths=None, cre_lines=None, reporter_lines=None, transgenic_lines=None, include_failed=False, simple=False*)

filter_experiments_and_containers (*self, objs, ids=None, targeted_structures=None, imaging_depths=None, cre_lines=None, reporter_lines=None, transgenic_lines=None, include_failed=False*)

filter_ophys_experiments (*self, experiments, ids=None, experiment_container_ids=None, targeted_structures=None, imaging_depths=None, cre_lines=None, reporter_lines=None, transgenic_lines=None, stimuli=None, session_types=None, include_failed=False, require_eye_tracking=False, simple=False*)

get_cell_metrics (*self, cell_specimen_ids=None, *args, **kwargs*)

Get cell metrics by id

Parameters

cell_metrics_ids [integer or list of integers, optional] only select specific cell metric records.

Returns

dict [cell metric metadata]

get_cell_specimen_id_mapping (*self, file_name, mapping_table_id=None*)

Download mapping table from old to new cell specimen IDs.

The mapping table is a CSV file that maps cell specimen ids that have changed between processing runs of the Brain Observatory pipeline.

Parameters

file_name [string] Filename to save locally.

mapping_table_id [integer] ID of the mapping table file. Defaults to the most recent mapping table.

Returns

pandas.DataFrame Mapping table as a DataFrame.

get_column_definitions (*self*, *api_class_name=None*)

Get column definitions

Parameters

api_class_names [string or list of strings, optional] only select specific column definition records.

Returns

dict [column definition metadata]

get_experiment_container_metrics (*self*, *experiment_container_metric_ids=None*)

Get experiment container metrics by id

Parameters

isi_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [isi experiment metadata]

get_experiment_containers (*self*, *experiment_container_ids=None*)

Get experiment container by id

Parameters

experiment_container_ids [integer or list of integers, optional] only select specific experiment containers.

Returns

dict [experiment container metadata]

get_isi_experiments (*self*, *isi_experiment_ids=None*)

Get ISI Experiments by id

Parameters

isi_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [isi experiment metadata]

get_ophys_experiments (*self*, *ophys_experiment_ids=None*)

Get OPhys Experiments by id

Parameters

ophys_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [ophys experiment metadata]

get_stimulus_mappings (*self*, *stimulus_mapping_ids=None*)

Get stimulus mappings by id

Parameters

stimulus_mapping_ids [integer or list of integers, optional] only select specific stimulus mapping records.

Returns

dict [stimulus mapping metadata]

list_column_definition_class_names (*self*)

Get column definitions

Returns

list [api class name strings]

list_isi_experiments (*self*, *isi_ids=None*)

List ISI experiments available through the Allen Institute API

Parameters

neuronal_model_ids [integer or list of integers, optional] only select specific isi experiments.

Returns

dict [neuronal model metadata]

rma_templates = {'brain_observatory_queries': [{'name': 'list_isi_experiments', 'des

save_ophys_experiment_analysis_data (*self*, *ophys_experiment_id*, *file_name*)

save_ophys_experiment_data (*self*, *ophys_experiment_id*, *file_name*)

save_ophys_experiment_event_data (*self*, *ophys_experiment_id*, *file_name*)

save_ophys_experiment_eye_gaze_data (*self*, *ophys_experiment_id*: int, *ophys_session_id*: int, *file_name*: str)

simplify_experiment_containers (*self*, *containers*)

simplify_ophys_experiments (*self*, *exps*)

allensdk.api.queries.brain_observatory_api.find_container_tags (*container*)

Custom logic for extracting tags from donor conditions. Filtering out tissueocyte tags.

allensdk.api.queries.brain_observatory_api.find_experiment_acquisition_age (*exp*)

allensdk.api.queries.brain_observatory_api.find_specimen_cre_line (*specimen*)

allensdk.api.queries.brain_observatory_api.find_specimen_reporter_line (*specimen*)

allensdk.api.queries.brain_observatory_api.find_specimen_transgenic_lines (*specimen*)

allensdk.api.queries.cell_types_api module

class allensdk.api.queries.cell_types_api.**CellTypesApi** (*base_uri=None*)

Bases: *allensdk.api.queries.rma_api.RmaApi*

HUMAN = 'Homo Sapiens'

MARKER_FILE_TYPE = '3DNeuronMarker'

MOUSE = 'Mus musculus'

```
NWB_FILE_TYPE = 'NWBDownload'
```

```
SWC_FILE_TYPE = '3DNeuronReconstruction'
```

```
filter_cells (self, cells, require_morphology, require_reconstruction, reporter_status, species)
```

Filter a list of cell specimens to those that optionally have morphologies or have morphological reconstructions.

Parameters

cells: list List of cell metadata dictionaries to be filtered

require_morphology: boolean Filter out cells that have no morphological images.

require_reconstruction: boolean Filter out cells that have no morphological reconstructions.

reporter_status: list Filter for cells that have a particular cell reporter status

species: list Filter for cells that belong to one or more species. If None, return all. Must be one of [CellTypesApi.MOUSE, CellTypesApi.HUMAN].

```
filter_cells_api (self, cells, require_morphology=False, require_reconstruction=False, reporter_status=None, species=None, simple=True)
```

```
get_cell (self, id)
```

Query the API for a one cells in the Cell Types Database.

Returns

list Meta data for one cell.

```
get_ephys_features (self)
```

Query the API for the full table of EphysFeatures for all cells.

```
get_ephys_sweeps (self, specimen_id)
```

Query the API for a list of sweeps for a particular cell in the Cell Types Database.

Parameters

specimen_id: int Specimen ID of a cell.

Returns

list: List of sweep dictionaries belonging to a cell

```
get_morphology_features (self)
```

Query the API for the full table of morphology features for all cells

Notes

by default the tags column is removed because it isn't useful

```
list_cells (self, id=None, require_morphology=False, require_reconstruction=False, reporter_status=None, species=None)
```

Query the API for a list of all cells in the Cell Types Database.

Parameters

id: int ID of a cell. If not provided returns all matching cells.

require_morphology: boolean Only return cells that have morphology images.

require_reconstruction: boolean Only return cells that have morphological reconstructions.

reporter_status: list Return cells that have a particular cell reporter status.

species: list Filter for cells that belong to one or more species. If None, return all. Must be one of [CellTypesApi.MOUSE, CellTypesApi.HUMAN].

Returns

list Meta data for all cells.

list_cells_api (*self*, *id=None*, *require_morphology=False*, *require_reconstruction=False*, *reporter_status=None*, *species=None*)

save_ephys_data (*self*, *specimen_id*, *file_name*)

Save the electrophysiology recordings for a cell as an NWB file.

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the NWB file.

save_reconstruction (*self*, *specimen_id*, *file_name*)

Save the morphological reconstruction of a cell as an SWC file.

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the SWC file.

save_reconstruction_markers (*self*, *specimen_id*, *file_name*)

Save the marker file for the morphological reconstruction of a cell. These are comma-delimited files indicating points of interest in a reconstruction (truncation points, early tracing termination, etc).

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the marker file.

simplify_cells_api (*self*, *cells*)

allensdk.api.queries.connected_services module

class allensdk.api.queries.connected_services.ConnectedServices

Bases: object

A class representing a schema of informatics web services.

Notes

See [Connected Services and Pipes](#) for a human-readable list of services and parameters.

The URL format is documented at [Service Pipelines](#).

Connected Services only include API services that are accessed via the RMA endpoint using an rma::services stage.

ARRAY = 'array'

```
BOOLEAN = 'boolean'
FLOAT = 'float'
INTEGER = 'integer'
STRING = 'string'
build_url (self, service_name, kwargs)
    Create a single stage RMA url from a service name and parameters.
classmethod schema ()
    Dictionary of service names and parameters.
```

Notes

See [Connected Services and Pipes](#) for a human-readable list of connected services and their parameters.

allensdk.api.queries.glif_api module

```
class allensdk.api.queries.glif_api.GlifApi (base_uri=None)
    Bases: allensdk.api.queries.rma_template.RmaTemplate
    GLIF_TYPES = [395310498, 395310469, 395310475, 395310479, 471355161]
    NWB_FILE_TYPE = None
    cache_stimulus_file (self, output_file_name)
        DEPRECATED Download the NWB file for the current neuronal model and save it to a file.

    Parameters
        output_file_name: string File name to store the NWB file.

    get_ephys_sweeps (self)
        DEPRECATED Retrieve ephys sweep information out of downloaded metadata for a neuronal model

    Returns
        list A list of sweeps metadata dictionaries

    get_neuron_config (self, output_file_name=None)
        DEPRECATED Retrieve a model configuration file from the API, optionally save it to disk, and return the
        contents of that file as a dictionary.

    Parameters
        output_file_name: string File name to store the neuron configuration (optional).

    get_neuron_configs (self, neuronal_model_ids=None)

    get_neuronal_model (self, neuronal_model_id)
        DEPRECATED Query the current RMA endpoint with a neuronal_model id to get the corresponding well
        known files and meta data.

    Returns
        dict A dictionary containing

    get_neuronal_model_templates (self)

    get_neuronal_models (self, ephys_experiment_ids=None)

    get_neuronal_models_by_id (self, neuronal_model_ids=None)
```


`list_neuronal_models(self)`

DEPRECATED Query the API for a list of all GLIF neuronal models.

Returns

list Meta data for all GLIF neuronal models.

```
rma_templates = {'glif_queries':  [{'name':  'neuronal_model_templates', 'description':
```

allensdk.api.queries.grid_data_api module

class `allensdk.api.queries.grid_data_api.GridDataApi` (*resolution=None*,
base_uri=None)

Bases: `allensdk.api.queries.rma_api.RmaApi`

HTTP Client for the Allen 3-D Expression Grid Data Service.

See: [Downloading 3-D Expression Grid Data](#)

DATA_MASK = 'data_mask'

DENSITY = 'density'

ENERGY = 'energy'

INJECTION_DENSITY = 'injection_density'

INJECTION_ENERGY = 'injection_energy'

INJECTION_FRACTION = 'injection_fraction'

INTENSITY = 'intensity'

PROJECTION_DENSITY = 'projection_density'

PROJECTION_ENERGY = 'projection_energy'

download_alignment3d (*self*, *section_data_set_id*, *num_rows='all'*, *count=False*, ***kwargs*)

Download the parameters of the 3D affine tranformation mapping this section data set's image-space stack to CCF-space (or vice-versa).

Parameters

section_data_set_id [int] download the parameters for this data set.

Returns

dict : parameters of this section data set's alignment3d

download_deformation_field (*self*, *section_data_set_id*, *header_path=None*,
voxel_path=None, *voxel_type='DeformationFieldVoxels'*,
header_type='DeformationFieldHeader')

Download the local alignment parameters for this dataset. This a 3D vector image (3 components) describing a deformable local mapping from CCF voxels to this section data set's affine-aligned image stack.

Parameters

section_data_set_id [int]

Download the deformation field for this data set

header_path [str, optional] If supplied, the deformation field header will be downloaded to this path.

voxel_path [str, optional] If supplied, the deformation field voxels will be downloaded to this path.

voxel_type [str] WellKnownFileType of this dataset's data file

header_type [str] WellKnownFileType of this dataset's header file

download_expression_grid_data (*self, section_data_set_id, include=None, path=None*)
Download in zipped metaimage format.

Parameters

section_data_set_id [integer] What to download.

include [list of strings, optional] Image volumes. 'energy' (default), 'density', 'intensity'.

path [string, optional] File name to save as.

Returns

file [3-D expression grid data packaged into a compressed archive file (.zip).]

download_gene_expression_grid_data (*self, section_data_set_id, volume_type, path*)
Download a metaimage file containing registered gene expression grid data

Parameters

section_data_set_id [int] Download data from this experiment

volume_type [str] Download this type of data (options are GridDataApi.ENERGY, GridDataApi.DENSITY, GridDataApi.INTENSITY)

path [str] Download to this path

download_projection_grid_data (*self, section_data_set_id, image=None, resolution=None, save_file_path=None*)

Download in NRRD format.

Parameters

section_data_set_id [integer] What to download.

image [list of strings, optional] Image volume. 'projection_density', 'projection_energy', 'injection_fraction', 'injection_density', 'injection_energy', 'data_mask'.

resolution [integer, optional] in microns. 10, 25, 50, or 100 (default).

save_file_path [string, optional] File name to save as.

Notes

See [Downloading 3-D Projection Grid Data](#) for additional documentation.

allensdk.api.queries.image_download_api module

class allensdk.api.queries.image_download_api.**ImageDownloadApi** (*base_uri=None*)
Bases: *allensdk.api.queries.rma_template.RmaTemplate*

HTTP Client to download whole or partial two-dimensional images from the Allen Institute with the Section-Image, AtlasImage and ProjectionImage Download Services.

See [Downloading an Image](#) for more documentation.

```
COLORMAPS = {'aba': 8, 'aibsmmap_alt': 9, 'blue': 6, 'colormap': 10, 'expression':
```

```
atlas_image_query(self, atlas_id, image_type_name=None)
```

List atlas images belonging to a specified atlas

Parameters

atlas_id [integer, optional] Find images from this atlas.

image_type_name [string, optional] Restrict response to images of this type. If not provided, the query will get it from the atlas id.

Returns

list of dict : Each element is an AtlasImage record.

Notes

See [Downloading Atlas Images and Graphics](#) for additional documentation. `allensdk.api.queries.ontologies_api.OntologiesApi.get_atlases()` can also be used to list atlases along with their ids.

```
download_atlas_image(self, atlas_image_id, file_path=None, **kwargs)
```

```
download_image(self, image_id, file_path=None, endpoint=None, **kwargs)
```

Download whole or partial two-dimensional images from the Allen Institute with the SectionImage or AtlasImage service.

Parameters

image_id [integer] SubImage to download.

file_path [string, optional] where to put it, defaults to image_id.jpg

downsample [int, optional] Number of times to downsample the original image.

quality [int, optional] jpeg quality of the returned image, 0 to 100 (default)

expression [boolean, optional] Request the expression mask for the SectionImage.

view [string, optional] 'expression', 'projection', 'tumor_feature_annotation' or 'tumor_feature_boundary'

top [int, optional] Index of the topmost row of the region of interest.

left:int, optional Index of the leftmost column of the region of interest.

width [int, optional] Number of columns in the output image.

height [int, optional] Number of rows in the output image.

range [list of ints, optional] Filter to specify the RGB channels. low,high,low,high,low,high

colormap [list of floats, optional] Filter to specify the RGB channels. [lower_threshold,colormap] gain 0-1, colormap id is a string from ImageDownloadApi.COLORMAPS

rgb [list of floats, optional] Filter to specify the RGB channels. [red,green,blue] 0-1

contrast [list of floats, optional] Filter to specify contrast parameters. [gain,bias] 0-1

annotation [boolean, optional] Request the annotated AtlasImage

atlas [int, optional] Specify the desired Atlas' annotations.

projection [boolean, optional] Request projection for the specified image.

downsample_dimensions [boolean, optional] Indicates if the width and height should be adjusted to account for downsampling.

Returns

None the file is downloaded and saved to the path.

Notes

By default, an unfiltered full-sized image with the highest quality is returned as a download if no parameters are provided.

'downsample=1' halves the number of pixels of the original image both horizontally and vertically. `range_list = kwargs.get('range', None)`

Specifying 'downsample=2' quarters the height and width values.

Quality must be an integer from 0, for the lowest quality, up to as high as 100. If it is not specified, it defaults to the highest quality.

Top is specified in full-resolution (largest tier) pixel coordinates. `SectionImage.y` is the default value.

Left is specified in full-resolution (largest tier) pixel coordinates. `SectionImage.x` is the default value.

Width is specified in tier-resolution (desired tier) pixel coordinates. `SectionImage.width` is the default value. It is automatically adjusted when downsampled.

Height is specified in tier-resolution (desired tier) pixel coordinates. `SectionImage.height` is the default value. It is automatically adjusted when downsampled.

The range parameter consists of 6 comma delimited integers that define the lower (0) and upper (4095) bound for each channel in red-green-blue order (i.e. "range=0,1500,0,1000,0,4095"). The default range values can be determined by referring to the following fields on the Equalization model associated with the `SectionDataSet`: `red_lower`, `red_upper`, `green_lower`, `green_upper`, `blue_lower`, `blue_upper`. For more information, see the [Image Controls](#) section of the Allen Mouse Brain Connectivity Atlas: [Projection Dataset](#) help topic. See: 'Image Download Service' <<http://help.brain-map.org/display/api/Downloading+an+Image>>_

download_projection_image (*self*, *projection_image_id*, *file_path=None*, ***kwargs*)

download_section_image (*self*, *section_image_id*, *file_path=None*, ***kwargs*)

get_section_data_sets_by_product (*self*, *product_ids*, *include_failed=False*, *num_rows='all'*, *count=False*, ***kwargs*)

List all of the section data sets produced as part of one or more products

Parameters

product_ids [list of int] Integer specifiers for Allen Institute products. A product is a set of related data.

include_failed [bool, optional] If True, find both failed and passed datasets. Default is False

num_rows [int, optional] how many records to retrieve. Default is 'all'.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

Returns

list of dict : Each returned element is a section data set record.

Notes

See <http://api.brain-map.org/api/v2/data/query.json?criteria=model::Product> for a list of products.

get_section_image_ranges (*self*, *section_image_ids*, *num_rows*='all', *count*=False, *as_lists*=True, ***kwargs*)

Section images from the Mouse Connectivity Atlas are displayed on connectivity.brain-map.org after having been linearly windowed and leveled. This method obtains parameters defining channelwise upper and lower bounds of the windows used for one or more images.

Parameters

section_image_ids [list of int] Each element is a unique identifier for a section image.

num_rows [int, optional] how many records to retrieve. Default is 'all'.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

as_lists [bool, optional] If True, return the window parameters in a list, rather than a dict (this is the format of the range parameter on ImageDownloadApi.download_image). Default is False.

Returns

list of dict or list of list : For each section image id provided, return the window bounds for each channel.

rma_templates = {'image_queries': [{'name': 'section_image_ranges', 'description':

section_image_query (*self*, *section_data_set_id*, *num_rows*='all', *count*=False, ***kwargs*)

List section images belonging to a specified section data set

Parameters

atlas_id [integer, optional] Find images from this section data set.

num_rows [int] how many records to retrieve. Default is 'all'

count [bool] If True, return a count of the lines found by the query.

Returns

list of dict : Each element is an SectionImage record.

Notes

The SectionDataSet model is used to represent single experiments which produce an array of images. This includes Mouse Connectivity and Mouse Brain Atlas experiments, among other projects. You may see references to the ids of experiments from those projects. These are the same as section data set ids.

allensdk.api.queries.mouse_atlas_api module

class allensdk.api.queries.mouse_atlas_api.**MouseAtlasApi** (*base_uri*=None)

Bases: *allensdk.api.queries.reference_space_api.ReferenceSpaceApi*, *allensdk.api.queries.grid_data_api.GridDataApi*

Downloads Mouse Brain Atlas grid data, reference volumes, and metadata.

DEVMOUSE_ATLAS_PRODUCTS = (3,)

HUMAN_ORGANISM = (1,)

```
MOUSE_ATLAS_PRODUCTS = (1,)
MOUSE_ORGANISM = (2,)
download_expression_density(self, path, experiment_id)
download_expression_energy(self, path, experiment_id)
download_expression_intensity(self, path, experiment_id)
get_genes(self, organism_ids=None, chromosome_ids=None, **kwargs)
    Download a list of genes
```

Parameters

organism_ids [list of int, optional] Filter genes to those appearing in these organisms. Defaults to mouse (2).

chromosome_ids [list of int, optional] Filter genes to those appearing on these chromosomes. Defaults to all.

Returns

list of dict: Each element is a gene record, with a nested chromosome record (also a dict).

```
get_section_data_sets(self, gene_ids=None, product_ids=None, **kwargs)
    Download a list of section data sets (experiments) from the Mouse Brain Atlas project.
```

Parameters

gene_ids [list of int, optional] Filter results based on the genes whose expression was characterized in each experiment. Default is all.

product_ids [list of int, optional] Filter results to a subset of products. Default is the Mouse Brain Atlas.

Returns

list of dict : Each element is a section data set record, with one or more gene records nested in a list.

`allensdk.api.queries.mouse_connectivity_api` module

```
class allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi (base_uri=None)
    Bases: allensdk.api.queries.reference_space_api.ReferenceSpaceApi, allensdk.api.queries.grid_data_api.GridDataApi
```

HTTP Client for the Allen Mouse Brain Connectivity Atlas.

See: [Mouse Connectivity API](#)

```
PRODUCT_IDS = [5, 31]
```

```
build_reference_aligned_image_channel_volumes_url(self, data_set_id)
    Construct url to download the red, green, and blue channels aligned to the 25um adult mouse brain reference space volume.
```

Parameters

data_set_id [integerallensdk.api.queries] aka attachable_id

Notes

See: [Reference-aligned Image Channel Volumes](#) for additional documentation.

calculate_injection_centroid (*self*, *injection_density*, *injection_fraction*, *resolution=25*)
 Compute the centroid of an injection site.

Parameters

injection_density: **np.ndarray** The injection density volume of an experiment

injection_fraction: **np.ndarray** The injection fraction volume of an experiment

download_data_mask (*self*, *path*, *experiment_id*, *resolution*)

download_injection_density (*self*, *path*, *experiment_id*, *resolution*)

download_injection_fraction (*self*, *path*, *experiment_id*, *resolution*)

download_projection_density (*self*, *path*, *experiment_id*, *resolution*)

download_reference_aligned_image_channel_volumes (*self*, *data_set_id*,
save_file_path=None)

Returns

The well known file is downloaded

experiment_correlation_search (*self*, ***kwargs*)

Select a seed experiment and a domain over which the similarity comparison is to be made.

Parameters

row [integer] SectionDataSet.id to correlate against.

structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

hemisphere [string, optional] Use 'right' or 'left'. Defaults to both hemispheres.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Correlation Search](#) and [service::mouse_connectivity_correlation](#).

experiment_injection_coordinate_search (*self*, ***kwargs*)

User specifies a seed location within the 3D reference space. The service returns a rank list of experiments by distance of its injection site to the specified seed location.

Parameters

seed_point [list of floats] The coordinates of a point in 3-D SectionDataSet space.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Injection Coordinate Search](#) and `service::mouse_connectivity_injection_coordinate`.

experiment_source_search (*self*, ***kwargs*)

Search over the whole projection signal statistics dataset to find experiments with specific projection profiles.

Parameters

injection_structures [list of integers or strings] Integer Structure.id or String Structure.acronym.

target_domain [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

injection_hemisphere [string, optional] ‘right’ or ‘left’, Defaults to both hemispheres.

target_hemisphere [string, optional] ‘right’ or ‘left’, Defaults to both hemispheres.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_domain [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Source Search](#), [Target Search](#), and `service::mouse_connectivity_injection_structure`.

experiment_spatial_search (*self*, ***kwargs*)

Displays all SectionDataSets with projection signal density ≥ 0.1 at the seed point. This service also returns the path along the most dense pixels from the seed point to the center of each injection site..

Parameters

seed_point [list of floats] The coordinates of a point in 3-D SectionDataSet space.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

section_data_sets [list of integers, optional] Ids to filter the results.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Spatial Search](#) and `service::mouse_connectivity_target_spatial`.

get_experiment_detail (*self*, *experiment_id*)

Retrieve the experiments data.

get_experiments (*self*, *structure_ids*, ***kwargs*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas.

Parameters

structure_ids [integer or list, optional] injection structure

Returns

url [string] The constructed URL

get_experiments_api (*self*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas via the ApiConnectivity table.

Returns

url [string] The constructed URL

get_manual_injection_summary (*self*, *experiment_id*)

Retrieve manual injection summary.

get_projection_image_info (*self*, *experiment_id*, *section_number*)

Fetch meta-information of one projection image.

Parameters

experiment_id [integer]

section_number [integer]

Notes

See: image examples under [Experimental Overview](#) and [Metadata](#) for additional documentation. Download the image using `allensdk.api.queries.image_download_api.ImageDownloadApi.download_section_image()`

get_reference_aligned_image_channel_volumes_url (*self*, *data_set_id*)

Retrieve the download link for a specific data set. Notes — See [Reference-aligned Image Channel Volumes](#) for additional documentation.

get_structure_unionizes (*self*, *experiment_ids*, *is_injection=None*, *structure_name=None*, *structure_ids=None*, *hemisphere_ids=None*, *normalized_projection_volume_limit=None*, *include=None*, *debug=None*, *order=None*)

allensdk.api.queries.ontologies_api module

class allensdk.api.queries.ontologies_api.**OntologiesApi** (*base_uri=None*)

Bases: *allensdk.api.queries.rma_template.RmaTemplate*

See: [Atlas Drawings and Ontologies](#)

get_atlases (*self*)

get_atlases_table (*self*, *atlas_ids=None*, *brief=True*)

List Atlases available through the API with associated ontologies and structure graphs.

Parameters

atlas_ids [integer or list of integers, optional] only select specific atlases

brief [boolean, optional] True (default) requests only name and id fields.

Returns

dict [atlas metadata]

Notes

This query is based on the [table of available Atlases](#). See also: [Class: Atlas](#)

get_structure_graphs (*self*)

get_structure_sets (*self*, *structure_set_ids=None*)

get_structures (*self*, *structure_graph_ids=None*, *structure_graph_names=None*,
structure_set_ids=None, *structure_set_names=None*, *or-
der=['structures.graph_order']*, *num_rows='all'*, *count=False*, ***kwargs*)

Retrieve data about anatomical structures.

Parameters

structure_graph_ids [int or list of ints, optional] database keys to get all structures in particular graphs

structure_graph_names [string or list of strings, optional] list of graph names to narrow the query

structure_set_ids [int or list of ints, optional] database keys to get all structures in a particular set

structure_set_names [string or list of strings, optional] list of set names to narrow the query.

order [list of strings] list of RMA order clauses for sorting

num_rows [int] how many records to retrieve

Returns

dict the parsed json response containing data from the API

Notes

Only one of the methods of limiting the query should be used at a time.

get_structures_with_sets (*self*, *structure_graph_ids*, *order=['structures.graph_order']*,
num_rows='all', *count=False*, ***kwargs*)

Download structures along with the sets to which they belong.

Parameters

structure_graph_ids [int or list of int] Only fetch structure records from these graphs.

order [list of strings] list of RMA order clauses for sorting

num_rows [int] how many records to retrieve

Returns

dict the parsed json response containing data from the API

rma_templates = {'ontology_queries': [{'name': 'structures_by_graph_ids', 'descripti

unpack_structure_set_ancestors (*self*, *structure_dataframe*)

Convert a slash-separated structure_id_path field to a list.

Parameters

structure_dataframe [DataFrame] structure data from the API

Returns

None A new column is added to the dataframe containing the ancestor list.

allensdk.api.queries.reference_space_api module

class allensdk.api.queries.reference_space_api.**ReferenceSpaceApi** (*base_uri=None*)

Bases: *allensdk.api.queries.rma_api.RmaApi*

ARA_NISSL = 'ara_nissl'

AVERAGE_TEMPLATE = 'average_template'

CCF_2015 = 'annotation/ccf_2015'

CCF_2016 = 'annotation/ccf_2016'

CCF_2017 = 'annotation/ccf_2017'

CCF_VERSION_DEFAULT = 'annotation/ccf_2017'

DEVMOUSE_2012 = 'annotation/devmouse_2012'

MOUSE_2011 = 'annotation/mouse_2011'

VOXEL_RESOLUTION_100_MICRONS = 100

VOXEL_RESOLUTION_10_MICRONS = 10

VOXEL_RESOLUTION_25_MICRONS = 25

VOXEL_RESOLUTION_50_MICRONS = 50

build_volumetric_data_download_url (*self*, *data_path*, *file_name*, *voxel_resolution=None*,
release=None, *coordinate_framework=None*)

Construct url to download 3D reference model in NRRD format.

Parameters

data_path [string] 'average_template', 'ara_nissl', 'annotation/ccf_{year}', 'annotation/mouse_2011', or 'annotation/devmouse_2012'

voxel_resolution [int] 10, 25, 50 or 100

coordinate_framework [string] 'mouse_ccf' (default) or 'mouse_annotation'

Notes

See: [3-D Reference Models](#) for additional documentation.

download_annotation_volume (*self, ccf_version, resolution, file_name*)

Download the annotation volume at a particular resolution.

Parameters

ccf_version: **string** Which reference space version to download. Defaults to “annotation/ccf_2017”

resolution: **int** Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name: **string** Where to save the annotation volume.

Note: the parameters must be used as positional parameters, not keywords

download_mouse_atlas_volume (*self, age, volume_type, file_name*)

Download a reference volume (annotation, grid annotation, atlas volume) from the mouse brain atlas project

Parameters

age [str] Specify a mouse age for which to download the reference volume

volume_type [str] Specify the type of volume to download

file_name [str] Specify the path to the downloaded volume

download_structure_mask (*self, structure_id, ccf_version, resolution, file_name*)

Download an indicator mask for a specific structure.

Parameters

structure_id [int] Unique identifier for the annotated structure

ccf_version [string] Which reference space version to download. Defaults to “annotation/ccf_2017”

resolution [int] Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name [string] Where to save the downloaded mask.

download_structure_mesh (*self, structure_id, ccf_version, file_name*)

Download a Wavefront obj file containing a triangulated 3d mesh built from an annotated structure.

Parameters

structure_id [int] Unique identifier for the annotated structure

ccf_version [string] Which reference space version to download. Defaults to “annotation/ccf_2017”

file_name [string] Where to save the downloaded mask.

download_template_volume (*self, resolution, file_name*)

Download the registration template volume at a particular resolution.

Parameters

resolution: **int** Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name: **string** Where to save the registration template volume.

```
download_volumetric_data (self, data_path, file_name, voxel_resolution=None,
                           save_file_path=None, release=None, coordinate_framework=None)
```

Download 3D reference model in NRRD format.

Parameters

data_path [string] 'average_template', 'ara_nissl', 'annotation/ccf_{year}', 'annotation/mouse_2011', or 'annotation/devmouse_2012'

file_name [string] server-side file name. 'annotation_10.nrrd' for example.

voxel_resolution [int] 10, 25, 50 or 100

coordinate_framework [string] 'mouse_ccf' (default) or 'mouse_annotation'

Notes

See: [3-D Reference Models](#) for additional documentation.

allensdk.api.queries.rma_api module

```
class allensdk.api.queries.rma_api.RmaApi (base_uri=None)
```

Bases: [allensdk.api.api.Api](#)

See: [RESTful Model Access \(RMA\)](#)

ALL = 'all'

COUNT = 'count'

CRITERIA = 'rma::criteria'

DEBUG = 'debug'

EQ = '\$eq'

EXCEPT = 'except'

EXCPT = 'excpt'

FALSE = 'false'

INCLUDE = 'rma::include'

IS = '\$is'

MODEL = 'model::'

NUM_ROWS = 'num_rows'

ONLY = 'only'

OPTIONS = 'rma::options'

ORDER = 'order'

PIPE = 'pipe::'

PREVIEW = 'preview'

SERVICE = 'service::'

START_ROW = 'start_row'

TABULAR = 'tabular'

TRUE = 'true'

build_query_url (*self*, *stage_clauses*, *fmt*='json')

Combine one or more RMA query stages into a single RMA query.

Parameters

stage_clauses [list of strings] subqueries

fmt [string, optional] json (default), xml, or csv

Returns

string complete RMA url

build_schema_query (*self*, *clazz*=None, *fmt*='json')

Build the URL that will fetch the data schema.

Parameters

clazz [string, optional] Name of a specific class or None (default).

fmt [string, optional] json (default) or xml

Returns

url [string] The constructed URL

Notes

If a class is specified, only the schema information for that class will be requested, otherwise the url requests the entire schema.

debug_clause (*self*, *debug_value*=None)

Construct a debug clause for use in an rma::options clause. Parameters ——— debug_value : string or boolean

True, False, None (default) or 'preview'

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

True will request debugging information in the response. False will request no debugging information. None will return an empty clause. 'preview' will request debugging information without the query being run.

filter (*self*, *key*, *value*)

serialize a single RMA query filter clause.

Parameters

key [string] keys for narrowing a query.

value [string] value for narrowing a query.

Returns

string a single filter clause for an RMA query string.

filters (*self*, *filters*)
serialize RMA query filter clauses.

Parameters

filters [dict] keys and values for narrowing a query.

Returns

string filter clause for an RMA query string.

get_schema (*self*, *clazz=None*)
Retrieve schema information.

model_query (*self*, **args*, ***kwargs*)
Construct and execute a model stage of an RMA query string.

Parameters

model [string] The top level data type

filters [dict] key, value comparisons applied to the top-level model to narrow the results.

criteria [string] raw RMA criteria clause to choose what object are returned

include [string] raw RMA include clause to return associated objects

only [list of strings, optional] to be joined into an rma::options only filter to limit what data is returned

except [list of strings, optional] to be joined into an rma::options except filter to limit what data is returned

excpt [list of strings, optional] synonym for except parameter to avoid a reserved word conflict.

tabular [list of string, optional] return columns as a tabular data structure rather than a nested tree.

count [boolean, optional] False to skip the extra database count query.

debug [string, optional] 'true', 'false' or 'preview'

num_rows [int or string, optional] how many database rows are returned (may not correspond directly to JSON tree structure)

start_row [int or string, optional] which database row is start of returned data (may not correspond directly to JSON tree structure)

Notes

See [RMA Path Syntax](#) for a brief overview of the normalized RMA syntax. Normalized RMA syntax differs from the legacy syntax used in much of the RMA documentation. Using the `&debug=true` option with an RMA URL will include debugging information in the response, including the normalized query.

model_stage (*self*, *model*, ***kwargs*)
Construct a model stage of an RMA query string.

Parameters

model [string] The top level data type

filters [dict] key, value comparisons applied to the top-level model to narrow the results.

criteria [string] raw RMA criteria clause to choose what object are returned

include [string] raw RMA include clause to return associated objects

only [list of strings, optional] to be joined into an `rma::options` only filter to limit what data is returned

except [list of strings, optional] to be joined into an `rma::options` except filter to limit what data is returned

tabular [list of string, optional] return columns as a tabular data structure rather than a nested tree.

count [boolean, optional] False to skip the extra database count query.

debug [string, optional] 'true', 'false' or 'preview'

num_rows [int or string, optional] how many database rows are returned (may not correspond directly to JSON tree structure)

start_row [int or string, optional] which database row is start of returned data (may not correspond directly to JSON tree structure)

Notes

See [RMA Path Syntax](#) for a brief overview of the normalized RMA syntax. Normalized RMA syntax differs from the legacy syntax used in much of the RMA documentation. Using the `&debug=true` option with an RMA URL will include debugging information in the response, including the normalized query.

only_except_tabular_clause (*self*, *filter_type*, *attribute_list*)

Construct a clause to filter which attributes are returned for use in an `rma::options` clause.

Parameters

filter_type [string] 'only', 'except', or 'tabular'

attribute_list [list of strings] for example ['acronym', 'products.name', 'structure.id']

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

The title of tabular columns can be set by adding '+as+<title>' to the attribute. The tabular filter type requests a response that is row-oriented rather than a nested structure. Because of this, the tabular option can mask the lazy query behavior of an `rma::include` clause. The tabular option does not mask the inner-join behavior of an `rma::include` clause. The tabular filter is required for .csv format RMA requests.

options_clause (*self*, ***kwargs*)

build `rma::options` clause.

Parameters

only [list of strings, optional]

except [list of strings, optional]

tabular [list of string, optional]

count [boolean, optional]

debug [string, optional] 'true', 'false' or 'preview'

num_rows [int or string, optional]

start_row [int or string, optional]

order_clause (*self*, *order_list=None*)

Construct a debug clause for use in an rma::options clause.

Parameters

order_list [list of strings] for example ['acronym', 'products.name+asc', 'structure.id+desc']

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

Optionally adding '+asc' (default) or '+desc' after an attribute will change the sort order.

pipe_stage (*self*, *pipe_name*, *parameters*)

Connect model and service stages via their JSON responses.

Notes

See: [Service Pipelines](#) and [Connected Services and Pipes](#)

quote_string (*self*, *the_string*)

Wrap a clause in single quotes.

Parameters

the_string [string] a clause to be included in an rma query that needs to be quoted

Returns

string input wrapped in single quotes

service_query (*self*, **args*, ***kwargs*)

Construct and Execute a single-stage RMA query to send a request to a connected service.

Parameters

service_name [string] Name of a documented connected service.

parameters [dict] key-value pairs as in the online documentation.

Notes

See: [Service Pipelines](#) and [Connected Services and Pipes](#)

service_stage (*self*, *service_name*, *parameters=None*)

Construct an RMA query fragment to send a request to a connected service.

Parameters

service_name [string] Name of a documented connected service.

parameters [dict] key-value pairs as in the online documentation.

Notes

See: [Service Pipelines and Connected Services and Pipes](#)

tuple_filters (*self*, *filters*)
Construct an RMA filter clause.

Notes

See [RMA Path Syntax - Square Brackets for Filters](#) for additional documentation.

allensdk.api.queries.rma_pager module

```
class allensdk.api.queries.rma_pager.RmaPager
    Bases: object

    static pager (fn, *args, **kwargs)

allensdk.api.queries.rma_pager.pageable (total_rows=None, num_rows=None)
```

allensdk.api.queries.rma_template module

```
class allensdk.api.queries.rma_template.RmaTemplate (base_uri=None,
                                                    query_manifest=None)
    Bases: allensdk.api.queries.rma_api.RmaApi

    See: Atlas Drawings and Ontologies

    template_query (self, template_name, entry_name, **kwargs)

    to_filter_rhs (self, rhs)
```

allensdk.api.queries.svg_api module

```
class allensdk.api.queries.svg_api.SvgApi (base_uri=None)
    Bases: allensdk.api.api.Api

    build_query (self, section_image_id, groups=None, download=False)
        Build the URL that will fetch meta data for the specified structure.

        Parameters
            section_image_id [integer] Key of the object to be retrieved.
            groups [array of integers] Keys of the group labels to filter the svg types that are returned.

        Returns
            url [string] The constructed URL

    download_svg (self, section_image_id, groups=None, file_path=None)
        Download the svg file

    get_svg (self, section_image_id, groups=None)
        Get the svg document.
```

allensdk.api.queries.synchronization_api module

class allensdk.api.queries.synchronization_api.**SynchronizationApi** (*base_uri=None*)
 Bases: *allensdk.api.api.Api*

HTTP client for image synchronization services uses the image alignment results from the Informatics Data Processing Pipeline. Note: all locations on SectionImages are reported in pixel coordinates and all locations in 3-D ReferenceSpaces are reported in microns.

See [Image to Image Synchronization](#) for additional documentation.

get_image_to_atlas (*self, section_image_id, x, y, atlas_id*)

For a specified Atlas, find the closest annotated SectionImage and (x,y) location as defined by a seed SectionImage and seed (x,y) location.

Parameters

section_image_id [integer] Seed for spatial sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

atlas_id [int] Target Atlas for image sync.

Returns

dict The parsed json response

get_image_to_image (*self, section_image_id, x, y, section_data_set_ids*)

For a list of target SectionDataSets, find the closest SectionImage and (x,y) location as defined by a seed SectionImage and seed (x,y) pixel location.

Parameters

section_image_id [integer] Seed for spatial sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

section_data_set_ids [list of integers] Target SectionDataSet IDs for image sync.

Returns

dict The parsed json response

get_image_to_image_2d (*self, section_image_id, x, y, section_image_ids*)

For a list of target SectionImages, find the closest (x,y) location as defined by a seed SectionImage and seed (x,y) location.

Parameters

section_image_id [integer] Seed for image sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

section_image_ids [list of ints] Target SectionImage IDs for image sync.

Returns

dict The parsed json response

get_image_to_reference (*self*, *section_image_id*, *x*, *y*)

For a specified SectionImage and (x,y) location, return the (x,y,z) location in the ReferenceSpace of the associated SectionDataSet.

Parameters

section_image_id [integer] Seed for image sync.

x [float] Pixel coordinate on the specified SectionImage.

y [float] Pixel coordinate on the specified SectionImage.

Returns

dict The parsed json response

get_reference_to_image (*self*, *reference_space_id*, *x*, *y*, *z*, *section_data_set_ids*)

For a list of target SectionDataSets, find the closest SectionImage and (x,y) location as defined by a (x,y,z) location in a specified ReferenceSpace.

Parameters

reference_space_id [integer] Seed for spatial sync.

x [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

y [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

z [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

section_data_set_ids [list of ints] Target SectionDataSets IDs for image sync.

Returns

dict The parsed json response

get_structure_to_image (*self*, *section_data_set_id*, *structure_ids*)

For a list of target structures, find the closest SectionImage and (x,y) location as defined by the centroid of each Structure.

Parameters

section_data_set_id [integer] primary key

structure_ids [list of integers] primary key

Returns

dict The parsed json response

allensdk.api.queries.tree_search_api module

class allensdk.api.queries.tree_search_api.**TreeSearchApi** (*base_uri=None*)

Bases: [allensdk.api.api.Api](#)

See [Searching a Specimen or Structure Tree](#) for additional documentation.

get_tree (*self*, *kind*, *db_id*, *ancestors=None*, *descendants=None*)

Fetch meta data for the specified structure or specimen.

Parameters

kind [string] 'Structure' or 'Specimen'

db_id [integer] The id of the structure or specimen to search.

ancestors [boolean, optional] whether to include ancestors in the response (defaults to False)

descendants [boolean, optional] whether to include descendants in the response (defaults to False)

Returns

dict parsed json response data

Module contents

allensdk.api.warehouse_cache package

Submodules

allensdk.api.warehouse_cache.cache module

```
class allensdk.api.warehouse_cache.cache.Cache (manifest=None, cache=True, version=None, **kwargs)
```

Bases: object

add_manifest_paths (*self, manifest_builder*)

Add cache-class specific paths to the manifest. In derived classes, should call super.

build_manifest (*self, file_name*)

Creation of default path specifications.

Parameters

file_name [string] where to save it

static **cache_csv** ()

static **cache_csv_dataframe** ()

static **cache_csv_json** ()

static **cache_json** ()

static **cache_json_dataframe** ()

static **cacher** (*fn, *args, **kwargs*)

make an rma query, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

strategy [string or None, optional] ‘create’ always generates the data, ‘file’ loads from disk, ‘lazy’ queries the server if no file exists, None generates the data and bypasses all caching behavior

pre [function] dfjson->dfjson, takes one data argument and returns filtered version, None for pass-through

post [function] dfjson->?, takes one data argument and returns Object

reader [function, optional] path -> data, default NOP

writer [function, optional] path, data -> None, default NOP

kwargs [objects] passed through to the query function

Returns

Object or None data type depends on fn, reader and/or post methods.

static csv_writer (*pth, gen*)

get_cache_path (*self, file_name, manifest_key, *args*)

Helper method for accessing path specs from manifest keys.

Parameters

file_name [string]

manifest_key [string]

args [ordered parameters]

Returns

string or None path

static json_remove_keys (*data, keys*)

static json_rename_columns (*data, new_old_name_tuples=None*)

Convenience method to rename columns in a pandas dataframe.

Parameters

data [dataframe] edited in place.

new_old_name_tuples [list of string tuples (new, old)]

load_csv (*self, path, rename=None, index=None*)

Read a csv file as a pandas dataframe.

Parameters

rename [list of string tuples (new old), optional] columns to rename

index [string, optional] post-rename column to use as the row label.

load_json (*self, path, rename=None, index=None*)

Read a json file as a pandas dataframe.

Parameters

rename [list of string tuples (new old), optional] columns to rename

index [string, optional] post-rename column to use as the row label.

load_manifest (*self, file_name, version=None*)

Read a keyed collection of path specifications.

Parameters

file_name [string] path to the manifest file

Returns

Manifest

manifest_dataframe (*self*)

Convenience method to view manifest as a pandas dataframe.

static nocache_dataframe ()

static nocache_json ()

static pathfinder (*file_name_position*, *secondary_file_name_position=None*,
path_keyword=None)
 helper method to find path argument in legacy methods written prior to the @cacheable decorator. Do not use for new @cacheable methods.

Parameters

file_name_position [integer] zero indexed position in the decorated method args where file path may be found.

secondary_file_name_position [integer] zero indexed position in the decorated method args where the file path may be found.

path_keyword [string] kwarg that may have the file path.

Notes

This method is only intended to provide backward-compatibility for some methods that otherwise do not follow the path conventions of the @cacheable decorator.

static remove_keys (*data*, *keys=None*)
 DataFrame version

static rename_columns (*data*, *new_old_name_tuples=None*)
 Convenience method to rename columns in a pandas dataframe.

Parameters

data [dataframe] edited in place.

new_old_name_tuples [list of string tuples (new, old)]

wrap (*self*, *fn*, *path*, *cache*, *save_as_json=True*, *return_dataframe=False*, *index=None*, *rename=None*,
***kwargs*)
 make an rma query, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

cache [boolean] True will make the query, False just loads from disk

save_as_json [boolean, optional] True (default) will save data as json, False as csv

return_dataframe [boolean, optional] True will cast the return value to a pandas dataframe,
 False (default) will not

index [string, optional] column to use as the pandas index

rename [list of string tuples, optional] (new, old) columns to rename

kwargs [objects] passed through to the query function

Returns

dict or DataFrame data type depends on return_dataframe option.

Notes

Column renaming happens after the file is reloaded for json

```
allensdk.api.warehouse_cache.cache.cacheable (strategy=None, pre=None,
                                                writer=None, reader=None, post=None,
                                                pathfinder=None)
```

decorator for rma queries, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

strategy [string or None, optional] 'create' always gets the data from the source (server or generated), 'file' loads from disk, 'lazy' creates the data and saves to file if no file exists, None queries the server and bypasses all caching behavior

pre [function] dfjson->dfjson, takes one data argument and returns filtered version, None for pass-through

post [function] dfjson->?, takes one data argument and returns Object

reader [function, optional] path -> data, default NOP

writer [function, optional] path, data -> None, default NOP

kwargs [objects] passed through to the query function

Returns

dict or DataFrame data type depends on dataframe option.

Notes

Column renaming happens after the file is reloaded for json

```
allensdk.api.warehouse_cache.cache.get_default_manifest_file (cache_name)
```

```
allensdk.api.warehouse_cache.cache.memoize (f)
```

Creates an unbound cache of function calls and results. Note that arguments of different types are not cached separately (so `f(3.0)` and `f(3)` are not treated as distinct calls)

Arguments to the cached function must be hashable.

View the cache size with `f.cache_size()`. Clear the cache with `f.cache_clear()`. Access the underlying function with `f.__wrapped__`.

allensdk.api.warehouse_cache.caching_utilities module

```
allensdk.api.warehouse_cache.caching_utilities.call_caching(fetch: Callable[[  
    ~Q], ~Q], write: Callable[[~Q],  
    NoneType], read: Union[Callable[[  
    ~P], NoneType] = None, pre_write:  
    Union[Callable[[~Q], ~Q], NoneType] =  
    None, cleanup: Union[Callable[[  
    NoneType], NoneType] = None,  
    lazy: bool = True, num_tries: int = 1,  
    failure_message: str = ") → Union[~P,  
    NoneType]
```

Access data, caching on a local store for future accesses.

Parameters

fetch : Function which pulls data from a remote/expensive source.

write : Function which stores data in a local/inexpensive store.

read : Function which pulls data from a local/inexpensive store.

pre_write : Function applied to obtained data after fetching, but before writing.

cleanup : Function for fixing a failed fetch. e.g. unlinking a partially downloaded file. Exceptions raised by cleanup are not themselves handled

lazy : If True, attempt to read the data from the local/inexpensive store before fetching it. If False, forcibly fetch from the remote/expensive store.

num_tries : How many fetches to attempt before (re)raising an exception. A fetch is failed if reading the result raises an exception.

failure_message : Provides additional context in the event of a failed download. Emitted when retrying, and when a fetch failure occurs after tries are exhausted

Returns

The result of calling read

```
allensdk.api.warehouse_cache.caching_utilities.one_file_call_caching(path:
                                                                    Union[pathlib.Path,
                                                                    str],
                                                                    fetch:
                                                                    Callable[[],
                                                                    ~Q],
                                                                    write:
                                                                    Callable[[Union[pathlib.Path,
                                                                    str],
                                                                    ~Q],
                                                                    None-
                                                                    Type],
                                                                    read:
                                                                    Union[Callable[[Union[pathlib
                                                                    str]],
                                                                    ~P],
                                                                    None-
                                                                    Type]
                                                                    =
                                                                    None,
                                                                    pre_write:
                                                                    Union[Callable[[~Q],
                                                                    ~Q],
                                                                    None-
                                                                    Type]
                                                                    =
                                                                    None,
                                                                    cleanup:
                                                                    Union[Callable[[],
                                                                    None-
                                                                    Type],
                                                                    None-
                                                                    Type]
                                                                    =
                                                                    None,
                                                                    lazy:
                                                                    bool =
                                                                    True,
                                                                    num_tries:
                                                                    int =
                                                                    1, failure_message:
                                                                    str =
                                                                    ") →
                                                                    Union[~P,
                                                                    None-
                                                                    Type]
```

A call_caching variant where the local store is a single file. See call_caching for complete documentation.

Parameters

path : Path at which the data will be stored

Module contents

Submodules

allensdk.api.api module

class allensdk.api.api.**Api** (*api_base_url_string=None*)

Bases: object

cleanup_truncated_file (*self, file_path*)

Helper for removing files.

Parameters

file_path [string] Absolute path including the file name to remove.

construct_well_known_file_download_url (*self, well_known_file_id*)

Join data api endpoint and id.

Parameters

well_known_file_id [integer or string representing an integer] well known file id

Returns

string the well-known-file download url for the current api server

See also:

[*retrieve_file_over_http*](#) Can be used to retrieve the file from the url.

default_api_url = 'http://api.brain-map.org'

do_query (*self, url_builder_fn, json_traversal_fn, *args, **kwargs*)

Bundle an query url construction function with a corresponding response json traversal function.

Parameters

url_builder_fn [function] A function that takes parameters and returns an rma url.

json_traversal_fn [function] A function that takes a json-parsed python data structure and returns data from it.

post [boolean, optional kwarg] True does an HTTP POST, False (default) does a GET

args [arguments] Arguments to be passed to the url builder function.

kwargs [keyword arguments] Keyword arguments to be passed to the rma builder function.

Returns

any type The data extracted from the json response.

Examples

A simple Api subclass example.

do_rma_query (*self, rma_builder_fn, json_traversal_fn, *args, **kwargs*)

Bundle an RMA query url construction function with a corresponding response json traversal function.

..note:: **Deprecated in AllenSDK 0.9.2** *do_rma_query* will be removed in AllenSDK 1.0, it is replaced by *do_query* because the latter is more general.

Parameters

- rma_builder_fn** [function] A function that takes parameters and returns an rma url.
- json_traversal_fn** [function] A function that takes a json-parsed python data structure and returns data from it.
- args** [arguments] Arguments to be passed to the rma builder function.
- kwargs** [keyword arguments] Keyword arguments to be passed to the rma builder function.

Returns

- any type** The data extracted from the json response.

Examples

A simple `Api` subclass example.

```
download_url = 'http://download.alleninstitute.org'
```

```
json_msg_query (self, url, dataframe=False)
```

Common case where the url is fully constructed and the response data is stored in the 'msg' field.

Parameters

- url** [string] Where to get the data in json form
- dataframe** [boolean] True converts to a pandas dataframe, False (default) doesn't

Returns

- dict or DataFrame** returned data; type depends on dataframe option

```
load_api_schema (self)
```

Download the RMA schema from the current RMA endpoint

Returns

- dict** the parsed json schema message

Notes

This information and other [Allen Brain Atlas Data Portal Data Model](#) documentation is also available as a [Class Hierarchy](#) and [Class List](#).

```
read_data (self, parsed_json)
```

Return the message data from the parsed query.

Parameters

- parsed_json** [dict] A python structure corresponding to the JSON data returned from the API.

Notes

See [API Response Formats - Response Envelope](#) for additional documentation.

```
retrieve_file_over_http (self, url, file_path, zipped=False)
```

Get a file from the data api and save it.

Parameters

url [string] Url[R099781a1d33c-1]_ from which to get the file.

file_path [string] Absolute path including the file name to save.

zipped [bool, optional] If true, assume that the response is a zipped directory and attempt to extract contained files into the directory containing file_path. Default is False.

See also:

[`construct_well_known_file_download_url`](#) Can be used to construct the url.

References

[1]

retrieve_parsed_json_over_http (*self, url, post=False*)

Get the document and put it in a Python data structure

Parameters

url [string] Full API query url.

post [boolean] True does an HTTP POST, False (default) encodes the URL and does a GET

Returns

dict Result document as parsed by the JSON library.

retrieve_xml_over_http (*self, url*)

Get the document and put it in a Python data structure

Parameters

url [string] Full API query url.

Returns

string Unparsed xml string.

set_api_urls (*self, api_base_url_string*)

Set the internal RMA and well known file download endpoint urls based on a api server endpoint.

Parameters

api_base_url_string [string] url of the api to point to

set_default_working_directory (*self, working_directory*)

Set the working directory where files will be saved.

Parameters

working_directory [string] the absolute path string of the working directory.

allensdk.api.api.stream_file_over_http (*url, file_path, timeout=(9.05, 31.1)*)

Supply an http get request and stream the response to a file.

Parameters

url [str] Send the request to this url

file_path [str] Stream the response to this path

timeout [float or tuple of float, optional] Specify a timeout for the request. If a tuple, specify separate connect and read timeouts.

`allensdk.api.api.stream_zip_directory_over_http` (*url*, *directory*, *members=None*, *time-out=(9.05, 31.1)*)

Supply an http get request and stream the response to a file.

Parameters

url [str] Send the request to this url

directory [str] Extract the response to this directory

members [list of str, optional] Extract only these files

timeout [float or tuple of float, optional] Specify a timeout for the request. If a tuple, specify separate connect and read timeouts.

Module contents

Subclasses of `allensdk.api.api.Api` to implement specific queries to the [Allen Brain Atlas Data Portal](#).

6.1.2 allensdk.brain_observatory package

Subpackages

`allensdk.brain_observatory.behavior` package

Subpackages

`allensdk.brain_observatory.behavior.behavior_project_cache` package

Subpackages

`allensdk.brain_observatory.behavior.behavior_project_cache.external` package

Submodules

`allensdk.brain_observatory.behavior.behavior_project_cache.external.behavior_project_metadata_writer` module

Module contents

`allensdk.brain_observatory.behavior.behavior_project_cache.tables` package

Subpackages

`allensdk.brain_observatory.behavior.behavior_project_cache.tables.util` package

Submodules

`allensdk.brain_observatory.behavior.behavior_project_cache.tables.util.prior_exposure_processing` module

Module contents

Submodules

`allensdk.brain_observatory.behavior.behavior_project_cache.tables.experiments_table` module

class `allensdk.brain_observatory.behavior.behavior_project_cache.tables.experiments_table.ExperimentsTable`

Bases: `allensdk.brain_observatory.behavior.behavior_project_cache.tables.project_table.ProjectTable`, `allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_mixin.OphysMixin`

Class for storing and manipulating project-level data at the behavior-ophys experiment level

postprocess_additional (*self*)
Additional postprocessing should be overridden by subclassess

`allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_mixin` module

class `allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_mixin.OphysMixin`
Bases: `object`

A mixin class for ophys project data

`allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_sessions_table` module

class `allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_sessions_table.OphysSessionsTable`

Bases: `allensdk.brain_observatory.behavior.behavior_project_cache.tables.project_table.ProjectTable`, `allensdk.brain_observatory.behavior.behavior_project_cache.tables.ophys_mixin.OphysMixin`

Class for storing and manipulating project-level data at the behavior-ophys session level

postprocess_additional (*self*)
Additional postprocessing should be overridden by subclassess

allensdk.brain_observatory.behavior.behavior_project_cache.tables.project_table module

class allensdk.brain_observatory.behavior.behavior_project_cache.tables.project_table.**Proj**

Bases: abc.ABC
Class for storing and manipulating project-level data
postprocess (*self*)
Postprocess loop
postprocess_additional (*self*)
Additional postprocessing should be overridden by subclassess
postprocess_base (*self*)
Postprocessing to apply to all project-level data
table

allensdk.brain_observatory.behavior.behavior_project_cache.tables.sessions_table module

Module contents

Submodules

allensdk.brain_observatory.behavior.behavior_project_cache.behavior_project_cache module

Module contents

allensdk.brain_observatory.behavior.metadata package

Submodules

allensdk.brain_observatory.behavior.metadata.behavior_metadata module

class allensdk.brain_observatory.behavior.metadata.behavior_metadata.**BehaviorMetadata** (extractor, stimulus, numpy, behavior_stimulus)

Bases: object

Container class for behavior metadata

age_in_days

Converts the age cod into a numeric days representation

behavior_session_id

behavior_session_uuid

Get the universally unique identifier (UUID)

cre_line

Parses cre_line from full_genotype

date_of_acquisition

Return the timestamp for when experiment was started in UTC

NOTE: This method will only get acquisition datetime from extractor (data from LIMS) methods. As a sanity check, it will also read the acquisition datetime from the behavior stimulus (.pkl) file and raise a warning if the date differs too much from the datetime obtained from the behavior stimulus (.pkl) file.

Return type datetime

driver_line

equipment_name

full_genotype

get_extractor (self)

indicator

Parses indicator from reporter

mouse_id

static parse_age_in_days (age: str, warn=False) → Union[int, NoneType]

Converts the age code into a numeric days representation

Parameters

age age code, ie P123

warn Whether to output warning if parsing fails

static parse_cre_line (full_genotype: str, warn=False) → Union[str, NoneType]

Parameters

full_genotype formatted from LIMS, e.g. Vip-IRES-Cre/wt;Ai148(TIT2L-GC6f-ICL-tTA2)/wt

warn Whether to output warning if parsing fails

Returns

cre_line just the Cre line, e.g. Vip-IRES-Cre, or None if not possible to parse

static parse_indicator (*reporter_line*: Union[str, NoneType], *warn*=False) → Union[str, NoneType]

Parses indicator from reporter

static parse_reporter_line (*reporter_line*: Union[List[str], NoneType], *warn*=False) → Union[str, NoneType]

There can be multiple reporter lines, so it is returned from LIMS as a list. But there shouldn't be more than 1 for behavior. This tries to convert to str

Parameters

reporter_line List of reporter line

warn Whether to output warnings if parsing fails

Returns

single reporter line, or None if not possible

reporter_line

session_type

sex

stimulus_frame_rate

to_dict (*self*) → dict

Returns dict representation of all properties in class

`allensdk.brain_observatory.behavior.metadata.behavior_metadata.get_expt_description` (*session_type*: str) → str

Determine a behavior ophys session's experiment description based on session type. Matches the regex patterns defined as the keys in `description_dict`

Parameters

session_type [str] A session description string (e.g. OPHYS_1_images_B)

Returns

str A description of the experiment based on the session_type.

Raises

RuntimeError Behavior ophys sessions should only have 6 different session types. Unknown session types (or malformed session_type strings) will raise an error.

`allensdk.brain_observatory.behavior.metadata.behavior_metadata.get_task_parameters` (*data*: Dict) → Dict

Read task_parameters metadata from the behavior stimulus pickle file.

Parameters

data: dict The nested dict read in from the behavior stimulus pickle file. All of the data expected by this method lives under data['items']['behavior']

Returns

dict A dict containing the task_parameters associated with this session.

allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata module

class allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata.**BehaviorOphysMet**

Bases: `allensdk.brain_observatory.behavior.metadata.behavior_metadata.BehaviorMetadata`

Container class for behavior ophys metadata

emission_lambda

excitation_lambda

experiment_container_id

field_of_view_height

field_of_view_width

imaging_depth

imaging_plane_group

imaging_plane_group_count

ophys_experiment_id

ophys_frame_rate

ophys_session_id

project_code

targeted_structure

to_dict (*self*) → dict

Returns dict representation of all properties in class

Module contents

allensdk.brain_observatory.behavior.session_apis package

Subpackages

`allensdk.brain_observatory.behavior.session_apis.abcs` package

Subpackages

`allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base` package

Submodules

`allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_data_extractor_base` module

class `allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_data_extractor_base`
Bases: `abc.ABC`

Abstract base class implementing required methods for extracting data (from LIMS or from JSON) that will be transformed or passed on to fill behavior session data.

get_age (*self*) → str
Get the age code of the subject (ie P123)

get_behavior_session_id (*self*) → int
Get the ID of the behavior session

get_behavior_stimulus_file (*self*) → str
Get the filepath to the StimulusPickle file for the session

get_date_of_acquisition (*self*) → `datetime.datetime`
Get the acquisition date of an experiment in UTC

get_driver_line (*self*) → List[str]
Get the (gene) driver line(s) for the subject associated with a behavior or behavior + ophys experiment

get_equipment_name (*self*) → str
Get the name of the experiment rig (ex: CAM2P.3)

get_foraging_id (*self*) → int
Get the foraging ID for the behavior session

get_full_genotype (*self*) → str
Get the full genotype of the subject associated with a behavior or behavior + ophys experiment

get_mouse_id (*self*) → int
Get the mouse id (LabTracks ID) for the subject associated with a behavior experiment

get_reporter_line (*self*) → List[str]
Get the (gene) reporter line(s) for the subject associated with a behavior or behavior + ophys experiment

get_sex (*self*) → str
Get the sex of the subject (ex: 'M', 'F', or 'unknown')

get_stimulus_name (*self*) → str
Get the name of the stimulus presented for a behavior or behavior + ophys experiment

allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base module

class allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base

Bases: `allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base.BehaviorDataExtractorBase`

Abstract base class implementing required methods for extracting data (from LIMS or from JSON) that will be transformed or passed on to fill behavior + ophys session data.

get_average_intensity_projection_image_file (*self*) → str
Get the avg intensity project image filepath associated with an ophys experiment

get_demix_file (*self*) → str
Get the filepath of the demixed traces file associated with an ophys experiment

get_dff_file (*self*) → str
Get the filepath of the dff trace file associated with an ophys experiment.

get_event_detection_filepath (*self*) → str
Get the filepath of the .h5 events file associated with an ophys experiment

get_eye_tracking_filepath (*self*) → dict
Get the eye tracking filepath containing ellipse fits

get_eye_tracking_rig_geometry (*self*) → dict
Get the eye tracking rig geometry associated with an ophys experiment

get_field_of_view_shape (*self*) → Dict[str, int]
Get a field of view dictionary for a given ophys experiment. ex: {"width": int, "height": int}

get_imaging_depth (*self*) → int
Get the imaging depth for an ophys experiment (ex: 400, 500, etc.)

get_imaging_plane_group (*self*) → Union[int, NoneType]
Get the imaging plane group number. This is a numeric index that indicates the order that the frames were acquired when there is more than one frame acquired concurrently. Relevant for mesoscope data timestamps, as the laser jumps between plane groups during the scan. Will be None for non-mesoscope data.

get_max_projection_file (*self*) → str
Get the filepath of the max projection image associated with the ophys experiment

get_ophys_cell_segmentation_run_id (*self*) → int
Get the ophys cell segmentation run id associated with an ophys experiment id

get_ophys_container_id (*self*) → int
Get the experiment container id associated with an ophys experiment

get_ophys_experiment_id (*self*) → int
Return the ophys experiment id (experiments are an internal alias for an imaging plane)

get_ophys_session_id (*self*) → int
Return the ophys session id

get_plane_group_count (*self*) → int
Gets the total number of plane groups in the session. This is required for resampling ophys timestamps for mesoscope data. Will be 0 if the scope did not capture multiple concurrent frames. See `get_imaging_plane_group` for more info.

get_project_code (*self*) → str
Get the project code.

get_raw_cell_specimen_table_dict (*self*) → dict
Get the cell_rois table from LIMS in dictionary form

get_rigid_motion_transform_file (*self*) → str
Get the filepath for the motion transform file (.csv) associated with an ophys experiment

get_surface_2p_pixel_size_um (*self*) → float
Get the pixel size for 2-photon movies in micrometers

get_sync_file (*self*) → str
Get the filepath of the sync timing file associated with the ophys experiment

get_targeted_structure (*self*) → str
Get the targeted structure (acronym) for an ophys experiment (ex: “Visp”)

Module contents

allensdk.brain_observatory.behavior.session_apis.abcs.session_base package

Submodules

allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_base module

class allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_base.**BehaviorBase**
Bases: abc.ABC

Abstract base class implementing required methods for interacting with behavior session data.

Child classes should be instantiated with a fetch API that implements these methods.

get_behavior_session_id (*self*) → int
Returns the behavior_session_id associated with this experiment, if applicable.

get_licks (*self*) → pandas.core.frame.DataFrame
Get lick data from pkl file.

Returns

pd.DataFrame A dataframe containing lick timestamps.

get_metadata (*self*) → Union[allensdk.brain_observatory.behavior.metadata.behavior_metadata.BehaviorMetadata, dict]
Get metadata for Session

Returns

dict if NWB

BehaviorMetadata otherwise

get_rewards (*self*) → pandas.core.frame.DataFrame
Get reward data from pkl file.

Returns

pd.DataFrame A dataframe containing timestamps of delivered rewards.

get_running_acquisition_df (*self*) → pandas.core.frame.DataFrame
Get running speed acquisition data from a behavior pickle file.

Returns

pd.DataFrame

Dataframe with an index of timestamps and the following columns: “speed”: computed running speed “dx”: angular change, computed during data collection “v_sig”: voltage signal from the encoder “v_in”: the theoretical maximum voltage that the encoder

will reach prior to “wrapping”. This should theoretically be 5V (after crossing 5V goes to 0V, or vice versa). In practice the encoder does not always reach this value before wrapping, which can cause transient spikes in speed at the voltage “wraps”.

get_running_speed (*self*) → pandas.core.frame.DataFrame

Get running speed using timestamps from self.get_stimulus_timestamps.

NOTE: Do not correct for monitor delay.

Returns**pd.DataFrame**

timestamps [np.ndarray] index consisting of timestamps of running speed data samples

speed [np.ndarray] Running speed of the experimental subject (in cm / s).

get_stimulus_presentations (*self*) → pandas.core.frame.DataFrame

Get stimulus presentation data.

NOTE: Uses timestamps that do not account for monitor delay.

Returns

pd.DataFrame Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_templates (*self*) → allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate

Get stimulus templates (movies, scenes) for behavior session.

Returns**StimulusTemplate**

get_stimulus_timestamps (*self*) → numpy.ndarray

Get stimulus timestamps from.pkl file.

NOTE: Located with behavior_session_id

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor that do no account for monitor delay.

get_task_parameters (*self*) → dict

Get task parameters from.pkl file.

Returns

dict A dictionary containing parameters used to define the task runtime behavior.

get_trials (*self*) → pandas.core.frame.DataFrame

Get trials from.pkl file

Returns

pd.DataFrame A dataframe containing behavioral trial start/stop times, and trial data

allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_ophys_base module

class allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_ophys_base

Bases: `allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_base.BehaviorBase`

Abstract base class implementing required methods for interacting with behavior + ophys session data.

Child classes should be instantiated with a fetch API that implements these methods.

get_average_projection (*self*) → allensdk.brain_observatory.behavior.image_api.Image

Get an image whose values are the average obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to avg projection image data and metadata.

get_cell_specimen_table (*self*) → pandas.core.frame.DataFrame

Get a cell specimen dataframe containing ROI information about cells identified in an ophys experiment.

Returns

pd.DataFrame Cell ROI information organized into a dataframe. Index is the cell ROI IDs.

get_corrected_fluorescence_traces (*self*) → pandas.core.frame.DataFrame

Get motion-corrected fluorescence traces.

Returns

pd.DataFrame Motion-corrected fluorescence traces organized into a dataframe. Index is the cell ROI IDs.

get_dff_traces (*self*) → pandas.core.frame.DataFrame

Get a table of delta fluorescence over fluorescence traces.

Returns

pd.DataFrame The traces of dff (normalized fluorescence) organized into a dataframe. Index is the cell ROI IDs.

get_events (*self*) → pandas.core.frame.DataFrame

Get event detection data

Returns

pd.DataFrame

index: cell_specimen_id: int

cell_roi_id: int events: np.array filtered_events: np.array lambdas: float64 noise_std: float64

get_eye_tracking (*self*) → Union[pandas.core.frame.DataFrame, NoneType]

Get eye tracking data from behavior + ophys session.

Returns

pd.DataFrame A refined eye tracking dataframe that contains information about eye tracking ellipse fits, frame times, eye areas, pupil areas, and frames with likely blinks/outliers.

get_eye_tracking_rig_geometry (*self*) → dict

Get eye tracking rig metadata from behavior + ophys session.

Returns

dict Includes geometry of monitor, camera, LED

get_max_projection (*self*) → allensdk.brain_observatory.behavior.image_api.Image

Get an image whose values are the maximum obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to max projection image data and metadata.

get_metadata (*self*) → Union[allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata.BehaviorOphysMeta
dict]

Get behavior+ophys session metadata.

Returns

dict if NWB

BehaviorOphysMetadata otherwise

get_motion_correction (*self*) → pandas.core.frame.DataFrame

Get motion correction trace data.

Returns

pd.DataFrame A dataframe containing trace data used during motion correction computation.

get_ophys_experiment_id (*self*) → Union[int, NoneType]

Returns the ophys_experiment_id for the instantiated BehaviorOphys Session (or BehaviorOphys data fetcher) if applicable.

get_ophys_session_id (*self*) → Union[int, NoneType]

Returns the behavior + ophys_session_id associated with this experiment, if applicable.

get_ophys_timestamps (*self*) → numpy.ndarray

Get optical physiology frame timestamps.

Returns

np.ndarray Timestamps associated with frames captured by the microscope.

get_stimulus_presentations (*self*) → pandas.core.frame.DataFrame

Get stimulus presentation data.

NOTE: Uses monitor delay corrected stimulus timestamps.

Returns

pd.DataFrame Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_timestamps (*self*) → numpy.ndarray

Get stimulus timestamps.

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor after accounting for monitor delay.

Module contents

Module contents

allensdk.brain_observatory.behavior.session_apis.data_io package

Submodules

allensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api module

class allensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.**BehaviorJSONAPI**

Bases: *allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_transforms.BehaviorDataTransforms*

A data fetching and processing class that serves processed data from a specified raw data source (extractor). Contains all methods needed to fill a BehaviorSession.

class allensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.**BehaviorJSONExtractor**

Bases: *allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_data_extractor_base.BehaviorDataExtractorBase*

A class which ‘extracts’ data from a json file. The extracted data is necessary (but not sufficient) for populating a ‘BehaviorSession’.

Most data provided by this extractor needs to be processed by BehaviorDataTransforms methods in order to be usable by ‘BehaviorSession’s.

This class is used by the write_nwb module for behavior sessions.

get_age (*self*) → str

Get the age code of the subject (ie P123)

get_behavior_session_id (*self*) → int

Get the ID of the behavior session

get_behavior_stimulus_file (*self*) → str

Get the filepath to the StimulusPickle file for the session

get_date_of_acquisition (*self*) → datetime.datetime

Get the acquisition date of an experiment (in UTC)

NOTE: LIMS writes to JSON in local time. Needs to be converted to UTC

get_driver_line (*self*) → str

Get the (gene) driver line for the subject associated with an experiment

get_equipment_name (*self*) → str

Get the name of the experiment rig (ex: CAM2P.3)

get_foraging_id (*self*) → int

Get the foraging ID for the behavior session

get_full_genotype (*self*) → str

Get the full genotype of the subject associated with an experiment

get_mouse_id (*self*) → int

Get the external specimen id (LabTracks ID) for the subject associated with a behavior experiment

get_reporter_line (*self*) → str

Get the (gene) reporter line for the subject associated with an experiment

get_sex(*self*) → str
Get the sex of the subject (ex: 'M', 'F', or 'unknown')

get_stimulus_name(*self*) → str
Get the name of the stimulus presented for the experiment

allensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api module

class allensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api.**BehaviorL**

Bases: *allensdk.brain_observatory.behavior.session_apis.data_transforms.
behavior_data_transforms.BehaviorDataTransforms, allensdk.core.
cache_method_utilities.CachedInstanceMethodMixin*

A data fetching and processing class that serves processed data from a specified raw data source (extractor).
Contains all methods needed to fill a BehaviorSession.

class allensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api.**BehaviorL**

Bases: *allensdk.brain_observatory.behavior.session_apis.abcs.
data_extractor_base.behavior_data_extractor_base.BehaviorDataExtractorBase*

A data fetching class that serves as an API for fetching 'raw' data from LIMS necessary (but not sufficient) for filling a 'BehaviorSession'.

Most ‘raw’ data provided by this API needs to be processed by BehaviorDataTransforms methods in order to usable by ‘BehaviorSession’s

```
classmethod from_foraging_id (foraging_id: Union[str, uuid.UUID, int], lims_credentials:
                                Union[allensdk.core.authentication.DbCredentials, None-
                                Type] = None) → 'BehaviorLimsApi'
```

Create a BehaviorLimsAPI instance from a foraging_id instead of a behavior_session_id.

NOTE: ‘foraging_id’ in the LIMS behavior_session table should be the same as the ‘behavior_session_uuid’ in mtrain which should also be the same as the ‘session_uuid’ field in the .pkl returned by ‘get_behavior_stimulus_file()’.

```
get_age (self) → str
```

Return the age code of the subject (ie P123) :rtype: str

```
get_behavior_session_id (self) → int
```

Getter to be consistent with BehaviorOphysLimsApi.

```
get_behavior_stimulus_file (self) → str
```

Return the path to the StimulusPickle file for a session. :rtype: str

```
get_birth_date (self) → datetime.datetime
```

Returns the birth date of the animal. :rtype: datetime.date

```
get_date_of_acquisition (self) → datetime.datetime
```

Get the acquisition date of a behavior_session in UTC :rtype: datetime

```
get_driver_line (self) → List[str]
```

Returns the genotype name(s) of the driver line(s). :rtype: list

```
get_equipment_name (self) → str
```

Returns the name of the experimental rig. :rtype: str

```
get_foraging_id (self) → int
```

Get the foraging ID for the behavior session

```
get_full_genotype (self) → str
```

Return the name of the subject’s genotype :rtype: str

```
get_mouse_id (self) → int
```

Returns the LabTracks ID :rtype: int

```
get_ophys_container_id (self) → Union[int, NoneType]
```

```
get_ophys_experiment_ids (self) → Union[List[int], NoneType]
```

```
get_ophys_session_id (self) → Union[int, NoneType]
```

```
get_reporter_line (self) → List[str]
```

Returns the genotype name(s) of the reporter line(s). :rtype: list

```
get_sex (self) → str
```

Returns sex of the animal (M/F) :rtype: str

```
get_stimulus_name (self) → str
```

Get the stimulus set used from the behavior session.pkl file :rtype: str

allensdk.brain_observatory.behavior.session_apis.data_io.behavior_nwb_api module

allensdk.brain_observatory.behavior.session_apis.data_io.behavior_ophys_json_api module

`allensdk.brain_observatory.behavior.session_apis.data_io.behavior_ophys_lims_api` module

`allensdk.brain_observatory.behavior.session_apis.data_io.behavior_ophys_nwb_api` module

`allensdk.brain_observatory.behavior.session_apis.data_io.ophys_lims_api` module

Module contents

`allensdk.brain_observatory.behavior.session_apis.data_transforms` package

Submodules

`allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_transforms` module

class `allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_trans`

Bases: `allensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_base.BehaviorBase`

This class provides methods that transform data extracted from LIMS or JSON data sources into final data products necessary for populating a BehaviorSession.

get_behavior_session_id (*self*)

Returns the behavior_session_id associated with this experiment, if applicable.

get_extended_trials (*self*) → `pandas.core.frame.DataFrame`

Get extended trials from pickle file

Returns

pd.DataFrame A dataframe containing extended behavior trial information.

get_licks (*self*) → `pandas.core.frame.DataFrame`

Get lick data from pickle file. This function assumes that the first sensor in the list of lick_sensors is the desired lick sensor. If this changes we need to update to get the proper line.

Since licks can occur outside of a trial context, the lick times are extracted from the vsyncs and the frame number in *lick_events*. Since we don't have a timestamp for when in "experiment time" the vsync stream starts (from *self.get_stimulus_timestamps*), we compute it by fitting a linear regression (frame number x time) for the *start_trial* and *end_trial* events in the *trial_log*, to true up these time streams.

Returns `pd.DataFrame` Two columns: "time", which contains the sync time of the licks that occurred in this session and "frame", the frame numbers of licks that occurred in this session

get_metadata (*self*) → `allensdk.brain_observatory.behavior.metadata.behavior_metadata.BehaviorMetadata`
Return metadata about the session. :rtype: BehaviorMetadata

get_monitor_delay (*self*) → float

Return monitor delay for behavior only sessions (in seconds)

get_rewards (*self*) → `pandas.core.frame.DataFrame`

Get reward data from pickle file, based on pickle file timestamps (not sync file).

Returns `pd.DataFrame` – A dataframe containing timestamps of delivered rewards.

get_running_acquisition_df (*self*, *lowpass=True*, *zscore_threshold=10.0*) → pandas.core.frame.DataFrame
Get running speed acquisition data from a behavior pickle file.

NOTE: Rebases timestamps with the `self.get_stimulus_timestamps()` method which varies between the BehaviorDataTransformer and the BehaviorOphysDataTransformer.

Parameters

lowpass: bool (default=True) Whether to apply a 10Hz low-pass filter to the running speed data.

zscore_threshold: float The threshold to use for removing outlier running speeds which might be noise and not true signal

Returns

pd.DataFrame

Dataframe with an index of timestamps and the following columns: “speed”: computed running speed “dx”: angular change, computed during data collection “v_sig”: voltage signal from the encoder “v_in”: the theoretical maximum voltage that the encoder

will reach prior to “wrapping”. This should theoretically be 5V (after crossing 5V goes to 0V, or vice versa). In practice the encoder does not always reach this value before wrapping, which can cause transient spikes in speed at the voltage “wraps”.

get_running_speed (*self*, *lowpass=True*) → pandas.core.frame.DataFrame
Get running speed using timestamps from `self.get_stimulus_timestamps`.

NOTE: Do not correct for monitor delay.

Returns pd.DataFrame index: timestamps speed : subject’s running speeds (in cm/s)

get_stimulus_frame_rate (*self*) → float

get_stimulus_presentations (*self*) → pandas.core.frame.DataFrame
Get stimulus presentation data.

NOTE: Uses timestamps that do not account for monitor delay.

Returns pd.DataFrame – Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_templates (*self*) → Union[allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.S
NoneType]
Get stimulus templates (movies, scenes) for behavior session.

Returns

StimulusTemplate or None if there are no images for the experiment

get_stimulus_timestamps (*self*) → numpy.ndarray
Get stimulus timestamps (vsyncs) from pkl file. Align to the (frame, time) points in the trial events.

NOTE: Located with `behavior_session_id`. Does not use the `sync_file` which requires `ophys_session_id`.

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor that do no account for monitor delay.

get_task_parameters (*self*) → dict
Get task parameters from pkl file.

Returns

dict A dictionary containing parameters used to define the task runtime behavior.

get_trials (*self*) → pandas.core.frame.DataFrame
Get trials from pkl file

Returns

pd.DataFrame A dataframe containing behavioral trial start/stop times, and trial data

allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_ophys_data_transforms module

class allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_ophys_data_transforms

Bases: *allensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_transforms.BehaviorDataTransforms*, *allensdk.brain_observatory.behavior.session_apis.abcs.session_base.BehaviorOphysBase*

This class provides methods that transform data extracted from LIMS or JSON data sources into final data products necessary for populating a BehaviorOphysExperiment

get_average_projection (*self*, *image_api=None*)
Get an image whose values are the average obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to avg projection image data and metadata.

get_cell_roi_ids (*self*)

get_cell_specimen_table (*self*)
Get a cell specimen dataframe containing ROI information about cells identified in an ophys experiment.

Returns

pd.DataFrame Cell ROI information organized into a dataframe. Index is the cell ROI IDs.

get_corrected_fluorescence_traces (*self*)
Get motion-corrected fluorescence traces.

Returns

pd.DataFrame Motion-corrected fluorescence traces organized into a dataframe. Index is the cell ROI IDs.

get_dff_traces (*self*)
Get a table of delta fluorescence over fluorescence traces.

Returns

pd.DataFrame The traces of dff (normalized fluorescence) organized into a dataframe. Index is the cell ROI IDs.

get_events (*self*, *filter_scale*: float = 2, *filter_n_time_steps*: int = 20) → pandas.core.frame.DataFrame
Returns events in dataframe format

Parameters

filter_scale: float See filter_events_array for description

filter_n_time_steps: int See filter_events_array for description

See behavior_ophys_experiment.events for return type

get_eye_tracking (*self*, *z_threshold*: float = 3.0, *dilation_frames*: int = 2) → Union[pandas.core.frame.DataFrame, NoneType]
Gets corneal, eye, and pupil ellipse fit data

Parameters

z_threshold [float, optional] The z-threshold when determining which frames likely contain outliers for eye or pupil areas. Influences which frames are considered ‘likely blinks’. By default 3.0

dilation_frames [int, optional]

Determines the number of additional adjacent frames to mark as ‘likely_blink’, by default 2.

Returns

Optional[pd.DataFrame] *_area *_center_x *_center_y *_height *_phi *_width
likely_blink

where “*” can be “corneal”, “pupil” or “eye”

Will return None if class attr _skip_eye_tracking is True.

get_eye_tracking_rig_geometry (*self*) → Union[dict, NoneType]
Get eye tracking rig metadata from behavior + ophys session.

Returns

dict Includes geometry of monitor, camera, LED

get_max_projection (*self*, *image_api*=None)
Get an image whose values are the maximum obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to max projection image data and metadata.

get_metadata (*self*) → allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata.BehaviorOphysMetadata
Return metadata about the session. :rtype: BehaviorOphysMetadata

get_monitor_delay (*self*)
Return the monitor delay (in seconds)

get_motion_correction (*self*)
Get motion correction trace data.

Returns

pd.DataFrame A dataframe containing trace data used during motion correction computation.

get_ophys_experiment_id (*self*)

Returns the ophys_experiment_id for the instantiated BehaviorOphys Session (or BehaviorOphys data fetcher) if applicable.

get_ophys_session_id (*self*)

Returns the behavior + ophys_session_id associated with this experiment, if applicable.

get_ophys_timestamps (*self*)

Get optical physiology frame timestamps.

Returns

np.ndarray Timestamps associated with frames captured by the microscope.

get_raw_dff_data (*self*)

get_rewards (*self*)

Get reward data from pkl file, based on pkl file timestamps (not sync file).

Returns pd.DataFrame – A dataframe containing timestamps of delivered rewards.

get_roi_masks_by_cell_roi_id (*self*, *cell_roi_ids*: Union[int, Iterable[int], NoneType] = None)

Obtains boolean masks indicating the location of one or more ROIs in this session.

Parameters

cell_roi_ids [array-like of int, optional] ROI masks for these rois will be returned. The default behavior is to return masks for all rois.

Returns

result [xr.DataArray]

dimensions are:

- **roi_id** : which roi is described by this mask?
- **row** : index within the underlying image
- **column** : index within the image

values are 1 where an ROI was present, otherwise 0.

Notes

This method helps Allen Institute scientists to look at sessions that have not yet had cell specimen ids assigned. You probably want to use get_roi_masks instead.

get_stimulus_timestamps (*self*)

Return a numpy array of stimulus timestamps uncorrected for monitor delay (in seconds)

get_sync_data (*self*)

get_sync_licks (*self*)

Module contents

Module contents

allensdk.brain_observatory.behavior.stimulus_processing package

Submodules

allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates module

class allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.**StimulusImage**

Bases: object

Container class for image stimuli

name

class allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.**StimulusImageFactory**

Bases: object

Factory for StimulusImage

static from_processed (*warped: numpy.ndarray, unwarped: numpy.ndarray, name: str*) → allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusImage
Creates a StimulusImage from processed input (usually nwb). Image has already been warped and pre-processed

from_unprocessed (*self, input_array: numpy.ndarray, name: str*) → allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusImage
Creates a StimulusImage from unprocessed input (usually pkl). Image needs to be warped and pre-processed

class allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.**StimulusImageSet**

Bases: object

Container class for a collection of image stimuli

image_names

image_set_name

images

items (*self*)

keys (*self*)

to_dataframe (*self*) → pandas.core.frame.DataFrame

values (*self*)

class allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.**StimulusTemplate**

Bases: object

Factory for StimulusTemplate

static from_processed (*image_set_name: str, image_attributes: List[dict], unwarped: List[numpy.ndarray], warped: List[numpy.ndarray]*) → allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate

Create StimulusTemplate from nwb or other processed input. Stimulus templates created this way DO NOT need to be processed to acquire unwarped versions of the images presented.

NOTE: The ordering of image_attributes, unwarped, and warped matter!

NOTE: Warped images display what was seen on a monitor by a subject. Unwarped images display a ‘diagnostic’ version of the stimuli to be presented.

Parameters

image_set_name [str]

The name of the image set. Example: Natural_Images_Lum_Matched_set_TRAINING_2017.07.14

image_attributes [List[dict]] A list of dictionaries containing image metadata. Must at least contain the key:

image_name

But will usually also contain: image_category, orientation, phase, spatial_frequency, image_index

unwarped [List[np.ndarray]] A list of unwarped image arrays

warped [List[np.ndarray]] A list of warped image arrays

Returns

StimulusTemplate A StimulusTemplate object

```
static from_unprocessed (image_set_name: str, image_attributes:
                        List[dict], images: List[np.ndarray]) → al-
                        lensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate
```

Create StimulusTemplate from pkl or unprocessed input. Stimulus templates created this way need to be processed to acquire unwarped versions of the images presented.

NOTE: The ordering of image_attributes and images matter!

NOTE: Warped images display what was seen on a monitor by a subject. Unwarped images display a ‘diagnostic’ version of the stimuli to be presented.

Parameters

image_set_name [str]

The name of the image set. Example: Natural_Images_Lum_Matched_set_TRAINING_2017.07.14

image_attributes [List[dict]] A list of dictionaries containing image metadata. Must at least contain the key:

image_name

But will usually also contain: image_category, orientation, phase, spatial_frequency, image_index

images [List[np.ndarray]] A list of image arrays

Returns

StimulusTemplate A StimulusTemplate object

allensdk.brain_observatory.behavior.stimulus_processing.util module

allensdk.brain_observatory.behavior.stimulus_processing.util.**convert_filepath_caseinsensitive**

allensdk.brain_observatory.behavior.stimulus_processing.util.**get_image_set_name** (*image_set_path*

str)

→

str

Strips the stem from the image_set filename

Module contents

allensdk.brain_observatory.behavior.stimulus_processing.**get_gratings_metadata** (*stimuli:*

Dict,

start_idx:

int

=

0)

→

pan-

das.core.frame.Dat

This function returns the metadata for each unique grating that was presented during the experiment. If no gratings were displayed during this experiment it returns an empty dataframe with the expected columns. Parameters

—— stimuli:

The stimuli field (pkl[‘items’][‘behavior’][‘stimuli’]) loaded from the experiment pkl file.

start_idx: The index to start index column

Returns

pd.DataFrame: DataFrame containing the unique stimuli presented during an experiment.

The columns contained in this DataFrame are ‘image_category’, ‘image_name’, ‘image_set’, ‘phase’, ‘spatial_frequency’, ‘orientation’, and ‘image_index’. This returns empty if no gratings were presented.

allensdk.brain_observatory.behavior.stimulus_processing.**get_images_dict** (*pkl*)

→

Dict

Gets the dictionary of images that were presented during an experiment along with image set metadata and the image specific metadata. This function uses the path to the image pkl file to read the images and their metadata from the pkl file and return this dictionary. Parameters ——— pkl: The pkl file containing the data for the stimuli presented during

experiment

Returns

Dict: A dictionary containing keys images, metadata, and image_attributes. These correspond to paths to image arrays presented, metadata on the whole set of images, and metadata on specific images, respectively.

allensdk.brain_observatory.behavior.stimulus_processing.**get_stimulus_metadata** (*pkl*)

→

pan-

das.core.frame.Dat

Gets the stimulus metadata for each type of stimulus presented during the experiment. The metadata is return for gratings, images, and omitted stimuli. Parameters ——— pkl: the pkl file containing the information about what stimuli were

presented during the experiment

Returns

pd.DataFrame: The dataframe containing a row for every stimulus that was presented during the experiment. The row contains the following data, image_category, image_name, image_set, phase, spatial_frequency, orientation, and image index.

```
allensdk.brain_observatory.behavior.stimulus_processing.get_stimulus_presentations (data,
stim-
u-
lus_timesta
→
pan-
das.core.fra
```

This function retrieves the stimulus presentation dataframe and renames the columns, adds a stop_time column, and set's index to stimulus_presentation_id before sorting and returning the dataframe. :param data: stimulus file associated with experiment id :param stimulus_timestamps: timestamps indicating when stimuli switched

during experiment

Returns stimulus_table: dataframe containing the stimuli metadata as well as what stimuli was presented

```
allensdk.brain_observatory.behavior.stimulus_processing.get_stimulus_templates (pkl:
dict,
grat-
ing_images_dict:
Union[dict,
None-
Type]
=
None)
→
Union[allensdk.b
None-
Type]
```

Gets images presented during experiments from the behavior stimulus file (*.pkl)

Parameters

pkl [dict] Loaded pkl dict containing data for the presented stimuli.

grating_images_dict [Optional[dict]] Because behavior pkl files do not contain image versions of grating stimuli, they must be obtained from an external source. The grating_images_dict is a nested dictionary where top level keys correspond to grating image names (e.g. 'gratings_0.0', 'gratings_270.0') as they would appear in table returned by get_gratings_metadata(). Sub-nested dicts are expected to have 'warped' and 'unwarped' keys where values are numpy image arrays of aforementioned warped or unwarped grating stimuli.

Returns

StimulusTemplate: StimulusTemplate object containing images that were presented during the experiment

```
allensdk.brain_observatory.behavior.stimulus_processing.get_visual_stimuli_df(data,  
                                                                           time)  
                                                                           →  
                                                                           pan-  
                                                                           das.core.frame.DataFrame
```

This function loads the stimuli and the omitted stimuli into a dataframe. These stimuli are loaded from the input data, where the `set_log` and `draw_log` contained within are used to calculate the epochs. These epochs are used as `start_frame` and `end_frame` and converted to times by input stimulus timestamps. The omitted stimuli do not have a `end_frame` by design though there duration is always 250ms. :param data: the behavior data file :param time: the stimulus timestamps indicating when each stimulus is

displayed

Returns df: a pandas dataframe containing the stimuli and omitted stimuli that were displayed with their frame, end_frame, start_time, and duration

```
allensdk.brain_observatory.behavior.stimulus_processing.is_change_event(stimulus_presentations:  
                                                                        pan-  
                                                                        das.core.frame.DataFrame)  
                                                                        →  
                                                                        pan-  
                                                                        das.core.series.Series
```

Returns whether a stimulus is a change stimulus A change stimulus is defined as the first presentation of a new image_name Omitted stimuli are ignored The first stimulus in the session is ignored

:param stimulus_presentations The stimulus presentations table

Returns is_change: pd.Series indicating whether a given stimulus is a change stimulus

```
allensdk.brain_observatory.behavior.stimulus_processing.load_pickle(pstream)
```

```
allensdk.brain_observatory.behavior.stimulus_processing.unpack_change_log(change)
```

allensdk.brain_observatory.behavior.sync package

Submodules

allensdk.brain_observatory.behavior.sync.process_sync module

```
allensdk.brain_observatory.behavior.sync.process_sync.calculate_delay(sync_data,  
                                                                      stim_vsync_fall,  
                                                                      sample  
                                                                      frequency)
```

```
allensdk.brain_observatory.behavior.sync.process_sync.filter_digital(rising,  
                                                                      falling,  
                                                                      threshold  
                                                                      old=0.0001)
```

Removes short transients from digital signal.

Rising and falling should be same length and units in seconds.

Kwargs: threshold (float): transient width

Module contents

Created on Sunday July 15 2018

@author: marinag

```
allensdk.brain_observatory.behavior.sync.get_behavior_monitoring(dataset: al-
                                                                    lensdk.brain_observatory.sync_dataset.Dataset,
                                                                    permissive:
                                                                    bool =
                                                                    False) →
                                                                    Union[numpy.ndarray,
                                                                    NoneType]
```

Report the timestamps of each frame of the behavior monitoring video

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If] False, a KeyError will be raised

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no behavior monitoring timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_eye_tracking(dataset: al-
                                                                    lensdk.brain_observatory.sync_dataset.Dataset,
                                                                    permissive:
                                                                    bool =
                                                                    False) →
                                                                    Union[numpy.ndarray,
                                                                    NoneType]
```

Report the timestamps of each frame of the eye tracking video

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If] False, a KeyError will be raised

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no eye tracking timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_lick_times(dataset: al-
                                                                    lensdk.brain_observatory.sync_dataset.Dataset,
                                                                    permissive:
                                                                    bool = False)
→ Union[numpy.ndarray,
NoneType]
```

Report the timestamps of each detected lick

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If] False, a KeyError will be raised

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no lick timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_ophys_frames(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool  
= False) →  
numpy.ndarray
```

Report the timestamps of each optical physiology video frame

Parameters

dataset [describes experiment timing]

Returns

array of timestamps (floating point; seconds; relative to experiment start).

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Notes

use rising edge for Scientifica, falling edge for Nikon <http://confluence.corp.alleninstitute.org/display/IT/Ophys+Time+Sync> This function uses rising edges

```
allensdk.brain_observatory.behavior.sync.get_raw_stimulus_frames(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool  
= False) →  
numpy.ndarray
```

Report the raw timestamps of each stimulus frame. This corresponds to the time at which the psychopy window's flip method returned, but not necessarily to the time at which the stimulus frame was displayed.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

array of timestamps (floating point; seconds; relative to experiment start).

```
allensdk.brain_observatory.behavior.sync.get_stim_photodiode(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool  
= False) →  
Union[List[float],  
NoneType]
```

Report the timestamps of each detected sync square transition (both black -> white and white -> black) in this experiment.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no photodiode timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_sync_data(sync_path: str, permissive:
                                                    bool = False) → Dict[str,
                                                    Union[List, numpy.ndarray,
                                                    NoneType]]
```

Convenience function for extracting several timestamp arrays from a sync file.

Parameters

sync_path [The hdf5 file here ought to be a Visual Behavior sync output] file. See `allensdk.brain_observatory.sync_dataset` for more details of this format.

permissive [If True, None will be returned if timestamps are not found. If False, a `KeyError` will be raised

Returns

A dictionary with the following keys. All timestamps in seconds: `ophys_frames` : timestamps of each optical physiology frame `lick_times` : timestamps of each detected lick `ophys_trigger` : The time at which ophys acquisition was started `eye_tracking` : timestamps of each eye tracking video frame `behavior_monitoring` : timestamps of behavior monitoring video frame `stim_photodiode` : timestamps of each photodiode transition `stimulus_times_no_delay` : raw stimulus frame timestamps

Some values may be None. This indicates that the corresponding timestamps were not located in this sync file.

```
allensdk.brain_observatory.behavior.sync.get_trigger(dataset: allensdk.brain_observatory.sync_dataset.Dataset,
                                                    permissive: bool = False) →
                                                    Union[numpy.ndarray, NoneType]
```

Returns (as a 1-element array) the time at which optical physiology acquisition was started.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a `KeyError` will be raised

Returns

timestamps (floating point; seconds; relative to experiment start) or None. If None, no timestamps were found in this sync dataset.

Notes

Ophys frame timestamps can be recorded before acquisition start when experimenters are setting up the recording session. These do not correspond to acquired ophys frames.

`allensdk.brain_observatory.behavior.write_behavior_nwb` package

Module contents

`allensdk.brain_observatory.behavior.write_nwb` package

Subpackages

`allensdk.brain_observatory.behavior.write_nwb.extensions` package

Subpackages

`allensdk.brain_observatory.behavior.write_nwb.extensions.event_detection` package

Submodules

`allensdk.brain_observatory.behavior.write_nwb.extensions.event_detection.extension_builder` module

`allensdk.brain_observatory.behavior.write_nwb.extensions.event_detection.ndx_ophys_events` module

Module contents

`allensdk.brain_observatory.behavior.write_nwb.extensions.stimulus_template` package

Submodules

`allensdk.brain_observatory.behavior.write_nwb.extensions.stimulus_template.extension_builder` module

`allensdk.brain_observatory.behavior.write_nwb.extensions.stimulus_template.ndx_stimulus_template` module

Module contents

Module contents

Module contents

Submodules

`allensdk.brain_observatory.behavior.behavior_ophys_analysis` module

`allensdk.brain_observatory.behavior.behavior_ophys_experiment` module

`allensdk.brain_observatory.behavior.behavior_ophys_session` module

allensdk.brain_observatory.behavior.behavior_session module**allensdk.brain_observatory.behavior.criteria module**

Functions for calculating mtrain state transitions. If criteria are met, return true. Otherwise, return false.

`allensdk.brain_observatory.behavior.criteria.consistency_is_key(session_summary)`
 need some way to judge consistency of various parameters

- dprime
- num trials
- hit rate
- fa rate
- lick timing

`allensdk.brain_observatory.behavior.criteria.consistent_behavior_within_session(session_summary)`
 need some way to measure consistent performance within a session

- compare peak to overall dprime?
- variance in rolling window dprime?

`allensdk.brain_observatory.behavior.criteria.meets_engagement_criteria(session_summary)`
 Returns true if engagement criteria were met for the past 3 days, else false. Args:

`session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak' and 'num_engaged_trials', ordered ascending by training day, for at least 3 days. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function) The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.mostly_useful(trials)`
 Returns True if fewer than half the trial time on the last day were aborted trials.

Args: `trials` (pd.DataFrame): Pandas dataframe with columns 'training_day', 'trial_type', and 'trial_length'.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.n_complete(threshold, count)`
 For compatibility with original API. If count >= threshold, return True. Otherwise return False. Args:

`threshold` (numeric): Threshold for the count to meet. `count` (numeric): The count to compare to the threshold.

Returns: True if count >= threshold, otherwise False.

`allensdk.brain_observatory.behavior.criteria.no_response_bias(session_summary)`
 the mouse meets this criterion if their last session exhibited a response bias between 10% and 90%

Args: `session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'response_bias', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior

(sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

```
allensdk.brain_observatory.behavior.criteria.summer_over(trials)
```

Returns true if the maximum value of 'training_day' in the trials dataframe is ≥ 40 , else false.

```
allensdk.brain_observatory.behavior.criteria.two_out_of_three_aint_bad(session_summary)
```

Returns true if 2 of the last 3 days showed a peak d-prime above 2.

Args: session_summary (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak', ordered ascending by training day, for at least the past 3 days. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

```
allensdk.brain_observatory.behavior.criteria.whole_lotta_trials(session_summary)
```

Mouse meets this criterion if the last session has more than 300 trials. Args:

session_summary (pd.DataFrame): Pandas dataframe with daily values for 'num_contingent_trials', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

```
allensdk.brain_observatory.behavior.criteria.yesterday_was_good(session_summary)
```

Returns true if the last day showed a peak d-prime above 2 Args:

session_summary (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

allensdk.brain_observatory.behavior.dprime module

```
allensdk.brain_observatory.behavior.dprime.get_catch_responses(correct_reject=None,
                                                                false_alarm=None,
                                                                aborted=None)
```

```
allensdk.brain_observatory.behavior.dprime.get_dprime(hit_rate, fa_rate, sliding_window=100)
```

calculates the d-prime for a given hit rate and false alarm rate https://en.wikipedia.org/wiki/Sensitivity_index

Parameters ——— hit_rate : float

rate of hits in the True class

fa_rate [float] rate of false alarms in the False class

limits [tuple, optional] limits on extreme values, which distort. default: (0.01,0.99)

d_prime

```
allensdk.brain_observatory.behavior.dprime.get_false_alarm_rate(correct_reject=None,
                                                                false_alarm=None,
                                                                aborted=None,
                                                                sliding_window=100)
```

```
allensdk.brain_observatory.behavior.dprime.get_go_responses(hit=None,
                                                            miss=None,
                                                            aborted=None)
```

```
allensdk.brain_observatory.behavior.dprime.get_hit_rate(hit=None, miss=None,
                                                         aborted=None, sliding_window=100)
```

```
allensdk.brain_observatory.behavior.dprime.get_rolling_dprime(rolling_hit_rate,
                                                                rolling_fa_rate,
                                                                sliding_window=100)
```

```
allensdk.brain_observatory.behavior.dprime.get_trial_count_corrected_false_alarm_rate(correct_reject=None,
                                                                                       false_alarm=None,
                                                                                       aborted=None,
                                                                                       sliding_window=100)
```

```
allensdk.brain_observatory.behavior.dprime.get_trial_count_corrected_hit_rate(hit=None,
                                                                                miss=None,
                                                                                aborted=None,
                                                                                sliding_window=100)
```

```
allensdk.brain_observatory.behavior.dprime.trial_number_limit(p, N)
```

allensdk.brain_observatory.behavior.event_detection module

```
allensdk.brain_observatory.behavior.event_detection.filter_events_array(arr:
                                                                    numpy.ndarray,
                                                                    scale:
                                                                    float
                                                                    =
                                                                    2,
                                                                    n_time_steps:
                                                                    int
                                                                    =
                                                                    20)
→
numpy.ndarray
```

Convolve the trace array with a 1d causal half-gaussian filter to smooth it for visualization

Uses a halfnorm distribution as weights to the filter

Modified from initial implementation by Nick Ponvert

Parameters

arr: **np.ndarray** Trace matrix of dimension n traces x n frames

scale: float std deviation of halfnorm distribution

n_time_steps: int number of time steps to use for the convolution operation

Returns

np.ndarray: Output of the convolution operation

allensdk.brain_observatory.behavior.eye_tracking_processing module

`allensdk.brain_observatory.behavior.eye_tracking_processing.compute_circular_area(df_row: pandas.core.series.Series) → float`

Calculate the area of the pupil as a circle using the max of the height/width as radius.

Note: This calculation assumes that the pupil is a perfect circle and any eccentricity is a result of the angle at which the pupil is being viewed.

Parameters

df_row [pd.Series] A row from an eye tracking dataframe containing only “pupil_width” and “pupil_height”.

Returns

float The circular area of the pupil in pixels².

`allensdk.brain_observatory.behavior.eye_tracking_processing.compute_elliptical_area(df_row: pandas.core.series.Series) → float`

Calculate the area of corneal reflection (cr) or eye ellipse fits using the ellipse formula.

Parameters

df_row [pd.Series] A row from an eye tracking dataframe containing either: “cr_width”, “cr_height” or “eye_width”, “eye_height”

Returns

float The elliptical area of the eye or cr in pixels²

`allensdk.brain_observatory.behavior.eye_tracking_processing.determine_likely_blinks` (*eye_areas*: `pan-das.core.Series`, *pupil_areas*: `pan-das.core.Series`, *outliers*: `pan-das.core.Series`, *dilation_frames*: `int`)
`= 2)`
`→`
`pan-das.core.Series`

Determine eye tracking frames which contain likely blinks or outliers

Parameters

eye_areas [`pd.Series`] A pandas series of eye areas.

pupil_areas [`pd.Series`] A pandas series of pupil areas.

outliers [`pd.Series`] A pandas series containing bool values of outlier rows.

dilation_frames [`int`, optional] Determines the number of additional adjacent frames to mark as 'likely_blink', by default 2.

Returns

pd.Series A pandas series of bool values that has the same length as the number of eye tracking dataframe rows (frames).

`allensdk.brain_observatory.behavior.eye_tracking_processing.determine_outliers` (*data_df*: `pan-das.core.DataFrame`, *z_threshold*: `float`)
`→`
`pan-das.core.Series`

Given a dataframe and some z-score threshold return a pandas boolean Series where each entry indicates whether a given row contains at least one outlier (where outliers are calculated along columns).

Parameters

data_df [`pd.DataFrame`] A dataframe containing only columns where outlier detection is desired. (e.g. "cr_area", "eye_area", "pupil_area")

z_threshold [`float`] z-score values higher than the `z_threshold` will be considered outliers.

Returns

pd.Series A pandas boolean Series whose length == `len(data_df.index)`. True denotes that a row in the `data_df` contains at least one outlier.

```
allensdk.brain_observatory.behavior.eye_tracking_processing.load_eye_tracking_hdf(eye_tracking_file,  
                                     pathlib.Path)  
→  
pan-  
das.core.frame
```

Load a DeepLabCut hdf5 file containing eye tracking data into a dataframe.

Note: The eye tracking hdf5 file contains 3 separate dataframes. One for corneal reflection (cr), eye, and pupil ellipse fits. This function loads and returns this data as a single dataframe.

Parameters

eye_tracking_file [Path] Path to an hdf5 file produced by the DeepLabCut eye tracking pipeline. The hdf5 file will contain the following keys: “cr”, “eye”, “pupil”. Each key has an associated dataframe with the following columns: “center_x”, “center_y”, “height”, “width”, “phi”.

Returns

pd.DataFrame A dataframe containing combined corneal reflection (cr), eyelid (eye), and pupil data. Column names for each field will be renamed by prepending the field name. (e.g. center_x -> eye_center_x)

```
allensdk.brain_observatory.behavior.eye_tracking_processing.process_eye_tracking_data(eye_data,  
                                          pan-  
das.co  
frame_  
pan-  
das.co  
z_thre  
float  
=  
3.0,  
di-  
la-  
tion_fr  
int  
=  
2),  
→  
pan-  
das.co
```

Processes and refines raw eye tracking data by adding additional computed feature columns.

Parameters

eye_data [pd.DataFrame] A ‘raw’ eye tracking dataframe produced by load_eye_tracking_hdf()

frame_times [pd.Series] A series of frame times acquired from a behavior + ophy session ‘sync file’.

z_threshold [float] z-score values higher than the z_threshold will be considered outliers, by default 3.0.

dilation_frames [int, optional] Determines the number of additional adjacent frames to mark as ‘likely_blink’, by default 2.

Returns

pd.DataFrame A refined eye tracking dataframe that contains additional information about frame times, eye areas, pupil areas, and frames with likely blinks/outliers.

Raises

RuntimeError If the number of sync file frame times does not match the number of eye tracking frames.

allensdk.brain_observatory.behavior.image_api module

```
class allensdk.brain_observatory.behavior.image_api.Image
    Bases: tuple
    Describes a 2D Image
    data [np.ndarray] Image data points
    spacing [tuple] Spacing describes the physical size of each pixel
    unit [str] Physical unit of the spacing (currently constrained to be isotropic)
    data
        Alias for field number 0
    spacing
        Alias for field number 1
    unit
        Alias for field number 2
class allensdk.brain_observatory.behavior.image_api.ImageApi
    Bases: object
    static deserialize (img)
    static serialize (data, spacing, unit)
```

allensdk.brain_observatory.behavior.mtrain module

```
class allensdk.brain_observatory.behavior.mtrain.ExtendedTrialSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
    Bases: marshmallow.schema.Schema
    This schema describes the edf core trial structure
    opts = <marshmallow.schema.SchemaOpts object>
class allensdk.brain_observatory.behavior.mtrain.FriendlyDate (format=None,
                                                                    **kwargs)
    Bases: marshmallow.fields.Date
```

class allensdk.brain_observatory.behavior.mtrain.**FriendlyDateTime** (*format=None*,
***kwargs*)

Bases: `marshmallow.fields.DateTime`

`allensdk.brain_observatory.behavior.mtrain.annotate_change_detect` (*trials*)
adds *change* and *detect* columns to dataframe

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

`allensdk.brain_observatory.behavior.mtrain.annotate_trials` (*trials*)
performs multiple annotations:

- `annotate_change_detect`
- `fix_change_time`
- `explode_response_window`

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

`allensdk.brain_observatory.behavior.mtrain.assign_session_id` (*trials*)

adds a column with a unique ID for the session defined as a combination of the mouse ID and startdatetime

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

`allensdk.brain_observatory.behavior.mtrain.explode_response_window` (*trials*)
explodes the *response_window* column in lower & upper columns

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

`allensdk.brain_observatory.behavior.mtrain.fix_change_time(trials)`
 forces *None* values in the *change_time* column to numpy NaN

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

allensdk.brain_observatory.behavior.rewards_processing module

`allensdk.brain_observatory.behavior.rewards_processing.get_rewards` (*data:*
Dict,
times-
tamps:
numpy.ndarray)
 → *pan-*
das.core.frame.DataFrame

Construct and return a pandas DataFrame containing reward data for this session

data: Dict The dict that results from reading the stimulus pickle file associated with the session

timestamps: np.ndarray[1d] A numpy array of timestamps associated with the stimulus frames in this session.
 timestamps[ii] is the clock time of the iiith frame.

Returns

pd.DataFrame containing the data associated with rewards given in this session

allensdk.brain_observatory.behavior.running_processing module

`allensdk.brain_observatory.behavior.running_processing.calc_deriv` (*x, time*)
`allensdk.brain_observatory.behavior.running_processing.deg_to_dist` (*angular_speed:*
numpy.ndarray)
 →
numpy.ndarray

Takes the angular speed (radians/s) at each step in radians, and computes the linear speed in cm/s.

Parameters

angular_speed: np.ndarray (1d) 1d array of angular speed in radians/s

Returns

—

np.ndarray (1d) Linear speed in cm/s at each time point.

```
allensdk.brain_observatory.behavior.running_processing.get_running_df(data,  
                                                                    time:  
                                                                    numpy.ndarray,  
                                                                    low-  
                                                                    pass:  
                                                                    bool  
                                                                    =  
                                                                    True,  
                                                                    zs-  
                                                                    core_threshold=10.0)
```

Given the data from the behavior ‘pkl’ file object and a 1d array of timestamps, compute the running speed. Returns a dataframe with the raw voltage data as well as the computed speed at each timestamp. By default, the running speed is filtered with a 10 Hz Butterworth lowpass filter to remove artifacts caused by the rotary encoder.

Parameters

data Deserialized ‘behavior pkl’ file data

time: np.ndarray (1d) Timestamps for running data measurements

lowpass: bool (default=True) Whether to apply a 10Hz low-pass filter to the running speed data.

zscore_threshold: float The threshold to use for removing outlier running speeds which might be noise and not true signal

Returns

pd.DataFrame Dataframe with an index of timestamps and the following columns:

“speed”: computed running speed “dx”: angular change, computed during data collection “v_sig”: voltage signal from the encoder “v_in”: the theoretical maximum voltage that the encoder

will reach prior to “wrapping”. This should theoretically be 5V (after crossing 5V goes to 0V, or vice versa). In practice the encoder does not always reach this value before wrapping, which can cause transient spikes in speed at the voltage “wraps”.

The raw data are provided so that the user may compute their own speed from source, if desired.

Notes

Though the angular change is available in the raw data (key=“dx”), this method recomputes the angular change from the voltage signal (key=“vsig”) due to very specific, low-level artifacts in the data caused by the encoder. See method docstrings for more detailed information. The raw data is included in the final output in case the end user wants to apply their own corrections and compute running speed from the raw source.

allensdk.brain_observatory.behavior.schemas module

```
class allensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.RaisingSchema*

This schema contains metadata pertaining to behavior.

```
neurodata_doc = 'Metadata for behavior and behavior + ophys experiments'
neurodata_skip = {'date_of_acquisition'}
neurodata_type = 'BehaviorMetadata'
neurodata_type_inc = 'LabMetaData'
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.RaisingSchema*

This schema encompasses task parameters used for behavior or ophys + behavior.

```
neurodata_doc = 'Metadata for behavior or behavior + ophys task parameters'
neurodata_type = 'BehaviorTaskParameters'
neurodata_type_inc = 'LabMetaData'
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.CompleteOphysBehaviorMetadataSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
par-
tial=False,
un-
known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema*,
allensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema

This schema combines fields from behavior, ophys, and subject schemas. Metadata info is passed by the behavior+ophys session in a combined lump containing all the field types.

opts = <marshmallow.schema.SchemaOpts object>

```
class allensdk.brain_observatory.behavior.schemas.EyeTrackingRigGeometry (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
par-
tial=False,
un-
known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.RaisingSchema*

Eye tracking rig geometry

opts = <marshmallow.schema.SchemaOpts object>

```
class allensdk.brain_observatory.behavior.schemas.NwbOphysMetadataSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
par-
tial=False,
un-
known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.RaisingSchema*

This schema contains fields that will be stored in pyNWB base classes pertaining to optical physiology.

opts = <marshmallow.schema.SchemaOpts object>

```
class allensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema (only=None,
                                                                 ex-
                                                                 clude=(),
                                                                 many=False,
                                                                 con-
                                                                 text=None,
                                                                 load_only=(),
                                                                 dump_only=(),
                                                                 par-
                                                                 tial=False,
                                                                 un-
                                                                 known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema`,
`allensdk.brain_observatory.behavior.schemas.OphysMetadataSchema`

This schema contains fields pertaining to ophys+behavior. It is used as a template for generating our custom NWB behavior + ophys extension.

```
neurodata_doc = 'Metadata for behavior + ophys experiments'
```

```
neurodata_skip = {'date_of_acquisition', 'emission_lambda', 'excitation_lambda', 'indi
```

```
neurodata_type = 'OphysBehaviorMetadata'
```

```
neurodata_type_inc = 'BehaviorMetadata'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.OphysEyeTrackingRigMetadataSchema (only=None,
                                                                 ex-
                                                                 clude=(),
                                                                 many=False,
                                                                 con-
                                                                 text=None,
                                                                 load_only=(),
                                                                 dump_only=(),
                                                                 par-
                                                                 tial=False,
                                                                 un-
                                                                 known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.RaisingSchema`

This schema encompasses metadata for ophys experiment rig

```
neurodata_doc = 'Metadata for ophys experiment rig'
```

```
neurodata_type = 'OphysEyeTrackingRigMetadata'
```

```
neurodata_type_inc = 'NWBDataInterface'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.OphysMetadataSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.NwbOphysMetadataSchema*

This schema contains metadata pertaining to optical physiology (ophys).

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.RaisingSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    context=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False, un-
                                                                    known=None)
```

Bases: *marshmallow.schema.Schema*

```
class Meta
```

Bases: *object*

```
unknown = 'raise'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: *allensdk.brain_observatory.behavior.schemas.RaisingSchema*

This schema contains metadata pertaining to a subject in either a behavior or behavior + ophys experiment.

```
neurodata_doc = 'Metadata for an AIBS behavior or behavior + ophys subject'
```

```
neurodata_skip = {'age_in_days', 'genotype', 'sex', 'subject_id'}
```

```
neurodata_type = 'BehaviorSubject'
```

```
neurodata_type_inc = 'Subject'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```


allensdk.brain_observatory.behavior.session_metrics module

`allensdk.brain_observatory.behavior.session_metrics.num_contingent_trials(session_trials)`

Returns the number of “go” and “catch” trials in a training session dataframe. Args:

session_trials (pandas.DataFrame): a pandas.DataFrame describing behavior training trials, with the string column “trial_type” describing the type of trial.

Returns (int): Number of “go” and “catch” trials

`allensdk.brain_observatory.behavior.session_metrics.response_bias(trials, detect_col, trial_types=('go', 'catch'))`

Calculate the response bias for a subset of trial types from a behavioral training dataframe. Args:

trials (pandas.DataFrame): Dataframe containing trial-level information from a behavioral training session. Required columns: “trial_type”, *detect_col*.

detect_col (str): Name of column containing boolean or numeric codings (0/1) for whether or not the mouse had a response.

trial_types (iterable<str>): Iterable containing string trial types to check for the response bias. Trials of types not included in this iterable will be ignored. Default=(“go”, “catch”)

Return: The response bias (or average value of the *detect_col*) for trials in *trial_types*.

allensdk.brain_observatory.behavior.trial_masks module

`allensdk.brain_observatory.behavior.trial_masks.contingent_trials(trials)`

GO & CATCH trials only

Parameters

trials [pandas DataFrame] dataframe of trials

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

`allensdk.brain_observatory.behavior.trial_masks.reward_rate(trials, thresh=2.0)`

masks trials where the reward rate (per minute) is below some threshold.

This de facto omits trials in which the animal was not licking for extended periods or periods when they were licking indiscriminantly.

Parameters

trials [pandas DataFrame] dataframe of trials

thresh [float, optional] threshold under which trials will not be included, default: 2.0

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

`allensdk.brain_observatory.behavior.trial_masks.trial_types(trials, trial_types)`

only include trials of certain trial types

Parameters

trials [pandas DataFrame] dataframe of trials

trial_types [list or other iterator]

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

allensdk.brain_observatory.behavior.trials_processing module

`allensdk.brain_observatory.behavior.trials_processing.calculate_response_latency_list` (*trials: pandas.core.series.Series, response_latency: float*)
→ List
per trial, detemines a response latency

Parameters

trials: pd.DataFrame contains columns “lick_times” and “change_times”

response_window_start: float [seconds] relative to the non-display-lag-compensated presentation of the change-image

Returns

response_latency_list: list len() = trials.shape[0] value is ‘inf’ if there are no valid licks in the trial

`allensdk.brain_observatory.behavior.trials_processing.calculate_reward_rate` (*response_latency=None, start_time=None, win_dow=0.75, trial_window=25, initial_trials=10*)
`allensdk.brain_observatory.behavior.trials_processing.calculate_reward_rate_fix_nans` (*trials: pandas.core.series.Series, response_latency: float*)
→ numpy.ndarray
per trial, detemines the reward rate, replacing infs with nans

Parameters

trials: pd.DataFrame contains columns “lick_times”, “change_times”, and “start_time”

response_window_start: float [seconds] relative to the non-display-lag-compensated presentation of the change-image

Returns

reward_rate: np.ndarray size = trials.shape[0] value is nan if calculate_reward_rate evaluates to ‘inf’

```
allensdk.brain_observatory.behavior.trials_processing.categorize_one_trial(tr)
```

```
allensdk.brain_observatory.behavior.trials_processing.colormap(trial_type,  
                                                                response_type)
```

```
allensdk.brain_observatory.behavior.trials_processing.construct_rolling_performance_df(trials  
pan-  
das.c  
re-  
spon-  
ses-  
sion_  
→  
pan-  
das.c
```

Return a DataFrame containing trial by trial behavior response performance metrics.

Parameters

trials: **pd.DataFrame** contains columns “lick_times”, “change_times”, and “start_time”

response_window_start: **float** [seconds] relative to the non-display-lag-compensated presentation of the change-image

session_type: **str** used to check if this was a passive session

Returns

pd.DataFrame

A pandas DataFrame containing:

trials_id [index]: Index of the trial. All trials, including aborted trials, are assigned an index starting at 0 for the first trial.

reward_rate: Rewards earned in the previous 25 trials, normalized by the elapsed time of the same 25 trials. Units are rewards/minute.

hit_rate_raw: Fraction of go trials where the mouse licked in the response window, calculated over the previous 100 non-aborted trials. Without trial count correction applied.

hit_rate: Fraction of go trials where the mouse licked in the response window, calculated over the previous 100 non-aborted trials. With trial count correction applied.

false_alarm_rate_raw: Fraction of catch trials where the mouse licked in the response window, calculated over the previous 100 non-aborted trials. Without trial count correction applied.

false_alarm_rate: Fraction of catch trials where the mouse licked in the response window, calculated over the previous 100 non-aborted trials. Without trial count correction applied.

rolling_dprime: d prime calculated using the rolling hit_rate and rolling false_alarm_rate.

```
allensdk.brain_observatory.behavior.trials_processing.create_extended_trials(trials=None,  
meta-  
data=None,  
time=None,  
licks=None)
```

```
allensdk.brain_observatory.behavior.trials_processing.data_to_licks(data,
                                                                    time)
allensdk.brain_observatory.behavior.trials_processing.data_to_metadata(data,
                                                                    time)
allensdk.brain_observatory.behavior.trials_processing.find_licks(reward_times,
                                                                    licks, win-
                                                                    dow=3.5)
allensdk.brain_observatory.behavior.trials_processing.get_change_time_frame_response_latency
allensdk.brain_observatory.behavior.trials_processing.get_even_sampling(data)
    Get status of even_sampling
```

Parameters

data: Mapping foraging2 experiment output data

Returns

bool: True if even_sampling is enabled

```
allensdk.brain_observatory.behavior.trials_processing.get_extended_trials(data,
                                                                    time=None)
allensdk.brain_observatory.behavior.trials_processing.get_image_info_from_trial(trial_log,
                                                                    ti)
allensdk.brain_observatory.behavior.trials_processing.get_mouse_id(exp_data)
allensdk.brain_observatory.behavior.trials_processing.get_ori_info_from_trial(trial_log,
                                                                    ti)
allensdk.brain_observatory.behavior.trials_processing.get_params(exp_data)
allensdk.brain_observatory.behavior.trials_processing.get_response_latency(change_event,
                                                                    trial)
allensdk.brain_observatory.behavior.trials_processing.get_response_type(trials)
allensdk.brain_observatory.behavior.trials_processing.get_stimulus_attr_changes(stim_dict,
                                                                    change_frame,
                                                                    first_frame,
                                                                    last_frame)
```

Notes

- assumes only two stimuli are ever shown
- converts attr_names to lowercase
- gets the net attr changes from the start of a trial to the end of a trial

```
allensdk.brain_observatory.behavior.trials_processing.get_time(exp_data)
allensdk.brain_observatory.behavior.trials_processing.get_trial_bounds(trial_log:
                                                                    List)
                                                                    →
                                                                    List
```

Adjust trial boundaries from a trial_log so that there is no dead time between trials.

Parameters

trial_log: list The trial_log read in from the well known behavior stimulus pickle file

Returns

list Each element in the list is a tuple of the form (start_frame, end_frame) so that the ith element of the list gives the start and end frames of the ith trial. The endframe of the last trial will be -1, indicating that it should map to the last timestamp in the session

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_image_names (trial,
                                                                              stimuli)
                                                                              →
                                                                              Dict[str,
                                                                              str]
```

Gets the name of the stimulus presented at the beginning of the trial and what is it changed to at the end of the trial. Parameters ——— trial: A trial in a behavior ophys session stimuli: The stimuli presentation log for the behavior session

Returns

A dictionary indicating the starting_stimulus and what the stimulus is changed to.

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_reward_time (rebased_reward_time,
                                                                              start_time,
                                                                              stop_time)
```

extract reward times in time range

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_timing (event_dict:
                                                                              dict,
                                                                              licks:
                                                                              List[float],
                                                                              go:
                                                                              bool,
                                                                              catch:
                                                                              bool,
                                                                              auto_rewarded:
                                                                              bool,
                                                                              hit:
                                                                              bool,
                                                                              false_alarm:
                                                                              bool,
                                                                              aborted:
                                                                              bool,
                                                                              times-
                                                                              tamps:
                                                                              numpy.ndarray,
                                                                              mon-
                                                                              i-
                                                                              tor_delay:
                                                                              float)
```

Extract a dictionary of trial timing data. See trial_data_from_log for a description of the trial types.

Parameters

event_dict: dict Dictionary of trial events in the well-known *pkl* file

licks: List[float] list of lick timestamps, from the *get_licks* response for the BehaviorOphysExperiment.api.

go: bool True if “go” trial, False otherwise. Mutually exclusive with *catch*.

catch: bool True if “catch” trial, False otherwise. Mutually exclusive with *go*.

auto_rewarded: bool True if “auto_rewarded” trial, False otherwise.

hit: bool True if “hit” trial, False otherwise

false_alarm: bool True if “false_alarm” trial, False otherwise

aborted: bool True if “aborted” trial, False otherwise

timestamps: np.ndarray[1d] Array of ground truth timestamps for the session (sync times, if available)

monitor_delay: float The monitor delay in seconds associated with the session

Returns

dict

start_time: float The time the trial started (in seconds elapsed from recording start)

stop_time: float The time the trial ended (in seconds elapsed from recording start)

trial_length: float Duration of the trial in seconds

response_time: float The response time, for non-aborted trials. This is equal to the first lick in the trial. For aborted trials or trials without licks, *response_time* is NaN.

change_frame: int The frame number that the stimulus changed

change_time: float The time in seconds that the stimulus changed

response_latency: float or None The time in seconds between the stimulus change and the animal’s lick response, if the trial is a “go”, “catch”, or “auto_rewarded” type. If the animal did not respond, return *float(“inf”)*. In all other cases, return None.

Notes

The following parameters are mutually exclusive (exactly one can be true):

hit, miss, false_alarm, aborted, auto_rewarded

`allensdk.brain_observatory.behavior.trials_processing.get_trials_from_data_transform(input_tr`
→
pan-
das.core

Create and return a pandas DataFrame containing data about the trials associated with this session

Parameters

input_transform: An instantiation of a class that inherits from either `BehaviorDataTransform` or `BehaviorOphysDataTransform`. This object will be used to get at the data needed by this method to create the trials dataframe.

Returns

pd.DataFrame A dataframe containing data pertaining to the trials that make up this session

Notes

The `input_transform` object must have the following methods:

input_transform.behavior_stimulus_file Which returns the dict resulting from reading in this session’s stimulus_data pickle file

input_transform.get_rewards Which returns a dataframe containing data about rewards given during this session, i.e. the output of `allensdk/brain_observatory/behavior/rewards_processing.get_rewards`

input_transform.get_licks Which returns a dataframe containing the columns *time* and *frame* denoting the time (in seconds) and frame number at which licks occurred during this session

input_transform.get_stimulus_timestamps Which returns a `numpy.ndarray` of timestamps (in seconds) associated with the frames presented in this session.

input_transform.get_monitor_delay Which returns the monitor delay (in seconds) associated with the experimental rig

```
allensdk.brain_observatory.behavior.trials_processing.get_trials_v0(data,
                                                                    time)
allensdk.brain_observatory.behavior.trials_processing.local_time(iso_timestamp,
                                                                    time-
                                                                    zone=None)
allensdk.brain_observatory.behavior.trials_processing.resolve_initial_image(stimuli,
                                                                              start_frame)
```

Attempts to resolve the initial image for a given `start_frame` for a trial

Parameters

stimuli: Mapping foraging2 shape stimuli mapping

start_frame: int start frame of the trial

Returns

initial_image_category_name: str stimulus category of initial image

initial_image_group: str group name of the initial image

initial_image_name: str name of the initial image

```
allensdk.brain_observatory.behavior.trials_processing.trial_data_from_log(trial)
```

Infer trial logic from trial log. Returns a dictionary.

- **reward volume:** volume of water delivered on the trial, in mL

Each of the following values is boolean:

Trial category values are mutually exclusive * **go:** trial was a go trial (trial with a stimulus change) * **catch:** trial was a catch trial (trial with a sham stimulus change)

stimulus_change/sham_change are mutually exclusive * **stimulus_change:** did the stimulus change (True on 'go' trials) * **sham_change:** stimulus did not change, but response was evaluated

(True on 'catch' trials)

Each trial can be one (and only one) of the following: * **hit** (stimulus changed, animal responded in response window) * **miss** (stimulus changed, animal did not respond in response window) * **false_alarm** (stimulus did not change,

animal responded in response window)

- **correct_reject** (stimulus did not change, animal did not respond in response window)
- **aborted** (animal responded before change time)
- **auto_rewarded** (reward was automatically delivered following the change. This will bias the animals choice and should not be categorized as hit/miss)

`allensdk.brain_observatory.behavior.trials_processing.validate_trial_condition_exclusivity`

ensure that only one of N possible mutually exclusive trial conditions is True

`allensdk.brain_observatory.behavior.validation` module

Module contents

`allensdk.brain_observatory.ecephys` package

Subpackages

`allensdk.brain_observatory.ecephys.align_timestamps` package

Submodules

`allensdk.brain_observatory.ecephys.align_timestamps.barcode` module

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.extract_barcodes_from_times` (*on_*

off_
in-
ter_
bar
bar
cod
nbin

Read barcodes from timestamped rising and falling edges.

Parameters

on_times [numpy.ndarray] Timestamps of rising edges on the barcode line

off_times [numpy.ndarray] Timestamps of falling edges on the barcode line

inter_barcode_interval [numeric, optional] Minimum duration of time between barcodes.

bar_duration [numeric, optional] A value slightly shorter than the expected duration of each bar

barcode_duration_ceiling [numeric, optional] The maximum duration of a single barcode

nbits [int, optional] The bit-depth of each barcode

Returns

barcode_start_times [list of numeric] For each detected barcode, the time at which that barcode started

barcodes [list of int] For each detected barcode, the value of that barcode as an integer.

Notes

ignores first code in prod (ok, but not intended) ignores first on pulse (intended - this is needed to identify that a barcode is starting)

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.find_matching_index` (*master_barcode*, *probe_barcode*, *alignment_type*=*'start'* or *'end'*)

Given a set of barcodes for the master clock and the probe clock, find the indices of a matching set, either starting from the beginning or the end of the list.

Parameters

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe line. One per barcode

alignment_type [string] *'start'* or *'end'*

Returns

master_barcode_index [int] matching index for master barcodes (None if not found)

probe_barcode_index [int] matching index for probe barcodes (None if not found)

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.get_probe_time_offset` (*master_times*, *master_barcodes*, *probe_times*, *probe_barcode_index*, *acq_start_index*, *local_probe_rate*)

Time offset between master clock and recording probes. For converting probe time to master clock.

Parameters

master_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

acq_start_index [int] sample index of probe acquisition start time

local_probe_rate [float] the probe's apparent sampling rate

Returns

total_time_shift [float] Time at which the probe started acquisition, assessed on the master clock. If < 0 , the probe started earlier than the master line.

probe_rate [float] The probe's sampling rate, assessed on the master clock

master_endpoints [iterable] Defines the start and end times of the sync interval on the master clock

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.linear_transform_from_intervals`

Find a scale and translation which aligns two 1d segments

Parameters

master [iterable] Pair of floats defining the master interval. Order is [start, end].

probe [iterable] Pair of floats defining the probe interval. Order is [start, end].

Returns

scale [float] Scale factor. If > 1.0 , the probe clock is running fast compared to the master clock. If < 1.0 , the probe clock is running slow.

translation [float] If > 0 , the probe clock started before the master clock. If < 0 , after.

Notes

solves $(\text{master} + \text{translation}) * \text{scale} = \text{probe}$

for scale and translation

```
allensdk.brain_observatory.ecephys.align_timestamps.barcode.match_barcodes(master_times,  
                                                                           mas-  
                                                                           ter_barcodes,  
                                                                           probe_times,  
                                                                           probe_barcodes)
```

Given sequences of barcode values and (local) times on a probe line and a master line, find the time points on each clock corresponding to the first and last shared barcode.

If there's only one probe barcode, only the first matching timepoint is returned.

Parameters

master_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

Returns

probe_interval [np.ndarray] Start and end times of the matched interval according to the probe_clock.

master_interval [np.ndarray] Start and end times of the matched interval according to the master clock

allensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset module

```
class allensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset.BarcodeSyncDataset
```

Bases: `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset`

barcode_line

Obtain the index of the barcode line for this dataset.

```
extract_barcodes(self, **barcode_kwargs)
```

Read barcodes and their times from this dataset's barcode line.

Parameters

****barcode_kwargs** : Will be passed to `.barcode.extract_barcodes_from_times`

Returns

times [np.ndarray] The start times of each detected barcode.

codes [np.ndarray] The values of each detected barcode

get_barcode_table (*self*, ***barcode_kwargs*)

A convenience method for getting barcode times and codes in a dictionary.

Notes

This method is deprecated!

allensdk.brain_observatory.ecephys.align_timestamps.channel_states module

allensdk.brain_observatory.ecephys.align_timestamps.channel_states.**extract_barcodes_from_st**

Obtain barcodes from timestamped rising/falling edges.

Parameters

channel_states [numpy.ndarray] Rising and falling edges, denoted 1 and -1

timestamps [numpy.ndarray] Sample index of each event.

sampling_rate [numeric] Samples / second

****barcode_kwargs** : Additional parameters describing the barcodes.

allensdk.brain_observatory.ecephys.align_timestamps.channel_states.**extract_splits_from_stat**

Obtain barcodes from timestamped rising/falling edges.

Parameters

channel_states [numpy.ndarray] Rising and falling edges, denoted 1 and -1

timestamps [numpy.ndarray] Sample index of each event.

sampling_rate [numeric] Samples / second

****barcode_kwargs** : Additional parameters describing the barcodes.

allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer module

class allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer.**ProbeSynchronizer**

Bases: object

classmethod compute (*master_barcode_times*, *master_barcodes*, *probe_barcode_times*,
probe_barcodes, *min_time*, *max_time*, *probe_start_index*, *local_probe_sampling_rate*)

Compute a transform from probe samples to master times by aligning barcodes.

Parameters

master_barcode_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_barcode_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

min_time [Float] time (in seconds) of first barcode to align

max_time [Float] time (in seconds) of last barcode to align

probe_start_index [int] sample index of probe acquisition start time

local_probe_sampling_rate [float] the probe's apparent sampling rate

Returns

ProbeSynchronizer : When called, applies the transform computed here to samples on the probe clock.

sampling_rate_scale

The ratio of the probe's sampling rate assessed on the global clock to the probe's locally assessed sampling rate.

Module contents

allensdk.brain_observatory.ecephys.copy_utility package

Module contents

allensdk.brain_observatory.ecephys.current_source_density package

Module contents

allensdk.brain_observatory.ecephys.ecephys_project_api package

Submodules

allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api module

class allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.**EcephysProjectApi**

Bases: object

get_channels (*self*, *channel_ids*: Union[~ArrayLike, NoneType] = None, *probe_ids*: Union[~ArrayLike, NoneType] = None, *session_ids*: Union[~ArrayLike, NoneType] = None, *published_at*: Union[str, NoneType] = None)

get_isi_experiments (*self*, *args, **kwargs)

get_natural_movie_template (*self*, *number*) → Iterable

get_natural_scene_template (*self*, *number*) → Iterable

get_probe_lfp_data (*self*, *probe_id*: int) → Iterable

get_probes (*self*, *probe_ids*: Union[~ArrayLike, NoneType] = None, *session_ids*: Union[~ArrayLike, NoneType] = None, *published_at*: Union[str, NoneType] = None)

get_session_data (*self*, *session_id*: int) → Iterable

get_sessions (*self*, *session_ids*: Union[~ArrayLike, NoneType] = None, *published_at*: Union[str, NoneType] = None)

get_unit_analysis_metrics (*self*, *unit_ids*: Union[~ArrayLike, NoneType] = None, *ecephys_session_ids*: Union[~ArrayLike, NoneType] = None, *session_types*: Union[~ArrayLike, NoneType] = None) → pandas.core.frame.DataFrame

get_units (*self*, *unit_ids*: Union[~ArrayLike, NoneType] = None, *channel_ids*: Union[~ArrayLike, NoneType] = None, *probe_ids*: Union[~ArrayLike, NoneType] = None, *session_ids*: Union[~ArrayLike, NoneType] = None, *published_at*: Union[str, NoneType] = None)

allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_fixed_api module

allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_lims_api module

allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_warehouse_api module

allensdk.brain_observatory.ecephys.ecephys_project_api.http_engine module

allensdk.brain_observatory.ecephys.ecephys_project_api.rma_engine module

allensdk.brain_observatory.ecephys.ecephys_project_api.utilities module

Module contents

allensdk.brain_observatory.ecephys.ecephys_session_api package

Submodules

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_nwb1_session_api` module

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_nwb_session_api` module

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api` module

class `allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi`

Bases: `object`

get_channels (*self*) → `pandas.core.frame.DataFrame`

get_ecephys_session_id (*self*) → `int`

get_invalid_times (*self*) → `pandas.core.frame.DataFrame`

get_lfp (*self*, *probe_id: int*) → `xarray.core.dataarray.DataArray`

get_mean_waveforms (*self*) → `Dict[int, numpy.ndarray]`

get_metadata (*self*)

get_optogenetic_stimulation (*self*) → `pandas.core.frame.DataFrame`

get_probes (*self*) → `pandas.core.frame.DataFrame`

get_pupil_data (*self*) → `Union[pandas.core.frame.DataFrame, NoneType]`

get_rig_metadata (*self*) → `Union[dict, NoneType]`

get_running_speed (*self*) → `allensdk.brain_observatory.running_speed.RunningSpeed`

get_screen_gaze_data (*self*, *include_filtered_data=False*) → `Union[pandas.core.frame.DataFrame, NoneType]`

get_session_start_time (*self*) → `datetime.datetime`

get_spike_amplitudes (*self*) → `Dict[int, numpy.ndarray]`

get_spike_times (*self*) → `Dict[int, numpy.ndarray]`

get_stimulus_presentations (*self*) → `pandas.core.frame.DataFrame`

get_units (*self*) → `pandas.core.frame.DataFrame`

session_na = `-1`

test (*self*) → `bool`

Module contents

`allensdk.brain_observatory.ecephys.file_io` package

Submodules

allensdk.brain_observatory.ecephys.file_io.continuous_file module

```
class allensdk.brain_observatory.ecephys.file_io.continuous_file.ContinuousFile(data_path,
                                     times-
                                     tamps_path,
                                     to-
                                     tal_num_chann
                                     dtype=<class
                                     'numpy.int16'>)
```

Bases: object

Represents a continuous (.dat) file, and its associated timestamps

```
get_lfp_channel_order(self)
    Returns the channel ordering for LFP data extracted from NPX files.

    None
```

```
load(self, memmap=False, memmap_thresh=10000000000.0)
    Reads lfp data and timestamps from the filesystem
```

memmap [bool, optional] If True, the returned data array will be a memory map of the file on disk. Default is True.

memmap_thresh [float, optional] Files above this size in bytes will be memory-mapped, regardless of memmap setting

allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset module

```
class allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset
    Bases: allensdk.brain_observatory.sync_dataset.Dataset
```

```
extract_frame_times(self, strategy, photodiode_cycle=60, frame_keys=('frames', 'stim_vsync'),
                    photodiode_keys=('photodiode', 'stim_photodiode'))
```

```
extract_frame_times_from_photodiode(self, photodiode_cycle=60, frame_keys=('frames',
                                'stim_vsync'), photodiode_keys=('photodiode',
                                'stim_photodiode'))
```

```
extract_frame_times_from_vsyncs(self, photodiode_cycle=60, frame_keys=('frames',
                                'stim_vsync'), photodiode_keys=('photodiode',
                                'stim_photodiode'))
```

```
extract_led_times(self, keys=('LED_sync', 'opto_trial'), fallback_line=18)
```

```
classmethod factory(path)
    Build a new SyncDataset.
```

Parameters

path [str] Filesystem path to the h5 file containing sync information to be loaded.

sample_frequency

allensdk.brain_observatory.ecephys.file_io.stim_file module

```
class allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile(data,
                                                                                       **kwargs)
```

Bases: object

angular_wheel_rotation

Extract the total rotation of the running wheel on each frame.

angular_wheel_velocity

Extract the mean angular velocity of the running wheel (degrees / s) for each frame.

classmethod factory (*path*, ***kwargs*)**frames_per_second**

Framerate of stimulus presentation

pre_blank_sec

Time (s) before initial stimulus presentation

stimuli

List of dictionaries containing information about individual stimuli

vin**vsig**

Running speed signal voltage

Module contents

allensdk.brain_observatory.ecephys.lfp_subsampling package

Submodules

allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling module

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.remove_lfp_noise` (*lfp*, *surface_channel*, *channel_numbers*, *channel_max=384*, *channel_limit=380*)

Subtract mean of channels out of brain to remove noise

lfp [numpy.ndarray] 2D array of LFP values (time x channels)

surface_channel [int] Surface channel (relative to original probe)

channel_numbers [numpy.ndarray] Channel numbers in 'lfp' array (relative to original probe)

Returns:

lfp_noise_removed [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.remove_lfp_offset` (*lfp*, *sampling_frequency*, *cut_off_frequency*, *filter_order*)

High-pass filters LFP data to remove offset

lfp [numpy.ndarray] 2D array of LFP values (time x channels)

sampling_frequency [float] Sampling frequency in Hz

cutoff_frequency [float] Cutoff frequency for highpass filter

filter_order [int] Butterworth filter order

Returns:

lfp_filtered [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.select_channels` (*total_channels*,
sur-
face_channel,
sur-
face_padding,
start_channel_offset,
chan-
nel_stride,
chan-
nel_order,
noisy_channels,
dtype=float64),
re-
move_noisy_channels,
ref-
er-
ence_channels,
dtype=float64),
re-
move_references)

Selects a subset of channels for spatial downsampling

total_channels [int] Number of channels in the original data file

surface_channel [int] Index of channel at brain surface

surface_padding [int] Number of channels above surface to save

start_channel_offset [int] First channel to save

channel_stride [int] Number of channels to skip in output

channel_order [np.ndarray] Actual order of LFP channels (needed to account for the bug in NPX extraction)

noisy_channels [numpy.ndarray] Array indicating noisy channels

remove_noisy_channels [bool] Flag to remove noisy channels

reference_channels [numpy.ndarray] Array indicating reference channels

remove_references [bool] Flag to remove reference channels

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.subsample_lfp` (*lfp_raw*,
se-
lected_channels,
sub-
sam-
pling_factor)

Subsamples LFP data

lfp_raw [numpy.ndarray] 2D array of LFP values (time x channels)

selected_channels [numpy.ndarray] Indices of channels to select (spatial subsampling)

downsampling_factor [int] Factor by which to subsample in time

Returns:

lfp_subsampled [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.subsample_timestamps` (*timestamp*
sub-
sam-
pling_fac)

Subsamples an array of timestamps

timestamps [numpy.ndarray] 1D array of timestamp values

downsampling_factor [int] Factor by which to subsample the timestamps

Returns:

timestamps_sub [numpy.ndarray] New 1D array of timestamps

Module contents

`allensdk.brain_observatory.ecephys.nwb` package

Submodules

`allensdk.brain_observatory.ecephys.nwb.ecephys_nwb_extension_builder` module

Module contents

`allensdk.brain_observatory.ecephys.optotagging_table` package

Module contents

`allensdk.brain_observatory.ecephys.stimulus_analysis` package

Submodules

`allensdk.brain_observatory.ecephys.stimulus_analysis.dot_motion` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.drifting_gratings` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.flashes` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.natural_movies` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.natural_scenes` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.receptive_field_mapping` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.static_gratings` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.stimulus_analysis` module

Module contents

`allensdk.brain_observatory.ecephys.stimulus_table` package

Subpackages

`allensdk.brain_observatory.ecephys.stimulus_table.visualization` package

Submodules

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks` module

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.build_colormap`

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.get_blocks` (*table*)

```
allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.main(table_csv_path)
allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.plot_blocks(block_id, color, norm)
```

Module contents

Submodules

allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes module

Created on Fri Dec 16 15:11:23 2016

@author: Xiaoxuan Jia

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.apply_display_sequence(sweep_frames_table, frame_display_sequence, frame_times)
    Adjust raw sweep frames for a stimulus based on the display sequence for that stimulus.
```

Parameters

sweep_frames_table [pd.DataFrame] Each row is a sweep. Has two columns, 'start' and 'end', which describe (in frames) when that sweep began and ended.

frame_display_sequence [np.ndarray] 2D array. Rows are display intervals. The 0th column is the start frame of that interval, the 1st the end frame.

Returns

sweep_frames_table [pd.DataFrame] As above, but start and end frames have been adjusted based on the display sequence.

Notes

The frame values in the raw sweep_frames_table are given in 0-indexed offsets from the start of display for this stimulus. This domain only takes into account frames which are part of a display interval for that stimulus, so the frame ids need to be adjusted to lie on the global frame sequence.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.apply_frame_times(stimulus_table, frame_times, sweep_frames_table, extra_frames_table, map_color='End'))
    Converts sweep times from frames to seconds.
```

Parameters

stimulus_table [pd.DataFrame] Rows are sweeps. Columns are stimulus parameters as well as start and end frames for each sweep.

frame_times [numpy.ndarray] Gives the time in seconds at which each frame (indices) began.

frames_per_second [numeric, optional] If provided, and `extra_frame_time` is True, will be used to calculate the `extra_frame_time`.

extra_frame_time [float, optional] If provided, an additional frame time will be appended. The time will be incremented by `extra_frame_time` from the previous last frame time, to denote the time at which the last frame ended. If False, no extra time will be appended. If None (default), the increment will be 1.0/fps.

map_columns [tuple of str, optional] Which columns to replace with times. Defaults to 'Start' and 'End'

Returns

stimulus_table [pd.DataFrame] As above, but with `map_columns` values converted to seconds from frames.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.assign_sweep_values(stim_
swee
on='
drop:
tmp_
```

Left joins a stimulus table to a sweep table in order to associate epochs in time with stimulus characteristics.

Parameters

stim_table [pd.DataFrame] Each row is a stimulus epoch, with start and end times and a foreign key onto a particular sweep.

sweep_table [pd.DataFrame] Each row is a sweep. Should have columns in common with the `stim_table` - the resulting table will use values from the `sweep_table`.

on [str, optional] Column on which to join.

drop [bool, optional] If True (default), the join column (argument `on`) will be dropped from the output.

tmp_suffix [str, optional] Will be used to identify overlapping columns. Should not appear in the name of any column in either dataframe.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.build_stimuluswise_table
```

Construct a table of sweeps, including their times on the experiment-global clock and the values of each relevant parameter.

Parameters

stimulus [dict] Describes presentation of a stimulus on a particular experiment. Has a number of fields, of which we are using:

stim_path [str] windows file path to the stimulus data

sweep_frames [list of lists] rows are sweeps, columns are start and end frames of that sweep (in the stimulus-specific frame domain). C-order.

sweep_order [list of int] indices are frames, values are the sweep on that frame

display_sequence [list of list]

rows are intervals in which the stimulus was displayed. Columns are start and end times (s, global) of the display. C-order.

dimnames [list of str] Names of parameters for this stimulus (such as “Contrast”)

sweep_table [list of tuple] Each element is a tuple of parameter values (1 per dimname) describing a single sweep.

seconds_to_frames [function] Converts experiment seconds to frames

start_key [str, optional] key to use for start frame indices. Defaults to ‘Start’

end_key [str, optional] key to use for end frame indices. Defaults to ‘End’

name_key [str, optional] key to use for stimulus name annotations. Defaults to ‘stimulus_name’

block_key [str, optional] key to use for the 0-index position of this stimulus block

get_stimulus_name [function | dict -> str, optional] extracts stimulus name from the stimulus dictionary. Default is read_stimulus_name_from_path

Returns

list of pandas.DataFrame : Each table corresponds to an entry in the display sequence. Rows are sweeps, columns are stimulus parameter values as well as “Start” and “End”.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.create_stim_table (stimuli, stim-  
u-  
lus_table, spon-  
ta-  
neous_activity_table, sort_key, block_key, in-  
dex_key)
```

Build a full stimulus table

Parameters

stimuli [list of dict] Each element is a stimulus dictionary, as provided by the stim.pkl file.

stimulus_table [function] A function which takes a single stimulus dictionary as its argument and returns a stimulus table dataframe.

spontaneous_activity_table [function] A function which takes a list of stimulus tables as arguments and returns a list of 0 or more tables describing spontaneous activity sweeps.

sort_key [str, optional] Sort the final stimulus table in ascending order by this key. Defaults to 'Start'.

Returns

stim_table_full [pandas.DataFrame] Each row is a sweep. Has columns describing (in frames) the start and end times of each sweep. Other columns describe the values of stimulus parameters on those sweeps.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.make_spontaneous_activity_table
```

Fills in frame gaps in a set of stimulus tables. Suitable for use as the spontaneous_activity_table in create_stim_table.

Parameters

stimulus_tables [list of pd.DataFrame] Input tables - should have start_key and end_key columns.

start_key [str, optional] Column name for the start of a sweep. Defaults to 'Start'.

end_key [str, optional] Column name for the end of a sweep. Defaults to 'End'.

duration_threshold [numeric or None] If not None (default is 0), remove spontaneous activity sweeps whose duration is less than this threshold.

Returns

list : Either empty, or contains a single pd.DataFrame. The rows of the dataframe are spontaneous activity sweeps.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.read_stimulus_name_from_filename
```

Obtains a human-readable stimulus name by looking at the filename of the 'stim_path' item.

Parameters

stimulus [dict] must contain a 'stim_path' item.

Returns

str : name of stimulus

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.split_column(table, column, new_columns, drop_old=True)
```

Divides a dataframe column into multiple columns.

Parameters

table [pandas.DataFrame] Columns will be drawn from and assigned to this dataframe. This dataframe will NOT be modified inplace.

column [str] This column will be split.

new_columns [dict, mapping strings to functions] Each key will be the name of a new column, while its value (a function) will be used to build the new column's values. The functions should map from a single value of the original column to a single value of the new column.

drop_old [bool, optional] If True, the original column will be dropped from the table.

Returns

table [pd.DataFrame] The modified table

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities module

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**add_number_to_shuffled_r**

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**collapse_columns**(table)
merge, where possible, columns that describe the same parameter. This is pretty conservative - it only matches columns by capitalization and it only overrides nans.

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**drop_empty_columns**(table)
Remove from the stimulus table columns whose values are all nan

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**map_column_names**(table,
name_ma
ig-
nore_cas

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**map_stimulus_names**(table,
name_
stim_c

Applies a mapping to the stimulus names in a stimulus table

Parameters

table [pd.DataFrame] the input stimulus table

name_map [dict, optional] rename the stimuli according to this mapping

stim_colname: str, optional look in this column for stimulus names


```
allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.standardize_movie_numbers
```

Natural movie stimuli in visual coding are numbered using words, like “natural_movie_two” rather than “natural_movie_2”. This function ensures that all of the natural movie stimuli in an experiment are named by that convention.

Parameters

- table** [pd.DataFrame] the incoming stimulus table
- movie_re** [re.Pattern, optional] regex that matches movie stimulus names
- numeral_re** [re.Pattern, optional] regex that extracts movie numbers from stimulus names
- digit_names** [dict, optional] map from numerals to english words
- stim_colname** [str, optional] the name of the dataframe column that contains stimulus names

Returns

- table** [pd.DataFrame] the stimulus table with movie numerals having been mapped to english words

allensdk.brain_observatory.ecephys.stimulus_table.output_validation module

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_epoch_duration
```

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_epoch_order (table,
time_start, time_end,
'EPOCH_DURATION')
```

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_max_spontaneous
```

allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction module

```
allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.extract_com
```

Parameters which are not set as sweep_params in the stimulus script (usually because they are not varied during the course of the session) are not output in an easily machine-readable format. This function attempts to recover them by parsing the string repr of the stimulus.

Parameters

stim_repr [str]

The repr of the camstim stimulus object. Served up per-stimulus in the stim pickle.

repr_params_re [re.Pattern] Extracts attributes as “=”-seperated strings

array_re [re.Pattern] Extracts list reprs from numpy array reprs.

Returns

repr_params [dict] dictionary of parameter keys and values extracted from the stim repr.
Where possible, the values are converted to native Python types.

```
allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.extract_st
```

```
allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.parse_stim
```

Read the string representation of a psychopy stimulus and extract stimulus parameters.

Parameters

stim_repr [str]

drop_params [tuple]

repr_params_re [re.Pattern]

array_re [re.Pattern]

Returns

dict : maps extracted parameter names to values

Module contents

allensdk.brain_observatory.ecephys.visualization package

Module contents

```
allensdk.brain_observatory.ecephys.visualization.plot_mean_waveforms(mean_waveforms,  
                                                                    unit_ids,  
                                                                    peak_channels)
```

Utility for plotting mean waveforms on each unit's peak channel

Parameters

mean_waveforms [dictionary] Maps unit ids to channelwise average spike waveforms for those units

unit_ids [array-like] unique integer identifiers for units to be included

```
allensdk.brain_observatory.ecephys.visualization.plot_spike_counts(data_array,  
                                                                    time_coords,  
                                                                    cbar_label,  
                                                                    title, xlabel=  
                                                                    'time  
                                                                    relative to  
                                                                    stimulus  
                                                                    onset  
                                                                    (s)', ylabel=  
                                                                    'unit',  
                                                                    xtick_step=20)
```

Utility for making a simple spike counts plot.

Parameters

data_array [xarray.DataArray] 2D data array unitwise values per time bin. See EcephysSession.sweepwise_spike_counts

```
allensdk.brain_observatory.ecephys.visualization.raster_plot(spike_times,  
                                                                figsize=(8, 8),  
                                                                cmap=<matplotlib.colors.ListedColormap  
                                                                object at  
                                                                0x7f5f14378750>,  
                                                                title='spike  
                                                                raster', cycle_  
                                                                colors=False)
```

allensdk.brain_observatory.ecephys.write_nwb package

Module contents

Submodules

allensdk.brain_observatory.ecephys.ecephys_project_cache module

allensdk.brain_observatory.ecephys.ecephys_session module

allensdk.brain_observatory.ecephys.stimulus_sync module

```
allensdk.brain_observatory.ecephys.stimulus_sync.allocate_by_vsync(vs_diff,  
                                                                    index,  
                                                                    starts,  
                                                                    ends,  
                                                                    frame_duration,  
                                                                    irregularity, cycle)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.assign_to_last(index,  
                                                                    starts, ends,  
                                                                    frame_duration,  
                                                                    irregularity,  
                                                                    cycle)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.compute_frame_times(photodiode_times,
                                                                    frame_duration,
                                                                    num_frames,
                                                                    cycle,
                                                                    irregular_interval_policy=<function>
                                                                    as-
                                                                    sign_to_last
                                                                    at
                                                                    0x7f5ee6c00050>)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.correct_on_off_effects(pd_times)
```

Notes

This cannot (without additional info) determine whether an assymmetric offset is odd-long or even-long.

```
allensdk.brain_observatory.ecephys.stimulus_sync.estimate_frame_duration(pd_times,
                                                                    cycle
                                                                    cle=60)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.fix_unexpected_edges(pd_times,
                                                                    ndevs=10,
                                                                    cycle
                                                                    cle=60,
                                                                    max_frame_offset=4)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.flag_unexpected_edges(pd_times,
                                                                    ndevs=10)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.trim_border_pulses(pd_times,
                                                                    vs_times,
                                                                    frame_interval=0.016666666666666666
                                                                    num_frames=5)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.trimmed_stats(data, pc-
                                                                    tiles=(10, 90))
```

Module contents

```
allensdk.brain_observatory.ecephys.get_unit_filter_value(key, pop=True, re-
                                                                    place_none=True,
                                                                    **source)
```

allensdk.brain_observatory.extract_running_speed package

Module contents

allensdk.brain_observatory.gaze_mapping package

Module contents

allensdk.brain_observatory.nwb package

Subpackages

`allensdk.brain_observatory.nwb.eye_tracking` package

Submodules

`allensdk.brain_observatory.nwb.eye_tracking.extension_builder` module

`allensdk.brain_observatory.nwb.eye_tracking.ndx_ellipse_eye_tracking` module

Module contents

Submodules

`allensdk.brain_observatory.nwb.behavior_ophys_nwb_extension_builder` module

`allensdk.brain_observatory.nwb.metadata` module

`allensdk.brain_observatory.nwb.nwb_api` module

`allensdk.brain_observatory.nwb.nwb_utils` module

`allensdk.brain_observatory.nwb.schemas` module

Module contents

`allensdk.brain_observatory.ophys` package

Subpackages

`allensdk.brain_observatory.ophys.trace_extraction` package

Module contents

Module contents

`allensdk.brain_observatory.receptive_field_analysis` package

Submodules

`allensdk.brain_observatory.receptive_field_analysis.chisquarperf` module

`allensdk.brain_observatory.receptive_field_analysis.chisquarperf.NLL_to_pvalue` (*NLLs*,
log_base=10.0)

`allensdk.brain_observatory.receptive_field_analysis.chisquareref.build_trial_matrix` (*LSN_template*, *num_trials*, *on_off_lum*, 0))

Construct indicator arrays for on/off pixels across trials.

Parameters

LSN_template [np.ndarray] Dimensions are (nTrials, nYPixels, nXPixels). Luminance values per pixel and trial. The size of the first dimension may be larger than the `num_trials` argument (in which case only the first `num_trials` slices will be used) but may not be smaller.

num_trials [int] The number of trials (left-justified) to build indicators for.

on_off_luminance [array-like, optional] The zeroth element is the luminance value of a pixel when on, the first when off. Defaults are [255, 0].

Returns

trial_mat [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}, nTrials). Boolean values indicate that a pixel was on/off on a particular trial.

`allensdk.brain_observatory.receptive_field_analysis.chisquareref.chi_square_binary` (*events*, *LSN_template*)

`allensdk.brain_observatory.receptive_field_analysis.chisquareref.chi_square_within_mask` (*exclusion_mask*, *events_per_pixel*, *trials_per_pixel*, *exclusion_mask*)

Determine if cells respond preferentially to on/off pixels in a mask using a chi2 test.

Parameters

exclusion_mask [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer indicator for INCLUSION (!) of a pixel within the testing region.

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Integer values are response counts by cell to on/off luminance at each pixel.

trials_per_pixel [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer values are counts of trials where a pixel is on/off.

Returns

p_vals [np.ndarray] One-dimensional, of length nCells. Float values are p-values for the hypothesis that a given cell has a receptive field within the exclusion mask.

chi [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Values (float) are squared residual event counts divided by expected event counts.

`allensdk.brain_observatory.receptive_field_analysis.chisquareref.deinterpolate_RF` (*rf_map*, *x_pts*, *y_pts*, *deg_per_pt*)

Downsample an image

Parameters

rf_map [np.ndarray] Input image

x_pts [np.ndarray] Count of sample points along the first (column) axis

y_pts [np.ndarray] Count of sample points along the zeroth (row) axis

deg_per_pnt [numeric] scale factor

Returns

sampled_yx [np.ndarray] Downsampled image

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.get_disc_masks` (*LSN_template*,
radius,
on_luminance=255,
off_luminance=0)

Obtain an indicator mask surrounding each pixel. The mask is a square, excluding pixels which are coactive on any trial with the main pixel.

Parameters

LSN_template [np.ndarray] Dimensions are (nTrials, nYPixels, nXPixels). Luminance values per pixel and trial.

radius [int] The base mask will be a box whose sides are $2 * \text{radius} + 1$ in length.

on_luminance [int, optional] The value of the luminance for on trials. Default is 255

off_luminance [int, optional] The value of the luminance for off trials. Default is 0

Returns

masks [np.ndarray] Dimensions are (nYPixels, nXPixels, nYPixels, nXPixels). The first 2 dimensions describe the pixel from which the mask was computed. The last 2 serve as the dimensions of the mask images themselves. Masks are binary arrays of type float, with 1 indicating inside, 0 outside.

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.get_events_per_pixel` (*responses*,
trial_matrix)

Obtain a matrix linking cellular responses to pixel activity.

Parameters

responses_np [np.ndarray] Dimensions are (nTrials, nCells). Boolean values indicate presence/absence of a response on a given trial.

trial_matrix [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}, nTrials). Boolean values indicate that a pixel was on/off on a particular trial.

Returns

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Values for each cell, pixel, and on/off state are the sum of events for that cell across all trials where the pixel was in the on/off state.

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.get_expected_events_by_pixel`

Calculate expected number of events per pixel

Parameters

exclusion_mask [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer indicator for INCLUSION (!) of a pixel within the testing region.

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Integer values are response counts by cell to on/off luminance at each pixel.

trials_per_pixel [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer values are counts of trials where a pixel is on/off.

Returns

np.ndarray : Dimensions (nCells, nYPixels, nXPixels, {on, off}). Float values are pixelwise counts of events expected if events are evenly distributed in mask across trials.

```
allensdk.brain_observatory.receptive_field_analysis.chisquareref.get_peak_significance(chi_sq,
                                                                                   LSN_t,
                                                                                   al-
                                                                                   pha=0)
```

allensdk.brain_observatory.receptive_field_analysis.chisquareref.interpolate_RF(*rf_map*,
deg_per_pnt)

Upsample an image

Parameters

rf_map [np.ndarray] Input image
deg_per_pnt [numeric] scale factor

Returns

interpolated [np.ndarray] Upsampled image

```
allensdk.brain_observatory.receptive_field_analysis.chisquareref.locate_median(y,
                                                                                   x)
```

allensdk.brain_observatory.receptive_field_analysis.chisquareref.pvalue_to_NLL(*p_values*,
max_NLL=10.0)

allensdk.brain_observatory.receptive_field_analysis.chisquareref.smooth_STA(*STA*,
gauss_std=0.75,
to-
tal_degrees=64)

Smooth an image by convolution with a gaussian kernel

Parameters

STA [np.ndarray] Input image
gauss_std [numeric, optional] Standard deviation of the gaussian kernel. Will be applied to the upsampled image, so units are visual degrees. Default is 0.75
total_degrees [int, optional] Size in visual degrees of the input image along its zeroth (row) axis. Used to set the scale factor for up/downsampling.

Returns

STA_smoothed [np.ndarray] Smoothed image

allensdk.brain_observatory.receptive_field_analysis.eventdetection module

```
allensdk.brain_observatory.receptive_field_analysis.eventdetection.detect_events(data,
                                                                                   cell_index,
                                                                                   stim-
                                                                                   u-
                                                                                   lus,
                                                                                   de-
                                                                                   bug_plots=Fa)
```

allensdk.brain_observatory.receptive_field_analysis.fit_parameters module

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.add_to_fit_parameters_d`

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.compute_distance` (*center_on*,
center_off)

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.compute_overlap` (*data_fitted*,
data_fitted)

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.get_gaussian_fit_single`

allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D module

exception `allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.GaussianFitError`
Bases: `RuntimeError`

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.fitgaussian2D` (*data*)
Fit a 2D gaussian to an image

Parameters

data [np.ndarray] input image

Returns

p2 [list] height row mean column mean row standard deviation column standard deviation
rotation

Notes

see `gaussian2D` for details about output values

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.gaussian2D` (*height*,
center_x,
center_y,
width_x,
width_y,
rotation)

Build a function which evaluates a scaled 2d gaussian pdf

Parameters

height [float] scale factor

center_x [float] first coordinate of mean

center_y [float] second coordinate of mean

width_x [float] standard deviation along x axis

width_y [float] standard deviation along y axis

rotation [float] degrees clockwise by which to rotate the gaussian

Returns

rotgauss: fn parameters are x and y positions (row/column semantics are set by your inputs to this function). Return value is the scaled gaussian pdf evaluated at the argued point.

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.moments2(data)`
Treating input image data as an independent multivariate gaussian, estimate mean and standard deviations

Parameters

data [np.ndarray] 2d numpy array.

Returns

height [float] The maximum observed value in the data

y [float] Mean row index

x [float] Mean column index

width_y [float] The standard deviation along the mean row

width_x [float] The standard deviation along the mean column

None : This function returns an instance of None.

Notes

uses original method from website for finding center

allensdk.brain_observatory.receptive_field_analysis.postprocessing module

`allensdk.brain_observatory.receptive_field_analysis.postprocessing.get_gaussian_fit(rf)`

`allensdk.brain_observatory.receptive_field_analysis.postprocessing.run_postprocessing(data, rf)`

allensdk.brain_observatory.receptive_field_analysis.receptive_field module

`allensdk.brain_observatory.receptive_field_analysis.receptive_field.compute_receptive_field`

`allensdk.brain_observatory.receptive_field_analysis.receptive_field.compute_receptive_field`

```
allensdk.brain_observatory.receptive_field_analysis.receptive_field.events_to_pvalues_no_f
```

```
allensdk.brain_observatory.receptive_field_analysis.receptive_field.get_attribute_dict(rf)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.print_summary(rf)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.read_h5_group(g)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.read_receptive_field_f
allensdk.brain_observatory.receptive_field_analysis.receptive_field.write_receptive_field_t
```

allensdk.brain_observatory.receptive_field_analysis.tools module

```
allensdk.brain_observatory.receptive_field_analysis.tools.dict_generator(indict,
                                                                           pre=None)
allensdk.brain_observatory.receptive_field_analysis.tools.list_of_dicts_to_dict_of_lists(lis
allensdk.brain_observatory.receptive_field_analysis.tools.read_h5_group(g)
```

allensdk.brain_observatory.receptive_field_analysis.utilities module

```
allensdk.brain_observatory.receptive_field_analysis.utilities.convolve(img,
                                                                           sigma=4)
    2D Gaussian convolution
allensdk.brain_observatory.receptive_field_analysis.utilities.get_A(data,
                                                                           stimu-
                                                                           lus)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_A_blur(data,
                                                                           stim-
                                                                           u-
                                                                           lus)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_attribute_dict(rf)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_components(receptive_field_data
allensdk.brain_observatory.receptive_field_analysis.utilities.get_shuffle_matrix(data,
                                                                           event_vector,
                                                                           A,
                                                                           num-
                                                                           ber_of_shuffle
                                                                           re-
                                                                           sponse_detect
```

```
allensdk.brain_observatory.receptive_field_analysis.utilities.get_sparse_noise_epoch_mask_
```

```
allensdk.brain_observatory.receptive_field_analysis.utilities.smooth(x, win-
                                                                    dow_len=11,
                                                                    win-
                                                                    dow='hanning',
                                                                    mode='valid')
```

smooth the data using a window with requested size.

This method is based on the convolution of a scaled window with the signal. The signal is prepared by introducing reflected copies of the signal (with the window size) in both ends so that transient parts are minimized in the beginning and end part of the output signal.

input: x: the input signal window_len: the dimension of the smoothing window; should be an odd integer
window: the type of window from 'flat', 'hanning', 'hamming', 'bartlett', 'blackman'

flat window will produce a moving average smoothing.

output: the smoothed signal

example:

```
t=linspace(-2,2,0.1) x=sin(t)+randn(len(t))*0.1 y=smooth(x)
```

see also:

numpy.hanning, numpy.hamming, numpy.bartlett, numpy.blackman, numpy.convolve scipy.signal.lfilter

TODO: the window parameter could be the window itself if an array instead of a string NOTE: length(output) != length(input), to correct this: return y[(window_len/2-1):-(window_len/2)] instead of just y.

```
allensdk.brain_observatory.receptive_field_analysis.utilities.upsample_image_to_degrees (img
```

allensdk.brain_observatory.receptive_field_analysis.visualization module

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_chi_square_summary (n
```

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_ellipses (gaussian_fit_di
                                                                    ax=None,
                                                                    show=True,
                                                                    close=True,
                                                                    save_file_name
                                                                    color='b')
```

Example Usage: oeid, cell_index, stimulus = 512176430, 12, 'locally_sparse_noise' brain_observatory_cache = BrainObservatoryCache() data_set = brain_observatory_cache.get_ophys_experiment_data(oeid) lsn = LocallySparseNoise(data_set, stimulus) result = compute_receptive_field_with_postprocessing(data_set, cell_index, stimulus, alpha=.05, number_of_shuffles=5000) plot_ellipses(result['off']['gaussian_fit'], color='r')

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_fields(on_data,
                                                                              off_data,
                                                                              on_axes,
                                                                              off_axes,
                                                                              cbar_axes=None,
                                                                              clim=None,
                                                                              cmap='magma')

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_gaussian_fit(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=
                                                                              cmap='m

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_mask(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=None,
                                                                              cmap='magma')

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_msr_summary(lsn,
                                                                              cell_index,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=N
                                                                              cmap=Non

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_p_values(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=None,
                                                                              cmap='magma

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_receptive_field_data

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_rts_blur_summary(rf_a
                                                                              ax_
                                                                              ax_
                                                                              ax_
                                                                              cma

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_rts_summary(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=N
                                                                              cmap='ma

allensdk.brain_observatory.receptive_field_analysis.visualization.pvalue_to_NLL(p_values,
                                                                              max_NLL=10.0
```

Module contents

allensdk.brain_observatory.sync_utilities package

Module contents

```
allensdk.brain_observatory.sync_utilities.get_synchronized_frame_times(session_sync_file:
                                                                    path-
                                                                    lib.Path,
                                                                    sync_line_label_keys:
                                                                    Tu-
                                                                    ple[str,
                                                                    ...],
                                                                    trim_after_spike:
                                                                    bool
                                                                    =
                                                                    True)
                                                                    →
                                                                    pan-
                                                                    das.core.series.Series
```

Get experimental frame times from an experiment session sync file.

Parameters

session_sync_file [Path] Path to an ephys session sync file. The sync file contains rising/falling edges from a daq system which indicates when certain events occur (so they can be related to each other).

sync_line_label_keys [Tuple[str, ...]] Line label keys to get times for. See class attributes of `allensdk.brain_observatory.sync_dataset.Dataset` for a listing of possible keys.

trim_after_spike [bool = True] If True, will call `trim_discontiguous_times` on the frame times before returning them, which will detect any spikes in the data and remove all elements for the list which come after the spike.

Returns

pd.Series An array of times when frames for the eye tracking camera were acquired.

```
allensdk.brain_observatory.sync_utilities.trim_discontiguous_times(times,
                                                                    thresh-
                                                                    old=100)
```

allensdk.brain_observatory.visualization package

Module contents

```
allensdk.brain_observatory.visualization.plot_running_speed(timestamps,      val-
                                                                    ues,
                                                                    start_index=0,
                                                                    stop_index=None,
                                                                    step=1,
                                                                    yla-
                                                                    bel='running speed
                                                                    (cm/s)', xlabel='time
                                                                    (s)', title=None)
```

Make a simple plot of a running speed trace

Parameters

timestamps [numpy.ndarray] Times at which running speed samples were collected

values [numpy.ndarray] Running speed values (by default: linear cm / s with negative values indicating backwards movement)

Submodules

allensdk.brain_observatory.argschema_utilities module

```
class allensdk.brain_observatory.argschema_utilities.ArgSchemaParserPlus (*args,  
                                                                    **kwargs)  
    Bases: argschema.argschema_parser.ArgSchemaParser
```

```
class allensdk.brain_observatory.argschema_utilities.InputFile (default=<marshmallow.missing>,  
                                                                    attribute=None,  
                                                                    data_key=None,  
                                                                    error=None,  
                                                                    validate=None,  
                                                                    required=False,  
                                                                    al-  
                                                                    low_none=None,  
                                                                    load_only=False,  
                                                                    dump_only=False,  
                                                                    miss-  
                                                                    ing=<marshmallow.missing>,  
                                                                    er-  
                                                                    ror_messages=None,  
                                                                    **metadata)  
    Bases: marshmallow.fields.String
```

A marshmallow String field subclass which deserializes json str fields that represent a desired input path to pathlib.Path. Also performs read access checking.

```
class allensdk.brain_observatory.argschema_utilities.OutputFile (default=<marshmallow.missing>,  
                                                                    at-  
                                                                    tribute=None,  
                                                                    data_key=None,  
                                                                    error=None,  
                                                                    vali-  
                                                                    date=None,  
                                                                    re-  
                                                                    quired=False,  
                                                                    al-  
                                                                    low_none=None,  
                                                                    load_only=False,  
                                                                    dump_only=False,  
                                                                    miss-  
                                                                    ing=<marshmallow.missing>,  
                                                                    er-  
                                                                    ror_messages=None,  
                                                                    **metadata)  
    Bases: marshmallow.fields.String
```

A marshmallow String field subclass which deserializes json str fields that represent a desired output file path to a pathlib.Path. Also performs write access checking.


```

class allensdk.brain_observatory.argschema_utilities.RaisingSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)

    Bases: argschema.schemas.DefaultSchema

    class Meta
        Bases: object

        unknown = 'raise'

        opts = <marshmallow.schema.SchemaOpts object>

allensdk.brain_observatory.argschema_utilities.check_read_access (path)
allensdk.brain_observatory.argschema_utilities.check_write_access (filepath, al-
                                                                    low_exists=False)
allensdk.brain_observatory.argschema_utilities.check_write_access_dir (dirpath)
allensdk.brain_observatory.argschema_utilities.check_write_access_overwrite (path)
allensdk.brain_observatory.argschema_utilities.optional_lims_inputs (argv, in-
                                                                    put_schema,
                                                                    out-
                                                                    put_schema,
                                                                    lims_input_getter)
allensdk.brain_observatory.argschema_utilities.write_or_print_outputs (data,
                                                                    parser)

```

allensdk.brain_observatory.brain_observatory_exceptions module

```

exception allensdk.brain_observatory.brain_observatory_exceptions.BrainObservatoryAnalysisException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.EpochSeparationException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.MissingStimulusException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.NoEyeTrackingException
    Bases: Exception

```

allensdk.brain_observatory.brain_observatory_plotting module

```

allensdk.brain_observatory.brain_observatory_plotting.plot_drifting_grating_traces (dg,
                                                                    save_dir)
    saves figures with a Ori X TF grid of mean responses

```

```
allensdk.brain_observatory.brain_observatory_plotting.plot_lsn_traces (lsn,  
                                                                    save_dir,  
                                                                    suf-  
                                                                    fix="")  
allensdk.brain_observatory.brain_observatory_plotting.plot_ns_traces (nsa,  
                                                                    save_dir)  
allensdk.brain_observatory.brain_observatory_plotting.plot_running_a (dg,  
                                                                    nm1,  
                                                                    nm3,  
                                                                    save_dir)  
allensdk.brain_observatory.brain_observatory_plotting.plot_sg_traces (sg,  
                                                                    save_dir)
```

allensdk.brain_observatory.chisquare_categorical module

Created on Wed Jun 5 15:52:22 2019

@author: dan

```
allensdk.brain_observatory.chisquare_categorical.advance_combination (curr_combination,  
                                                                    op-  
                                                                    tions_per_column)  
allensdk.brain_observatory.chisquare_categorical.chisq_from_stim_table (stim_table,  
                                                                    columns,  
                                                                    mean_sweep_events,  
                                                                    num_shuffles=1000,  
                                                                    ver-  
                                                                    bose=False)  
allensdk.brain_observatory.chisquare_categorical.compute_chi (observed,      ex-  
                                                                    pected)  
allensdk.brain_observatory.chisquare_categorical.compute_chi_shuffle (mean_sweep_events,  
                                                                    sweep_categories,  
                                                                    num_shuffles=1000)  
allensdk.brain_observatory.chisquare_categorical.compute_expected (mean_sweep_events,  
                                                                    sweep_conditions)  
allensdk.brain_observatory.chisquare_categorical.compute_observed (mean_sweep_events,  
                                                                    sweep_conditions)  
allensdk.brain_observatory.chisquare_categorical.make_category_dummy (sweep_categories)  
allensdk.brain_observatory.chisquare_categorical.stim_table_to_categories (stim_table,  
                                                                    columns,  
                                                                    ver-  
                                                                    bose=False)
```

allensdk.brain_observatory.circle_plots module

```
class allensdk.brain_observatory.circle_plots.CoronaPlotter (angle_start=270,  
                                                                    plot_scale=1.2,  
                                                                    inner_radius=0.3,  
                                                                    *args, **kwargs)
```

Bases: `allensdk.brain_observatory.circle_plots.PolarPlotter`

```

infer_dims (self, category_data)

plot (self, category_data, data=None, clim=None, cmap=<matplotlib.colors.LinearSegmentedColormap
      object at 0x7f5ee6ec3690>)

set_dims (self, categories)

show_arrow (self, color=None)

show_circle (self, color=None)

class allensdk.brain_observatory.circle_plots.FanPlotter (group_scale=0.9, *args,
                                                         **kwargs)
    Bases: allensdk.brain_observatory.circle_plots.PolarPlotter

    static for_drifting_gratings ()

    static for_static_gratings ()

    infer_dims (self, r_data, angle_data, group_data)

    plot (self, r_data, angle_data, group_data=None, data=None,
          cmap=<matplotlib.colors.LinearSegmentedColormap object at 0x7f5ee6ec3690>, clim=None,
          rmap=None, rlim=None, axis_color=None, label_color=None)

    set_dims (self, rs, angles, groups)

    show_angle_labels (self, angles=None, labels=None, color=None, offset=0.05, fontdict=None)

    show_axes (self, angles=None, radii=None, closed=False, color=None)

    show_group_labels (self, groups=None, color=None, fontdict=None)

    show_r_labels (self, radii=None, labels=None, color=None, offset=0.1, fontdict=None)

class allensdk.brain_observatory.circle_plots.PolarPlotter (direction=-1, an-
                                                         gle_start=0, cir-
                                                         cle_scale=1.1, in-
                                                         ner_radius=None,
                                                         plot_center=(0.0,
                                                         0.0), plot_scale=0.9)

    Bases: object

    DIR_CCW = 1

    DIR_CW = -1

    finalize (self)

class allensdk.brain_observatory.circle_plots.TrackPlotter (direction=-1, an-
                                                         gle_start=270.0,
                                                         inner_radius=0.45,
                                                         ring_length=None,
                                                         *args, **kwargs)

    Bases: allensdk.brain_observatory.circle_plots.PolarPlotter

    plot (self, data, clim=None, cmap=<matplotlib.colors.LinearSegmentedColormap object at
          0x7f5ee6ec3690>, mean_cmap=<matplotlib.colors.LinearSegmentedColormap object at
          0x7f5ee73ac190>, norm=None)

    show_arrow (self, color=None)

allensdk.brain_observatory.circle_plots.add_angle_labels (ax, angles, labels, ra-
                                                         dius, color=None, font-
                                                         dict=None, offset=0.05)

```

```
allensdk.brain_observatory.circle_plots.add_arrow(ax, radius, start_angle, end_angle,
                                                  color=None, width=18.0)
allensdk.brain_observatory.circle_plots.angle_lines(angles, inner_radius,
                                                    outer_radius)
allensdk.brain_observatory.circle_plots.build_hex_pack(n)
allensdk.brain_observatory.circle_plots.hex_pack(radius, n)
allensdk.brain_observatory.circle_plots.make_pincushion_plot(data, trials, on,
                                                            nrows, ncols,
                                                            clim=None,
                                                            color_map=None,
                                                            radius=None)
allensdk.brain_observatory.circle_plots.polar_line_circles(radii, theta,
                                                           start_r=0)
allensdk.brain_observatory.circle_plots.polar_linspace(radius, start_angle,
                                                         stop_angle, num, end-
                                                         point=False, degrees=True)
    Evenly distributed list of x,y coordinates from an input range of angles and a radius in polar coordinates.
allensdk.brain_observatory.circle_plots.polar_to_xy(angles, radius)
    Convert an array of angles (in radians) and a radius in polar coordinates to an array of x,y coordinates.
allensdk.brain_observatory.circle_plots.radial_arcs(rs, start_theta, end_theta)
allensdk.brain_observatory.circle_plots.radial_circles(rs)
allensdk.brain_observatory.circle_plots.reset_hex_pack()
allensdk.brain_observatory.circle_plots.rings_in_hex_pack(ct)
allensdk.brain_observatory.circle_plots.spiral_trials(radii, x=0.0, y=0.0)
allensdk.brain_observatory.circle_plots.spiral_trials_polar(r, theta, radii, off-
                                                            set=None)
allensdk.brain_observatory.circle_plots.wedge_ring(N, inner_radius, outer_radius,
                                                    start=0, stop=360)
```

allensdk.brain_observatory.demixer module

```
allensdk.brain_observatory.demixer.demix_time_dep_masks(raw_traces:
                                                         numpy.ndarray,
                                                         stack: numpy.ndarray,
                                                         masks: numpy.ndarray,
                                                         max_block_size: int
                                                         = 1000) → Tuple[
                                                         numpy.ndarray, list]
```

Demix traces of potentially overlapping masks extracted from a single 2p recording.

Parameters

- **raw_traces** – 2d array of traces for each mask, of dimensions (n, t), where *t* is the number of time points and *n* is the number of masks.
- **stack** – 3d array representing a 1p recording movie, of dimensions (t, H, W) or corresponding hdf5 dataset.
- **masks** – 3d array of binary roi masks, of shape (n, H, W), where *n* is the number of masks, and HW are the dimensions of an individual frame in the movie *stack*.

Max_block_size int representing maximum number of movie frames to read at a time (-1 for full length *t* of *stack*) (the default is 1000)

Returns Tuple of demixed traces and whether each frame was skipped in the demixing calculation.

```
allensdk.brain_observatory.demixer.find_negative_baselines(trace)
allensdk.brain_observatory.demixer.find_negative_transients_threshold(trace,
                                                                    win-
                                                                    dow=500,
                                                                    length=10,
                                                                    std_devs=3)
allensdk.brain_observatory.demixer.find_zero_baselines(traces)
allensdk.brain_observatory.demixer.identify_valid_masks(mask_array)
allensdk.brain_observatory.demixer.plot_negative_baselines(raw_traces,
                                                         demix_traces,
                                                         mask_array,
                                                         roi_ids_mask,
                                                         plot_dir, ext='png')
allensdk.brain_observatory.demixer.plot_negative_transients(raw_traces,
                                                           demix_traces,
                                                           valid_roi,
                                                           mask_array,
                                                           roi_ids_mask,
                                                           plot_dir, ext='png')
allensdk.brain_observatory.demixer.plot_overlap_masks_lengthOne(roi_ind, masks,
                                                                savefile=None,
                                                                weighted=False)
allensdk.brain_observatory.demixer.plot_traces(raw_trace, demix_trace, roi_id, roi_ind,
                                                save_file)
allensdk.brain_observatory.demixer.plot_transients(roi_ind, t_trans, masks, traces,
                                                    demix_traces, savefile)
allensdk.brain_observatory.demixer.rolling_window(trace, window=500)
```

Parameters

- **trace** –
- **window** –

Returns

allensdk.brain_observatory.dff module

```
allensdk.brain_observatory.dff.calculate_dff(traces, dff_computation_cb=None,
                                              save_plot_dir=None)
```

Apply dF/F computation to a set of traces.

The default computation method is `compute_dff_windowed_median()` using default window parameters.

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

dff_computation_cb [function] Function that takes traces as an argument and returns an array of the same shape that is the calculated dF/F.

save_plot_dir [str] Directory to save dF/F plots to. By default no plots are saved.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.compute_dff_windowed_median(traces, median_kernel_long=5401, median_kernel_short=101, noise_stds=None, n_small_baseline_frames=None, **kwargs)
```

Compute dF/F of a set of traces with median filter detrending.

The operation is basically:

```
T_long = windowed_median(T) # long timescale kernel
T_dff1 = (T - T_long) / elementwise_max(T_long, noise_std(T))
T_short = windowed_median(T_dff1) # short timescale kernel
T_dff = T_dff1 - elementwise_min(T_short, 2.5*noise_std(T_dff1))
```

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

median_kernel_long [int] Window size to use for long timescale median detrending.

median_kernel_short [int] Window size to use for short timescale median detrending.

noise_stds [list] List that will contain noise_std(T_dff1) for each trace. The value for each trace will be appended to the list if provided.

n_small_baseline_frames [list] List that will contain the number of frames for each trace where the long-timescale median window is less than noise_std(T). The value for each trace will be appended to the list if provided.

kwargs: Additional keyword arguments are passed to `noise_std()`.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.compute_dff_windowed_mode(traces, mode_kernelsize=5400, mean_kernelsize=3000)
```

Compute dF/F of a set of traces using a low-pass windowed-mode operator.

The operation is basically:

```
T_mm = windowed_mean(windowed_mode(T))
T_dff = (T - T_mm) / T_mm
```

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

mode_kernelsize [int] Window size to use for windowed_mode.

mean_kernelsize [int] Window size to use for windowed_mean.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.main()
```

```
allensdk.brain_observatory.dff.movingaverage(x, kernelsize, y)
```

Compute the windowed average of an array.

Parameters

x [np.ndarray] Array to be analyzed

kernelsize [int] Size of the moving window

y [np.ndarray] Output array to store the results

```
allensdk.brain_observatory.dff.movingmode_fast(x, kernelsize, y)
```

Compute the windowed mode of an array. A running mode is initialized with a histogram of values over the initial kernelsize/2 values. The mode is then updated as the kernel moves by adding and subtracting values from the histogram.

Parameters

x [np.ndarray] Array to be analyzed

kernelsize [int] Size of the moving window

y [np.ndarray] Output array to store the results

```
allensdk.brain_observatory.dff.noise_std(x, noise_kernel_length=31, positive_peak_scale=1.5, outlier_std_scale=2.5)
```

Robust estimate of the standard deviation of the trace noise.

```
allensdk.brain_observatory.dff.plot_onetrace(dff, fc)
```

Debug plotting function

```
allensdk.brain_observatory.dff.robust_std(x)
```

Robust estimate of standard deviation.

Estimate of the standard deviation using the median absolute deviation of x.

allensdk.brain_observatory.drifting_gratings module

```
class allensdk.brain_observatory.drifting_gratings.DriftingGratings(data_set, **kwargs)
```

Bases: `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`

Perform tuning analysis specific to drifting gratings stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

static from_analysis_file (data_set, analysis_file)

get_noise_correlation (self, corr='spearman')

get_peak (self)

Computes metrics related to each cell's peak response condition.

Returns

Pandas data frame containing the following columns (_dg suffix is for drifting grating):

- ori_dg (orientation)
- tf_dg (temporal frequency)
- reliability_dg
- osi_dg (orientation selectivity index)
- dsi_dg (direction selectivity index)
- peak_dff_dg (peak dF/F)
- ptest_dg
- p_run_dg
- run_modulation_dg
- cv_dg (circular variance)

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# orientations, # temporal frequencies, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant response ($p < 0.05$) to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

number_ori

number_tf

open_star_plot (*self*, *cell_specimen_id*=None, *include_labels*=False, *cell_index*=None)

orivals

plot_direction_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_orientation_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_direction (*self*, *include_labels*=False, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_temporal_frequency (*self*, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim x repetition for noise correlations

tfvals

allensdk.brain_observatory.findlevel module

`allensdk.brain_observatory.findlevel.findlevel` (*inwave, threshold, direction='both'*)

allensdk.brain_observatory.locally_sparse_noise module

class `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` (*data_set, stimulus=None, **kwargs*)

Bases: `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`

Perform tuning analysis specific to the locally sparse noise stimulus.

Parameters

data_set: `BrainObservatoryNwbDataSet` object

stimulus: `string` Name of locally sparse noise stimulus. See `brain_observatory.stimulus_info`.

nrows: `int` Number of rows in the stimulus template

ncol: `int` Number of columns in the stimulus template

LSN

LSN_GREY = 127

LSN_OFF = 0

LSN_OFF_SCREEN = 64

LSN_ON = 255

LSN_mask

cell_index_receptive_field_analysis_data

extralength

static from_analysis_file (*data_set, analysis_file, stimulus*)

get_mean_response (*self*)

get_peak (*self*)

Implemented by subclasses.

get_receptive_field (*self*)

Calculates receptive fields for each cell

get_receptive_field_analysis_data (*self*)

Calculates receptive fields for each cell

get_receptive_field_attribute_df (*self*)

interlength

mean_response

static merge_mean_response (*rc1, rc2*)

Move out of this class, to session analysis

open_pincushion_plot (*self, on, cell_specimen_id=None, color_map=None, cell_index=None*)

```
plot_cell_receptive_field(self, on, cell_specimen_id=None, color_map=None, clim=None,
                           mask=None, cell_index=None, scalebar=True)
plot_population_receptive_field(self, color_map='RdPu', clim=None, mask=None,
                                scalebar=True)
plot_receptive_field_analysis_data(self, cell_index, **kwargs)
populate_stimulus_table(self)
    Implemented by subclasses.
static read_cell_index_receptive_field_analysis(file_handle, prefix, path=None)
receptive_field
static save_cell_index_receptive_field_analysis(cell_index_receptive_field_analysis_data,
                                                new_nwb, prefix)

sort_trials(self)
sweeplength
```

allensdk.brain_observatory.natural_movie module

```
class allensdk.brain_observatory.natural_movie.NaturalMovie(data_set,
                                                             movie_name,
                                                             **kwargs)
Bases: allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
```

Perform tuning analysis specific to natural movie stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

movie_name: string

one of [stimulus_info.NATURAL_MOVIE_ONE, stimulus_info.NATURAL_MOVIE_TWO,
stimulus_info.NATURAL_MOVIE_THREE]

```
static from_analysis_file(data_set, analysis_file, movie_name)
```

```
get_peak(self)
```

Computes properties of the peak response condition for each cell.

Returns

Pandas data frame with the below fields. A suffix of “nm1”, “nm2” or “nm3” is appended to the field name

on which of three movie clips was presented.

- peak_nm1 (frame with peak response)
- response_variability_nm1

```
get_sweep_response(self)
```

Returns the dF/F response for each cell

Returns

Numpy array

```
open_track_plot(self, cell_specimen_id=None, cell_index=None)
```

populate_stimulus_table (*self*)

Implemented by subclasses.

sweep_response

sweeplength

allensdk.brain_observatory.natural_scenes module

class allensdk.brain_observatory.natural_scenes.**NaturalScenes** (*data_set*,
***kwargs*)

Bases: *allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis*

Perform tuning analysis specific to natural scenes stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

extralength

static from_analysis_file (*data_set*, *analysis_file*)

get_noise_correlation (*self*, *corr*='spearman')

get_peak (*self*)

Computes metrics about peak response condition for each cell.

Returns

Pandas data frame with the following fields ('_ns' suffix is for natural scene):

- scene_ns (scene number)
- reliability_ns
- peak_dff_ns (peak dF/F)
- ptest_ns
- p_run_ns
- run_modulation_ns
- time_to_peak_ns

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# scenes, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant ($p < 0.05$) response to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

interlength

number_scenes

open_corona_plot (*self*, *cell_specimen_id*=None, *cell_index*=None)

plot_time_to_peak (*self*, *p_value_max*=0.05, *color_map*=<matplotlib.colors.LinearSegmentedColormap object at 0x7f5ee73e2410>)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim x repetition for noise correlations

sweeplength

allensdk.brain_observatory.observatory_plots module

```
class allensdk.brain_observatory.observatory_plots.DimensionPatchHandler (vals,  
                                                                    start_color,  
                                                                    end_color,  
                                                                    *args,  
                                                                    **kwargs)  
  
    Bases: object  
  
    dim_color (self, index)  
  
    legend_artist (self, legend, orig_handle, fontsize, handlebox)  
  
allensdk.brain_observatory.observatory_plots.figure_in_px (w, h, file_name,  
                                                            dpi=96.0, transparent=False)  
  
allensdk.brain_observatory.observatory_plots.finalize_no_axes (pad=0.0)  
allensdk.brain_observatory.observatory_plots.finalize_no_labels (pad=0.3, legend=False)  
allensdk.brain_observatory.observatory_plots.finalize_with_axes (pad=0.3)  
allensdk.brain_observatory.observatory_plots.float_label (n)  
allensdk.brain_observatory.observatory_plots.plot_cell_correlation (sig_corrs,  
                                                                    labels,  
                                                                    colors,  
                                                                    scale=15)  
  
allensdk.brain_observatory.observatory_plots.plot_combined_speed (binmed_resp_vis,  
                                                                    binmed_dx_vis,  
                                                                    binmed_resp_sp,  
                                                                    binmed_dx_sp,  
                                                                    evoked_color,  
                                                                    spont_color)  
  
allensdk.brain_observatory.observatory_plots.plot_condition_histogram (vals,  
                                                                    bins,  
                                                                    color='#ccccdd')  
  
allensdk.brain_observatory.observatory_plots.plot_mask_outline (mask, ax,  
                                                                    color='k')
```

```

allensdk.brain_observatory.observatory_plots.plot_pupil_location(xy_deg, s=1,
                                                                    c=None,
                                                                    cmap=<matplotlib.colors.LinearSegmentedColormap:
                                                                    object at
                                                                    0x7f5ee6f83f90>,
                                                                    edge-
                                                                    color="", in-
                                                                    clude_labels=True)

allensdk.brain_observatory.observatory_plots.plot_radial_histogram(angles,
                                                                    counts,
                                                                    all_angles=None,
                                                                    in-
                                                                    clude_labels=False,
                                                                    off-
                                                                    set=180.0,
                                                                    direction=-
                                                                    1,
                                                                    closed=False,
                                                                    color='#ccccdd')

allensdk.brain_observatory.observatory_plots.plot_receptive_field(rf,
                                                                    color_map=None,
                                                                    clim=None,
                                                                    mask=None,
                                                                    out-
                                                                    line_color='#cccccc',
                                                                    scale-
                                                                    bar=True)

allensdk.brain_observatory.observatory_plots.plot_representational_similarity(rs,
                                                                    dims=None,
                                                                    dim_labels=None,
                                                                    col-
                                                                    ors=None,
                                                                    dim_order=None,
                                                                    la-
                                                                    bels=True)

allensdk.brain_observatory.observatory_plots.plot_selectivity_cumulative_histogram(sis,
                                                                    xla-
                                                                    bel,
                                                                    si_range=[
                                                                    1.5],
                                                                    n_hist_bins,
                                                                    color='#cccccc')

allensdk.brain_observatory.observatory_plots.plot_speed(binned_resp, binned_dx,
                                                                    num_bins, color)

allensdk.brain_observatory.observatory_plots.plot_time_to_peak(msrs, ttps,
                                                                    t_start, t_end,
                                                                    stim_start,
                                                                    stim_end,
                                                                    cmap)

```

```
allensdk.brain_observatory.observatory_plots.population_correlation_scatter(sig_corrs,  
                                                                           noise_corrs,  
                                                                           la-  
                                                                           bels,  
                                                                           col-  
                                                                           ors,  
                                                                           scale=15)
```

allensdk.brain_observatory.r_neuropil module

```
class allensdk.brain_observatory.r_neuropil.NeuropilSubtract(lam=0.05, dt=1.0,  
                                                             folds=4)
```

Bases: object

TODO: docs

```
estimate_error(self, r)
```

Estimate error values for a given r for each fold and return the mean.

```
fit(self, r_range=[0.0, 2.0], iterations=3, dr=0.1, dr_factor=0.1)
```

Estimate error values for a range of r values. Identify a new r range around the minimum error values and repeat multiple times. TODO: docs

```
fit_block_coordinate_desc(self, r_init=5.0, min_delta_r=1e-08)
```

```
set_F(self, F_M, F_N)
```

Break the F_M and F_N traces into the number of folds specified in the class constructor and normalize each fold of F_M and R_N relative to F_N.

```
allensdk.brain_observatory.r_neuropil.ab_from_T(T, lam, dt)
```

```
allensdk.brain_observatory.r_neuropil.ab_from_diagonals(mat_dict)
```

Constructs value for scipy.linalg.solve_banded

Parameters

mat_dict: dictionary of diagonals keyed by offsets

Returns

ab: value for scipy.linalg.solve_banded

```
allensdk.brain_observatory.r_neuropil.alpha_filter(A=1.0, alpha=0.05, beta=0.25,  
                                                    T=100)
```

```
allensdk.brain_observatory.r_neuropil.error_calc(F_M, F_N, F_C, r)
```

```
allensdk.brain_observatory.r_neuropil.error_calc_outlier(F_M, F_N, F_C, r)
```

```
allensdk.brain_observatory.r_neuropil.estimate_contamination_ratios(F_M,  
                                                                       F_N,  
                                                                       lam=0.05,  
                                                                       folds=4,  
                                                                       itera-  
                                                                       tions=3,  
                                                                       r_range=[0.0,  
                                                                       2.0],  
                                                                       dr=0.1,  
                                                                       dr_factor=0.1)
```

Calculates neuropil contamination of ROI

Parameters

F_M: ROI trace F_N: Neuropil trace

Returns

dictionary: key-value pairs

- 'r': the contamination ratio – corrected trace = $M - r \cdot N$
- 'err': RMS error
- 'min_error': minimum error
- 'bounds_error': boolean. True if error or R are outside tolerance

`allensdk.brain_observatory.r_neuropil.get_diagonals_from_sparse(mat)`

Returns a dictionary of diagonals keyed by offsets

Parameters

mat: `scipy.sparse` matrix

Returns

dictionary: diagonals keyed by offsets

`allensdk.brain_observatory.r_neuropil.normalize_F(F_M, F_N)`

`allensdk.brain_observatory.r_neuropil.synthesize_F(T, af1, af2, p1=0.05, p2=0.1)`

Build a synthetic F_C, F_M, F_N, and r of length T TODO: docs

`allensdk.brain_observatory.r_neuropil.validate_with_synthetic_F(T, N)`

Compute N synthetic traces of length T with known values of r, then estimate r. TODO: docs

`allensdk.brain_observatory.roi_masks` module

class `allensdk.brain_observatory.roi_masks.Mask`(*image_w*, *image_h*, *label*,
mask_group)

Bases: `object`

Abstract class to represent image segmentation mask. Its two main subclasses are `RoiMask` and `NeuropilMask`. The former represents the mask of a region of interest (ROI), such as a cell observed in 2-photon imaging. The latter represents the neuropil around that cell, and is useful when subtracting the neuropil signal from the measured ROI signal.

This class should not be instantiated directly.

Parameters

image_w: integer Width of image that ROI resides in

image_h: integer Height of image that ROI resides in

label: text User-defined text label to identify mask

mask_group: integer User-defined number to help put masks into different categories

get_mask_plane(*self*)

Returns mask content on full-size image plane

Returns

numpy 2D array [*img_rows*][*img_cols*]

init_by_pixels(*self*, *border*, *pix_list*)

Initialize mask using a list of mask pixels

Parameters

border: float[4] Coordinates defining useable area of image. See `create_roi_mask()`

pix_list: integer[][2] List of pixel coordinates (x,y) that define the mask

overlaps_motion_border

class `allensdk.brain_observatory.roi_masks.NeuropilMask` (*w, h, label, mask_group*)

Bases: `allensdk.brain_observatory.roi_masks.Mask`

init_by_mask (*self, border, array*)

Initialize mask using spatial mask

Parameters

border: float[4] Border widths on the [right, left, down, up] sides. The resulting neuropil mask will not include pixels falling into a border.

array: integer[image height][image width] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

class `allensdk.brain_observatory.roi_masks.RoiMask` (*image_w, image_h, label, mask_group*)

Bases: `allensdk.brain_observatory.roi_masks.Mask`

init_by_mask (*self, border, array*)

Initialize mask using spatial mask

Parameters

border: float[4] Coordinates defining useable area of image. See `create_roi_mask()`.

roi_mask: integer[image height][image width] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

`allensdk.brain_observatory.roi_masks.calculate_roi_and_neuropil_traces` (*movie_h5, roi_mask_list, motion_border*)

get roi and neuropil masks

`allensdk.brain_observatory.roi_masks.calculate_traces` (*stack, mask_list, block_size=1000*)

Calculates the average response of the specified masks in the image stack

Parameters

stack: float[image height][image width] Image stack that masks are applied to

mask_list: list<Mask> List of masks

Returns

float[number masks][number frames] This is the average response for each Mask in each image frame

`allensdk.brain_observatory.roi_masks.create_neuropil_mask` (*roi, border, combined_binary_mask, label=None*)

Convenience function to create and initialize a Neuropil mask. Neuropil masks are defined as the region around an ROI, up to 13 pixels out, that does not include other ROIs

Parameters

roi: RoiMask object The ROI that the neuropil masks will be based on

border: float[4] Border widths on the [right, left, down, up] sides. The resulting neuropil mask will not include pixels falling into a border.

combined_binary_mask List of pixel coordinates (x,y) that define the mask

combined_binary_mask: integer[image_h][image_w] Image-sized array that shows the position of all ROIs in the image. ROI masks should have a value of one. Background pixels must be zero. In other words, the combined_binary_mask is a bitmap union of all ROI masks

label: text User-defined text label to identify the mask

Returns

NeuropilMask object

```
allensdk.brain_observatory.roi_masks.create_roi_mask(image_w, image_h, border,
                                                    pix_list=None,
                                                    roi_mask=None, label=None,
                                                    mask_group=-1)
```

Convenience function to create and initialize an RoiMask

Parameters

image_w: integer Width of image that ROI resides in

image_h: integer Height of image that ROI resides in

border: float[4] Coordinates defining useable area of image. If the entire image is usable, and masks are valid anywhere in the image, this should be [0, 0, 0, 0]. The following constants help describe the array order:

RIGHT_SHIFT = 0

LEFT_SHIFT = 1

DOWN_SHIFT = 2

UP_SHIFT = 3

When parts of the image are unusable, for example due to motion correction shifting of different image frames, the border array should store the usable image area

pix_list: integer[][2] List of pixel coordinates (x,y) that define the mask

roi_mask: integer[image_h][image_w] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

label: text User-defined text label to identify mask

mask_group: integer User-defined number to help put masks into different categories

Returns

RoiMask object

```
allensdk.brain_observatory.roi_masks.create_roi_mask_array(rois)
```

Create full image mask array from list of RoiMasks.

Parameters

rois: list<RoiMask> List of roi masks.

Returns

np.ndarray: NxWxH array Boolean array of len(rois) image masks.

`allensdk.brain_observatory.roi_masks.validate_mask(mask)`
Check a given roi or neuropil mask for (a subset of) disqualifying problems.

`allensdk.brain_observatory.running_speed` module

class `allensdk.brain_observatory.running_speed.RunningSpeed`
Bases: `tuple`
Describes the rate at which an experimental subject ran during a session.
values `[np.ndarray]` running speed (cm/s) at each sample point
timestamps `[np.ndarray]` The time at which each sample was collected (s).
timestamps
Alias for field number 0
values
Alias for field number 1

`allensdk.brain_observatory.session_analysis` module

class `allensdk.brain_observatory.session_analysis.SessionAnalysis` (*nwb_path*,
save_path)
Bases: `object`

Run all of the stimulus-specific analyses associated with a single experiment session.

Parameters

nwb_path: string, path to NWB file

save_path: string, path to HDF5 file to store outputs. Recommended NOT to modify the NWB file.

append_experiment_metrics (*self*, *metrics*)
Extract stimulus-agnostic metrics from an experiment into a dictionary

append_metadata (*self*, *df*)
Append the metadata fields from the NWB file as columns to a `pd.DataFrame`

append_metrics_drifting_grating (*self*, *metrics*, *dg*)
Extract metrics from the DriftingGratings peak response table into a dictionary.

append_metrics_locally_sparse_noise (*self*, *metrics*, *lsn*)
Extract metrics from the LocallySparseNoise peak response table into a dictionary.

append_metrics_natural_movie_one (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_ONE) peak response table into a dictionary.

append_metrics_natural_movie_three (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_THREE) peak response table into a dictionary.

append_metrics_natural_movie_two (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_TWO) peak response table into a dictionary.

append_metrics_natural_scene (*self*, *metrics*, *ns*)
Extract metrics from the NaturalScenes peak response table into a dictionary.

append_metrics_static_grating (*self, metrics, sg*)

Extract metrics from the StaticGratings peak response table into a dictionary.

save_session_a (*self, dg, nm1, nm3, peak*)

Save the output of session A analysis to self.save_path.

Parameters

dg: DriftingGratings instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm3: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_THREE

peak: pd.DataFrame The combined peak response property table created in self.session_a().

save_session_b (*self, sg, nm1, ns, peak*)

Save the output of session B analysis to self.save_path.

Parameters

sg: StaticGratings instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

ns: NaturalScenes instance

peak: pd.DataFrame The combined peak response property table created in self.session_b().

save_session_c (*self, lsn, nm1, nm2, peak*)

Save the output of session C analysis to self.save_path.

Parameters

lsn: LocallySparseNoise instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm2: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_TWO

peak: pd.DataFrame The combined peak response property table created in self.session_c().

save_session_c2 (*self, lsn4, lsn8, nm1, nm2, peak*)

Save the output of session C2 analysis to self.save_path.

Parameters

lsn4: LocallySparseNoise instance This LocallySparseNoise instance should have been created with self.stimulus = stimulus_info.LOCALLY_SPARSE_NOISE_4DEG.

lsn8: LocallySparseNoise instance This LocallySparseNoise instance should have been created with self.stimulus = stimulus_info.LOCALLY_SPARSE_NOISE_8DEG.

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm2: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_TWO

peak: `pd.DataFrame` The combined peak response property table created in `self.session_c2()`.

session_a (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural movie one, natural movie three, and drifting gratings. The input NWB be for a `stimulus_info.THREE_SESSION_A` experiment.

Parameters

plot_flag: `bool` Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: `bool` Whether to save the output of analysis to `self.save_path` upon completion.

session_b (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural scenes, static gratings, and natural movie one. The input NWB be for a `stimulus_info.THREE_SESSION_B` experiment.

Parameters

plot_flag: `bool` Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: `bool` Whether to save the output of analysis to `self.save_path` upon completion.

session_c (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural movie one, natural movie two, and locally sparse noise. The input NWB be for a `stimulus_info.THREE_SESSION_C` experiment.

Parameters

plot_flag: `bool` Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: `bool` Whether to save the output of analysis to `self.save_path` upon completion.

session_c2 (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for locally sparse noise (4 deg.), locally sparse noise (8 deg.), natural movie one, and natural movie two. The input NWB be for a `stimulus_info.THREE_SESSION_C2` experiment.

Parameters

plot_flag: `bool` Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: `bool` Whether to save the output of analysis to `self.save_path` upon completion.

verify_roi_lists_equal (*self*, *roi1*, *roi2*)

TODO: replace this with simpler numpy comparisons

`allensdk.brain_observatory.session_analysis.main()`

`allensdk.brain_observatory.session_analysis.multi_dataframe_merge(dfs)`

merge a number of `pd.DataFrames` into a single dataframe on their index columns. If any columns are duplicated, prefer the first occuring instance of the column

```
allensdk.brain_observatory.session_analysis.run_session_analysis(nwb_path,
                                                                save_path,
                                                                plot_flag=False,
                                                                save_flag=True)
```

Inspect an NWB file to determine which experiment session was run and compute all stimulus-specific analyses.

Parameters

nwb_path: **string** Path to NWB file.

save_path: **string** path to save results. Recommended NOT to use NWB file.

plot_flag: **bool** Whether to save brain_observatory_plotting work plots.

save_flag: **bool** Whether to save results to save_path.

allensdk.brain_observatory.session_api_utils module

```
class allensdk.brain_observatory.session_api_utils.ParamsMixin(ignore: set =
                                                                {'api'})
```

Bases: object

This mixin adds parameter management functionality to the class it is mixed into.

This mixin expects that the class it is mixed into will have an `__init__` with type annotated parameters. It also expects for the class to have semi-private attributes of the `__init__` type annotated parameters.

Example:

SomeClassWhereParamManagementIsDesired(ParamsMixin):

```
# Managed params should be typed (with simple types if possible)!
def __init__(self,
             param_to_ignore, a_param_1: int, a_param_2: float,
             b_param_1: list):
    # Parameters can be ignored by the mixin super().__init__(ignore={'param_to_ignore'})
    # Pay attention to the naming scheme! self._a_param_1 = a_param_1 self._a_param_2
    = a_param_2 self._b_param_1 = b_param_1
    ...
```

After being mixed in, methods like `'get_params'`, `'set_params'`, `'needs_data_refresh'`, and `'clear_updated_params'` will be available.

clear_updated_params (*self*, *data_params: set*)

This method clears 'updated params' whose data have been updated

get_params (*self*) → Dict[str, Any]

Get managed params and their values

needs_data_refresh (*self*, *data_params: set*) → bool

Check if specific params have been updated via `set_params()`

set_params (*self*, ***params*)

Set managed params

```
allensdk.brain_observatory.session_api_utils.compare_session_fields(x1: Any,
                                                                    x2: Any,
                                                                    err_msg="")
```

Helper function to compare if two fields (attributes) from a Session object are equal to one another.

Parameters

x1 [Any] The field from the first session to compare

x2 [Any] The corresponding field from the second session to compare

err_msg [str, optional] The error message to display if two compared fields do not equal one another, by default "" (an empty string)

`allensdk.brain_observatory.session_api_utils.is_equal(a: Any, b: Any) → bool`

Function to deal with checking if two variables of possibly mixed types have the same value.

`allensdk.brain_observatory.session_api_utils.sessions_are_equal(A, B, reraise=False) → bool`

Check if two Session objects are equal (have same methods and attributes).

Parameters

A [Session A] The first session to compare

B [Session B] The second session to compare

reraise [bool, optional] Whether to reraise when encountering an Assertion or AttributeError, by default False

Returns

bool Whether the two sessions are equal to one another.

`allensdk.brain_observatory.static_gratings` module

class `allensdk.brain_observatory.static_gratings.StaticGratings` (*data_set, **kwargs*)

Bases: `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`

Perform tuning analysis specific to static gratings stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

extralength

static_from_analysis_file (*data_set, analysis_file*)

get_noise_correlation (*self, corr='spearman'*)

get_peak (*self*)

Computes metrics related to each cell's peak response condition.

Returns

Panda data frame with the following fields (_sg suffix is for static grating):

- `ori_sg` (orientation)
- `sf_sg` (spatial frequency)
- `phase_sg`
- `response_variability_sg`
- `osi_sg` (orientation selectivity index)
- `peak_dff_sg` (peak dF/F)

- ptest_sg
- time_to_peak_sg

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# orientations, # spatial frequencies, # phasees, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant response ($p < 0.05$) to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

interlength

number_ori

number_phase

number_sf

open_fan_plot (*self*, *cell_specimen_id*=None, *include_labels*=False, *cell_index*=None)

orivals

phasevals

plot_orientation_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_orientation (*self*, *include_labels*=False, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_spatial_frequency (*self*, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_time_to_peak (*self*, *p_value_max*=0.05, *color_map*=<matplotlib.colors.LinearSegmentedColormap
object at 0x7f5ee73e2410>)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim conditions x repetition for noise correlations

this is a re-organization of the mean sweep response table

sfvals

sweeplength

allensdk.brain_observatory.stimulus_analysis module

class allensdk.brain_observatory.stimulus_analysis.**StimulusAnalysis** (*data_set*)

Bases: object

Base class for all response analysis code. Subclasses are responsible for computing metrics and traces relevant to a particular stimulus. The base class contains methods for organizing sweep responses row of a stimulus stable (get_sweep_response). Subclasses implement the get_response method, computes the mean sweep response to all sweeps for a each stimulus condition.

Parameters

data_set: BrainObservatoryNwbDataSet instance

speed_tuning: boolean, deprecated Whether or not to compute speed tuning histograms

acquisition_rate

binned_cells_sp

binned_cells_vis

binned_dx_sp

binned_dx_vis

cell_id

celltraces

dfftraces

dxcm

dxtime

get_fluorescence (*self*)

get_peak (*self*)

Implemented by subclasses.

get_response (*self*)

Implemented by subclasses.

get_speed_tuning (*self*, *binsize*)

Calculates speed tuning, spontaneous versus visually driven. The return is a 5-tuple of speed and dF/F histograms.

binned_dx_sp: (bins,2) np.ndarray of running speeds binned during spontaneous activity stimulus. The first bin contains all speeds below 1 cm/s. Dimension 0 is mean running speed in the bin. Dimension 1 is the standard error of the mean.

binned_cells_sp: (bins,2) np.ndarray of fluorescence during spontaneous activity stimulus. First bin contains all data for speeds below 1 cm/s. Dimension 0 is mean fluorescence in the bin. Dimension 1 is the standard error of the mean.

binned_dx_vis: (bins,2) np.ndarray of running speeds outside of spontaneous activity stimulus. The first bin contains all speeds below 1 cm/s. Dimension 0 is mean running speed in the bin. Dimension 1 is the standard error of the mean.

binned_cells_vis: np.ndarray of fluorescence outside of spontaneous activity stimulus. First bin contains all data for speeds below 1 cm/s. Dimension 0 is mean fluorescence in the bin. Dimension 1 is the standard error of the mean.

peak_run: pd.DataFrame of speed-related properties of a cell.

Returns

tuple: **binned_dx_sp**, **binned_cells_sp**, **binned_dx_vis**, **binned_cells_vis**, **peak_run**

get_sweep_response (*self*)

Calculates the response to each sweep in the stimulus table for each cell and the mean response. The return is a 3-tuple of:

- **sweep_response:** pd.DataFrame of response dF/F traces organized by cell (column) and sweep (row)

- `mean_sweep_response`: mean values of the traces returned in `sweep_response`
- `pval`: p value from 1-way ANOVA comparing response during sweep to response prior to sweep

Returns

3-tuple: `sweep_response`, `mean_sweep_response`, `pval`

`mean_sweep_response`

`numbercells`

`peak`

`peak_run`

`plot_representational_similarity` (*self*, *repsim*, *stimulus=False*)

`plot_running_speed_histogram` (*self*, *xlim=None*, *nbins=None*)

`plot_speed_tuning` (*self*, *cell_specimen_id=None*, *cell_index=None*, *evoked_color='b30000'*, *spontaneous_color='#0000b3'*)

`populate_stimulus_table` (*self*)

Implemented by subclasses.

`pval`

`response`

`roi_id`

`row_from_cell_id` (*self*, *csid=None*, *idx=None*)

`stim_table`

`sweep_response`

`timestamps`

`allensdk.brain_observatory.stimulus_analysis.nonraising_ks_2samp` (*data1*, *data2*, ***kwargs*)

`scipy.stats.ks_2samp` now raises a `ValueError` if one of the input arrays is of length 0. Previously it signaled this case by returning nans. This function restores the prior behavior.

`allensdk.brain_observatory.stimulus_info` module

class `allensdk.brain_observatory.stimulus_info.BinaryIntervalSearchTree` (*search_list*)
Bases: `object`

add (*self*, *input_list*, *tmap=None*)

static from_df (*input_df*)

search (*self*, *fi*, *tmap=None*)

class `allensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor` (*experiment_geometry=None*)
Bases: `allensdk.brain_observatory.stimulus_info.Monitor`

http://help.brain-map.org/display/observatory/Documentation?preview=/10616846/10813485/VisualCoding_VisualStimuli.pdf <https://www.cnet.com/products/asus-pa248q/specs/>

grating_to_screen (*self*, *phase*, *spatial_frequency*, *orientation*, ***kwargs*)

lsn_image_to_screen (*self*, *img*, ***kwargs*)

```
pixels_to_visual_degrees (self, n, **kwargs)
visual_degrees_to_pixels (self, vd, **kwargs)
warp_image (self, img, **kwargs)
class allensdk.brain_observatory.stimulus_info.ExperimentGeometry (distance,
                                                                    mon_height_cm,
                                                                    mon_width_cm,
                                                                    mon_res,
                                                                    eyepoint)
    Bases: object
    generate_warp_coordinates (self)
    warp_coordinates
class allensdk.brain_observatory.stimulus_info.Monitor (n_pixels_r,      n_pixels_c,
                                                         panel_size, spatial_unit)
    Bases: object
    aspect_ratio
    get_mask (self)
    grating_to_screen (self, phase, spatial_frequency, orientation, distance_from_monitor,
                       p2p_amp=256, baseline=127, translation=(0, 0))
    height
    lsn_image_to_screen (self, img, stimulus_type, origin='lower', background_color=127, translation=(0, 0))
    map_stimulus (self, source_stimulus_coordinate, source_stimulus_type, target_stimulus_type)
    mask
    natural_movie_image_to_screen (self, img, origin='lower', translation=(0, 0))
    natural_scene_image_to_screen (self, img, origin='lower', translation=(0, 0))
    panel_size
    pixel_size
    pixels_to_visual_degrees (self, n, distance_from_monitor, small_angle_approximation=True)
    set_spatial_unit (self, new_unit)
    show_image (self, img, ax=None, show=True, mask=False, warp=False, origin='lower')
    spatial_frequency_to_pix_per_cycle (self, spatial_frequency, distance_from_monitor)
    visual_degrees_to_pixels (self,          vd,          distance_from_monitor,
                              small_angle_approximation=True)
    width
class allensdk.brain_observatory.stimulus_info.StimulusSearch (nwb_dataset)
    Bases: object
    search (self, fi)
allensdk.brain_observatory.stimulus_info.all_stimuli ()
    Return a list of all stimuli in the data set
```

```

allensdk.brain_observatory.stimulus_info.get_spatial_grating(height=None, aspect_ratio=None,
                                                             ori=None,
                                                             pix_per_cycle=None,
                                                             phase=None,
                                                             p2p_amp=2,
                                                             baseline=0)

allensdk.brain_observatory.stimulus_info.get_spatio_temporal_grating(t,
                                                                      tempo-
                                                                      ral_frequency=None,
                                                                      **kwargs)

allensdk.brain_observatory.stimulus_info.lsn_coordinate_to_monitor_coordinate(lsn_coordinate,
                                                                              mon-
                                                                              i-
                                                                              tor_shape,
                                                                              stim-
                                                                              u-
                                                                              lus_type)

allensdk.brain_observatory.stimulus_info.make_display_mask(display_shape=(1920,
                                                            1200))
    Build a display-shaped mask that indicates which pixels are on screen after warping the stimulus.

allensdk.brain_observatory.stimulus_info.map_monitor_coordinate_to_stimulus_coordinate(moni-
                                                                                      mon-
                                                                                      i-
                                                                                      tor_s
                                                                                      stim-
                                                                                      u-
                                                                                      lus_t)

allensdk.brain_observatory.stimulus_info.map_monitor_coordinate_to_template_coordinate(moni-
                                                                                      mon-
                                                                                      i-
                                                                                      tor_s
                                                                                      tem-
                                                                                      plate)

allensdk.brain_observatory.stimulus_info.map_stimulus(source_stimulus_coordinate,
                                                       source_stimulus_type, target_stimulus_type,
                                                       monitor_shape)

allensdk.brain_observatory.stimulus_info.map_stimulus_coordinate_to_monitor_coordinate(tempo-
                                                                                      mon-
                                                                                      i-
                                                                                      tor_s
                                                                                      stim-
                                                                                      u-
                                                                                      lus_t)

allensdk.brain_observatory.stimulus_info.map_template_coordinate_to_monitor_coordinate(tempo-
                                                                                      mon-
                                                                                      i-
                                                                                      tor_s
                                                                                      tem-
                                                                                      plate)

```

```
allensdk.brain_observatory.stimulus_info.mask_stimulus_template(template_display_coords,  
                                                                tem-  
                                                                plate_shape,  
                                                                dis-  
                                                                play_mask=None,  
                                                                thresh-  
                                                                old=1.0)
```

Build a mask for a stimulus template of a given shape and display coordinates that indicates which part of the template is on screen after warping.

Parameters

template_display_coords: **list** list of (x,y) display coordinates

template_shape: **tuple** (width,height) of the display template

display_mask: **np.ndarray** boolean 2D mask indicating which display coordinates are on screen after warping.

threshold: **float** Fraction of pixels associated with a template display coordinate that should remain on screen to count as belonging to the mask.

Returns

tuple: (template mask, pixel fraction)

```
allensdk.brain_observatory.stimulus_info.monitor_coordinate_to_lsn_coordinate(monitor_coordinate,  
                                                                mon-  
                                                                i-  
                                                                tor_shape,  
                                                                stim-  
                                                                u-  
                                                                lus_type)
```

```
allensdk.brain_observatory.stimulus_info.monitor_coordinate_to_natural_movie_coordinate(monitor_coordinate,  
                                                                mon-  
                                                                i-  
                                                                tor_shape,  
                                                                stim-  
                                                                u-  
                                                                lus_type)
```

```
allensdk.brain_observatory.stimulus_info.natural_movie_coordinate_to_monitor_coordinate(natural_movie_coordinate,  
                                                                mon-  
                                                                i-  
                                                                tor_shape,  
                                                                stim-  
                                                                u-  
                                                                lus_type)
```

```
allensdk.brain_observatory.stimulus_info.natural_scene_coordinate_to_monitor_coordinate(natural_scene_coordinate,  
                                                                mon-  
                                                                i-  
                                                                tor_shape,  
                                                                stim-  
                                                                u-  
                                                                lus_type)
```

```
allensdk.brain_observatory.stimulus_info.rotate(X, Y, theta)
```

```
allensdk.brain_observatory.stimulus_info.sessions_with_stimulus(stimulus)
```

Return the names of the sessions that contain a given stimulus.

```
allensdk.brain_observatory.stimulus_info.stimuli_in_session(session, allow_unknown=True)
```

Return a list what stimuli are available in a given session.

Parameters

session: **string** Must be one of: [stimulus_info.THREE_SESSION_A, stimulus_info.THREE_SESSION_B, stimulus_info.THREE_SESSION_C, stimulus_info.THREE_SESSION_C2]

```
allensdk.brain_observatory.stimulus_info.translate_image_and_fill(img, translation=(0,
                                                                    0))

allensdk.brain_observatory.stimulus_info.warp_stimulus_coords(vertices, distance=15.0,
                                                                mon_height_cm=32.5,
                                                                mon_width_cm=51.0,
                                                                mon_res=(1920,
                                                                1200), eyepoint=(0.5, 0.5))
```

For a list of screen vertices, provides a corresponding list of texture coordinates.

Parameters

vertices: `numpy.ndarray` `[[x0,y0], [x1,y1], ...]` A set of vertices to convert to texture positions.

distance: `float` distance from the monitor in cm.

mon_height_cm: `float` monitor height in cm

mon_width_cm: `float` monitor width in cm

mon_res: `tuple` monitor resolution (x,y)

eyepoint: `tuple`

Returns

np.ndarray x,y coordinates shaped like the input that describe what pixel coordinates are displayed an the input coordinates after warping the stimulus.

allensdk.brain_observatory.sync_dataset module

dataset.py

Dataset object for loading and unpacking an HDF5 dataset generated by sync.py

@author: derricw

Allen Institute for Brain Science

Dependencies

numpy <http://www.numpy.org/> h5py <http://www.h5py.org/>

class allensdk.brain_observatory.sync_dataset.**Dataset** (*path*)
Bases: object

A sync dataset. Contains methods for loading and parsing the binary data.

Parameters

path [str] Path to HDF5 file.

Examples

```
>>> dset = Dataset('my_h5_file.h5')
>>> logger.info(dset.meta_data)
>>> dset.stats()
>>> dset.close()
```

```
>>> with Dataset('my_h5_file.h5') as d:
...     logger.info(dset.meta_data)
...     dset.stats()
```

The sync file documentation from MPE can be found at [sharepoint > Instrumentation > Shared Documents > Sync_line_labels_discussion_2020-01-27-.xlsx](https://alleninstitute.sharepoint.com/:x/s/Instrumentation/ES2bi1xJ3E9NupX-zQeXTIYBS2mVVySycfbCQhsD_jPMUw?e=Z9jCwH) # NOQA E501 Direct link: https://alleninstitute.sharepoint.com/:x/s/Instrumentation/ES2bi1xJ3E9NupX-zQeXTIYBS2mVVySycfbCQhsD_jPMUw?e=Z9jCwH

```
BEHAVIOR_TRACKING_KEYS = ('beh_frame_received', 'cam1_exposure', 'behavior_monitoring')
```

```
DEPRECATED_KEYS = {}
```

```
EYE_TRACKING_KEYS = ('eye_frame_received', 'cam2_exposure', 'eyetracking', 'eye_tracking')
```

```
FRAME_KEYS = ('frames', 'stim_vsync')
```

```
OPTOGENETIC_STIMULATION_KEYS = ('LED_sync', 'opto_trial')
```

```
PHOTODIODE_KEYS = ('photodiode', 'stim_photodiode')
```

```
analog_meta_data
```

```
close(self)
```

Closes the dataset.

```
duty_cycle(self, line)
```

Doesn't work right now. Freezes python for some reason.

Returns the duty cycle of a line.

```
frequency(self, line, edge='rising')
```

Returns the average frequency of a line.

```
get_all_bits(self)
```

Returns the data for all bits.

```
get_all_events(self)
```

Returns all counter values and their cooresponding IO state.

```
get_all_times(self, units='samples')
```

Returns all counter values.

Parameters

units [str] Return times in 'samples' or 'seconds'

```
get_analog_channel(self, channel, start_time=0.0, stop_time=None, downsample=1)
```

Returns the data from the specified analog channel between the timepoints.

Args: channel (int, str): desired channel index or label start_time (Optional[float]): start time in seconds

stop_time (Optional[float]): stop time in seconds downsample (Optional[int]): downsample factor

Returns: ndarray: slice of data for specified channel

Raises: KeyError: no analog data present

get_analog_meta (*self*)

Returns the metadata for the analog data.

get_bit (*self*, *bit*)

Returns the values for a specific bit.

Parameters

bit [int] Bit to return.

get_bit_changes (*self*, *bit*)

Returns the first derivative of a specific bit. Data points are 1 on rising edges and 255 on falling edges.

Parameters

bit [int] Bit for which to return changes.

get_edges (*self*, *kind*: str, *keys*: Union[str, Sequence[str]], *units*: str = 'seconds', *permissive*: bool = False) → Union[numpy.ndarray, NoneType]

Utility function for extracting edge times from a line

Parameters

kind [One of “rising”, “falling”, or “all”. Should this method return] timestamps for rising, falling or both edges on the appropriate line

keys [These will be checked in sequence. Timestamps will be returned] for the first which is present in the line labels

units [one of “seconds”, “samples”, or “indices”. The returned] “time”stamps will be given in these units.

raise_missing [If True and no matching line is found, a KeyError will] be raised

Returns

An array of edge times. If raise_missing is False and none of the keys were found, returns None.

Raises

KeyError [none of the provided keys were found among this dataset’s] line labels

get_events_by_bit (*self*, *bit*, *units*=‘samples’)

Returns all counter values for transitions (both rising and falling) for a specific bit.

Parameters

bit [int] Bit for which to return events.

get_events_by_line (*self*, *line*, *units*=‘samples’)

Returns all counter values for transitions (both rising and falling) for a specific line.

Parameters

line [str] Line for which to return events.

get_falling_edges (*self*, *line*, *units*=‘samples’)

Returns the counter values for the falling edges for a specific bit or line.

Parameters

line [str] Line for which to return edges.

get_line (*self*, *line*)

Returns the values for a specific line.

Parameters

line [str] Line to return.

get_line_changes (*self*, *line*)

Returns the first derivative of a specific line. Data points are 1 on rising edges and 255 on falling edges.

Parameters

line [(str)] Line name for which to return changes.

get_nearest (*self*, *source*, *target*, *source_edge*='rising', *target_edge*='rising', *direction*='previous', *units*='indices')

For all values of the source line, finds the nearest edge from the target line.

By default, returns the indices of the target edges.

Args: *source* (str, int): desired source line *target* (str, int): desired target line *source_edge* [Optional(str)]: “rising” or “falling” source edges *target_edge* [Optional(str)]: “rising” or “falling” target edges *direction* (str): “previous” or “next”. Whether to prefer the previous edge or the following edge.
units (str): “indices”

get_rising_edges (*self*, *line*, *units*='samples')

Returns the counter values for the rizing edges for a specific bit or line.

Parameters

line [str] Line for which to return edges.

line_stats (*self*, *line*, *print_results*=True)

Quick-and-dirty analysis of a bit.

##TODO: Split this up into smaller functions.

load (*self*, *path*)

Loads an hdf5 sync dataset.

Parameters

path [str] Path to hdf5 file.

period (*self*, *line*, *edge*='rising')

Returns a dictionary with avg, min, max, and st of period for a line.

plot_all (*self*, *start_time*, *stop_time*, *auto_show*=True)

Plot all active bits.

Yikes. Come up with a better way to show this.

plot_bit (*self*, *bit*, *start_time*=0.0, *end_time*=None, *auto_show*=True, *axes*=None, *name*="")

Plots a specific bit at a specific time period.

plot_bits (*self*, *bits*, *start_time=0.0*, *end_time=None*, *auto_show=True*)
Plots a list of bits.

plot_line (*self*, *line*, *start_time=0.0*, *end_time=None*, *auto_show=True*)
Plots a specific line at a specific time period.

plot_lines (*self*, *lines*, *start_time=0.0*, *end_time=None*, *auto_show=True*)
Plots specific lines at a specific time period.

sample_freq

stats (*self*)

Quick-and-dirty analysis of all bits. Prints a few things about each bit where events are found.

`allensdk.brain_observatory.sync_dataset.get_bit (uint_array, bit)`
Returns a bool array for a specific bit in a uint ndarray.

Parameters

uint_array [(numpy.ndarray)] The array to extract bits from.

bit [(int)] The bit to extract.

`allensdk.brain_observatory.sync_dataset.unpack_uint32 (uint32_array, endian='L')`
Unpacks an array of 32-bit unsigned integers into bits.

Default is least significant bit first.

***Not currently used by sync dataset because get_bit is better and does** basically the same thing. I'm just leaving it in because it could potentially account for endianness and possibly have other uses in the future.

Module contents

class `allensdk.brain_observatory.JSONEncoder` (*, *skipkeys=False*, *ensure_ascii=True*,
check_circular=True, *allow_nan=True*,
sort_keys=False, *indent=None*, *separators=None*, *default=None*)

Bases: `json.encoder.JSONEncoder`

default (*self*, *o*)

Implement this method in a subclass such that it returns a serializable object for *o*, or calls the base implementation (to raise a `TypeError`).

For example, to support arbitrary iterators, you could implement default like this:

```
def default(self, o):
    try:
        iterable = iter(o)
    except TypeError:
        pass
    else:
        return list(iterable)
    # Let the base class default method raise the TypeError
    return JSONEncoder.default(self, o)
```

`allensdk.brain_observatory.dict_to_indexed_array (dc, order=None)`

Given a dictionary and an ordered arr, build a concatenation of the dictionary's values and an index describing how that concatenation can be unpacked

`allensdk.brain_observatory.hook (json_dict)`

6.1.3 allensdk.config package

Subpackages

allensdk.config.app package

Submodules

allensdk.config.app.application_config module

```
class allensdk.config.app.application_config.ApplicationConfig(defaults,
                                                                name='app',
                                                                help='Run ap-
                                                                plication.', de-
                                                                fault_log_config=None)
```

Bases: object

Convenience class that handles of application configuration from environment variables, .conf files and the command line using Python standard libraries and formats.

apply_configuration_from_command_line (*self*, *parsed_args*)

Read application configuration variables from the command line.

Unassigned variables are left unchanged if previously assigned, set to their default values, or None if no default is specified at init time. Assigned variables will overwrite the previous value.

see: <https://docs.python.org/2/howto/argparse.html>

Parameters

parsed_args [dict] the arguments as parsed from the command line.

apply_configuration_from_environment (*self*)

Read application configuration variables from the environment.

The variable names are upper case and have a prefix defined by the application.

See: <https://docs.python.org/2/library/os.html>

apply_configuration_from_file (*self*, *config_file_path*)

Read application configuration variables from a .conf file.

Unassigned variables are set to their default values or None if no default is specified at init time. The variables are found in a section named by the application.

Parameters

config_file_path [string] path to to an INI (.conf) or JSON format application config file.

Returns

see: <https://docs.python.org/2/library/configparser.html>

create_argparser (*self*)

Initialization for the command-line parsing stage.

An application specific prefix is applied to argument names.

Parameters

prog [string] Application specific prefix for argument names.

description [string] A brief 'help' description of the application.

Returns

argParse.ArgumentParser The initialized argument parser object.

Notes

Defaults are set at the first environment reading. Command line args only override them when present

from_json_file (*self*, *json_path*)

Read an application configuration from a JSON format file.

Parameters

json_path [string] Path to the JSON file.

Returns

string An application configuration in INI format

from_json_string (*self*, *json_string*)

Read a configuration from a JSON format string.

Parameters

json_string [string] A JSON-formatted string containing an application configuration.

Returns

string An application configuration in INI format

load (*self*, *command_line_args*, *disable_existing_loggers=True*)

Load application configuration options, first from the environment, then from the configuration file, then from the command line.

Each stage of loading can override the previous stage.

Parameters

command_line_args [dict] Parameters passed to the application.

disable_existing_loggers [boolean] Reset the logging system or not.

Returns

fileConfig Configuration object with all levels applied

parse_command_line_args (*self*, *args*)

Simply call the internal argparser object.

Parameters

args [array] Parameters passed to the application.

Returns

Namespace Parsed parameters.

to_config_string (*self*, *description*)

Create a configuration string from a dict.

Parameters

description [dict] Configuration options for an application.

Returns

string Equivalent configuration as an INI format string

Notes

The Python configparser library natively supports this functionality in Python 3.

Module contents

`allensdk.config.app` is a package that assists in configuring application software, as opposed to domain-specific configuration.

`allensdk.config.model` package

Subpackages

`allensdk.config.model.formats` package

Submodules

`allensdk.config.model.formats.hdf5_util` module

```
class allensdk.config.model.formats.hdf5_util.Hdf5Util
    Bases: object
    read (self, file_path)
    write (self, file_path, m)
```

`allensdk.config.model.formats.json_description_parser` module

```
class allensdk.config.model.formats.json_description_parser.JsonDescriptionParser
    Bases: allensdk.config.model.description_parser.DescriptionParser
    log = <Logger allensdk.config.model.formats.json_description_parser (WARNING)>
    read (self, file_path, description=None, section=None, **kwargs)
        Parse a complete or partial configuration.
```

Parameters

- json_string** [string] Input to parse.
- description** [Description, optional] Where to put the parsed configuration. If None a new one is created.
- section** [string, optional] Where to put the parsed configuration within the description.

Returns

- Description** The input description with parsed configuration added.
- Section is only specified for “bare” objects that are to be added to a section array.**

```
read_string (self, json_string, description=None, section=None, **kwargs)
    Parse a complete or partial configuration.
```

Parameters

- json_string** [string] Input to parse.

description [Description, optional] Where to put the parsed configuration. If None a new one is created.

section [string, optional] Where to put the parsed configuration within the description.

Returns

Description The input description with parsed configuration added.

Section is only specified for “bare” objects that are to be added to a section array.

write (*self*, *filename*, *description*)

Write the description to a JSON file.

Parameters

description [Description] Object to write.

write_string (*self*, *description*)

Write the description to a JSON string.

Parameters

description [Description] Object to write.

Returns

string JSON serialization of the input.

allensdk.config.model.formats.pycfg_description_parser module

class allensdk.config.model.formats.pycfg_description_parser.**PycfgDescriptionParser**
Bases: *allensdk.config.model.description_parser.DescriptionParser*

log = <Logger allensdk.config.model.formats.pycfg_description_parser (WARNING)>

read (*self*, *pycfg_file_path*, *description=None*, *section=None*, ***kwargs*)

Read a serialized description from a Python (.pycfg) file.

Parameters

filename [string] Name of the .pycfg file.

Returns

Description Configuration object.

read_string (*self*, *python_string*, *description=None*, *section=None*, ***kwargs*)

Read a serialized description from a Python (.pycfg) string.

Parameters

python_string [string] Python string with a serialized description.

Returns

Description Configuration object.

write (*self*, *filename*, *description*)

Write the description to a Python (.pycfg) file.

Parameters

filename [string] Name of the file to write.

write_string (*self*, *description*)

Write the description to a pretty-printed Python string.

Parameters

description [Description] Configuration object to write.

Module contents

Submodules

allensdk.config.model.description module

class allensdk.config.model.description.**Description**

Bases: object

fix_unary_sections (*self*, *section_names=None*)

Wrap section contents that don't have the proper array surrounding them in an array.

Parameters

section_names [list of strings, optional] Keys of sections that might not be in array form.

is_empty (*self*)

Check if anything is in the object.

Returns

boolean true if self.data is missing or empty

unpack (*self*, *data*, *section=None*)

Read the manifest and other stand-alone configuration structure, or insert a configuration object into a section of an existing configuration.

Parameters

data [dict] A configuration object including top level sections, or an configuration object to be placed within a section.

section [string, optional.] If this is present, place data within an existing section array.

unpack_manifest (*self*, *data*)

Pull the manifest configuration section into a separate place.

Parameters

data [dict] A configuration structure that still has a manifest section.

update_data (*self*, *data*, *section=None*)

Merge configuration data possibly from multiple files.

Parameters

data [dict] Configuration structure to add.

section [string, optional] What configuration section to read it into if the file does not specify.

allensdk.config.model.description_parser module

```
class allensdk.config.model.description_parser.DescriptionParser
```

```
Bases: object
```

```
log = <Logger allensdk.config.model.description_parser (WARNING)>
```

```
parser_for_extension (self, filename)
```

```
    Choose a subclass that can read the format.
```

Parameters

```
    filename [string] For the extension.
```

Returns

```
    DescriptionParser Appropriate subclass.
```

```
read (self, file_path, description=None, section=None, **kwargs)
```

```
    Parse data needed for a simulation.
```

Parameters

```
    description [dict] Configuration from parsing previous files.
```

```
    section [string, optional] What configuration section to read it into if the file does not
    specify.
```

```
read_string (self, data_string, description=None, section=None, header=None)
```

```
    Parse data needed for a simulation from a string.
```

```
write (self, filename, description)
```

```
    Save the configuration.
```

Parameters

```
    filename [string] Name of the file to write.
```

Module contents**Submodules****allensdk.config.manifest module**

```
class allensdk.config.manifest.Manifest (config=None,      relative_base_dir='.',      ver-
```

```
sion=None)
```

```
Bases: object
```

```
    Manages the location of external files referenced in an Allen SDK configuration
```

```
DIR = 'dir'
```

```
DIRNAME = 'dir_name'
```

```
FILE = 'file'
```

```
VERSION = 'manifest_version'
```

```
add_file (self, file_key, file_name, dir_key=None, path_format=None)
```

```
    Insert a new file entry.
```

Parameters

file_key [string] Reference to the entry.

file_name [string] Substitutions of the %s, %d style allowed.

dir_key [string] Reference to the parent directory entry.

path_format [string, optional] File type for further parsing.

add_path (*self*, *key*, *path*, *path_type*=*'dir'*, *absolute*=*True*, *path_format*=*None*, *parent_key*=*None*)
Insert a new entry.

Parameters

key [string] Identifier for referencing the entry.

path [string] Specification for a path using %s, %d style substitution.

path_type [string enumeration] 'dir' (default) or 'file'

absolute [boolean] Is the spec relative to the process current directory.

path_format [string, optional] Indicate a known file type for further parsing.

parent_key [string] Refer to another entry.

add_paths (*self*, *path_info*)
add information about paths stored in the manifest.

Parameters

path_info [dict] Information about the new paths

as_dataframe (*self*)

check_dir (*self*, *path_key*, *do_exit*=*False*)
Verify a directories existence or optionally exit.

Parameters

path_key [string] Reference to the entry.

do_exit [boolean] What to do if the directory is not present.

create_dir (*self*, *path_key*)
Make a directory for an entry.

Parameters

path_key [string] Reference to the entry.

get_format (*self*, *path_key*)
Retrieve the type of a path entry.

Parameters

path_key [string] reference to the entry

Returns

string File type.

get_path (*self*, *path_key*, **args*)
Retrieve an entry with substitutions.

Parameters

path_key [string] Refer to the entry to retrieve.

args [any types, optional] arguments to be substituted into the path spec for %s, %d, etc.

Returns

string Path with parent structure and substitutions applied.

load_config (*self*, *config*, *version=None*)

Load paths into the manifest from an Allen SDK config section.

Parameters

config [Config] Manifest section of an Allen SDK config.

log = <Logger allensdk.config.manifest (WARNING)>

resolve_paths (*self*, *description_dict*, *suffix='_key'*)

Walk input items and expand those that refer to a manifest entry.

Parameters

description_dict [dict] Any entries with key names ending in suffix will be expanded.

suffix [string] Indicates the entries to be expanded.

classmethod safe_make_parent_dirs (*file_name*)

Create a parent directories for file.

Parameters

file_name [string]

Returns

leftmost [string] most rootward directory created

classmethod safe_mkdir (*directory*)

Create path if not already there.

Parameters

directory [string] create it if it doesn't exist

Returns

leftmost [string] most rootward directory created

exception allensdk.config.manifest.**ManifestVersionError** (*message*, *version*,
found_version)

Bases: Exception

outdated

allensdk.config.manifest_builder module

class allensdk.config.manifest_builder.**ManifestBuilder**

Bases: object

add_path (*self*, *key*, *spec*, *typename='dir'*, *parent_key=None*, *format=None*)

add_section (*self*, *name*, *contents*)

as_dataframe (*self*)

df_columns = ['key', 'parent_key', 'spec', 'type', 'format']

from_dataframe (*self*, *df*)

get_config (*self*)

```
get_manifest (self)
set_version (self, value)
write_json_file (self, path, overwrite=False)
write_json_string (self)
```

Module contents

`allensdk.config.enable_console_log (level=None)`
configure allensdk logging to output to the console.

Parameters

level [int] logging level 0-50 (logging.INFO, logging.DEBUG, etc.)

Notes

See: [Logging Cookbook](#)

6.1.4 allensdk.core package

Subpackages

`allensdk.core.lazy_property` package

Submodules

`allensdk.core.lazy_property.lazy_property` module

```
class allensdk.core.lazy_property.lazy_property.LazyProperty (api_method:
                                                                Callable, wrap-
                                                                pers:      Iter-
                                                                able[T_co] =
                                                                (), settable: bool
                                                                = False, *args,
                                                                **kwargs)
```

Bases: object

calculate (self)

`allensdk.core.lazy_property.lazy_property_mixin` module

```
class allensdk.core.lazy_property.lazy_property_mixin.LazyPropertyMixin
```

Bases: object

LazyProperty

Module contents

Submodules

allensdk.core.auth_config module

allensdk.core.authentication module

class allensdk.core.authentication.CredentialProvider

Bases: abc.ABC

METHOD = 'custom'

provide (*self*, *credential*)

class allensdk.core.authentication.DbCredentials (*dbname*, *user*, *host*, *port*, *password*)

Bases: tuple

dbname

Alias for field number 0

host

Alias for field number 2

password

Alias for field number 4

port

Alias for field number 3

user

Alias for field number 1

class allensdk.core.authentication.EnvCredentialProvider (*environ*: *Optional*[*Dict*[*str*, *Any*]] = *None*)

Bases: *allensdk.core.authentication.CredentialProvider*

Provides credentials from environment variables for variables listed in CREDENTIAL_KEYS.

METHOD = 'env'

provide (*self*, *credential*)

allensdk.core.authentication.credential_injector (*credential_map*: *Dict*[*str*, *Any*], *provider*: *Union*[*allensdk.core.authentication.CredentialProvider*, *NoneType*] = *None*)

Decorator used to inject credentials from another source if not explicitly provided in the function call. This function will only supply values for keyword arguments. All keys defined in *credential_map* must correspond to keyword arguments in the function signature.

Parameters

credential_map: *Dict*[*Str*: *Any*] Dictionary where the keys are the keyword of a credential kwarg passed to the decorated function, and the values are the name of the credential in the credential provider (see CREDENTIAL_KEYS).

Example of *credential_map* for PostgresQueryMixin connecting to LIMS database:

```
{ "dbname": "LIMS_DBNAME", "user": "LIMS_USER", "host":
  "LIMS_HOST", "password": "LIMS_PASSWORD", "port": "LIMS_PORT" }
```

```
}
```

provider: Optional[CredentialProvider] Subclass of CredentialProvider to provide credentials to the wrapped function. If left unspecified, will default to EnvCredentialProvider, which provides credentials from environment variables.

```
allensdk.core.authentication.get_credential_provider()
```

```
allensdk.core.authentication.set_credential_provider(provider)
```

allensdk.core.brain_observatory_cache module

allensdk.core.brain_observatory_nwb_data_set module

class allensdk.core.brain_observatory_nwb_data_set.**BrainObservatoryNwbDataSet** (*nwb_file*)
Bases: object

FILE_METADATA_MAPPING = {'age': 'general/subject/age', 'device_string': 'general/dev

MOTION_CORRECTION_DATASETS = ['MotionCorrection/2p_image_series/xy_translations', 'Mot

PIPELINE_DATASET = 'brain_observatory_pipeline'

STIMULUS_TABLE_TYPES = {'abstract_feature_series': ['drifting_gratings', 'static_grat

SUPPORTED_PIPELINE_VERSION = '3.0'

get_cell_specimen_ids (*self*)

Returns an array of cell IDs for all cells in the file

Returns

cell specimen IDs: list

get_cell_specimen_indices (*self, cell_specimen_ids*)

Given a list of cell specimen ids, return their index based on their order in this file.

Parameters

cell_specimen_ids: list of cell specimen ids

get_corrected_fluorescence_traces (*self, cell_specimen_ids=None*)

Returns an array of demixed and neuropil-corrected fluorescence traces for all ROIs and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Corrected fluorescence traces for each cell

get_demixed_traces (*self, cell_specimen_ids=None*)

Returns an array of demixed fluorescence traces for all ROIs and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Demixed fluorescence traces for each cell

get_dff_traces (*self*, *cell_specimen_ids=None*)

Returns an array of dF/F traces for all ROIs and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return data for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

dF/F: 2D numpy array dF/F values for each cell

get_fluorescence_timestamps (*self*)

Returns an array of timestamps in seconds for the fluorescence traces

get_fluorescence_traces (*self*, *cell_specimen_ids=None*)

Returns an array of fluorescence traces for all ROI and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Fluorescence traces for each cell

get_locally_sparse_noise_stimulus_template (*self*, *stimulus*, *mask_off_screen=True*)

Return an array of the stimulus template for the specified stimulus.

Parameters

stimulus: string

Which locally sparse noise stimulus to retrieve. Must be one of:

stimulus_info.LOCALLY_SPARSE_NOISE	stimu-
lus_info.LOCALLY_SPARSE_NOISE_4DEG	stimu-
lus_info.LOCALLY_SPARSE_NOISE_8DEG	

mask_off_screen: boolean Set off-screen regions of the stimulus to LocallySparseNoise.LSN_OFF_SCREEN.

Returns

tuple: (template, off-screen mask)

get_max_projection (*self*)

Returns the maximum projection image for the 2P movie.

Returns

max projection: np.ndarray

get_metadata (*self*)

Returns a dictionary of meta data associated with each experiment, including Cre line, specimen number, visual area imaged, imaging depth

Returns

metadata: dictionary

get_motion_correction (*self*)

Returns a Panda DataFrame containing the x- and y- translation of each image used for image alignment

get_neuropil_r (*self*, *cell_specimen_ids=None*)

Returns a scalar value of r for neuropil correction of fluorescence traces

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then results for all are returned

Returns

r: 1D numpy array, len(r)=len(cell_specimen_ids) Scalar for neuropil subtraction for each cell

get_neuropil_traces (*self*, *cell_specimen_ids=None*)

Returns an array of neuropil fluorescence traces for all ROIs and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Neuropil fluorescence traces for each cell

get_pupil_location (*self*, *as_spherical=True*)

Returns the x, y pupil location.

Parameters

as_spherical [bool] Whether to return the location as spherical (default) or not. If true, the result is altitude and azimuth in degrees, otherwise it is x, y in centimeters. (0,0) is the center of the monitor.

Returns

(timestamps, location) Timestamps is an (Nx1) array of timestamps in seconds. Location is an (Nx2) array of spatial location.

get_pupil_size (*self*)

Returns the pupil area in pixels.

Returns

(timestamps, areas) Timestamps is an (Nx1) array of timestamps in seconds. Areas is an (Nx1) array of pupil areas in pixels.

get_roi_ids (*self*)

Returns an array of IDs for all ROIs in the file

Returns

ROI IDs: list

get_roi_mask (*self*, *cell_specimen_ids=None*)

Returns an array of all the ROI masks

Parameters

cell specimen IDs: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

List of ROI_Mask objects

get_roi_mask_array (*self*, *cell_specimen_ids=None*)

Return a numpy array containing all of the ROI masks for requested cells. If *cell_specimen_ids* is omitted, return all masks.

Parameters

cell_specimen_ids: list List of cell specimen ids. Default None.

Returns

np.ndarray: NxWxH array, where N is number of cells

get_running_speed (*self*)

Returns the mouse running speed in cm/s

get_session_type (*self*)

Returns the type of experimental session, presently one of the following: *three_session_A*, *three_session_B*, *three_session_C*

Returns

session type: string

get_stimulus (*self*, *frame_ind*)

get_stimulus_epoch_table (*self*)

Returns a pandas dataframe that summarizes the stimulus epoch duration for each acquisition time index in the experiment

Parameters

None

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Fluorescence traces for each cell

get_stimulus_table (*self*, *stimulus_name*)

Return a stimulus table given a stimulus name

Notes

For more information, see: http://help.brain-map.org/display/observatory/Documentation?preview=/10616846/10813485/VisualCoding_VisualStimuli.pdf

get_stimulus_template (*self*, *stimulus_name*)

Return an array of the stimulus template for the specified stimulus.

Parameters

stimulus_name: string Must be one of the strings returned by *list_stimuli()*.

Returns

stimulus table: pd.DataFrame

list_stimuli (*self*)

Return a list of the stimuli presented in the experiment.

Returns

stimuli: list of strings

number_of_cells

Number of cells in the experiment

save_analysis_arrays (*self*, **datasets*)

save_analysis_dataframes (*self*, **tables*)

stimulus_search

`allensdk.core.brain_observatory_nwb_data_set.align_running_speed(dxcm, dx-`
`time, times-`
`tamps)`

If running speed timestamps differ from fluorescence timestamps, adjust by inserting NaNs to running speed.

Returns

tuple: dxcm, dxtime

`allensdk.core.brain_observatory_nwb_data_set.get_epoch_mask_list(st, threshold,`
`max_cuts=2)`

Convenience function to cut a stim table into multiple epochs

Parameters

- **st** – input stimtable
- **threshold** – threshold on the max duration of a subepoch
- **max_cuts** – maximum number of allowed epochs to cut into

Returns epoch_mask_list, a list of indices that define the start and end of sub-epochs

allensdk.core.cache_method_utilities module

class `allensdk.core.cache_method_utilities.CachedInstanceMethodMixin`

Bases: `object`

cache_clear (*self*)

Calls `cache_clear` method on all bound methods in this instance (where valid). Intended to clear calls cached with the `memoize` decorator. Note that this will also clear functions decorated with `lru_cache` and `lru_cache` in this class (or any other function with `cache_clear` attribute).

allensdk.core.cell_types_cache module

class `allensdk.core.cell_types_cache.CellTypesCache` (*cache=True*, *manifest_file=None*, *base_uri=None*)

Bases: `allensdk.api.warehouse_cache.cache.Cache`

Cache class for storing and accessing data from the Cell Types Database. By default, this class will cache any downloaded metadata or files in well known locations defined in a manifest file. This behavior can be disabled.

Parameters

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. `get_ephys_data(file_name='file.nwb')`).

manifest_file: string File name of the manifest to be read. Default is "cell_types_manifest.json".

Attributes

api: CellTypesApi instance The object used for making API queries related to the Cell Types Database

CELLS_KEY = 'CELLS'

EPHYS_DATA_KEY = 'EPHYS_DATA'

EPHYS_FEATURES_KEY = 'EPHYS_FEATURES'

EPHYS_SWEEPS_KEY = 'EPHYS_SWEEPS'

MANIFEST_VERSION = '1.1'

MARKER_KEY = 'MARKER'

MORPHOLOGY_FEATURES_KEY = 'MORPHOLOGY_FEATURES'

RECONSTRUCTION_KEY = 'RECONSTRUCTION'

build_manifest (*self, file_name*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

get_all_features (*self, dataframe=False, require_reconstruction=True*)

Download morphology and electrophysiology features for all cells and merge them into a single table.

Parameters

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

require_reconstruction: boolean Only return ephys and morphology features for cells that have reconstructions. Default True.

get_cells (*self, file_name=None, require_morphology=False, require_reconstruction=False, reporter_status=None, species=None, simple=True*)

Download metadata for all cells in the database and optionally return a subset filtered by whether or not they have a morphology or reconstruction.

Parameters

file_name: string File name to save/read the cell metadata as JSON. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

require_morphology: boolean Filter out cells that have no morphological images.

require_reconstruction: boolean Filter out cells that have no morphological reconstructions.

reporter_status: list Filter for cells that have one or more cell reporter statuses.

species: list Filter for cells that belong to one or more species. If None, return all. Must be one of [CellTypesApi.MOUSE, CellTypesApi.HUMAN].

get_ephys_data (*self, specimen_id, file_name=None*)

Download electrophysiology traces for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

NwbDataSet A class instance with helper methods for retrieving stimulus and response traces out of an NWB file.

get_ephys_features (*self, dataframe=False, file_name=None*)

Download electrophysiology features for all cells in the database.

Parameters

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

get_ephys_sweeps (*self, specimen_id, file_name=None*)

Download sweep metadata for a single cell specimen.

Parameters

specimen_id: int ID of a cell.

get_morphology_features (*self, dataframe=False, file_name=None*)

Download morphology features for all cells with reconstructions in the database.

Parameters

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

get_reconstruction (*self, specimen_id, file_name=None*)

Download and open a reconstruction for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the reconstruction SWC. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

Morphology A class instance with methods for accessing morphology compartments.

get_reconstruction_markers (*self, specimen_id, file_name=None*)

Download and open a reconstruction marker file for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the reconstruction marker. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

Morphology A class instance with methods for accessing morphology compartments.

class allensdk.core.cell_types_cache.**ReporterStatus**

Bases: object

Valid strings for filtering by cell reporter status.

INDETERMINATE = None

NA = None

NEGATIVE = 'negative'

POSITIVE = 'positive'

allensdk.core.dat_utilities module

class allensdk.core.dat_utilities.**DatUtilities**

Bases: object

classmethod **save_voltage** (*output_path*, *v*, *t*)

Save a single voltage output result into a simple text format.

The output file is one t v pair per line.

Parameters

output_path [string] file name for output

v [numpy array] voltage

t [numpy array] time

allensdk.core.exceptions module

exception allensdk.core.exceptions.**DataFrameIndexError** (*msg*,
caught_exception=None)

Bases: LookupError

More verbose method for accessing invalid rows or columns in a dataframe. Should be used when an index error is thrown on a dataframe.

exception allensdk.core.exceptions.**DataFrameKeyError** (*msg*, *caught_exception=None*)

Bases: LookupError

More verbose method for accessing invalid rows or columns in a dataframe. Should be used when a keyerror is thrown on a dataframe.

exception allensdk.core.exceptions.**MissingDataError**

Bases: ValueError

allensdk.core.h5_utilities module

`allensdk.core.h5_utilities.decode_bytes` (*bytes_dataset*, *encoding='UTF-8'*)

Convert the elements of a dataset of bytes to str

`allensdk.core.h5_utilities.h5_object_matcher_relname_in` (*relnames*,
h5_object_name,
h5_object)

Asks if an h5 object's relative name (the final section of its absolute name) is contained within a provided array

Parameters

relnames [array-like] Relative names against which to match

h5_object_name [str] Full name (path from origin) of h5 object

h5_object [h5py.Group, h5py.Dataset] Check this object's relative name

Returns

bool : whether the match succeeded

h5_object [h5py.group, h5py.Dataset] the argued object

`allensdk.core.h5_utilities.keyed_locate_h5_objects` (*matcher_cbs*, *h5_file*,
start_node=None)

Traverse an h5 file and build up a dictionary mapping supplied keys to located objects

`allensdk.core.h5_utilities.load_datasets_by_relnames` (*relnames*, *h5_file*, *start_node*)

A convenience function for finding and loading into memory one or more datasets from an h5 file

`allensdk.core.h5_utilities.locate_h5_objects` (*matcher_cb*, *h5_file*, *start_node=None*)

Traverse an h5 file and return objects matching supplied criteria

`allensdk.core.h5_utilities.traverse_h5_file` (*callback*, *h5_file*, *start_node=None*)

Traverse an h5 file and apply a callback to each node

allensdk.core.json_utilities module

class `allensdk.core.json_utilities.JsonComments`

Bases: object

classmethod `read_file` (*file_name*)

classmethod `read_string` (*json_string*)

classmethod `remove_comments` (*json_string*)

Strip single and multiline javascript-style comments.

Parameters

json [string] Json string with javascript-style comments.

Returns

string Copy of the input with comments removed.

Note: A JSON decoder MAY accept and ignore comments.

classmethod `remove_multiline_comments` (*json_string*)

Rebuild input without substrings matching */.../*.

Parameters

json_string [string] may or may not contain multiline comments.

Returns

string Copy of the input without the comments.

`allensdk.core.json_utilities.json_handler(obj)`

Used by write_json convert a few non-standard types to things that the json package can handle.

`allensdk.core.json_utilities.read(file_name)`

Shortcut reading JSON from a file.

`allensdk.core.json_utilities.read_url(url, method='POST')`

`allensdk.core.json_utilities.read_url_get(url)`

Transform a JSON contained in a file into an equivalent nested python dict.

Parameters

url [string] where to get the json.

Returns

dict Python version of the input

Note: if the input is a bare array or literal, for example,

the output will be of the corresponding type.

`allensdk.core.json_utilities.read_url_post(url)`

Transform a JSON contained in a file into an equivalent nested python dict.

Parameters

url [string] where to get the json.

Returns

dict Python version of the input

Note: if the input is a bare array or literal, for example,

the output will be of the corresponding type.

`allensdk.core.json_utilities.write(file_name, obj)`

Shortcut for writing JSON to a file. This also takes care of serializing numpy and data types.

`allensdk.core.json_utilities.write_string(obj)`

Shortcut for writing JSON to a string. This also takes care of serializing numpy and data types.

allensdk.core.mouse_connectivity_cache module

```
class allensdk.core.mouse_connectivity_cache.MouseConnectivityCache(resolution=None,
                                                                    cache=True,
                                                                    mani-
                                                                    fest_file=None,
                                                                    ccf_version=None,
                                                                    base_uri=None,
                                                                    ver-
                                                                    sion=None)
```

Bases: `allensdk.core.reference_space_cache.ReferenceSpaceCache`

Cache class for storing and accessing data related to the adult mouse Connectivity Atlas. By default, this class will cache any downloaded metadata or files in well known locations defined in a manifest file. This behavior can be disabled.

Parameters

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

ccf_version: string Desired version of the Common Coordinate Framework. This affects the annotation volume (`get_annotation_volume`) and structure masks (`get_structure_mask`). Must be one of (`MouseConnectivityApi.CCF_2015`, `MouseConnectivityApi.CCF_2016`). Default: `MouseConnectivityApi.CCF_2016`

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. `get_projection_density(file_name='file.nrrd')`).

manifest_file: string File name of the manifest to be read. Default is "mouse_connectivity_manifest.json".

Attributes

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

api: MouseConnectivityApi instance Used internally to make API queries.

```
ALIGNMENT3D_KEY = 'ALIGNMENT3D'
DATA_MASK_KEY = 'DATA_MASK'
DEFAULT_STRUCTURE_SET_IDS = (167587189,)
DEFORMATION_FIELD_HEADER_KEY = 'DEFORMATION_FIELD_HEADER'
DEFORMATION_FIELD_VOXEL_KEY = 'DEFORMATION_FIELD_VOXELS'
DFMFLD_RESOLUTIONS = (25,)
EXPERIMENTS_KEY = 'EXPERIMENTS'
INJECTION_DENSITY_KEY = 'INJECTION_DENSITY'
INJECTION_FRACTION_KEY = 'INJECTION_FRACTION'
MANIFEST_VERSION = 1.3
PROJECTION_DENSITY_KEY = 'PROJECTION_DENSITY'
STRUCTURE_UNIONIZES_KEY = 'STRUCTURE_UNIONIZES'
SUMMARY_STRUCTURE_SET_ID = 167587189
add_manifest_paths(self, manifest_builder)
    Construct a manifest for this Cache class and save it in a file.
```

Parameters

file_name: string File location to save the manifest.

default_structure_ids

filter_experiments (*self, experiments, cre=None, injection_structure_ids=None*)

Take a list of experiments and filter them by cre status and injection structure.

Parameters

cre: boolean or list If True, return only cre-positive experiments. If False, return only cre-negative experiments. If None, return all experiments. If list, return all experiments with cre line names in the supplied list. Default None.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If None, return all experiments. Default None.

filter_structure_unionizes (*self, unionizes, is_injection=None, structure_ids=None, include_descendants=False, hemisphere_ids=None*)

Take a list of unionizes and return a subset of records filtered by injection status, structure, and hemisphere.

Parameters

is_injection: boolean If True, only return unionize records that disregard non-injection pixels. If False, only return unionize records that disregard injection pixels. If None, return all records. Default None.

structure_ids: list Only return unionize records for a set of structures. If None, return all records. Default None.

include_descendants: boolean Include all descendant records for specified structures. Default False.

hemisphere_ids: list Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If None, include all records [1, 2, 3]. Default None.

get_affine_parameters (*self, section_data_set_id, direction='trv', file_name=None*)

Extract the parameters of the 3D affine tranformation mapping this section data set's image-space stack to CCF-space (or vice-versa).

Parameters

section_data_set_id [int] download the parameters for this data set.

direction [str, optional]

Valid options are:

trv ["transform from reference to volume". Maps CCF points to image space points. If you are] resampling data into CCF, this is the direction you want.

tvr : "transform from volume to reference". Maps image space points to CCF points.

file_name [str] If provided, store the downloaded file here.

Returns

alignment [numpy.ndarray]

4 X 3 matrix. In order to transform a point [X_1, X_2, X_3] run `np.dot([X_1, X_2, X_3, 1], alignment)`. In

to build a SimpleITK affine transform run: `transform = sitk.AffineTransform(3)`
`transform.SetParameters(alignment.flatten())`

get_data_mask (*self, experiment_id, file_name=None*)

Read a data mask volume for a single experiment. Download it first if it doesn't exist. Data mask is a binary mask of voxels that have valid data. Only use valid data in analysis!

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to `section_data_set_id` in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If `file_name` is `None`, the `file_name` will be pulled out of the manifest. Default is `None`.

get_deformation_field (*self*, *section_data_set_id*, *header_path=None*, *voxel_path=None*)

Extract the local alignment parameters for this dataset. This a 3D vector image (3 components) describing a deformable local mapping from CCF voxels to this section data set's affine-aligned image stack.

Parameters

section_data_set_id [int]

Download the deformation field for this data set

header_path [str, optional] If supplied, the deformation field header will be downloaded to this path.

voxel_path [str, optional] If supplied, the deformation field voxels will be downloaded to this path.

Returns

numpy.ndarray : 3D X 3 component vector array (origin 0, 0, 0; 25-micron isotropic resolution) defining a deformable transformation from CCF-space to affine-transformed image space.

get_experiment_structure_unionizes (*self*, *experiment_id*, *file_name=None*, *is_injection=None*, *structure_ids=None*, *include_descendants=False*, *hemisphere_ids=None*)

Retrieve the structure unionize data for a specific experiment. Filter by structure, injection status, and hemisphere.

Parameters

experiment_id: int ID of the experiment of interest. Corresponds to `section_data_set_id` in the API.

file_name: string File name to save/read the experiments list. If `file_name` is `None`, the `file_name` will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is `None`.

is_injection: boolean If `True`, only return unionize records that disregard non-injection pixels. If `False`, only return unionize records that disregard injection pixels. If `None`, return all records. Default `None`.

structure_ids: list Only return unionize records for a specific set of structures. If `None`, return all records. Default `None`.

include_descendants: boolean Include all descendant records for specified structures. Default `False`.

hemisphere_ids: list Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If `None`, include all records [1, 2, 3]. Default `None`.

get_experiments (*self*, *dataframe=False*, *file_name=None*, *cre=None*, *injection_structure_ids=None*)

Read a list of experiments that match certain criteria. If caching is enabled, this will save the whole (unfiltered) list of experiments to a file.

Parameters

dataframe: boolean Return the list of experiments as a Pandas DataFrame. If False, return a list of dictionaries. Default False.

file_name: string File name to save/read the structures table. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

cre: boolean or list If True, return only cre-positive experiments. If False, return only cre-negative experiments. If None, return all experiments. If list, return all experiments with cre line names in the supplied list. Default None.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If None, return all experiments. Default None.

get_injection_density (*self*, *experiment_id*, *file_name=None*)

Read an injection density volume for a single experiment. Download it first if it doesn't exist. Injection density is the proportion of projecting pixels in a grid voxel only including pixels that are part of the injection site in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_injection_fraction (*self*, *experiment_id*, *file_name=None*)

Read an injection fraction volume for a single experiment. Download it first if it doesn't exist. Injection fraction is the proportion of pixels in the injection site in a grid voxel in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_projection_density (*self*, *experiment_id*, *file_name=None*)

Read a projection density volume for a single experiment. Download it first if it doesn't exist. Projection density is the proportion of projecting pixels in a grid voxel in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_projection_matrix (*self*, *experiment_ids*, *projection_structure_ids=None*, *hemisphere_ids=None*, *parameter='projection_volume'*, *dataframe=False*)

```
get_structure_unionizes (self, experiment_ids, is_injection=None, structure_ids=None, include_descendants=False, hemisphere_ids=None)
```

Get structure unionizes for a set of experiment IDs. Filter the results by injection status, structure, and hemisphere.

Parameters

experiment_ids: list List of experiment IDs. Corresponds to `section_data_set_id` in the API.

is_injection: boolean If True, only return unionize records that disregard non-injection pixels. If False, only return unionize records that disregard injection pixels. If None, return all records. Default None.

structure_ids: list Only return unionize records for a specific set of structures. If None, return all records. Default None.

include_descendants: boolean Include all descendant records for specified structures. Default False.

hemisphere_ids: list Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If None, include all records [1, 2, 3]. Default None.

```
rank_structures (self, experiment_ids, is_injection, structure_ids=None, hemisphere_ids=None,  
                 rank_on='normalized_projection_volume', n=5, threshold=0.01)
```

Produces one or more (per experiment) ranked lists of brain structures, using a specified data field.

Parameters

experiment_ids [list of int] Obtain injection_structures for these experiments.

is_injection [boolean] Use data from only injection (or non-injection) unionizes.

structure_ids [list of int, optional] Consider only these structures. It is a good idea to make sure that these structures are not spatially overlapping; otherwise your results will contain redundant information. Defaults to the summary structures - a brain-wide list of nonoverlapping mid-level structures.

hemisphere_ids [list of int, optional] Consider only these hemispheres (1: left, 2: right, 3: both). Like with structures, you might get redundant results if you select overlapping options. Defaults to [1, 2].

rank_on [str, optional] Rank unionize data using this field (descending). Defaults to `normalized_projection_volume`.

n [int, optional] Return only the top n structures.

threshold [float, optional] Consider only records whose data value - specified by the `rank_on` parameter - exceeds this value.

Returns

list : Each element (1 for each input experiment) is a list of dictionaries. The dictionaries describe the top injection structures in descending order. They are specified by their structure and hemisphere id fields and additionally report the value specified by the `rank_on` parameter.

`allensdk.core.nwb_data_set` module

```
class allensdk.core.nwb_data_set.NwbDataSet (file_name, spike_time_key=None)
```

Bases: object

A very simple interface for extracting electrophysiology data from an NWB file.

DEPRECATED_SPIKE_TIMES = 'aibs_spike_times'

SPIKE_TIMES = 'spike_times'

fill_sweep_responses (*self*, *fill_value*=0.0, *sweep_numbers*=None, *extend_experiment*=False)

Fill sweep response arrays with a single value.

Parameters

fill_value: float Value used to fill sweep response array

sweep_numbers: list List of integer sweep numbers to be filled (default all sweeps)

extend_experiment: bool If True, extend experiment epoch length to the end of the sweep (undo any truncation)

get_experiment_sweep_numbers (*self*)

Get all of the sweep numbers for experiment epochs in the file, not including test sweeps.

get_pipeline_version (*self*)

Returns the AI pipeline version number, stored in the metadata field 'generated_by'. If that field is missing, version 0.0 is returned.

Returns

int tuple: (major, minor)

get_spike_times (*self*, *sweep_number*, *key*=None)

Return any spike times stored in the NWB file for a sweep.

Parameters

sweep_number: int index to access

key [string] label where the spike times are stored (default NwbDataSet.SPIKE_TIMES)

Returns

list list of spike times in seconds relative to the start of the sweep

get_sweep (*self*, *sweep_number*)

Retrieve the stimulus, response, index_range, and sampling rate for a particular sweep. This method hides the NWB file's distinction between a "Sweep" and an "Experiment". An experiment is a subset of a sweep that excludes the initial test pulse. It also excludes any erroneous response data at the end of the sweep (usually for ramp sweeps, where recording was terminated mid-stimulus).

Some sweeps do not have an experiment, so full data arrays are returned. Sweeps that have an experiment return full data arrays (include the test pulse) with any erroneous data trimmed from the back of the sweep.

Parameters

sweep_number: int

Returns

dict A dictionary with 'stimulus', 'response', 'index_range', and 'sampling_rate' elements. The index range is a 2-tuple where the first element indicates the end of the test pulse and the second index is the end of valid response data.

get_sweep_metadata (*self*, *sweep_number*)

Retrieve the sweep level metadata associated with each sweep. Includes information on stimulus parameters like its name and amplitude as well as recording quality metadata, like access resistance and seal quality.

Parameters**sweep_number: int****Returns****dict** A dictionary with 'aibs_stimulus_amplitude_pa', 'aibs_stimulus_name', 'gain', 'initial_access_resistance', 'seal' elements. These specific fields are ones encoded in the original AIBS in vitro .nwb files.**get_sweep_numbers** (*self*)

Get all of the sweep numbers in the file, including test sweeps.

set_spike_times (*self*, *sweep_number*, *spike_times*, *key=None*)

Set or overwrite the spikes times for a sweep.

Parameters**sweep_number** [int] index to access**key** [string] where the times are stored (default NwbDataSet.SPIKE_TIME)**spike_times: np.array** array of spike times in seconds**set_sweep** (*self*, *sweep_number*, *stimulus*, *response*)

Overwrite the stimulus or response of an NWB file. If the supplied arrays are shorter than stored arrays, they are padded with zeros to match the original data size.

Parameters**sweep_number: int****stimulus: np.array** Overwrite the stimulus with this array. If None, stimulus is unchanged.**response: np.array** Overwrite the response with this array. If None, response is unchanged.**allensdk.core.obj_utilities module**`allensdk.core.obj_utilities.parse_obj` (*lines*)

Parse a wavefront obj file into a triplet of vertices, normals, and faces. This parser is specific to obj files generated from our annotation volumes

Parameters**lines** [list of str] Lines of input obj file**Returns****vertices** [np.ndarray] Dimensions are (nSamples, nCoordinates=3). Locations in the reference space of vertices**vertex_normals** [np.ndarray] Dimensions are (nSample, nElements=3). Vectors normal to vertices.**face_vertices** [np.ndarray] Dimensions are (sample, nVertices=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertices that make up each face.**face_normals** [np.ndarray] Dimensions are (sample, nNormals=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertex normals that make up each face.

Notes

This parser is specialized to the obj files that the Allen Institute for Brain Science generates from our own structure annotations.

```
allensdk.core.obj_utilities.read_obj(path)
```

allensdk.core.ontology module

```
class allensdk.core.ontology.Ontology(df)
    Bases: object
```

Note: Deprecated from 0.12.5 *Ontology* has been replaced by *StructureTree*.

```
get_child_ids(self, structure_ids)
```

Find the set of ids that are immediate children of one or more structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

set Set of child structure ids

```
get_children(self, structure_ids)
```

Find the set of structures that are immediate children of one or more structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

pandas.DataFrame Set of child structures

```
get_descendant_ids(self, structure_ids)
```

Find the set of the ids of structures that are descendants of one or more structures. The returned set will include the input structure ids.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

set Set of descendant structure ids.

```
get_descendants(self, structure_ids)
```

Find the set of structures that are descendants of one or more structures. The returned set will include the input structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

pandas.DataFrame Set of descendant structures.

structure_descends_from (*self*, *child_id*, *parent_id*)

Return whether one structure id is a descendant of another structure id.

allensdk.core.ophys_experiment_session_id_mapping module

allensdk.core.reference_space module

class `allensdk.core.reference_space.ReferenceSpace` (*structure_tree*, *annotation*, *resolution*)

Bases: `object`

static check_and_write (*base_dir*, *structure_id*, *fn*)

A many_structure_masks callback that writes the mask to a nrrd file if the file does not already exist.

check_coverage (*self*, *structure_ids*, *domain_mask*)

Determines whether a spatial domain is completely covered by structures in a set.

Parameters

structure_ids [list of int] Specifies the set of structures to check.

domain_mask [numpy ndarray] Same shape as annotation. 1 inside the mask, 0 out. Specifies spatial domain.

Returns

numpy ndarray : 1 where voxels are missing from the candidate, 0 where the candidate exceeds the domain

direct_voxel_counts (*self*)

Determines the number of voxels directly assigned to one or more structures.

Returns

dict : Keys are structure ids, values are the number of voxels directly assigned to those structures.

direct_voxel_map

downsample (*self*, *target_resolution*)

Obtain a smaller reference space by downsampling

Parameters

target_resolution [tuple of numeric] Resolution in microns of the output space.

interpolator [string] Method used to interpolate the volume. Currently only 'nearest' is supported

Returns

ReferenceSpace : A new ReferenceSpace with the same structure tree and a downsampled annotation.

export_itksnap_labels (*self*, *id_type*=<class 'numpy.uint16'>, *label_description_kwargs*=None)

Produces itksnap labels, remapping large ids if needed.

Parameters

id_type [np.integer, optional] Used to determine the type of the output annotation and whether ids need to be remapped to smaller values.

label_description_kwargs [dict, optional] Keyword arguments passed to StructureTree.export_label_description

Returns

np.ndarray : Annotation volume, remapped if needed

pd.DataFrame label_description dataframe

get_slice_image (*self*, *axis*, *position*, *cmap=None*)

Produce a AxBy3 RGB image from a slice in the annotation

Parameters

axis [int] Along which to slice the annotation volume. 0 is coronal, 1 is horizontal, and 2 is sagittal.

position [int] In microns. Take the slice from this far along the specified axis.

cmap [dict, optional] Keys are structure ids, values are rgb triplets. Defaults to structure rgb_triplets.

Returns

np.ndarray : RGB image array.

Notes

If you assign a custom colormap, make sure that you take care of the background in addition to the structures.

make_structure_mask (*self*, *structure_ids*, *direct_only=False*)

Return an indicator array for one or more structures

Parameters

structure_ids [list of int] Make a mask that indicates the union of these structures' voxels

direct_only [bool, optional] If True, only include voxels directly assigned to a structure in the mask. Otherwise include voxels assigned to descendants.

Returns

numpy ndarray : Same shape as annotation. 1 inside mask, 0 outside.

many_structure_masks (*self*, *structure_ids*, *output_cb=None*, *direct_only=False*)

Build one or more structure masks and do something with them

Parameters

structure_ids [list of int] Specify structures to be masked

output_cb [function, optional] Must have the following signature: `output_cb(structure_id, fn)`. On each requested id, `fn` will be carried to make a mask for that id. Defaults to returning the structure id and mask.

direct_only [bool, optional] If True, only include voxels directly assigned to a structure in the mask. Otherwise include voxels assigned to descendants.

Yields

Return values of **output_cb** called on each **structure_id**, **structure_mask** pair.

Notes

output_cb is called on every yield, so any side-effects (such as writing to a file) will be carried out regardless of what you do with the return values. You do actually have to iterate through the output, though.

remove_unassigned (*self*, *update_self=True*)

Obtains a structure tree consisting only of structures that have at least one voxel in the annotation.

Parameters

update_self [bool, optional] If True, the contained structure tree will be replaced,

Returns

list of dict : elements are filtered structures

static return_mask_cb (*structure_id*, *fn*)

A basic callback for many_structure_masks

total_voxel_counts (*self*)

Determines the number of voxels assigned to a structure or its descendants

Returns

dict : Keys are structure ids, values are the number of voxels assigned to structures' descendants.

total_voxel_map

validate_structures (*self*, *structure_ids*, *domain_mask*)

Determines whether a set of structures produces an exact and nonoverlapping tiling of a spatial domain

Parameters

structure_ids [list of int] Specifies the set of structures to check.

domain_mask [numpy ndarray] Same shape as annotation. 1 inside the mask, 0 out. Specifies spatial domain.

Returns

set : Ids of structures that are the ancestors of other structures in the supplied set.

numpy ndarray : Indicator for missing voxels.

write_itksnap_labels (*self*, *annotation_path*, *label_path*, ***kwargs*)

Generate a label file (nrrd) and a label_description file (csv) for use with ITKSnap

Parameters

annotation_path [str] write generated label file here

label_path [str] write generated label_description file here

****kwargs** : will be passed to self.export_itksnap_labels

allensdk.core.reference_space_cache module

```
class allensdk.core.reference_space_cache.ReferenceSpaceCache (resolution, reference_space_key,
                                                                **kwargs)
```

Bases: `allensdk.api.warehouse_cache.cache.Cache`

ANNOTATION_KEY = 'ANNOTATION'


```

MANIFEST_VERSION = 1.2
REFERENCE_SPACE_VERSION_KEY = 'REFERENCE_SPACE_VERSION'
STRUCTURES_KEY = 'STRUCTURES'
STRUCTURE_MASK_KEY = 'STRUCTURE_MASK'
STRUCTURE_MESH_KEY = 'STRUCTURE_MESH'
STRUCTURE_TREE_KEY = 'STRUCTURE_TREE'
TEMPLATE_KEY = 'TEMPLATE'

```

add_manifest_paths (*self*, *manifest_builder*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

get_annotation_volume (*self*, *file_name=None*)

Read the annotation volume. Download it first if it doesn't exist.

Parameters

file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If *file_name* is None, the *file_name* will be pulled out of the manifest. Default is None.

get_reference_space (*self*, *structure_file_name=None*, *annotation_file_name=None*)

Build a ReferenceSpace from this cache's annotation volume and structure tree. The ReferenceSpace does operations that relate brain structures to spatial domains.

Parameters

structure_file_name: string File name to save/read the structures table. If *file_name* is None, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

annotation_file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If *file_name* is None, the *file_name* will be pulled out of the manifest. Default is None.

get_structure_mask (*self*, *structure_id*, *file_name=None*, *annotation_file_name=None*)

Read a 3D numpy array shaped like the annotation volume that has non-zero values where voxels belong to a particular structure. This will take care of identifying substructures.

Parameters

structure_id: int ID of a structure.

file_name: string File name to store the structure mask. If it already exists, it will be read from this file. If *file_name* is None, the *file_name* will be pulled out of the manifest. Default is None.

annotation_file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If *file_name* is None, the *file_name* will be pulled out of the manifest. Default is None.

Notes

This method downloads structure masks from the Allen Institute. To make your own locally, see `ReferenceSpace.many_structure_masks`.

get_structure_mesh (*self, structure_id, file_name=None*)

Obtain a 3D mesh specifying the surface of an annotated structure.

Parameters

structure_id: int ID of a structure.

file_name: string File name to store the structure mesh. If it already exists, it will be read from this file. If *file_name* is *None*, the *file_name* will be pulled out of the manifest. Default is *None*.

Returns

vertices [np.ndarray] Dimensions are (nSamples, nCoordinates=3). Locations in the reference space of vertices

vertex_normals [np.ndarray] Dimensions are (nSample, nElements=3). Vectors normal to vertices.

face_vertices [np.ndarray] Dimensions are (sample, nVertices=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertices that make up each face.

face_normals [np.ndarray] Dimensions are (sample, nNormals=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertex normals that make up each face.

Notes

These meshes are meant for 3D visualization and as such have been smoothed. If you are interested in performing quantitative analyses, we recommend that you use the structure masks instead.

get_structure_tree (*self, file_name=None, structure_graph_id=1*)

Read the list of adult mouse structures and return an *StructureTree* instance.

Parameters

file_name: string File name to save/read the structures table. If *file_name* is *None*, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is *None*.

structure_graph_id: int Build a tree using structure only from the identified structure graph.

get_template_volume (*self, file_name=None*)

Read the template volume. Download it first if it doesn't exist.

Parameters

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If *file_name* is *None*, the *file_name* will be pulled out of the manifest. Default is *None*.

classmethod validate_structure_id (*structure_id*)

classmethod validate_structure_ids (*structure_ids*)

allensdk.core.simple_tree module

class allensdk.core.simple_tree.**SimpleTree** (*nodes, node_id_cb, parent_id_cb*)

Bases: *object*

ancestor_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' ancestors

Parameters

node_ids [list of hashable] Items are ids of nodes whose ancestors you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' ancestors' ids.

Notes

Given the tree: A -> B -> C

‘-> D

The ancestors of C are [C, B, A]. The ancestors of A are [A]. The ancestors of D are [D, A]

ancestors (*self*, *node_ids*)

Get one or mode nodes' ancestor nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose ancestors will be found.

Returns

list of list of dict : Items are lists of ancestor nodes corresponding to argued ids.

child_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' children

Parameters

node_ids [list of hashable] Items are ids of nodes whose children you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' children's ids.

children (*self*, *node_ids*)

Get one or mode nodes' child nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose children will be found.

Returns

list of list of dict : Items are lists of child nodes corresponding to argued ids.

descendant_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' descendants

Parameters

node_ids [list of hashable] Items are ids of nodes whose descendants you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' descendants' ids.

Notes

Given the tree: A -> B -> C

‘-> D

The descendants of A are [B, C, D]. The descendants of C are [].

descendants (*self*, *node_ids*)

Get one or more nodes' descendant nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose descendants will be found.

Returns

list of list of dict : Items are lists of descendant nodes corresponding to argued ids.

filter_nodes (*self*, *criterion*)

Obtain a list of nodes filtered by some criterion

Parameters

criterion [function | node dict => bool] Only nodes for which criterion returns true will be returned.

Returns

list of dict : Items are node dictionaries that passed the filter.

node (*self*, *node_ids=None*)

node_ids (*self*)

Obtain the node ids of each node in the tree

Returns

list : elements are node ids

nodes (*self*, *node_ids=None*)

Get one or more nodes' full dictionaries from their ids.

Parameters

node_ids [list of hashable] Items are ids of nodes to be returned. Default is all.

Returns

list of dict : Items are nodes corresponding to argued ids.

nodes_by_property (*self*, *key*, *values*, *to_fn=None*)

Get nodes by a specified property

Parameters

key [hashable or function] The property used for lookup. Should be unique. If a function, will be invoked on each node.

values [list] Select matching elements from the lookup.

to_fn [function, optional] Defines the outputs, on a per-node basis. Defaults to returning the whole node.

Returns

list : outputs, 1 for each input value.

parent (*self*, *node_ids*)

parent_id (*self*, *node_ids*)

parent_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' parents

Parameters

node_ids [list of hashable] Items are ids of nodes whose parents you wish to find.

Returns

list of hashable : Items are ids of input nodes' parents in order.

parents (*self*, *node_ids*)

Get one or more nodes' parent nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose parents will be found.

Returns

list of dict : Items are parents of nodes corresponding to argued ids.

value_map (*self*, *from_fn*, *to_fn*)

Obtain a look-up table relating a pair of node properties across nodes

Parameters

from_fn [function | node dict => hashable value] The keys of the output dictionary will be obtained by calling from_fn on each node. Should be unique.

to_fn [function | node_dict => value] The values of the output function will be obtained by calling to_fn on each node.

Returns

dict : Maps the node property defined by from_fn to the node property defined by to_fn across nodes.

allensdk.core.sitk_utilities module

`allensdk.core.sitk_utilities.fix_array_dimensions` (*array*, *ncomponents=1*)

Convenience function that reorders ndarray dimensions for io with SimpleITK

Parameters

array [np.ndarray] The array to be reordered

ncomponents [int, optional] Number of components per pixel, default 1.

Returns

np.ndarray : Reordered array

`allensdk.core.sitk_utilities.get_sitk_image_information` (*image*)

Extract information about a SimpleITK image

Parameters

image [sitk.Image] Extract information about this image.

Returns

dict : Extracted information. Includes spacing, origin, size, direction, and number of components per pixel

`allensdk.core.sitk_utilities.read_ndarray_with_sitk(path)`

Read a numpy array from a file using SimpleITK

Parameters

path [str] Read from this path

Returns

image [np.ndarray] Obtained array

information [dict] Additional information about the array

`allensdk.core.sitk_utilities.set_sitk_image_information(image, information)`

Set information on a SimpleITK image

Parameters

image [sitk.Image] Set information on this image.

information [dict] Stores information to be set. Supports spacing, origin, direction. Also checks (but cannot set) size and number of components per pixel

`allensdk.core.sitk_utilities.write_ndarray_with_sitk(array, path, **information)`

Write a numpy array to a file using SimpleITK

Parameters

array [np.ndarray] Array to be written.

path [str] Write to here

****information** [dict] Contains additional information to be stored in the image file. See `set_sitk_image_information` for more information.

allensdk.core.structure_tree module

`class allensdk.core.structure_tree.StructureTree(nodes)`

Bases: `allensdk.core.simple_tree.SimpleTree`

static clean_structures (structures, whitelist=None, data_transforms=None, renames=None)

Convert structures_with_sets query results into a form that can be used to construct a StructureTree

Parameters

structures [list of dict] Each element describes a structure. Should have a structure id path field (str values) and a structure_sets field (list of dict).

whitelist [list of str, optional] Only these fields will be included in the final structure record. Default is the output of `StructureTree.whitelist`.

data_transforms [dict, optional] Keys are str field names. Values are functions which will be applied to the data associated with those fields. Default is to map colors from hex to rgb and convert the structure id path to a list of int.

renames [dict, optional] Controls the field names that appear in the output structure records. Default is to map 'color_hex_triplet' to 'rgb_triplet'.

Returns

list of dict : structures, after conversion of structure_id_path and structure_sets

static collect_sets (*structure*)

Structure sets may be specified by full records or id. This method collects all of the structure set records/ids in a structure record and replaces them with a single list of id records.

static data_transforms ()

export_label_description (*self*, *alphas=None*, *exclude_label_vis=None*, *exclude_mesh_vis=None*, *label_key='acronym'*)

Produces an itksnap label_description table from this structure tree

Parameters

alphas [dict, optional] Maps structure ids to alpha levels. Optional - will only use provided ids.

exclude_label_vis [list, optional] The structures denoted by these ids will not be visible in ITKSnap.

exclude_mesh_vis [list, optional] The structures denoted by these ids will not have visible meshes in ITKSnap.

label_key: str, optional Use this column for display labels.

Returns

pd.DataFrame : Contains data needed for loading as an ITKSnap label description file.

get_ancestor_id_map (*self*)

Get a dictionary mapping structure ids to ancestor ids across all nodes.

Returns

dict : Keys are structure ids. Values are lists of ancestor ids.

get_colormap (*self*)

Get a dictionary mapping structure ids to colors across all nodes.

Returns

dict : Keys are structure ids. Values are RGB lists of integers.

get_id_acronym_map (*self*)

Get a dictionary mapping structure acronyms to ids across all nodes.

Returns

dict : Keys are structure acronyms. Values are structure ids.

get_name_map (*self*)

Get a dictionary mapping structure ids to names across all nodes.

Returns

dict : Keys are structure ids. Values are structure name strings.

get_structure_sets (*self*)

Lists all unique structure sets that are assigned to at least one structure in the tree.

Returns

list of int : Elements are ids of structure sets.

get_structures_by_acronym (*self*, *acronyms*)

Obtain a list of brain structures from their acronyms

Parameters

names [list of str] Get structures corresponding to these acronyms.

Returns

list of dict : Each item describes a structure.

get_structures_by_id (*self*, *structure_ids*)

Obtain a list of brain structures from their structure ids

Parameters

structure_ids [list of int] Get structures corresponding to these ids.

Returns

list of dict : Each item describes a structure.

get_structures_by_name (*self*, *names*)

Obtain a list of brain structures from their names,

Parameters

names [list of str] Get structures corresponding to these names.

Returns

list of dict : Each item describes a structure.

get_structures_by_set_id (*self*, *structure_set_ids*)

Obtain a list of brain structures from by the sets that contain them.

Parameters

structure_set_ids [list of int] Get structures belonging to these structure sets.

Returns

list of dict : Each item describes a structure.

has_overlaps (*self*, *structure_ids*)

Determine if a list of structures contains structures along with their ancestors

Parameters

structure_ids [list of int] Check this set of structures for overlaps

Returns

set : Ids of structures that are the ancestors of other structures in the supplied set.

static hex_to_rgb (*hex_color*)

Convert a hexadecimal color string to a uint8 triplet

Parameters

hex_color [string] Must be 6 characters long, unless it is 7 long and the first character is #. If hex_color is a triplet of int, it will be returned unchanged.

Returns

list of int : 3 characters long - 1 per two characters in the input string.

static path_to_list (*path*)

Structure id paths are sometimes formatted as "/"-seperated strings. This method converts them to a list of integers, if needed.

static renames ()

structure_descends_from (*self*, *child_id*, *parent_id*)

Tests whether one structure descends from another.

Parameters

child_id [int] Id of the putative child structure.

parent_id [int] Id of the putative parent structure.

Returns

bool : True if the structure specified by **child_id** is a descendant of the one specified by **parent_id**. Otherwise False.

static whitelist()

allensdk.core.swc module

class allensdk.core.swc.**Compartment** (*args, **kwargs)

Bases: dict

A dictionary class storing information about a single morphology node

print_node (self)

print out compartment information with field names

class allensdk.core.swc.**Marker** (*args, **kwargs)

Bases: dict

Simple dictionary class for handling reconstruction marker objects.

CUT_DENDRITE = 10

NO_RECONSTRUCTION = 20

SPACING = [0.1144, 0.1144, 0.28]

class allensdk.core.swc.**Morphology** (compartment_list=None, compartment_index=None)

Bases: object

Keep track of the list of compartments in a morphology and provide a few helper methods (soma, tree information, pruning, etc).

APICAL_DENDRITE = 4

AXON = 2

BASAL_DENDRITE = 3

DENDRITE = 3

NODE_TYPES = [1, 2, 3, 3, 4]

SOMA = 1

append (self, node_list)

Add additional nodes to this Morphology. Those nodes must originate from another morphology object.

Parameters

node_list: list of Morphology nodes

apply_affine (self, aff, scale=None)

Apply an affine transform to all compartments in this morphology. Node radius is adjusted as well.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

change_parent (*self, child, parent*)

Change the parent of a node. The child node is adjusted to point to the new parent, the child is taken off of the previous parent's child list, and it is added to the new parent's child list.

Parameters

child: integer or Morphology Object The ID of the child node, or the child node itself

parent: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

Nothing

children_of (*self, seg*)

Returns a list of the children of the specified node

Parameters

seg: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

A list of the child morphology objects. If the ID of the parent node is invalid, None is returned.

compartment_index

Return the compartment index. This is a property to ensure that the compartment list and compartment index are in sync.

compartment_index_by_type (*self, compartment_type*)

Return an dictionary of compartments indexed by id that all have a particular compartment type.

Parameters

compartment_type: int Desired compartment type

Returns

A dictionary of Morphology Objects, indexed by ID

compartment_list

Return the compartment list. This is a property to ensure that the compartment list and compartment index are in sync.

compartment_list_by_type (*self, compartment_type*)

Return an list of all compartments having the specified compartment type.

Parameters

compartment_type: int Desired compartment type

Returns

A list of of Morphology Objects

convert_type (*self*, *old_type*, *new_type*)

Converts all compartments from one type to another. Nodes of the original type are not affected so this procedure can also be used as a merge procedure.

Parameters

old_type: enum The compartment type to be changed. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

new_type: enum The target compartment type. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

delete_tree (*self*, *n*)

Delete tree, and all of its compartments, from the morphology.

Parameters

n: Integer The tree number to delete

find (*self*, *x*, *y*, *z*, *dist*, *node_type=None*)

Returns a list of Morphology Objects located within ‘dist’ of coordinate (x,y,z). If node_type is specified, the search will be constrained to return only nodes of that type.

Parameters

x, y, z: float The x,y,z coordinates from which to search around

dist: float The search radius

node_type: enum (optional) One of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE or APICAL_DENDRITE

Returns

A list of all Morphology Objects matching the search criteria

node (*self*, *n*)

Returns the morphology node having the specified ID.

Parameters

n: integer ID of desired node

Returns

A morphology object having the specified ID, or None if such a node doesn’t exist

num_nodes

Return the number of compartments in the morphology.

num_trees

Return the number of trees in the morphology. A tree is defined as everything following from a single root compartment.

parent_of (*self*, *seg*)

Returns parent of the specified node.

Parameters

seg: integer or Morphology Object The ID of the child node, or the child node itself

Returns

A morphology object, or None if no parent exists or if the specified node ID doesn't exist

root

[deprecated] Returns root node of soma, if present. Use 'soma' instead of 'root'

save (*self*, *file_name*)

Write this morphology out to an SWC file

Parameters

file_name: string desired name of your SWC file

soma

Returns root node of soma, if present

sparsify (*self*, *modulo*, *compress_ids=False*)

Return a new Morphology object that has a given number of non-leaf, non-root nodes removed. IDs can be reassigned so as to be continuous.

Parameters

modulo: int keep 1 out of every modulo nodes.

compress_ids: boolean Reassign ids so that ids are continuous (no missing id numbers).

Returns

Morphology A new morphology instance

strip_all_other_types (*self*, *node_type*, *keep_soma=True*)

Strips everything from the morphology except for the specified type. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The compartment type to keep in the morphology. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

keep_soma: Boolean (optional) True (default) if soma nodes should remain in the morphology, and False if the soma should also be stripped

strip_type (*self*, *node_type*)

Strips all compartments of the specified type from the morphology. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The compartment type to strip from the morphology. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

stumpify_axon (*self*, *count=10*)

Remove all axon compartments except the first 'count' nodes, as counted from the connected axon root.

Parameters

count: Integer The length of the axon 'stump', in number of compartments

tree (*self*, *n*)

Returns a list of all Morphology Nodes within the specified tree. A tree is defined as a fully connected graph of nodes. Each tree has exactly one root.

Parameters

n: integer ID of desired tree

Returns

A list of all morphology objects in the specified tree, or None
if the tree doesn't exist

write (*self*, *file_name*)

`allensdk.core.swc.read_marker_file` (*file_name*)

read in a marker file and return a list of dictionaries

`allensdk.core.swc.read_swc` (*file_name*, *columns*='NOT_USED', *numeric_columns*='NOT_USED')

Read in an SWC file and return a Morphology object.

Parameters

file_name: string SWC file name.

Returns

Morphology A Morphology instance.

allensdk.core.typing module

class `allensdk.core.typing.SupportsStr`

Bases: `typing._Protocol`

Classes that support the `__str__` method

Module contents

6.1.5 allensdk.ephys package

Submodules

allensdk.ephys.ephys_extractor module

class `allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` (*ramps_ext*,
short_squares_ext,
long_squares_ext,
subthresh_min_amp=-100)

Bases: `object`

SAG_TARGET = -100.0

SUBTHRESH_MAX_AMP = 0

as_dict (*self*)

Create dict of cell features.

cell_features (*self*)

long_squares_features (*self*, *option*=None)

long_squares_stim_amps (*self*, *option*=None)

process (*self*, *keys=None*)

Processes features. Can take a specific key (or set of keys) to do a subset of processing.

ramps_features (*self*, *all=False*)

short_squares_features (*self*)

```
class allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor (t=None,  
                                                             v=None,  
                                                             i=None,  
                                                             start=None,  
                                                             end=None,  
                                                             filter=10.0,  
                                                             dv_cutoff=20.0,  
                                                             max_interval=0.005,  
                                                             min_height=2.0,  
                                                             min_peak=-  
                                                             30.0,  
                                                             thresh_frac=0.05,  
                                                             base-  
                                                             line_interval=0.1,  
                                                             base-  
                                                             line_detect_thresh=0.3,  
                                                             id=None)
```

Bases: object

Feature calculation for a sweep (voltage and/or current time series).

as_dict (*self*)

Create dict of features and spikes.

burst_metrics (*self*)

Find bursts and return max “burstiness” index (normalized max rate in burst vs out).

Returns

max_burstiness_index [max “burstiness” index across detected bursts]

num_bursts [number of bursts detected]

delay_metrics (*self*)

Calculates ratio of latency to dominant time constant of rise before spike

Returns

delay_ratio [ratio of latency to tau (higher means more delay)]

tau [dominant time constant of rise before spike]

estimate_sag (*self*, *peak_width=0.005*)

Calculate the sag in a hyperpolarizing voltage response.

Parameters

peak_width [window width to get more robust peak estimate in sec (default 0.005)]

Returns

sag [fraction that membrane potential relaxes back to baseline]

estimate_time_constant (*self*)

Calculate the membrane time constant by fitting the voltage response with a single exponential.

Returns

tau [membrane time constant in seconds]

is_spike_feature_affected_by_clipping (*self*, *key*)

pause_metrics (*self*)

Estimate average number of pauses and average fraction of time spent in a pause

Attempts to detect pauses with a variety of conditions and averages results together.

Pauses that are consistently detected contribute more to estimates.

Returns

avg_n_pauses [average number of pauses detected across conditions]

avg_pause_frac [average fraction of interval (between start and end) spent in a pause]

max_reliability [max fraction of times most reliable pause was detected given weights tested]

n_max_rel_pauses [number of pauses detected with *max_reliability*]

process_new_spike_feature (*self*, *feature_name*, *feature_func*, *affected_by_clipping=False*)

Add new spike-level feature calculation function

The function should take this sweep extractor as its argument. Its results can be accessed by calling the method `spike_feature(<feature_name>)`.

process_new_sweep_feature (*self*, *feature_name*, *feature_func*)

Add new sweep-level feature calculation function

The function should take this sweep extractor as its argument. Its results can be accessed by calling the method `sweep_feature(<feature_name>)`.

process_spikes (*self*)

Perform spike-related feature analysis

set_stimulus_amplitude_calculator (*self*, *function*)

spike_feature (*self*, *key*, *include_clipped=False*, *force_exclude_clipped=False*)

Get specified feature for every spike.

Parameters

key [feature name]

include_clipped: return values for every identified spike, even when clipping means they will be incorrect/

Returns

spike_feature_values [ndarray of features for each spike]

spike_feature_keys (*self*)

Get list of every available spike feature.

spikes (*self*)

Get all features for each spike as a list of records.

stimulus_amplitude (*self*)

sweep_feature (*self*, *key*, *allow_missing=False*)

Get sweep-level feature (*key*).

Parameters

key [name of sweep-level feature]

allow_missing [return np.nan if key is missing for sweep (default False)]

Returns

sweep_feature [sweep-level feature value]

sweep_feature_keys (*self*)

Get list of every available sweep-level feature.

voltage_deflection (*self*, *deflect_type=None*)

Measure deflection (min or max, between start and end if specified).

Parameters

deflect_type [measure minimal ('min') or maximal ('max') voltage deflection] If not specified, it will check to see if the current (i) is positive or negative between start and end, then choose 'max' or 'min', respectively If the current is not defined, it will default to 'min'.

Returns

deflect_v [peak]

deflect_index [index of peak deflection]

```
class allensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor (t_set=None,  
                                                                v_set=None,  
                                                                i_set=None,  
                                                                start=None,  
                                                                end=None,  
                                                                filter=10.0,  
                                                                dv_cutoff=20.0,  
                                                                max_interval=0.005,  
                                                                min_height=2.0,  
                                                                min_peak=-30.0,  
                                                                thresh_frac=0.05,  
                                                                base_line_interval=0.1,  
                                                                base_line_detect_thresh=0.3,  
                                                                id_set=None)
```

Bases: object

classmethod from_sweeps (*sweep_list*)

Initialize EphysSweepSetFeatureExtractor object with a list of pre-existing sweep feature extractor objects.

process_spikes (*self*)

Analyze spike features for all sweeps.

spike_feature_averages (*self*, *key*)

Get nparray of average spike-level feature (*key*) for all sweeps

sweep_features (*self*, *key*, *allow_missing=False*)

Get nparray of sweep-level feature (*key*) for all sweeps

Parameters

key [name of sweep-level feature]

allow_missing [return np.nan if key is missing for sweep (default False)]

Returns**sweep_feature** [ndarray of sweep-level feature values]**sweeps** (*self*)

Get list of EphysSweepFeatureExtractor objects.

```
allensdk.ephys.ephys_extractor.cell_extractor_for_nwb(dataset, ramps,
                                                    short_squares, long_squares,
                                                    subthresh_min_amp=-100)
```

Initialize EphysCellFeatureExtractor object from NWB data set

Parameters**dataset** [NwbDataSet]**ramps** [list of sweep numbers of ramp sweeps]**short_squares** [list of sweep numbers of short square sweeps]**long_squares** [list of sweep numbers of long square sweeps]

```
allensdk.ephys.ephys_extractor.extractor_for_nwb_sweeps(dataset, sweep_numbers,
                                                         fixed_start=None,
                                                         fixed_end=None,
                                                         dv_cutoff=20.0,
                                                         thresh_frac=0.05)
```

```
allensdk.ephys.ephys_extractor.fit_fit_slope(ext)
```

Fit the rate and stimulus amplitude to a line and return the slope of the fit.

```
allensdk.ephys.ephys_extractor.input_resistance(ext)
```

Estimate input resistance in MOhms, assuming all sweeps in passed extractor are hyperpolarizing responses.

```
allensdk.ephys.ephys_extractor.membrane_time_constant(ext)
```

Average the membrane time constant values estimated from each sweep in passed extractor.

```
allensdk.ephys.ephys_extractor.reset_long_squares_start(when)
```

allensdk.ephys.ephys_features module**exception** allensdk.ephys.ephys_features.FeatureError

Bases: Exception

Generic Python-exception-derived object raised by feature detection functions.

```
allensdk.ephys.ephys_features.adaptation_index(isis)
```

Calculate adaptation index of *isis*.

```
allensdk.ephys.ephys_features.analyze_trough_details(v, t, spike_indexes,
                                                       peak_indexes, clipped=None,
                                                       end=None, filter=10.0,
                                                       heavy_filter=1.0,
                                                       term_frac=0.01,
                                                       adp_thresh=0.5, tol=0.5,
                                                       flat_interval=0.002,
                                                       adp_max_delta_t=0.005,
                                                       adp_max_delta_v=10.0,
                                                       dvd_t=None)
```

Analyze trough to determine if an ADP exists and whether the reset is a 'detour' or 'direct'

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
spike_indexes [numpy array of spike indexes]
peak_indexes [numpy array of spike peak indexes]
end [end of time window (optional)]
filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (default 1)]
heavy_filter [lower cutoff frequency for 4-pole low-pass Bessel filter in kHz (default 1)]
thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]
adp_thresh: minimum dV/dt in V/s to exceed to be considered to have an ADP (optional, default 1.5)

tol [tolerance for evaluating whether Vm drops appreciably further after end of spike (default 1.0 mV)]
flat_interval: if the trace is flat for this duration, stop looking for an ADP (default 0.002 s)

adp_max_delta_t: max possible ADP delta t (default 0.005 s)
adp_max_delta_v: max possible ADP delta v (default 10 mV)
dvdv [pre-calculated time-derivative of voltage (optional)]

Returns

isi_types [numpy array of isi reset types (direct or detour)]
fast_trough_indexes [numpy array of indexes at the start of the trough (i.e. end of the spike)]
adp_indexes [numpy array of adp indexes (np.nan if there was no ADP in that ISI)]
slow_trough_indexes [numpy array of indexes at the minimum of the slow phase of the trough] (if there wasn't just a fast phase)

`allensdk.ephys.ephys_features.average_rate(t, spikes, start, end)`
Calculate average firing rate during interval between *start* and *end*.

Parameters

t [numpy array of times in seconds]
spikes [numpy array of spike indexes]
start [start of time window for spike detection]
end [end of time window for spike detection]

Returns

avg_rate [average firing rate in spikes/sec]

`allensdk.ephys.ephys_features.average_voltage(v, t, start=None, end=None)`
Calculate average voltage between start and end.

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
start [start of time window for spike detection (optional, default None)]

end [end of time window for spike detection (optional, default None)]

Returns

v_avg [average voltage]

`allensdk.ephys.ephys_features.calculate_dvdt(v, t, filter=None)`

Low-pass filters (if requested) and differentiates voltage by time.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default None)]

Returns

dvdt [numpy array of time-derivative of voltage (V/s = mV/ms)]

`allensdk.ephys.ephys_features.check_thresholds_and_peaks(v, t, spike_indexes, peak_indexes, upstroke_indexes, end=None, max_interval=0.005, thresh_frac=0.05, filter=10.0, dvdt=None, tol=1.0)`

Validate thresholds and peaks for set of spikes

Check that peaks and thresholds for consecutive spikes do not overlap Spikes with overlapping thresholds and peaks will be merged.

Check that peaks and thresholds for a given spike are not too far apart.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]

peak_indexes [numpy array of indexes of spike peaks]

upstroke_indexes [numpy array of indexes of spike upstrokes]

max_interval [maximum allowed time between start of spike and time of peak in sec (default 0.005)]

thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

tol [tolerance for returning to threshold in mV (optional, default 1)]

Returns

spike_indexes [numpy array of modified spike indexes]

peak_indexes [numpy array of modified spike peak indexes]

upstroke_indexes [numpy array of modified spike upstroke indexes]

clipped [numpy array of clipped status of spikes]

`allensdk.ephys.ephys_features.detect_bursts` (*isis, isi_types, fast_tr_v, fast_tr_t, slow_tr_v, slow_tr_t, thr_v, tol=0.5, pause_cost=1.0*)

Detect bursts in spike train.

Parameters

isis [numpy array of n interspike intervals]

isi_types [numpy array of n interspike interval types]

fast_tr_v [numpy array of fast trough voltages for the n + 1 spikes of the train]

fast_tr_t [numpy array of fast trough times for the n + 1 spikes of the train]

slow_tr_v [numpy array of slow trough voltages for the n + 1 spikes of the train]

slow_tr_t [numpy array of slow trough times for the n + 1 spikes of the train]

thr_v [numpy array of threshold voltages for the n + 1 spikes of the train]

tol [tolerance for the difference in slow trough voltages and thresholds (default 0.5 mV)] Used to identify “delay” interspike intervals that occur within a burst

Returns

bursts [list of bursts] Each item in list is a tuple of the form (burst_index, start, end) where *burst_index* is a comparison index between the highest instantaneous rate within the burst vs the highest instantaneous rate outside the burst. *start* is the index of the first ISI of the burst, and *end* is the ISI index immediately following the burst.

`allensdk.ephys.ephys_features.detect_pauses` (*isis, isi_types, cost_weight=1.0*)

Determine which ISIs are “pauses” in ongoing firing.

Pauses are unusually long ISIs with a “detour reset” among “direct resets”.

Parameters

isis [numpy array of interspike intervals]

isi_types [numpy array of interspike interval types (‘direct’ or ‘detour’)]

cost_weight [weight for cost function for calling an ISI a pause] Higher cost weights lead to fewer ISIs identified as pauses. The cost function also depends on the difference between the duration of the “pause” ISIs and the average duration and standard deviation of “non-pause” ISIs.

Returns

pauses [numpy array of indices corresponding to pauses in *isis*]

`allensdk.ephys.ephys_features.detect_putative_spikes` (*v, t, start=None, end=None, filter=10.0, dv_cutoff=20.0*)

Perform initial detection of spikes and return their indexes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

start [start of time window for spike detection (optional)]

end [end of time window for spike detection (optional)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dv_cutoff [minimum dV/dt to qualify as a spike in V/s (optional, default 20)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

putative_spikes [numpy array of preliminary spike indexes]

```
allensdk.ephys.ephys_features.estimate_adjusted_detection_parameters(v_set,
                                                                    t_set,
                                                                    inter-
                                                                    val_start,
                                                                    inter-
                                                                    val_end,
                                                                    fil-
                                                                    ter=10)
```

Estimate adjusted values for spike detection by analyzing a period when the voltage changes quickly but passively (due to strong current stimulation), which can result in spurious spike detection results.

Parameters

v_set [list of numpy arrays of voltage time series in mV]

t_set [list of numpy arrays of times in seconds]

interval_start [start of analysis interval (sec)]

interval_end [end of analysis interval (sec)]

Returns

new_dv_cutoff [adjusted dv/dt cutoff (V/s)]

new_thresh_frac [adjusted fraction of avg upstroke to find threshold]

```
allensdk.ephys.ephys_features.filter_putative_spikes(v, t, spike_indexes,
                                                    peak_indexes, min_height=2.0,
                                                    min_peak=-30.0, filter=10.0,
                                                    dvdt=None)
```

Filter out events that are unlikely to be spikes based on:

- Voltage failing to go down between peak and the next spike's threshold
- Height (threshold to peak)
- Absolute peak level

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of preliminary spike indexes]

peak_indexes [numpy array of indexes of spike peaks]

min_height [minimum acceptable height from threshold to peak in mV (optional, default 2)]

min_peak [minimum acceptable absolute peak level in mV (optional, default -30)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

spike_indexes [numpy array of threshold indexes]

peak_indexes [numpy array of peak indexes]

```
allensdk.ephys.ephys_features.find_downstroke_indexes(v, t, peak_indexes,
                                                         trough_indexes,
                                                         clipped=None, filter=10.0,
                                                         dvdt=None)
```

Find indexes of minimum voltage (troughs) between spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

peak_indexes [numpy array of spike peak indexes]

trough_indexes [numpy array of threshold indexes]

clipped: boolean array - False if spike not clipped by edge of window

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

downstroke_indexes [numpy array of downstroke indexes]

```
allensdk.ephys.ephys_features.find_peak_indexes(v, t, spike_indexes, end=None)
```

Find indexes of spike peaks.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of preliminary spike indexes]

end [end of time window for spike detection (optional)]

```
allensdk.ephys.ephys_features.find_time_index(t, t_0)
```

Find the index value of a given time (t_0) in a time series (t).

```
allensdk.ephys.ephys_features.find_trough_indexes(v, t, spike_indexes, peak_indexes,
                                                         clipped=None, end=None)
```

Find indexes of minimum voltage (trough) between spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]

peak_indexes [numpy array of spike peak indexes]

end [end of time window (optional)]

Returns

trough_indexes [numpy array of threshold indexes]

```
allensdk.ephys.ephys_features.find_upstroke_indexes(v, t, spike_indexes, peak_indexes,
                                                         filter=10.0, dvdt=None)
```

Find indexes of maximum upstroke of spike.

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
spike_indexes [numpy array of preliminary spike indexes]
peak_indexes [numpy array of indexes of spike peaks]
filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]
dvdv [pre-calculated time-derivative of voltage (optional)]

Returns

upstroke_indexes [numpy array of upstroke indexes]

`allensdk.ephys.ephys_features.find_widths(v, t, spike_indexes, peak_indexes, trough_indexes, clipped=None)`

Find widths at half-height for spikes.

Widths are only returned when heights are defined

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
spike_indexes [numpy array of spike indexes]
peak_indexes [numpy array of spike peak indexes]
trough_indexes [numpy array of trough indexes]

Returns

widths [numpy array of spike widths in sec]

`allensdk.ephys.ephys_features.fit_membrane_time_constant(v, t, start, end, min_rsme=0.0001)`

Fit an exponential to estimate membrane time constant between start and end

Parameters

v [numpy array of voltages in mV]
t [numpy array of times in seconds]
start [start of time window for exponential fit]
end [end of time window for exponential fit]
min_rsme: minimal acceptable root mean square error (default 1e-4)

Returns

a, inv_tau, y0 [Coefficients of equation $y_0 + a * \exp(-\text{inv_tau} * x)$]
returns np.nan for values if fit fails

`allensdk.ephys.ephys_features.fit_prespike_time_constant(v, t, start, spike_time, dv_limit=-0.001, tau_limit=0.3)`

Finds the dominant time constant of the pre-spike rise in voltage

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

start [start of voltage rise (seconds)]

spike_time [time of first spike (seconds)]

dv_limit [dV/dt cutoff (default -0.001)] Shortens fit window if rate of voltage drop exceeds this limit

tau_limit [upper bound for slow time constant (seconds, default 0.3)] If the slower time constant of a double-exponential fit is twice that of the faster and exceeds this limit, the faster one will be considered the dominant one

Returns

tau [dominant time constant (seconds)]

`allensdk.ephys.ephys_features.get_isis(t, spikes)`
Find interspike intervals in sec between spikes (as indexes).

`allensdk.ephys.ephys_features.has_fixed_dt(t)`
Check that all time intervals are identical.

`allensdk.ephys.ephys_features.latency(t, spikes, start)`
Calculate time to the first spike.

`allensdk.ephys.ephys_features.norm_diff(a)`
Calculate average of $(a[i] - a[i+1]) / (a[i] + a[i+1])$.

`allensdk.ephys.ephys_features.norm_sq_diff(a)`
Calculate average of $(a[i] - a[i+1])^2 / (a[i] + a[i+1])^2$.

`allensdk.ephys.ephys_features.refine_threshold_indexes(v, t, upstroke_indexes,`
thresh_frac=0.05, filter=10.0, dvdt=None

Refine threshold detection of previously-found spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

upstroke_indexes [numpy array of indexes of spike upstrokes (for threshold target calculation)]

thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

threshold_indexes [numpy array of threshold indexes]

`allensdk.ephys.extract_cell_features` module

`allensdk.ephys.extract_cell_features.extract_cell_features(data_set,`
ramp_sweep_numbers,
short_square_sweep_numbers,
long_square_sweep_numbers,
sub-
thresh_min_amp=None


```
allensdk.ephys.extract_cell_features.extract_sweep_features (data_set,
                                                             sweeps_by_type)
allensdk.ephys.extract_cell_features.get_ramp_stim_characteristics (i, t)
    Identify the start time and start index of a ramp sweep.
allensdk.ephys.extract_cell_features.get_square_stim_characteristics (i, t,
                                                                    no_test_pulse=False)
    Identify the start time, duration, amplitude, start index, and end index of a square stimulus. This assumes that
    there is a test pulse followed by the stimulus square.
allensdk.ephys.extract_cell_features.get_stim_characteristics (i, t,
                                                            no_test_pulse=False)
    Identify the start time, duration, amplitude, start index, and end index of a general stimulus. This assumes that
    there is a test pulse followed by the stimulus square.
allensdk.ephys.extract_cell_features.mean_features_spike_zero (sweeps)
    Compute mean feature values for the first spike in list of extractors
```

allensdk.ephys.feature_extractor module

```
class allensdk.ephys.feature_extractor.EphysFeatureExtractor
    Bases: object
    adaptation_index (self, spikes, stim_end)
    calculate_trough (self, spike, v, curr, t, next_idx)
    isicv (self, spikes)
    process_instance (self, name, v, curr, t, onset, dur, stim_name)
    push_summary (self, new_summary)
    score_feature_set (self, set_num)
    summarize (self, summary)
class allensdk.ephys.feature_extractor.EphysFeatures (name)
    Bases: object
    clone (self, param_dict)
    print_out (self)
```

Module contents

6.1.6 allensdk.internal package

Subpackages

allensdk.internal.api package

Subpackages

allensdk.internal.api.queries package

Submodules

`allensdk.internal.api.queries.biophysical_module_api` module

class `allensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi` (*base_uri=None*)

Bases: `allensdk.api.queries.rma_template.RmaTemplate`

get_neuronal_model_runs (*self, neuronal_model_run_ids=None*)

List Neuronal Model Rusn available through LIMS with associated info needed to run in NEURON.

Parameters

neuronal_model_run_ids [integer or list of integers, optional] only select specific neuronal_model_runs.

Returns

dict [neuronal model run metadata]

get_neuronal_models (*self, neuronal_model_ids=None*)

List Neuronal Models available through LIMS with associated info needed to run in NEURON.

Parameters

neuronal_model_ids [integer or list of integers, optional] only select specific neuronal_models.

Returns

dict [neuronal model metadata]

rma_templates = {'biophysical_lims_queries': [{'name': 'neuronal_model_runs_by_ids',

`allensdk.internal.api.queries.biophysical_module_reader` module

class `allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader`

Bases: `object`

MOD_FILE_TYPE_ID = 292178729

MORPHOLOGY_TYPE_ID = 303941301

STIMULUS_CONTENT_TYPE = None

fit_parameters_file_entries (*self*)

read the fit_parameter file path from the lims result corresponding to the stimulus file :return: well_known_file entries :rtype: array of dicts

fit_parameters_path (*self*)

Get the path to the fit parameters file from the lims result. :return: path to file :rtype: string

lims_working_directory (*self*)

While this is the same directory as the neuronal_model_run directory, it can be mocked out for testing if the other directory is read only.

mod_file_entries (*self*)

read the NERUON .mod file entries from the lims result corresponding to the NeuronModel :return: well known file entries :rtype: array of dicts

mod_file_paths (*self*)

Get the paths to the mod files from the lims result. :return: paths to mod files :rtype: array of strings

model_type (*self*)
 TODO: comment

morphology_file_entries (*self*)
 read the well known file paths from the lims result corresponding to the morphology

Returns

array of dicts: well known file entries

morphology_path (*self*)
 Get the path to the morphology file from the lims result. :return: path to morphology file :rtype: string

neuronal_model_run_dir (*self*)
 read the directory path where output goes from the lims optimization config json

Returns

string: directory path

read_json (*self*, *path*)

read_json_string (*self*, *json_string*)

read_lims_file (*self*, *lims_path*)

read_lims_message (*self*, *message*, *lims_path*)

set_workflow_state (*self*, *state*)

stimulus_file_entries (*self*)
 read the well known file path from the lims result corresponding to the stimulus file :return:
 well_known_file entries :rtype: array of dicts

stimulus_path (*self*)
 Get the path to the stimulus file from the lims result. :return: path to stimulus file :rtype: string

sweep_entries (*self*)
 read the sweep entries from the lims result corresponding to the stimulus :return: stimulus sweep entries
 :rtype: array of dicts

sweep_numbers (*self*)
 Get the stimulus sweep numbers from the lims result :return: list of sweep numbers :rtype: array of ints

sweep_numbers_by_type (*self*)

to_manifest (*self*, *manifest_path*=None)

update_well_known_file (*self*, *path*, *well_known_file_type_id*=None)

write_file (*self*, *path*)

allensdk.internal.api.queries.grid_data_api_prerelease module

class allensdk.internal.api.queries.grid_data_api_prerelease.**GridDataApiPrerelease** (*storage_dir*,
res-
o-
lu-
tion=None,
base_uri=None)

Bases: *allensdk.api.queries.grid_data_api.GridDataApi*

Client for retrieving prereleased mouse connectivity data from lims.

Parameters

base_uri [string, optional] Does not affect pulling from lims.

file_name [string, optional] File name to save/read storage_directories dict. Passed to GridDataApiPrerelease constructor.

GRID_DATA_DIRECTORY = 'grid'

download_projection_grid_data (*self, path, experiment_id, file_name*)

Copy data from path to file_name.

Parameters

path [string] path to file in shared directory (copy source)

experiment_id [int] image series id.

file_name [string] path to file destination (copy target)

classmethod from_file_name (*file_name, cache=True, **kwargs*)

Alternative constructor using cache path file_name.

Parameters

file_name [string] Path where storage_directories will be saved.

****kwargs** Keyword arguments to be supplied to `__init__`

Returns

cls [instance of GridDataApiPrerelease]

allensdk.internal.api.queries.mouse_connectivity_api_prerelease module

class allensdk.internal.api.queries.mouse_connectivity_api_prerelease.**MouseConnectivityApi**

Bases: *allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi*

Client for retrieving prereleased mouse connectivity data from lims.

Parameters

base_uri [string, optional] Does not affect pulling from lims.

file_name [string, optional] File name to save/read storage_directories dict. Passed to GridDataApiPrerelease constructor.

download_data_mask (*self, path, experiment_id, resolution*)

download_injection_density (*self, path, experiment_id, resolution*)

download_injection_fraction (*self, path, experiment_id, resolution*)

download_projection_density (*self, path, experiment_id, resolution*)

get_experiments (*self*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas.

Parameters

structure_ids [integer or list, optional] injection structure

Returns

url [string] The constructed URL

```
get_structure_unionizes(self)
```

allensdk.internal.api.queries.optimize_config_reader module

```
class allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader
```

```
Bases: object
```

```
MOD_FILE_TYPE_ID = 292178729
```

```
MORPHOLOGY_TYPE_ID = 303941301
```

```
NEURONAL_MODEL_PARAMETERS = 329230374
```

```
STIMULUS_CONTENT_TYPE = None
```

```
build_manifest(self, manifest_path=None)
```

```
lims_working_directory(self)
```

While this is the same directory as the optimize directory, it can be mocked out for testing if the optimize directory is write only.

```
mod_file_entries(self)
```

read the NERUON .mod file entries from the lims result corresponding to the NeuronModel :return: well known file entries :rtype: array of dicts

```
mod_file_paths(self)
```

Get the paths to the mod files from the lims result. :return: paths to mod files :rtype: array of strings

```
morphology_file_entries(self)
```

read the well known file paths from the lims result corresponding to the morphology

Returns

array of dicts: well known file entries

```
morphology_path(self)
```

Get the path to the morphology file from the lims result. :return: path to morphology file :rtype: string

```
neuronal_model_optimize_dir(self)
```

read the directory path where output goes from the lims optimization config json

Returns

string: directory path

```
output_directory(self)
```

```
read_json(self, path)
```

```
read_json_string(self, json_string)
```

```
read_lims_file(self, lims_path)
```

```
read_lims_message(self, message, lims_path)
```

```
stimulus_file_entries(self)
```

read the well known file path from the lims result corresponding to the stimulus file :return: well-known_file entries :rtype: array of dicts

```
stimulus_path(self)
```

Get the path to the stimulus file from the lims result. :return: path to stimulus file :rtype: string

sweep_entries (*self*)
read the sweep entries from the lims result corresponding to the stimulus :return: stimulus sweep entries
:rtype: array of dicts

sweep_numbers (*self*)
Get the stimulus sweep numbers from the lims result :return: list of sweep numbers :rtype: array of ints

to_manifest (*self*, *manifest_path=None*)

update_well_known_file (*self*, *path*, *well_known_file_type_id=None*)

write_file (*self*, *path*)

allensdk.internal.api.queries.pre_release module

Module contents

Submodules

allensdk.internal.api.api_prerelease module

class allensdk.internal.api.api_prerelease.**ApiPrerelease** (*api_base_url_string=None*)
Bases: *allensdk.api.api.Api*

Extends allensdk.api.api to copy files ‘locally’ from shared storage.

retrieve_file_from_storage (*self*, *storage_path*, *save_file_path*)
Copy data from path to file_name.

Parameters

storage_path [string] path to file in shared directory (copy source)

save_file_name [string] path to file destination (copy target)

allensdk.internal.api.lims_api module

class allensdk.internal.api.lims_api.**LimsApi** (*lims_credentials:* *Optional[allensdk.core.authentication.DbCredentials]* *= None*)
Bases: object

get_behavior_tracking_video_filepath_df (*self*)

get_experiment_id (*self*)

get_eye_tracking_video_filepath_df (*self*)

allensdk.internal.api.mtrain_api module

Module contents

exception allensdk.internal.api.**OneOrMoreResultExpectedError**
Bases: RuntimeError

```

class allensdk.internal.api.PostgresQueryMixin (*, dbname, user, host, password, port)
    Bases: object

    fetchall (self, query, strict=True)
    fetchone (self, query, strict=True)
    get_connection (self)
    get_cursor (self)
    select (self, query)
    select_one (self, query)

allensdk.internal.api.db_connection_creator (credentials: Union[allensdk.core.authentication.DbCredentials,
    NoneType] = None, fallback_credentials:
    Union[dict, NoneType] = None) → al-
    lensdk.internal.api.PostgresQueryMixin

```

Create a db connection using credentials. If credentials are not provided then use fallback credentials (which attempt to read from shell environment variables).

Note: Must provide one of either ‘credentials’ or ‘fallback_credentials’. If both are provided, ‘credentials’ will take precedence.

Parameters

credentials [Optional[DbCredentials], optional] User specified credentials, by default None

fallback_credentials [dict] Fallback credentials to use for creating the DB connection in the case that no ‘credentials’ are provided, by default None.

Fallback credentials will attempt to get db connection info from shell environment variables.

Some examples of environment variables that fallback credentials will try to read from can be found in `allensdk.core.auth_config`.

Returns

PostgresQueryMixin A DB connection instance which can execute queries to the DB specified by credentials or fallback_credentials.

Raises

RuntimeError If neither ‘credentials’ nor ‘fallback_credentials’ were provided.

```
allensdk.internal.api.psycpg2_select (query, database, host, port, username, password)
```

allensdk.internal.brain_observatory package

Subpackages

allensdk.internal.brain_observatory.resources package

Module contents

Submodules

allensdk.internal.brain_observatory.annotated_region_metrics module

Module for calculating annotated region metrics from ISI data

```
allensdk.internal.brain_observatory.annotated_region_metrics.create_region_mask(image_shape,  
                                                                                  x,  
                                                                                  y,  
                                                                                  width,  
                                                                                  height,  
                                                                                  mask)
```

Create mask for region on retinotopic map

Parameters

image_shape [tuple] (height, width) of retinotopic map
x [int] x offset of region mask within retinotopic map
y [int] y offset of region mask within retinotopic map
width [int] width of region mask
height [int] height of region mask
mask [list] region mask as a list of lists

Returns

numpy.ndarray Region mask

```
allensdk.internal.brain_observatory.annotated_region_metrics.eccentricity(az,  
                                                                              alt,  
                                                                              az_center,  
                                                                              alt_center)
```

Compute eccentricity

Parameters

az [numpy.ndarray] Azimuth retinotopic map
alt [numpy.ndarray] Altitude retinotopic map
az_center [float] Azimuth value to use as center of eccentricity map
alt_center [float] Altitude value to use as center of eccentricity map

Returns

numpy.ndarray Eccentricity map

```
allensdk.internal.brain_observatory.annotated_region_metrics.get_metrics(altitude_phase,  
                                                                              az-  
                                                                              imuth_phase,  
                                                                              x=None,  
                                                                              y=None,  
                                                                              width=None,  
                                                                              height=None,  
                                                                              mask=None,  
                                                                              al-  
                                                                              ti-  
                                                                              tude_scale=0.322,  
                                                                              az-  
                                                                              imuth_scale=0.383)
```

Calculate annotated region metrics

`allensdk.internal.brain_observatory.annotated_region_metrics.retinotopy_metric` (*mask*,
isi_map)

Compute retinotopic metrics for a responding area

Parameters

mask [numpy.ndarray] Mask representing the area over which to calculate metrics

isi_map [numpy.ndarray] Retinotopic map

Returns

(float, float, float, float) tuple min, max, range, bias of retinotopic map over masked region

`allensdk.internal.brain_observatory.demix_report` module

`allensdk.internal.brain_observatory.demix_report.background_trace` (*trace*,
save_dir,
data_set=None)

`allensdk.internal.brain_observatory.demix_report.compute_correlations` (*dm*,
movie_path,
movie_dataset)

`allensdk.internal.brain_observatory.demix_report.compute_correlations_without_masks` (*dm*)

`allensdk.internal.brain_observatory.demix_report.compute_non_overlap_masks` (*dm*)

`allensdk.internal.brain_observatory.demix_report.compute_non_overlap_traces` (*dm*,
movie_path,
movie_dataset)

`allensdk.internal.brain_observatory.demix_report.correlation_report` (*dm*,
save_dir,
with-
out_masks=True)

parameters: *dm*: [DeMix object] *without_masks*: boolean

`allensdk.internal.brain_observatory.demix_report.plot_masks` (*dm*, *save_dir*,
movie_file,
movie_dataset,
window=150,
add_background=True)

`allensdk.internal.brain_observatory.demixer` module

`allensdk.internal.brain_observatory.demixer.demix_time_dep_masks` (*raw_traces*,
stack,
masks)

Parameters

- **raw_traces** – extracted traces
- **stack** – movie (same length as traces)
- **masks** – binary roi masks

Returns

demixed traces

`allensdk.internal.brain_observatory.demixer.find_negative_baselines` (*trace*)

```
allensdk.internal.brain_observatory.demixer.find_negative_transients_threshold(trace,
                                                                              win-
                                                                              dow=500,
                                                                              length=10,
                                                                              std_devs=3)

allensdk.internal.brain_observatory.demixer.find_zero_baselines(traces)

allensdk.internal.brain_observatory.demixer.identify_valid_masks(mask_array)

allensdk.internal.brain_observatory.demixer.plot_negative_baselines(raw_traces,
                                                                    demix_traces,
                                                                    mask_array,
                                                                    roi_ids_mask,
                                                                    plot_dir,
                                                                    ext='png')

allensdk.internal.brain_observatory.demixer.plot_negative_transients(raw_traces,
                                                                    demix_traces,
                                                                    valid_roi,
                                                                    mask_array,
                                                                    roi_ids_mask,
                                                                    plot_dir,
                                                                    ext='png')

allensdk.internal.brain_observatory.demixer.plot_overlap_masks_lengthOne(roi_ind,
                                                                            masks,
                                                                            save-
                                                                            file=None,
                                                                            weighted=False)

allensdk.internal.brain_observatory.demixer.plot_traces(raw_trace,    demix_trace,
                                                         roi_id, roi_ind, save_file)

allensdk.internal.brain_observatory.demixer.plot_transients(roi_ind,    t_trans,
                                                             masks,    traces,
                                                             demix_traces, save-
                                                             file)

allensdk.internal.brain_observatory.demixer.rolling_window(trace, window=500)
```

Parameters

- **trace** –
- **window** –

Returns

allensdk.internal.brain_observatory.eye_calibration module

```
class allensdk.internal.brain_observatory.eye_calibration.EyeCalibration (monitor_position=array([
    8.62,
    3.16]),
    mon-
    i-
    tor_rotations=array([0.,
    0.,
    0.]),
    led_position=array([25.8,
    -
    6.12,
    3.21]),
    cam-
    era_position=array([13.,
    0.,
    0.]),
    cam-
    era_rotations=array([0.,
    0.,
    0.22863813]),
    eye_radius=0.1682,
    cm_per_pixel=0.0010199
```

Bases: `object`

Class for performing eye-tracking calibration.

Provides methods for estimating the position of the pupil in 3D space and projecting the gaze onto the monitor in both 3D space and monitor space given the experimental geometry.

Parameters

monitor_position [numpy.ndarray] [x,y,z] position of monitor in cm.

monitor_rotations [numpy.ndarray] [x,y,z] rotations of monitor in radians.

led_position [numpy.ndarray] [x,y,z] position of LED in cm.

camera_position [numpy.ndarray] [x,y,z] position of camera in cm.

camera_rotations [numpy.ndarray] [x,y,z] rotations for camera in radians. X and Y must be 0.

eye_radius [float] Radius of the eye in cm.

cm_per_pixel [float] Pixel size of eye-tracking camera.

compute_area (*self*, *pupil_parameters*)

Compute the area of the pupil.

Assume the pupil is a circle, and that as it moves off-axis with the camera the observed ellipse major axis remains the diameter of the circle.

Parameters

pupil_parameters [numpy.ndarray] [nx5] array of pupil parameters.

Returns

numpy.ndarray [nx1] array of pupil areas in estimated pixels.

static cr_position_in_mouse_eye_coordinates (*led_position, eye_radius*)

Determine the 3D position of the corneal reflection.

The eye is modeled as a spherical mirror, so the reflection appears to be half the radius of the eye from the origin along the eye-LED axis.

Parameters

led_position [numpy.ndarray] [x,y,z] position of the LED in eye coordinates.

eye_radius [float] Radius of the eye in centimeters.

Returns

numpy.ndarray [x,y,z] location of the corneal reflection in eye coordinates.

pupil_position_in_mouse_eye_coordinates (*self, pupil_parameters, cr_parameters*)

Compute the 3D pupil position in mouse eye coordinates.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_parameters [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimates in eye coordinates.

pupil_position_on_monitor_in_cm (*self, pupil_parameters, cr_parameters*)

Compute the pupil position on the monitor in cm.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_parameters [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimates in eye coordinates.

pupil_position_on_monitor_in_degrees (*self, pupil_parameters, cr_parameters*)

Get pupil position on monitor measured in visual degrees.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_parameters [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimate in visual degrees.

allensdk.internal.brain_observatory.eye_calibration.base_object_to_eye_rotation_matrix (*object*)

Rotation matrix to rotate base object frame to eye coordinates.

By convention, any other object's coordinate frame before rotations is set with positive Z pointing from the object's position back to the origin of the eye coordinate system, with X parallel to the eye X-Y plane.

Parameters

object_position [np.ndarray] [x, y, z] position of object in eye coordinates.

Returns

numpy.ndarray [3x3] rotation matrix.

```
allensdk.internal.brain_observatory.eye_calibration.object_norm_eye_coordinates(object_position,
                                                                              x_rotation,
                                                                              y_rotation,
                                                                              z_rotation)
```

Get the normal vector for the object plane in eye coordinates.

Parameters

object_position [numpy.ndarray] [x, y, z] location of the object in eye coordinates.

x_rotation [float] Rotation about the x-axis in radians.

y_rotation [float] Rotation about the y-axis in radians.

z_rotation [float] Rotation about the z-axis in radians.

Returns

numpy.ndarray Endpoint of the object plane vector in eye coordinates.

```
allensdk.internal.brain_observatory.eye_calibration.object_rotation_matrix(x_rotation,
                                                                              y_rotation,
                                                                              z_rotation)
```

Rotation matrix in object coordinate frame.

The rotation matrix for rotating the object coordinate frame from the initial position. This is done by rotating around x, then around y', then around z".

Parameters

x_rotation [float] Rotation about x axis in radians.

y_rotation [float] Rotation about y axis in radians.

z_rotation [float] Rotation about z axis in radians.

Returns

numpy.ndarray [3x3] rotation matrix.

```
allensdk.internal.brain_observatory.eye_calibration.project_to_plane(plane_normal,
                                                                      plane_point,
                                                                      points)
```

Project from the origin through points onto a plane.

Parameters

plane_normal [numpy.ndarray] [x, y, z] normal unit vector to the plane.

plane_point [numpy.ndarray] [x, y, z] point on the plane.

points [numpy.ndarray] [nx3] points in space through which to project.

Returns

numpy.ndarray [nx3] points projected on the plane.

allensdk.internal.brain_observatory.fit_ellipse module

```
class allensdk.internal.brain_observatory.fit_ellipse.FitEllipse(min_points,  
                                                                max_iter,  
                                                                threshold,  
                                                                num_close)
```

Bases: object

```
choose_inliers (self, candidate_points)
```

```
fit_ellipse (self, inlier_points)
```

```
outlier_cost (self, outlier_points, params)
```

```
ransac_fit (self, candidate_points)
```

```
allensdk.internal.brain_observatory.fit_ellipse.ellipse_angle_of_rotation (a)
```

```
allensdk.internal.brain_observatory.fit_ellipse.ellipse_angle_of_rotation2 (a)
```

```
allensdk.internal.brain_observatory.fit_ellipse.ellipse_axis_length (a)
```

```
allensdk.internal.brain_observatory.fit_ellipse.ellipse_center (a)
```

```
allensdk.internal.brain_observatory.fit_ellipse.fit_ellipse (candidate_points)
```

```
allensdk.internal.brain_observatory.fit_ellipse.rotate_vector (y, x, theta)
```

```
allensdk.internal.brain_observatory.fit_ellipse.test_fit ()
```

allensdk.internal.brain_observatory.frame_stream module

```
class allensdk.internal.brain_observatory.frame_stream.CvInputStream(movie_path,  
                                                                num_frames=None,  
                                                                block_size=1,  
                                                                cache_frames=False)
```

Bases: object

```
close (self)
```

```
open (self)
```

```
class allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream(movie_path,  
                                                                frame_shape,  
                                                                ffm-  
                                                                peg_bin='ffmpeg',  
                                                                num_frames=None,  
                                                                block_size=1,  
                                                                cache_frames=False,  
                                                                pro-  
                                                                cess_frame_cb=None)
```

Bases: `allensdk.internal.brain_observatory.frame_stream.FrameInputStream`

```
close (self)
```

```
create_images (self, output_directory, image_type)
```

```
open (self)
```

```

class allensdk.internal.brain_observatory.frame_stream.FfmpegOutputStream (frame_shape,
                                                                    ffm-
                                                                    peg_bin='ffmpeg',
                                                                    block_size=1)

    Bases: allensdk.internal.brain_observatory.frame_stream.FrameOutputStream

    close (self)

    open (self, movie_path)

class allensdk.internal.brain_observatory.frame_stream.FrameInputStream (movie_path,
                                                                    num_frames=None,
                                                                    block_size=1,
                                                                    cache_frames=False,
                                                                    process_frame_cb=None)

    Bases: object

    close (self)

    create_images (self, output_directory, image_type)

    open (self)

class allensdk.internal.brain_observatory.frame_stream.FrameOutputStream (block_size=1)

    Bases: object

    close (self)

    open (self, movie_path)

    write (self, frame)

class allensdk.internal.brain_observatory.frame_stream.ImageOutputStream (block_size=1)

    Bases: allensdk.internal.brain_observatory.frame_stream.FrameOutputStream

```

allensdk.internal.brain_observatory.itracker module

allensdk.internal.brain_observatory.itracker_utils module

```

allensdk.internal.brain_observatory.itracker_utils.default_ray (n)

allensdk.internal.brain_observatory.itracker_utils.eccentricity (a1, a2)

allensdk.internal.brain_observatory.itracker_utils.filter_bad_params (params,
                                                                    frame_width,
                                                                    frame_height)

    Replace positions outside image with nan

allensdk.internal.brain_observatory.itracker_utils.generate_rays (image_array,
                                                                    seed_pixel)

allensdk.internal.brain_observatory.itracker_utils.initial_cr_point (image_array,
                                                                    bbox=None)

    bbox is a tuple of (xmin, xmax, ymin, ymax)

allensdk.internal.brain_observatory.itracker_utils.initial_pupil_point (image_array,
                                                                    bbox=None)

    bbox is a tuple of (xmin, xmax, ymin, ymax)

allensdk.internal.brain_observatory.itracker_utils.medfilt_custom (x,      kernel_size=3)

    This median filter returns 'nan' whenever any value in the kernel width is 'nan' and the median otherwise

```

```
allensdk.internal.brain_observatory.itracker_utils.median_absolute_deviation(a,
```

```
consistency_constant=1.4
```

Calculate the median absolute deviation of a univariate dataset.

Parameters

a [numpy.ndarray] Sample data.

consistency_constant [float] Constant to make the MAD a consistent estimator of the population standard deviation (1.4826 for a normal distribution).

Returns

float Median absolute deviation of the data.

```
allensdk.internal.brain_observatory.itracker_utils.post_process_cr(cr_params)
```

This will replace questionable values of the CR x and y position with 'nan'

- 1) threshold ellipse area by 99th percentile area distribution
- 2) median filter using custom median filter
- 3) remove deviations from discontinuous jumps

The 'nan' values likely represent obscured CRs, secondary reflections, merges with the secondary reflection, or visual distortions due to the whisker or deformations of the eye

```
allensdk.internal.brain_observatory.itracker_utils.post_process_pupil(pupil_params)
```

Filter pupil parameters to replace outliers with nan

Parameters

pupil_params [numpy.ndarray] (Nx5) array of pupil parameters [x, y, angle, axis1, axis2].

Returns

numpy.ndarray Pupil parameters with outliers replaced with nan

```
allensdk.internal.brain_observatory.itracker_utils.rotate_ray(ray, theta)
```

```
allensdk.internal.brain_observatory.itracker_utils.sobel_grad(image_array)
```

allensdk.internal.brain_observatory.mask_set module

```
class allensdk.internal.brain_observatory.mask_set.MaskSet(masks)
```

Bases: object

close (self, mask_idx, max_dist)

close_sets (self, set_size, max_dist)

count

detect_duplicates (self, overlap_threshold)

detect_unions (self, set_size=2, max_dist=10, threshold=0.7)

distance (self, mask_idx)

intersection (self, mask_idx)

intersection_size (self, mask_idx)

mask (self, mask_idx)

mask_is_union_of_set (*self*, *mask_idx*, *set_idx*s, *threshold*)

overlap_fraction (*self*, *idx0*, *idx1*)

size (*self*, *mask_idx*)

union (*self*, *mask_idx*s)

union_size (*self*, *mask_idx*s)

`allensdk.internal.brain_observatory.mask_set.bb_dist` (*bbs*)

`allensdk.internal.brain_observatory.mask_set.make_bbs` (*masks*)

allensdk.internal.brain_observatory.ophys_session_decomposition module

`allensdk.internal.brain_observatory.ophys_session_decomposition.export_frame_to_hdf5` (*raw_filename*, *data_hdf5_filename*, *auxiliary_hdf5_filename*, *channel_descriptions*, *frame_number*, *compression*, *use_gzip*, *use_compression*, *use_presession_ops*)

Export a frame from raw to hdf5.

Data with the channel_description *data* is stored in the *data_hdf5_filename*, while any other data is stored in the *auxiliary_hdf5_filename*

`allensdk.internal.brain_observatory.ophys_session_decomposition.load_frame` (*raw_filename*, *json_meta*, *use_memmap=False*)

Load a frame of a multi-frame raw file.

`allensdk.internal.brain_observatory.ophys_session_decomposition.open_view_on_binary` (*file_like*, *dtype=<dtype>*, *numpy.dtype*, *mode='r'*, *offset=0*, *shape=None*, *order='C'*, *strides=None*)

Open a view into a memory-mapped binary file.

Parameters

file_like [{string, file object}] File to open.

dtype [numpy.dtype] Numpy dtype to open the memory-mapped array as.

mode [string] Mode to open the file in.

offset [integer] Offset (in bytes) into the file at which to start the memory map.

shape [(tuple, list)] Shape of the array.

order [{"C", "F"}] C or Fortran ordering.

strides [(tuple, list)] Strides along each axis for reading the array.

Returns

numpy.memmap Strided view into memory-mapped array.

`allensdk.internal.brain_observatory.ophys_session_decomposition.read_strided(filename, dtype, off-set, shape, strides)`

Load a frame without memory-mapping.

`allensdk.internal.brain_observatory.roi_filter` module

`allensdk.internal.brain_observatory.roi_filter_utils` module

`allensdk.internal.brain_observatory.roi_filter_utils.CRITERIA()`

class `allensdk.internal.brain_observatory.roi_filter_utils.TrainingLabelClassifier(criteria)`
Bases: `object`

Very basic threshold_based classifier.

Has a decision function that is just the number of distinct criteria met by the classifier. Criteria are defined as a list of strings used with `pandas.DataFrame.eval`.

Parameters

criteria [list] List of evaluation strings.

decision_function (*self*, *X*)

Get the distance from the decision boundary.

Parameters

X [array-like] Features for each ROI.

Returns

T [array-like] Distance for each sample from the decision boundary.

class `allensdk.internal.brain_observatory.roi_filter_utils.TrainingMultiLabelClassifier(criteria)`
Bases: `object`

Multilabel classifier using groups of `TrainingLabelClassifiers`.

This was used to generate labeling for training the original SVM for classification.

Parameters

criteria [dictionary] Label names and criteria for each label.

get_eXcluded (*self*, *X*)

Get the calculated value of the eXcluded column.

This is useful for comparison with the original classifier implementation.

Parameters

X [`pandas.DataFrame`] Object features from the object list file.

Returns

numpy.ndarray Calculated eXcluded score from the classifier.

label_data (*self, X, as_columns=True*)
Generate labels for each row in X.

Parameters

X [pandas.DataFrame] Object features from the object list file.

Returns

numpy.ndarray Array of label codes representing the combination of labels found for each row.

`allensdk.internal.brain_observatory.roi_filter_utils.calculate_max_border` (*motion_df, max_shift*)

Calculate motion boundary from frame offsets.

When the motion correction algorithm fails to find sufficient matches, it generates very large frame offsets. The use of *max_shift* avoids filtering too many cells due to the large offsets, with the tradeoff that those frames will be noise.

Parameters

motion_df [pandas.DataFrame] Dataframe containing the x, y offsets from motion correction.

max_shift [float] Maximum shift to allow when considering motion correction. Any larger shifts are considered outliers.

Returns

list [right_shift, left_shift, down_shift, up_shift]

`allensdk.internal.brain_observatory.roi_filter_utils.get_indices_by_distance` (*object_list_points, mask_points*)

Find indices of nearest neighbor matches.

Require a distance of 0 (perfect match) and a unique match between masks and object_list entries.

`allensdk.internal.brain_observatory.roi_filter_utils.get_rois` (*segmentation_stack, border=None*)

Extract a list of rois from the segmentation data array.

Parameters

segmentation_stack [numpy.ndarray] The array from the maxInt_masks file showing the object masks.

border [list] [right_shift, left_shift, down_shift, up_shift] bounding box determined from motion correction.

Returns

list List of RoiMask objects.

`allensdk.internal.brain_observatory.roi_filter_utils.order_rois_by_object_list` (*object_data, rois*)

Reorder rois by matching bounding boxes to object list.

Parameters

object_data [pandas.DataFrame] Object list data.

rois [list] List of RoiMasks.

Returns

list The list of rois reordered to index the same as object_data.

allensdk.internal.brain_observatory.run_itracker module**allensdk.internal.brain_observatory.time_sync module**

```
class allensdk.internal.brain_observatory.time_sync.OphysTimeAligner(sync_file,  
                                                                    scan-  
                                                                    ner=None,  
                                                                    dff_file=None,  
                                                                    stimu-  
                                                                    lus_pkl=None,  
                                                                    eye_video=None,  
                                                                    behav-  
                                                                    ior_video=None,  
                                                                    long_stim_threshold=0.2)
```

Bases: object

behavior_video_timestamps

clipped_stim_timestamps

Return the stimulus timestamps with the erroneous initial spike removed (if relevant)

Returns

timestamps: np.ndarray An array of stimulus timestamps in seconds with the monitor delay added

delta: int Difference between the length of timestamps and the number of frames reported in the stimulus pickle file, i.e. `len(timestamps) - len(pkl_file['items']['behavior']['intervalsms'])`

corrected_behavior_video_timestamps

corrected_eye_video_timestamps

corrected_ophys_timestamps

corrected_stim_timestamps

The stimulus timestamps corrected for monitor delay

Returns

timestamps: np.ndarray An array of stimulus timestamps in seconds with the monitor delay added

delta: int Difference between the length of timestamps and the number of frames reported in the stimulus pickle file, i.e. `len(timestamps) - len(pkl_file['items']['behavior']['intervalsms'])`

delay: float The monitor delay in seconds

dataset

eye_video_timestamps

monitor_delay

The monitor delay (in seconds) associated with the session

ophys_timestamps

Get the timestamps for the ophys data.

stim_timestamps

```
allensdk.internal.brain_observatory.time_sync.calculate_monitor_delay(sync_dset,
                                                                    stim_times,
                                                                    pho-
                                                                    todi-
                                                                    ode_key,
                                                                    tran-
                                                                    si-
                                                                    tion_frame_interval=60,
                                                                    max_monitor_delay=0.07)
```

Calculate monitor delay.

```
allensdk.internal.brain_observatory.time_sync.corrected_video_timestamps(video_name,
                                                                    times-
                                                                    tamps,
                                                                    data_length)
```

```
allensdk.internal.brain_observatory.time_sync.get_alignment_array(ref, other,
                                                                    int_method=<ufunc
                                                                    'floor'>)
```

Generate an alignment array

```
allensdk.internal.brain_observatory.time_sync.get_keys(sync_dset:          al-
                                                                    lensdk.brain_observatory.sync_dataset.Dataset)
                                                                    → dict
```

Gets the correct keys for the sync file by searching the sync file line labels. Removes key from the dictionary if it is not in the sync dataset line labels. Args:

sync_dset: The sync dataset to search for keys within

Returns:

key_dict: dictionary of key value pairs for finding data in the sync file

```
allensdk.internal.brain_observatory.time_sync.get_ophys_data_length(filename)
allensdk.internal.brain_observatory.time_sync.get_photodiode_events(sync_dset,
                                                                    photodi-
                                                                    ode_key)
```

Returns the photodiode events with the start/stop indicators and the window init flash stripped off. These transitions occur roughly ~1.0s apart, since the sync square changes state every N frames (where N = 60, and frame rate is 60 Hz). Because there are no markers for when the first transition of this type started, we estimate based on the event intervals. For the first valid event, find the first two events that both meet the following criteria:

The next event occurs ~1.0s later

First the last valid event, find the first two events that both meet the following criteria:

The last valid event occurred ~1.0s before

```
allensdk.internal.brain_observatory.time_sync.get_real_photodiode_events(sync_dset,
                                                                    pho-
                                                                    to-
                                                                    di-
                                                                    ode_key,
                                                                    anomaly_threshold=0.5)
```

Gets the photodiode events with the anomalies removed.

```
allensdk.internal.brain_observatory.time_sync.get_stim_data_length(filename:
                                                                    str) → int
```

Get stimulus data length from .pkl file.

Parameters

filename [str] Path of stimulus data .pkl file.

Returns

int Stimulus data length.

```
allensdk.internal.brain_observatory.time_sync.get_video_length(filename)
```

Module contents

allensdk.internal.core package

Submodules

allensdk.internal.core.lims_pipeline_module module

```
class allensdk.internal.core.lims_pipeline_module.PipelineModule(description="",
                                                                    parser=None)
```

Bases: object

args

input_data (*self*)

write_output_data (*self*, *data*)

```
allensdk.internal.core.lims_pipeline_module.default_argument_parser(description="")
```

```
allensdk.internal.core.lims_pipeline_module.run_module(module,          input_data,
                                                         storage_directory,    op-
                                                         tional_args=None,
                                                         python='/shared/utils.x86_64/python-
                                                         2.7/bin/python',
                                                         sdk_path='/shared/bioapps/infoapps/lims2_modules/
                                                         local=False, pbs=None)
```

allensdk.internal.core.lims_utilities module

```
allensdk.internal.core.lims_utilities.append_well_known_file(wkfs,          path,
                                                             wkf_type_id=None,
                                                             con-
                                                             tent_type=None)
```

```
allensdk.internal.core.lims_utilities.connect(user='limsreader',          host='limsdb2',
                                              database='lims2',      password='limsro',
                                              port=5432)
```

```
allensdk.internal.core.lims_utilities.convert_from_titan_linux(file_name)
```

```
allensdk.internal.core.lims_utilities.get_input_json(object_id, object_class, strat-
                                                         egy_class,          host='lims2',
                                                         **kwargs)
```

```

allensdk.internal.core.lims_utilities.get_well_known_file_by_name(wkfs, file-
                                                                    name)
allensdk.internal.core.lims_utilities.get_well_known_file_by_type(wkfs,
                                                                    wkf_type_id)
allensdk.internal.core.lims_utilities.get_well_known_files_by_name(wkfs, file-
                                                                    name)
allensdk.internal.core.lims_utilities.get_well_known_files_by_type(wkfs,
                                                                    wkf_type_id)
allensdk.internal.core.lims_utilities.linux_to_windows(file_name)
allensdk.internal.core.lims_utilities.query(query, user='limsreader', host='limsdb2',
                                             database='lims2', password='limsro',
                                             port=5432)
allensdk.internal.core.lims_utilities.safe_system_path(file_name)
allensdk.internal.core.lims_utilities.select(cursor, query)

```

allensdk.internal.core.mouse_connectivity_cache_prerelease module

```
class allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePre
```

Bases: `allensdk.core.mouse_connectivity_cache.MouseConnectivityCache`

Extends MouseConnectivityCache to use prereleased data from lims.

Parameters

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

ccf_version: string Desired version of the Common Coordinate Framework. This affects the annotation volume (`get_annotation_volume`) and structure masks (`get_structure_mask`). Must be one of (`MouseConnectivityApi.CCF_2015`, `MouseConnectivityApi.CCF_2016`). Default: `MouseConnectivityApi.CCF_2016`

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. `get_projection_density(file_name='file.nrrd')`).

manifest_file: string File name of the manifest to be read. Default is "mouse_connectivity_manifest.json".

Attributes

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

api: MouseConnectivityApiPrerelease instance Used internally to make API queries.

EXPERIMENTS_PRERELEASE_KEY = 'EXPERIMENTS_PRERELEASE'

STORAGE_DIRECTORIES_PRERELEASE_KEY = 'STORAGE_DIRECTORIES_PRERELEASE'

add_manifest_paths (*self, manifest_builder*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

filter_experiments (*self, experiments, cre=None, injection_structure_ids=None, age=None, gender=None, workflow_state=None, workflows=None, project_code=None*)

Take a list of experiments and filter them by cre status and injection structure.

Parameters

cre: boolean or list If True, return only cre-positive experiments. If False, return only cre-negative experiments. If None, return all experiments. If list, return all experiments with cre line names in the supplied list. Default None.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If None, return all experiments. Default None.

age [list] Only return experiments with specimens with ages provided here. If None, return all experiments. Default None.

get_experiments (*self, dataframe=False, file_name=None, cre=None, injection_structure_ids=None, age=None, gender=None, workflow_state=None, workflows=None, project_code=None*)

Read a list of experiments.

If caching is enabled, this will save the whole (unfiltered) list of experiments to a file.

Parameters

dataframe: boolean Return the list of experiments as a Pandas DataFrame. If False, return a list of dictionaries. Default False.

file_name: string File name to save/read the structures table. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

allensdk.internal.core.simpletree module

class allensdk.internal.core.simpletree.**SimpleTree** (*nodes, node_id_cb, parent_id_cb*)

Bases: object

ancestor_ids (*self, nid*)

ancestors (*self, nid*)

child_ids (*self, nid*)

children (*self, nid*)

descendant_ids (*self, nid*)


```

descendants (self, nid)
node (self, nid)
node_ids (self)
nodes (self, nids=None)
parent (self, nid)
parent_id (self, nid)

```

allensdk.internal.core.swc module

```

class allensdk.internal.core.swc.Marker (*args, **kwargs)
    Bases: dict

    Simple dictionary class for handling reconstruction marker objects.

    CUT_DENDRITE = 10
    NO_RECONSTRUCTION = 20
    SPACING = [0.1144, 0.1144, 0.28]

allensdk.internal.core.swc.read_marker_file (file_name)
    read in a marker file and return a list of dictionaries

allensdk.internal.core.swc.read_swc (file_name)
    Read in an SWC file and return a Morphology object.

```

Parameters

file_name: string SWC file name.

Returns

Morphology A Morphology instance.

Module contents

allensdk.internal.ephys package

Submodules

allensdk.internal.ephys.core_feature_extract module

```

allensdk.internal.ephys.core_feature_extract.extract_data (data, nwb_file)
allensdk.internal.ephys.core_feature_extract.filter_sweeps (sweeps, types=None,
                                                             passed_only=True,
                                                             iclamp_only=True)
allensdk.internal.ephys.core_feature_extract.filtered_sweep_numbers (sweeps,
                                                                        types=None,
                                                                        passed_only=True,
                                                                        iclamp_only=True)

```

`allensdk.internal.ephys.core_feature_extract.find_coarse_long_square_amp_delta` (*sweeps, dec-i-mals=0*)

Find the delta between amplitudes of coarse long square sweeps. Includes failed sweeps.

`allensdk.internal.ephys.core_feature_extract.find_stim_start` (*stim, idx0=0*)

Find the index of the first nonzero positive or negative jump in an array.

Parameters

stim: `np.ndarray` Array to be searched

idx0: `int` Start searching with this index (default: 0).

Returns

`int`

`allensdk.internal.ephys.core_feature_extract.find_sweep_stim_start` (*data_set, sweep_number*)

`allensdk.internal.ephys.core_feature_extract.generate_output_cell_features` (*cell_features, sweep_features, sweep_index*)

`allensdk.internal.ephys.core_feature_extract.nan_get` (*obj, key*)

Return a value from a dictionary. If it does not exist, return None. If it is NaN, return None

`allensdk.internal.ephys.core_feature_extract.save_qc_figures` (*qc_fig_dir, nwb_file, out-put_data, plot_cell_figures*)

`allensdk.internal.ephys.core_feature_extract.update_output_sweep_features` (*cell_features, sweep_features, sweep_index*)

`allensdk.internal.ephys.plot_qc_figures` module

`allensdk.internal.ephys.plot_qc_figures.exp_curve` (*x, a, inv_tau, y0*)

Function used for tau curve fitting

`allensdk.internal.ephys.plot_qc_figures.get_features` (*sweep_features, sweep_number*)

`allensdk.internal.ephys.plot_qc_figures.get_spikes` (*sweep_features, sweep_number*)

`allensdk.internal.ephys.plot_qc_figures.get_time_string` ()

`allensdk.internal.ephys.plot_qc_figures.load_experiment` (*file_name, sweep_number*)

`allensdk.internal.ephys.plot_qc_figures.main` ()

`allensdk.internal.ephys.plot_qc_figures.make_cell_html` (*image_files, ephys_roi_result, file_name, relative_sweep_link*)

`allensdk.internal.ephys.plot_qc_figures.make_cell_page` (*nwb_file, ephys_roi_result, working_dir, save_cell_plots=True*)

`allensdk.internal.ephys.plot_qc_figures.make_sweep_html` (*sweep_files, file_name*)

```

allensdk.internal.ephys.plot_qc_figures.make_sweep_page(nwb_file, ephys_roi_result,
                                                         working_dir)

allensdk.internal.ephys.plot_qc_figures.mask_nulls(data)

allensdk.internal.ephys.plot_qc_figures.plot_cell_figures(nwb_file,
                                                         ephys_roi_result,
                                                         image_dir, sizes)

allensdk.internal.ephys.plot_qc_figures.plot_fi_curve_figures(nwb_file,
                                                             cell_features,
                                                             lims_features,
                                                             sweep_features,
                                                             image_dir, sizes,
                                                             cell_image_files)

allensdk.internal.ephys.plot_qc_figures.plot_hero_figures(nwb_file, cell_features,
                                                         lims_features,
                                                         sweep_features,
                                                         image_dir,      sizes,
                                                         cell_image_files)

allensdk.internal.ephys.plot_qc_figures.plot_images(ephys_roi_result,      image_dir,
                                                    sizes, image_sets)

allensdk.internal.ephys.plot_qc_figures.plot_instantaneous_threshold_thumbnail(nwb_file,
                                                                              sweep_numbers,
                                                                              cell_features,
                                                                              lims_features,
                                                                              sweep_features,
                                                                              color='red')

allensdk.internal.ephys.plot_qc_figures.plot_long_square_summary(nwb_file,
                                                                cell_features,
                                                                lims_features,
                                                                sweep_features)

allensdk.internal.ephys.plot_qc_figures.plot_ramp_figures(nwb_file,
                                                         cell_specimen,
                                                         cell_features,
                                                         lims_features,
                                                         sweep_features,
                                                         image_dir,      sizes,
                                                         cell_image_files)

allensdk.internal.ephys.plot_qc_figures.plot_rheo_figures(nwb_file, cell_features,
                                                         lims_features,
                                                         sweep_features,
                                                         image_dir,      sizes,
                                                         cell_image_files)

allensdk.internal.ephys.plot_qc_figures.plot_sag_figures(nwb_file, cell_features,
                                                         lims_features,
                                                         sweep_features,
                                                         image_dir,      sizes,
                                                         cell_image_files)

```

```
allensdk.internal.ephys.plot_qc_figures.plot_short_square_figures(nwb_file,  
                                                                    cell_features,  
                                                                    lims_features,  
                                                                    sweep_features,  
                                                                    image_dir,  
                                                                    sizes,  
                                                                    cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_single_ap_values(nwb_file,  
                                                                sweep_numbers,  
                                                                lims_features,  
                                                                sweep_features,  
                                                                cell_features,  
                                                                type_name)  
  
allensdk.internal.ephys.plot_qc_figures.plot_subthreshold_long_square_figures(nwb_file,  
                                                                                cell_features,  
                                                                                lims_features,  
                                                                                sweep_features,  
                                                                                im-  
                                                                                age_dir,  
                                                                                sizes,  
                                                                                cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_sweep_figures(nwb_file,  
                                                            ephys_roi_result,  
                                                            image_dir, sizes)  
  
allensdk.internal.ephys.plot_qc_figures.plot_sweep_set_summary(nwb_file, high-  
                                                                light_sweep_number,  
                                                                sweep_numbers,  
                                                                high-  
                                                                light_color='#0779BE',  
                                                                back-  
                                                                ground_color='#dddddd')
```

allensdk.internal.ephys.plot_qc_figures3 module

```
allensdk.internal.ephys.plot_qc_figures3.exp_curve(x, a, inv_tau, y0)  
    Function used for tau curve fitting  
  
allensdk.internal.ephys.plot_qc_figures3.get_features(sweep_features,  
                                                       sweep_number)  
  
allensdk.internal.ephys.plot_qc_figures3.get_spikes(sweep_features, sweep_number)  
  
allensdk.internal.ephys.plot_qc_figures3.get_time_string()
```

```

allensdk.internal.ephys.plot_qc_figures3.load_experiment (file_name,
                                                         sweep_number)
allensdk.internal.ephys.plot_qc_figures3.make_cell_html (image_files, file_name,
                                                         relative_sweep_link,
                                                         specimen_info, fields)
allensdk.internal.ephys.plot_qc_figures3.make_cell_page (nwb_file, cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         sweep_info,
                                                         well_known_files,
                                                         specimen_info, working_dir,
                                                         fields_to_show,
                                                         save_cell_plots=True)

nwb_file: name of nwb file (string)
cell_features:
rheo_features: dict containing extracted features from rheobase sweep
sweep_features:
sweep_info:
well_known_files: LIMS-output information containing graphics file names
working_dir:
save_cell_plots:

allensdk.internal.ephys.plot_qc_figures3.make_sweep_html (sweep_files, file_name)
allensdk.internal.ephys.plot_qc_figures3.make_sweep_page (nwb_file, working_dir,
                                                         sweep_data)
allensdk.internal.ephys.plot_qc_figures3.mask_nulls (data)
allensdk.internal.ephys.plot_qc_figures3.plot_cell_figures (nwb_file,
                                                         cell_features,
                                                         sweep_features,
                                                         rheo_features, image_dir,
                                                         sweep_info,
                                                         sizes)
allensdk.internal.ephys.plot_qc_figures3.plot_fi_curve_figures (nwb_file,
                                                         cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         image_dir,
                                                         sizes,
                                                         cell_image_files)
allensdk.internal.ephys.plot_qc_figures3.plot_hero_figures (nwb_file,
                                                         cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         image_dir, sizes,
                                                         cell_image_files)
allensdk.internal.ephys.plot_qc_figures3.plot_images (well_known_files, image_dir,
                                                         sizes, image_sets)

```

```
allensdk.internal.ephys.plot_qc_figures3.plot_instantaneous_threshold_thumbnail(nwb_file,
                                                                              sweep_numbers,
                                                                              cell_features,
                                                                              rheo_features,
                                                                              sweep_features,
                                                                              color='red')

allensdk.internal.ephys.plot_qc_figures3.plot_long_square_summary(nwb_file,
                                                                    cell_features,
                                                                    rheo_features,
                                                                    sweep_features)

allensdk.internal.ephys.plot_qc_figures3.plot_ramp_figures(nwb_file, sweep_info,
                                                            cell_features,
                                                            rheo_features,
                                                            sweep_features,
                                                            image_dir,      sizes,
                                                            cell_image_files)

allensdk.internal.ephys.plot_qc_figures3.plot_rheo_figures(nwb_file,
                                                            cell_features,
                                                            rheo_features,
                                                            sweep_features,
                                                            image_dir,      sizes,
                                                            cell_image_files)

allensdk.internal.ephys.plot_qc_figures3.plot_sag_figures(nwb_file, cell_features,
                                                            rheo_features,
                                                            sweep_features,
                                                            image_dir,      sizes,
                                                            cell_image_files)

allensdk.internal.ephys.plot_qc_figures3.plot_short_square_figures(nwb_file,
                                                                    cell_features,
                                                                    rheo_features,
                                                                    sweep_features,
                                                                    im-
                                                                    age_dir,
                                                                    sizes,
                                                                    cell_image_files)

allensdk.internal.ephys.plot_qc_figures3.plot_single_ap_values(nwb_file,
                                                                sweep_numbers,
                                                                rheo_features,
                                                                sweep_features,
                                                                cell_features,
                                                                type_name)

allensdk.internal.ephys.plot_qc_figures3.plot_subthreshold_long_square_figures(nwb_file,
                                                                                  cell_features,
                                                                                  rheo_features,
                                                                                  sweep_features,
                                                                                  im-
                                                                                  age_dir,
                                                                                  sizes,
                                                                                  cell_image_files)

allensdk.internal.ephys.plot_qc_figures3.plot_sweep_figures(nwb_file,
                                                              sweep_data,    im-
                                                              age_dir, sizes)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_sweep_set_summary (nwb_file, high-  
light_sweep_number,  
sweep_numbers,  
high-  
light_color='#0779BE',  
back-  
ground_color='#dddddd')  
  
allensdk.internal.ephys.plot_qc_figures3.plot_sweep_value_figures (sweep_info,  
image_dir,  
sizes,  
cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.save_figure (fig, image_name, im-  
age_set_name, image_dir,  
sizes, image_sets, scalew=1,  
scaleh=1, ext='jpg')
```

Module contents

allensdk.internal.model package

Subpackages

allensdk.internal.model.biophysical package

Subpackages

allensdk.internal.model.biophysical.fits package

Subpackages

allensdk.internal.model.biophysical.fits.fit_styles package

Module contents

Module contents

allensdk.internal.model.biophysical.passive_fitting package

Subpackages

allensdk.internal.model.biophysical.passive_fitting.passive package

Module contents

Submodules

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.arg_parser()
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.main()
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.process_inputs(parser)
```

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit2 module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit2.main()
```

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit_elec module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit_elec.main()
```

allensdk.internal.model.biophysical.passive_fitting.neuron_utils module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.get_h()
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.load_morphology(filename)
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.parse_neuron_output(output_s)
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.read_neuron_fit_stdout(functi
```

allensdk.internal.model.biophysical.passive_fitting.output_grabber module

```
class allensdk.internal.model.biophysical.passive_fitting.output_grabber.OutputGrabber(stream, thread)
    Bases: object
    Class used to grab standard output or another stream.
    escape_char = '\x08'
    readOutput(self)
        Read the stream data (one byte at a time) and save the text in capturedtext.
    start(self)
        Start capturing the stream data.
    stop(self)
        Stop capturing the stream data and save the text in capturedtext.
```

allensdk.internal.model.biophysical.passive_fitting.preprocess module

```
allensdk.internal.model.biophysical.passive_fitting.preprocess.get_cap_check_indices(i)
allensdk.internal.model.biophysical.passive_fitting.preprocess.get_passive_fit_data(cap_checks, data_set)
allensdk.internal.model.biophysical.passive_fitting.preprocess.main()
```


Module contents

Submodules

allensdk.internal.model.biophysical.biophysical_archiver module

```

class allensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver (archive_dir)
    Bases: object
    archive_cell (self, ephys_result_id, specimen_id, template, neuronal_model_id)
    get_cells (self)
    get_neuronal_models (self, specimen_ids)
    get_stimulus_file (self, neuronal_model_id)
    get_template_names (self)

```

allensdk.internal.model.biophysical.check_fi_shift module

```

allensdk.internal.model.biophysical.check_fi_shift.calculate_fi_curves (data_set,
                                                                    sweeps)
allensdk.internal.model.biophysical.check_fi_shift.estimate_fi_shift (data_set,
                                                                    sweeps)

```

allensdk.internal.model.biophysical.deap_utils module

```

class allensdk.internal.model.biophysical.deap_utils.Utils (description)
    Bases: allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils
    actual_parameters_from_normalized (self, params)
    calculate_feature_errors (self, t_ms, v, i)
    generate_morphology (self, morph_filename)
    insert_iclamp (self)
    load_cell_parameters (self)
    normalize_actual_parameters (self, params)
    record_values (self)
    set_actual_parameters (self, params)
    set_iclamp_params (self, amp, delay, dur)
    set_normalized_parameters (self, params)

```

allensdk.internal.model.biophysical.ephys_utils module

```

allensdk.internal.model.biophysical.ephys_utils.get_step_stim_characteristics (i,
                                                                    t)
allensdk.internal.model.biophysical.ephys_utils.get_sweep_v_i_t_from_set (data_set,
                                                                    sweep_number)

```

```
allensdk.internal.model.biophysical.ephys_utils.get_sweeps_of_type(sweep_type, sweeps)
```

allensdk.internal.model.biophysical.fit_stage_1 module

allensdk.internal.model.biophysical.fit_stage_2 module

allensdk.internal.model.biophysical.make_deap_fit_json module

```
class allensdk.internal.model.biophysical.make_deap_fit_json.Report(top_level_description, fit_type)
    Bases: object
    best_fit_value(self)
    check_org_selections_for_noise_block(self)
    gather_from_seeds(self)
    generate_fit_file(self)
    make_fit_json_file(self)
    setup_model(self)
```

allensdk.internal.model.biophysical.neuron_parallel module

allensdk.internal.model.biophysical.optimize module

allensdk.internal.model.biophysical.run_optimize module

```
class allensdk.internal.model.biophysical.run_optimize.RunOptimize(input_json, out-put_json)
    Bases: object
    copy_local(self)
    generate_manifest_lims(self, lims_json_path, manifest_path)
    generate_manifest_rma(self, neuronal_model_id, manifest_path, api_url=None)
    info(self, lims_json_path)
        return a string that a bash script can use to find the working directory, etc. to clean up.
    load_manifest(self)
    make_fit(self)
    nrnivmodl(self)
    start_specimen(self)
```

```
allensdk.internal.model.biophysical.run_optimize.main(command, input_json, out-put_json)
```

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.internal.model.biophysical.run_optimize_workflow module**allensdk.internal.model.biophysical.run_passive_fit module**

`allensdk.internal.model.biophysical.run_passive_fit.main` (*limit, manifest_path*)

`allensdk.internal.model.biophysical.run_passive_fit.run_passive_fit` (*description*)

allensdk.internal.model.biophysical.run_simulate_lims module

class `allensdk.internal.model.biophysical.run_simulate_lims.RunSimulateLims` (*input_json, out-put_json*)

Bases: `allensdk.model.biophysical.run_simulate.RunSimulate`

copy_local (*self*)

generate_manifest_lims (*self, lims_data_path, manifest_path*)

generate_manifest_rma (*self, neuronal_model_run_id, manifest_path, api_url=None*)

`allensdk.internal.model.biophysical.run_simulate_lims.main` (*command, lims_strategy_json, lims_response_json*)

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.internal.model.biophysical.run_simulate_workflow module**Module contents****allensdk.internal.model.glif package****Submodules****allensdk.internal.model.glif.ASGLM module****allensdk.internal.model.glif.MLIN module**

`allensdk.internal.model.glif.MLIN.MLIN` (*voltage, current, res, cap, dt, MAKE_PLOT=False, SHOW_PLOT=False, BLOCK=False, PUBLICATION_PLOT=False*)

voltage, current input:

voltage: numpy array of voltage with test pulse cut out
current: numpy array of stimulus with test pulse cut out

`allensdk.internal.model.glif.MLIN.autocorr` (*x*)

`allensdk.internal.model.glif.MLIN.exp_decay` (*time, amp, tau*)

`allensdk.internal.model.glif.MLIN.expsymm_cdf` (*v, dv*)

`allensdk.internal.model.glif.MLIN.expsymm_pdf` (*v, dv*)

```
allensdk.internal.model.glif.MLIN.find_bin_center(edges)
```

allensdk.internal.model.glif.are_two_lists_of_arrays_the_same module

```
allensdk.internal.model.glif.are_two_lists_of_arrays_the_same.are_two_lists_of_arrays_the_same
```

returns False if to lists of arrays are different. otherwise the function returns True.

allensdk.internal.model.glif.configure_model module

allensdk.internal.model.glif.error_functions module

```
allensdk.internal.model.glif.error_functions.MLIN_list_error(param_guess,
                                                             experiment,    in-
                                                             put_data)
```

allensdk.internal.model.glif.find_spikes module

```
allensdk.internal.model.glif.find_spikes.align_and_cut_spikes(voltage_list,
                                                             current_list,    dt,
                                                             spike_window=None)
```

This function aligns the spikes to some criteria and returns a current and voltage trace of of the spike over a time window. Also returns zero crossing, and threshold in reference to the aligned spikes.

```
allensdk.internal.model.glif.find_spikes.find_spikes_list(voltage_list, dt)
allensdk.internal.model.glif.find_spikes.find_spikes_list_old(voltage_list, dt)
allensdk.internal.model.glif.find_spikes.find_spikes_old(v, dt)
allensdk.internal.model.glif.find_spikes.find_spikes_ssq_list(voltage_list,
                                                             dt,    dv_cutoff,
                                                             thresh_frac)
```

allensdk.internal.model.glif.find_sweeps module

exception allensdk.internal.model.glif.find_sweeps.MissingSweepException

Bases: Exception

```
allensdk.internal.model.glif.find_sweeps.find_long_square_sweeps(sweeps)
```

```
allensdk.internal.model.glif.find_sweeps.find_noise_sweeps(sweeps)
```

Find 1) the noise1 sweeps 2) the noise2 sweeps 4) all noise sweeps

```
allensdk.internal.model.glif.find_sweeps.find_ramp_sweeps(sweeps)
```

Find 1) all ramp sweeps

2) all subthreshold ramps

3) all superthreshold ramps

```
allensdk.internal.model.glif.find_sweeps.find_ramp_to_rheo_sweeps(sweeps)
```

```
allensdk.internal.model.glif.find_sweeps.find_ranked_sweep(sweep_list, key, re-
                                                             verse=False)
```

```
allensdk.internal.model.glif.find_sweeps.find_short_square_sweeps(sweeps)
```

Find 1) all of the subthreshold short square sweeps

2) all of the superthreshold short square sweeps

3) the subthresholds short square sweep with maximum stimulus amplitude

```
allensdk.internal.model.glif.find_sweeps.find_sweeps(sweep_list)
```

```
allensdk.internal.model.glif.find_sweeps.get_sweep_numbers(sweep_list)
```

```
allensdk.internal.model.glif.find_sweeps.get_sweeps_by_name(sweeps,
                                                             sweep_type)
```

```
allensdk.internal.model.glif.find_sweeps.main()
```

```
allensdk.internal.model.glif.find_sweeps.organize_sweeps_by_name(sweeps,
                                                                  name)
```

```
allensdk.internal.model.glif.find_sweeps.parse_arguments()
```

allensdk.internal.model.glif.glif_experiment module

```
class allensdk.internal.model.glif.glif_experiment.GlifExperiment(neuron, dt,
                                                                  stim_list,
                                                                  resp_list,
                                                                  spike_time_steps,
                                                                  grid_spike_times,
                                                                  grid_spike_voltages,
                                                                  param_fit_names,
                                                                  **kwargs)
```

Bases: object

neuron_parameter_count (*self*)

run (*self*, *param_guess*)

This code will run the loaded neuron model in reference to the target neuron spikes. inputs:

self: is the instance of the neuron model and parameters alone with the values of the target spikes.

NOTE the values in each array of the self.gridSpikeIndexTarge_list and the self.interpolated_spike_times are in reference to the time start of of the stim in each induvidual array (not the universal time)

param_guess: array of scalars of the values that will be inserted into the mapping function below.

returns:

voltage_list: list of array of voltage values. NOTE: IF THE MODEL NEURON SPIKES BEFORE THE TARG
NOT BE CALCULATED THEREFORE THE RESULTING VECTOR WILL NOT BE AS
LONG AS THE TARGET AND ALSO WILL NOT MAKE SENSE WITH THE STIMULUS
UNLESS YOU CUT IT AND OUTPUT IT TOO.

grid_spike_times_list: interpolated_spike_time_list: an array of the actual times of the spikes.
NOTE: THESE TIMES ARE CALCULATED BY ADDING THE

TIME OF THE INDIVIDUAL SPIKE TO THE TIME OF THE LAST SPIKE.

gridISIFromLastTargSpike_list: list of arrays of spike times of the model in reference to the last target (biologi
spike (not in reference to sweep start)

interpolatedISIFromLastTargSpike_list: list of arrays of spike times of the model in reference to the last target spike (not in reference to sweep start)

voltageOfModelAtGridBioSpike_list: list of arrays of scalars that contain the voltage of the model neuron when the target or bio neuron spikes. **theshOfModelAtGridBioSpike_list:** list of arrays of scalars that contain the threshold of the model neuron when the target or bio neuron spikes.

run_base_model (*self*, *param_guess*)

This code will run the loaded neuron model. inputs:

self: is the instance of the neuron model and parameters alone with the values of the target spikes.

NOTE the values in each array of the `self.gridSpikeIndexTarge_list` and the `self.interpolated_spike_times` are in reference to the time start of of the stim in each induvidual array (not the universal time)

`param_guess`: array of scalars of the values that will be inserted into the mapping function below.

returns:

voltage_list: list of array of voltage values. **NOTE: IF THE MODEL NEURON SPIKES BEFORE THE TARG**
NOT BE CALCULATED THEREFORE THE RESULTING VECTOR WILL NOT BE AS
LONG AS THE TARGET AND ALSO WILL NOT MAKE SENSE WITH THE STIMULUS
UNLESS YOU CUT IT AND OUTPUT IT TOO.

`gridTime_list`: `interpolatedTime_list`: an array of the actual times of the spikes. NOTE: THESE
TIMES ARE CALCULATED BY ADDING THE

TIME OF THE INDIVIDUAL SPIKE TO THE TIME OF THE LAST SPIKE.

grid_ISI_list: list of arrays of spike times of the model in reference to the last target (biological)
spike (not in reference to sweep start)

interpolated_ISI_list: list of arrays of spike times of the model in reference to the last target (biological)
spike (not in reference to sweep start)

`grid_spike_voltage_list`: list of arrays of scalars that contain the voltage of the model neuron when
the target or bio neuron spikes. `grid_spike_threshold_list`: list of arrays of scalars that contain the
threshold of the model neuron when the target or bio neuron spikes.

set_neuron_parameters (*self*, *param_guess*)

Maps the parameter guesses to the coefficients of the model. input:

`param_guess` is vector of values. It is assumed that the length will be

allensdk.internal.model.glif.glif_optimizer module

```
class allensdk.internal.model.glif.glif_optimizer.GlifOptimizer (experiment, dt,
                                                             outer_iterations,
                                                             in-
                                                             ner_iterations,
                                                             sigma_outer,
                                                             sigma_inner,
                                                             param_fit_names,
                                                             stim,    xtol,
                                                             ftol,    inter-
                                                             nal_iterations,
                                                             bessel,    er-
                                                             ror_function=None,
                                                             er-
                                                             ror_function_data=None,
                                                             init_params=None)
```

Bases: object

evaluate (*self, x, dt_multiplier=100*)

initiate_unique_seed (*self, seed=None*)

randomize_parameter_values (*self, values, sigma*)

run_many (*self, iteration_finished_callback=None, seed=None*)

run_once (*self, param0*)

@param param0: a list of the initial guesses for the optimizer @return: tuple including parameters that optimize function and value - see fmin docs

run_once_bound (*self, low_bound, high_bound*)

@param low_bound: a scalar initial guess for the optimizer @param high_bound: a scalar high bound for the optimizer @return: tuple including parameters that optimize function and value - see fmin docs

to_dict (*self*)

allensdk.internal.model.glif.glif_optimizer_neuron module

```
exception allensdk.internal.model.glif.glif_optimizer_neuron.GlifBadInitializationException
```

Bases: Exception

Exception raised when voltage is above threshold at the beginning of a sweep. i.e. probably caused by the optimizer.

```
exception allensdk.internal.model.glif.glif_optimizer_neuron.GlifNeuronException (message,
                                                                              data)
```

Bases: Exception

Exception for catching simulation errors and reporting intermediate data.

```
class allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron (*args,
                                                                              **kwargs)
```

Bases: *allensdk.model.glif.glif_neuron.GlifNeuron*

Contains methods for running the neuron model in a “forced-spike” paradigm used during optimization.

TYPE = 'GLIF'

classmethod `from_dict(d)`

classmethod `from_dict_legacy(d)`

run_until_biological_spike (*self, voltage_t0, threshold_t0, AScurrents_t0, stimulus, response, start_index, after_end_index, bio_spike_time_steps*)

Run the neuron simulation over a segment of a stimulus given initial conditions for use in the “forced spike” optimization paradigm. [Note: the section of stimulus is meant to be between two biological neuron spikes. Thus the stimulus is during the interspike interval (ISI)]. The model is simulated until either the model spikes or the end of the segment is reached. If the model does not spike, a spike time is extrapolated past the end of the simulation segment.

This function also returns the initial conditions for the subsequent stimulus segment. In the forced spike paradigm there are several ways

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

stimulus [np.ndarray] the full stimulus array (not just the segment of data being simulated)

response [np.ndarray] the full response array (not just the segment of data being simulated)

start_index [int] index of global stimulus at which to start simulation

after_end_index [int] index of global stimulus *after* the last index to be simulated

bio_spike_time_steps [list] time steps of input spikes

Returns

dict

a dictionary containing: ‘voltage’: simulated voltage value ‘threshold’: simulated threshold values ‘AScurrent_matrix’: afterspike current values during the simulation ‘grid_model_spike_time’: model spike time (in units of dt) ‘interpolated_model_spike_time’: model spike time (in units of dt) interpolated between time steps ‘voltage_t0’: reset voltage value to be used in subsequent simulation interval ‘threshold_t0’: reset threshold value to be used in subsequent simulation interval ‘AScurrents_t0’: reset afterspike current value to be used in subsequent simulation interval ‘grid_bio_spike_model_voltage’: model voltage at the time of the input spike ‘grid_bio_spike_model_threshold’: model threshold at the time of the input spike

run_with_biological_spikes (*self, stimulus, response, bio_spike_time_steps*)

Run the neuron simulation over a stimulus, but do not allow the model to spike on its own. Rather, force the simulation to spike and reset at a given set of spike indices. Dynamics rules are applied between spikes regardless of the simulated voltage and threshold values. Reset rules are applied only at input spike times. This is used during optimization to force the model to follow the spikes of biological data. The model is optimized in this way so that history effects due to spiking can be adequately modeled. For example, every time the model spikes a new set of afterspike currents will be initiated. To ensure that afterspike currents can be optimized, we force them to be initiated at the time of the biological spike.

Parameters

stimulus [np.ndarray] vector of scalar current values

responses [np.ndarray] vector of scalar voltage values

bio_spike_time_steps [list] spike time step indices

Returns

dict

a dictionary containing: 'voltage': simulated voltage values, 'threshold': simulated threshold values, 'AScurrent_matrix': afterspike currents during the simulation, 'grid_model_spike_times': spike times of the model aligned to the simulation grid (when it would have spiked), 'interpolated_model_spike_times': spike times of the model linearly interpolated between time steps, 'grid_ISI': interspike interval between grid model spike times, 'interpolated_ISI': interspike interval between interpolated model spike times, 'grid_bio_spike_model_voltage': voltage of the model at biological/input spike times, 'grid_bio_spike_model_threshold': voltage of the model at biological/input spike times interpolated between time steps

to_dict (*self*)

Convert the neuron to a serializable dictionary.

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_model_spike_from_endpoints (n
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_model_spike_from_endpoints_s
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_spike_time (dt,
                                                                           num_time_steps,
                                                                           thresh-
                                                                           old_t0,
                                                                           thresh-
                                                                           old_t1,
                                                                           volt-
                                                                           age_t0,
                                                                           volt-
                                                                           age_t1)
```

Given two voltage and threshold values and an interval between them, extrapolate a spike time by intersecting lines the thresholds and voltages.

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_spike_voltage(dt,  
                                                                           num_time_steps,  
                                                                           thresh-  
                                                                           old_t0,  
                                                                           thresh-  
                                                                           old_t1,  
                                                                           volt-  
                                                                           age_t0,  
                                                                           volt-  
                                                                           age_t1)
```

Given two voltage and threshold values and an interval between them, extrapolate a spike time by intersecting lines the thresholds and voltages.

```
allensdk.internal.model.glif.glif_optimizer_neuron.find_first_model_spike(voltage,  
                                                                           thresh-  
                                                                           old,  
                                                                           volt-  
                                                                           age_t1,  
                                                                           thresh-  
                                                                           old_t1,  
                                                                           dt)
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.interpolate_spike_voltage(dt,  
                                                                           time_step,  
                                                                           thresh-  
                                                                           old_t0,  
                                                                           thresh-  
                                                                           old_t1,  
                                                                           volt-  
                                                                           age_t0,  
                                                                           volt-  
                                                                           age_t1)
```

Given two voltage and threshold values, the dt between them and the initial time step, interpolate a spike time within the dt interval by intersecting the two lines.

allensdk.internal.model.glif.optimize_neuron module

allensdk.internal.model.glif.plotting module

Written by Corinne Teeter 3-31-14

```
allensdk.internal.model.glif.plotting.checkPreprocess(originalStim_list,    pro-  
                                                         cessedStim_list,    orig-  
                                                         inalVoltage_list,    pro-  
                                                         cessedVoltage_list,  config,  
                                                         blockME=False)
```

```
allensdk.internal.model.glif.plotting.checkSpikeCutting(originalStim_list,    cut-  
                                                         Stim_list,    originalVolt-  
                                                         age_list,    cutVoltage_list,  
                                                         allindOfNonSpiking_list,  
                                                         config, blockME=False)
```

```
allensdk.internal.model.glif.plotting.plotLineRegress1(slope, intercept, r, xlim)
```

```
allensdk.internal.model.glif.plotting.plotLineRegressRed(slope, intercept, r, xlim)
```

```
allensdk.internal.model.glif.plotting.plotSpikes(voltage_list, spike_ind_list, dt,  
                                                blockME=False, method=False)
```

allensdk.internal.model.glif.preprocess_neuron module

allensdk.internal.model.glif.rc module

```
allensdk.internal.model.glif.rc.least_squares_RCEl_calc_tested(voltage_list,  
                                                                current_list, dt)
```

Calculate resistance, capacitance and resting potential by performing least squares on current and voltage.

Parameters

voltage_list: list of arrays voltage responses for several sweep repeats

current_list: list of arrays current injections for several sweep repeats

dt: float time step size in voltage and current traces

Returns

r_list: list of floats each value corresponds to the resistance of a sweep

c_list: list of floats each value corresponds to the capacitance of a sweep

el_list: list of floats each value corresponds to the resting potential of a sweep

allensdk.internal.model.glif.spike_cutting module

```
allensdk.internal.model.glif.spike_cutting.calc_spike_cut_and_v_reset_via_expvar_residuals
```

This function calculates where the spike should be cut based on explained variance. The goal is to find a model where the voltage after a spike maximally explains the voltage before a spike. This will also specify the voltage reset rule inputs:

spike_determination_method: string specifying the method used to find threshold
all_current_list: list of current (list of current traces injected into neuron)
all_voltage_list: list of voltages (list of voltage trace)

The change is that if the slope is greater than one or intercept is greater than zero it forces it. Regardless of required force the residuals are used.

```
allensdk.internal.model.glif.spike_cutting.plotLineRegress1(slope, intercept, r,  
                                                            xlim)
```

```
allensdk.internal.model.glif.spike_cutting.plotLineRegressRed(slope, intercept, r,  
                                                            xlim)
```

allensdk.internal.model.glif.threshold_adaptation module

`allensdk.internal.model.glif.threshold_adaptation.calc_spike_component_of_threshold_from_m`

Calculate the spike components of the threshold by fitting a decaying exponential function to data to threshold versus time since last spike in the multiblip data. The exponential is forced to decay to the local `th_inf` (calculated as the mean all of the threshold values of the first spikes in each individual triblip stimulus). For each multiblip stimulus in a stimulus set if there is more than one spike the difference in voltages from the first and second spike are plotted versus the separation in time. Note that this algorithm should only be implemented on multiblips sweeps where the neuron spike on the first and second blip. Since there is no easy way to do this, this erroneous data should not be provided to this algorithm (i.e is should be visually checked and eliminated the preprocessor should hold back this data manually for now.)

#TODO: check to see if this is still true. Notes: The standard SDK spike detection algorithm does not work with the multiblip stimulus due to artifacts when the stimulus turns on and off. Please see the `find_multiblip_spikes` module for more information.

Input:

multi_SS: dictionary contains multiblip information such as current and stimulus

dt: float time step in seconds

Returns:

const_to_add_to_thresh_for_reset: float amplitude of the exponential fit otherwise known as `a_spike`. Note that this is without any spike cutting

decay_const: float decay constant of exponential. Note the function fit is a negative exponential which will mean this value will either have to be negated when it is used or the functions used will have to have to include the negative.

`thresh_inf: float`

`allensdk.internal.model.glif.threshold_adaptation.exp_fit_c(t, a1, k1, const)`

`allensdk.internal.model.glif.threshold_adaptation.exp_force_c(t_const, a1, k1)`

`allensdk.internal.model.glif.threshold_adaptation.fit_avoltage_bvoltage(x,`

`v_trace_list,`
`El_list,`
`spike_cut_length,`
`all_spikeInd_list,`
`th_inf,`
`dt,`
`a_spike,`
`b_spike,`
`fake=False)`

This is a version of `fit_avoltage_bvoltage_debug` that does not require the `th_trace,`

`v_component_of_thresh_trace`, and `spike_component_of_thresh_trace` needed for debugging. A test should be run to make sure the same output comes out from this and the debug function

This function returns the squared error for the difference between the ‘known’ voltage component of the threshold obtained from the biological neuron and the voltage component of the threshold of the model obtained with the input parameters (so that the minimum can be searched for via `fmin`). The overall threshold is the sum of threshold infinity the spike component of the threshold and the voltage component of the threshold. Therefore threshold infinity and the spike component of the threshold must be subtracted from the threshold of the neuron in order to isolate the voltage component of the threshold. In the evaluation of the model the actual voltage of the neuron is used so that any errors in the other components of the model will not influence the fits here (for example, if a afterspike current was estimated incorrectly)

Notes: * The spike component of the threshold is subtracted from the

voltage which means that the voltage component of the threshold should only be added to rules.

- **b_spike was fit using a negative value in the function therefore the negative is placed in the equation.**
- **values in this function are in ‘real’ voltage as opposed to voltage** relative to resting potential.
- **current injection during the spike is not taken into account. This seems reasonable as the** ion channels are open during this time and injected current may not greatly influence the neuron.

x: numpy array `x[0]`=a_voltage input, `x[1]` is b_voltage_input, `x[2]` is th_inf

v_trace_list: list of numpy arrays voltage traces (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

El_list: list of floats reversal potential (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

spike_cut_length: int number of indices removed after initiation of a spike

all_spikeInd_list: list of numpy arrays indices of spike trains

th_inf: float threshold infinity (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

dt: float size of time step (SI units)

a_spike: float amplitude of spike component of threshold.

b_spike: float decay constant in spike component of the threshold

fake: Boolean if True makes uses the voltage value of spike step-1 because there is not a voltage value at the spike step because it is set to nan in the simulator.

```
allensdk.internal.model.glif.threshold_adaptation.fit_avoltage_bvoltage_th(x,
                                                                           v_trace_list,
                                                                           El_list,
                                                                           spike_cut_length,
                                                                           all_spikeInd_list,
                                                                           dt,
                                                                           a_spike,
                                                                           b_spike,
                                                                           fake=False)
```

This is a version of `fit_avoltage_bvoltage_th_debug` that does not require the `th_trace`, `v_component_of_thresh_trace`, and `spike_component_of_thresh_trace` needed for debugging. A test should be run to make sure the same output comes out from this and the debug function

This function returns the squared error for the difference between the ‘known’ voltage component of the threshold obtained from the biological neuron and the voltage component of the threshold of the model obtained with the input parameters (so that the minimum can be searched for via `fmin`). The overall threshold is the sum of

threshold infinity the spike component of the threshold and the voltage component of the threshold. Therefore threshold infinity and the spike component of the threshold must be subtracted from the threshold of the neuron in order to isolate the voltage component of the threshold. In the evaluation of the model the actual voltage of the neuron is used so that any errors in the other components of the model will not influence the fits here (for example, if a afterspike current was estimated incorrectly)

Notes: * The spike component of the threshold is subtracted from the

voltage which means that the voltage component of the threshold should only be added to rules.

- **b_spike** was fit using a negative value in the function therefore the negative is placed in the equation.
- **values in this function are in ‘real’ voltage as opposed to voltage** relative to resting potential.
- **current injection during the spike is not taken into account. This seems reasonable as the** ion channels are open during this time and injected current may not greatly influence the neuron.

x: numpy array x[0]=a_voltage input, x[1] is b_voltage_input, x[2] is th_inf

v_trace_list: list of numpy arrays voltage traces (v_trace, El, and th_inf must be in the same frame of reference)

El_list: list of floats reversal potential (v_trace, El, and th_inf must be in the same frame of reference)

spike_cut_length: int number of indicies removed after initiation of a spike

all_spikeInd_list: list of numpy arrays indicies of spike trains

dt: float size of time step (SI units)

a_spike: float amplitude of spike component of threshold.

b_spike: float decay constant in spike component of the threshold

fake: Boolean if True makes uses the voltage value of spike step-1 because there is not a voltage value at the spike step because it is set to nan in the simulator.

`allensdk.internal.model.glif.threshold_adaptation.get_peaks(voltage, above-Value=0)`

This function was written by Corinne Teeter and calculates the action potential peaks of a voltage equation” inputs

voltage: numpy array of voltages aboveValue: scalar voltage value over which voltage is considered a spike.

outputs: peakInd: array of indicies of peaks

Module contents

Submodules

`allensdk.internal.model.AIC` module

`allensdk.internal.model.AIC.AIC(RSS, k, n)`

Computes the Akaike Information Criterion.

RSS-residual sum of squares of the fitting errors. k - number of fitted parameters. n - number of observations.

```
allensdk.internal.model.AIC.AICc (RSS, k, n)
```

Corrected AIC. formula from Wikipedia.

```
allensdk.internal.model.AIC.BIC (RSS, k, n)
```

Bayesian information criterion or Schwartz information criterion. Formula from wikipedia.

allensdk.internal.model.GLM module

```
allensdk.internal.model.GLM.create_basis_IPSP (neye, ncos, kpeaks, ks, DTsim, t0, I_stim,  
                                                nkt, flag_exp, npcut)
```

```
allensdk.internal.model.GLM.ff (x, c, dc)
```

```
allensdk.internal.model.GLM.invnl (x)
```

```
allensdk.internal.model.GLM.makeBasis_StimKernel (kbasprs, nkt)
```

```
allensdk.internal.model.GLM.makeBasis_StimKernel_exp (kbasprs, nkt)
```

```
allensdk.internal.model.GLM.makeFitStruct_GLM (dtsim, kbasprs, nkt, flag_exp)
```

```
allensdk.internal.model.GLM.nlin (x)
```

```
allensdk.internal.model.GLM.normalizecols (A)
```

```
allensdk.internal.model.GLM.sameconv (A, B)
```

allensdk.internal.model.data_access module

```
allensdk.internal.model.data_access.load_sweep (file_name, sweep_number, de-  
                                                sired_dt=None, cut=0, bessel=False)
```

load a data sweep and do specified data processing. Inputs:

file_name: string name of .nwb data file

sweep_number: number specifying the sweep to be loaded

desired_dt: the size of the time step the data should be subsampled to

cut: indicie of which to start reporting data (i.e. cut off data before this indicie)

bessel: dictionary contains parameters 'N' and 'Wn' to implement standard python bessel filtering

Returns:

dictionary containing voltage: array current: array dt: time step of the returned data start_idx: the index at which the first stimulus starts (excluding the test pulse)

```
allensdk.internal.model.data_access.load_sweeps (file_name, sweep_numbers, dt=None,  
                                                cut=0, bessel=False)
```

load sweeps and do specified data processing. Inputs:

file_name: string name of .nwb data file

sweep_numbers: sweep numbers to be loaded

desired_dt: the size of the time step the data should be subsampled to

cut: indicie of which to start reporting data (i.e. cut off data before this indicie)

bessel: dictionary contains parameters 'N' and 'Wn' to implement standard python bessel filtering

Returns:

dictionary containing voltage: list of voltage trace arrays current: list of current trace arrays dt: list of time step corresponding to each array of the returned data start_idx: list of the indices at which the first stimulus starts (excluding the test pulse) in each returned sweep

```
allensdk.internal.model.data_access.subsample_data(data, method, present_time_step,
                                                    desired_time_step)
```

Module contents

allensdk.internal.morphology package

Submodules

allensdk.internal.morphology.compartment module

```
class allensdk.internal.morphology.compartment.Compartment (node1, node2)
    Bases: object
```

allensdk.internal.morphology.morphology module

```
class allensdk.internal.morphology.morphology.Morphology (node_list=None)
    Bases: object
```

Keep track of the list of nodes in a morphology and provide a few helper methods (soma, tree information, pruning, etc).

APICAL_DENDRITE = 4

AXON = 2

BASAL_DENDRITE = 3

NODE_TYPES = [1, 2, 3, 4]

SOMA = 1

append (*self*, *nodes*)

Add additional nodes to this Morphology. Those nodes must originate from another morphology object.

Parameters

nodes: list of Morphology nodes

apply_affine (*self*, *aff*, *scale=None*)

Apply an affine transform to all nodes in this morphology. Compartment radius is adjusted as well.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

apply_affine_only_rotation (*self*, *aff*)

Apply an affine transform to all nodes in this morphology. Only the rotation element of the transform is performed (i.e., although the entire transformation and translation matrix is supplied, only the rotation element is used). The morphology is translated to the point where the soma root is at 0,0,0.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

change_parent (*self*, *child*, *parent*)

Change the parent of a node. The child node is adjusted to point to the new parent, the child is taken off of the previous parent's child list, and it is added to the new parent's child list.

Parameters

child: integer or Morphology Object The ID of the child node, or the child node itself

parent: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

Nothing

children_of (*self*, *seg*)

Returns a list of the children of the specified node

Parameters

seg: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

A list of the child morphology objects. If the ID of the parent node is invalid, None is returned.

clone (*self*)

Create a clone (deep copy) of this morphology

compartment (*self*, *n*)

Returns the morphology Compartment having the specified ID.

Parameters

n: integer ID of desired compartment

Returns

A morphology object having the specified ID, or None if such a

node doesn't exist

compartment_list

convert_type (*self*, *from_type*, *to_type*)

Convert all nodes in morphology from one type to another

Parameters

from_type: enum The node type that will be eliminated and replaced. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

to_type: enum The new type that will replace it. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

delete_tree (*self*, *n*)

Delete tree, and all of its nodes, from the morphology.

Parameters

n: Integer The tree number to delete

find (*self*, *x*, *y*, *z*, *dist*, *node_type=None*)

Returns a list of Morphology Objects located within 'dist' of coordinate (x,y,z). If node_type is specified, the search will be constrained to return only nodes of that type.

Parameters

x, y, z: float The x,y,z coordinates from which to search around

dist: float The search radius

node_type: enum (optional) One of the following constants: SOMA, AXON, BASAL_DENDRITE or APICAL_DENDRITE

Returns

A list of all Morphology Objects matching the search criteria

get_dimensions (*self*)

Returns tuple of overall width, height and depth of morphology. WARNING: if locations of nodes in morphology are manipulated then this value can become incorrect. It can be reset and recalculated by programmitically setting self.dims to None.

Returns

3 real arrays: [width, height, depth], [min_x, min_y, min_z],

[max_x, max_y, max_z]

node (*self*, *n*)

Returns the morphology node having the specified ID.

Parameters

n: integer ID of desired node

Returns

A morphology node having the specified ID, or None if such a node doesn't exist

node_list

Return the node list. This is a property to ensure that the node list and node index are in sync.

node_list_by_type (*self*, *node_type*)

Return an list of all nodes having the specified node type.

Parameters

node_type: int Desired node type

Returns

A list of of Morphology Objects

num_nodes

Return the number of nodes in the morphology.

num_trees

Return the number of trees in the morphology. A tree is defined as everything following from a single root node.

parent_of (*self*, *seg*)

Returns parent of the specified node.

Parameters

seg: integer or Morphology Object The ID of the child node, or the child node itself

Returns

A morphology object, or None if no parent exists or if the specified node ID doesn't exist

save (*self*, *file_name*)

Write this morphology out to an SWC file

Parameters

file_name: string desired name of your SWC file

soma_root (*self*)

Returns root node of soma, if present

sparsify (*self*, *modulo*)

Return a new Morphology object that has a given number of non-leaf, non-root nodes removed.

Parameters

modulo: int keep 1 out of every modulo nodes.

Returns

Morphology A new morphology instance

strip_all_other_types (*self*, *node_type*, *keep_soma=True*)

Strips everything from the morphology except for the specified type. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The node type to keep in the morphology. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

keep_soma: Boolean (optional) True (default) if soma nodes should remain in the morphology, and False if the soma should also be stripped

strip_type (*self*, *node_type*)

Strips all nodes of the specified type from the morphology. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The node type to strip from the morphology. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

stumpify_axon (*self*, *count=10*)

Remove all axon nodes except the first ‘count’ nodes, as counted from the connected axon root.

Parameters

count: Integer The length of the axon ‘stump’, in number of nodes

to_dict (*self*)

Returns a dictionary of Node objects. These Nodes are a copy of the Morphology. Modifying them will not modify anything in the Morphology itself.

tree (*self*, *n*)

Returns a list of all Morphology nodes within the specified tree. A tree is defined as a fully connected graph of nodes. Each tree has exactly one root.

Parameters

n: integer ID of desired tree

Returns

**A list of all morphology objects in the specified tree, or None
if the tree doesn’t exist**

write (*self*, *file_name*)

allensdk.internal.morphology.morphvis module

class allensdk.internal.morphology.morphvis.**MorphologyColors**

Bases: object

set_apical_color (*self*, *r*, *g*, *b*)

set_axon_color (*self*, *r*, *g*, *b*)

set_basal_color (*self*, *r*, *g*, *b*)

set_soma_color (*self*, *r*, *g*, *b*)

allensdk.internal.morphology.morphvis.**calculate_scale** (*morph*, *pix_width*,
pix_height)

Calculates scaling factor and x,y insets required to auto-scale and center morphology into box with specified numbers of pixels

Parameters

morph: AISDK Morphology object

pix_width: int

Number of image pixels on X axis

pix_height: int

Number of image pixels on Y axis

Returns

real, real, real

First return value is the scaling factor. Second is the number of pixels needed to adjust x-coordinates so that the morphology is horizontally centered. Third is the number of pixels needed to adjust the y-coordinates so that the morphology is vertically centered.

```
allensdk.internal.morphology.morphvis.create_image(w, h, color=None, alpha=False)
allensdk.internal.morphology.morphvis.draw_density_hist(img, morph, vert_scale,
                                                         inset_left=0, inset_right=0, inset_top=0,
                                                         inset_bottom=0,
                                                         num_bins=None, colors=None)
```

Draws density histogram onto image When no scaling is applied, and no insets are provided, the coordinates of the morphology are used directly – i.e., 100 in morphology coordinates is equal to 100 pixels.

The scale factor is multiplied to morphology coordinates before being drawn. If scale_factor=2 then 50 in morphology coordinates is 100 pixels. Left and top insets shift the coordinate axes for drawing. E.g., if left=10 and top=5 then 0,0 in morphology coordinates is 10,5 in pixel space. Bottom and right insets are ignored.

If scale_to_fit is set then scale factor is ignored. The morphology is scaled to be the maximum size that fits in the image, taking into account insets. In a 100x100 image, if all insets=10, then the image is scaled to fit into the center 80x80 pixel area, and nothing is drawn in the inset border areas.

Axons are drawn before soma and dendrite compartments.

Parameters

img: PIL image object

morph: AISDK Morphology object

vert_scale: real

This is the amount required to multiply to a morphology

y-coordinate to convert it to relative cortical depth (on [0,1]).

This is the inverse of the cortical thickness.

inset_*: real

This is the number of pixels to use as border on top/bottom/

right/left. If scale_to_fit is false then only the top/left

values are used, as the scale_factor will determine how

large the morphology is (it can be drawn beyond insets and even beyond image boundaries)

num_bins: int

The number of bins in the histogram

colors: MorphologyColors object

This is the color scheme used to draw the morphology. If

colors=None then default coloring is used

Returns

Histogram arrays: [hist, hist2, hist3, hist4]

**where hist is the histogram of all neurites, and hist[234] are
the histograms of SWC types 2,3,4**

```
allensdk.internal.morphology.morphvis.draw_morphology(img, morph, inset_left=0,  
                                                         inset_right=0, inset_top=0, inset_bottom=0,  
                                                         scale_to_fit=False,  
                                                         scale_factor=1.0, colors=None)
```

Draws morphology onto image When no scaling is applied, and no insets are provided, the coordinates of the morphology are used directly – i.e., 100 in morphology coordinates is equal to 100 pixels.

The scale factor is multiplied to morphology coordinates before being drawn. If `scale_factor=2` then 50 in morphology coordinates is 100 pixels. Left and top insets shift the coordinate axes for drawing. E.g., if `left=10` and `top=5` then 0,0 in morphology coordinates is 10,5 in pixel space. Bottom and right insets are ignored.

If `scale_to_fit` is set then scale factor is ignored. The morphology is scaled to be the maximum size that fits in the image, taking into account insets. In a 100x100 image, if all insets=10, then the image is scaled to fit into the center 80x80 pixel area, and nothing is drawn in the inset border areas.

Axons are drawn before soma and dendrite compartments.

Parameters

img: PIL image object

morph: AISDK Morphology object

inset_*: real

**This is the number of pixels to use as border on top/bottom/
right/left. If `scale_to_fit` is false then only the top/left
values are used, as the `scale_factor` will determine how
large the morphology is (it can be drawn beyond insets and even
beyond image boundaries)**

scale_to_fit: boolean

**If true then morphology is scaled to the inset area of the
image and `scale_factor` is ignored. Morphology is centered
in the image in the sense that the top/bottom and left/right
edges of the morphology are equidistant from image borders.**

scale_factor: real

**A scalar amount that is multiplied to morphology coordinates
before drawing**

colors: MorphologyColors object

**This is the color scheme used to draw the morphology. If
`colors=None` then default coloring is used**

Returns

2-dimensional array, the pixel coordinates of the soma root [x,y]

allensdk.internal.morphology.node module

class allensdk.internal.morphology.node.**Node** (*n, t, x, y, z, r, pn, **kwargs*)

Bases: object

Represents node in SWC morphology file

classmethod **from_dict** (*d*)

short_string (*self*)

create string with node information in succinct, single-line form

to_dict (*self*)

Convert the node into a serializable dictionary

allensdk.internal.morphology.node.**euclidean_distance** (*node1, node2*)

allensdk.internal.morphology.node.**midpoint** (*node1, node2*)

allensdk.internal.morphology.validate_swc module

class allensdk.internal.morphology.validate_swc.**TestNode** (*n, t, x, y, z, r, pn*)

Bases: object

allensdk.internal.morphology.validate_swc.**main** ()

allensdk.internal.morphology.validate_swc.**resave_swc** (*orig_swc, new_file*)

Reads SWC file into AllenSDK Morphology object and resaves it. This can fix some problems in an SWC file that may disrupt other software tools reading the file (e.g., NEURON)

Parameters

orig_swc: string Name of SWC file to read

new_file: string Name of output SWC file

allensdk.internal.morphology.validate_swc.**validate_swc** (*swc_file*)

Tests SWC files for compatibility with AllenSDK

To be compatible with NEURON, SWC files must have the following properties:

- 1) a single root node with parent ID '-1'
- 2) sequentially increasing ID numbers
- 3) immediate children of the soma cannot branch

To be compatible with feature analysis, SWC files can only have node types in the range 1-4:

1 = soma 2 = axon 3 = [basal] dendrite 4 = apical dendrite

Module contents

allensdk.internal.mouse_connectivity package

Subpackages

allensdk.internal.mouse_connectivity.interval_unionize package

Submodules

allensdk.internal.mouse_connectivity.interval_unionize.cav_unionize module

allensdk.internal.mouse_connectivity.interval_unionize.cav_unionizer module

allensdk.internal.mouse_connectivity.interval_unionize.data_utilities module

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_cav_density(cav_data)`

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_injection_data(injection_data)`
Read nrrd files containing injection signal data

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_projection_data(projection_data)`
Read nrrd files containing global signal data

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_sum_pixel_intensities(sum_pixels)`
`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.load_annotation(annotation_data)`
Read data files segmenting the reference space into regions of valid and invalid data, then further among brain structures

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.read(path)`

allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer module

class `allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer.IntervalUnionizer`
Bases: `object`

direct_unionize (*self*, *data_arrays*, *pre_sorted=False*, ***kwargs*)
Obtain unionize records from directly annotated regions.

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened arrays.

sorted [bool, optional] If False, data arrays will be sorted.

extract_data (*self*, *data_arrays*, *low*, *high*, ***kwargs*)
Given flattened data arrays and a specified interval, generate summary data

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened, sorted arrays.

low [int] Index at which interval of interest begins. Inclusive.

high [int] Index at which interval of interest ends. Exclusive.

postprocess_unionizes (*self*, *raw_unionizes*, ***kwargs*)

Carry out additional calculations/formatting derivative of core unionization.

Parameters

raw_unionizes [list of unionizes] Each entry is a unionize record.

classmethod propagate_record (*child_record*, *ancestor_record*, *copy_all=False*)

Updates one unionize corresponding to a rootward structure with information from a unionize corresponding to a leafward structure

Parameters

child_record [unionize] Data will be drawn from this record

ancestor_record [unionize] This record will be updated

classmethod propagate_to_bilateral (*lateral_unionizes*)

classmethod propagate_unionizes (*direct_unionizes*, *ancestor_id_map*)

Structures are arranged in a tree, whose leafward-oriented edges indicate physical containment. This method updates rootward unionize records with information from leafward ones.

Parameters

direct_unionizes [list of unionizes] Each entry is a unionize record produced from a collection of directly labeled voxels in the segmentation volume.

ancestor_id_map [dict] Keys are structure ids. Values are ids of all structures rootward in

the tree, including the key node

Returns

output_unionizes [list of unionizes] Contains completed unionize records at all depths in the structure tree

classmethod record_cb ()

setup_interval_map (*self*, *annotation*)

Build a map from structure ids to intervals in the sorted flattened reference space.

Parameters

annotation [np.ndarray] Segmentation label array.

sort_data_arrays (*self*, *data_arrays*)

Apply the precomputed sort to flattened data arrays

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened, unsorted arrays.

Returns

dict : As input, but values are sorted

allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_cav module

allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic module

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.get`

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.get`

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.run`

allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record module

class `allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.Ti`

Bases: `allensdk.internal.mouse_connectivity.interval_unionize.unionize_record.Unionize`

direct_sum_projection_pixels

max_voxel_density

max_voxel_index

output (*self*, *output_spacing_iso*, *volume_scale*, *target_shape*, *sort*)
Generate derived data for this unionize

Parameters

output_spacing_iso [numeric] Isometric spacing of reference space in microns

volume_scale [numeric] Scale factor mapping pixels to microns³

target_shape [array-like of numeric] Shape of reference space

projection_density

projection_energy

projection_intensity

propagate (*self*, *ancestor*, *copy_all=False*)
Update a rootward unionize with data from this unionize record

Parameters

ancestor [TissuecyteBaseUnionize] will be updated

Returns

ancestor [TissuecyteBaseUnionize]

set_max_voxel (*self*, *density_array*, *low*)
Find the voxel of greatest density in this unionizes spatial domain

Parameters

density_array [ndarray] Float values are densities per voxel

low [int] index in full flattened, sorted array of starting voxel

sum_pixel_intensity

sum_pixels

sum_projection_pixel_intensity

sum_projection_pixels

```
class allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteUnionizeRecord
    Bases: allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteBaseUnionize
```

```
    calculate (self, low, high, data_arrays)
```

```
class allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteUnionizeRecord
    Bases: allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteBaseUnionize
```

```
    calculate (self, low, high, data_arrays, ij_record)
```

allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer module

```
class allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer.TissuecyteUnionizer
    Bases: allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer.IntervalUnionizer
```

A specialization of the IntervalUnionizer set up for unionizing Tissuecyte-derived projection data.

```
    extract_data (self, data_arrays, low, high)
        As parent
```

```
    postprocess_unionizes (self, raw_unionizes, image_series_id, output_spacing_iso, volume_scale, target_shape, sort)
        As parent
```

```
    classmethod propagate_record (child_record, ancestor_record, copy_all=False)
        As parent
```

```
    classmethod record_cb ()
```

allensdk.internal.mouse_connectivity.interval_unionize.unionize_record module

```
class allensdk.internal.mouse_connectivity.interval_unionize.unionize_record.UnionizeRecord (*args, **kwargs)
```

Bases: object

Abstract base class for unionize records.

```
    calculate (self, *args, **kwargs)
```

```
    output (self, *args, **kwargs)
```

```
    propagate (self, ancestor, copy_all, *args, **kwargs)
```

```
    slice_arrays (self, low, high, data_arrays)
        Extract a slice from several aligned arrays
```

Parameters

low [int] start of slice, inclusive

high [int] end of slice, exclusive

data_arrays [dict] keys are varieties of data. values are sorted, flattened data arrays

Module contents

`allensdk.internal.mouse_connectivity.projection_thumbnail` package

Submodules

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip` module

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.apply_c`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.blend_v`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.do_blur`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.handle`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.max_cb`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.run (vol`
`imin`
`ima`
`ro-`
`ta-`
`tion`
`col-`
`orn`

`allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.simple`

`allensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet` module

```
class allensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet.ImageSheet
    Bases: object
    append (self, new_cell)
    apply (self, fn, *args, **kwargs)
```

```

static build_from_image (image, n, axis)
copy (self)
get_output (self, axis)

```

allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions module

```

allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.convert_axis
allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.max_project

```

```

allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.template_pr

```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities module

```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.blend (ima
                                                                    wei

```

Parameters

image_stack :: list of np.ndarray The images to be blended. Shapes cannot differ

weight_stack :: list of np.ndarray The weight of each image at each pixel. Will be normal-
ized.

```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.convert_c

```

Generates a matplotlib continuous colormap on [0, 1] from a discrete colormap at N evenly spaced points.

Parameters

data [list of list] Sublists are [r, g, b].

Returns

matplotlib.colors.LinearSegmentedColormap Gamma is 1. Output space is 3 X [0, 1]

```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.minmax_n
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.normalize

```

```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.sitk_safe

```

allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector module

```
class allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector.VolumePro
    Bases: object

    build_rotation_transform (self, from_axis, to_axis, angle)

    extract (self, cb, volume=None)

    classmethod fixed_factory (volume, size)

    rotate (self, from_axis, to_axis, angle)

    rotate_and_extract (self, from_axes, to_axes, angles, cb)

    classmethod safe_factory (volume)
```

allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities module

```
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_center
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_diagonal
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_image_p
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_size_p
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_paste_into
```

Module contents**allensdk.internal.mouse_connectivity.tissuecyte_stitching package****Submodules****allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher module**

```
class allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.Stitcher (image_dimensions,
                                         tiles,
                                         average_tiles,
                                         channels)

    Bases: object

    run (self, cb=<built-in function array>)

    stitch (self, slice_image, stitched_indicator, tile, cb=<built-in function array>)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.blend_component_from_pos
```

Obtains a normalized component of the blend, which describes depth of overlap along a specified axis in a specified direction

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_blend(indicator_region,
                                                                           stap,
                                                                           cb=<built-
                                                                           in
                                                                           func-
                                                                           tion
                                                                           ar-
                                                                           ray>)
```

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_blend_component(indic
                                                                                       lg,
                                                                                       axis,
                                                                                       mesh)
```

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_indicator_bound_point
```

Finds the index of first change in a binary mask along a specified axis in a specified direction

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_overall_blend(indicator
                                                                                       meshes)
```

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.initialize_image(dimension
                                                                                       nchan-
                                                                                       nels,
                                                                                       dtype,
                                                                                       or-
                                                                                       der='C')
```

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.initialize_images(dimension
                                                                                       nchan-
                                                                                       nels)
```

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.make_blended_tile(blend,
                                                                                       tile,
                                                                                       cur-
                                                                                       rent_reg)
```

allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile module

```
class allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile(index,
                                                                           im-
                                                                           age,
                                                                           is_missing,
                                                                           bounds,
                                                                           chan-
                                                                           nel,
                                                                           size,
                                                                           mar-
                                                                           gins,
                                                                           *args,
                                                                           **kwargs)
```

Bases: object

```
apply_average_tile(self, average_tile)
```

```
apply_average_tile_to_self(self, average_tile)
```

```
average_tile_is_untrimmed(self, average_tile)
```

```
get_image_region (self)
get_missing_path (self)
initialize_image (self)
trim (self, image)
trim_self (self)
```

Module contents

Module contents

allensdk.internal.pipeline_modules package

Subpackages

allensdk.internal.pipeline_modules.gbm package

Submodules

allensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_run_records module

```
allensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_run_records.main (analysis_records_  
db_host,  
db_port,  
db_name,  
db_user,  
db_passwd)
```

allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap module

```
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_gene_fpkm_table (analysis_ru  
Creates a a matrix ("rows x columns = genes x samples") of fpkm gene expression values for each particular  
(gene, sample) pair. Rows are sorted by entrez_id and columns are by rna_well_id  
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_genes_for_transcripts (an  
Creates a list that contains the associated gene for each transcript sorted alphabetically  
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_sample_metadata (sample_me  
Creates a table of sample metadata sorted by rna_well_id  
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_transcript_fpkm_table (an  
Creates a a matrix ("rows x columns = transcripts x samples") of fpkm gene expression values for each particular  
(transcript, sample) pair. Rows are sorted by transcript id and columns are by rna_well_id  
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_transcripts_for_genes (an  
Creates a list that contains the associated transcript for each gene sorted by entrez_id  
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.main ()
```


allensdk.internal.pipeline_modules.gbm.generate_gbm_sample_metadata module

```
allensdk.internal.pipeline_modules.gbm.generate_gbm_sample_metadata.main(sample_metadata_json_lo-  
db_host,  
db_port,  
db_name,  
db_user,  
db_passwd)
```

Module contents

Submodules

allensdk.internal.pipeline_modules.run_annotated_region_metrics module

Run annotated region metrics calculations

```
allensdk.internal.pipeline_modules.run_annotated_region_metrics.debug(region_id,  
stor-  
age_directory='./',  
lo-  
cal=True,  
sdk_path='/data/informatics/C  
script_path='/data/informatics  
lims_host='lims2')  
  
allensdk.internal.pipeline_modules.run_annotated_region_metrics.load_arrays(h5_file)  
allensdk.internal.pipeline_modules.run_annotated_region_metrics.main()
```

allensdk.internal.pipeline_modules.run_demixing module

```
allensdk.internal.pipeline_modules.run_demixing.assert_exists(file_name)  
allensdk.internal.pipeline_modules.run_demixing.debug(experiment_id, local=False)  
allensdk.internal.pipeline_modules.run_demixing.get_path(obj, key, check_exists)  
allensdk.internal.pipeline_modules.run_demixing.main()  
allensdk.internal.pipeline_modules.run_demixing.parse_input(data, ex-  
clude_labels)
```

allensdk.internal.pipeline_modules.run_dff_computation module

```
allensdk.internal.pipeline_modules.run_dff_computation.main()  
allensdk.internal.pipeline_modules.run_dff_computation.parse_input(data)
```

allensdk.internal.pipeline_modules.run_eye_tracking module

allensdk.internal.pipeline_modules.run_neuropil_correction module

```
allensdk.internal.pipeline_modules.run_neuropil_correction.adjust_r_for_negativity(r,
                                                                                   F_C,
                                                                                   F_M,
                                                                                   F_N)

allensdk.internal.pipeline_modules.run_neuropil_correction.debug(experiment_id,
                                                                    local=False)

allensdk.internal.pipeline_modules.run_neuropil_correction.debug_plot(file_name,
                                                                           roi_trace,
                                                                           neu-
                                                                           ropil_trace,
                                                                           cor-
                                                                           rected_trace,
                                                                           r,
                                                                           r_vals=None,
                                                                           err_vals=None)

allensdk.internal.pipeline_modules.run_neuropil_correction.main()
```

allensdk.internal.pipeline_modules.run_observatory_analysis module

```
allensdk.internal.pipeline_modules.run_observatory_analysis.debug(experiment_ids,
                                                                    lo-
                                                                    cal=False,
                                                                    OUT-
                                                                    PUT_DIR='/data/informatics/CAM/
                                                                    SDK_PATH='/data/informatics/CAM
                                                                    wall-
                                                                    time='10:00:00',
                                                                    python='/shared/utils.x86_64/python
                                                                    2.7/bin/python',
                                                                    queue='braintv')

allensdk.internal.pipeline_modules.run_observatory_analysis.get_experiment_nwb_file(experiment_id)

allensdk.internal.pipeline_modules.run_observatory_analysis.get_experiment_session(experiment_id)

allensdk.internal.pipeline_modules.run_observatory_analysis.main()
```

allensdk.internal.pipeline_modules.run_observatory_container_thumbnails module

allensdk.internal.pipeline_modules.run_observatory_thumbnails module

```

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_cell_plots(cell_specimen_id,
                                         pre-
                                         fix,
                                         as-
                                         spect,
                                         con-
                                         figs,
                                         out-
                                         put_dir,
                                         axes=None,
                                         trans-
                                         par-
                                         ent=False)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_correlation_plots(data_series,
                                              anal-
                                              y-
                                              sis_file,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_drifting_gratings(dga,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_experiment_thumbnails(data_dir,
                                              a-
                                              n-
                                              a-
                                              l-
                                              y-
                                              s-
                                              i-
                                              s-
                                              f-
                                              i-
                                              l-
                                              e)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_eye_tracking_plots(data_dir,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_locally_sparse_noise(local_noise_dir,
                                              co-
                                              n-
                                              figs,
                                              ou-
                                              t-
                                              pu-
                                              t-
                                              o-
                                              n)

```

```
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_natural_movie(nma,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir,  
                                          name)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_natural_scenes(nsa,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_plots(prefix,  
                                   as-  
                                   pect,  
                                   con-  
                                   figs,  
                                   out-  
                                   put_dir,  
                                   axes=None,  
                                   trans-  
                                   par-  
                                   ent=False)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_receptive_field(lsna,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_speed_tuning(analysis,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_static_gratings(sga,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_type(nwb_file,  
                                   data_file,  
                                   con-  
                                   figs,  
                                   out-  
                                   put_dir,  
                                   type_name)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.debug(experiment_id,  
                              plots=None,  
                              lo-  
                              cal=False)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_analysis_file  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_files(experiment  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_nwb_file(experim
```

```
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_input_data(experiment_id)
allensdk.internal.pipeline_modules.run_observatory_thumbnails.lsna_check_hvas(data_set,
                                                                              data_file)
allensdk.internal.pipeline_modules.run_observatory_thumbnails.main()
allensdk.internal.pipeline_modules.run_observatory_thumbnails.parse_input(data)
```

allensdk.internal.pipeline_modules.run_ophys_eye_calibration module

```
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.debug(experiment_id,
                                                                    lo-
                                                                    cal=False)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.get_wkf(wkf_type,
                                                                        experi-
                                                                        ment_id)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.main()
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.parse_input_data(data)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.write_output(filename,
                                                                              po-
                                                                              si-
                                                                              tion_degrees,
                                                                              po-
                                                                              si-
                                                                              tion_cm,
                                                                              ar-
                                                                              eas)
```

allensdk.internal.pipeline_modules.run_ophys_session_decomposition module

```
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.convert_frame(conversion_de
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.create_fake_metadata(exp
raw
cha
nel
wid
heig
item
size
n_p

allensdk.internal.pipeline_modules.run_ophys_session_decomposition.debug(experiment_id,
                                                                    lo-
                                                                    cal=False,
                                                                    raw_path=None)
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.main()
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.parse_input(data)
    Load all input data from the input json.
```

allensdk.internal.pipeline_modules.run_ophys_time_sync module

```
class allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs
    Bases: tuple

    Schema for synchronization outputs

    behavior_alignment
        Alias for field number 12

    behavior_delta
        Alias for field number 5

    behavior_times
        Alias for field number 9

    experiment_id
        Alias for field number 0

    eye_alignment
        Alias for field number 11

    eye_delta
        Alias for field number 4

    eye_times
        Alias for field number 8

    ophys_delta
        Alias for field number 2

    ophys_times
        Alias for field number 6

    stimulus_alignment
        Alias for field number 10

    stimulus_delay
        Alias for field number 1

    stimulus_delta
        Alias for field number 3

    stimulus_times
        Alias for field number 7

class allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncWriter(output_h5_path:
                                                                    str,
                                                                    out-
                                                                    put_json_path:
                                                                    Op-
                                                                    tional[str]
                                                                    =
                                                                    None)

    Bases: object

    validate_paths(self)
        Determines whether we can actually write to the specified paths, allowing for creation of intermediate
        directories. It is a good idea to run this before doing any heavy calculations!

    write(self, outputs: allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs)
        Convenience for writing both an output h5 and (if applicable) an output json.
```

Parameters**outputs** [the data to be written]**write_output_h5** (*self*, *outputs*)

Write (mainly) heaviweight data to an h5 file.

Parameters**outputs** [the data to be written]**write_output_json** (*self*, *outputs*)

Write lightweight data to a json

Parameters**outputs** [the data to be written]

```
allensdk.internal.pipeline_modules.run_ophys_time_sync.check_stimulus_delay(obt_delay:
float,
min_delay:
float,
max_delay:
float)
```

Raise an exception if the monitor delay is not within specified bounds

Parameters**obt_delay** [obtained monitor delay (s)]**min_delay** [lower threshold (s)]**max_delay** [upper threshold (s)]

```
allensdk.internal.pipeline_modules.run_ophys_time_sync.main()
```

```
allensdk.internal.pipeline_modules.run_ophys_time_sync.run_ophys_time_sync(aligner:
al-
lensdk.internal.brain_d
ex-
per-
i-
ment_id:
int,
min_stimulus_delay:
float,
max_stimulus_delay:
float)
→
al-
lensdk.internal.pipeline
```

Carry out synchronization of timestamps across the data streams of an ophys experiment.

Parameters**aligner** [drives alignment. See OphysTimeAligner for details of the] attributes and properties that must be implemented.**experiment_id** [unique identifier for the experiment being aligned]**min_stimulus_delay** [reject alignment run (raise a ValueError) if the] calculated monitor delay is below this value (s).

max_stimulus_delay [reject alignment run (raise a ValueError) if the] calculated monitor delay is above this value (s).

Returns

A **TimeSyncOutputs** (see definition for more information) of **output** parameters and arrays of aligned timestamps.

`allensdk.internal.pipeline_modules.run_roi_filter` module

`allensdk.internal.pipeline_modules.run_tissuecyte_projection_thumbnail_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_stitching_classic` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_cav_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_counts_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_from_json.main()`

Module contents

Module contents

6.1.7 allensdk.model package

Subpackages

`allensdk.model.biophys_sim` package

Subpackages

`allensdk.model.biophys_sim.neuron` package

Submodules

`allensdk.model.biophys_sim.neuron.hoc_utils` module

class `allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils` (*description*)

Bases: `object`

A helper class for containing references to NEUORN.

Attributes

h [object] The NEURON hoc object.

nrn [object] The NEURON python object.

neuron [module] The NEURON module.


```

h = None
initialize_hoc (self)
    Basic setup for NEURON.

neuron = None
nrn = None

```

Module contents

allensdk.model.biophys_sim.scripts package

Module contents

Submodules

allensdk.model.biophys_sim.bps_command module

```

allensdk.model.biophys_sim.bps_command.choose_bps_command (command='bps_simple',
                                                             conf_file=None)

allensdk.model.biophys_sim.bps_command.run_module (description, module_name, func-
                                                    tion_name)

```

allensdk.model.biophys_sim.config module

class allensdk.model.biophys_sim.config.**Config**

Bases: *allensdk.config.app.application_config.ApplicationConfig*

load (*self*, *config_path*, *disable_existing_logs*=*False*)

Parse the application configuration then immediately load the model configuration files.

Parameters

disable_existing_logs [boolean, optional] If false (default) leave existing logs after configuration.

read_model_description (*self*)

parse the model_file field of the application configuration and read the files.

The model_file field of the application configuration is first split at commas, since it may list more than one file.

The files may be uris of the form *file:filename?section=name*, in which case a bare configuration object is read from filename into the configuration section with key 'name'.

A simple filename without a section option is treated as a standard multi-section configuration file.

Returns

description [Description] Configuration object.

Module contents

allensdk.model.biophysical package

Submodules

allensdk.model.biophysical.run_simulate module

class allensdk.model.biophysical.run_simulate.**RunSimulate** (*input_json*, *out-put_json*)

Bases: object

load_manifest (*self*)

nrnivmodl (*self*)

simulate (*self*)

allensdk.model.biophysical.run_simulate.**main** (*command*, *lims_strategy_json*, *lims_response_json*)

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.model.biophysical.runner module

allensdk.model.biophysical.runner.**load_description** (*args_dict*)

Read configurations.

Parameters

args_dict [dict] Parsed arguments dictionary with following keys.

manifest_file [string] .json file with containing the experiment configuration

axon_type [string] Axon handling for the all-active models

Returns

Config Object with all information needed to run the experiment.

allensdk.model.biophysical.runner.**prepare_nwb_output** (*nwb_stimulus_path*, *nwb_result_path*)

Copy the stimulus file, zero out the recorded voltages and spike times.

Parameters

nwb_stimulus_path [string] NWB file name

nwb_result_path [string] NWB file name

allensdk.model.biophysical.runner.**run** (*args*, *sweeps=None*, *procs=6*)

Main function for simulating sweeps in a biophysical experiment.

Parameters

args [dict] Parsed arguments to run the experiment.

procs [int] number of sweeps to simulate simultaneously.

sweeps [list] list of experiment sweep numbers to simulate. If None, simulate all sweeps.

`allensdk.model.biophysical.runner.run_sync` (*description*, *sweeps=None*)

Single-process main function for simulating sweeps in a biophysical experiment.

Parameters

description [Config] All information needed to run the experiment.

sweeps [list] list of experiment sweep numbers to simulate. If None, simulate all sweeps.

`allensdk.model.biophysical.runner.save_nwb` (*output_path*, *v*, *sweep*, *sweeps_by_type*)

Save a single voltage output result into an existing sweep in a NWB file. This is intended to overwrite a recorded trace with a simulated voltage.

Parameters

output_path [string] file name of a pre-existing NWB file.

v [numpy array] voltage

sweep [integer] which entry to overwrite in the file.

allensdk.model.biophysical.utils module

class `allensdk.model.biophysical.utils.AllActiveUtils` (*description*, *axon_type*)

Bases: `allensdk.model.biophysical.utils.Utils`

generate_morphology (*self*, *morph_filename*)

Load a neuroLucida or swc-format cell morphology file.

Parameters

morph_filename [string] Path to morphology.

load_cell_parameters (*self*)

Configure a neuron after the cell morphology has been loaded.

class `allensdk.model.biophysical.utils.Utils` (*description*)

Bases: `allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils`

A helper class for NEURON functionality needed for biophysical simulations.

Attributes

h [object] The NEURON hoc object.

nrn [object] The NEURON python object.

neuron [module] The NEURON module.

generate_morphology (*self*, *morph_filename*)

Load a swc-format cell morphology file.

Parameters

morph_filename [string] Path to swc.

get_recorded_data (*self*, *vec*)

Extract recorded voltages and timestamps given the recorded Vector instance. If `self.stimulus_sampling_rate` is smaller than `self.simulation_sampling_rate`, resample to `self.stimulus_sampling_rate`.

Parameters

vec [neuron.Vector] constructed by `self.record_values`

Returns

dict with two keys: 'v' = numpy.ndarray with voltages, 't' = numpy.ndarray with timestamps

load_cell_parameters (*self*)

Configure a neuron after the cell morphology has been loaded.

static nearest_neuron_sampling_rate (*hz*, *target_hz=40000*)

read_stimulus (*self*, *stimulus_path*, *sweep=0*)

Load current values for a specific experiment sweep and setup simulation and stimulus sampling rates.

NOTE: NEURON only allows simulation timestamps of multiples of 40KHz. To avoid aliasing, we set the simulation sampling rate to the least common multiple of the stimulus sampling rate and 40KHz.

Parameters

stimulus_path [string] NWB file name

sweep [integer, optional] sweep index

record_values (*self*)

Set up output voltage recording.

setup_iclamp (*self*, *stimulus_path*, *sweep=0*)

Assign a current waveform as input stimulus.

Parameters

stimulus_path [string] NWB file name

update_default_cell_hoc (*self*, *description*, *default_cell_hoc='cell.hoc'*)

replace the default 'cell.hoc' path in the manifest with 'cell.hoc' packaged within AllenSDK if it does not exist

`allensdk.model.biophysical.utils.create_utils` (*description*, *model_type=None*)

Factory method to create a Utils subclass.

Parameters

description [Config instance] used to initialize Utils subclass

model_type [string] Must be one of [PERISOMATIC_TYPE, ALL_ACTIVE_TYPE]. If none, defaults to PERISOMATIC_TYPE

Returns

Utils instance

Module contents**allensdk.model.glif package****Submodules****allensdk.model.glif.glif_neuron module**

exception `allensdk.model.glif.glif_neuron.GlifBadResetException` (*message*, *dv*)

Bases: `Exception`

Exception raised when voltage is still above threshold after a reset rule is applied.

```
class allensdk.model.glif.glif_neuron.GlifNeuron (El, dt, asc_tau_array, R_input, C,  

asc_amp_array, spike_cut_length,  

th_inf, th_adapt, coeffs, AScurrent_dynamics_method, voltage_dynamics_method,  

old_dynamics_method, threshold_dynamics_method,  

current_reset_method, voltage_reset_method,  

old_reset_method, init_voltage,  

init_threshold, init_AScurrents,  

**kwargs)
```

Bases: object

Implements the current-based Mihalas Neiber GLIF neuron. Simulations model the voltage, threshold, and afterspike currents of a neuron given an input stimulus. A set of modular dynamics rules are applied until voltage crosses threshold, at which point a set of modular reset rules are applied. See `glif_neuron_methods.py` for a list of what options there are for voltage, threshold, and afterspike current dynamics and reset rules.

Parameters

El [float]
resting potential

dt [float] duration between time steps

asc_tau_array: np.ndarray TODO

R_input [float] input resistance

C [float] capacitance

asc_amp_array [np.ndarray] afterspike current vector. one element per element of `asc_tau_array`.

spike_cut_length [int] how many time steps to replace with NaNs when a spike occurs.

th_inf [float] instantaneous threshold

coeffs [dict] dictionary coefficients premultiplied to neuron properties during simulation. used for optimization.

AScurrent_dynamics_method [dict] dictionary containing the ‘name’ of the afterspike current dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

voltage_dynamics_method [dict] dictionary containing the ‘name’ of the voltage dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

threshold_dynamics_method [dict] dictionary containing the ‘name’ of the threshold dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

AScurrent_reset_method [dict] dictionary containing the ‘name’ of the afterspike current dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

voltage_reset_method [dict] dictionary containing the ‘name’ of the voltage dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

threshold_reset_method [dict] dictionary containing the ‘name’ of the threshold dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

init_voltage [float] initial voltage value

init_threshold [float] initial spike threshold value

init_AScurrents [np.ndarray] initial afterspike current vector. one element per element of asc_tau_array.

TYPE = 'GLIF'

append_threshold_components (*self, spike, voltage*)

static configure_library_method (*method_type, params*)

Create a GlifNeuronMethod instance out of a library of functions organized by type name. This refers to the METHOD_LIBRARY in glif_neuron_methods.py, which lays out the available functions that can be used for dynamics and reset rules.

Parameters

method_type [string] the name of a function category (e.g. 'AScurrent_dynamics_method' for the afterspike current dynamics methods)

params [dict] a dictionary with two members. 'name': the string name of function you want, and 'params': parameters you want to pass to that function

Returns

GlifNeuronMethod a GlifNeuronMethod instance

static configure_method (*method_name, method, method_params*)

Create a GlifNeuronMethod instance given a name, a function, and function parameters. This is just a shortcut to the GlifNeuronMethod constructor.

Parameters

method_name [string] name for referring to this method later

method [function] a python function

method_parameters [dict] function arguments whose values should be fixed

Returns

GlifNeuronMethod a GlifNeuronMethod instance

dynamics (*self, voltage_t0, threshold_t0, AScurrents_t0, inj, time_step, spike_time_steps*)

Update the voltage, threshold, and afterspike currents of the neuron for a single time step.

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

inj [float] the current value of the current injection into the neuron

time_step [int] the current time step of the neuron simulation

spike_time_steps [list] a list of all of the time steps of spikes in the neuron

Returns

tuple voltage_t1 (voltage at next time step), threshold_t1 (threshold at next time step), AScurrents_t1 (afterspike currents at next time step)

classmethod from_dict (*d*)

reset (*self*, *voltage_t0*, *threshold_t0*, *AScurrents_t0*)

Apply reset rules to the neuron's voltage, threshold, and afterspike currents assuming a spike has occurred (voltage is above threshold).

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

Returns

tuple *voltage_t1* (voltage at next time step), *threshold_t1* (threshold at next time step), *AScurrents_t1* (afterspike currents at next time step)

run (*self*, *stim*)

Run neuron simulation over a given stimulus. This steps through the stimulus applying dynamics equations. After each step it checks if voltage is above threshold. If so, *self.spike_cut_length* NaNs are inserted into the output voltages, reset rules are applied to the voltage, threshold, and afterspike currents, and the simulation resumes.

Parameters

stim [np.ndarray] vector of scalar current values

Returns

dict

a dictionary containing: 'voltage': simulated voltage values, 'threshold': threshold values during the simulation, 'AScurrents': afterspike current values during the simulation, 'grid_spike_times': spike times (in units of *self.dt*) aligned to simulation time steps, 'interpolated_spike_times': spike times (in units of *self.dt*) linearly interpolated between time steps, 'spike_time_steps': the indices of grid spike times, 'interpolated_spike_voltage': voltage of the simulation at interpolated spike times, 'interpolated_spike_threshold': threshold of the simulation at interpolated spike times

tau_m

to_dict (*self*)

Convert the neuron to a serializable dictionary.

`allensdk.model.glif.glif_neuron.interpolate_spike_time` (*dt*, *time_step*, *threshold_t0*, *threshold_t1*, *voltage_t0*, *voltage_t1*)

Given two voltage and threshold values, the *dt* between them and the initial time step, interpolate a spike time within the *dt* interval by intersecting the two lines.

`allensdk.model.glif.glif_neuron.interpolate_spike_value` (*dt*, *interpolated_spike_time_offset*, *v0*, *v1*)

Take a value at two adjacent time steps and linearly interpolate what the value would be at an offset between the two time steps.

`allensdk.model.glif.glif_neuron.line_crossing_x` (*dx*, *a0*, *a1*, *b0*, *b1*)

Find the *x* value of the intersection of two lines.

`allensdk.model.glif.glif_neuron.line_crossing_y` (*dx*, *a0*, *a1*, *b0*, *b1*)

Find the *y* value of the intersection of two lines.

allensdk.model.glif.glif_neuron_methods module

The methods in this module are used for configuring dynamics and reset rules for the GlifNeuron. For more details on how to use these methods, see [Generalized LIF Models](#).

```
class allensdk.model.glif.glif_neuron_methods.GlifNeuronMethod(method_name,  
                                                             method,  
                                                             method_params)
```

Bases: object

A simple class to keep track of the name and parameters associated with a neuron method. This class is initialized with a name, function, and parameters to pass to the function. The function then has those passed parameters fixed to a partial function using `functools.partial`. This class then mimics a function itself using the `__call__` convention. Parameters that are not fixed in this way are assumed to be passed into the method when it is called. If the passed parameters contain an argument that is not part of the function signature, an exception will be raised.

Parameters

method_name [string] A shorthand name that will be used to reference this method in the *GlifNeuron*.

method [function] A python function to be called when this instance is called.

method_params [dict] A dictionary mapping function arguments to values for values that should be fixed.

```
modify_parameter(self, param, operator)
```

Modify a function parameter needs to be modified after initialization.

Parameters

param [string] the name of the parameter to modify

operator [callable] a function or lambda that returns the desired modified value

Returns

type the new value of the variable that was just modified.

```
to_dict(self)
```

```
allensdk.model.glif.glif_neuron_methods.dynamics_ACurrent_exp(neuron, AS-  
                                                                currents_t0,  
                                                                time_step,  
                                                                spike_time_steps)
```

Exponential afterspike current dynamics method takes a current at t0 and returns the current at a time step later.

```
allensdk.model.glif.glif_neuron_methods.dynamics_ACurrent_none(neuron, AS-  
                                                                currents_t0,  
                                                                time_step,  
                                                                spike_time_steps)
```

This method always returns zeros for the afterspike currents, regardless of input.

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_inf(neuron, thresh-  
                                                                old_t0,    volt-  
                                                                age_t0,    AS-  
                                                                currents_t0,  
                                                                inj)
```

Set threshold to the neuron's instantaneous threshold.

Parameters

neuron [class]
threshold_t0 [not used here]
voltage_t0 [not used here]
AScurrents_t0 [not used here]
inj [not used here]
AScurrents_t0 [not used here]
inj [not used here]

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_spike_component(neuron,
                                                                           thresh-
                                                                           old_t0,
                                                                           volt-
                                                                           age_t0,
                                                                           AS-
                                                                           cur-
                                                                           rents_t0,
                                                                           inj,
                                                                           a_spike,
                                                                           b_spike,
                                                                           a_voltage,
                                                                           b_voltage)
```

Analytical solution for spike component of threshold. The threshold will adapt via a component initiated by a spike which decays as an exponential. The component is in reference to threshold infinity and are recorded in the neuron's threshold components. The voltage component of the threshold is set to zero in the threshold components because it is zero here. The third component refers to `th_inf` which is added separately as opposed to being included in the voltage component of the threshold as is done in equation 2.1 of Mihalas and Nieber 2009. Threshold infinity is removed for simple optimization.

Parameters

neuron [class]
threshold_t0 [float] threshold input to function
voltage_t0 [float] voltage input to function
AScurrents_t0 [vector] values of after spike currents
inj [float] current injected into the neuron

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_three_components_exact(neuron,
                                                                           thresh-
                                                                           old_t0,
                                                                           volt-
                                                                           age_t0,
                                                                           AS-
                                                                           cur-
                                                                           rents_t0,
                                                                           inj,
                                                                           a_spike,
                                                                           b_spike,
                                                                           a_voltage,
                                                                           b_voltage)
```

Analytical solution for threshold dynamics. The threshold will adapt via two mechanisms: 1. a voltage dependent adaptation. 2. a component initiated by a spike which decays as an exponential. These two component are

in reference to threshold infinity and are recorded in the neuron's threshold components. The third component refers to `th_inf` which is added separately as opposed to being included in the voltage component of the threshold as is done in equation 2.1 of Mihalas and Nieber 2009. Threshold infinity is removed for simple optimization.

Parameters

neuron [class]
threshold_t0 [float] threshold input to function
voltage_t0 [float] voltage input to function
AScurrents_t0 [vector] values of after spike currents
inj [float] current injected into the neuron

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_exact(neuron,  
                                                                       volt-  
                                                                       age_t0,  
                                                                       AS-  
                                                                       cur-  
                                                                       rents_t0,  
                                                                       inj)
```

(TODO) Linear voltage dynamics.

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_forward_euler(neuron,  
                                                                               volt-  
                                                                               age_t0,  
                                                                               AS-  
                                                                               cur-  
                                                                               rents_t0,  
                                                                               inj)
```

(TODO) Linear voltage dynamics.

```
allensdk.model.glif.glif_neuron_methods.max_of_line_and_const(x, b, c, d)  
Find the maximum of a value and a position on a line
```

Parameters

x: float x position on line 1
c: float slope of line 1
d: float y-intercept of line 1
b: float y-intercept of line 2

Returns

float the max of a line value and a constant

```
allensdk.model.glif.glif_neuron_methods.min_of_line_and_zero(x, c, d)  
Find the minimum of a value and a position on a line
```

Parameters

x: float x position on line 1
c: float slope of line 1
d: float y-intercept of line 1
b: float y-intercept of line 2

Returns

float the max of a line value and a constant

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_none(neuron, AScurrents_t0)
```

Reset afterspike currents to zero.

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_sum(neuron, AScurrents_t0, r)
```

Reset afterspike currents by adding summed exponentials. Left over currents from last spikes as well as newly initiated currents from current spike. Currents amplitudes in `neuron.asc_amp_array` need to be the amplitudes advanced though the spike cutting. I.e. In the preprocessor if the after spike currents are calculated via the GLM from spike initiation the amplitude at the time after the spike cutting needs to be calculated and `neuron.asc_amp_array` needs to be set to this value.

Parameters

r [np.ndarray] a coefficient vector applied to the afterspike currents

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_inf(neuron, threshold_old_t0, voltage_v1)
```

Reset the threshold to instantaneous threshold.

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_three_components(neuron, threshold_old_t0, voltage_v1, a_spike, b_spike)
```

This method calculates the two components of the threshold: a spike (fast) component and a voltage (slow) component. The `threshold_components` vectors are then updated so that the traces match the voltage, current, and total threshold traces. The spike component of the threshold decays via an exponential fit specified by the amplitude `a_spike` and the time constant `b_spike` fit via the multiblip data. The voltage component does not change during the duration of the spike. The spike component are threshold component are summed along with threshold infinity to return the total threshold. Note that in the current implementation `a_spike` is added to the last value of the `threshold_components` which means that `a_spike` is the amplitude after spike cutting (if there is any).

Inputs:

neuron: class contains attributes of the neuron

threshold_t0, voltage_t0: float are not used but are here for consistency with other methods

a_spike: float amplitude of the exponential decay of spike component of threshold after spike cutting has been implemented.

b_spike: float amplitude of the exponential decay of spike component of threshold

Outputs:

Returns: float the total threshold which is the sum of the spike component of threshold, the voltage component of threshold and threshold infinity (with it's corresponding coefficient)

neuron.threshold_components: dictionary containing

a_spike: list vector of spiking component of threshold that corresponds to the voltage, current, and total threshold traces

b_spike: list

vector of voltage component of threshold that corresponds to the voltage, current, and total threshold traces.

Note that this function can be changed to use `a_spike` at the time of the spike and then have the the spike component plus the residual decay thought the spike. There are benefits and drawbacks to this. This potential change would be beneficial as it perhaps makes more biological sense for the threshold to go up at the time of spike if the traces are ever used. Also this would mean that `a_spike` would not have to be adjusted thought the spike cutting after the multiblip fit. However the current implementation makes sense in that it is similar to how afterspike currents are implemented.

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_v_before(neuron, voltage_t0, a, b)
```

Reset voltage to the previous value with a scale and offset applied.

Parameters

a [float] voltage scale constant

b [float] voltage offset constant

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_zero(neuron, voltage_t0)
```

Reset voltage to zero.

```
allensdk.model.glif.glif_neuron_methods.spike_component_of_threshold_exact(th0, b_spike, t)
```

Spike component of threshold modeled as an exponential decay. Implemented here as exact analytical solution.

Parameters

th0 [float] threshold input to function

b_spike [float] decay constant of exponential

t [float or array] time step if used in an Euler setup time if used analytically

```
allensdk.model.glif.glif_neuron_methods.spike_component_of_threshold_forward_euler(th_t0, b_spike, dt)
```

Spike component of threshold modeled as an exponential decay. Implemented here for forward Euler

Parameters

th_t0 [float] threshold input to function

b_spike [float] decay constant of exponential

dt [float] time step

```
allensdk.model.glif.glif_neuron_methods.voltage_component_of_threshold_exact(th0, v0, I, t, a_voltage, b_voltage, C, g, EI)
```

Note this function is the exact formulation; however, `dt` is used because `t0` is the initial time and `dt` is the time the function is exactly evaluated at. Note: that here, this equation is in reference to `th_inf`. Therefore `th0` is the total threshold-`thr_inf` (`threshold_inf` replaced with 0 in the equation to be verbose). This is done so that `th_inf` can be optimized without affecting this function.

Parameters

th0 [float] threshold input to function

v0 [float] voltage input to function

I [float] total current entering neuron (note if there are after spike currents these must be included in this value)

t [float or array] time step if used in an Euler setup time if used analytically

a_voltage [float] constant a

b_voltage [float] constant b

C [float] capacitance

g [float] conductance (1/resistance)

El [float] reversal potential

```
allensdk.model.glif.glif_neuron_methods.voltage_component_of_threshold_forward_euler(th_t0,
                                                                                       v_t0,
                                                                                       dt,
                                                                                       a_voltage,
                                                                                       b_voltage,
                                                                                       El)
```

Equation 2.1 of Mihalas and Nieber, 2009 implemented for use in forward Euler. Note here all variables are in reference to threshold infinity. Therefore thr_inf is zero here (replaced threshold_inf with 0 in the equation to be verbose). This is done so that th_inf can be optimized without affecting this function.

Parameters

th_t0 [float] threshold input to function

v_t0 [float] voltage input to function

dt [float] time step

a_voltage [float] constant a

b_voltage [float] constant b

El [float] reversal potential

allensdk.model.glif.simulate_neuron module

```
allensdk.model.glif.simulate_neuron.load_sweep(file_name, sweep_number)
```

Load the stimulus for a sweep from file.

```
allensdk.model.glif.simulate_neuron.main()
```

```
allensdk.model.glif.simulate_neuron.parse_arguments()
```

Use argparse to get required arguments from the command line

```
allensdk.model.glif.simulate_neuron.simulate_neuron(neuron, sweep_numbers,
                                                    input_file_name, out-
                                                    put_file_name, spike_cut_value)
```

```
allensdk.model.glif.simulate_neuron.simulate_sweep(neuron, stimulus,
                                                    spike_cut_value)
```

Simulate a neuron given a stimulus and initial conditions.

```
allensdk.model.glif.simulate_neuron.simulate_sweep_from_file(neuron,  
                                                             sweep_number,  
                                                             input_file_name,  
                                                             output_file_name,  
                                                             spike_cut_value)
```

Load a sweep stimulus, simulate the response, and write it out.

```
allensdk.model.glif.simulate_neuron.write_sweep_response(file_name,  
                                                         sweep_number, response,  
                                                         spike_times)
```

Overwrite the response in a file.

Module contents

A Generalized Linear Integrate and Fire (GLIF) neuron modeling package. Use this code to run the GLIF models available in the Allen Cell Types Atlas. See [Generalized LIF Models](#) for more details.

Module contents

6.1.8 allensdk.morphology package

Submodules

allensdk.morphology.validate_swc module

```
allensdk.morphology.validate_swc.main()  
allensdk.morphology.validate_swc.validate_swc(swc_file)
```

To be compatible with NEURON, SWC files must have the following properties:

- 1) a single root node with parent ID '-1'
- 2) sequentially increasing ID numbers
- 3) immediate children of the soma cannot branch

Module contents

6.1.9 allensdk.mouse_connectivity package

Subpackages

allensdk.mouse_connectivity.grid package

Subpackages

allensdk.mouse_connectivity.grid.subimage package

Submodules

allensdk.mouse_connectivity.grid.subimage.base_subimage module

```

class allensdk.mouse_connectivity.grid.subimage.base_subimage.IntensitySubImage (reduce_level,
                                          in_dims,
                                          in_spacing,
                                          coarse_spacing,
                                          in-
                                          ten-
                                          sity_paths,
                                          *args,
                                          **kwargs)

    Bases: allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage

    get_intensity (self)

    required_intensities = []

    setup_images (self)

class allensdk.mouse_connectivity.grid.subimage.base_subimage.PolygonSubImage (reduce_level,
                                          in_dims,
                                          in_spacing,
                                          coarse_spacing,
                                          poly-
                                          gon_info,
                                          *args,
                                          **kwargs)

    Bases: allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage

    get_polygons (self)

    optional_polys = []

    required_polys = []

    setup_images (self)

class allensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentationSubImage (reduce_level,
                                          in_dims,
                                          in_spacing,
                                          coarse_spac-
                                          seg-
                                          men-
                                          ta-
                                          tion_paths,
                                          *args,
                                          **kwargs)

    Bases: allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage

    extract_injection_from_segmentation (self, segmentation_name='segmentation', injec-
                                          tion_name='injection')

```

Notes

Currently, the segmentation uses a series of codes to map 8-bit values onto meaningful classifications. The code for signal pixels is a 1 in at least one of the 5 rightmost bits.

```

extract_signal_from_segmentation (self, segmentation_name='segmentation', sig-
                                  nal_name='signal')

```

Notes

Currently, the segmentation uses a series of codes to map 8-bit values onto meaningful classifications. The code for signal pixels is a 1 in the leftmost bit.

In some cases, bit 5 indicates that the pixel was not removed in a posfiltering process. Optionally, this postfilter can be applied in gridding.

```
get_segmentation (self)
```

```
process_segmentation (self)
```

```
read_segmentation_image (self, segmentation_name='segmentation')
```

Notes

We downsample in memory rather than using the jp2 pyramid because the segmentation is a label image.

```
required_segmentations = []
```

```
setup_images (self)
```

```
class allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage (reduce_level,  
                                                                    in_dims,  
                                                                    in_spacing,  
                                                                    coarse_spacing,  
                                                                    *args,  
                                                                    **kwargs)
```

Bases: object

```
apply_mask (self, image_name, mask_name, positive=True)
```

```
apply_pixel_counter (self, accumulator_name, image)
```

```
binarize (self, image_name)
```

```
compute_coarse_planes (self)
```

```
make_pixel_counter (self)
```

```
pixel_counter
```

```
setup_images (self)
```

```
allensdk.mouse_connectivity.grid.subimage.base_subimage.run_subimage (input_data)
```

allensdk.mouse_connectivity.grid.subimage.cav_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage (reduce_level,  
                                                                    in_dims,  
                                                                    in_spacing,  
                                                                    coarse_spacing,  
                                                                    polygon_info,  
                                                                    *args,  
                                                                    **kwargs)
```

Bases: `allensdk.mouse_connectivity.grid.subimage.base_subimage.`

`PolygonSubImage`

```
compute_coarse_planes (self)
```



```
required_polys = ['missing_tile', 'cav_tracer']
```

allensdk.mouse_connectivity.grid.subimage.classic_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage(reduce_level,
in_dims,
in_spacing,
coarse_spacing,
poly-
gon_info,
seg-
men-
ta-
tion_paths,
in-
ten-
sity_paths,
in-
jec-
tion_polygon,
*args,
**kwargs)
```

Bases: `allensdk.mouse_connectivity.grid.subimage.base_subimage.IntensitySubImage`, `allensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentationSubImage`, `allensdk.mouse_connectivity.grid.subimage.base_subimage.PolygonSubImage`

```
compute_coarse_planes(self)
```

```
compute_injection(self)
```

```
compute_intensity(self)
```

```
compute_projection(self)
```

```
compute_sum_pixels(self)
```

```
optional_polys = ['aav_tracer']
```

```
process_segmentation(self)
```

```
required_intensities = ['green']
```

```
required_polys = ['missing_tile', 'no_signal', 'aav_exclusion']
```

```
required_segmentations = ['segmentation']
```

allensdk.mouse_connectivity.grid.subimage.count_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage(reduce_level,  
                                                                           in_dims,  
                                                                           in_spacing,  
                                                                           coarse_spacing,  
                                                                           poly-  
                                                                           gon_info,  
                                                                           seg-  
                                                                           men-  
                                                                           ta-  
                                                                           tion_paths,  
                                                                           in-  
                                                                           jec-  
                                                                           tion_polygon_key=''  
                                                                           *args,  
                                                                           **kwargs)
```

```
Bases:      allensdk.mouse_connectivity.grid.subimage.base_subimage.  
            SegmentationSubImage,      allensdk.mouse_connectivity.grid.subimage.  
            base_subimage.PolygonSubImage
```

```
compute_coarse_planes(self)  
compute_injection(self)  
compute_projection(self)  
compute_sum_pixels(self)  
process_segmentation(self)  
required_polys = ['missing_tile', 'no_signal', 'aav_exclusion']  
required_segmentations = ['segmentation']
```

Module contents

```
allensdk.mouse_connectivity.grid.subimage.run_subimage(input_data)
```

allensdk.mouse_connectivity.grid.utilities package

Submodules

allensdk.mouse_connectivity.grid.utilities.downsampling_utilities module

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.apply_divisions(image,  
                                                                           win-  
                                                                           dow_size)  
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.block_average(volume,  
                                                                           fac-  
                                                                           tor)
```

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.conv(image,
                                                                    fac-
                                                                    tor,
                                                                    win-
                                                                    dow_size)
```

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.downsample_average(volume,
                                                                                      cur-
                                                                                      rent_spa-
                                                                                      tar-
                                                                                      get_spa-
```

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.extract(image,
                                                                    fac-
                                                                    tor,
                                                                    win-
                                                                    dow_size,
                                                                    win-
                                                                    dow_step,
                                                                    out-
                                                                    put_shape)
```

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.window_average(volume,
                                                                                      fac-
                                                                                      tor)
```

allensdk.mouse_connectivity.grid.utilities.image_utilities module

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.block_apply(in_image,
                                                                    out_shape,
                                                                    dtype,
                                                                    blocks,
                                                                    fn)
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.build_affine_transform(aff_params)
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.build_composite_transform(dfmfield,
                                                                                      aff_para-
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.compute_coarse_parameters(in_dims,
                                                                                      in_spa-
                                                                                      ce,
                                                                                      out_spa-
                                                                                      ce,
                                                                                      re-
                                                                                      duce_le-
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.grid_image_blocks(in_shape,
                                                                    in_spacing,
                                                                    out_spacing)
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.image_from_array(array,
                                                                    spac-
                                                                    ing,
                                                                    ori-
                                                                    gin=True)
```

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.new_image(dims,  
                                                                    spacing,  
                                                                    dtype,  
                                                                    origin=True)  
  
allensdk.mouse_connectivity.grid.utilities.image_utilities.np_sitk_convert(np_type)  
allensdk.mouse_connectivity.grid.utilities.image_utilities.rasterize_polygons(shape,  
                                                                    scale,  
                                                                    polys)  
  
allensdk.mouse_connectivity.grid.utilities.image_utilities.read_intensity_image(path)  
allensdk.mouse_connectivity.grid.utilities.image_utilities.read_segmentation_image(path)  
allensdk.mouse_connectivity.grid.utilities.image_utilities.resample_into_volume(image,  
                                                                    transform,  
                                                                    z,  
                                                                    vol,  
                                                                    dtype=8)  
  
allensdk.mouse_connectivity.grid.utilities.image_utilities.resample_volume(volume,  
                                                                    dims,  
                                                                    spacing,  
                                                                    interpolator=None,  
                                                                    transform=None)  
  
allensdk.mouse_connectivity.grid.utilities.image_utilities.set_image_spacing(image,  
                                                                    spacing,  
                                                                    origin=True)  
  
allensdk.mouse_connectivity.grid.utilities.image_utilities.sitk_np_convert(sitk_type)  
allensdk.mouse_connectivity.grid.utilities.image_utilities.write_volume(volume,  
                                                                    name,  
                                                                    prefix=None,  
                                                                    specify_resolution=None,  
                                                                    extension='.nrrd',  
                                                                    paths=None)
```

Module contents

allensdk.mouse_connectivity.grid.writers package

Module contents

`allensdk.mouse_connectivity.grid.writers.cav_writer`(*gridder, grid_prefix, accumulator_prefix, **kwargs*)

`allensdk.mouse_connectivity.grid.writers.classic_writer`(*gridder, grid_prefix, accumulator_prefix, target_spacings, **kwargs*)

`allensdk.mouse_connectivity.grid.writers.count_writer`(*gridder, grid_prefix, accumulator_prefix, target_spacings, **kwargs*)

`allensdk.mouse_connectivity.grid.writers.handle_pyramid`(*isg, key, target_spacings, prefix, paths*)

`allensdk.mouse_connectivity.grid.writers.ratio_and_pyramid`(*isg, num, den, out, accumulator_prefix, grid_prefix, target_spacings, paths*)

Submodules

allensdk.mouse_connectivity.grid.image_series_gridder module

class `allensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`(*in_dims, in_spacing, out_dims, out_spacing, reduce_level, subimages, subimage_kwargs, nprocesses, affine_params, dfm fld_path*)

Bases: `object`

accumulator_to_numpy(*self, key, cb*)

build_coarse_grids(*self*)

consume_volume(*self, key, cb*)

initialize_coarse_volume(*self, key, dtype*)

make_ratio_volume(*self, num_key, den_key, ratio_key*)
assume parents numpified

paste_slice(*self, key, index, slice_array*)

```
paste_subimage (self, index, output)  
    Inserts planar accumulators into coarse grid volumes  
resample_volume (self, key)  
set_coarse_grid_parameters (self)  
setup_subimages (self)  
transform
```

Module contents

Module contents

6.1.10 allensdk.test_utilities package

Submodules

allensdk.test_utilities.custom_comparators module

```
class allensdk.test_utilities.custom_comparators.WhitespaceStrippedString (string:  
                                                                    str,  
                                                                    whites-  
                                                                    pace_chars:  
                                                                    str  
                                                                    =  
                                                                    '\s',  
                                                                    ASCII:  
                                                                    bool  
                                                                    =  
                                                                    False)
```

Bases: object

Comparator class to compare strings that have been stripped of whitespace. By default removes any unicode whitespace character that matches the regex `s`, (which includes `[` `],` and other unicode whitespace characters).

allensdk.test_utilities.regression_fixture module

```
allensdk.test_utilities.regression_fixture.get_list_of_path_dict()
```

allensdk.test_utilities.temp_dir module

```
allensdk.test_utilities.temp_dir.temp_dir (request)
```

Module contents

6.2 Submodules

6.2.1 allensdk.deprecated module

`allensdk.deprecated.class_deprecated(message=None)`

`allensdk.deprecated.deprecated(message=None)`

`allensdk.deprecated.legacy(message=None)`

6.3 Module contents

exception `allensdk.OneResultExpectedError`

Bases: `RuntimeError`

`allensdk.one(x)`

The Allen Software Development Kit houses source code for reading and processing Allen Brain Atlas data. The Allen SDK focuses on the Allen Brain Observatory, Cell Types Database, and Mouse Brain Connectivity Atlas.

Attention: As of October 2019, we have dropped Python 2 support and any files with a py2 dependency (for example analysis files) have been updated.



Allen Brain Observatory


The Allen Brain Observatory is a collection of data resources for understanding sensory processing in the mouse visual cortex. These resources systematically measure visual responses in multiple cortical areas and layers using two-photon calcium imaging or high-density extracellular electrophysiology (Neuropixels) probes. Recordings are performed on mice passively viewing visual stimuli or trained to actively perform an image change detection task.

Behavior	Modality	Resource	Initial Release
Passive	Optical physiology	Visual Coding - Optical Physiology	June 2016
Passive	Extracellular electrophysiology	Visual Coding - Neuropixels	October 2019
Active	Optical physiology	Visual Behavior - Optical Physiology	March 2021
Active	Extracellular electrophysiology	Visual Behavior - Neuropixels	Coming soon...

Experiment and stimulus data are provided in [Neurodata Without Borders](#) (NWB) files. The AllenSDK provides code to:

- download and organize experiment data according to cortical area, imaging depth, and Cre line
- access experiment metadata and data streams
- transform and analyze data

More information about each study is provided in the linked pages. A web-based entry point to the Visual Coding – Optical physiology data is available at <http://observatory.brain-map.org/visualcoding>.

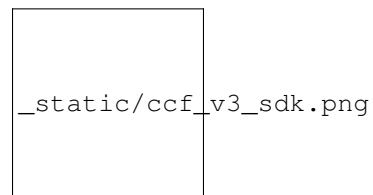

Allen Cell Types Database

The [Allen Cell Types Database](#) contains electrophysiological and morphological characterizations of individual neurons in the mouse primary visual cortex. The Allen SDK provides Python code for accessing electrophysiology measurements ([NWB files](#)) for all neurons and morphological reconstructions ([SWC files](#)) for a subset of neurons.

The Database also contains two classes of models fit to this data set: biophysical models produced using the NEURON simulator and generalized leaky integrate and fire models (GLIFs) produced using custom Python code provided with this toolkit.

The Allen SDK provides sample code demonstrating how to download neuronal model parameters from the Allen Brain Atlas API and run your own simulations using stimuli from the Allen Cell Types Database or custom current injections:

- *Biophysical Models*
- *Generalized LIF Models*



Allen Mouse Brain Connectivity Atlas

The [Allen Mouse Brain Connectivity Atlas](#) is a high-resolution map of neural connections in the mouse brain. Built on an array of transgenic mice genetically engineered to target specific cell types, the Atlas comprises a unique compendium of projections from selected neuronal populations throughout the brain. The primary data of the Atlas consists of high-resolution images of axonal projections targeting different anatomic regions or various cell types using Cre-dependent specimens. Each data set is processed through an informatics data analysis pipeline to obtain spatially mapped quantified projection information.

The Allen SDK provides Python code for accessing experimental metadata along with projection signal volumes registered to a common coordinate framework. This framework has structural annotations, which allows users to compute structure-level signal statistics.

See the [mouse connectivity section](#) for more details.

CHAPTER 10

What's New - 2.11.0

- Refactored the schema for the Ecephys copy utility to avoid raising an error when a previous output file already exists.

CHAPTER 11

What's New - 2.11.0

- python 3.8 compatibility
- CloudCache (the class supporting cloud-based data releases) is now smart enough to construct symlinks between files that are identical across dataset versions (rather than downloading duplicate copies of files).
- VisualBehaviorOphysProjectCache supports user-controlled switching between dataset versions.

CHAPTER 12

What's New - 2.10.3

- Adds restriction to require hdmf version to be strictly less than 2.5.0 which accidentally introduced a major version breaking change

CHAPTER 13

What's New - 2.10.2

- This version marks the release of Visual Behavior Optical Physiology data! For more details please visit the: [Visual Behavior - Optical Physiology Project Page](#)
- Update documentation to support visual behavior data release
- Fixes a bug with the dictionary returned by BehaviorSession get get_performance_metrics() method
- Adds docstrings to the BehaviorSession get_performance_metrics(), get_rolling_performance_df(), and get_reward_rate() methods

CHAPTER 14

What's New - 2.10.1

- Changes name of BehaviorProjectCache to VisualBehaviorOphysProjectCache
- Changes VisualBehaviorOphysProjectCache method `get_session_table()` to `get_ophys_session_table()`
- Changes VisualBehaviorOphysProjectCache method `get_experiment_table()` to `get_ophys_experiment_table()`
- VisualBehaviorOphysProjectCache is enabled to instantiate `from_s3_cache()` and `from_local_cache()`
- Improvements to BehaviorProjectCache
- Adds project metadata writer

CHAPTER 15

What's New - 2.9.0

- Updates to Session metadata; refactors implementation to use class rather than dict internally
- Fixes a bug that was preventing Allen Institute Windows users from accessing gratings images

CHAPTER 16

What's New - 2.8.0

- Created lookup table to get monitor_delay for cases where calculation from data fails
- If sync timestamp file has more timestamps than eye tracking moving has frame, trim excess timestamps (up to 15)
- Session API returns both warped and unwarped stimulus images, and both are written to NWB

CHAPTER 17

What's New - 2.7.0

- Refactored behavior and ophys session and data APIs to remove a circular inheritance issue
- Fixed segmentation mask and roi_mask misregistration in 'BehaviorOphysSession'
- Replaces BehaviorOphysSession.get_roi_masks() method with roi_masks property
- Fixes bug which prevented the SDK from loading stimuli dataframes for static gratings
- Return event detection data through session API
- Read/write event detection data from/to NWB
- Time stamps for events in trial_log are set to the exact sync timestamp of the corresponding frame.
- For behavior-only sessions, sync-like timestamp of the first frame is set to zero.
- Refactored BehaviorOphysSession to inherit methods and properties from BehaviorSession
- Fixed a test for checking that Behavior and BehaviorOphysSessions contain the same data regardless of which API (LIMS/JSON/NWB) is used. Also fixed resulting failure cases.

CHAPTER 18

What's New - 2.6.0

- Adds ability to write and read behavior only experiments
- Adds eye tracking ellipse fits and metadata as new NWB data stream
- OPhys Behavior data retrieval methods no longer depend on ROIs being ordered identically in different files.

CHAPTER 19

What's New - 2.5.0 (January 29, 2021)

- Adds unfiltered running speed as new data stream
- `run_demixing` gracefully ignores any ROIs that are not in the input trace file

What's New - 2.4.0 (December 21, 2020)

As of the 2.4.0 release: - When running `raster_plot` on a `spike_times` dataframe, the spike times from each unit are plotted twice. (thank you @dgmurx) - improvements and fixes to behavior ophys NWB files. - improvements and fixes to BehaviorProjectCache tables including new column “donor_id” - implemented a timeout to obtaining an ecephys session. (thank you @wesley-jones) - big overhaul of how Behavior and BehaviorOphys classes are structured for the visual behavior project. See <https://github.com/AllenInstitute/AllenSDK/pull/1789>

CHAPTER 21

What's New - 2.3.2 (October 19, 2020)

As of the 2.3.2 release:

- (Internal) Fixed a `running_processing` bug for behavior ophys experiments when the input data would have one more encoder entry than timestamp. The behavior of the code now matches what the warning says.

CHAPTER 22

What's New - 2.3.1 (October 13, 2020)

As of the 2.3.1 release:

- (Internal) Fixed a write_nwb bug for behavior ophys experiments involving the BehaviorOphysJsonApi expecting a mesoscope-specific method.

What's New - 2.3.0 (October 9, 2020)

As of the 2.3.0 release:

- Visual behavior running speed is now low-pass filtered at 10Hz. The raw running speed data is still available. The running speed is corrected for encoder threshold crossing artifacts.
- Support for stimulus gratings for visual behavior.
- Fixed an eye-tracking sync problem.
- Updates to some visual behavior pynwb implementations.
- Adds load sync data for individual plane sets to relate accurate event timings to mesoscope data.
- Adds public API method to access the `behavior_session_id` from an instance of `BehaviorOphysSession`.

What's New - 2.2.0 (September 3, 2020)

As of the 2.2.0 release:

- AllenSDK HTTP engine streaming requests now include a progress bar
- *ImportError: cannot import name 'MultiContainerInterface' from 'hdmf.container'* errors should now be resolved (by removing explicit version bounds on the *hdmf* package).
- The optical physiology 2-photon trace demixer has been modified to be more memory friendly and should no longer result in out of memory errors when trying to demix very large movie stacks.

What's New - 2.1.0 (July 16, 2020)

As of the 2.1.0 release:

- behavior ophys nwb files can now be written using updated pynwb and hdmf
- A warning has been added if you are using AllenSDK with outdated NWB files
- A new documentation file has been added which will contain Visual Behavior specific terms for quick lookup

What's New - 2.0.0 (June 11, 2020)

As of the 2.0.0 release:

- pynwb and hdmf version requirements have been made less strict
- The organization of data for ecephys neuropixels Neurodata Without Borders (NWB) files has been significantly changed to conform with NWB specifications and best practices
- CCF locations for ecephys neuropixels electrodes are now written to NWB files
- Examples for accessing eye tracking ellipse fit and screen gaze location data have been added to ecephys example notebooks

Important Note: Due to newer versions of pynwb/hdmf having issues reading previously released Visual Coding Neuropixels NWB files and due to the significant reorganization of their NWB file contents, this release contains breaking changes that necessitate a major version revision. NWB files released prior to 6/11/2020 are not guaranteed to work with the 2.0.0 version of AllenSDK. If you cannot or choose not to re-download the updated NWB files, you can continue using a prior version of AllenSDK (< 2.0.0) to access them. However, no further features or bugfixes for AllenSDK (< 2.0.0) are planned. Data released for other projects (Cell Types, Mouse Connectivity, etc.) are *NOT* affected and will *NOT* need to be re-downloaded

CHAPTER 27

Previous Release Notes

- 1.8.0
- 1.7.1
- 1.7.0
- 1.6.0
- 1.5.0
- 1.4.0
- 1.3.0
- 1.2.0
- 1.1.1
- 1.1.0
- 1.0.2
- 0.16.3
- 0.16.2
- 0.16.1
- 0.16.0
- 0.14.5
- 0.14.4
- 0.14.3
- 0.14.2
- 0.13.2
- 0.13.1
- 0.13.0

- [0.12.4](#)

Bibliography

[1] Allen Brain Atlas Data Portal: [Downloading a WellKnownFile](#).

a

allensdk, 371
allensdk.api, 106
allensdk.api.api, 103
allensdk.api.cloud_cache, 67
allensdk.api.cloud_cache.file_attributes, 65
allensdk.api.cloud_cache.manifest, 66
allensdk.api.cloud_cache.utils, 67
allensdk.api.queries, 97
allensdk.api.queries.annotated_section_data_sets_api, 68
allensdk.api.queries.biophysical_api, 69
allensdk.api.queries.brain_observatory_api, 70
allensdk.api.queries.cell_types_api, 73
allensdk.api.queries.connected_services, 75
allensdk.api.queries.glif_api, 76
allensdk.api.queries.grid_data_api, 77
allensdk.api.queries.image_download_api, 78
allensdk.api.queries.mouse_atlas_api, 81
allensdk.api.queries.mouse_connectivity_api, 82
allensdk.api.queries.ontologies_api, 86
allensdk.api.queries.reference_space_api, 87
allensdk.api.queries.rma_api, 89
allensdk.api.queries.rma_pager, 94
allensdk.api.queries.rma_template, 94
allensdk.api.queries.svg_api, 94
allensdk.api.queries.synchronization_api, 95
allensdk.api.queries.tree_search_api, 96
allensdk.api.warehouse_cache, 103
allensdk.api.warehouse_cache.cache, 97
allensdk.api.warehouse_cache.caching_utilities, 101

b

allensdk.brain_observatory, 221
allensdk.brain_observatory.argschema_utilities, 188
allensdk.brain_observatory.behavior, 156
allensdk.brain_observatory.behavior.behavior_project, 108
allensdk.brain_observatory.behavior.behavior_project_api, 107
allensdk.brain_observatory.behavior.behavior_project_data, 107
allensdk.brain_observatory.behavior.behavior_project_data_api, 107
allensdk.brain_observatory.behavior.behavior_project_data_api, 108
allensdk.brain_observatory.behavior.behavior_project_data_api, 107
allensdk.brain_observatory.behavior.behavior_project_data_api, 107
allensdk.brain_observatory.behavior.behavior_project_data_api, 108
allensdk.brain_observatory.behavior.behavior_project_data_api, 107
allensdk.brain_observatory.behavior.behavior_project_data_api, 135
allensdk.brain_observatory.behavior.dprime, 136
allensdk.brain_observatory.behavior.event_detection, 137
allensdk.brain_observatory.behavior.eye_tracking_project, 138
allensdk.brain_observatory.behavior.image_api, 141
allensdk.brain_observatory.behavior.metadata, 111
allensdk.brain_observatory.behavior.metadata.behavior_project, 109
allensdk.brain_observatory.behavior.metadata.behavior_project, 111
allensdk.brain_observatory.behavior.mtrain, 141

168 allensdk.brain_observatory.roi_masks,
 allensdk.brain_observatory.ecephys.stimulus_table, 203
 172 allensdk.brain_observatory.running_speed,
 allensdk.brain_observatory.ecephys.stimulus_table, 206
 173 allensdk.brain_observatory.output_validation,
 allensdk.brain_observatory.session_analysis,
 allensdk.brain_observatory.ecephys.stimulus_table, 206
 174 allensdk.brain_observatory.stimulus_parameter_extraction,
 allensdk.brain_observatory.session_api_utils,
 allensdk.brain_observatory.ecephys.stimulus_table, 209
 168 allensdk.brain_observatory.visualization,
 allensdk.brain_observatory.static_gratings,
 allensdk.brain_observatory.ecephys.stimulus_table, 210
 167 allensdk.brain_observatory.visualization.view_blocks,
 allensdk.brain_observatory.stimulus_analysis,
 allensdk.brain_observatory.ecephys.visualization, 211
 175 allensdk.brain_observatory.stimulus_info,
 allensdk.brain_observatory.ecephys.write_nwb, 213
 176 allensdk.brain_observatory.sync_dataset,
 allensdk.brain_observatory.extract_running_speed, 217
 177 allensdk.brain_observatory.sync_utilities,
 allensdk.brain_observatory.findlevel, 187
 197 allensdk.brain_observatory.visualization,
 allensdk.brain_observatory.gaze_mapping, 187
 177
 allensdk.brain_observatory.locally_sparse_noise,
 197 allensdk.config, 230
 allensdk.brain_observatory.natural_movie, allensdk.config.app, 224
 198 allensdk.config.app.application_config,
 allensdk.brain_observatory.natural_scenes, 222
 199 allensdk.config.manifest, 227
 allensdk.brain_observatory.observatory_plots, allensdk.config.manifest_builder, 229
 200 allensdk.config.model, 227
 allensdk.brain_observatory.ophys, 178 allensdk.config.model.description, 226
 allensdk.brain_observatory.ophys.trace_extraction, allensdk.config.model.description_parser,
 178 227
 allensdk.brain_observatory.r_neuropil, allensdk.config.model.formats, 226
 202 allensdk.config.model.formats.hdf5_util,
 allensdk.brain_observatory.receptive_field_analysis, 224
 186 allensdk.config.model.formats.json_description_parser,
 allensdk.brain_observatory.receptive_field_analysis, 224
 178 allensdk.config.model.formats.pycfg_description_parser,
 allensdk.brain_observatory.receptive_field_analysis, 225
 181 allensdk.eventdetection,
 allensdk.brain_observatory.receptive_field_analysis, 265
 182 allensdk.core, 231
 allensdk.brain_observatory.receptive_field_analysis, 231
 182 allensdk.core.authentication, 231
 allensdk.brain_observatory.receptive_field_analysis, 232
 182 allensdk.brain_observatory.brain_observatory_nwb_data_set,
 allensdk.brain_observatory.receptive_field_analysis, 232
 183 allensdk.core.postprocessing_utilities,
 allensdk.brain_observatory.receptive_field_analysis, 236
 183 allensdk.core.receptive_field_cache, 236
 allensdk.brain_observatory.receptive_field_analysis, 239
 184 allensdk.core.dat_utilities, 239
 allensdk.brain_observatory.receptive_field_analysis, 239
 184 allensdk.core.exceptions, 239
 allensdk.brain_observatory.receptive_field_analysis, 240
 184 allensdk.core.h5_utilities, 240
 allensdk.brain_observatory.receptive_field_analysis, 240
 184 allensdk.core.json_utilities, 240
 allensdk.brain_observatory.receptive_field_analysis, 231
 185 allensdk.core.lazy_property, 231
 allensdk.brain_observatory.receptive_field_analysis, 235
 185 allensdk.core.visualization_utilities, 235

[allensdk.core.lazy_property.lazy_property, 230](#)
[allensdk.core.mouse_connectivity_cache, 241](#)
[allensdk.core.nwb_data_set, 246](#)
[allensdk.core.obj_utilities, 248](#)
[allensdk.core.ontology, 249](#)
[allensdk.core.ophys_experiment_session_id_mapper, 250](#)
[allensdk.core.reference_space, 250](#)
[allensdk.core.reference_space_cache, 252](#)
[allensdk.core.simple_tree, 254](#)
[allensdk.core.sitk_utilities, 257](#)
[allensdk.core.structure_tree, 258](#)
[allensdk.core.swc, 261](#)
[allensdk.core.typing, 265](#)
d
[allensdk.deprecated, 371](#)
e
[allensdk.ephys, 277](#)
[allensdk.ephys.ephys_extractor, 265](#)
[allensdk.ephys.ephys_features, 269](#)
[allensdk.ephys.extract_cell_features, 276](#)
[allensdk.ephys.feature_extractor, 277](#)
i
[allensdk.internal, 348](#)
[allensdk.internal.api, 282](#)
[allensdk.internal.api.api_prerelease, 282](#)
[allensdk.internal.api.lims_api, 282](#)
[allensdk.internal.api.queries, 282](#)
[allensdk.internal.api.queries.biophysical_module_api, 278](#)
[allensdk.internal.api.queries.biophysical_module_reader, 278](#)
[allensdk.internal.api.queries.grid_data_api_prerelease, 279](#)
[allensdk.internal.api.queries.mouse_connectivity_api_prerelease, 280](#)
[allensdk.internal.api.queries.optimize_config_reader, 281](#)
[allensdk.internal.brain_observatory, 298](#)
[allensdk.internal.brain_observatory.annotated_region_metrics, 284](#)
[allensdk.internal.brain_observatory.demix_report, 285](#)
[allensdk.internal.brain_observatory.demixer, 285](#)
[allensdk.internal.brain_observatory.eye_calibration, 287](#)
[allensdk.internal.brain_observatory.fit_ellipse, 290](#)
[allensdk.internal.brain_observatory.frame_stream, 290](#)
[allensdk.internal.brain_observatory.itracker_utils, 291](#)
[allensdk.internal.brain_observatory.mask_set, 293](#)
[allensdk.internal.brain_observatory.ophys_session_id_mapper, 293](#)
[allensdk.internal.brain_observatory.resources, 283](#)
[allensdk.internal.brain_observatory.roi_filter_utilities, 294](#)
[allensdk.internal.brain_observatory.time_sync, 296](#)
[allensdk.internal.core, 301](#)
[allensdk.internal.core.lims_pipeline_module, 298](#)
[allensdk.internal.core.lims_utilities, 298](#)
[allensdk.internal.core.mouse_connectivity_cache_prerelease, 299](#)
[allensdk.internal.core.simpletree, 300](#)
[allensdk.internal.core.swc, 301](#)
[allensdk.internal.ephys, 307](#)
[allensdk.internal.ephys.core_feature_extract, 301](#)
[allensdk.internal.ephys.plot_qc_figures, 302](#)
[allensdk.internal.ephys.plot_qc_figures3, 304](#)
[allensdk.internal.model, 324](#)
[allensdk.internal.model.AIC, 322](#)
[allensdk.internal.model.biophysical, 311](#)
[allensdk.internal.model.biophysical.biophysical_array_api, 309](#)
[allensdk.internal.model.biophysical.check_fit_shift, 309](#)
[allensdk.internal.model.biophysical.deap_utils, 309](#)
[allensdk.internal.model.biophysical.ephys_utils, 309](#)
[allensdk.internal.model.biophysical.fits, 307](#)
[allensdk.internal.model.biophysical.fits.fit_styles, 307](#)
[allensdk.internal.model.biophysical.make_deap_fit, 310](#)
[allensdk.internal.model.biophysical.passive_fitting_report, 309](#)
[allensdk.internal.model.biophysical.passive_fitting_utilities, 309](#)
[allensdk.internal.model.biophysical.passive_fitting_utilities, 308](#)
[allensdk.internal.model.biophysical.passive_fitting_utilities, 308](#)

308 allensdk.internal.morphology.validate_swc,
 allensdk.internal.model.biophysical.passive_fit_elec, 331
 308 allensdk.internal.mouse_connectivity,
 allensdk.internal.model.biophysical.passive_fit_neuron_utils, 340
 308 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.passive_fit_output_grabber, 336
 308 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.passive_fit_preprocess, 332
 307 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.passive_fit_preprocess, 332
 308 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.run_optimize, 234
 310 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.run_optimize_workflow, 234
 311 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.run_passive_fit, 336
 311 allensdk.internal.mouse_connectivity.interval_union,
 allensdk.internal.model.biophysical.run_simulate_limits, 335
 311 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.biophysical.run_simulate_workflow, 338
 311 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.data_access, 323 336
 allensdk.internal.model.glif, 322 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.glif.are_two_lists_of_days_the_same, 336
 312 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.glif.error_functions, 337
 312 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.glif.find_spikes, 337
 312 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.glif.find_sweeps, 338
 312 allensdk.internal.mouse_connectivity.projection_threshold,
 allensdk.internal.model.glif.glif_experiment, 338
 313 allensdk.internal.mouse_connectivity.tissuecyte_stimulus,
 allensdk.internal.model.glif.glif_optimizer, 340
 315 allensdk.internal.mouse_connectivity.tissuecyte_stimulus,
 allensdk.internal.model.glif.glif_optimizer_neuron, 338
 315 allensdk.internal.mouse_connectivity.tissuecyte_stimulus,
 allensdk.internal.model.glif.MLIN, 311 339
 allensdk.internal.model.glif.plotting, allensdk.internal.pipeline_modules, 348
 318 allensdk.internal.pipeline_modules.gbm,
 allensdk.internal.model.glif.rc, 319 341
 allensdk.internal.model.glif.spike_cutting, allensdk.internal.pipeline_modules.gbm.generate_gbm, 340
 319 allensdk.internal.pipeline_modules.gbm.generate_gbm,
 allensdk.internal.model.glif.threshold_adaptation, allensdk.internal.pipeline_modules.gbm.generate_gbm, 340
 320 allensdk.internal.pipeline_modules.gbm.generate_gbm,
 allensdk.internal.model.GLM, 323 allensdk.internal.pipeline_modules.gbm.generate_gbm,
 allensdk.internal.morphology, 331 341
 allensdk.internal.morphology.compartment allensdk.internal.pipeline_modules.run_annotated_reconstruction, 341
 324 allensdk.internal.pipeline_modules.run_demixing,
 allensdk.internal.morphology.morphology, allensdk.internal.pipeline_modules.run_demixing, 341
 324 allensdk.internal.pipeline_modules.run_dff_computation,
 allensdk.internal.morphology.morphvis, allensdk.internal.pipeline_modules.run_dff_computation, 341
 328 allensdk.internal.pipeline_modules.run_neuropil_computation,
 allensdk.internal.morphology.node, 331 allensdk.internal.pipeline_modules.run_neuropil_computation, 341

[342](#) `allensdk.mouse_connectivity.grid.utilities.image_utilities,`
`allensdk.internal.pipeline_modules.run_observatory_analysis,`
[342](#) `allensdk.mouse_connectivity.grid.writers,`
`allensdk.internal.pipeline_modules.run_observatory_thumbnails,`
[343](#)
`allensdk.internal.pipeline_modules.runophys_eye_calibration,`
[345](#) `allensdk.test_utilities, 371
allensdk.internal.pipeline_modules.run_ophysestest_decomposition,
345 allensdk.test_utilities.custom_comparators,
345 allensdk.test_utilities.regression_fixture,
346 allensdk.test_utilities, 370
allensdk.internal.pipeline_modules.run_tassensdktest_utilizeitclassitemfromj350,
348`

m

`allensdk.model,` [362](#)
`allensdk.model.biophys_sim,` [350](#)
`allensdk.model.biophys_sim.bps_command,`
[349](#)
`allensdk.model.biophys_sim.config,` [349](#)
`allensdk.model.biophys_sim.neuron,` [349](#)
`allensdk.model.biophys_sim.neuron.hoc_utils,`
[348](#)
`allensdk.model.biophys_sim.scripts,` [349](#)
`allensdk.model.biophysical,` [352](#)
`allensdk.model.biophysical.run_simulate,`
[350](#)
`allensdk.model.biophysical.runner,` [350](#)
`allensdk.model.biophysical.utils,` [351](#)
`allensdk.model.glif,` [362](#)
`allensdk.model.glif.glif_neuron,` [352](#)
`allensdk.model.glif.glif_neuron_methods,`
[356](#)
`allensdk.model.glif.simulate_neuron,` [361](#)
`allensdk.morphology,` [362](#)
`allensdk.morphology.validate_swc,` [362](#)
`allensdk.mouse_connectivity,` [370](#)
`allensdk.mouse_connectivity.grid,` [370](#)
`allensdk.mouse_connectivity.grid.image_series_gridder,`
[369](#)
`allensdk.mouse_connectivity.grid.subimage,`
[366](#)
`allensdk.mouse_connectivity.grid.subimage.base_subimage,`
[363](#)
`allensdk.mouse_connectivity.grid.subimage.cav_subimage,`
[364](#)
`allensdk.mouse_connectivity.grid.subimage.classic_subimage,`
[365](#)
`allensdk.mouse_connectivity.grid.subimage.count_subimage,`
[366](#)
`allensdk.mouse_connectivity.grid.utilities,`
[368](#)
`allensdk.mouse_connectivity.grid.utilities.downsampling_utilities,`
[366](#)

A

`ab_from_diagonals()` (in module *allensdk.brain_observatory.r_neuropil*), 202
`ab_from_T()` (in module *allensdk.brain_observatory.r_neuropil*), 202
`accumulator_to_numpy()` (*allensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder* method), 369
`acquisition_rate` (*allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis* attribute), 212
`actual_parameters_from_normalized()` (*allensdk.internal.model.biophysical.deap_utils.Utils* method), 309
`adaptation_index()` (*allensdk.ephys.feature_extractor.EphysFeatureExtractor* method), 277
`adaptation_index()` (in module *allensdk.ephys.ephys_features*), 269
`add()` (*allensdk.brain_observatory.stimulus_info.BinaryIntervalSearchFree* method), 213
`add_angle_labels()` (in module *allensdk.brain_observatory.circle_plots*), 191
`add_arrow()` (in module *allensdk.brain_observatory.circle_plots*), 191
`add_file()` (*allensdk.config.manifest.Manifest* method), 227
`add_manifest_paths()` (*allensdk.api.warehouse_cache.cache.Cache* method), 97
`add_manifest_paths()` (*allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* method), 242
`add_manifest_paths()` (*allensdk.core.reference_space_cache.ReferenceSpaceCache* method), 253
`add_manifest_paths()` (*allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease* method), 300
`add_number_to_shuffled_movie()` (in module *allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities*), 172
`add_path()` (*allensdk.config.manifest.Manifest* method), 228
`add_path()` (*allensdk.config.manifest_builder.ManifestBuilder* method), 229
`add_paths()` (*allensdk.config.manifest.Manifest* method), 228
`add_section()` (*allensdk.config.manifest_builder.ManifestBuilder* method), 229
`add_to_fit_parameters_dict_single()` (in module *allensdk.brain_observatory.receptive_field_analysis.fit_parameters*), 182
`adjust_r_for_negativity()` (in module *allensdk.internal.pipeline_modules.run_neuropil_correction*), 342
`advance_combination()` (in module *allensdk.brain_observatory.chisquare_categorical*), 190
`age_in_days` (*allensdk.brain_observatory.behavior.metadata.behavior_metadata* attribute), 109
`AIC()` (in module *allensdk.internal.model.AIC*), 322
`AICc()` (in module *allensdk.internal.model.AIC*), 322
`align_and_cut_spikes()` (in module *allensdk.internal.model.glif.find_spikes*), 312
`align_running_speed()` (in module *allensdk.core.brain_observatory_nwb_data_set*), 236
`ALIGNMENT3D_KEY` (*allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* attribute), 242
`ALL` (*allensdk.api.queries.rma_api.RmaApi* attribute), 89
`all_stimuli()` (in module *allensdk.brain_observatory.stimulus_info*), 214
`AllActiveUtils` (class in *allensdk*)

(module), 118
 allensdk.brain_observatory.behavior.session_handlers (module), 119
 allensdk.brain_observatory.behavior.session_handlers.brain_observatory.ecephys.align_timestamps (module), 125
 allensdk.brain_observatory.behavior.session_handlers.brain_observatory.ecephys.align_timestamps (module), 121
 allensdk.brain_observatory.behavior.session_handlers.brain_observatory.ecephys.align_timestamps (module), 123
 allensdk.brain_observatory.behavior.session_handlers.brain_observatory.ecephys.current_source_model (module), 149
 allensdk.brain_observatory.behavior.stimulus_handlers.brain_observatory.ecephys.ecephys_project (module), 128
 allensdk.brain_observatory.behavior.stimulus_handlers.brain_observatory.ecephys.ecephys_session (module), 126
 allensdk.brain_observatory.behavior.stimulus_handlers.brain_observatory.ecephys.file_io (module), 128
 allensdk.brain_observatory.behavior.sync_handlers.brain_observatory.ecephys.file_io.continuous (module), 131
 allensdk.brain_observatory.behavior.sync_handlers.brain_observatory.ecephys.file_io.ecephys (module), 130
 allensdk.brain_observatory.behavior.trials_handlers.brain_observatory.ecephys.file_io.stim_files (module), 149
 allensdk.brain_observatory.behavior.trials_handlers.brain_observatory.ecephys.lfp_subsampling (module), 150
 allensdk.brain_observatory.behavior.write_handlers.brain_observatory.ecephys.lfp_subsampling (module), 134
 allensdk.brain_observatory.behavior.write_handlers.brain_observatory.ecephys.optotagging_tables (module), 134
 allensdk.brain_observatory.behavior.write_handlers.brain_observatory.ecephys.stimulus_sync (module), 134
 allensdk.brain_observatory.behavior.write_handlers.brain_observatory.ecephys.stimulus_table (module), 134
 allensdk.brain_observatory.behavior.write_handlers.brain_observatory.ecephys.stimulus_table.ecephys (module), 134
 allensdk.brain_observatory.brain_observatory_handlers.brain_observatory.ecephys.stimulus_table.ecephys (module), 189
 allensdk.brain_observatory.brain_observatory_handlers.brain_observatory.ecephys.stimulus_table.ecephys (module), 189
 allensdk.brain_observatory.chisquare_category (module), 190
 allensdk.brain_observatory.circle_plots (module), 190
 allensdk.brain_observatory.demixer (module), 192
 allensdk.brain_observatory.dff (module), 193
 allensdk.brain_observatory.drifting_gratings (module), 195
 allensdk.brain_observatory.ecephys (module), 177
 allensdk.brain_observatory.ecephys.align_timestamps (module), 160
 allensdk.brain_observatory.ecephys.align_timestamps.brain_observatory.gaze_mapping (module), 156
 allensdk.brain_observatory.ecephys.align_timestamps.brain_observatory.gaze_mapping (module), 158
 allensdk.brain_observatory.ecephys.align_timestamps.brain_observatory.gaze_mapping (module), 159
 allensdk.brain_observatory.ecephys.align_timestamps.brain_observatory.gaze_mapping (module), 160
 allensdk.brain_observatory.ecephys.align_timestamps.brain_observatory.gaze_mapping (module), 160
 allensdk.brain_observatory.ecephys.current_source_model (module), 160
 allensdk.brain_observatory.ecephys.ecephys_project (module), 161
 allensdk.brain_observatory.ecephys.ecephys_session (module), 162
 allensdk.brain_observatory.ecephys.file_io (module), 164
 allensdk.brain_observatory.ecephys.file_io.continuous (module), 163
 allensdk.brain_observatory.ecephys.file_io.ecephys (module), 163
 allensdk.brain_observatory.ecephys.file_io.stim_files (module), 163
 allensdk.brain_observatory.ecephys.lfp_subsampling (module), 166
 allensdk.brain_observatory.ecephys.lfp_subsampling (module), 164
 allensdk.brain_observatory.ecephys.optotagging_tables (module), 166
 allensdk.brain_observatory.ecephys.stimulus_sync (module), 176
 allensdk.brain_observatory.ecephys.stimulus_table (module), 175
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 168
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 172
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 173
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 174
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 168
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys (module), 167
 allensdk.brain_observatory.ecephys.visualization (module), 175
 allensdk.brain_observatory.ecephys.write_nwb (module), 176
 allensdk.brain_observatory.extract_running_speed (module), 177
 allensdk.brain_observatory.findlevel (module), 197
 allensdk.brain_observatory.gaze_mapping (module), 197

(module), 177
 allensdk.brain_observatory.locally_sparse_embeddings (module), 197
 allensdk.brain_observatory.natural_movie_embeddings (module), 198
 allensdk.brain_observatory.natural_scenes_embeddings (module), 199
 allensdk.brain_observatory.observatory_plots (module), 200
 allensdk.brain_observatory.ophys (module), 178
 allensdk.brain_observatory.ophys.trace_embeddings (module), 178
 allensdk.brain_observatory.r_neuropil (module), 202
 allensdk.brain_observatory.receptive_field_analysis (module), 186
 allensdk.brain_observatory.receptive_field_analysis_high_resolution (module), 178
 allensdk.brain_observatory.receptive_field_analysis_low_resolution (module), 181
 allensdk.brain_observatory.receptive_field_analysis_tools (module), 182
 allensdk.brain_observatory.receptive_field_analysis_utilities (module), 182
 allensdk.brain_observatory.receptive_field_analysis_processing (module), 183
 allensdk.brain_observatory.receptive_field_analysis_receptive_field (module), 183
 allensdk.brain_observatory.receptive_field_analysis_tools (module), 184
 allensdk.brain_observatory.receptive_field_analysis_utilities (module), 184
 allensdk.brain_observatory.receptive_field_analysis_utilities (module), 185
 allensdk.brain_observatory.roi_masks (module), 203
 allensdk.brain_observatory.running_speed (module), 206
 allensdk.brain_observatory.session_analysis (module), 206
 allensdk.brain_observatory.session_api_utilities (module), 209
 allensdk.brain_observatory.static_gratings (module), 210
 allensdk.brain_observatory.stimulus_analysis (module), 211
 allensdk.brain_observatory.stimulus_info (module), 213
 allensdk.brain_observatory.sync_dataset (module), 217
 allensdk.brain_observatory.sync_utilities (module), 187
 allensdk.brain_observatory.visualization (module), 187
 allensdk.config (module), 230
 allensdk.config.app (module), 224
 allensdk.config.app.application_config (module), 222
 allensdk.config.manifest (module), 227
 allensdk.config.manifest_builder (module), 229
 allensdk.config.model (module), 227
 allensdk.config.model.description (module), 226
 allensdk.config.model.description_parser (module), 227
 allensdk.config.model.formats (module), 226
 allensdk.config.model.formats.hdf5_util (module), 224
 allensdk.config.model.formats.json_description_parser (module), 224
 allensdk.config.model.formats.pycfg_description_parser (module), 225
 allensdk.core (module), 231
 allensdk.core.auth_config (module), 231
 allensdk.core.authentication_and_session (module), 231
 allensdk.core.brain_observatory_nwb_data_set (module), 232
 allensdk.core.cache_method_utilities (module), 236
 allensdk.core.cell_types_cache (module), 236
 allensdk.core.dat_utilities (module), 239
 allensdk.core.exceptions (module), 239
 allensdk.core.h5_utilities (module), 240
 allensdk.core.initialization_utilities (module), 240
 allensdk.core.lazy_property (module), 231
 allensdk.core.lazy_property.lazy_property (module), 230
 allensdk.core.lazy_property.lazy_property_mixin (module), 230
 allensdk.core.mouse_connectivity_cache (module), 241
 allensdk.core.nwb_data_set (module), 246
 allensdk.core.obj_utilities (module), 248
 allensdk.core.ontology (module), 249
 allensdk.core.ophys_experiment_session_id_mapping (module), 250
 allensdk.core.reference_space (module), 250
 allensdk.core.reference_space_cache (module), 252
 allensdk.core.simple_tree (module), 254
 allensdk.core.sitk_utilities (module), 257
 allensdk.core.structure_tree (module), 258
 allensdk.core.swc (module), 261

427

allensdk.internal.model.biophysical.run_optimize (module), 332
 (module), 311
 allensdk.internal.model.biophysical.run_passive (module), 334
 (module), 311
 allensdk.internal.model.biophysical.run_simulate (module), 334
 (module), 311
 allensdk.internal.model.biophysical.run_simulate (module), 335
 (module), 311
 allensdk.internal.model.data_access (module), 335
 (module), 323
 allensdk.internal.model.glif (module), 322 (module), 338
 allensdk.internal.model.glif.are_two_lists_of_neurons_the_same (module), 336
 (module), 312
 allensdk.internal.model.glif.error_function (module), 336
 (module), 312
 allensdk.internal.model.glif.find_spikes (module), 337
 (module), 312
 allensdk.internal.model.glif.find_sweeps (module), 337
 (module), 312
 allensdk.internal.model.glif.glif_experiment (module), 338
 (module), 313
 allensdk.internal.model.glif.glif_optimize (module), 338
 (module), 315
 allensdk.internal.model.glif.glif_optimize (module), 340
 (module), 315
 allensdk.internal.model.glif.MLIN (module), 311 (module), 338
 allensdk.internal.model.glif.plotting (module), 339
 (module), 318
 allensdk.internal.model.glif.rc (module), 319 (module), 348
 allensdk.internal.model.glif.spike_cutting (module), 341
 (module), 319
 allensdk.internal.model.glif.threshold_adaptation (module), 340
 (module), 320
 allensdk.internal.model.GLM (module), 323 (module), 340
 allensdk.internal.morphology (module), 331 (module), 341
 allensdk.internal.morphology.compartmental (module), 341
 (module), 324
 allensdk.internal.morphology.morphology (module), 341
 (module), 324
 allensdk.internal.morphology.morphvis (module), 341
 (module), 328
 allensdk.internal.morphology.node (module), 341
 (module), 331
 allensdk.internal.morphology.validate_swap (module), 342
 (module), 331
 allensdk.internal.mouse_connectivity (module), 342
 (module), 340
 allensdk.internal.mouse_connectivity.interval_union (module), 343
 (module), 336
 allensdk.internal.mouse_connectivity.interval_union_endpoints (module), 345
 (module), 332
 allensdk.internal.mouse_connectivity.interval_union_endpoints (module), 345
 (module), 332
 allensdk.internal.mouse_connectivity.interval_union_endpoints (module), 345
 (module), 332
 allensdk.internal.pipeline_modules (module), 348
 allensdk.internal.pipeline_modules.gbm (module), 341
 allensdk.internal.pipeline_modules.gbm.generate_gbm (module), 340
 allensdk.internal.pipeline_modules.gbm.generate_gbm (module), 340
 allensdk.internal.pipeline_modules.gbm.generate_gbm (module), 341
 allensdk.internal.pipeline_modules.run_annotated_re (module), 341
 allensdk.internal.pipeline_modules.run_demixing (module), 341
 allensdk.internal.pipeline_modules.run_dff_computat (module), 341
 allensdk.internal.pipeline_modules.run_neuropil_cor (module), 342
 allensdk.internal.pipeline_modules.run_observatory (module), 342
 allensdk.internal.pipeline_modules.run_observatory (module), 343
 allensdk.internal.pipeline_modules.run_ophys_eye_ca (module), 345
 allensdk.internal.pipeline_modules.run_ophys_sessio (module), 345

(*module*), 345

allensdk.internal.pipeline_modules.run_ophyensdkemysyn_connectivity.grid.writers
(*module*), 346

allensdk.internal.pipeline_modules.run_tassensdketanonizeitclass(*module*), 371

allensdk.model (*module*), 362

allensdk.model.biophys_sim (*module*), 350

allensdk.model.biophys_sim.bps_command
(*module*), 349

allensdk.model.biophys_sim.config (*module*), 349

allensdk.model.biophys_sim.neuron (*module*), 349

allensdk.model.biophys_sim.neuron.hoc_util_alpha_filter() (in *module* *allensdk.brain_observatory.r_neuropil*), 202

allensdk.model.biophys_sim.scripts (*module*), 349

allensdk.model.biophysical (*module*), 352

allensdk.model.biophysical.run_simulate
(*module*), 350

allensdk.model.biophysical.runner (*module*), 350

allensdk.model.biophysical.utils (*module*), 351

allensdk.model.glif (*module*), 362

allensdk.model.glif.glif_neuron (*module*), 352

allensdk.model.glif.glif_neuron_methods
(*module*), 356

allensdk.model.glif.simulate_neuron
(*module*), 361

allensdk.morphology (*module*), 362

allensdk.morphology.validate_swc (*module*), 362

allensdk.mouse_connectivity (*module*), 370

allensdk.mouse_connectivity.grid (*module*), 370

allensdk.mouse_connectivity.grid.image_series_*attribute*, 164
(*module*), 369

allensdk.mouse_connectivity.grid.subimage
(*module*), 366

allensdk.mouse_connectivity.grid.subimage_enhanced_subimage() (in *module* *allensdk.brain_observatory.behavior.mtrain*), 142
(*module*), 363

allensdk.mouse_connectivity.grid.subimage.cav_subimage
(*module*), 364

allensdk.mouse_connectivity.grid.subimage.class_*attribute*, 68
(*module*), 365

allensdk.mouse_connectivity.grid.subimage_annotation_key (*module*), 366

allensdk.mouse_connectivity.grid.utilities
(*module*), 368

allensdk.mouse_connectivity.grid.utilities.abicallsempring (*module*), 366

allensdk.mouse_connectivity.grid.utilities.abicallgendrites (*module*), 367

allensdk.mouse_connectivity.grid.utilities.abicallsempring (*module*), 369

allensdk.test_utilities.custom_comparators
(*module*), 370

allensdk.test_utilities.regression_fixture
(*module*), 370

allensdk.test_utilities.temp_dir (*module*), 370

allocate_by_vsync() (in *module* *allensdk.brain_observatory.ecephys.stimulus_sync*), 176

alpha_filter() (in *module* *allensdk.brain_observatory.r_neuropil*), 202

analog_meta_data (*allensdk.brain_observatory.sync_dataset.Dataset*
attribute), 218

analyze_trough_details() (in *module* *allensdk.ephys.ephys_features*), 269

ancestor_ids() (*allensdk.core.simple_tree.SimpleTree*
method), 254

ancestor_ids() (*allensdk.internal.core.simpletree.SimpleTree*
method), 300

ancestors() (*allensdk.core.simple_tree.SimpleTree*
method), 255

ancestors() (*allensdk.internal.core.simpletree.SimpleTree*
method), 300

angle_lines() (in *module* *allensdk.brain_observatory.circle_plots*), 192

angular_wheel_rotation (*allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOneP*
attribute), 163

angular_wheel_velocity (*allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOneP*
attribute), 164

annotate_change_detect() (in *module* *allensdk.brain_observatory.behavior.mtrain*), 142

annotated_section_data_sets_api (*class* in *allensdk.api.query*
allensdk.api.query.annotated_section_data_sets_api), 68

ANNOTATION_KEY (*allensdk.core.reference_space_cache.ReferenceSpaceCache*
attribute), 252

Api (*class* in *allensdk.api.api*), 103

ABICALLSEMPRING (*allensdk.core.swc.Morphology*
attribute), 261

ABICALLGENDRITES (*allensdk.core.swc.Morphology*
attribute), 261

[lensdk.internal.morphology.morphology.Morphology.apply_average_tile\(\)](#) (al-
[attribute](#)), 324 [lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile](#)
[lensdk.internal.api.api_prerelease](#) (class in al-
[lensdk.internal.api.api_prerelease](#)), 282 [method](#)), 339
[append\(\)](#) ([allensdk.core.swc.Morphology](#) [method](#)), [lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile](#)
261 [method](#)), 339
[append\(\)](#) ([lensdk.internal.morphology.morphology.Morphology.colormap\(\)](#) (in module al-
[method](#)), 324 [lensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet](#)
[append\(\)](#) ([lensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet](#)
[method](#)), 336 [apply_configuration_from_command_line\(\)](#)
[append_experiment_metrics\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 222
[method](#)), 206 [apply_configuration_from_environment\(\)](#)
[append_metadata\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 222
[method](#)), 206 [apply_configuration_from_file\(\)](#) (al-
[append_metrics_drifting_grating\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 222
[method](#)), 206 [lensdk.config.app.application_config.ApplicationConfig](#)
[append_metrics_locally_sparse_noise\(\)](#) ([allensdk.brain_observatory.session_analysis.SessionAnalysis](#)
[method](#)), 206 [apply_display_sequence\(\)](#) (in module al-
[lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes](#)
[method](#)), 168
[append_metrics_natural_movie_one\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 206
[method](#)), 206 [lensdk.mouse_connectivity.grid.utilities.downsampling_utilities](#),
366
[append_metrics_natural_movie_three\(\)](#) ([allensdk.brain_observatory.session_analysis.SessionAnalysis](#)
[method](#)), 206 [apply_frame_times\(\)](#) (in module al-
[lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes](#)
[method](#)), 168
[append_metrics_natural_movie_two\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 206
[method](#)), 206 [apply_mask\(\)](#) ([allensdk.mouse_connectivity.grid.subimage.base_subimage](#)
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 364
[method](#)), 206 [apply_pixel_counter\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 364
[method](#)), 206 [lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage](#)
[append_metrics_natural_scene\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 206
[method](#)), 206 [apply_ssl\(\)](#) ([allensdk.api.queries.reference_space_api.ReferenceSpaceApi](#)
[attribute](#)), 87
[append_metrics_static_grating\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 207
[method](#)), 207 [archive_cell\(\)](#) (al-
[lensdk.brain_observatory.session_analysis.SessionAnalysis](#) [method](#)), 309
[method](#)), 309
[append_threshold_components\(\)](#) (al-
[lensdk.model.glif.glif_neuron.GlifNeuron](#) [method](#)), 354
[method](#)), 354 (in module al-
[lensdk.internal.model.glif.are_two_lists_of_arrays_the_same](#)),
312
[append_well_known_file\(\)](#) (in module al-
[lensdk.internal.core.lims_utilities](#)), 298 [arg_parser\(\)](#) (in module al-
[lensdk.config.app.application_config](#)), 222 [lensdk.internal.model.biophysical.passive_fitting.neuron_passive](#)
308
[apply\(\)](#) ([allensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet](#)
[method](#)), 336 [lensdk.internal.pipeline_module.PipelineModule](#)
[attribute](#)), 298
[apply_affine\(\)](#) ([allensdk.core.swc.Morphology](#) [ArgSchemaParserPlus](#) (class in al-
[method](#)), 261 [lensdk.brain_observatory.argschema_utilities](#)),
188
[apply_affine\(\)](#) (al-
[lensdk.internal.morphology.morphology.Morphology](#) [method](#)), 324
[method](#)), 324 [ARRAY](#) ([allensdk.api.queries.connected_services.ConnectedServices](#)
[attribute](#)), 75
[apply_affine_only_rotation\(\)](#) (al-
[lensdk.internal.morphology.morphology.Morphology](#) [method](#)), 228
[method](#)), 325 [as_dataframe\(\)](#) ([allensdk.config.manifest.Manifest](#)
[method](#)), 228
[as_dataframe\(\)](#) (al-

[lensdk.config.manifest_builder.ManifestBuilder](#) (in module [al-lensdk.internal.brain_observatory.eye_calibration](#)),
[method](#)), 229
[as_dict\(\)](#) ([allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor](#)
[method](#)), 265
[as_dict\(\)](#) ([allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor](#)
[method](#)), 266
[aspect_ratio\(\)](#) ([allensdk.brain_observatory.stimulus_info.Monitor](#)
[attribute](#)), 214
[assert_exists\(\)](#) (in module [al-lensdk.internal.pipeline_modules.run_demixing](#)), [behavior_delta](#)
341
[assign_session_id\(\)](#) (in module [al-lensdk.brain_observatory.behavior.mtrain](#)), [behavior_session_id](#)
142
[assign_sweep_values\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.stimulus_table.EphysVireospires](#)
169
[assign_to_last\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.stimulus_sync](#)), [behavior_times](#)
176
[atlas_image_query\(\)](#) ([al-lensdk.api.queries.image_download_api.ImageDownloadAPI](#)
[method](#)), 79
[autocorr\(\)](#) (in module [al-lensdk.internal.model.glif.MLIN](#)), 311
[average_rate\(\)](#) (in module [al-lensdk.ephys.ephys_features](#)), 270
AVERAGE_TEMPLATE ([al-lensdk.api.queries.reference_space_api.ReferenceSpaceApi](#)
[attribute](#)), 87
[average_tile_is_untrimmed\(\)](#) ([al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile_tile](#)
[method](#)), 339
[average_voltage\(\)](#) (in module [al-lensdk.ephys.ephys_features](#)), 270
AXON ([allensdk.core.swc.Morphology](#) [attribute](#)), 261
AXON ([allensdk.internal.morphology.morphology.Morphology](#)
[attribute](#)), 324
B
[background_trace\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)),
285
[barcode_line\(\)](#) ([allensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset.BarcodeSyncDataset](#)
[attribute](#)), 158
BarcodeSyncDataset (class in [al-lensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset](#)),
158
BASAL_DENDRITE ([allensdk.core.swc.Morphology](#) [attribute](#)), 261
BASAL_DENDRITE ([al-lensdk.internal.morphology.morphology.Morphology](#)
[attribute](#)), 324
[base_object_to_eye_rotation_matrix\(\)](#)

BehaviorOphysBase (class in al- 337
 lensdk.brain_observatory.behavior.session_apis.abstraction_base.behavior_ophys_base), (in module al-
 116 lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher),
 BehaviorOphysDataExtractorBase (class in al- 338
 lensdk.brain_observatory.behavior.session_apis.abstraction_base.behavior_ophys_data_extractor_base),
 113 lensdk.internal.mouse_connectivity.projection_thumbnail.generator,
 BehaviorOphysDataTransforms (class in al- 336
 lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_ophys_data_transforms), al-
 123 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 BehaviorOphysMetadata (class in al- 367
 lensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata), (in module al-
 111 lensdk.mouse_connectivity.grid.utilities.downsampling_utilities),
 BehaviorOphysSessionsTable (class in al- 366
 lensdk.brain_observatory.behavior.behavior_projects_backend.behavior_ophys_sessions_table), (in module al-
 107 lensdk.brain_observatory.behavior.projects_backend.behavior_ophys_sessions_table), (in module al-
 75 lensdk.brain_observatory.projects_backend.behavior_ophys_sessions_table), (in module al-
 BrainObservatoryAnalysisException, 189
 BehaviorTaskParametersSchema (class in al- BrainObservatoryApi (class in al-
 lensdk.brain_observatory.behavior.schemas), 145
 lensdk.api.queries.brain_observatory_api),
 best_fit_value() (al- 70
 lensdk.internal.model.biophysical.make_deep_fit_for_report),
 BIC() (in module allensdk.internal.model.AIC), 323 213
 binarize() (allensdk.mouse_connectivity.grid.subimage_base.subimage_subimage_webDataSet (class in al-
 method), 364 lensdk.core.brain_observatory_nwb_data_set),
 BinaryIntervalSearchTree (class in al- 232
 lensdk.brain_observatory.stimulus_info), bucket_name_from_url() (in module al-
 213 lensdk.api.cloud_cache.utils), 67
 binned_cells_sp (al- build_affine_transform() (in module al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis, 367
 attribute), 212 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 binned_cells_vis (al- build_cell_plots() (in module al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis, 343
 attribute), 212 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 binned_dx_sp (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis, 343
 attribute), 212 lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder, 369
 binned_dx_vis (al- build_lormap() (in module al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis, 343
 attribute), 212 lensdk.brain_observatory.ecephys.stimulus_table.visualization.visualization_utilities),
 BIOPHYSICAL_MODEL_TYPE_IDS (al- 167
 lensdk.api.queries.biophysical_api.BiophysicalApi, 69
 attribute), 69 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 BiophysicalApi (class in al- 367
 lensdk.api.queries.biophysical_api), 69 build_correlation_plots() (in module al-
 BiophysicalArchiver (class in al- lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 lensdk.internal.model.biophysical.biophysical_archiver), 343
 309 build_drifting_gratings() (in module al-
 BiophysicalModuleApi (class in al- lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 lensdk.internal.api.queries.biophysical_module_api), 343
 278 build_experiment_thumbnails()
 BiophysicalModuleReader (class in al- (in module al-
 lensdk.internal.api.queries.biophysical_module_reader), lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 278 343
 blend() (in module al- build_eye_tracking_plots() (in module al-
 lensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities), lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 343

[343](#)
 build_from_image() (al- [build_trial_matrix\(\)](#) (in module al-
 lensdk.internal.mouse_connectivity.projection_thumbnail_image_shoot.
 static method), [336](#)
[169](#)
 build_hex_pack() (in module al- [build_type\(\)](#) (in module al-
 lensdk.brain_observatory.circle_plots), [192](#) lensdk.internal.pipeline_modules.run_observatory_thumbnails),
[178](#)
 build_locally_sparse_noise() (in module al- [344](#)
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), (allensdk.api.queries.connected_services.ConnectedService
[343](#) method), [76](#)
 build_manifest() (al- [build_volumetric_data_download_url\(\)](#)
 lensdk.api.warehouse_cache.cache.Cache (allensdk.api.queries.reference_space_api.ReferenceSpaceApi
 method), [97](#) method), [87](#)
 build_manifest() (al- [burst_metrics\(\)](#) (al-
 lensdk.core.cell_types_cache.CellTypesCache lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor
 method), [237](#) method), [266](#)
 build_manifest() (al- [OptimizeConfigReader](#)
 lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader
 method), [281](#) Cache (class in allensdk.api.warehouse_cache.cache),
[97](#)
 build_natural_movie() (in module al- [97](#)
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), [cache_clear\(\)](#) (al-
[343](#) lensdk.core.cache_method_utilities.CachedInstanceMethodMixin
 method), [236](#)
 build_natural_scenes() (in module al- [344](#)
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), (allensdk.api.warehouse_cache.cache.Cache
 static method), [97](#)
 build_plots() (in module al- [cache_csv_dataframe\(\)](#) (al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), lensdk.api.warehouse_cache.cache.Cache
[344](#) static method), [97](#)
 build_query() (allensdk.api.queries.svg_api.SvgApi [cache_csv_json\(\)](#) (al-
 method), [94](#) lensdk.api.warehouse_cache.cache.Cache
 static method), [97](#)
 build_query_url() (al- [cache_data\(\)](#) (allensdk.api.queries.biophysical_api.BiophysicalApi
 lensdk.api.queries.rma_api.RmaApi method), [90](#) method), [69](#)
 build_receptive_field() (in module al- [cache_json\(\)](#) (allensdk.api.warehouse_cache.cache.Cache
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), static method), [97](#)
[344](#) [cache_json_dataframe\(\)](#) (al-
 build_reference_aligned_image_channel_volumes_url() lensdk.api.warehouse_cache.cache.Cache
 (allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi static method), [97](#)
 method), [82](#) [cache_stimulus_file\(\)](#) (al-
 build_rma() (allensdk.api.queries.biophysical_api.BiophysicalApi lensdk.api.queries.glif_api.GlifApi method),
 method), [69](#) [76](#)
 build_rotation_transform() (al- [cacheable\(\)](#) (in module al-
 lensdk.internal.mouse_connectivity.projection_thumbnail.volumetric_projector.VolumetricProjector, [99](#)
 method), [338](#) lensdk.api.warehouse_cache.cache.Cache),
 CachedInstanceMethodMixin (class in al-
 build_schema_query() (al- [lensdk.core.cache_method_utilities\)](#), [236](#)
 lensdk.api.queries.rma_api.RmaApi method), CacheFileAttributes (class in al-
[90](#) lensdk.api.cloud_cache.file_attributes), [65](#)
 build_speed_tuning() (in module al- [cacher\(\)](#) (allensdk.api.warehouse_cache.cache.Cache
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), static method), [97](#)
[344](#)
 build_static_gratings() (in module al- [calc_deriv\(\)](#) (in module al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails), lensdk.brain_observatory.behavior.running_processing),
[344](#) [43](#)
 build_stimuluswise_table() (in module al- [calc_spike_component_of_threshold_from_multiblip\(\)](#)
 lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes), (in module al-
 lensdk.internal.model.glif.threshold_adaptation),

[320](#)
`calc_spike_cut_and_v_reset_via_expvar_residuals()` (in module `lensdk.internal.model.glif.spike_cutting`), [319](#)
`calculate()` (`allensdk.core.lazy_property.lazy_property.LazyProperty` method), [230](#)
`calculate()` (`allensdk.internal.mouse_connectivity.interval_unionize_tissuecyte_unionize_record.TissuecyteInjectionUnionize` method), [335](#)
`calculate()` (`allensdk.internal.mouse_connectivity.interval_unionize_tissuecyte_unionize_record.TissuecyteInjectionUnionize` method), [335](#)
`calculate()` (`allensdk.internal.mouse_connectivity.interval_unionize_tissuecyte_unionize_record.TissuecyteInjectionUnionize` method), [335](#)
`calculate_delay()` (in module `allensdk.brain_observatory.behavior.sync.process_sync`), [130](#)
`calculate_dff()` (in module `allensdk.brain_observatory.dff`), [193](#)
`calculate_dvdt()` (in module `allensdk.ephys.ephys_features`), [271](#)
`calculate_feature_errors()` (`allensdk.internal.model.biophysical.deap_utils.Utils` method), [309](#)
`calculate_fi_curves()` (in module `allensdk.internal.model.biophysical.check_fi_shift`), [309](#)
`calculate_injection_centroid()` (`allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), [83](#)
`calculate_max_border()` (in module `allensdk.internal.brain_observatory.roi_filter_utils`), [295](#)
`calculate_monitor_delay()` (in module `allensdk.internal.brain_observatory.time_sync`), [297](#)
`calculate_response_latency_list()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [150](#)
`calculate_reward_rate()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [150](#)
`calculate_reward_rate_fix_nans()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [150](#)
`calculate_roi_and_neuropil_traces()` (in module `allensdk.brain_observatory.roi_masks`), [204](#)
`calculate_scale()` (in module `allensdk.internal.morphology.morphvis`), [328](#)
`calculate_traces()` (in module `allensdk.brain_observatory.roi_masks`), [204](#)
`calculate_trough()` (in module `allensdk.brain_observatory.roi_masks`), [204](#)
`lensdk.ephys.feature_extractor.EphysFeatureExtractor` method), [277](#)
`call_caching()` (in module `allensdk.api.warehouse_cache.caching_utilities`), [101](#)
`CanPickleStimFile` (class in `allensdk.brain_observatory.ecephys.file_io.stim_file`), [162](#)
`categorize_one_trial()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [150](#)
`cell_extractor_for_nwb()` (in module `allensdk.ephys.ephys_extractor`), [269](#)
`cell_features()` (`allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` method), [265](#)
`cell_id()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` attribute), [212](#)
`cell_index_receptive_field_analysis_data()` (`allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` attribute), [197](#)
`CELL_MAPPING_ID` (`allensdk.brain_observatory.behavior.trials_processing` attribute), [70](#)
`CELLS_KEY` (`allensdk.core.cell_types_cache.CellTypesCache` attribute), [237](#)
`celltraces()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` attribute), [212](#)
`CellTypesApi` (class in `allensdk.api.queries.cell_types_api`), [73](#)
`CellTypesCache` (class in `allensdk.core.cell_types_cache`), [236](#)
`change_parent()` (`allensdk.core.swc.Morphology` method), [262](#)
`change_parent()` (`allensdk.internal.morphology.morphology.Morphology` method), [325](#)
`check_and_write()` (`allensdk.core.reference_space.ReferenceSpace` attribute), [277](#)

static method), 250

check_coverage() (allensdk.core.reference_space.ReferenceSpace method), 250

check_dir() (allensdk.config.manifest.Manifest method), 228

check_org_selections_for_noise_block() (allensdk.internal.model.biophysical.make_deap_fit_report method), 310

check_read_access() (in module allensdk.brain_observatory.argschema_utilities), 189

check_stimulus_delay() (in module allensdk.internal.pipeline_modules.run_ophys_time_sync method), 347

check_thresholds_and_peaks() (in module allensdk.ephys.ephys_features), 271

check_write_access() (in module allensdk.brain_observatory.argschema_utilities), 189

check_write_access_dir() (in module allensdk.brain_observatory.argschema_utilities), 189

check_write_access_overwrite() (in module allensdk.brain_observatory.argschema_utilities), 189

checkPreprocess() (in module allensdk.internal.model.glif.plotting), 318

checkSpikeCutting() (in module allensdk.internal.model.glif.plotting), 318

chi_square_binary() (in module allensdk.brain_observatory.receptive_field_analysis.chisquare method), 179

chi_square_within_mask() (in module allensdk.brain_observatory.receptive_field_analysis.chisquare method), 179

child_ids() (allensdk.core.simple_tree.SimpleTree method), 255

child_ids() (allensdk.internal.core.simpletree.SimpleTree method), 300

children() (allensdk.core.simple_tree.SimpleTree method), 255

children() (allensdk.internal.core.simpletree.SimpleTree method), 300

children_of() (allensdk.core.swc.Morphology method), 262

children_of() (allensdk.internal.morphology.morphology.Morphology method), 325

chisq_from_stim_table() (in module allensdk.brain_observatory.chisquare_categorical), 190

choose_bps_command() (in module allensdk.model.biophys_sim.bps_command), 349

choose_inliers() (allensdk.internal.brain_observatory.fit_ellipse.FitEllipse method), 290

class_deprecated() (in module allensdk.deprecated), 371

class_report_iter() (in module allensdk.mouse_connectivity.grid.writers), 369

ClassicSubImage (class in allensdk.mouse_connectivity.grid.subimage.classic_subimage), 365

clone_structures() (allensdk.core.structure_tree.StructureTree static method), 258

cleanup_truncated_file() (allensdk.api.api.Api method), 103

clear_updated_params() (allensdk.brain_observatory.session_api_utils.ParamsMixin method), 209

clipped_stim_timestamps (allensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute), 296

clone() (allensdk.ephys.feature_extractor.EphysFeatures method), 277

clone() (allensdk.internal.morphology.morphology.Morphology method), 325

close() (allensdk.brain_observatory.sync_dataset.Dataset method), 218

close() (allensdk.internal.brain_observatory.frame_stream.CvInputStream method), 290

close() (allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream method), 290

close() (allensdk.internal.brain_observatory.frame_stream.FfmpegOutputStream method), 291

close() (allensdk.internal.brain_observatory.frame_stream.FrameInputStream method), 291

close() (allensdk.internal.brain_observatory.frame_stream.FrameOutputStream method), 291

close() (allensdk.internal.brain_observatory.mask_set.MaskSet method), 292

close_sets() (allensdk.internal.brain_observatory.mask_set.MaskSet method), 292

collapse_columns() (in module allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities), 172

collect_sets() (allensdk.core.structure_tree.StructureTree static method), 259

colormap() (in module allensdk.brain_observatory.behavior.trials_processing), 151

COLORMAPS (allensdk.api.queries.image_download_api.ImageDownloadAPI

[attribute](#)), 78
[compare_session_fields\(\)](#) (in module [al-lensdk.brain_observatory.session_api_utils](#)), 209
[Compartment](#) (class in [allensdk.core.swc](#)), 261
[Compartment](#) (class in [al-lensdk.internal.morphology.compartment](#)), 324
[compartment\(\)](#) ([al-lensdk.internal.morphology.morphology.Morphology](#) method), 325
[compartment_index](#) ([allensdk.core.swc.Morphology](#) attribute), 262
[compartment_index_by_type\(\)](#) ([al-lensdk.core.swc.Morphology](#) method), 262
[compartment_list](#) ([allensdk.core.swc.Morphology](#) attribute), 262
[compartment_list](#) ([al-lensdk.internal.morphology.morphology.Morphology](#) attribute), 326
[compartment_list_by_type\(\)](#) ([al-lensdk.core.swc.Morphology](#) method), 262
[CompleteOphysBehaviorMetadataSchema](#) (class in [al-lensdk.brain_observatory.behavior.schemas](#)), 145
[compute\(\)](#) ([allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer.ProbeSynchronizer](#) class method), 160
[compute_area\(\)](#) ([al-lensdk.internal.brain_observatory.eye_calibration.EyeCalibrator](#) method), 287
[compute_chi\(\)](#) (in module [al-lensdk.brain_observatory.chisquare_categorical](#)), 190
[compute_chi_shuffle\(\)](#) (in module [al-lensdk.brain_observatory.chisquare_categorical](#)), 190
[compute_circular_area\(\)](#) (in module [al-lensdk.brain_observatory.behavior.eye_tracking_processing](#)), 138
[compute_coarse_parameters\(\)](#) (in module [al-lensdk.mouse_connectivity.grid.utilities.image_utilities](#)), 367
[compute_coarse_planes\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage](#) method), 364
[compute_coarse_planes\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage](#) method), 364
[compute_coarse_planes\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage](#) method), 365
[compute_coarse_planes\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage](#) method), 365
[compute_correlations\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)), 285
[compute_correlations_without_masks\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)), 285
[compute_dff_windowed_median\(\)](#) (in module [allensdk.brain_observatory.dff](#)), 194
[compute_dff_windowed_mode\(\)](#) (in module [al-lensdk.brain_observatory.dff](#)), 194
[compute_distance\(\)](#) (in module [al-lensdk.brain_observatory.receptive_field_analysis.fit_parameters](#)), 182
[compute_elliptical_area\(\)](#) (in module [al-lensdk.brain_observatory.behavior.eye_tracking_processing](#)), 138
[compute_expected\(\)](#) (in module [al-lensdk.brain_observatory.chisquare_categorical](#)), 190
[compute_frame_times\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.stimulus_sync](#)), 176
[compute_injection\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage](#) method), 366
[compute_injection\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage](#) method), 366
[compute_intensity\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage](#) method), 365
[compute_non_overlap_masks\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)), 285
[compute_non_overlap_traces\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)), 285
[compute_observed\(\)](#) (in module [al-lensdk.brain_observatory.chisquare_categorical](#)), 190
[compute_overlap\(\)](#) (in module [al-lensdk.brain_observatory.receptive_field_analysis.fit_parameters](#)), 182
[compute_projection\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage](#) method), 365
[compute_projection\(\)](#) ([al-lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage](#) method), 365
[compute_receptive_field\(\)](#) (in module [al-lensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis](#)), 182

compute_receptive_field_with_postprocessing() (in module al-
lensdk.brain_observatory.behavior.stimulus_processing.util),
183

compute_sum_pixels() (al-
lensdk.mouse_connectivity.grid.subimage.classic_subimage.
ClassicSubImage method), 365

compute_sum_pixels() (al-
lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage
method), 366

Config (class in allensdk.model.biophys_sim.config), 349

configure_library_method() (al-
lensdk.model.glif.glif_neuron.GlifNeuron
static method), 354

configure_method() (al-
lensdk.model.glif.glif_neuron.GlifNeuron
static method), 354

connect() (in module al-
lensdk.internal.core.lims_utilities), 298

ConnectedServices (class in al-
lensdk.api.queries.connected_services), 75

consistency_is_key() (in module al-
lensdk.brain_observatory.behavior.criteria),
135

consistent_behavior_within_session() (in module al-
lensdk.brain_observatory.behavior.criteria),
135

construct_rolling_performance_df() (in module al-
lensdk.brain_observatory.behavior.trials_processing),
151

construct_well_known_file_download_url() (allensdk.api.api.Api method), 103

consume_volume() (al-
lensdk.mouse_connectivity.grid.image_series_griddler.ImageSeriesGriddler
method), 369

contingent_trials() (in module al-
lensdk.brain_observatory.behavior.trial_masks),
149

ContinuousFile (class in al-
lensdk.brain_observatory.ecephys.file_io.continuous_file.
ContinuousFile), 163

conv() (in module al-
lensdk.mouse_connectivity.grid.utilities.downsampling_utils.
DownsamplingUtils), 366

convert_axis() (in module al-
lensdk.internal.mouse_connectivity.projection_thumbnail.projection_utilities),
337

convert_discrete_colormap() (in module al-
lensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities),
337

convert_filepath_caseinsensitive() (al-
lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage
(class in al-
lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage
method), 263

convert_frame() (in module al-
lensdk.internal.pipeline_modules.run_ophys_session_decomposition.
RunOphysSessionDecomposition), 295

convert_from_titan_linux() (in module al-
lensdk.internal.core.lims_utilities), 298

convert_type() (al-
lensdk.internal.morphology.morphology.Morphology
method), 326

convolve() (in module al-
lensdk.brain_observatory.receptive_field_analysis.utilities),
184

copy() (allensdk.internal.mouse_connectivity.projection_thumbnail.image_utilities.
ImageUtilities), 337

copy_local() (allensdk.internal.model.biophysical.run_optimize.RunOptimize.
RunOptimize), 310

copy_local() (allensdk.internal.model.biophysical.run_simulate_lims.
RunSimulateLims), 311

CoronaPlotter (class in al-
lensdk.brain_observatory.circle_plots), 190

correct_on_off_effects() (in module al-
lensdk.brain_observatory.ecephys.stimulus_sync),
177

corrected_behavior_video_timestamps (al-
lensdk.internal.brain_observatory.time_sync.OphysTimeAligner
attribute), 296

corrected_eye_video_timestamps (al-
lensdk.internal.brain_observatory.time_sync.OphysTimeAligner
attribute), 296

corrected_ophys_timestamps (al-
lensdk.internal.brain_observatory.time_sync.OphysTimeAligner
attribute), 296

corrected_video_timestamps (al-
lensdk.internal.brain_observatory.time_sync.OphysTimeAligner
attribute), 296

corrected_video_timestamps() (in module al-
lensdk.internal.brain_observatory.time_sync),
297

correlation_report() (in module al-
lensdk.internal.brain_observatory.demix_report),
285

count (allensdk.internal.brain_observatory.mask_set.MaskSet
function), 292

count_writer() (in module al-
lensdk.mouse_connectivity.grid.writers,
366

CountSubImage (class in al-
lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage),
263

366
 cr_position_in_mouse_eye_coordinates() (allensdk.internal.brain_observatory.eye_calibration. EyeCalibration static method), 287
 cre_line (allensdk.brain_observatory.behavior.metadata.behavior_metadata. BehaviorMetadata attribute), 109
 create_argparser() (allensdk.config.app.application_config.ApplicationConfig method), 222
 create_basis_IPSP() (in module allensdk.internal.model.GLM), 323
 create_dir() (allensdk.config.manifest.Manifest method), 228
 create_extended_trials() (in module allensdk.brain_observatory.behavior.trials_processing), 151
 create_fake_metadata() (in module allensdk.internal.pipeline_modules.run_ophys_sessions_by_decomposition), 345
 create_gene_fpkm_table() (in module allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap), 340
 create_genes_for_transcripts() (in module allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap), 340
 create_image() (in module allensdk.internal.morphology.morphvis), 329
 create_images() (allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream method), 290
 create_images() (allensdk.internal.brain_observatory.frame_stream.FrameInputStream method), 291
 create_manifest() (allensdk.api.queries.biophysical_api.BiophysicalApi method), 69
 create_neuropil_mask() (in module allensdk.brain_observatory.roi_masks), 204
 create_region_mask() (in module allensdk.internal.brain_observatory.annotated_region_metrics), 284
 create_roi_mask() (in module allensdk.brain_observatory.roi_masks), 205
 create_roi_mask_array() (in module allensdk.brain_observatory.roi_masks), 205
 create_sample_metadata() (in module allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap), 340
 create_stim_table() (in module allensdk.brain_observatory.ecephys.stimulus_table.ephys_stimulus_table), 170
 create_transcript_fpkm_table() (in module allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap), 340
 create_utils() (in module allensdk.model.biophysical.utils), 352
 credential_injector() (in module allensdk.core.authentication), 231
 CredentialProvider (class in allensdk.core.authentication), 231
 CRITERIA (allensdk.api.queries.rma_api.RmaApi attribute), 89
 CRITERIA() (in module allensdk.internal.brain_observatory.roi_filter_utils), 294
 CUT_DENDRITE (allensdk.core.swc.Marker attribute), 301
 CUT_DENDRITE (allensdk.internal.core.swc.Marker attribute), 301
 CvInputStream (class in allensdk.internal.brain_observatory.frame_stream), 290
 D
 data (allensdk.brain_observatory.behavior.image_api.Image attribute), 141
 data_file_attributes() (allensdk.api.cloud_cache.manifest.Manifest method), 66
 DATA_MASK (allensdk.api.queries.grid_data_api.GridDataApi attribute), 77
 DATA_MASK_KEY (allensdk.core.mouse_connectivity_cache.MouseConnectivityCache attribute), 242
 data_to_licks() (in module allensdk.brain_observatory.behavior.trials_processing), 151
 data_to_metadata() (in module allensdk.brain_observatory.behavior.trials_processing), 152
 data_transforms() (allensdk.core.structure_tree.StructureTree static method), 259
 data_query() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 70
 data_query_string() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 71
 DataFrameIndexError, 239

[DataFrameKeyError](#), 239
[dataset](#) (`allensdk.internal.brain_observatory.time_sync.OphysTimeAlignment.alignment_parser()` (in module `allensdk.internal.core.lims_pipeline_module`), attribute), 296
[Dataset](#) (class in `allensdk.brain_observatory.sync_dataset`), 217
[date_of_acquisition](#) (`allensdk.brain_observatory.behavior.metadata.behavior_metadata.BehaviorMetadata` (attribute), 109
[DatUtilities](#) (class in `allensdk.core.dat_utilities`), 239
[db_connection_creator\(\)](#) (in module `allensdk.internal.api`), 283
[DbCredentials](#) (class in `allensdk.core.authentication`), 231
[dbname](#) (`allensdk.core.authentication.DbCredentials` attribute), 231
[DEBUG](#) (`allensdk.api.queries.rma_api.RmaApi` attribute), 89
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_annotated_region_demixing`), 341
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_demixing`), 341
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_neuropil_correction`), 342
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_observatory_analysis_method`), 342
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnail_method`), 344
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_ophys_eye_calibration`), 345
[debug\(\)](#) (in module `allensdk.internal.pipeline_modules.run_ophys_session_decomposition`), 345
[debug_clause\(\)](#) (`allensdk.api.queries.rma_api.RmaApi` method), 90
[debug_plot\(\)](#) (in module `allensdk.internal.pipeline_modules.run_neuropil_correction`), 342
[decision_function\(\)](#) (`allensdk.internal.brain_observatory.roi_filter_utils.TrainingLabelClassifier.nwb_data_set.NwbDataSet` attribute), 294
[decode_bytes\(\)](#) (in module `allensdk.core.h5_utilities`), 240
[default\(\)](#) (`allensdk.brain_observatory.JSONEncoder` method), 221
[default_api_url](#) (`allensdk.api.api.Api` attribute), 103
[default_ray\(\)](#) (in module `allensdk.internal.brain_observatory.itracker_utils`), 298
[DEFAULT_STRUCTURE_SET_IDS](#) (`allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), 242
[DEFORMATION_FIELD_HEADER_KEY](#) (`allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), 242
[DEFORMATION_FIELD_VOXEL_KEY](#) (`allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), 242
[deg_to_dist\(\)](#) (in module `allensdk.brain_observatory.behavior.running_processing`), 143
[deinterpolate_RF\(\)](#) (in module `allensdk.brain_observatory.receptive_field_analysis.chisquarerf`), 179
[delay_metrics\(\)](#) (`allensdk.epphys.epphys_extractor.EphysSweepFeatureExtractor` method), 266
[delete_tree\(\)](#) (`allensdk.core.swc.Morphology` method), 263
[delete_tree\(\)](#) (`allensdk.internal.morphology.morphology.Morphology` method), 326
[demix_time_dep_masks\(\)](#) (in module `allensdk.brain_observatory.demixer`), 192
[demix_time_dep_masks\(\)](#) (in module `allensdk.internal.brain_observatory.demixer`), 285
[DENSITY](#) (`allensdk.api.queries.grid_data_api.GridDataApi` attribute), 77
[deprecated\(\)](#) (in module `allensdk.deprecated`), 371
[DEPRECATED_KEYS](#) (`allensdk.brain_observatory.sync_dataset.Dataset` attribute), 218
[DEPRECATED_SPIKE_TIMES](#) (`allensdk.brain_observatory.sync_dataset.Dataset` attribute), 218
[descendant_ids\(\)](#) (`allensdk.core.simple_tree.SimpleTree` method), 255
[descendant_ids\(\)](#) (`allensdk.internal.core.simpletree.SimpleTree` method), 255

`method`), 300
`descendants()` (`al-`
`lensdk.core.simple_tree.SimpleTree` `method`),
256
`descendants()` (`al-`
`lensdk.internal.core.simpletree.SimpleTree`
`method`), 300
`Description` (`class` `in` `al-`
`lensdk.config.model.description`), 226
`DescriptionParser` (`class` `in` `al-`
`lensdk.config.model.description_parser`),
227
`deserialize()` (`al-`
`lensdk.brain_observatory.behavior.image_api.ImageApi`
`static method`), 141
`detect_bursts()` (`in` `module` `al-`
`lensdk.ephys.ephys_features`), 271
`detect_duplicates()` (`al-`
`lensdk.internal.brain_observatory.mask_set.MaskSet`
`method`), 292
`detect_events()` (`in` `module` `al-`
`lensdk.brain_observatory.receptive_field_analysis.eventdetector`
181
`detect_pauses()` (`in` `module` `al-`
`lensdk.ephys.ephys_features`), 272
`detect_putative_spikes()` (`in` `module` `al-`
`lensdk.ephys.ephys_features`), 272
`detect_unions()` (`al-`
`lensdk.internal.brain_observatory.mask_set.MaskSet`
`method`), 292
`determine_likely_blinks()` (`in` `module` `al-`
`lensdk.brain_observatory.behavior.eye_tracking_processing`
138
`determine_outliers()` (`in` `module` `al-`
`lensdk.brain_observatory.behavior.eye_tracking_processing`
139
`DEVMOUSE_2012` (`al-`
`lensdk.api.queries.reference_space_api.ReferenceSpaceApi`
`attribute`), 87
`DEVMOUSE_ATLAS_PRODUCTS` (`al-`
`lensdk.api.queries.mouse_atlas_api.MouseAtlasApi`
`attribute`), 81
`df_columns` (`allensdk.config.manifest_builder.ManifestBuilder`
`attribute`), 229
`dfftraces` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`attribute`), 212
`DFMFLD_RESOLUTIONS` (`al-`
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache`
`attribute`), 242
`dict_generator()` (`in` `module` `al-`
`lensdk.brain_observatory.receptive_field_analysis.tools`),
184
`dict_to_indexed_array()` (`in` `module` `al-`
`lensdk.brain_observatory`), 221
`dim_color()` (`allensdk.brain_observatory.observatory_plots.Dimension`
`method`), 200
`DimensionPatchHandler` (`class` `in` `al-`
`lensdk.brain_observatory.observatory_plots`),
200
`DIR` (`allensdk.config.manifest.Manifest` `attribute`), 227
`DIR_CCW` (`allensdk.brain_observatory.circle_plots.PolarPlotter`
`attribute`), 191
`DIR_CW` (`allensdk.brain_observatory.circle_plots.PolarPlotter`
`attribute`), 191
`direct_sum_projection_pixels` (`al-`
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
`attribute`), 334
`direct_unionize()` (`al-`
`lensdk.internal.mouse_connectivity.interval_unionize.interval_un`
`method`), 332
`direct_voxel_counts()` (`al-`
`lensdk.core.reference_space.ReferenceSpace`
`method`), 250
`direct_voxel_map` (`al-`
`lensdk.core.reference_space.ReferenceSpace`
`attribute`), 250
`DIRNAME` (`allensdk.config.manifest.Manifest` `attribute`),
227
`distance()` (`allensdk.internal.brain_observatory.mask_set.MaskSet`
`method`), 292
`do_blur()` (`in` `module` `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.generat`
336
`do_query()` (`allensdk.api.api.Api` `method`), 103
`do_rma_query()` (`allensdk.api.api.Api` `method`), 103
`download_alignment3d()` (`al-`
`lensdk.api.queries.grid_data_api.GridDataApi`
`method`), 77
`download_annotation_volume()` (`al-`
`lensdk.api.queries.reference_space_api.ReferenceSpaceApi`
`method`), 88
`download_atlas_image()` (`al-`
`lensdk.api.queries.image_download_api.ImageDownloadApi`
`method`), 79
`download_data_mask()` (`al-`
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityAp`
`method`), 83
`download_data_mask()` (`al-`
`lensdk.internal.api.queries.mouse_connectivity_api_prerelease.M`
`method`), 280
`download_deformation_field()` (`al-`
`lensdk.api.queries.grid_data_api.GridDataApi`
`method`), 77
`download_expression_density()` (`al-`
`lensdk.api.queries.mouse_atlas_api.MouseAtlasApi`
`method`), 82
`download_expression_energy()` (`al-`
`lensdk.api.queries.mouse_atlas_api.MouseAtlasApi`

method), 82

download_expression_grid_data() (allensdk.api.queries.grid_data_api.GridDataApi method), 78

download_expression_intensity() (allensdk.api.queries.mouse_atlas_api.MouseAtlasApi method), 82

download_gene_expression_grid_data() (allensdk.api.queries.grid_data_api.GridDataApi method), 78

download_image() (allensdk.api.queries.image_download_api.ImageDownloadApi method), 79

download_injection_density() (allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi method), 83

download_injection_density() (allensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease method), 280

download_injection_fraction() (allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi method), 83

download_injection_fraction() (allensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease method), 280

download_mouse_atlas_volume() (allensdk.api.queries.reference_space_api.ReferenceSpaceApi method), 88

download_projection_density() (allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi method), 83

download_projection_density() (allensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease method), 280

download_projection_grid_data() (allensdk.api.queries.grid_data_api.GridDataApi method), 78

download_projection_grid_data() (allensdk.internal.api.queries.grid_data_api_prerelease.GridDataApiPrerelease method), 280

download_projection_image() (allensdk.api.queries.image_download_api.ImageDownloadApi method), 80

download_reference_aligned_image_channel_volume() (allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi method), 83

download_section_image() (allensdk.api.queries.image_download_api.ImageDownloadApi method), 80

download_structure_mask() (allensdk.api.queries.reference_space_api.ReferenceSpaceApi method), 88

download_structure_mesh() (allensdk.api.queries.reference_space_api.ReferenceSpaceApi method), 88

download_svg() (allensdk.api.queries.svg_api.SvgApi method), 94

download_template_volume() (allensdk.api.queries.reference_space_api.ReferenceSpaceApi method), 88

download_url (allensdk.api.api.Api attribute), 104

download_volumetric_data() (allensdk.api.queries.reference_space_api.ReferenceSpaceApi method), 88

download_volume() (allensdk.core.reference_space.ReferenceSpace method), 250

downsample_average() (in module allensdk.api.mouse_connectivity.grid.utilities.downsampling_utilities), 367

draw_density_hist() (in module allensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease), 329

draw_morphology() (in module allensdk.internal.morphology.morphvis), 330

driftings (class in allensdk.brain_observatory.drifting_gratings), 195

drop_empty_columns() (in module allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities), 172

duty_cycle() (allensdk.brain_observatory.sync_dataset.Dataset method), 218

dxcm (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis attribute), 212

dxcm (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis attribute), 212

dynamics() (allensdk.model.glif.glif_neuron.GlifNeuron method), 354

dynamics_ACurrent_exp() (in module allensdk.model.glif.glif_neuron_methods), 356

dynamics_ACurrent_exp() (in module allensdk.model.glif.glif_neuron_methods), 356

dynamics_ACurrent_exp() (in module allensdk.model.glif.glif_neuron_methods), 356

dynamics_threshold_inf() (in module allensdk.model.glif.glif_neuron_methods), 356

dynamics_threshold_spike_component() (in module allensdk.model.glif.glif_neuron_methods), 357

dynamics_threshold_three_components_exact() (in module allensdk.model.glif.glif_neuron_methods), 357

dynamics_voltage_linear_exact() (in module allensdk.model.glif.glif_neuron_methods), 358

dynamics_voltage_linear_forward_euler() (in module allensdk.model.glif.glif_neuron_methods), 358

(in module *allensdk.model.glif.glif_neuron_methods*), 358

E

eccentricity() (in module *allensdk.internal.brain_observatory.annotated_region_methods*), 284

eccentricity() (in module *allensdk.internal.brain_observatory.itracker_utils*), 291

EcephysProjectApi (class in *allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api*), 161

EcephysSessionApi (class in *allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api*), 162

EcephysSyncDataset (class in *allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset*), 163

ellipse_angle_of_rotation() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 290

ellipse_angle_of_rotation2() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 290

ellipse_axis_length() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 290

ellipse_center() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 290

emission_lambda (*allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata* attribute), 111

enable_console_log() (in module *allensdk.config*), 230

ENERGY (*allensdk.api.queries.grid_data_api.GridDataApi* attribute), 77

EnvCredentialProvider (class in *allensdk.core.authentication*), 231

EPHYS_DATA_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 237

EPHYS_FEATURES_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 237

EPHYS_SWEEPS_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 237

EphysCellFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 265

EphysFeatureExtractor (class in *allensdk.ephys.feature_extractor*), 277

EphysFeatures (class in *allensdk.ephys.feature_extractor*), 277

EphysSweepFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 266

EphysSweepSetFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 268

EpochSeq (class in *allensdk.internal.model.separation_exception*), 189

EQ (*allensdk.api.queries.rma_api.RmaApi* attribute), 89

equipment_name (*allensdk.brain_observatory.behavior.metadata.behavior_metadata* attribute), 109

error_calc() (in module *allensdk.brain_observatory.r_neuropil*), 202

error_calc_outlier() (in module *allensdk.brain_observatory.r_neuropil*), 202

escape_ophys_session() (*allensdk.internal.model.biophysical.passive_fitting.output* attribute), 308

estimate_adjusted_detection_parameters() (in module *allensdk.ephys.ephys_features*), 273

estimate_contamination_ratios() (in module *allensdk.brain_observatory.r_neuropil*), 202

estimate_error() (*allensdk.brain_observatory.r_neuropil.NeuropilSubtract* method), 202

estimate_fi_shift() (in module *allensdk.internal.model.biophysical.check_fi_shift*), 309

estimate_frame_duration() (in module *allensdk.brain_observatory.ecephys.stimulus_sync*), 177

estimate_sag() (*allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor* method), 266

estimate_time_constant() (*allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor* method), 266

euclidean_distance() (in module *allensdk.internal.morphology.node*), 331

evaluate() (*allensdk.internal.model.glif.glif_optimizer.GlifOptimizer* method), 315

events_to_pvalues_no_fdr_correction() (in module *allensdk.brain_observatory.receptive_field_analysis.receptive_field*), 183

EXCEPT (*allensdk.api.queries.rma_api.RmaApi* attribute), 89

excitation_lambda (*allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata* attribute), 111

EXCPT (*allensdk.api.queries.rma_api.RmaApi* attribute), 89

exp_curve() (in module *allensdk.internal.ephys.plot_qc_figures*), 302

`exp_curve()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 304
`exp_decay()` (in module `allensdk.internal.model.glif.MLIN`), 311
`exp_fit_c()` (in module `allensdk.internal.model.glif.threshold_adaptation`), 320
`exp_force_c()` (in module `allensdk.internal.model.glif.threshold_adaptation`), 320
`experiment_container_id` (al- `allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata.align_timestamps.barcode_syn` attribute), 111
`experiment_correlation_search()` (al- `allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 83
`experiment_id` (al- `allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs` attribute), 346
`experiment_injection_coordinate_search()` (al- `allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 83
`experiment_source_search()` (al- `allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 84
`experiment_spatial_search()` (al- `allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 84
`ExperimentGeometry` (class in al- `allensdk.brain_observatory.stimulus_info`), 214
`EXPERIMENTS_KEY` (al- `allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), 242
`EXPERIMENTS_PRERELEASE_KEY` (al- `allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease` attribute), 300
`ExperimentsTable` (class in al- `allensdk.brain_observatory.behavior.behavior_project_cache.tables.ExperimentsTable`), 107
`explode_response_window()` (in module al- `allensdk.brain_observatory.behavior.mtrain`), 142
`export_frame_to_hdf5()` (in module al- `allensdk.internal.brain_observatory.ophys_session_decomposition`), 293
`export_itksnap_labels()` (al- `allensdk.core.reference_space.ReferenceSpace` method), 250
`export_label_description()` (al- `allensdk.core.structure_tree.StructureTree` method), 259
`expsymm_cdf()` (in module al- `allensdk.internal.model.glif.MLIN`), 311
`expsymm_pdf()` (in module al- `allensdk.internal.model.glif.MLIN`), 311
`ExtendedTrialSchema` (class in al- `allensdk.brain_observatory.behavior.mtrain`), 141
`extract()` (al- `allensdk.internal.mouse_connectivity.projection_thumbnail.v` method), 338
`extract()` (in module al- `allensdk.mouse_connectivity.grid.utilities.downsampling_utilities`), 367
`extract_barcodes()` (al- `allensdk.brain_observatory.ecephys.align_timestamps.barcode_syn` method), 158
`extract_barcodes_from_states()` (in module al- `allensdk.brain_observatory.ecephys.align_timestamps.channel_stat`), 159
`extract_cell_features()` (in module al- `allensdk.ephys.extract_cell_features`), 276
`extract_params_from_stim_repr()` (in module al- `allensdk.brain_observatory.ecephys.stimulus_table.stimulus_param`), 156
`extract_data()` (al- `allensdk.internal.mouse_connectivity.interval_unionize.interval_un` method), 332
`extract_data()` (al- `allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u` method), 335
`extract_data()` (in module al- `allensdk.internal.ephys.core_feature_extract`), 367
`extract_frame_times()` (al- `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E` method), 163
`extract_frame_times_from_photodiode()` (al- `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E` method), 163
`extract_frame_times_from_vsyncs()` (al- `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E` method), 163
`extract_injection_from_segmentation()` (al- `allensdk.mouse_connectivity.grid.subimage.base_subimage.Segmen` method), 363
`extract_led_times()` (al- `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E` method), 163
`extract_signal_from_segmentation()` (al- `allensdk.mouse_connectivity.grid.subimage.base_subimage.Segmen` method), 363

method), 363
 extract_splits_from_states() (in module `allensdk.brain_observatory.ecephys.align_timestamps.channel_states`), 159
 extract_stim_class_from_repr() (in module `allensdk.brain_observatory.ecephys.stimulus_table.stimulus_table`), 174
 extract_sweep_features() (in module `allensdk.ephys.extract_cell_features`), 277
 extractor_for_nwb_sweeps() (in module `allensdk.ephys.ephys_extractor`), 269
 extralength(`allensdk.brain_observatory.locally_sparse_fids`), 197
 extralength(`allensdk.brain_observatory.natural_scenes.NaturalScenes`), 199
 extralength(`allensdk.brain_observatory.static_gratings.StaticGratings`), 210
 extrapolate_model_spike_from_endpoints() (in module `allensdk.internal.model.glif.glif_optimizer_neuron`), 317
 extrapolate_model_spike_from_endpoints_single_buffer() (in module `allensdk.internal.model.glif.glif_optimizer_neuron`), 317
 extrapolate_spike_time() (in module `allensdk.internal.model.glif.glif_optimizer_neuron`), 317
 extrapolate_spike_voltage() (in module `allensdk.internal.model.glif.glif_optimizer_neuron`), 317
 eye_alignment (in module `allensdk.internal.pipeline_modules.run_ophys_time_sync`), 346
 eye_delta(`allensdk.internal.pipeline_modules.run_ophys_time_sync`), 346
 eye_times(`allensdk.internal.pipeline_modules.run_ophys_time_sync`), 346
 EYE_TRACKING_KEYS (in module `allensdk.brain_observatory.sync_dataset.Dataset`), 218
 eye_video_timestamps (in module `allensdk.internal.brain_observatory.time_sync.OphysTimeAlignment`), 296
 EyeCalibration (class in `allensdk.internal.brain_observatory.eye_calibration`), 287
 EyeTrackingRigGeometry (class in `allensdk.brain_observatory.behavior.schemas`), 146
 F
 factory() (`allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset`), 163
 factory() (`allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimFile`), 164
 FALSE (`allensdk.api.queries.rma_api.RmaApi` attribute), 89
 FanPlotter (class in `allensdk.brain_observatory.circle_plots`), 191
 FeatureError, 269
 fetchall() (`allensdk.internal.api.PostgresQueryMixin` method), 283
 fetchone() (`allensdk.internal.api.PostgresQueryMixin` method), 283
 fix_id(`allensdk.internal.model.GLM`), 323
 FfmpegInputStream (class in `allensdk.internal.brain_observatory.frame_stream`), 290
 FfmpegOutputStream (class in `allensdk.internal.brain_observatory.frame_stream`), 290
 field_of_view_height (in `allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata`), 111
 field_of_view_width (in `allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata`), 111
 figure_in_px() (in module `allensdk.brain_observatory.observatory_plots`), 200
 FILE (`allensdk.config.manifest.Manifest` attribute), 227
 file_hash(`allensdk.api.cloud_cache.file_attributes.CacheFileAttributes` attribute), 66
 file_hash_from_path() (in module `allensdk.api.cloud_cache.utils`), 67
 file_id_values (in module `allensdk.api.cloud_cache.manifest.Manifest` attribute), 66
 FILE_METADATA_MAPPING (in module `allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet`), 232
 fill_sweep_responses() (in module `allensdk.core.nwb_data_set.NwbDataSet`), 247
 filter() (`allensdk.api.queries.rma_api.RmaApi` method), 90
 filter_bad_params() (in module `allensdk.internal.brain_observatory.itracker_utils`), 291
 filter_cell_specimens() (in module `allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 71
 factory() (`allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset`), 163

<code>lensdk.api.queries.cell_types_api.CellTypesApi</code> <code>method</code>), 74	<code>find()</code> (<code>allensdk.internal.morphology.morphology.Morphology</code> <code>method</code>), 326
<code>filter_cells_api()</code> (<code>allensdk.api.queries.cell_types_api.CellTypesApi</code> <code>method</code>), 74	<code>find_bin_center()</code> (<code>in module allensdk.internal.model.glif.MLIN</code>), 311
<code>filter_digital()</code> (<code>in module allensdk.brain_observatory.behavior.sync.process_sync</code>), 130	<code>find_coarse_long_square_amp_delta()</code> (<code>in module allensdk.internal.ephys.core.feature_extract</code>), 301
<code>filter_events_array()</code> (<code>in module allensdk.brain_observatory.behavior.event_detection</code>), 137	<code>find_container_tags()</code> (<code>in module allensdk.api.queries.brain_observatory_api</code>), 73
<code>filter_experiment_containers()</code> (<code>allensdk.api.queries.brain_observatory_api.BrainObservatoryApi</code> <code>method</code>), 71	<code>find_downstroke_indexes()</code> (<code>in module allensdk.ephys.ephys_features</code>), 274
<code>filter_experiments()</code> (<code>allensdk.core.mouse_connectivity_cache.MouseConnectivityCache</code> <code>method</code>), 242	<code>find_experiment_acquisition_age()</code> (<code>in module allensdk.api.queries.brain_observatory_api</code>), 73
<code>filter_experiments()</code> (<code>allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease</code> <code>method</code>), 300	<code>find_first_model_spike()</code> (<code>in module allensdk.interactive_calglif.PyLifeSpanAnalyzer_neuron</code>), 318
<code>filter_experiments_and_containers()</code> (<code>allensdk.api.queries.brain_observatory_api.BrainObservatoryApi</code> <code>method</code>), 71	<code>find_licks()</code> (<code>in module allensdk.brain_observatory.behavior.trials_processing</code>), 152
<code>filter_nodes()</code> (<code>allensdk.core.simple_tree.SimpleTree</code> <code>method</code>), 256	<code>find_long_square_sweeps()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
<code>filter_ophys_experiments()</code> (<code>allensdk.api.queries.brain_observatory_api.BrainObservatoryApi</code> <code>method</code>), 71	<code>find_matching_index()</code> (<code>in module allensdk.brain_observatory.ecephys.align_timestamps.barcode</code>), 156
<code>filter_putative_spikes()</code> (<code>in module allensdk.ephys.ephys_features</code>), 273	<code>find_negative_baselines()</code> (<code>in module allensdk.brain_observatory.demixer</code>), 193
<code>filter_structure_unionizes()</code> (<code>allensdk.core.mouse_connectivity_cache.MouseConnectivityCache</code> <code>method</code>), 243	<code>find_negative_baselines()</code> (<code>in module allensdk.internal.brain_observatory.demixer</code>), 285
<code>filter_sweeps()</code> (<code>in module allensdk.internal.ephys.core.feature_extract</code>), 301	<code>find_negative_transients_threshold()</code> (<code>in module allensdk.brain_observatory.demixer</code>), 193
<code>filtered_sweep_numbers()</code> (<code>in module allensdk.internal.ephys.core.feature_extract</code>), 301	<code>find_negative_transients_threshold()</code> (<code>in module allensdk.internal.brain_observatory.demixer</code>), 285
<code>filters()</code> (<code>allensdk.api.queries.rma_api.RmaApi</code> <code>method</code>), 90	<code>find_noise_sweeps()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
<code>finalize()</code> (<code>allensdk.brain_observatory.circle_plots.PolarPlot</code> <code>method</code>), 191	<code>find_peak_indexes()</code> (<code>in module allensdk.ephys.ephys_features</code>), 274
<code>finalize_no_axes()</code> (<code>in module allensdk.brain_observatory.observatory_plots</code>), 200	<code>find_ramp_sweeps()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
<code>finalize_no_labels()</code> (<code>in module allensdk.brain_observatory.observatory_plots</code>), 200	<code>find_ramp_to_rheo_sweeps()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
<code>finalize_with_axes()</code> (<code>in module allensdk.brain_observatory.observatory_plots</code>), 200	<code>find_ranked_sweep()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
<code>find()</code> (<code>allensdk.core.swc.Morphology</code> <code>method</code>), 263	<code>find_short_square_sweeps()</code> (<code>in module allensdk.internal.model.glif.find_sweeps</code>), 312
	<code>find_specimen_cre_line()</code> (<code>in module allensdk.api.queries.brain_observatory_api</code>),

73
[find_specimen_reporter_line\(\)](#) (in module [allensdk.api.queries.brain_observatory_api](#)), 73
[find_specimen_transgenic_lines\(\)](#) (in module [allensdk.api.queries.brain_observatory_api](#)), 73
[find_spikes_list\(\)](#) (in module [allensdk.internal.model.glif.find_spikes](#)), 312
[find_spikes_list_old\(\)](#) (in module [allensdk.internal.model.glif.find_spikes](#)), 312
[find_spikes_old\(\)](#) (in module [allensdk.internal.model.glif.find_spikes](#)), 312
[find_spikes_ssq_list\(\)](#) (in module [allensdk.internal.model.glif.find_spikes](#)), 312
[find_stim_start\(\)](#) (in module [allensdk.internal.ephys.core.feature_extract](#)), 302
[find_sweep_stim_start\(\)](#) (in module [allensdk.internal.ephys.core.feature_extract](#)), 302
[find_sweeps\(\)](#) (in module [allensdk.internal.model.glif.find_sweeps](#)), 313
[find_time_index\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 274
[find_trough_indexes\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 274
[find_upstroke_indexes\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 274
[find_widths\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 275
[find_zero_baselines\(\)](#) (in module [allensdk.brain_observatory.demixer](#)), 193
[find_zero_baselines\(\)](#) (in module [allensdk.internal.brain_observatory.demixer](#)), 286
[findlevel\(\)](#) (in module [allensdk.brain_observatory.findlevel](#)), 197
[fit\(\)](#) ([allensdk.brain_observatory.r_neuropil.NeuropilSubtract](#) method), 202
[fit_avoltage_bvoltage\(\)](#) (in module [allensdk.internal.model.glif.threshold_adaptation](#)), 320
[fit_avoltage_bvoltage_th\(\)](#) (in module [allensdk.internal.model.glif.threshold_adaptation](#)), 321
[fit_block_coordinate_desc\(\)](#) ([allensdk.brain_observatory.r_neuropil.NeuropilSubtract](#) method), 202
[fit_ellipse\(\)](#) ([allensdk.internal.brain_observatory.fit_ellipse.FitEllipse](#) method), 290
[fit_ellipse\(\)](#) (in module [allensdk.internal.brain_observatory.fit_ellipse](#)), 290
[fit_fit_slope\(\)](#) (in module [allensdk.ephys.ephys_extractor](#)), 269
[fit_membrane_time_constant\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 275
[fit_parameters_file_entries\(\)](#) ([allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 278
[fit_parameters_path\(\)](#) ([allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 278
[fit_prespike_time_constant\(\)](#) (in module [allensdk.ephys.ephys_features](#)), 275
[FitEllipse](#) (class in [allensdk.internal.brain_observatory.fit_ellipse](#)), 290
[fitgaussian2D\(\)](#) (in module [allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D](#)), 182
[fix_array_dimensions\(\)](#) (in module [allensdk.core.sitk_utilities](#)), 257
[fix_change_time\(\)](#) (in module [allensdk.brain_observatory.behavior.mtrain](#)), 142
[fix_unary_sections\(\)](#) ([allensdk.config.model.description.Description](#) method), 226
[fix_unexpected_edges\(\)](#) (in module [allensdk.brain_observatory.ecephys.stimulus_sync](#)), 177
[fixed_factory\(\)](#) ([allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projection_thumbnail](#) class method), 338
[flag_unexpected_edges\(\)](#) (in module [allensdk.brain_observatory.ecephys.stimulus_sync](#)), 177
[FLOAT](#) ([allensdk.api.queries.connected_services.ConnectedServices](#) attribute), 76
[float_label\(\)](#) (in module [allensdk.brain_observatory.observatory_plots](#)), 200
[for_drifting_gratings\(\)](#) ([allensdk.brain_observatory.circle_plots.FanPlotter](#) static method), 191
[for_static_gratings\(\)](#) ([allensdk.brain_observatory.circle_plots.FanPlotter](#) static method), 191
[FRAME_KEYS](#) ([allensdk.brain_observatory.sync_dataset.Dataset](#) attribute), 218
[FrameInputStream](#) (class in [allensdk.internal.brain_observatory.frame_stream](#)), 291
[FrameOutputStream](#) (class in [allensdk.internal.brain_observatory.frame_stream](#)), 291

`lensdk.internal.brain_observatory.frame_stream`), `from_processed()` (al-
291 `lensdk.brain_observatory.behavior.stimulus_processing.stimulus_`
`frames_per_second` (al- `static method`), 126
`lensdk.brain_observatory.ecephys.file_io.stim_file.CompatiblePickleStimFile` (al-
`attribute`), 164 `lensdk.brain_observatory.behavior.stimulus_processing.stimulus_`
`frequency()` (allensdk.brain_observatory.sync_dataset.Dataset `static method`), 126
`method`), 218 `from_sweeps()` (al-
`FriendlyDate` (class in al- `lensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor`
`lensdk.brain_observatory.behavior.mtrain`), `class method`), 268
141 `from_unprocessed()` (al-
`FriendlyDateTime` (class in al- `lensdk.brain_observatory.behavior.stimulus_processing.stimulus_`
`lensdk.brain_observatory.behavior.mtrain`), `method`), 126
141 `from_unprocessed()` (al-
`from_analysis_file()` (al- `lensdk.brain_observatory.behavior.stimulus_processing.stimulus_`
`lensdk.brain_observatory.drifting_gratings.DriftingGratings` `static method`), 127
`static method`), 195 `full_genotype` (al-
`from_analysis_file()` (al- `lensdk.brain_observatory.behavior.metadata.behavior_metadata.`
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` `static method`), 109
`static method`), 197
`from_analysis_file()` (al- **G**
`lensdk.brain_observatory.natural_movie.NaturalMovie` `gather_from_seeds()` (al-
`static method`), 198 `lensdk.internal.model.biophysical.make_deap_fit_json.Report`
`from_analysis_file()` (al- `method`), 310
`lensdk.brain_observatory.natural_scenes.NaturalScenes` `Gaussian2D()` (in module al-
`static method`), 199 `lensdk.brain_observatory.receptive_field_analysis.fitgaussian2D`)
`from_analysis_file()` (al- 182
`lensdk.brain_observatory.static_gratings.StaticGratings` `GaussianFitError`, 182
`static method`), 210 `generate_fit_file()` (al-
`from_dataframe()` (al- `lensdk.internal.model.biophysical.make_deap_fit_json.Report`
`lensdk.config.manifest_builder.ManifestBuilder` `method`), 310
`method`), 229 `generate_manifest_lims()` (al-
`from_df()` (allensdk.brain_observatory.stimulus_info.BinaryIntervalSearchTree `lensdk.internal.model.biophysical.run_optimize.RunOptimize`
`static method`), 213 `method`), 310
`from_dict()` (allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron `generate_manifest_lims()` (al-
`class method`), 315 `lensdk.internal.model.biophysical.run_simulate_lims.RunSimulate`
`from_dict()` (allensdk.internal.morphology.node.Node `method`), 311
`class method`), 331 `generate_manifest_rma()` (al-
`from_dict()` (allensdk.model.glif.glif_neuron.GlifNeuron `lensdk.internal.model.biophysical.run_optimize.RunOptimize`
`class method`), 354 `method`), 310
`from_dict_legacy()` (al- `generate_manifest_rma()` (al-
`lensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron` `lensdk.internal.model.biophysical.run_simulate_lims.RunSimulate`
`class method`), 316 `method`), 311
`from_file_name()` (al- `generate_morphology()` (al-
`lensdk.internal.api.queries.grid_data_api_prerelease.GridDataApiPrerelease` `lensdk.internal.model.biophysical.deap_utils.Utils`
`class method`), 280 `method`), 309
`from_foraging_id()` (al- `generate_morphology()` (al-
`lensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api.BehaviorLimsExtractor` `lensdk.internal.model.biophysical.deap_utils.Utils`
`class method`), 120 `method`), 351
`from_json_file()` (al- `generate_morphology()` (al-
`lensdk.config.app.application_config.ApplicationConfig` `lensdk.model.biophysical.utils.Utils` `method`),
`method`), 223 351
`from_json_string()` (al- `generate_output_cell_features()`
`lensdk.config.app.application_config.ApplicationConfig` (in module al-
`method`), 223 `lensdk.internal.ephys.core_feature_extract`),

302
generate_rays() (in module al-
lensdk.internal.brain_observatory.itracker_utils),
291
generate_warp_coordinates() (al-
lensdk.brain_observatory.stimulus_info.ExperimentGeometry), 214
get_A() (in module al-
lensdk.brain_observatory.receptive_field_analysis.utilities),
184
get_A_blur() (in module al-
lensdk.brain_observatory.receptive_field_analysis.utilities),
184
get_affine_parameters() (al-
lensdk.core.mouse_connectivity_cache.MouseConnectivityCache),
method), 243
get_age() (allensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorDataExtractorBase), 112
get_age() (allensdk.brain_observatory.behavior.session_apis.data_transforms.json_api.BehaviorJsonExtractor), 118
get_age() (allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_base.BehaviorSessionBase), 120
get_alignment_array() (in module al-
lensdk.internal.brain_observatory.time_sync),
297
get_all_bits() (al-
lensdk.brain_observatory.sync_dataset.Dataset),
method), 218
get_all_events() (al-
lensdk.brain_observatory.sync_dataset.Dataset),
method), 218
get_all_features() (al-
lensdk.core.cell_types_cache.CellTypesCache),
method), 237
get_all_times() (al-
lensdk.brain_observatory.sync_dataset.Dataset),
method), 218
get_analog_channel() (al-
lensdk.brain_observatory.sync_dataset.Dataset),
method), 218
get_analog_meta() (al-
lensdk.brain_observatory.sync_dataset.Dataset),
method), 218
get_ancestor_id_map() (al-
lensdk.core.structure_tree.StructureTree),
method), 259
get_ancestor_id_map() (in module al-
lensdk.internal.mouse_connectivity.interval_unionize.run_tissot_classifier),
334
get_annotated_section_data_sets() (al-
lensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi),
method), 68
get_annotated_section_data_sets_via_rma() (al-
lensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi),
method), 68
get_annotation_volume() (al-
lensdk.core.reference_space_cache.ReferenceSpaceCache),
method), 253
get_atlases() (al-
lensdk.api.queries.ontologies_api.OntologiesApi),
method), 86
get_atlases_table() (al-
lensdk.api.queries.ontologies_api.OntologiesApi),
method), 86
get_attribute_dict() (in module al-
lensdk.brain_observatory.receptive_field_analysis.receptive_field_utilities),
184
get_attribute_dict() (in module al-
lensdk.brain_observatory.receptive_field_analysis.utilities),
184
get_average_projection() (al-
lensdk.brain_observatory.behavior.session_apis.data_transforms.average_projection_api.AverageProjectionApi),
method), 123
get_barcode_table() (al-
lensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_api.BarcodeSyncApi),
method), 159
get_behavior_monitoring() (in module al-
lensdk.brain_observatory.behavior.sync),
131
get_behavior_session_id() (al-
lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorDataExtractorBase),
method), 112
get_behavior_session_id() (al-
lensdk.brain_observatory.behavior.session_apis.abcs.session_base.BehaviorSessionBase),
method), 114
get_behavior_session_id() (al-
lensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_base.BehaviorSessionBase),
method), 118
get_behavior_session_id() (al-
lensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_base.BehaviorSessionBase),
method), 120
get_behavior_session_id() (al-
lensdk.brain_observatory.behavior.session_apis.data_transforms.average_projection_api.AverageProjectionApi),
method), 121
get_behavior_stimulus_file() (al-
lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorDataExtractorBase),
method), 112
get_behavior_stimulus_file() (al-
lensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_base.BehaviorSessionBase),
method), 118
get_behavior_stimulus_file() (al-
lensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_base.BehaviorSessionBase),
method), 118

method), 120

get_behavior_tracking_video_filepath_df() (allensdk.internal.api.lims_api.LimsApi method), 282

get_birth_date() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 120

get_bit() (allensdk.brain_observatory.sync_dataset.Dataset method), 219

get_bit() (in module allensdk.brain_observatory.sync_dataset), 221

get_bit_changes() (allensdk.brain_observatory.sync_dataset.Dataset method), 219

get_blend() (in module allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher), 338

get_blend_component() (in module allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher), 339

get_blocks() (in module allensdk.brain_observatory.ecephys.stimulus_table.visualization.nwb_blocks), 167

get_cache_path() (allensdk.api.warehouse_cache.cache.Cache method), 98

get_cap_check_indices() (in module allensdk.internal.model.biophysical.passive_fitting.preprocessor), 308

get_catch_responses() (in module allensdk.brain_observatory.behavior.dprime), 136

get_cav_density() (in module allensdk.internal.mouse_connectivity.interval_unionize.data_utilities), 332

get_cell() (allensdk.api.queries.cell_types_api.CellTypesApi method), 74

get_cell_metrics() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 71

get_cell_roi_ids() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 123

get_cell_specimen_id_mapping() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 71

get_cell_specimen_ids() (allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet method), 232

get_cell_specimen_indices() (allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet method), 232

get_cell_specimen_table() (allensdk.brain_observatory.behavior.session_apis.abcs.session_base method), 116

get_cell_specimen_table() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 123

get_cells() (allensdk.internal.model.biophysical.biophysical_archiver method), 309

get_change_time_frame_response_latency() (in module allensdk.brain_observatory.behavior.trials_processing), 152

get_channels() (allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api method), 161

get_channels() (allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api method), 162

get_children() (allensdk.core.ontology.Ontology method), 249

get_colormap() (allensdk.core.structure_tree.StructureTree method), 259

get_column_definitions() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 72

get_components() (in module allensdk.brain_observatory.receptive_field_analysis.utilities), 184

get_compound_annotated_section_data_sets() (allensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi method), 68

get_config() (allensdk.config.manifest_builder.ManifestBuilder method), 229

get_connection() (allensdk.internal.api.PostgresQueryMixin method), 283

get_corrected_fluorescence_traces() (allensdk.brain_observatory.behavior.session_apis.abcs.session_base method), 116

get_corrected_fluorescence_traces() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 123

get_corrected_fluorescence_traces() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 123

get_corrected_fluorescence_traces() (allensdk.brain_observatory.behavior.session_apis.data_transforms.BehaviorOphysDataTransform method), 123

get_credential_provider() (in module allensdk.core.authentication), 232

get_data_mask() (allensdk.internal.api.PostgresQueryMixin method), 283

`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache_line()` (al-
 method), 243 `lensdk.brain_observatory.behavior.session_apis.data_io.behavior`
`get_date_of_acquisition()` (al- method), 118
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_data_extractor_base` (al-
 method), 112 `lensdk.brain_observatory.behavior.session_apis.data_io.behavior`
`get_date_of_acquisition()` (al- method), 120
`lensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.BehaviorJsonExtractor` (al-
 method), 118 `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_`
`get_date_of_acquisition()` (al- method), 162
`lensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.BehaviorJsonExtractor` (al-
 method), 120 `lensdk.brain_observatory.nwb_data_set.BrainObservatoryNwbDataSet`
`get_default_manifest_file()` (in module al- `get_ephys_data()` (al-
`lensdk.api.warehouse_cache.cache`), 100 `lensdk.core.cell_types_cache.CellTypesCache`
`get_deformation_field()` (al- method), 238
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` (al-
 method), 244 `lensdk.api.queries.cell_types_api.CellTypesApi`
`get_demix_file()` (al- method), 74
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base` (al-
 method), 113 `lensdk.core.cell_types_cache.CellTypesCache`
`get_demixed_traces()` (al- method), 238
`lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (al-
 method), 232 `lensdk.api.queries.cell_types_api.CellTypesApi`
`get_descendant_ids()` (al- method), 74
`lensdk.core.ontology.Ontology` (al- method), 249
`get_descendants()` (al- method), 249
`lensdk.core.ontology.Ontology` (al- method), 249
`get_dff_file()` (al- method), 238
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base` (al-
 method), 113 `lensdk.core.brain_observatory_nwb_data_set`,
`get_dff_traces()` (al- 236
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base` (al-
 method), 116 `lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base`
`get_dff_traces()` (al- method), 112
`lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_ophys_data_transforms.BehaviorOphysDataTran.`
 method), 123 `lensdk.brain_observatory.behavior.session_apis.data_io.behavior`
`get_dff_traces()` (al- method), 118
`lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (al-
 method), 233 `lensdk.brain_observatory.behavior.session_apis.data_io.behavior`
`get_diagonals_from_sparse()` (in module al- method), 120
`lensdk.brain_observatory.r_neuropil`), 203 `get_even_sampling()` (in module al-
`get_dimensions()` (al- `lensdk.brain_observatory.behavior.trials_processing`),
`lensdk.internal.morphology.morphology.Morphology` 152
 method), 326 `get_event_detection_filepath()` (al-
`get_disc_masks()` (in module al- `lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor`
`lensdk.brain_observatory.receptive_field_analysis.chisquare_method`), 113
 180 `get_events()` (allensdk.brain_observatory.behavior.session_apis.abcs.s
`get_dprime()` (in module al- method), 116
`lensdk.brain_observatory.behavior.dprime`), `get_events()` (allensdk.brain_observatory.behavior.session_apis.data_
 136 method), 123
`get_driver_line()` (al- `get_events_by_bit()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base`
 method), 112 method), 219

`get_events_by_line()` (al- `lensdk.internal.api.queries.mouse_connectivity_api_prerelease.M`
`lensdk.brain_observatory.sync_dataset.Dataset` method), 280
`method`), 219
`get_events_per_pixel()` (in module al- `lensdk.internal.core.mouse_connectivity_cache_prerelease.Mous`
`lensdk.brain_observatory.receptive_field_analysis.chisquare` method), 300
180
`get_excluded()` (al- `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityAp`
`lensdk.internal.brain_observatory.roi_filter_utils.TrainingModelClassifier`
`method`), 294
`get_expected_events_by_pixel()` (in module al- `lensdk.brain_observatory.behavior.metadata.behavior_metadata`)
110
`lensdk.brain_observatory.receptive_field_analysis.chisquare` method), 121
180
`get_experiment_analysis_file()` (in module al- `lensdk.brain_observatory.behavior.session_apis.data_transforms`
`method`), 121
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` method), 152
344
`get_experiment_container_metrics()` (al- `get_extractor()` (al-
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.brain_observatory.behavior.metadata.behavior_metadata`
`method`), 72
`method`), 109
`get_experiment_containers()` (al- `get_eye_tracking()` (al-
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.brain_observatory.behavior.session_apis.abcs.session_bas`
`method`), 72
`method`), 116
`get_experiment_detail()` (al- `get_eye_tracking()` (al-
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` `lensdk.brain_observatory.behavior.session_apis.data_transforms`
`method`), 85
`method`), 124
`get_experiment_files()` (in module al- `get_eye_tracking()` (in module al-
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` `lensdk.brain_observatory.behavior.sync`),
344
131
`get_experiment_id()` (al- `get_eye_tracking_filepath()` (al-
`lensdk.internal.api.lims_api.LimsApi` method), `lensdk.brain_observatory.behavior.session_apis.abcs.data_extrac`
282
`method`), 113
`get_experiment_nwb_file()` (in module al- `get_eye_tracking_rig_geometry()` (al-
`lensdk.internal.pipeline_modules.run_observatory_analysis` `lensdk.brain_observatory.behavior.session_apis.abcs.data_extrac`
342
`method`), 113
`get_experiment_nwb_file()` (in module al- `get_eye_tracking_rig_geometry()` (al-
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` `lensdk.brain_observatory.behavior.session_apis.abcs.session_bas`
344
`method`), 116
`get_experiment_session()` (in module al- `get_eye_tracking_rig_geometry()` (al-
`lensdk.internal.pipeline_modules.run_observatory_analysis` `lensdk.brain_observatory.behavior.session_apis.data_transforms`
342
`method`), 124
`get_experiment_structure_unionizes()` `get_eye_tracking_video_filepath_df()`
(al- `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` `lensdk.internal.api.lims_api.LimsApi`
`method`), 244
`method`), 282
`get_experiment_sweep_numbers()` (al- `get_falling_edges()` (al-
`lensdk.core.nwb_data_set.NwbDataSet` `lensdk.brain_observatory.sync_dataset.Dataset`
`method`), 247
`method`), 219
`get_experiments()` (al- `get_false_alarm_rate()` (in module al-
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` `lensdk.brain_observatory.behavior.dprime`),
`method`), 85
137
`get_experiments()` (al- `get_features()` (in module al-
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` `lensdk.internal.ephys.plot_qc_figures`), 302
`method`), 244
`method`), 244
`get_experiments()` (al- `get_features()` (in module al-
`lensdk.internal.ephys.plot_qc_figures3`),
3)

304
 get_field_of_view_shape() (al- lensdk.core.structure_tree.StructureTree
 lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorOp
 method), 113
 get_fluorescence() (al- lensdk.brain_observatory.behavior.trials_processing),
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 method), 212
 get_fluorescence_timestamps() (al- lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile
 lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet
 method), 233
 get_fluorescence_traces() (al- lensdk.brain_observatory.behavior.stimulus_processing.util),
 lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet
 method), 233
 get_foraging_id() (al- lensdk.api.queries.synchronization_api.SynchronizationApi
 lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorDataExtra
 method), 112
 get_foraging_id() (al- lensdk.api.queries.synchronization_api.SynchronizationApi
 lensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.BehaviorJsonExtractor
 method), 118
 get_foraging_id() (al- lensdk.api.queries.synchronization_api.SynchronizationApi
 lensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api.BehaviorLimsExtractor
 method), 120
 get_format() (allensdk.config.manifest.Manifest lensdk.api.queries.synchronization_api.SynchronizationApi
 method), 228 method), 95
 get_full_genotype() (al- get_images_dict() (in module al-
 lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.BehaviorDataExtra
 method), 112 128
 get_full_genotype() (al- get_imaging_depth() (al-
 lensdk.brain_observatory.behavior.session_apis.data_io.behavior_json_api.BehaviorJsonExtractor
 method), 118 method), 113
 get_full_genotype() (al- get_imaging_plane_group() (al-
 lensdk.brain_observatory.behavior.session_apis.data_io.behavior_lims_api.BehaviorLimsExtractor
 method), 120 method), 113
 get_gaussian_fit() (in module al- get_indicator_bound_point() (in module al-
 lensdk.brain_observatory.receptive_field_analysis.postprocessing), lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher),
 183 339
 get_gaussian_fit_single_channel() get_indices_by_distance() (in module al-
 (in module al- lensdk.internal.brain_observatory.roi_filter_utils),
 lensdk.brain_observatory.receptive_field_analysis.fit_parameters),
 182
 get_genes() (allensdk.api.queries.mouse_atlas_api.MouseAtlasApi lensdk.internal.mouse_connectivity.interval_unionize.data_utilitie
 method), 82 332
 get_go_responses() (in module al- get_injection_density() (al-
 lensdk.brain_observatory.behavior.dprime), lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
 137 method), 245
 get_gratings_metadata() (in module al- get_injection_fraction() (al-
 lensdk.brain_observatory.behavior.stimulus_processing), lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
 128 method), 245
 get_h() (in module al- get_input_data() (in module al-
 lensdk.internal.model.biophysical.passive_fitting.neuron_utils), lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 308 345
 get_hit_rate() (in module al- get_input_json() (in module al-
 lensdk.brain_observatory.behavior.dprime), lensdk.internal.core.lims_utilities), 298
 137
 get_intensity() (al-

`lensdk.mouse_connectivity.grid.subimage.base_subimage.ImageSubImage`
`method`), 363 `get_max_projection()` (al-
`get_invalid_times()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet`
`lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi`
`method`), 162 `get_max_projection_file()` (al-
`get_isi_experiments()` (al- `lensdk.brain_observatory.behavior.session_apis.abcs.data_transforms.ABCSessionBase`
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi`), 113
`method`), 72 `get_mean_response()` (al-
`get_isi_experiments()` (al- `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
`lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi`
`method`), 161 `get_mean_waveforms()` (al-
`get_isis()` (in module al- `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi`
`lensdk.ephys.ephys_features`), 276 `method`), 162
`get_keys()` (in module al- `get_metadata()` (al-
`lensdk.internal.brain_observatory.time_sync`), `lensdk.brain_observatory.behavior.session_apis.abcs.session_base.ABCSessionBase`
297 `method`), 114
`get_lfp()` (allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi (al-
`method`), 162 `lensdk.brain_observatory.behavior.session_apis.abcs.session_base.ABCSessionBase`
`get_lfp_channel_order()` (al- `method`), 117
`lensdk.brain_observatory.ecephys.file_io.continuous_file_io.ContinuousFileIO` (al-
`method`), 163 `lensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms`
`get_click_times()` (in module al- `method`), 121
`lensdk.brain_observatory.behavior.sync`), `get_metadata()` (al-
131 `lensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms`
`get_clicks()` (allensdk.brain_observatory.behavior.session_apis.abcs.session_base.ABCSessionBase.behavior_base.BehaviorBase
`method`), 114 `get_metadata()` (al-
`get_clicks()` (allensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms.behavior_data_transforms.BehaviorDataTransforms
`method`), 121 `method`), 162
`get_line()` (allensdk.brain_observatory.sync_dataset.Dataset) `get_metadata()` (al-
`method`), 220 `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet`
`get_line_changes()` (al- `method`), 233
`lensdk.brain_observatory.sync_dataset.Dataset` `get_metrics()` (in module al-
`method`), 220 `lensdk.internal.brain_observatory.annotated_region_metrics`),
284
`get_list_of_path_dict()` (in module al- `get_missing_path()` (al-
`lensdk.test_utilities.regression_fixture`), 370 `lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile`
`get_locally_sparse_noise_stimulus_template()` (allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet
`method`), 233 `get_monitor_delay()` (al-
`get_manifest()` (al- `lensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms`
`lensdk.config.manifest_builder.ManifestBuilder` `method`), 121
`method`), 229 `get_monitor_delay()` (al-
`get_manual_injection_summary()` (al- `lensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms`
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`), 124
`method`), 85 `get_morphology_features()` (al-
`get_mask()` (allensdk.brain_observatory.stimulus_info.Monitor `lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 214 `method`), 74
`get_mask_plane()` (al- `get_morphology_features()` (al-
`lensdk.brain_observatory.roi_masks.Mask` `lensdk.core.cell_types_cache.CellTypesCache`
`method`), 203 `method`), 238
`get_max_projection()` (al- `get_motion_correction()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.session_base.ABCSessionBase.behavior_observation_base.BehaviorObservationBase`
`method`), 117 `method`), 117
`get_max_projection()` (al- `get_motion_correction()` (al-
`lensdk.brain_observatory.behavior.session_apis.data_transforms.DataTransforms.behavior_observation_data_transforms.BehaviorObservationDataTransforms`
`method`), 117 `method`), 117

method), 124		method), 309
get_motion_correction()	(al- lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet. method), 234	(al- get_neuronal_models_by_id() lensdk.api.queries.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet. method), 76
get_mouse_id()	(al- lensdk.brain_observatory.behavior.session_apis.abcs.data_extraction_backend.behavior.session_apis.abcs.data_extraction_backend. method), 112	(al- get_neuropil_r() lensdk.brain_observatory.behavior.session_apis.abcs.data_extraction_backend.behavior.session_apis.abcs.data_extraction_backend. method), 234
get_mouse_id()	(al- lensdk.brain_observatory.behavior.session_apis.data_io.behavior.session_apis.data_io. method), 118	(al- get_neuropil_traces() lensdk.brain_observatory.behavior.session_apis.data_io.behavior.session_apis.data_io. method), 234
get_mouse_id()	(al- lensdk.brain_observatory.behavior.session_apis.data_io.behavior.session_apis.data_io. method), 120	(al- get_noise_correlation() lensdk.brain_observatory.behavior.session_apis.data_io.behavior.session_apis.data_io. method), 195
get_mouse_id() (in module lensdk.brain_observatory.behavior.trials_processing), 152	(al- get_noise_correlation() lensdk.brain_observatory.natural_scenes.NaturalScenes. method), 199	(al- get_noise_correlation() lensdk.brain_observatory.natural_scenes.NaturalScenes. method), 199
get_name_map()	(al- lensdk.core.structure_tree.StructureTree. method), 259	(al- get_noise_correlation() lensdk.brain_observatory.static_gratings.StaticGratings. method), 210
get_natural_movie_template()	(al- lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api. method), 161	(al- get_ophys_cell_segmentation_run_id() lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api. method), 113
get_natural_scene_template()	(al- lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api. method), 161	(al- get_ophys_container_id() lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api. method), 113
get_nearest()	(al- lensdk.brain_observatory.sync_dataset.Dataset. method), 220	(al- get_ophys_container_id() lensdk.brain_observatory.behavior.session_apis.data_io.behavior. method), 120
get_neuron_config()	(al- lensdk.api.queries.glif_api.GlifApi. method), 76	(al- get_ophys_data_length() (in module al- lensdk.internal.brain_observatory.time_sync), 297
get_neuron_configs()	(al- lensdk.api.queries.glif_api.GlifApi. method), 76	(al- get_ophys_experiment_id() lensdk.brain_observatory.behavior.session_apis.abcs.data_extraction_backend. method), 113
get_neuronal_model()	(al- lensdk.api.queries.glif_api.GlifApi. method), 76	(al- get_ophys_experiment_id() lensdk.brain_observatory.behavior.session_apis.abcs.session_base. method), 117
get_neuronal_model_runs()	(al- lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi. method), 278	(al- get_ophys_experiment_id() lensdk.brain_observatory.behavior.session_apis.data_transforms. method), 124
get_neuronal_model_templates()	(al- lensdk.api.queries.glif_api.GlifApi. method), 76	(al- get_ophys_experiment_ids() lensdk.brain_observatory.behavior.session_apis.data_io.behavior. method), 120
get_neuronal_models()	(al- lensdk.api.queries.biophysical_api.BiophysicalApi. method), 69	(al- get_ophys_experiments() lensdk.api.queries.brain_observatory_api.BrainObservatoryApi. method), 72
get_neuronal_models()	(al- lensdk.api.queries.glif_api.GlifApi. method), 76	(al- get_ophys_frames() (in module al- lensdk.brain_observatory.behavior.sync), 132
get_neuronal_models()	(al- lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi. method), 278	(al- get_ophys_session_id() lensdk.brain_observatory.behavior.session_apis.abcs.data_extraction_backend. method), 113
get_neuronal_models()	(al- lensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver. method), 278	(al- get_ophys_session_id() lensdk.brain_observatory.behavior.session_apis.abcs.session_base. method), 113

method), 117

get_ophys_session_id() (al- get_photodiode_events() (in module al-
lensdk.brain_observatory.behavior.session_apis.data_io.behavior.session_apis.data_io.
method), 120

get_ophys_session_id() (al- get_pipeline_version() (al-
lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior.session_apis.data_transforms.
method), 125

get_ophys_timestamps() (al- get_plane_group_count() (al-
lensdk.brain_observatory.behavior.session_apis.abcs.session_apis.abcs.data_extraction.
method), 117

get_ophys_timestamps() (al- get_polygons() (al-
lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior.session_apis.data_transforms.
method), 125

get_optogenetic_stimulation() (al- get_probe_lfp_data() (al-
lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.ecephys_session_api.ecephys_session_api.
method), 162

get_ori_info_from_trial() (in module al- get_probe_time_offset() (in module al-
lensdk.brain_observatory.behavior.trials_processing), lensdk.brain_observatory.ecephys.align_timestamps.barcode),
152

get_output() (allensdk.internal.mouse_connectivity.projection_image_info.
method), 337

get_overall_blend() (in module al- get_probes() (allensdk.brain_observatory.ecephys.ecephys_session_apis.
lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitching.
method), 162
339

get_params() (allensdk.brain_observatory.session_api_utils.ParallelSessionApi.
method), 209

get_params() (in module al- get_projection_data() (in module al-
lensdk.brain_observatory.behavior.trials_processing), lensdk.internal.mouse_connectivity.interval_unionize.data_utilitie
152

get_passive_fit_data() (in module al- get_projection_density() (al-
lensdk.internal.model.biophysical.passive_fitting.preprocess.
method), 245
308

get_path() (allensdk.config.manifest.Manifest get_projection_image_info() (al-
method), 228

get_path() (in module al-
lensdk.internal.pipeline_modules.run_demixing), get_projection_matrix() (al-
method), 85
341

get_peak() (allensdk.brain_observatory.drifting_gratings.DriftingGratings.
method), 195

get_peak() (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise.
method), 197

get_peak() (allensdk.brain_observatory.natural_movie.NaturalMovie.
method), 198

get_peak() (allensdk.brain_observatory.natural_scenes.NaturalScenes.
method), 199

get_peak() (allensdk.brain_observatory.static_gratings.StaticGratings.
method), 210

get_peak() (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis.
method), 212

get_peak_significance() (in module al- get_raw_cell_specimen_table_dict() (al-
lensdk.brain_observatory.receptive_field_analysis.chisquare.
method), 181

get_peaks() (in module al- get_raw_dff_data() (al-
lensdk.internal.model.glif.threshold_adaptation), lensdk.brain_observatory.behavior.session_apis.data_transforms.
method), 113

method), 125

get_raw_stimulus_frames() (in module `allensdk.brain_observatory.behavior.sync`), 132

get_real_photodiode_events() (in module `allensdk.internal.brain_observatory.time_sync`), 297

get_receptive_field() (in module `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`), 197

get_receptive_field_analysis_data() (in module `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`), 197

get_receptive_field_attribute_df() (in module `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`), 197

get_reconstruction() (in module `allensdk.core.cell_types_cache.CellTypesCache`), 238

get_reconstruction_markers() (in module `allensdk.core.cell_types_cache.CellTypesCache`), 238

get_recorded_data() (in module `allensdk.model.biophysical.utils.Utils`), 351

get_reference_aligned_image_channel_volumes() (in module `allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`), 85

get_reference_space() (in module `allensdk.core.reference_space_cache.ReferenceSpaceCache`), 253

get_reference_to_image() (in module `allensdk.api.queries.synchronization_api.SynchronizationApi`), 96

get_reporter_line() (in module `allensdk.brain_observatory.behavior.session_apis.abcs.data_extraction.Datasets`), 112

get_reporter_line() (in module `allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_apis.data_set.BrainObservatoryNWBDataIO`), 118

get_reporter_line() (in module `allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_apis.data_set.BrainObservatoryNWBDataIO`), 120

get_representational_similarity() (in module `allensdk.brain_observatory.drifting_gratings.DriftingGratings`), 196

get_representational_similarity() (in module `allensdk.brain_observatory.natural_scenes.NaturalScenes`), 199

get_representational_similarity() (in module `allensdk.brain_observatory.static_gratings.StaticGratings`), 211

get_response() (in module `allensdk.brain_observatory.drifting_gratings.DriftingGratings`), 196

get_response() (in module `allensdk.brain_observatory.natural_scenes.NaturalScenes`), 199

get_response() (in module `allensdk.brain_observatory.static_gratings.StaticGratings`), 211

get_response() (in module `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`), 212

get_response_latency() (in module `allensdk.brain_observatory.behavior.trials_processing`), 152

get_response_type() (in module `allensdk.brain_observatory.behavior.trials_processing`), 152

get_rewards() (in module `allensdk.brain_observatory.behavior.session_apis.abcs.session_apis`), 114

get_rewards() (in module `allensdk.brain_observatory.behavior.session_apis.data_transforms`), 121

get_rewards() (in module `allensdk.brain_observatory.behavior.session_apis.data_transforms`), 125

get_rewards() (in module `allensdk.brain_observatory.behavior.rewards_processing`), 143

get_rig_metadata() (in module `allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api`), 162

get_rigid_motion_transform_file() (in module `allensdk.brain_observatory.behavior.session_apis.abcs.data_extraction.Datasets`), 114

get_rising_edges() (in module `allensdk.brain_observatory.behavior.session_apis.abcs.data_extraction.Datasets`), 220

get_roi_ids() (in module `allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_apis.data_set.BrainObservatoryNWBDataIO`), 234

get_roi_mask() (in module `allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_apis.data_set.BrainObservatoryNWBDataIO`), 234

get_roi_mask_array() (in module `allensdk.brain_observatory.drifting_gratings.DriftingGratings`), 235

get_roi_masks_by_cell_roi_id() (in module `allensdk.brain_observatory.behavior.session_apis.data_transforms`), 125

get_rois() (in module `allensdk.internal.brain_observatory.roi_filter_utils`), 295

get_rolling_dprime() (in module `allensdk.brain_observatory.behavior.dprime`), 196

137 method), 118

get_running_acquisition_df() (al- get_sex() (allensdk.brain_observatory.behavior.session_apis.data_io.be
lensdk.brain_observatory.behavior.session_apis.abcs.session_apis.behavior_base.BehaviorBase
method), 114 get_shuffle_matrix() (in module al-

get_running_acquisition_df() (al- lensdk.brain_observatory.receptive_field_analysis.utilities),
lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_transforms.BehaviorDataTransforms
method), 121 get_signal_correlation() (al-

get_running_df() (in module al- lensdk.brain_observatory.drifting_gratings.DriftingGratings
lensdk.brain_observatory.behavior.running_processing), method), 196

143 get_signal_correlation() (al-

get_running_speed() (al- lensdk.brain_observatory.natural_scenes.NaturalScenes
lensdk.brain_observatory.behavior.session_apis.abcs.session_apis.behavior_base.BehaviorBase
method), 115 get_signal_correlation() (al-

get_running_speed() (al- lensdk.brain_observatory.static_gratings.StaticGratings
lensdk.brain_observatory.behavior.session_apis.data_transforms.behavior_data_transforms.BehaviorDataTransforms
method), 122 get_sitk_image_information() (in module al-

get_running_speed() (al- lensdk.core.sitk_utilities), 257

lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi (al-
method), 162 lensdk.core.reference_space.ReferenceSpace

get_running_speed() (al- method), 251

lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet epoch_mask_list() (in module al-

method), 235

get_schema() (allensdk.api.queries.rma_api.RmaApi lensdk.brain_observatory.receptive_field_analysis.utilities),
method), 91 184

get_screen_gaze_data() (al- get_spatial_grating() (in module al-
lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi),
method), 162 214

get_section_data_sets() (al- get_spatio_temporal_grating() (in module
lensdk.api.queries.mouse_atlas_api.MouseAtlasApi allensdk.brain_observatory.stimulus_info), 215
method), 82 get_speed_tuning() (al-

get_section_data_sets_by_product() (al- lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
lensdk.api.queries.image_download_api.ImageDownloadApi), method), 212

method), 80 get_spike_amplitudes() (al-

get_section_image_ranges() (al- lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_
lensdk.api.queries.image_download_api.ImageDownloadApi), method), 162

method), 81 get_spike_times() (al-

get_segmentation() (al- lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_
lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage get_subimage() (in module al-
method), 364 get_spike_times() (al-

get_session_data() (al- lensdk.core.nwb_data_set.NwbDataSet
lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi
method), 161 get_spikes() (in module al-

get_session_start_time() (al- lensdk.internal.ephys.plot_qc_figures), 302

lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi (al-
method), 162 lensdk.internal.ephys.plot_qc_figures3),

get_session_type() (al- 304

lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet characteristics() (in
method), 235 module allensdk.ephys.extract_cell_features),

get_sessions() (al- 277

lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi
method), 161 (in module al-

get_sex() (allensdk.brain_observatory.behavior.session_apis.abcs.session_apis.behavior_base.BehaviorDataTransforms
method), 112 309

get_sex() (allensdk.brain_observatory.behavior.session_apis.data_io.behavior_session_api.BehaviorSessionApi (in module al-

<code>lensdk.ephys.extract_cell_features()</code> , 277	<code>method)</code> , 235
<code>get_stim_data_length()</code> (in module <code>al-</code> <code>lensdk.internal.brain_observatory.time_sync</code>), 297	<code>get_stimulus_template()</code> (al- <code>lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet</code> <code>method)</code> , 235
<code>get_stim_photodiode()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.sync</code>), 132	<code>get_stimulus_templates()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 115
<code>get_stimulus()</code> (al- <code>lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet</code> <code>method)</code> , 235	<code>get_stimulus_templates()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 122
<code>get_stimulus_attr_changes()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.trials_processing</code>), 152	<code>get_stimulus_templates()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.stimulus_processing</code>), 129
<code>get_stimulus_epoch_table()</code> (al- <code>lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet</code> <code>method)</code> , 235	<code>get_stimulus_timestamps()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 115
<code>get_stimulus_file()</code> (al- <code>lensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver</code> <code>method)</code> , 309	<code>get_stimulus_timestamps()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 117
<code>get_stimulus_frame_rate()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 122	<code>get_stimulus_timestamps()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 122
<code>get_stimulus_mappings()</code> (al- <code>lensdk.api.queries.brain_observatory_api.BrainObservatoryApi</code> <code>method)</code> , 72	<code>get_stimulus_timestamps()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 125
<code>get_stimulus_metadata()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.stimulus_processing</code>), 128	<code>get_structure_graphs()</code> (al- <code>lensdk.api.queries.ontologies_api.OntologiesApi</code> <code>method)</code> , 86
<code>get_stimulus_name()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.data_io.behavior_session_base.SessionBase</code> <code>method)</code> , 112	<code>get_structure_mask()</code> (al- <code>lensdk.api.queries.behavior_session_base.SessionBase</code> <code>method)</code> , 253
<code>get_stimulus_name()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_io.behavior_session_base.SessionBase</code> <code>method)</code> , 119	<code>get_structure_mesh()</code> (al- <code>lensdk.api.queries.behavior_session_base.SessionBase</code> <code>method)</code> , 254
<code>get_stimulus_name()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_io.behavior_session_base.SessionBase</code> <code>method)</code> , 120	<code>get_structure_sets()</code> (al- <code>lensdk.api.queries.behavior_session_base.SessionBase</code> <code>method)</code> , 86
<code>get_stimulus_presentations()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 115	<code>get_structure_sets()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 259
<code>get_stimulus_presentations()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 117	<code>get_structure_to_image()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.abcs.session_base.SessionBase</code> <code>method)</code> , 96
<code>get_stimulus_presentations()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 122	<code>get_structure_tree()</code> (al- <code>lensdk.brain_observatory.behavior.session_api.data_transforms.DataTransforms</code> <code>method)</code> , 254
<code>get_stimulus_presentations()</code> (al- <code>lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.MouseConnectivityApi</code> <code>method)</code> , 162	<code>get_structure_unionizes()</code> (al- <code>lensdk.api.queries.ecephys_session_api.ecephys_session_api.MouseConnectivityApi</code> <code>method)</code> , 85
<code>get_stimulus_presentations()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.stimulus_processing</code>), 129	<code>get_structure_unionizes()</code> (al- <code>lensdk.core.mouse_connectivity_cache.MouseConnectivityCache</code> <code>method)</code> , 245
<code>get_stimulus_table()</code> (al- <code>lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet</code> <code>method)</code> , 235	<code>get_structure_unionizes()</code> (al- <code>lensdk.api.queries.mouse_connectivity_api.prerelease.MouseConnectivityApi</code> <code>method)</code> , 245

`method`), 281
`get_structures()` (al-
`lensdk.api.queries.ontologies_api.OntologiesApi`
`method`), 86
`get_structures_by_acronym()` (al-
`lensdk.core.structure_tree.StructureTree`
`method`), 259
`get_structures_by_id()` (al-
`lensdk.core.structure_tree.StructureTree`
`method`), 260
`get_structures_by_name()` (al-
`lensdk.core.structure_tree.StructureTree`
`method`), 260
`get_structures_by_set_id()` (al-
`lensdk.core.structure_tree.StructureTree`
`method`), 260
`get_structures_with_sets()` (al-
`lensdk.api.queries.ontologies_api.OntologiesApi`
`method`), 86
`get_sum_pixel_intensities()` (in module al-
`lensdk.internal.mouse_connectivity.interval_union_utilities`), 332
`get_sum_pixels()` (in module al-
`lensdk.internal.mouse_connectivity.interval_union_utilities`), 332
`get_surface_2p_pixel_size_um()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base.BehaviorOphysDataExtractorBase`
`method`), 114
`get_svg()` (allensdk.api.queries.svg_api.SvgApi
`method`), 94
`get_sweep()` (allensdk.core.nwb_data_set.NwbDataSet
`method`), 247
`get_sweep_metadata()` (al-
`lensdk.core.nwb_data_set.NwbDataSet`
`method`), 247
`get_sweep_numbers()` (al-
`lensdk.core.nwb_data_set.NwbDataSet`
`method`), 248
`get_sweep_numbers()` (in module al-
`lensdk.internal.model.glif.find_sweeps`), 313
`get_sweep_response()` (al-
`lensdk.brain_observatory.natural_movie.NaturalMovie`
`method`), 198
`get_sweep_response()` (al-
`lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`method`), 212
`get_sweep_v_i_t_from_set()` (in module al-
`lensdk.internal.model.biophysical.ephys_utils`), 309
`get_sweeps_by_name()` (in module al-
`lensdk.internal.model.glif.find_sweeps`), 313
`get_sweeps_of_type()` (in module al-
`lensdk.internal.model.biophysical.ephys_utils`), 309
`get_sync_data()` (al-
`lensdk.brain_observatory.behavior.session_apis.data_transforms`
`method`), 125
`get_sync_data()` (in module al-
`lensdk.brain_observatory.behavior.sync`), 133
`get_sync_file()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base.BehaviorOphysDataExtractorBase`
`method`), 114
`get_sync_ticks()` (al-
`lensdk.brain_observatory.behavior.session_apis.data_transforms`
`method`), 125
`get_synchronized_frame_times()` (in module
`allensdk.brain_observatory.sync_utilities`), 187
`get_targeted_structure()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.data_extractor_base.behavior_ophys_data_extractor_base.BehaviorOphysDataExtractorBase`
`method`), 114
`get_task_parameters()` (al-
`lensdk.brain_observatory.behavior.session_apis.abcs.session_base.behavior_ophys_data_extractor_base.BehaviorOphysDataExtractorBase`
`method`), 115
`get_template_utilities_parameters()` (al-
`lensdk.brain_observatory.behavior.session_apis.data_transforms`
`method`), 122
`get_template_utilities_parameters()` (in module al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata_utilities`), 110
`get_template_utilities_parameters()` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata_utilities`
`method`), 309
`get_template_volume()` (al-
`lensdk.core.reference_space_cache.ReferenceSpaceCache`
`method`), 254
`get_time()` (in module al-
`lensdk.brain_observatory.behavior.trials_processing`), 152
`get_time_string()` (in module al-
`lensdk.internal.ephys.plot_qc_figures`), 302
`get_time_string()` (in module al-
`lensdk.internal.ephys.plot_qc_figures3`), 304
`get_tree()` (allensdk.api.queries.tree_search_api.TreeSearchApi
`method`), 96
`get_trial_bounds()` (in module al-
`lensdk.brain_observatory.behavior.trials_processing`), 152
`get_trial_count_corrected_false_alarm_rate()`
(in module al-
`lensdk.brain_observatory.behavior.dprime`), 137
`get_trial_count_corrected_hit_rate()`
(in module al-
`lensdk.brain_observatory.behavior.dprime`), 137
`get_trial_image_names()` (in module al-

`lensdk.brain_observatory.behavior.trials_processing.GlifTypes (allensdk.api.queries.glif_api.GlifApi attribute), 76`
`153`
`get_trial_reward_time() (in module al- GlifApi (class in allensdk.api.queries.glif_api), 76`
`lensdk.brain_observatory.behavior.trials_processing), 315`
`153`
`GlifBadInitializationException, 315`
`GlifBadResetException, 352`
`get_trial_timing() (in module al- GlifExperiment (class in al-`
`lensdk.brain_observatory.behavior.trials_processing), lensdk.internal.model.glif.glif_experiment),`
`153`
`313`
`get_trials() (allensdk.brain_observatory.behavior.session_info.session_info_base.allensdk.model.glif.GlifExperimentBase`
`method), 115`
`352`
`get_trials() (allensdk.brain_observatory.behavior.session_info.session_info_base.allensdk.model.glif.GlifExperimentBase`
`method), 123`
`GlifNeuronMethod (class in al-`
`get_trials_from_data_transform() lensdk.model.glif.glif_neuron_methods),`
`(in module al- 356`
`lensdk.brain_observatory.behavior.trials_processing), 315`
`154`
`GlifOptimizer (class in al-`
`get_trials_v0() (in module al- 315`
`lensdk.brain_observatory.behavior.trials_processing), 315`
`155`
`GlifOptimizerNeuron (class in al-`
`get_trigger() (in module al- 315`
`lensdk.brain_observatory.behavior.sync), 315`
`133`
`grating_to_screen() (al-`
`get_unit_analysis_metrics() (al- method), 213`
`lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor`
`method), 161`
`lensdk.brain_observatory.ecephys.ecephys_project_api.EcephysProjectApi (al-`
`get_unit_filter_value() (in module al- method), 214`
`lensdk.brain_observatory.ecephys), 177`
`GRID_DATA_DIRECTORY (al-`
`get_units() (allensdk.brain_observatory.ecephys.ecephys_project_api.EcephysProjectApi (al-`
`method), 161`
`lensdk.brain_observatory.stimulus_info.Monitor`
`attribute), 280`
`get_units() (allensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi (al-`
`method), 162`
`lensdk.mouse_connectivity.grid.utilities.image_utilities),`
`367`
`get_video_length() (in module al- 367`
`lensdk.internal.brain_observatory.time_sync), 367`
`298`
`GridDataApi (class in al-`
`get_visual_stimuli_df() (in module al- GridDataApi (class in al-`
`lensdk.brain_observatory.behavior.stimulus_processing), lensdk.internal.api.queries.grid_data_api), 77`
`130`
`279`
`GridDataApiPrerelease (class in al-`
`get_volume_scale() (in module al- 279`
`lensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic),`
`334`
`h (allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils`
`get_well_known_file_by_name() (in module attribute), 348`
`allensdk.internal.core.lims_utilities), 298`
`h5_object_matcher_relname_in() (in module`
`get_well_known_file_by_type() (in module allensdk.core.h5_utilities), 240`
`allensdk.internal.core.lims_utilities), 299`
`handle_output_image() (in module al-`
`get_well_known_file_ids() (al- lensdk.internal.mouse_connectivity.projection_thumbnail.generator`
`lensdk.api.queries.biophysical_api.BiophysicalApi 336`
`method), 70`
`handle_pyramid() (in module al-`
`get_well_known_files_by_name() (in module lensdk.mouse_connectivity.grid.writers),`
`allensdk.internal.core.lims_utilities), 299`
`369`
`has_fixed_dt() (in module al-`
`get_well_known_files_by_type() (in module allensdk.ephys.ephys_features), 276`
`allensdk.internal.core.lims_utilities), 299`
`has_overlaps() (al-`
`get_wkf() (in module al- has_overlaps() (al-`
`lensdk.internal.pipeline_modules.run_ophys_eye_calibration), 260`
`345`
`lensdk.core.structure_tree.StructureTree`
`method), 260`

Hdf5Util (class in allensdk.config.model.formats.hdf5_util), 224

height (allensdk.brain_observatory.stimulus_info.Monitor attribute), 214

hex_pack() (in module allensdk.brain_observatory.circle_plots), 192

hex_to_rgb() (allensdk.core.structure_tree.StructureTree static method), 260

HocUtils (class in allensdk.model.biophys_sim.neuron.hoc_utils), 348

hook() (in module allensdk.brain_observatory), 221

host (allensdk.core.authentication.DbCredentials attribute), 231

HUMAN (allensdk.api.queries.cell_types_api.CellTypesApi attribute), 73

HUMAN_ORGANISM (allensdk.api.queries.mouse_atlas_api.MouseAtlasApi attribute), 81

imaging_depth (allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata attribute), 111

imaging_plane_group (allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata attribute), 111

imaging_plane_group_count (allensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata attribute), 111

INCLUDE (allensdk.api.queries.rma_api.RmaApi attribute), 89

INDETERMINATE

indicator (allensdk.brain_observatory.behavior.metadata.behavior_metadata attribute), 109

infer_dims() (allensdk.brain_observatory.circle_plots.CoronaPlotter method), 190

infer_dims() (allensdk.brain_observatory.circle_plots.FanPlotter method), 191

info() (allensdk.internal.model.biophysical.run_optimize.RunOptimize method), 310

init_by_mask() (allensdk.brain_observatory.roi_masks.NeuropilMask method), 204

init_by_mask() (allensdk.brain_observatory.roi_masks.RoiMask method), 204

init_by_pixels() (allensdk.brain_observatory.roi_masks.Mask method), 203

initial_processing_templates() (in module allensdk.internal.brain_observatory.itracker_utils), 291

initial_processing_templates() (in module allensdk.internal.brain_observatory.itracker_utils), 291

initialize_coarse_volume() (allensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder method), 369

initialize_hoc() (allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils method), 349

initialize_image() (allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile method), 340

initialize_image() (in module allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher), 339

initialize_images() (in module allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher), 339

initiate_unique_seed() (allensdk.internal.model.glif.glif_optimizer.GlifOptimizer

identify_valid_masks() (in module allensdk.brain_observatory.demixer), 193

identify_valid_masks() (in module allensdk.internal.brain_observatory.demixer), 286

Image (class in allensdk.brain_observatory.behavior.image_api), 141

image_from_array() (in module allensdk.mouse_connectivity.grid.utilities.image_utilities), 367

image_names (allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate attribute), 126

image_set_name (allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate attribute), 126

ImageApi (class in allensdk.brain_observatory.behavior.image_api), 141

ImageDownloadApi (class in allensdk.api.queries.image_download_api), 78

ImageOutputStream (class in allensdk.internal.brain_observatory.frame_stream), 291

images (allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate attribute), 126

ImageSeriesGridder (class in allensdk.mouse_connectivity.grid.image_series_gridder), 369

ImageSheet (class in allensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet), 336

method), 315

INJECTION_DENSITY (al-
lensdk.api.queries.grid_data_api.GridDataApi
attribute), 77

INJECTION_DENSITY_KEY (al-
lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
attribute), 242

INJECTION_ENERGY (al-
lensdk.api.queries.grid_data_api.GridDataApi
attribute), 77

INJECTION_FRACTION (al-
lensdk.api.queries.grid_data_api.GridDataApi
attribute), 77

INJECTION_FRACTION_KEY (al-
lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
attribute), 242

input_data() (allensdk.internal.core.lims_pipeline_module.PipelineModule
method), 298

input_resistance() (in module al-
lensdk.ephys.ephys_extractor), 269

InputFile (class in al-
lensdk.brain_observatory.argschema_utilities),
188

insert_iclamp() (al-
lensdk.internal.model.biophysical.deap_utils.Utils
method), 309

INTEGER (allensdk.api.queries.connected_services.ConnectedServices
attribute), 76

INTENSITY (allensdk.api.queries.grid_data_api.GridDataApi
attribute), 77

IntensitySubImage (class in al-
lensdk.mouse_connectivity.grid.subimage.base_subimage),
363

interlength (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise
attribute), 197

interlength (allensdk.brain_observatory.natural_scenes.NaturalScenes
attribute), 199

interlength (allensdk.brain_observatory.static_gratings.StaticGratings
attribute), 211

interpolate_RF() (in module al-
lensdk.brain_observatory.receptive_field_analysis.chisquare_tf),
181

interpolate_spike_time() (in module al-
lensdk.model.glif.glif_neuron), 355

interpolate_spike_value() (in module al-
lensdk.model.glif.glif_neuron), 355

interpolate_spike_voltage() (in module al-
lensdk.internal.model.glif.glif_optimizer_neuron),
318

intersection() (al-
lensdk.internal.brain_observatory.mask_set.MaskSet
method), 292

intersection_size() (al-
lensdk.internal.brain_observatory.mask_set.MaskSet
method), 292

IntervalUnionizer (class in al-
lensdk.internal.mouse_connectivity.interval_unionize.interval_uni-
332

invnl() (in module allensdk.internal.model.GLM), 323

is_change_event() (in module al-
lensdk.brain_observatory.behavior.stimulus_processing),
130

is_empty() (allensdk.config.model.description.Description
method), 226

is_equal() (in module al-
lensdk.brain_observatory.session_api_utils),
210

is_feature_affected_by_clipping() (allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor
method), 267

is_well_known_file_type() (al-
lensdk.api.queries.biophysical_api.BiophysicalApi
method), 70

isicv() (allensdk.ephys.feature_extractor.EphysFeatureExtractor
method), 277

items() (allensdk.brain_observatory.behavior.stimulus_processing.stimu-
method), 126

J

json_handler() (in module al-
lensdk.core.json_utilities), 241

json_msg_query() (allensdk.api.api.Api method),
104

json_remove_keys() (al-
lensdk.api.warehouse_cache.cache.Cache
static method), 98

json_remove_keys() (al-
lensdk.api.warehouse_cache.cache.Cache
static method), 98

JsonComments (class in allensdk.core.json_utilities),
440

JsonDescriptionParser (class in al-
lensdk.config.model.formats.json_description_parser),
27

JSONEncoder (class in allensdk.brain_observatory),
221

K

keyed_locate_h5_objects() (in module al-
lensdk.core.h5_utilities), 240

keys() (allensdk.brain_observatory.behavior.stimulus_processing.stimu-
method), 126

label_data() (allensdk.internal.brain_observatory.roi_filter_utils.Train-
method), 295

latency() (in module `allensdk.ephys.ephys_features`), 276
 LazyProperty (allensdk.core.lazy_property.lazy_property_mixin.LazyPropertyMixin attribute), 230
 LazyProperty (class in `allensdk.core.lazy_property.lazy_property`), 230
 LazyPropertyMixin (class in `allensdk.core.lazy_property.lazy_property_mixin`), 230
 least_squares_RCEl_calc_tested() (in module `allensdk.internal.model.glif.rc`), 319
 legacy() (in module `allensdk.deprecated`), 371
 legend_artist() (allensdk.brain_observatory.observatory_plots.DimensionPatchHandler method), 200
 lims_working_directory() (allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader method), 278
 lims_working_directory() (allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader method), 281
 LimsApi (class in `allensdk.internal.api.lims_api`), 282
 line_crossing_x() (in module `allensdk.model.glif.glif_neuron`), 355
 line_crossing_y() (in module `allensdk.model.glif.glif_neuron`), 355
 line_stats() (allensdk.brain_observatory.sync_dataset.Dataset method), 220
 linear_transform_from_intervals() (in module `allensdk.brain_observatory.ecephys.align_timestamps_barcode`), 157
 linux_to_windows() (in module `allensdk.internal.core.lims_utilities`), 299
 list_cells() (allensdk.api.queries.cell_types_api.CellTypesApi method), 74
 list_cells_api() (allensdk.api.queries.cell_types_api.CellTypesApi method), 75
 list_column_definition_class_names() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 73
 list_isi_experiments() (allensdk.api.queries.brain_observatory_api.BrainObservatoryApi method), 73
 list_neuronal_models() (allensdk.api.queries.glif_api.GlifApi method), 76
 list_of_dicts_to_dict_of_lists() (in module `allensdk.brain_observatory.receptive_field_analysis_tools`), 184
 list_stimuli() (allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataTestMethod), 235
 load() (allensdk.brain_observatory.sync_dataset.Dataset method), 220
 load() (allensdk.config.app.application_config.ApplicationConfig method), 223
 load() (allensdk.model.biophys_sim.config.Config method), 349
 load_annotation() (in module `allensdk.internal.mouse_connectivity.interval_unionize.data_utilities`), 332
 load_api_schema() (allensdk.api.api.Api method), 244
 load_arrays() (in module `allensdk.internal.pipeline_modules.run_annotated_region_metrics`), 311
 load_cell_parameters() (allensdk.internal.model.biophysical.deap_utils.Utils method), 300
 load_cell_parameters() (allensdk.model.biophysical.utils.AllActiveUtils method), 351
 load_cell_parameters() (allensdk.model.biophysical.utils.Utils method), 352
 load_config() (allensdk.config.manifest.Manifest method), 229
 load_csv() (allensdk.api.warehouse_cache.cache.Cache method), 98
 load_datasets_by_relnames() (in module `allensdk.core.h5_utilities`), 240
 load_description() (in module `allensdk.model.biophysical.runner`), 350
 load_experiment() (in module `allensdk.internal.ephys.plot_qc_figures`), 302
 load_experiment() (in module `allensdk.internal.ephys.plot_qc_figures3`), 304
 load_eye_tracking_hdf() (in module `allensdk.brain_observatory.behavior.eye_tracking_processing`), 139
 load_frame() (in module `allensdk.internal.brain_observatory.ophys_session_decomposition`), 293
 load_json() (allensdk.api.warehouse_cache.cache.Cache method), 98
 load_manifest() (allensdk.api.warehouse_cache.cache.Cache method), 98
 load_manifest() (allensdk.internal.model.biophysical.run_optimize.RunOptimize method), 310

`load_manifest()` (al- method), 214
`lensdk.model.biophysical.run_simulate.RunSimulate` (al- LSN_mask (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute), 197
method), 350
`load_morphology()` (in module al- LSN_OFF (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute), 197
`lensdk.internal.model.biophysical.passive_fitting.neuron_utilities`), 308
LSN_OFF_SCREEN (al-
`load_pickle()` (in module al- `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` attribute), 197
`lensdk.brain_observatory.behavior.stimulus_processing`), 130
LSN_ON (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute), 197
`load_sweep()` (in module al-
`lensdk.internal.model.data_access`), 323
`lsna_check_hvas()` (in module al-
`load_sweep()` (in module al- `lensdk.internal.pipeline_modules.run_observatory_thumbnails`), 345
`lensdk.model.glif.simulate_neuron`), 361
`load_sweeps()` (in module al-
`lensdk.internal.model.data_access`), 323
M
`local_path(allensdk.api.cloud_cache.file_attributes.CacheFileAttributes module allensdk.brain_observatory.dff)`, 195
attribute), 66
`main()` (in module al-
`local_time()` (in module al- `lensdk.brain_observatory.ecephys.stimulus_table.visualization.visualization`), 167
`lensdk.brain_observatory.behavior.trials_processing`), 155
`main()` (in module al-
`LocallySparseNoise` (class in al- `lensdk.brain_observatory.session_analysis`), 208
`lensdk.brain_observatory.locally_sparse_noise`), 197
`main()` (in module al-
`locate_h5_objects()` (in module al- `lensdk.internal.ephys.plot_qc_figures`), 302
`lensdk.core.h5_utilities`), 240
`main()` (in module al-
`locate_median()` (in module al- `lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit`), 181
`lensdk.brain_observatory.receptive_field_analysis.chisquare`), 181
`main()` (in module al-
`log(allensdk.config.manifest.Manifest attribute)`, 229
`lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit`), 308
`log(allensdk.config.model.description_parser.DescriptionParser attribute)`, 227
`main()` (in module al-
`log(allensdk.config.model.formats.json_description_parser.JsonDescriptionParser attribute)`, 224
`lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit`), 308
`log(allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser attribute)`, 225
`lensdk.internal.model.biophysical.passive_fitting.preprocess`), 308
`long_squares_features()` (al-
`lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` attribute), 265
`method`), 265
`long_squares_stim_amps()` (al-
`lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` attribute), 265
`method`), 265
`lensdk.internal.model.biophysical.run_optimize`), 310
`lensdk.internal.model.biophysical.run_passive_fit`),
`LSN(allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute)`, 197
`main()` (in module al-
`lsn_coordinate_to_monitor_coordinate()` (in module al- `lensdk.internal.model.biophysical.run_simulate_lims`), 311
`lensdk.brain_observatory.stimulus_info`), 215
`main()` (in module al-
`LSN_GREY(allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute)`, 197
`lensdk.internal.model.glif.find_sweeps`), 313
`in module al-
lsn_image_to_screen()` (al-
`lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor` attribute), 213
`method`), 213
`lensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_run`), 340
`lsn_image_to_screen()` (al-
`lensdk.brain_observatory.stimulus_info.Monitor` attribute), 213
`main()` (in module al-

```

lensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap(),
340 lensdk.internal.ephys.plot_qc_figures3),
305
main() (in module al- make_cell_page() (in module al-
lensdk.internal.pipeline_modules.gbm.generate_gbm_sample_heatmap(),
341 lensdk.internal.ephys.plot_qc_figures3), 302
341 make_cell_page() (in module al-
main() (in module al- lensdk.internal.ephys.plot_qc_figures3),
lensdk.internal.pipeline_modules.run_annotated_region_masks(),
341 make_display_mask() (in module al-
main() (in module al- lensdk.brain_observatory.stimulus_info),
lensdk.internal.pipeline_modules.run_demixing), 215
341 make_fit() (allensdk.internal.model.biophysical.run_optimize.RunOptimize
main() (in module al- method), 310
lensdk.internal.pipeline_modules.run_dff_computation), fit_json_file() (al-
341 lensdk.internal.model.biophysical.make_deep_fit_json.Report
main() (in module al- method), 310
lensdk.internal.pipeline_modules.run_neuropil_correction), incushion_plot() (in module al-
342 lensdk.brain_observatory.circle_plots), 192
main() (in module al- make_pixel_counter() (al-
lensdk.internal.pipeline_modules.run_observatory_analysis), lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage
342 method), 364
main() (in module al- make_ratio_volume() (al-
lensdk.internal.pipeline_modules.run_observatory_thumbnail), lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeries
345 method), 369
main() (in module al- make_spontaneous_activity_tables()
lensdk.internal.pipeline_modules.run_ophys_eye_calibration), (in module al-
345 lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes
main() (in module al- 171
lensdk.internal.pipeline_modules.run_ophys_session_decomposition), make_mask() (al-
345 lensdk.core.reference_space.ReferenceSpace
main() (in module al- method), 251
lensdk.internal.pipeline_modules.run_ophys_time_course), sweep_html() (in module al-
347 lensdk.internal.ephys.plot_qc_figures3), 302
main() (in module al- make_sweep_html() (in module al-
lensdk.internal.pipeline_modules.run_tissuecyte_unionize_cells), lensdk.internal.ephys.plot_qc_figures3),
348 305
main() (in module al- make_sweep_page() (in module al-
lensdk.model.biophysical.run_simulate), lensdk.internal.ephys.plot_qc_figures3), 302
350 make_sweep_page() (in module al-
main() (in module al- lensdk.internal.ephys.plot_qc_figures3),
lensdk.model.glif.simulate_neuron), 361 305
main() (in module allensdk.morphology.validate_swc), makeBasis_StimKernel() (in module al-
362 lensdk.internal.model.GLM), 323
make_bbs() (in module al- makeBasis_StimKernel_exp() (in module al-
lensdk.internal.brain_observatory.mask_set), lensdk.internal.model.GLM), 323
293 makeFitStruct_GLM() (in module al-
make_blended_tile() (in module al- lensdk.internal.model.GLM), 323
lensdk.internal.mouse_connectivity.tissuecyte_stitching_stitcher), Manifest (class in allensdk.api.cloud_cache.manifest),
339 66
make_category_dummy() (in module al- Manifest (class in allensdk.config.manifest), 227
lensdk.brain_observatory.chisquare_categorical), manifest_dataframe() (al-
190 lensdk.api.warehouse_cache.cache.Cache
make_cell_html() (in module al- method), 98
lensdk.internal.ephys.plot_qc_figures3), 302 MANIFEST_VERSION (al-
make_cell_html() (in module al- lensdk.core.cell_types_cache.CellTypesCache

```


`attribute`), 237
MANIFEST_VERSION (al-
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` (method), 292
`attribute`), 242
MANIFEST_VERSION (al-
`lensdk.core.reference_space_cache.ReferenceSpaceCache` (method), 252
`attribute`), 252
ManifestBuilder (class in al-
`lensdk.config.manifest_builder`), 229
ManifestVersionError, 229
many_structure_masks() (al-
`lensdk.core.reference_space.ReferenceSpace` (method), 251
map_column_names() (in module al-
`lensdk.brain_observatory.ecephys.stimulus_table.naming_utils`), 172
map_monitor_coordinate_to_stimulus_coordinate() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
map_monitor_coordinate_to_template_coordinate() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
map_stimulus() (al-
`lensdk.brain_observatory.stimulus_info.Monitor` (method), 214
map_stimulus() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
map_stimulus_coordinate_to_monitor_coordinate() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
map_stimulus_names() (in module al-
`lensdk.brain_observatory.ecephys.stimulus_table.naming_utils`), 172
map_template_coordinate_to_monitor_coordinate() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
Marker (class in `allensdk.core.swc`), 261
Marker (class in `allensdk.internal.core.swc`), 301
MARKER_FILE_TYPE (al-
`lensdk.api.queries.cell_types_api.CellTypesApi` (attribute), 73
MARKER_KEY (`allensdk.core.cell_types_cache.CellTypesCache` (attribute), 237
mask (`allensdk.brain_observatory.stimulus_info.Monitor` (attribute), 214
Mask (class in `allensdk.brain_observatory.roi_masks`), 203
mask() (`allensdk.internal.brain_observatory.mask_set.MaskSet` (method), 292
mask_is_union_of_set() (al-
`lensdk.internal.brain_observatory.mask_set.MaskSet` (method), 292
mask_nulls() (in module al-
`lensdk.internal.ephys.plot_qc_figures`), 303
mask_nulls() (in module al-
`lensdk.internal.ephys.plot_qc_figures3`), 305
mask_stimulus_template() (in module al-
`lensdk.brain_observatory.stimulus_info`), 215
MaskSet (class in al-
`lensdk.internal.brain_observatory.mask_set`), 292
match_barcodes() (in module al-
`lensdk.brain_observatory.ecephys.align_timestamps.barcode`), 158
match_barcodes() (in module al-
`lensdk.brain_observatory.ecephys.align_timestamps.barcode`), 158
max_of_line_and_const() (in module al-
`lensdk.model.glif.glif_neuron_methods`), 358
max_projection() (in module al-
`lensdk.internal.mouse_connectivity.projection_thumbnail.projection_thumbnail`), 337
max_voxel_density (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_utils` (attribute), 334
max_voxel_index (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_utils` (attribute), 334
mean_features_spike_zero() (in module al-
`lensdk.ephys.extract_cell_features`), 277
mean_response (al-
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` (attribute), 197
mean_response (al-
`lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` (attribute), 213
medfilt_custom() (in module al-
`lensdk.internal.brain_observatory.itracker_utils`), 291
median_absolute_deviation() (in module al-
`lensdk.internal.brain_observatory.itracker_utils`), 292
meets_engagement_criteria() (in module al-
`allensdk.brain_observatory.behavior.criteria`), 135
membrane_time_constant() (in module al-
`lensdk.ephys.ephys_extractor`), 269
memoize() (in module al-
`lensdk.api.warehouse_cache.cache`), 100
merge_mean_response() (al-
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` (static method), 197

[metadata_file_attributes\(\)](#) (alias [lensdk.api.cloud_cache.manifest.Manifest](#) method), 66
[metadata_file_names](#) (alias [lensdk.api.cloud_cache.manifest.Manifest](#) attribute), 67
[METHOD](#) ([allensdk.core.authentication.CredentialProvider](#) attribute), 231
[METHOD](#) ([allensdk.core.authentication.EnvCredentialProvider](#) attribute), 231
[midpoint\(\)](#) (in module [lensdk.internal.morphology.node](#)), 331
[min_of_line_and_zero\(\)](#) (in module [lensdk.model.glif.glif_neuron_methods](#)), 358
[minmax_norm\(\)](#) (in module [lensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities](#)), 337
[MissingDataError](#), 239
[MissingStimulusException](#), 189
[MissingSweepException](#), 312
[MLIN\(\)](#) (in module [allensdk.internal.model.glif.MLIN](#)), 311
[MLIN_list_error\(\)](#) (in module [lensdk.internal.model.glif.error_functions](#)), 312
[mod_file_entries\(\)](#) (alias [lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 278
[mod_file_entries\(\)](#) (alias [lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) method), 281
[mod_file_paths\(\)](#) (alias [lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 278
[mod_file_paths\(\)](#) (alias [lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) method), 281
[MOD_FILE_TYPE_ID](#) (alias [lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) attribute), 278
[MOD_FILE_TYPE_ID](#) (alias [lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) attribute), 281
[MODEL](#) ([allensdk.api.queries.rma_api.RmaApi](#) attribute), 89
[model_query\(\)](#) (alias [lensdk.api.queries.rma_api.RmaApi](#) method), 91
[model_stage\(\)](#) (alias [lensdk.api.queries.rma_api.RmaApi](#) method), 91
[model_type\(\)](#) ([allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 278
[modify_parameter\(\)](#) (alias [lensdk.model.glif.glif_neuron_methods.GlifNeuronMethod](#) method), 356
[moments2\(\)](#) (in module [lensdk.brain_observatory.receptive_field_analysis.fitgaussian2D](#)), 183
[Monitor](#) (class in [lensdk.brain_observatory.stimulus_info](#)), 214
[monitor_coordinate_to_lsn_coordinate\(\)](#) (in module [lensdk.brain_observatory.stimulus_info](#)), 216
[monitor_coordinate_to_natural_movie_coordinate\(\)](#) (in module [lensdk.brain_observatory.stimulus_info](#)), 216
[monitor_delay](#) (alias [lensdk.internal.brain_observatory.time_sync.OphysTimeAligner](#) attribute), 296
[Morphology](#) (class in [allensdk.core.swc](#)), 261
[Morphology](#) (class in [lensdk.internal.morphology.morphology](#)), 324
[MORPHOLOGY_FEATURES_KEY](#) (alias [lensdk.core.cell_types_cache.CellTypesCache](#) attribute), 237
[morphology_entries\(\)](#) (alias [lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 279
[morphology_entries\(\)](#) (alias [lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) method), 281
[morphology_paths\(\)](#) (alias [lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) method), 279
[morphology_paths\(\)](#) (alias [lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) method), 281
[MORPHOLOGY_COLORS](#) (class in [lensdk.internal.morphology.morphvis](#)), 328
[mostly_useful\(\)](#) (in module [lensdk.brain_observatory.behavior.criteria](#)), 135
[MOTION_CORRECTION_DATASETS](#) (alias [lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNWBDataSet](#) attribute), 237
[MOUSE](#) ([allensdk.api.queries.cell_types_api.CellTypesApi](#) attribute), 237

attribute), 73
 MOUSE_2011 (allensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 87
 MOUSE_ATLAS_PRODUCTS (al-
 lensdk.api.queries.mouse_atlas_api.MouseAtlasApi
 attribute), 81
 mouse_id (allensdk.brain_observatory.behavior.metadata.behavior_metadata.natural_movie),
 attribute), 109
 MOUSE_ORGANISM (al-
 lensdk.api.queries.mouse_atlas_api.MouseAtlasApi
 attribute), 82
 MouseAtlasApi (class in al-
 lensdk.api.queries.mouse_atlas_api), 81
 MouseConnectivityApi (class in al-
 lensdk.api.queries.mouse_connectivity_api),
 82
 MouseConnectivityApiPrerelease (class in al-
 lensdk.internal.api.queries.mouse_connectivity_api_prerelease),
 280
 MouseConnectivityCache (class in al-
 lensdk.core.mouse_connectivity_cache),
 241
 MouseConnectivityCachePrerelease
 (class in al-
 lensdk.internal.core.mouse_connectivity_cache_prerelease),
 299
 movingaverage() (in module al-
 lensdk.brain_observatory.dff), 195
 movingmode_fast() (in module al-
 lensdk.brain_observatory.dff), 195
 multi_dataframe_merge() (in module al-
 lensdk.brain_observatory.session_analysis),
 208
N
 n_complete() (in module al-
 lensdk.brain_observatory.behavior.criteria),
 135
 NA (allensdk.core.cell_types_cache.ReporterStatus at-
 tribute), 239
 name (allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusImage
 attribute), 126
 nan_get() (in module al-
 lensdk.internal.ephys.core.feature_extract),
 302
 natural_movie_coordinate_to_monitor_coordinate (in module al-
 lensdk.brain_observatory.stimulus_info),
 216
 natural_movie_image_to_screen() (al-
 lensdk.brain_observatory.stimulus_info.Monitor
 method), 214
 natural_scene_coordinate_to_monitor_coordinate (in module al-
 lensdk.brain_observatory.stimulus_info),
 216
 natural_scene_image_to_screen() (al-
 lensdk.brain_observatory.stimulus_info.Monitor
 method), 214
 NaturalMovie (class in al-
 lensdk.brain_observatory.natural_movies),
 198
 NaturalScenes (class in al-
 lensdk.brain_observatory.natural_scenes),
 199
 nearest_neuron_sampling_rate() (al-
 lensdk.model.biophysical.utils.Utils
 static
 method), 352
 needs_data_refresh() (al-
 lensdk.brain_observatory.session_api_utils.ParamsMixin
 method), 209
 NEURON_RELEASE_CELL_TYPES_CACHE (allensdk.core.cell_types_cache.ReporterStatus
 attribute), 239
 neurodata_doc (al-
 lensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema
 attribute), 145
 neurodata_doc (al-
 lensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema
 attribute), 145
 neurodata_doc (al-
 lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema
 attribute), 147
 neurodata_doc (al-
 lensdk.brain_observatory.behavior.schemas.OphysEyeTrackingRastersSchema
 attribute), 147
 neurodata_doc (al-
 lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema
 attribute), 148
 neurodata_skip (al-
 lensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema
 attribute), 145
 neurodata_skip (al-
 lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema
 attribute), 147
 neurodata_skip (al-
 lensdk.brain_observatory.behavior.schemas.OphysEyeTrackingRastersSchema
 attribute), 148
 neurodata_type (al-
 lensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema
 attribute), 145
 neurodata_type (al-
 lensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema
 attribute), 145
 neurodata_type (al-
 lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema
 attribute), 147
 neurodata_type (al-
 lensdk.brain_observatory.behavior.schemas.OphysEyeTrackingRastersSchema
 attribute), 147

attribute), 147
 neurodata_type (al- nocache_json() (al-
 lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema), 148
 attribute), 148
 neurodata_type_inc (al- Node (class in allensdk.internal.morphology.node), 331
 lensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema (allensdk.core.simple_tree.SimpleTree method),
 attribute), 145
 neurodata_type_inc (al- node() (allensdk.core.swc.Morphology method), 263
 lensdk.brain_observatory.behavior.schemas.BehaviorTaskParameterSchema (allensdk.internal.core.simpletree.SimpleTree
 attribute), 145
 neurodata_type_inc (al- node() (allensdk.internal.morphology.morphology.Morphology
 lensdk.brain_observatory.behavior.schemas.OphysBehaviorMethodSchema
 attribute), 147
 neurodata_type_inc (al- node_ids() (allensdk.core.simple_tree.SimpleTree
 lensdk.brain_observatory.behavior.schemas.OphysEyeTrackingRigMetadataSchema (allensdk.internal.core.simpletree.SimpleTree
 attribute), 147
 neurodata_type_inc (al- node_list (allensdk.internal.morphology.morphology.Morphology
 lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema), 326
 attribute), 148
 neuron (allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils lensdk.internal.morphology.morphology.Morphology
 attribute), 349
 neuron_parameter_count() (al- NODE_TYPES (allensdk.core.swc.Morphology attribute),
 lensdk.internal.model.glif.glif_experiment.GlifExperiment 261
 method), 313
 neuronal_model_optimize_dir() (al- attribute), 324
 lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader (allensdk.core.simple_tree.SimpleTree
 method), 281
 NEURONAL_MODEL_PARAMETERS (al- nodes() (allensdk.internal.core.simpletree.SimpleTree
 lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader
 attribute), 281
 neuronal_model_run_dir() (al- lensdk.core.simple_tree.SimpleTree method),
 lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader
 method), 279
 NeuropilMask (class in al- noise_std() (in module al-
 lensdk.brain_observatory.roi_masks), 204
 NeuropilSubtract (class in al- nonraising_ks_2samp() (in module al-
 lensdk.brain_observatory.r_neuropil), 202
 new_image() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities), diff() (in module al-
 367
 nlin() (in module allensdk.internal.model.GLM), 323
 nll_to_pvalue() (in module al- lensdk.ephys.ephys_features), 276
 lensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis.actual_parameters() (al-
 178
 NO_RECONSTRUCTION (allensdk.core.swc.Marker at-
 tribute), 261
 NO_RECONSTRUCTION (al- normalize_F() (in module al-
 lensdk.internal.core.swc.Marker attribute),
 301
 no_response_bias() (in module al-
 lensdk.brain_observatory.behavior.criteria),
 135
 nocache_dataframe() (al- np_sitk_convert() (in module al-
 lensdk.api.warehouse_cache.cache.Cache
 lensdk.mouse_connectivity.grid.utilities.image_utilities),

368
 nrn (*allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils* (in module *al-*
attribute), 349 *lensdk.internal.brain_observatory.eye_calibration*),
 289
 nrnivmodl () (*allensdk.internal.model.biophysical.run_optimize.RunOptimize* (in module *al-*
method), 310 *lensdk.internal.brain_observatory.eye_calibration*),
 289
 nrnivmodl () (*allensdk.model.biophysical.run_simulate.RunSimulate* (in module *al-*
method), 350 *one ()* (in module *allensdk*), 371
 num_contingent_trials () (in module *al-* *one_file_call_caching ()* (in module *al-*
lensdk.brain_observatory.behavior.session_metrics), *lensdk.api.warehouse_cache.caching_utilities*),
 149 101
 num_nodes (*allensdk.core.swc.Morphology* *attribute*), *OneOrMoreResultExpectedError*, 282
 263 *OneResultExpectedError*, 371
 num_nodes (*allensdk.internal.morphology.morphology.Morphology* (in module *al-* *lensdk.api.queries.rma_api.RmaApi* *attribute*),
attribute), 327 89
 NUM_ROWS (*allensdk.api.queries.rma_api.RmaApi* *attribute*), *only_except_tabular_clause ()* (*al-*
attribute), 89 *lensdk.api.queries.rma_api.RmaApi* *method*),
 num_trees (*allensdk.core.swc.Morphology* *attribute*), 92
 263 *OntologiesApi* (class in *al-*
 num_trees (*allensdk.internal.morphology.morphology.Morphology* (in module *al-* *lensdk.api.queries.ontologies_api*), 86
attribute), 327 *Ontology* (class in *allensdk.core.ontology*), 249
 number_of_cells (in module *al-* *open ()* (*allensdk.internal.brain_observatory.frame_stream.CvInputStream*
lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet
attribute), 236 *method*), 290
 number_ori (*allensdk.brain_observatory.drifting_gratings.DriftingGratings* (in module *al-* *open ()* (*allensdk.internal.brain_observatory.frame_stream.FfmpegInputS*
attribute), 196 *method*), 291
 number_ori (*allensdk.brain_observatory.static_gratings.StaticGratings* (in module *al-* *open ()* (*allensdk.internal.brain_observatory.frame_stream.FfmpegOutputS*
attribute), 211 *method*), 291
 number_phase (*allensdk.brain_observatory.static_gratings.StaticGratings* (in module *al-* *open ()* (*allensdk.internal.brain_observatory.frame_stream.FrameInputStr*
attribute), 211 *method*), 291
 number_scenes (in module *al-* *open ()* (*allensdk.internal.brain_observatory.frame_stream.FrameOutputS*
lensdk.brain_observatory.natural_scenes.NaturalScenes (in module *al-* *method*), 291
attribute), 199 *Scenes_corona_plot ()* (*al-*
lensdk.brain_observatory.natural_scenes.NaturalScenes
method), 199
 number_sf (*allensdk.brain_observatory.static_gratings.StaticGratings* (in module *al-* *open_fan_plot ()* (*al-*
attribute), 211 *method*), 211
 number_tf (*allensdk.brain_observatory.drifting_gratings.DriftingGratings* (in module *al-* *lensdk.brain_observatory.static_gratings.StaticGratings*
attribute), 196 *method*), 211
 numbercells (*allensdk.brain_observatory.stimulus_analysis.stimulus_analysis.StimulusAnalysis* (in module *al-* *open_stimulus_analysis_plot ()* (*al-*
attribute), 213 *lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise*
method), 197
 NWB_FILE_TYPE (in module *al-* *lensdk.api.queries.brain_observatory_api.BrainObservatoryApi* *plot ()* (*al-*
lensdk.api.queries.brain_observatory_api.BrainObservatoryApi *attribute*), 70 *lensdk.brain_observatory.drifting_gratings.DriftingGratings*
method), 196
 NWB_FILE_TYPE (in module *al-* *open_track_plot ()* (*al-*
lensdk.api.queries.cell_types_api.CellTypesApi *method*), 198
attribute), 73 *lensdk.brain_observatory.natural_movie.NaturalMovie*
 NWB_FILE_TYPE (*allensdk.api.queries.glif_api.GlifApi* (in module *al-* *method*), 198
attribute), 76 *open_view_on_binary ()* (in module *al-*
lensdk.internal.brain_observatory.ophys_session_decomposition
method), 293
 NwbDataSet (class in *allensdk.core.nwb_data_set*),
 246
 NwbOphysMetadataSchema (class in *al-* *OPHYS_ANALYSIS_FILE_TYPE* (in module *al-*
lensdk.brain_observatory.behavior.schemas), *lensdk.api.queries.brain_observatory_api.BrainObservatoryApi*
attribute), 70
 146 *ophys_delta* (*allensdk.internal.pipeline_modules.run_ophys_time_sync*
attribute), 346
 object_norm_eye_coordinates () *OPHYS_EVENTS_FILE_TYPE* (in module *al-*

`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` (al-
 attribute), 70
`ophys_experiment_id` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata` (al-
 attribute), 111
`ophys_frame_rate` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata` (al-
 attribute), 111
`ophys_session_id` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_ophys_metadata` (al-
 attribute), 111
`ophys_times` (`allensdk.internal.pipeline_modules.run_ophys_time_sync` (al-
 attribute), 346
`ophys_timestamps` (al-
`lensdk.internal.brain_observatory.time_sync.OphysTimeAligner` (al-
 attribute), 296
`OphysBehaviorMetadataSchema` (class in al-
`lensdk.brain_observatory.behavior.schemas`),
 146
`OphysEyeTrackingRigMetadataSchema`
 (class in al-
`lensdk.brain_observatory.behavior.schemas`),
 147
`OphysMetadataSchema` (class in al-
`lensdk.brain_observatory.behavior.schemas`),
 147
`OphysMixin` (class in al-
`lensdk.brain_observatory.behavior.behavior_project_cacheable_ophys_mixin`),
 107
`OphysTimeAligner` (class in al-
`lensdk.internal.brain_observatory.time_sync`),
 296
`OptimizeConfigReader` (class in al-
`lensdk.internal.api.queries.optimize_config_reader`),
 281
`optional_lims_inputs()` (in module al-
`lensdk.brain_observatory.argschema_utilities`),
 189
`optional_polys` (al-
`lensdk.mouse_connectivity.grid.subimage.base_subimage.PolygonsSubImage` (al-
 attribute), 363
`optional_polys` (al-
`lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClausoidSubImage` (al-
 attribute), 365
`OPTIONS` (`allensdk.api.queries.rma_api.RmaApi` at-
 tribute), 89
`options_clause()` (al-
`lensdk.api.queries.rma_api.RmaApi` method),
 92
`OPTOGENETIC_STIMULATION_KEYS` (al-
`lensdk.brain_observatory.sync_dataset.Dataset`
 attribute), 218
`opts` (`allensdk.brain_observatory.argschema_utilities.RaisingSchema` (al-
 attribute), 189

`overlap_fraction()` (al- 341
`lensdk.internal.brain_observatory.mask_set.MaskSet`
`method`), 293
`overlaps_motion_border` (al- 345
`lensdk.brain_observatory.roi_masks.Mask`
`attribute`), 204
P
`pageable()` (in module al-
`lensdk.api.queries.rma_pager`), 94
`pager()` (`allensdk.api.queries.rma_pager.RmaPager`
`static method`), 94
`panel_size` (`allensdk.brain_observatory.stimulus_info.Monitor`
`attribute`), 214
`ParamsMixin` (class in al-
`lensdk.brain_observatory.session_api_utils`),
209
`parent()` (`allensdk.core.simple_tree.SimpleTree`
`method`), 256
`parent()` (`allensdk.internal.core.simpletree.SimpleTree`
`method`), 301
`parent_id()` (`allensdk.core.simple_tree.SimpleTree`
`method`), 257
`parent_id()` (`allensdk.internal.core.simpletree.SimpleTree`
`method`), 301
`parent_ids()` (`allensdk.core.simple_tree.SimpleTree`
`method`), 257
`parent_of()` (`allensdk.core.swc.Morphology`
`method`), 263
`parent_of()` (`allensdk.internal.morphology.morphology.Morphology`
`method`), 327
`parents()` (`allensdk.core.simple_tree.SimpleTree`
`method`), 257
`parse_age_in_days()` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata`
`static method`), 109
`parse_arguments()` (in module al-
`lensdk.internal.model.glif.find_sweeps`), 313
`parse_arguments()` (in module al-
`lensdk.model.glif.simulate_neuron`), 361
`parse_command_line_args()` (al-
`lensdk.config.app.application_config.ApplicationConfig`
`method`), 223
`parse_cre_line()` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata`
`static method`), 109
`parse_indicator()` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata`
`static method`), 110
`parse_input()` (in module al-
`lensdk.internal.pipeline_modules.run_demixing`),
341
`parse_input()` (in module al-
`lensdk.internal.pipeline_modules.run_dff_computation`),
341
`parse_input()` (in module al-
`lensdk.internal.pipeline_modules.run_observatory_thumbnails`),
345
`parse_input()` (in module al-
`lensdk.internal.pipeline_modules.run_ophys_session_decomposition`),
345
`parse_input_data()` (in module al-
`lensdk.internal.pipeline_modules.run_ophys_eye_calibration`),
345
`parse_neuron_output()` (in module al-
`lensdk.internal.model.biophysical.passive_fitting.neuron_utils`),
308
`parse_obj()` (in module `allensdk.core.obj_utilities`),
248
`parse_reporter_line()` (al-
`lensdk.brain_observatory.behavior.metadata.behavior_metadata`
`static method`), 110
`parse_stim_repr()` (in module al-
`lensdk.brain_observatory.ecephys.stimulus_table.stimulus_params`),
174
`parser_for_extension()` (al-
`lensdk.config.model.description_parser.DescriptionParser`
`method`), 227
`password` (`allensdk.core.authentication.DbCredentials`
`attribute`), 231
`paste_slice()` (al-
`lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`
`method`), 369
`path_to_image()` (al-
`lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`
`method`), 369
`path_to_list()` (al-
`lensdk.core.structure_tree.StructureTree` `static`
`method`), 260
`pathfinder()` (`allensdk.api.warehouse_cache.cache.Cache`
`static method`), 98
`pause_metrics()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
`method`), 267
`peak` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`attribute`), 213
`peak_run` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`attribute`), 213
`period()` (`allensdk.brain_observatory.sync_dataset.Dataset`
`method`), 220
`phasevals` (`allensdk.brain_observatory.static_gratings.StaticGratings`
`attribute`), 211
`PHOTODIODE_KEYS` (al-
`lensdk.brain_observatory.sync_dataset.Dataset`
`attribute`), 218
`PIPE` (`allensdk.api.queries.rma_api.RmaApi` `attribute`),
89
`pipe_stage()` (`allensdk.api.queries.rma_api.RmaApi`
`method`), 89

method), 93

PIPELINE_DATASET (al-
lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet
attribute), 232

PipelineModule (class in al-
lensdk.internal.core.lims_pipeline_module),
298

pixel_counter (al-
lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage
attribute), 364

pixel_size (allensdk.brain_observatory.stimulus_info.Monitor
attribute), 214

pixels_to_visual_degrees() (al-
lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor
method), 213

pixels_to_visual_degrees() (al-
lensdk.brain_observatory.stimulus_info.Monitor
method), 214

plot() (allensdk.brain_observatory.circle_plots.CoronaPlotter
method), 191

plot() (allensdk.brain_observatory.circle_plots.FanPlotter
method), 191

plot() (allensdk.brain_observatory.circle_plots.TrackPlotter
method), 191

plot_all() (allensdk.brain_observatory.sync_dataset.Dataset
method), 220

plot_bit() (allensdk.brain_observatory.sync_dataset.Dataset
method), 220

plot_bits() (allensdk.brain_observatory.sync_dataset.Dataset
method), 220

plot_blocks() (in module al-
lensdk.brain_observatory.ecephys.stimulus_table.visualization), 168

plot_cell_correlation() (in module al-
lensdk.brain_observatory.observatory_plots),
200

plot_cell_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_cell_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures3),
305

plot_cell_receptive_field() (al-
lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise
method), 197

plot_chi_square_summary() (in module al-
lensdk.brain_observatory.receptive_field_analysis.visualization),
185

plot_combined_speed() (in module al-
lensdk.brain_observatory.observatory_plots),
200

plot_condition_histogram() (in module al-
lensdk.brain_observatory.observatory_plots),
200

plot_direction_selectivity() (al-
lensdk.brain_observatory.drifting_gratings.DriftingGratings
method), 196

plot_drifting_traces() (in module al-
lensdk.brain_observatory.brain_observatory_plotting),
189

plot_ellipses() (in module al-
lensdk.brain_observatory.receptive_field_analysis.visualization),
185

plot_fi_curve_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_fi_curve_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures3), 305

plot_gaussian_fit() (in module al-
lensdk.brain_observatory.receptive_field_analysis.visualization),
186

plot_hero_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_hero_figures() (in module al-
lensdk.internal.ephys.plot_qc_figures3),
305

plot_images() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_images() (in module al-
lensdk.internal.ephys.plot_qc_figures3),
305

plot_instantaneous_threshold_thumbnail() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_instantaneous_threshold_thumbnail() (in module al-
lensdk.internal.ephys.plot_qc_figures3),
305

plot_line() (allensdk.brain_observatory.sync_dataset.Dataset
method), 221

plot_lines() (allensdk.brain_observatory.sync_dataset.Dataset
method), 221

plot_long_square_summary() (in module al-
lensdk.internal.ephys.plot_qc_figures), 303

plot_long_square_summary() (in module al-
lensdk.internal.ephys.plot_qc_figures3), 306

plot_lsn_traces() (in module al-
lensdk.brain_observatory.brain_observatory_plotting),
189

plot_mask() (in module al-
lensdk.brain_observatory.receptive_field_analysis.visualization),
186

plot_mask_outline() (in module al-
lensdk.brain_observatory.observatory_plots),
200

plot_masks() (in module al-

`lensdk.internal.brain_observatory.demix_report()`, 200
`285`
`plot_mean_waveforms()` (in module `al-`
`lensdk.brain_observatory.ecephys.visualization`), 201
`175`
`plot_msr_summary()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 303
`186`
`plot_negative_baselines()` (in module `al-`
`lensdk.brain_observatory.demixer`), 193
`plot_negative_baselines()` (in module `al-`
`lensdk.internal.brain_observatory.demixer`), 201
`286`
`plot_negative_transients()` (in module `al-`
`lensdk.brain_observatory.demixer`), 193
`plot_negative_transients()` (in module `al-`
`lensdk.internal.brain_observatory.demixer`), 201
`286`
`plot_ns_traces()` (in module `al-`
`lensdk.brain_observatory.brain_observatory_plotting`), 213
`190`
`plot_onetrace()` (in module `al-`
`lensdk.brain_observatory.dff`), 195
`plot_orientation_selectivity()` (al-
`lensdk.brain_observatory.drifting_gratings.DriftingGratings`
`method`), 196
`plot_orientation_selectivity()` (al-
`lensdk.brain_observatory.static_gratings.StaticGratings`
`method`), 211
`plot_overlap_masks_lengthOne()` (in module
`allensdk.brain_observatory.demixer`), 193
`plot_overlap_masks_lengthOne()` (in module
`allensdk.internal.brain_observatory.demixer`), 201
`286`
`plot_p_values()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 303
`186`
`plot_population_receptive_field()` (al-
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
`method`), 198
`plot_preferred_direction()` (al-
`lensdk.brain_observatory.drifting_gratings.DriftingGratings`
`method`), 196
`plot_preferred_orientation()` (al-
`lensdk.brain_observatory.static_gratings.StaticGratings`
`method`), 211
`plot_preferred_spatial_frequency()` (al-
`lensdk.brain_observatory.static_gratings.StaticGratings`
`method`), 211
`plot_preferred_temporal_frequency()` (al-
`lensdk.brain_observatory.drifting_gratings.DriftingGratings`
`method`), 196
`plot_pupil_location()` (in module `al-`
`lensdk.brain_observatory.observatory_plots`), 201
`plot_radial_histogram()` (in module `al-`
`lensdk.brain_observatory.observatory_plots`), 201
`plot_ramp_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures`), 303
`plot_ramp_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures3`), 306
`plot_receptive_field()` (in module `al-`
`lensdk.brain_observatory.observatory_plots`), 201
`plot_receptive_field_analysis_data()`
`(allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
`method`), 198
`plot_receptive_field_data()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 306
`186`
`plot_representational_similarity()` (al-
`lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`method`), 213
`plot_representational_similarity()`
`(in module`
`lensdk.brain_observatory.observatory_plots`), 201
`plot_rheo_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures`), 303
`plot_rheo_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures3`), 306
`plot_rts_blur_summary()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 306
`186`
`plot_rts_summary()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 306
`186`
`plot_running_a()` (in module `al-`
`lensdk.brain_observatory.brain_observatory_plotting`), 213
`plot_running_speed()` (in module `al-`
`lensdk.brain_observatory.visualization`), 187
`plot_running_speed_histogram()` (al-
`lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`method`), 213
`plot_sag_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures`), 303
`plot_sag_figures()` (in module `al-`
`lensdk.internal.ephys.plot_qc_figures3`), 306
`plot_orientation_selectivity_cumulative_histogram()`
`(in module`
`lensdk.brain_observatory.observatory_plots`), 201

`plot_sg_traces()` (in module `allensdk.brain_observatory.brain_observatory_plotting`), 190
`plot_short_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 303
`plot_short_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 306
`plot_single_ap_values()` (in module `allensdk.internal.ephys.plot_qc_figures`), 304
`plot_single_ap_values()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 306
`plot_speed()` (in module `allensdk.brain_observatory.observatory_plots`), 201
`plot_speed_tuning()` (in module `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` method), 213
`plot_spike_counts()` (in module `allensdk.brain_observatory.ecephys.visualization`), 175
`plot_subthreshold_long_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 304
`plot_subthreshold_long_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 306
`plot_sweep_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 304
`plot_sweep_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 306
`plot_sweep_set_summary()` (in module `allensdk.internal.ephys.plot_qc_figures`), 304
`plot_sweep_set_summary()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 307
`plot_sweep_value_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 304
`plot_sweep_value_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 307
`plot_time_to_peak()` (in module `allensdk.brain_observatory.natural_scenes.NaturalScenes` method), 200
`plot_time_to_peak()` (in module `allensdk.brain_observatory.static_gratings.StaticGratings` method), 211
`plot_time_to_peak()` (in module `allensdk.brain_observatory.observatory_plots`), 201
`plot_traces()` (in module `allensdk.brain_observatory.demixer`), 193
`plot_traces()` (in module `allensdk.internal.brain_observatory.demixer`), 286
`plot_transients()` (in module `allensdk.brain_observatory.demixer`), 193
`plot_transients()` (in module `allensdk.internal.brain_observatory.demixer`), 286
`plotLineRegress1()` (in module `allensdk.internal.model.glif.plotting`), 318
`plotLineRegress1()` (in module `allensdk.internal.model.glif.spike_cutting`), 319
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.plotting`), 318
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.spike_cutting`), 319
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.plotting`), 318
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.spike_cutting`), 319
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.plotting`), 318
`polar_line_circles()` (in module `allensdk.brain_observatory.circle_plots`), 192
`polar_linspace()` (in module `allensdk.brain_observatory.circle_plots`), 192
`polar_to_xy()` (in module `allensdk.brain_observatory.circle_plots`), 192
`PolarPlotter` (class in `allensdk.brain_observatory.circle_plots`), 191
`PolygonSubImage` (class in `allensdk.mouse_connectivity.grid.subimage.base_subimage`), 363
`populate_stimulus_table()` (in module `allensdk.brain_observatory.drifting_gratings.DriftingGratings` method), 196
`populate_stimulus_table()` (in module `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 198
`populate_stimulus_table()` (in module `allensdk.brain_observatory.natural_movie.NaturalMovie` method), 198
`populate_stimulus_table()` (in module `allensdk.brain_observatory.natural_scenes.NaturalScenes` method), 200
`populate_stimulus_table()` (in module `allensdk.brain_observatory.static_gratings.StaticGratings` method), 211
`populate_stimulus_table()` (in module `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` method), 213
`population_correlation_scatter()` (in module `allensdk.brain_observatory.observatory_plots`), 201
`port` (`allensdk.core.authentication.DbCredentials` attribute), 231
`POSITIVE` (`allensdk.core.cell_types_cache.ReporterStatus` attribute), 239

`post_process_cr()` (in module `allensdk.internal.brain_observatory.itracker_utils`), 292
`post_process_pupil()` (in module `allensdk.internal.brain_observatory.itracker_utils`), 292
`PostgresQueryMixin` (class in `allensdk.internal.api`), 282
`postprocess()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ProjectTable`), 108
`postprocess_additional()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ExperimentsTable`), 107
`postprocess_additional()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ExperimentsTable`), 365
`postprocess_additional()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ExperimentsTable`), 107
`postprocess_additional()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ExperimentsTable`), 366
`postprocess_base()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ProjectTable`), 108
`postprocess_unionizes()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ProjectTable`), 268
`postprocess_unionizes()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ProjectTable`), 333
`postprocess_unionizes()` (in module `allensdk.brain_observatory.behavior.behavior_project_cache.tables.projects_table.ProjectTable`), 335
`pre_blank_sec` (in module `allensdk.brain_observatory.ecephys.file_io.stim_files_comp_stim_files_pickle_stim_files`), 164
`prepare_nwb_output()` (in module `allensdk.model.biophysical.runner`), 350
`PREVIEW` (in module `allensdk.api.queries.rma_api.RmaApi`), 89
`print_node()` (in module `allensdk.core.swc.Compartment`), 261
`print_out()` (in module `allensdk.ephys.feature_extractor.EphysFeatures`), 277
`print_summary()` (in module `allensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis`), 242
`ProbeSynchronizer` (class in `allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer`), 160
`process()` (in module `allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor`), 265
`process_eye_tracking_data()` (in module `allensdk.brain_observatory.behavior.eye_tracking_processing`), 140
`process_inputs()` (in module `allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fitting`), 308

propagate() (allensdk.internal.mouse_connectivity.interval_unionize.tissuocyte_unionize_record.TissuocyteBaseUnionize method), 334

propagate() (allensdk.internal.mouse_connectivity.interval_unionize_record.Unionize method), 335

propagate_record() (allensdk.brain_observatory.circle_plots), 192

propagate_record() (al-lensdk.brain_observatory.circle_plots), 192

propagate_record() (al-lensdk.internal.mouse_connectivity.interval_unionize.IntervalUnionizer class method), 333

propagate_record() (al-lensdk.brain_observatory.argschema_utilities), 188

propagate_record() (al-lensdk.internal.mouse_connectivity.interval_unionize.tissuocyte_unionize.TissuocyteUnionizer class method), 335

propagate_to_bilateral() (al-lensdk.brain_observatory.behavior.schemas), 148

propagate_to_bilateral() (al-lensdk.internal.mouse_connectivity.interval_unionize.IntervalUnionizer class method), 333

propagate_unionizes() (al-lensdk.brain_observatory.argschema_utilities), 189

propagate_unionizes() (al-lensdk.internal.mouse_connectivity.interval_unionize.IntervalUnionizer class method), 333

provide() (allensdk.core.authentication.CredentialProvider method), 231

provide() (allensdk.core.authentication.EnvCredentialProvider method), 231

psycpg2_select() (in module al-lensdk.internal.api), 283

pupil_position_in_mouse_eye_coordinates() (al-lensdk.brain_observatory.eye_calibration.EyeCalibration method), 288

pupil_position_in_mouse_eye_coordinates() (allensdk.internal.brain_observatory.eye_calibration.EyeCalibration method), 288

pupil_position_on_monitor_in_cm() (al-lensdk.internal.brain_observatory.eye_calibration.EyeCalibration method), 288

pupil_position_on_monitor_in_degrees() (al-lensdk.internal.brain_observatory.eye_calibration.EyeCalibration method), 288

pupil_position_on_monitor_in_degrees() (allensdk.internal.brain_observatory.eye_calibration.EyeCalibration method), 288

push_summary() (al-lensdk.ephys.feature_extractor.EphysFeatureExtractor method), 277

pval (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis attribute), 213

pvalue_to_NLL() (in module al-lensdk.brain_observatory.receptive_field_analysis.chisquare), 181

pvalue_to_NLL() (in module al-lensdk.brain_observatory.receptive_field_analysis.visualization), 186

PycfgDescriptionParser (class in al-lensdk.config.model.formats.pycfg_description_parser), 225

Q

query() (in module al-lensdk.internal.core.lims_utilities), 299

quote_string() (al-lensdk.api.queries.rma_api.RmaApi method), 93

radial_arcs() (in module al-lensdk.brain_observatory.circle_plots), 192

radial_circles() (in module al-lensdk.brain_observatory.circle_plots), 192

RaisingSchema (class in al-lensdk.brain_observatory.argschema_utilities), 188

RaisingSchema (class in al-lensdk.brain_observatory.behavior.schemas), 148

RaisingSchema (class in al-lensdk.brain_observatory.argschema_utilities), 189

RaisingSchema (class in al-lensdk.brain_observatory.behavior.schemas), 148

ramps_features() (al-lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor method), 266

randomize_parameter_values() (al-lensdk.internal.model.glif.glif_optimizer.GlifOptimizer method), 315

rank_structures() (al-lensdk.core.mouse_connectivity_cache.MouseConnectivityCache method), 246

ransac_fit() (allensdk.internal.brain_observatory.fit_ellipse.FitEllipse method), 290

raster_plot() (in module al-lensdk.brain_observatory.ecephys.visualization), 176

rasterize_polygons() (in module al-lensdk.mouse_connectivity.grid.utilities.image_utilities), 368

ratio_and_pyramid() (in module al-lensdk.mouse_connectivity.grid.writers), 369

read() (allensdk.config.model.description_parser.DescriptionParser method), 227

read() (allensdk.config.model.formats.hdf5_util.Hdf5Util method), 224

read() (allensdk.config.model.formats.json_description_parser.JsonDescriptionParser method), 224

read() (allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser method), 225

read() (in module allensdk.core.json_utilities), 241

read() (in module al-lensdk.internal.mouse_connectivity.interval_unionize.data_utilities), 332

read_cell_index_receptive_field_analysis() (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise static method), 198

read_data() (allensdk.api.api.Api method), 104

read_file() (allensdk.core.json_utilities.JsonComments

[illegible]

<code>lensdk.core.reference_space_cache.ReferenceSpaceCache</code>	<code>lensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage</code>
<code>attribute</code>), 253	<code>attribute</code>), 364
<code>ReferenceSpace</code> (class in <code>lensdk.core.reference_space</code>), 250	<code>required_polys</code> (al-
<code>ReferenceSpaceApi</code> (class in <code>lensdk.api.queries.reference_space_api</code>), 87	<code>lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage</code>
<code>ReferenceSpaceCache</code> (class in <code>lensdk.core.reference_space_cache</code>), 252	<code>attribute</code>), 365
<code>refine_threshold_indexes()</code> (in module <code>lensdk.ephys.ephys_features</code>), 276	<code>required_polys</code> (al-
<code>relative_path_from_url()</code> (in module <code>lensdk.api.cloud_cache.utils</code>), 67	<code>lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage</code>
<code>remove_comments()</code> (al-	<code>attribute</code>), 366
<code>lensdk.core.json_utilities.JsonComments</code>	<code>required_segmentations</code> (al-
<code>class method</code>), 240	<code>lensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentSubImage</code>
<code>remove_keys()</code> (al-	<code>attribute</code>), 364
<code>lensdk.api.warehouse_cache.cache.Cache</code>	<code>required_segmentations</code> (al-
<code>static method</code>), 99	<code>attribute</code>), 365
<code>remove_lfp_noise()</code> (in module <code>lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling</code>), 164	<code>required_segmentations</code> (al-
<code>lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling</code>	<code>lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage</code>
<code>remove_lfp_offset()</code> (in module <code>lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling</code>), 164	<code>attribute</code>), 366
<code>lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling</code>	<code>resample_into_volume()</code> (in module <code>lensdk.mouse_connectivity.grid.utilities.image_utilities</code>), 368
<code>remove_multiline_comments()</code> (al-	<code>resample_into_volume()</code> (al-
<code>lensdk.core.json_utilities.JsonComments</code>	<code>lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder</code>
<code>class method</code>), 240	<code>method</code>), 370
<code>remove_unassigned()</code> (al-	<code>resample_into_volume()</code> (in module <code>lensdk.mouse_connectivity.grid.utilities.image_utilities</code>), 368
<code>lensdk.core.reference_space.ReferenceSpace</code>	<code>resave_swc()</code> (in module <code>lensdk.internal.morphology.validate_swc</code>), 331
<code>method</code>), 252	<code>reset()</code> (<code>allensdk.model.glif.glif_neuron.GlifNeuron</code>
<code>rename_columns()</code> (al-	<code>method</code>), 354
<code>lensdk.api.warehouse_cache.cache.Cache</code>	<code>reset_AScurrent_none()</code> (in module <code>lensdk.model.glif.glif_neuron_methods</code>), 359
<code>static method</code>), 99	<code>reset_AScurrent_sum()</code> (in module <code>lensdk.model.glif.glif_neuron_methods</code>), 359
<code>renames()</code> (<code>allensdk.core.structure_tree.StructureTree</code>	<code>reset_hex_pack()</code> (in module <code>lensdk.brain_observatory.circle_plots</code>), 192
<code>static method</code>), 260	<code>reset_long_squares_start()</code> (in module <code>lensdk.brain_observatory.circle_plots</code>), 192
<code>Report</code> (class in <code>lensdk.internal.model.biophysical.make_deap_fit_job</code>), 310	<code>reset_threshold_inf()</code> (in module <code>lensdk.model.glif.glif_neuron_methods</code>), 359
<code>reporter_line</code> (al-	<code>reset_threshold_three_components()</code> (in module <code>lensdk.model.glif.glif_neuron_methods</code>), 359
<code>lensdk.brain_observatory.behavior.metadata.behavior_metadata</code>	<code>reset_voltage_zero()</code> (in module <code>lensdk.model.glif.glif_neuron_methods</code>), 360
<code>attribute</code>), 110	<code>reset_voltage_zero()</code> (al-
<code>ReporterStatus</code> (class in <code>lensdk.core.cell_types_cache</code>), 239	<code>reshape_response_array()</code> (al-
<code>required_intensities</code> (al-	
<code>lensdk.mouse_connectivity.grid.subimage.base_subimage.IntensitySubImage</code>	
<code>attribute</code>), 363	
<code>required_intensities</code> (al-	
<code>lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage</code>	
<code>attribute</code>), 365	
<code>required_polys</code> (al-	
<code>lensdk.mouse_connectivity.grid.subimage.base_subimage.PolySubImage</code>	
<code>attribute</code>), 363	
<code>required_polys</code> (al-	

`lensdk.brain_observatory.drifting_gratings.DriftingGratings` (attribute), 278
`method`), 196
`reshape_response_array()` (al- `RmaApi` (class in `allensdk.api.queries.rma_api`), 89
`lensdk.brain_observatory.natural_scenes.NaturalScenes` 94
`method`), 200
`reshape_response_array()` (al- `RmaTemplate` (class in al-
`lensdk.api.queries.rma_template`), 94
`lensdk.brain_observatory.static_gratings.StaticGratings` (in module al-
`method`), 211
`resolve_initial_image()` (in module al- `roi_id` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`lensdk.brain_observatory.behavior.trials_processing`), attribute), 213
155
`resolve_paths()` (`allensdk.config.manifest.Manifest` `lensdk.brain_observatory.roi_masks`), 204
`method`), 229
`response` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` (class in al-
`attribute`), 213
`rolling_window()` (in module al-
`response_bias()` (in module al- `lensdk.brain_observatory.demixer`), 193
`lensdk.brain_observatory.behavior.session_metrics`), `rolling_window()` (in module al-
149 286
`retinotopy_metric()` (in module al- `root` (`allensdk.core.swc.Morphology` attribute), 264
`lensdk.internal.brain_observatory.annotated_region_metrics`), `rotate()` (`allensdk.internal.mouse_connectivity.projection_thumbnail.volume`
284 338
`retrieve_file_from_storage()` (al- `rotate()` (in module al-
`lensdk.internal.api.api_prerelease.ApiPrerelease` 216
`method`), 282
`retrieve_file_over_http()` (al- `rotate_and_extract()` (al-
`lensdk.api.api.Api` method), 104
`retrieve_parsed_json_over_http()` (al- `lensdk.internal.mouse_connectivity.projection_thumbnail.volume`
`lensdk.api.api.Api` method), 105
`retrieve_xml_over_http()` (`allensdk.api.api.Api` `method`), 338
`method`), 105
`return_mask_cb()` (al- `rotate_ray()` (in module al-
`lensdk.core.reference_space.ReferenceSpace` `lensdk.internal.brain_observatory.itracker_utils`),
`static method`), 252 292
`reward_rate()` (in module al- `rotate_vector()` (in module al-
`lensdk.brain_observatory.behavior.trial_masks`), `lensdk.internal.brain_observatory.fit_ellipse`),
149 290
`rings_in_hex_pack()` (in module al- `row_from_cell_id()` (al-
`lensdk.brain_observatory.circle_plots`), 192
`rma_templates` (al- `lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`lensdk.api.queries.biophysical_api.BiophysicalApi` (in module al-
`attribute`), 70 `method`), 213
`rma_templates` (al- `run()` (`allensdk.internal.model.glif.glif_experiment.GlifExperiment`
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `method`), 313
`attribute`), 73 `run()` (`allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.S`
`method`), 338
`rma_templates` (`allensdk.api.queries.glif_api.GlifApi` `run()` (`allensdk.model.glif.glif_neuron.GlifNeuron`
`attribute`), 77 `method`), 355
`rma_templates` (al- `run()` (in module al-
`lensdk.api.queries.image_download_api.ImageDownloadApi` `lensdk.internal.mouse_connectivity.interval_unionize.run_tissuec`
`attribute`), 81 334
`rma_templates` (al- `run()` (in module al-
`lensdk.api.queries.ontologies_api.OntologiesApi` `lensdk.internal.mouse_connectivity.projection_thumbnail.generat`
`attribute`), 87 `method`), 336
`rma_templates` (al- `run_base_model()` (al-
`lensdk.api.queries.biophysical_module_api.BiophysicalModuleApi` `lensdk.internal.model.glif.glif_experiment.GlifExperiment`
`attribute`), 87 `method`), 314
`rma_templates` (al- `run_many()` (`allensdk.internal.model.glif.glif_optimizer.GlifOptimizer`
`lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi` `method`), 314

[run_module\(\)](#) (in module `allensdk.internal.core.lims_pipeline_module`), 298
[run_module\(\)](#) (in module `allensdk.model.biophys_sim.bps_command`), 349
[run_once\(\)](#) (`allensdk.internal.model.glif.glif_optimizer.GlifOptimizer` method), 315
[run_once_bound\(\)](#) (`allensdk.internal.model.glif.glif_optimizer.GlifOptimizer` method), 315
[run_ophys_time_sync\(\)](#) (in module `allensdk.internal.pipeline_modules.run_ophys_time_sync`), 347
[run_passive_fit\(\)](#) (in module `allensdk.internal.model.biophysical.run_passive_fit`), 311
[run_postprocessing\(\)](#) (in module `allensdk.brain_observatory.receptive_field_analysis.postprocessing`), 183
[run_session_analysis\(\)](#) (in module `allensdk.brain_observatory.session_analysis`), 208
[run_subimage\(\)](#) (in module `allensdk.mouse_connectivity.grid.subimage`), 366
[run_subimage\(\)](#) (in module `allensdk.mouse_connectivity.grid.subimage.base_subimage`), 364
[run_sync\(\)](#) (in module `allensdk.model.biophysical.runner`), 350
[run_until_biological_spike\(\)](#) (`allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron` method), 316
[run_with_biological_spikes\(\)](#) (`allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron` method), 316
[RunningSpeed](#) (class in `allensdk.brain_observatory.running_speed`), 206
[RunOptimize](#) (class in `allensdk.internal.model.biophysical.run_optimize`), 310
[RunSimulate](#) (class in `allensdk.model.biophysical.run_simulate`), 350
[RunSimulateLims](#) (class in `allensdk.internal.model.biophysical.run_simulate_lims`), 311
S
[safe_factory\(\)](#) (`allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector.VolumeProjector` class method), 338
[safe_make_parent_dirs\(\)](#) (`allensdk.config.manifest.Manifest` class method), 229
[safe_mkdir\(\)](#) (`allensdk.config.manifest.Manifest` class method), 229
[safe_system_path\(\)](#) (in module `allensdk.internal.core.lims_utilities`), 299
[SAG_TARGET](#) (`allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` attribute), 265
[save_reconv\(\)](#) (in module `allensdk.internal.model.GLM`), 323
[sample_freq](#) (`allensdk.brain_observatory.sync_dataset.Dataset` attribute), 221
[sample_frequency](#) (`allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset` attribute), 163
[sampling_rate_scale](#) (`allensdk.brain_observatory.ecephys.align_timestamps.probe_sync_utils` attribute), 160
[save\(\)](#) (`allensdk.core.swc.Morphology` method), 264
[save\(\)](#) (`allensdk.internal.morphology.morphology.Morphology` method), 327
[save_analysis_arrays\(\)](#) (`allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 236
[save_analysis_dataframes\(\)](#) (`allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 236
[save_cell_index_receptive_field_analysis\(\)](#) (`allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` static method), 198
[save_ephys_data\(\)](#) (`allensdk.api.queries.cell_types_api.CellTypesApi` method), 75
[save_figure\(\)](#) (in module `allensdk.internal.ephys.plot_qc_figures`), 304
[save_figure\(\)](#) (in module `allensdk.internal.ephys.plot_qc_figures3`), 307
[save_nwb\(\)](#) (in module `allensdk.model.biophysical.runner`), 351
[save_ophys_experiment_analysis_data\(\)](#) (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 73
[save_ophys_experiment_data\(\)](#) (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 73
[save_ophys_experiment_event_data\(\)](#) (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 73
[save_ophys_experiment_eye_gaze_data\(\)](#) (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 73
[save_qc_figures\(\)](#) (in module `allensdk.internal.ephys.plot_qc_figures`), 304

`lensdk.internal.ephys.core_feature_extract()`, 302
`save_reconstruction()` (al-
`lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 75
`save_reconstruction_markers()` (al-
`lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 75
`save_session_a()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 207
`save_session_b()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 207
`save_session_c()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 207
`save_session_c2()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 207
`save_voltage()` (al-
`lensdk.core.dat_utilities.DatUtilities`
`method`), 239
`schema()` (al-
`lensdk.api.queries.connected_services.ConnectedServices`
`class method`), 76
`score_feature_set()` (al-
`lensdk.ephys.feature_extractor.EphysFeatureExtractor`
`method`), 277
`search()` (al-
`lensdk.brain_observatory.stimulus_info.BinaryIntervalSearch`
`method`), 213
`search()` (al-
`lensdk.brain_observatory.stimulus_info.StimulusSearch`
`method`), 214
`section_image_query()` (al-
`lensdk.api.queries.image_download_api.ImageDownloadApi`
`method`), 81
`SegmentationSubImage` (class in al-
`lensdk.mouse_connectivity.grid.subimage.base_subimage`), 363
`select()` (al-
`lensdk.internal.api.PostgresQueryMixin`
`method`), 283
`select()` (in module al-
`lensdk.internal.core.lims_utilities`), 299
`select_channels()` (in module al-
`lensdk.brain_observatory.ecephys.lfp_subsampling_subsampling`), 165
`select_one()` (al-
`lensdk.internal.api.PostgresQueryMixin`
`method`), 283
`serialize()` (al-
`lensdk.brain_observatory.behavior.image_api.ImageApi`
`static method`), 141
`SERVICE` (al-
`lensdk.api.queries.rma_api.RmaApi`
`attribute`), 89
`service_query()` (al-
`lensdk.api.queries.rma_api.RmaApi`
`method`), 93
`service_stage()` (al-
`lensdk.api.queries.rma_api.RmaApi`
`method`), 93
`session_a()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 208
`session_b()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 208
`session_c()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 208
`session_c2()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 208
`session_na` (al-
`lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api`
`attribute`), 162
`session_type` (al-
`lensdk.brain_observatory.behavior.metadata.behavior.metadata`
`attribute`), 110
`SessionAnalysis` (class in al-
`lensdk.brain_observatory.session_analysis`), 206
`sessions_are_equal()` (in module al-
`lensdk.brain_observatory.session_api_utils`), 210
`sessions_with_stimulus()` (in module al-
`lensdk.brain_observatory.stimulus_info`), 216
`set_all_parameters()` (al-
`lensdk.internal.model.biophysical.deap_utils.Utils`
`method`), 309
`set_api_urls()` (al-
`lensdk.api.api.Api`
`method`), 105
`set_apical_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 328
`set_apical_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 328
`set_apical_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 328
`set_apical_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 328
`set_parse_grid_parameters()` (al-
`lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`
`method`), 370
`set_credential_provider()` (in module al-
`lensdk.core.authentication`), 232
`set_default_working_directory()` (al-
`lensdk.api.api.Api`
`method`), 105
`set_subsampling()` (al-
`lensdk.brain_observatory.circle_plots.CoronaPlotter`
`method`), 191
`set_dims()` (al-
`lensdk.brain_observatory.circle_plots.FanPlotter`
`method`), 191
`set_image_api()` (al-
`lensdk.brain_observatory.r_neuropil.NeuropilSubtract`
`method`), 202
`set_iclamp_params()` (al-
`lensdk.internal.model.biophysical.deap_utils.Utils`
`method`), 309
`set_image_spacing()` (in module al-
`lensdk.mouse_connectivity.grid.utilities.image_utilities`),

`lensdk.model.glif.simulate_neuron`), 361
`simulate_sweep_from_file()` (in module `al-`
`lensdk.model.glif.simulate_neuron`), 361
`sitk_get_center()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
338
`sitk_get_diagonal_length()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
338
`sitk_get_image_parameters()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
338
`sitk_get_size_parity()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
338
`sitk_np_convert()` (in module `al-`
`lensdk.mouse_connectivity.grid.utilities.image_utilities`),
368
`sitk_paste_into_center()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
338
`sitk_safe_ln()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
337
`size()` (`allensdk.internal.brain_observatory.mask_set.MaskSet`
method), 293
`slice_arrays()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize`
method), 335
`smooth()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.utilities`),
185
`smooth_STA()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.utilities`),
181
`sobel_grad()` (in module `al-`
`lensdk.internal.brain_observatory.itracker_utils`),
292
`SOMA` (`allensdk.core.swc.Morphology` attribute), 261
`soma` (`allensdk.core.swc.Morphology` attribute), 264
`SOMA` (`allensdk.internal.morphology.morphology.Morphology`
attribute), 324
`soma_root()` (`allensdk.internal.morphology.morphology.Morphology`
method), 327
`sort_data_arrays()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize`
method), 333
`sort_trials()` (al-
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
method), 198
`spacing` (`allensdk.brain_observatory.behavior.image_api.Image`
attribute), 141
`SPACING` (`allensdk.core.swc.Marker` attribute), 261
`SPACING` (`allensdk.internal.core.swc.Marker` attribute), 261
`sparsify()` (`allensdk.core.swc.Morphology` method),
264
`sparsify()` (`allensdk.internal.morphology.morphology.Morphology`
method), 324
`spatial_frequency_to_pix_per_cycle()`
(`allensdk.brain_observatory.stimulus_info.Monitor`
method), 214
`spike_component_of_threshold_exact()`
(in module `al-`
`lensdk.model.glif.glif_neuron_methods`),
360
`spike_component_of_threshold_forward_euler()`
(in module `al-`
`lensdk.model.glif.glif_neuron_methods`),
360
`spike_feature()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 267
`spike_visualization_utilities.pages()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor`
method), 268
`spike_visualization_utilities` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 267
`SPIKE_TIMES` (`allensdk.core.nwb_data_set.NwbDataSet`
attribute), 247
`size_in_bytes()` (`allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 267
`spiral_trials()` (in module `al-`
`lensdk.brain_observatory.circle_plots`), 192
`spiral_trials_polar()` (in module `al-`
`lensdk.brain_observatory.circle_plots`), 192
`square_column()` (in module `al-`
`lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spike`
method), 171
`standardize_movie_numbers()` (in module `al-`
`lensdk.brain_observatory.ecephys.stimulus_table.naming_utilities`
method), 172
`start()` (`allensdk.internal.model.biophysical.passive_fitting.output_grab`
method), 308
`START_ROW` (`allensdk.api.queries.rma_api.RmaApi` at-
tribute), 89
`start_specimen()` (al-
`lensdk.internal.model.biophysical.run_optimize.RunOptimize`
method), 260
`StaticGratings` (class in `al-`
`lensdk.brain_observatory.static_gratings`),
324
`stats()` (`allensdk.brain_observatory.sync_dataset.Dataset`
method), 221
`stim_table` (`allensdk.brain_observatory.stimulus_analysis.StimulusAna`
attribute), 213
`stim_table_to_categories()` (in module `al-`

[lensdk.brain_observatory.chisquare_categorical\(\)](#), 211
[lensdk.brain_observatory.chisquare_categorical](#), 190
[StimulusImage](#) (class in [al-](#)
[lensdk.brain_observatory.behavior.stimulus_processing.stimulus_](#)
[lensdk.internal.brain_observatory.time_sync.OphysTimeAlign](#), 226
[attribute](#)), 297
[StimulusImageFactory](#) (class in [al-](#)
[lensdk.brain_observatory.ecephys.file_io.stim_file.CamStimulusPickleSaveFile](#), 126
[attribute](#)), 164
[stimuli_in_session\(\)](#) (in module [al-](#) [StimulusSearch](#) (class in [al-](#)
[lensdk.brain_observatory.stimulus_info](#)), [lensdk.brain_observatory.stimulus_info](#)),
216, 214
[stimulus_alignment](#) (al- [StimulusTemplate](#) (class in [al-](#)
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#), 126
[attribute](#)), 346
[stimulus_amplitude\(\)](#) (al- [StimulusTemplateFactory](#) (class in [al-](#)
[lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor](#) [lensdk.brain_observatory.behavior.stimulus_processing.stimulus_](#)
[method](#)), 267, 126
[STIMULUS_CONTENT_TYPE](#) (al- [stitch\(\)](#) ([allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#)
[attribute](#)), 278, 278
[Stitcher](#) (class in [al-](#)
[STIMULUS_CONTENT_TYPE](#) (al- [lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)),
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#)
[attribute](#)), 281
[stop\(\)](#) ([allensdk.internal.model.biophysical.passive_fitting.output_grabber](#)
[method](#)), 308
[stimulus_delay](#) (al- [method](#)), 308
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#), 346
[attribute](#)), 346
[STIMULUS_PRERELEASE_KEY](#) (al-
[lensdk.internal.core.mouse_connectivity_cache_prerelease.Mouse](#)
[attribute](#)), 300
[stimulus_delta](#) (al- [method](#)), 300
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#), 346
[attribute](#)), 346
[stream_zip_directory_over_http\(\)](#) (in module [al-](#)
[lensdk.api.api](#)), 105
[stimulus_file_entries\(\)](#) (al- [stream_zip_directory_over_http\(\)](#) (in mod-
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#), 105
[method](#)), 279
[STRING](#) ([allensdk.api.queries.connected_services.ConnectedServices](#)
[attribute](#)), 76
[stimulus_file_entries\(\)](#) (al- [method](#)), 281
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) types () (al-
[method](#)), 281, [lensdk.core.swc.Morphology](#) method), 264
[stimulus_frame_rate](#) (al- [strip_all_other_types\(\)](#) (al-
[lensdk.brain_observatory.behavior.metadata.behavior_metadata_extractor.Morphology](#), 327
[attribute](#)), 110, [lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#) method), 327
[stimulus_path\(\)](#) (al- [strip_type\(\)](#) ([allensdk.core.swc.Morphology](#)
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#)
[method](#)), 279, [strip_type\(\)](#) ([allensdk.internal.morphology.morphology.Morphology](#)
[method](#)), 327
[stimulus_path\(\)](#) (al- [method](#)), 327
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) reads_from () (al-
[method](#)), 281, [lensdk.core.ontology.Ontology](#) method),
[stimulus_search](#) (al- 250
[lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet](#) from () (al-
[attribute](#)), 236, [lensdk.core.structure_tree.StructureTree](#)
[method](#)), 260
[STIMULUS_TABLE_TYPES](#) (al- [method](#)), 260
[lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet](#) (al-
[attribute](#)), 232, [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[attribute](#)), 253
[stimulus_times](#) (al- [method](#)), 253
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#) KEY (al-
[attribute](#)), 346, [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[attribute](#)), 253
[StimulusAnalysis](#) (class in [al-](#) [attribute](#)), 253
[lensdk.brain_observatory.stimulus_analysis](#)), [STRUCTURE_TREE_KEY](#) (al-

`lensdk.core.reference_space_cache.ReferenceSpaceCache` (class in `allensdk.core.typing`), 265
`attribute`), 253
`SvgApi` (class in `allensdk.api.queries.svg_api`), 94
STRUCTURE_UNIONIZES_KEY (alias for `SWC_FILE_TYPE`) (alias for `sweep_entries()`) (alias for `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` (class in `allensdk.core.cache`), 242
`attribute`), 242
`CellTypesApi` (class in `allensdk.api.queries.cell_types_api`), 74
STRUCTURES_KEY (alias for `sweep_entries()`) (alias for `lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader` (class in `allensdk.internal.module_reader`), 279
`method`), 279
StructureTree (class in `allensdk.core.structure_tree`), 258
`sweep_entries()` (alias for `lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader` (class in `allensdk.internal.module_reader`), 281
`method`), 281
stumpify_axon() (alias for `allensdk.core.swc.Morphology` (class in `allensdk.core.swc`), 264
`method`), 264
`sweep_feature()` (alias for `lensdk.epphys.epphys_extractor.EphysSweepFeatureExtractor` (class in `allensdk.epphys.feature_extractor`), 267
`method`), 267
stumpify_axon() (alias for `lensdk.internal.morphology.morphology.Morphology` (class in `allensdk.internal.morphology`), 328
`method`), 328
`sweep_feature_keys()` (alias for `lensdk.epphys.epphys_extractor.EphysSweepFeatureExtractor` (class in `allensdk.epphys.feature_extractor`), 268
`method`), 268
SubImage (class in `allensdk.mouse_connectivity.grid.subimage.base_subimage`), 364
`sweep_features()` (alias for `lensdk.epphys.epphys_extractor.EphysSweepSetFeatureExtractor` (class in `allensdk.epphys.set_extractor`), 268
`method`), 268
SubjectMetadataSchema (class in `allensdk.brain_observatory.behavior.schemas`), 148
`sweep_numbers()` (alias for `lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader` (class in `allensdk.internal.module_reader`), 279
`method`), 279
subsampling_data() (in module `allensdk.internal.model.data_access`), 324
`sweep_numbers()` (alias for `lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader` (class in `allensdk.internal.module_reader`), 282
`method`), 282
subsampling_lfp() (in module `allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling_lfp`), 165
`sweep_numbers_by_type()` (alias for `lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader` (class in `allensdk.internal.module_reader`), 279
`method`), 279
subsampling_timestamps() (in module `allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling_timestamps`), 166
`method`), 279
SUBTHRESH_MAX_AMP (alias for `lensdk.epphys.epphys_extractor.EphysCellFeatureExtractor` (class in `allensdk.epphys.cell_extractor`), 199
`attribute`), 199
sum_pixel_intensity (alias for `lensdk.brain_observatory.natural_movie.NaturalMovie` (class in `allensdk.brain_observatory.natural_movie`), 213
`attribute`), 213
sum_pixels (alias for `lensdk.brain_observatory.localizer.Localizer` (class in `allensdk.brain_observatory.localizer`), 198
`attribute`), 198
sum_projection_pixel_intensity (alias for `lensdk.brain_observatory.natural_movie.NaturalMovie` (class in `allensdk.brain_observatory.natural_movie`), 213
`attribute`), 213
sum_projection_pixels (alias for `lensdk.brain_observatory.natural_scenes.NaturalScenes` (class in `allensdk.brain_observatory.natural_scenes`), 200
`attribute`), 200
sum_statistics (alias for `lensdk.brain_observatory.static_gratings.StaticGratings` (class in `allensdk.brain_observatory.static_gratings`), 211
`attribute`), 211
summarize() (alias for `lensdk.epphys.epphys_extractor.EphysSweepSetFeatureExtractor` (class in `allensdk.epphys.set_extractor`), 269
`method`), 269
SUMMARY_STRUCTURE_SET_ID (alias for `SynchronizationApi` (class in `allensdk.core.synchronization`), 95
`attribute`), 242
summer_over() (in module `allensdk.brain_observatory.behavior.criteria`), 136
`synthesize_F()` (in module `allensdk.brain_observatory.r_neuropil`), 203
SUPPORTED_PIPELINE_VERSION (alias for `T` (class in `allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (class in `allensdk.core.brain_observatory_nwb_data_set`), 232
`attribute`), 232
`table` (alias for `allensdk.brain_observatory.behavior.behavior_project_cache.table` (class in `allensdk.brain_observatory.behavior.project_cache`), 108
`attribute`), 108

TABULAR	(allensdk.api.queries.rma_api.RmaApi	at-	to_dataframe()	(al-
tribute), 89			lensdk.brain_observatory.behavior.stimulus_processing.stimulus_	
targeted_structure	(al-		method), 126	
	lensdk.brain_observatory.behavior.metadata.behavior_ophys		method), 110	
attribute), 111			method), 110	
tau_m	(allensdk.model.glif.glif_neuron.GlifNeuron	at-	to_dict()	(allensdk.brain_observatory.behavior.metadata.behavior_ophys
tribute), 355			method), 111	
temp_dir()	(in module	al-	to_dict()	(allensdk.internal.model.glif.glif_optimizer.GlifOptimizer
	lensdk.test_utilities.temp_dir), 370		method), 315	
TEMPLATE_KEY	(allensdk.core.reference_space_cache.ReferenceSpaceCache	al-	to_dict()	(allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptim
attribute), 253			method), 317	
template_projection()	(in module	al-	to_dict()	(allensdk.internal.morphology.morphology.Morphology
	lensdk.internal.mouse_connectivity.projection_thumbnail.projection_thumbnail), 337		to_dict()	(allensdk.internal.morphology.node.Node
template_query()	(al-		method), 331	
	lensdk.api.queries.rma_template.RmaTemplate	to-	to_dict()	(allensdk.model.glif.glif_neuron.GlifNeuron
method), 94			method), 355	
test()	(allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.GlifNeuronMethods.GlifNeuronMethod		method), 162	
			method), 356	
test_fit()	(in module	al-	to_filter_rhs()	(al-
	lensdk.internal.brain_observatory.fit_ellipse),		lensdk.api.queries.rma_template.RmaTemplate	
290			method), 94	
TestNode	(class in	al-	to_manifest()	(al-
	lensdk.internal.morphology.validate_swc),		lensdk.internal.api.queries.biophysical_module_reader.Biophysic	
331			method), 279	
tfvals	(allensdk.brain_observatory.drifting_gratings.DriftingGratings		test()	(al-
attribute), 196			lensdk.internal.api.queries.optimize_config_reader.OptimizeConf	
Tile	(class in allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile), 339		total_voxel_counts()	(al-
timestamps	(allensdk.brain_observatory.running_speed.RunningSpeed		lensdk.core.reference_space.ReferenceSpace	
attribute), 206			method), 252	
timestamps	(allensdk.brain_observatory.stimulus_analysis.stimulus_analysis		lensdk.core.reference_space.ReferenceSpace	(al-
attribute), 213			attribute), 252	
TimeSyncOutputs	(class in	al-	TimeSyncPlotter	(class in
	lensdk.internal.pipeline_modules.run_ophys_time_sync),		lensdk.brain_observatory.circle_plots), 191	al-
346				
TimeSyncWriter	(class in	al-	TrainingLabelClassifier	(class in
	lensdk.internal.pipeline_modules.run_ophys_time_sync),		lensdk.internal.brain_observatory.roi_filter_utils),	al-
346			294	
TissuecyteBaseUnionize	(class in	al-	TrainingMultiLabelClassifier	(class in
	lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_base_unionize), 334		lensdk.internal.brain_observatory.roi_filter_utils),	al-
334			294	
TissuecyteInjectionUnionize	(class in	al-	transform	(allensdk.mouse_connectivity.grid.image_series_gridder.Image
	lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_injection_unionize), 335		translate_image_and_fill()	(in module
335			lensdk.brain_observatory.stimulus_info), 216	al-
TissuecyteProjectionUnionize	(class in	al-	lensdk.core.h5_utilities), 240	
	lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_projection_unionize), 335			
335				
TissuecyteUnionizer	(class in	al-	tree()	(allensdk.core.swc.Morphology
	lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer), 335		method), 264	
335			method), 328	
to_config_string()	(al-		TreeSearchApi	(class in
	lensdk.config.app.application_config.ApplicationConfig		lensdk.api.queries.tree_search_api), 96	al-
method), 223			trial_data_from_log()	(in module
				al-

[lensdk.brain_observatory.behavior.trials_processing](#), [method](#)), 226
[155](#) [unpack_structure_set_ancestors\(\)](#) ([al-](#)
[trial_number_limit\(\)](#) ([in module al-](#)
[lensdk.brain_observatory.behavior.dprime](#)), [method](#)), 87
[137](#) [unpack_uint32\(\)](#) ([in module al-](#)
[trial_types\(\)](#) ([in module al-](#)
[lensdk.brain_observatory.behavior.trial_masks](#)), [221](#)
[149](#) [update_data\(\)](#) ([al-](#)
[trim\(\)](#) ([allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile](#)
[method](#)), 340 [lensdk.config.model.description.Description](#)
[method](#)), 226
[trim_border_pulses\(\)](#) ([in module al-](#) [update_default_cell_hoc\(\)](#) ([al-](#)
[lensdk.brain_observatory.ecephys.stimulus_sync](#)), [lensdk.model.biophysical.utils.Utils](#) [method](#)),
[177](#) [352](#)
[trim_discontiguous_times\(\)](#) ([in module al-](#) [update_output_sweep_features\(\)](#) ([in module](#)
[lensdk.brain_observatory.sync_utilities](#)), 187 [allensdk.internal.ephys.core_feature_extract](#)),
[trim_self\(\)](#) ([allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile](#)
[method](#)), 340 [402](#)
[trimmed_stats\(\)](#) ([in module al-](#) [update_well_known_file\(\)](#) ([al-](#)
[lensdk.brain_observatory.ecephys.stimulus_sync](#)), [lensdk.internal.api.queries.biophysical_module_reader.Biophysic](#)
[177](#) [method](#)), 279
[TRUE](#) ([allensdk.api.queries.rma_api.RmaApi](#) [attribute](#)), [lensdk.internal.api.queries.optimize_config_reader.OptimizeConf](#)
[90](#) [method](#)), 282
[tuple_filters\(\)](#) ([al-](#) [upsample_image_to_degrees\(\)](#) ([in module al-](#)
[lensdk.api.queries.rma_api.RmaApi](#) [method](#)), [lensdk.brain_observatory.receptive_field_analysis.utilities](#)),
[94](#) [185](#)
[two_out_of_three_aint_bad\(\)](#) ([in module](#) [url](#) ([allensdk.api.cloud_cache.file_attributes.CacheFileAttributes](#)
[allensdk.brain_observatory.behavior.criteria](#)), [attribute](#)), 66
[136](#) [user](#) ([allensdk.core.authentication.DbCredentials](#) [at-](#)
[TYPE](#) ([allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron](#) [attribute](#)), 231
[attribute](#)), 315
[TYPE](#) ([allensdk.model.glif.glif_neuron.GlifNeuron](#)
[attribute](#)), 354 [Utils](#) ([class in allensdk.internal.model.biophysical.deap_utils](#)),
[309](#)
[Utils](#) ([class in allensdk.model.biophysical.utils](#)), 351

U

V

[union\(\)](#) ([allensdk.internal.brain_observatory.mask_set.MaskSet](#)
[method](#)), 293 [validate_epoch_durations\(\)](#) ([in module al-](#)
[union_size\(\)](#) ([allensdk.internal.brain_observatory.mask_set.MaskSet](#)
[method](#)), 293 [lensdk.brain_observatory.ecephys.stimulus_table.output_validation](#)
[Unionize](#) ([class in al-](#) [validate_epoch_order\(\)](#) ([in module al-](#)
[lensdk.internal.mouse_connectivity.interval_unionize.unionize](#) [record](#)),
[335](#) [validate_mask\(\)](#) ([in module al-](#)
[unit](#) ([allensdk.brain_observatory.behavior.image_api.Image](#) [lensdk.brain_observatory.roi_masks](#)), 206
[attribute](#)), 141 [validate_max_spontaneous_epoch_duration\(\)](#)
[unknown](#) ([allensdk.brain_observatory.argschema_utilities.RaisingSchema.Meta](#) [module](#) [al-](#)
[attribute](#)), 189 [lensdk.brain_observatory.ecephys.stimulus_table.output_validation](#)
[unknown](#) ([allensdk.brain_observatory.behavior.schemas.RaisingSchema.Meta](#)
[attribute](#)), 148 [73](#)
[unpack\(\)](#) ([allensdk.config.model.description.Description](#) [lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync](#)
[method](#)), 226 [method](#)), 346
[unpack_change_log\(\)](#) ([in module al-](#) [validate_structure_id\(\)](#) ([al-](#)
[lensdk.brain_observatory.behavior.stimulus_processing](#)), [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[130](#) [class method](#)), 254
[unpack_manifest\(\)](#) ([al-](#) [validate_structure_ids\(\)](#) ([al-](#)
[lensdk.config.model.description.Description](#) [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)

class method), 254
 validate_structures() (al-
 lensdk.core.reference_space.ReferenceSpace
 method), 252
 validate_swc() (in module al-
 lensdk.internal.morphology.validate_swc),
 331
 validate_swc() (in module al-
 lensdk.morphology.validate_swc), 362
 validate_trial_condition_exclusivity()
 (in module al-
 lensdk.brain_observatory.behavior.trials_processing),
 155
 validate_with_synthetic_F() (in module al-
 lensdk.brain_observatory.r_neuropil), 203
 value_map() (allensdk.core.simple_tree.SimpleTree
 method), 257
 values (allensdk.brain_observatory.running_speed.RunningSpeed
 attribute), 206
 values() (allensdk.brain_observatory.behavior.stimulus_processing.stimulus_templates.StimulusTemplate
 method), 126
 verify_roi_lists_equal() (al-
 lensdk.brain_observatory.session_analysis.SessionAnalysis
 method), 208
 version (allensdk.api.cloud_cache.manifest.Manifest
 attribute), 67
 VERSION (allensdk.config.manifest.Manifest attribute),
 227
 version_id (allensdk.api.cloud_cache.file_attributes.CacheFileAttributes
 attribute), 66
 vin (allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile
 attribute), 164
 visual_degrees_to_pixels() (al-
 lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor
 method), 214
 visual_degrees_to_pixels() (al-
 lensdk.brain_observatory.stimulus_info.Monitor
 method), 214
 voltage_component_of_threshold_exact()
 (in module al-
 lensdk.model.glif.glif_neuron_methods),
 360
 voltage_component_of_threshold_forward_euler()
 (in module al-
 lensdk.model.glif.glif_neuron_methods),
 361
 voltage_deflection() (al-
 lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor
 method), 268
 VolumeProjector (class in al-
 lensdk.internal.mouse_connectivity.projection_thumbnails.volume_utilities.volume_utilities.Morphology
 method), 338
 VOXEL_RESOLUTION_100_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 87
 VOXEL_RESOLUTION_10_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 87
 VOXEL_RESOLUTION_25_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 87
 VOXEL_RESOLUTION_50_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 87
 vsig (allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile
 attribute), 164
 W
 warp_coordinates (al-
 lensdk.brain_observatory.stimulus_info.ExperimentGeometry
 attribute), 214
 warp_image() (allensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor
 method), 214
 warp_stimulus_coordinates() (in module al-
 lensdk.brain_observatory.stimulus_info),
 217
 wedge_ring() (in module al-
 lensdk.brain_observatory.circle_plots), 192
 whitelist() (allensdk.core.structure_tree.StructureTree
 static method), 261
 WhitespaceStrippedString (class in al-
 lensdk.test_utilities.custom_comparators),
 214
 whole_lotta_trials() (in module al-
 lensdk.brain_observatory.behavior.criteria),
 136
 width (allensdk.brain_observatory.stimulus_info.Monitor
 attribute), 214
 window_average() (in module al-
 lensdk.mouse_connectivity.grid.utilities.downsampling_utilities),
 367
 wrap() (allensdk.api.warehouse_cache.cache.Cache
 method), 99
 write() (allensdk.config.model.description_parser.DescriptionParser
 method), 227
 write() (allensdk.config.model.formats.hdf5_util.Hdf5Util
 method), 224
 write() (allensdk.config.model.formats.json_description_parser.JsonDescriptionParser
 method), 225
 write() (allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser
 method), 225
 write() (allensdk.core.swc.Morphology method), 265
 write() (allensdk.internal.brain_observatory.frame_stream.FrameOutputWriter
 method), 291
 write_volume_utilities() (in module al-
 lensdk.internal.morphology.morphology.Morphology
 method), 328
 write() (allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync
 method), 346

`write()` (in module `allensdk.core.json_utilities`), 241

`write_file()` (`allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader` method), 279

`write_file()` (`allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader` method), 282

`write_itksnap_labels()` (`allensdk.core.reference_space.ReferenceSpace` method), 252

`write_json_file()` (`allensdk.config.manifest_builder.ManifestBuilder` method), 230

`write_json_string()` (`allensdk.config.manifest_builder.ManifestBuilder` method), 230

`write_ndarray_with_sitk()` (in module `allensdk.core.sitk_utilities`), 258

`write_or_print_outputs()` (in module `allensdk.brain_observatory.argschema_utilities`), 189

`write_output()` (in module `allensdk.internal.pipeline_modules.run_ophys_eye_calibration`), 345

`write_output_data()` (`allensdk.internal.core.lims_pipeline_module.PipelineModule` method), 298

`write_output_h5()` (`allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncWriter` method), 347

`write_output_json()` (`allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncWriter` method), 347

`write_receptive_field_to_h5()` (in module `allensdk.brain_observatory.receptive_field_analysis.receptive_field`), 184

`write_string()` (`allensdk.config.model.formats.json_description_parser.JsonDescriptionParser` method), 225

`write_string()` (`allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser` method), 225

`write_string()` (in module `allensdk.core.json_utilities`), 241

`write_sweep_response()` (in module `allensdk.model.glif.simulate_neuron`), 362

`write_volume()` (in module `allensdk.mouse_connectivity.grid.utilities.image_utilities`), 368

Y

`yesterday_was_good()` (in module `allensdk.brain_observatory.behavior.criteria`), 136