
Allen SDK Documentation

Release dev

Allen Institute for Brain Science

Oct 14, 2020

Contents

1	Install Guide	1
1.1	Quick Start Using Anaconda	1
1.2	Quick Start Using Pip	1
1.3	Other Distribution Formats	2
1.4	Installation with Docker (Optional)	2
2	Data Resources	3
2.1	Brain Observatory	3
2.2	Cell Types	6
2.3	Mouse Connectivity	9
2.4	Reference Space	11
2.5	API Access	11
2.6	Visual Coding – Neuropixels	14
3	Models	23
3.1	Generalized LIF Models	23
3.2	Biophysical Models	29
4	Examples	39
5	Authors	41
6	allensdk package	43
6.1	Subpackages	43
6.2	Submodules	338
6.3	Module contents	338
7	Allen Brain Observatory	339
8	Allen Cell Types Database	341
9	Allen Mouse Brain Connectivity Atlas	343
10	What’s New - 2.3.1 (October 13, 2020)	345
11	What’s New - 2.3.0 (October 9, 2020)	347
12	What’s New - 2.2.0 (September 3, 2020)	349

13 What's New - 2.1.0 (July 16, 2020)	351
14 What's New - 2.0.0 (June 11, 2020)	353
15 Previous Release Notes	355
Bibliography	357
Python Module Index	359
Index	365

CHAPTER 1

Install Guide

This guide is a resource for using the Allen SDK package. It is maintained by the [Allen Institute for Brain Science](#).

Attention: As of October 2019, we have dropped Python 2 support. The Allen SDK is developed and tested with Python 3.6 and 3.7. We do not guarantee consistent behavior with other Python versions.

1.1 Quick Start Using Anaconda

1. From the [Anaconda downloads page](#), download the Python 3.7 version for your operating system and run the installer.
2. After the installation is complete, open up a terminal (in Windows open Anaconda3 Command Prompt).
3. Install the AllenSDK using PIP:

```
pip install allensdk
```

4. Download one of our many [Jupyter Notebook examples](#) to a new folder.
5. In your terminal, navigate to the directory where you downloaded the Jupyter Notebook example and run the following command:

```
jupyter notebook
```

6. Your browser should open and you should see the Jupyter Notebook example. Enjoy using the Allen SDK!

1.2 Quick Start Using Pip

First ensure you have [pip](#) installed. It is included with the Anaconda distribution.

```
pip install allensdk
```

To uninstall the SDK:

```
pip uninstall allensdk
```

1.3 Other Distribution Formats

The Allen SDK is also available from the Github source repository.

1.4 Installation with Docker (Optional)

[Docker](#) is an open-source technology for building and deploying applications with a consistent environment including required dependencies. The AllenSDK is not distributed as a Docker image, but example Dockerfiles are available.

1. Ensure you have Docker installed.
2. Use Docker to build the image:

```
docker pull alleninstitute/allensdk
```

Other docker configurations are also available under docker directory in the source repository.

3. Run the docker image:

```
docker run -i -t -p 8888:8888 alleninstitute/allensdk /bin/bash
```

4. Run the SDK tests:

```
cd allensdk
make test
```

5. Start a Jupyter Notebook:

```
cd allensdk/doc_template/examples_root/examples/nb
jupyter notebook --ip=* --no-browser --allow-root
```


Using the browser on your host machine, navigate to the path provided by the output from the jupyter notebook command.

The Allen SDK features Python code to support data and model access for the Allen Cell Types Database. Resources for other Allen Brain Atlas data resources will come in future updates.

2.1 Brain Observatory

The [Allen Brain Observatory](#) is a database of the visually-evoked functional responses of neurons in mouse visual cortex based on 2-photon fluorescence imaging. Characterized responses include orientation tuning, spatial and temporal frequency tuning, temporal dynamics, and spatial receptive field structure.

The data is organized into experiments and experiment containers. An experiment container represents a group of experiments with the same targeted imaging area, imaging depth, and Cre line. The individual experiments within an experiment container have different stimulus protocols, but cover the same imaging field of view.

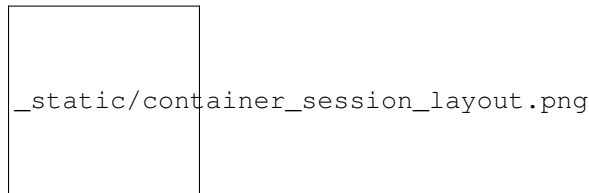

Note: Version 1.3 of scipy fixed an error in its 2 sample Kolmogorov-Smirnoff test implementation. The new version produces more accurate p values for small and medium-sized samples. This change impacts speed tuning analysis p values (as returned by *StimulusAnalysis.get_speed_tuning*). If you access precalculated analysis results via *BrainObservatoryCache.get_ophys_experiment_analysis*, you will see values calculated using an older version of scipy's *ks_2samp*. To access values calculated from the new version, install `scipy>=1.3.0` in your environment and construct a *StimulusAnalysis* object from a *BrainObservatoryNwbDataSet* (as returned by *BrainObservatoryCache.get_ophys_experiment_data*).

Note: Data collected after September 2016 uses a new session C stimulus designed to better-characterize spatial receptive fields in higher visual areas. The original locally sparse noise stimulus used 4.65 visual degree pixels. Session C2 broke that stimulus into two separate stimulus blocks: one with 4.65 degree pixels and one with 9.3 degree pixels.

Note that the `stimulus_info` module refers to these as `locally_sparse_noise_4deg` and `locally_sparse_noise_8deg`, respectively.

For more information on experimental design and a data overview, please visit the [Allen Brain Observatory data portal](#).

2.1.1 Data Processing

For all data in Allen Brain Observatory, we perform the following processing:

1. Segment cell masks from each experiment's 2-photon fluorescence video
2. Associate cells from experiments belonging to the same experiment container and assign unique IDs
3. Extract each cell's mean fluorescence trace
4. Extract mean fluorescence traces from each cell's surrounding neuropil
5. Demix traces from overlapping ROIs
6. Estimate neuropil-corrected fluorescence traces
7. Compute dF/F
8. Compute stimulus-specific tuning metrics

All traces and masks for segmented cells in an experiment are stored in a Neurodata Without Borders (NWB) file. Stored traces include the raw fluorescence trace, neuropil trace, demixed trace, and dF/F trace. Code for extracting neuropil-corrected fluorescence traces, computing dF/F, and computing tuning metrics is available in the SDK.

New in June 2017: Trace demixing is a new addition as of June 2017. All past data was reprocessed using the new demixing algorithm. We have also developed a new module to better characterize a cell's receptive field. Take a look at the [receptive field analysis example notebook](#)

For more information about data processing, please [read the technical whitepapers](#).

2.1.2 Getting Started

The Brain Observatory [Jupyter notebook](#) has many code samples to help get started with the available data:

- Download experimental metadata by visual area, imaging depth, and Cre line
- Find cells with specific response properties, like direction tuning
- Download data for an experiment
- Plot raw fluorescence traces, neuropil-corrected traces, and dF/F
- Find the ROI mask for a given cell
- Run neuropil correction
- Get pupil location and size

The code used to analyze and visualize data in the [Allen Brain Observatory data portal](#) is available as part of the SDK. Take a look at this [Jupyter notebook](#) to find out how to:

- Plot cell's response to its preferred stimulus condition
- Compute a cell's on/off receptive field based on the locally sparse noise stimulus

More detailed documentation is available demonstrating how to:

- Read and visualize the stimulus presentation tables in the NWB files
- Understand the layout of Brain Observatory NWB files

- Map previous cell specimen IDs to current cell specimen IDs

2.1.3 Precomputed Cell Metrics

A large table of precomputed metrics are available for download to support population analysis and filtering. The table below describes all of the metrics in the table. The `get_cell_specimens()` method will download this table as a list of dictionaries which can be converted to a pandas DataFrame as shown in this [Jupyter notebook](#).

Stimulus	Metric	Field Name
drifting gratings	orientation selectivity	osi_dg
	direction selectivity	dsi_dg
	preferred direction	pref_dir_dg
	preferred temporal frequency	pref_tf_dg
	response p value	p_dg
	global ori. selectivity	g_osi_dg
	global dir. selectivity	g_dsi_dg
	response reliability	reliability_dg
	running modulation	run_mod_dg
	running modulation p value	p_run_mod_dg
	pref. condition mean df/f	peak_dff_dg
	TF discrimination index	tfdi_dg
static gratings	orientation selectivity	osi_sg
	preferred orientation	pref_ori_sg
	preferred spatial frequency	pref_sf_sg
	preferred phase	pref_phase_sg
	mean time to peak response	time_to_peak_sg
	response p value	p_sg
	global ori. selectivity	g_osi_sg
	reponse reliability	reliability_sg
	running modulation	run_mod_sg
	running modulation p value	p_run_mod_sg
	pref. condition mean df/f	peak_dff_ns
	SF discrimination index	sfdi_sg
natural scenes	mean time to peak response	time_to_peak_ns
	preferred scene index	pref_scene_ns
	response p value	p_ns
	image selectivity	image_sel_ns
	running modulation	run_mod_ns
	running modulation p value	p_run_mod_ns
	pref. condition mean df/f	peak_dff_ns
natural movie 1	response reliability (session A)	reliability_nm1_a
	response reliability (session B)	reliability_nm1_b
	response reliability (session C)	reliability_nm1_c
natural movie 2	response reliability	reliability_nm2
natural movie 3	response reliability	reliability_nm3
locally sparse noise	RF area (on subunit)	rf_area_on_lsn
	RF area (off subunit)	rf_area_off_lsn
	RF center (on subunit)	rf_center_on_x, rf_center_on_y
	RF center (off subunit)	rf_center_off_x, rf_center_off_y
	RF chi^2	rf_chi2_lsn
	RF on-off subunit distance	rf_distance_lsn
	RF on-off subunit overlap index	rf_overlap_lsn

2.2 Cell Types

The Allen Cell Types data set is a database of mouse and human neuronal cell types based on multimodal characterization of single cells to enable data-driven approaches to classification and is fully integrated with other Allen Brain Atlas resources. The database currently includes:

- **electrophysiology**: whole cell current clamp recordings made from Cre-positive neurons
- **morphology**: 3D bright-field images of the complete structure of neurons from the visual cortex

This page describes how the SDK can be used to access data in the Cell Types Database. For more information, please visit the Cell Types Database [home page](#) and the [API documentation](#).

2.2.1 Examples

The Cell Types [Jupyter notebook](#) has many code samples to help get started with analysis:

- Download and plot stimuli and responses from an NWB file for a cell
- Download and plot a cell's morphological reconstruction
- Download and plot precomputed electrophysiology features
- Download precomputed morphology features to a table
- Compute electrophysiology features for a single sweep

2.2.2 Cell Types Cache

The `CellTypesCache` class provides a Python interface for downloading data in the Allen Cell Types Database into well known locations so that you don't have to think about file names and directories. The following example demonstrates how to download meta data for all cells with 3D reconstructions, then download the reconstruction and electrophysiology recordings for one of those cells:

```
from allensdk.core.cell_types_cache import CellTypesCache

ctc = CellTypesCache(manifest_file='cell_types/manifest.json')

# a list of cell metadata for cells with reconstructions, download if necessary
cells = ctc.get_cells(require_reconstruction=True)

# open the electrophysiology data of one cell, download if necessary
data_set = ctc.get_ephys_data(cells[0]['id'])

# read the reconstruction, download if necessary
reconstruction = ctc.get_reconstruction(cells[0]['id'])
```

`CellTypesCache` takes care of knowing if you've already downloaded some files and reads them from disk instead of downloading them again. All data is stored in the same directory as the `manifest_file` argument to the constructor.

2.2.3 Feature Extraction

The `EphysFeatureExtractor` class calculates electrophysiology features from cell recordings. `extract_cell_features()` can be used to extract the precise feature values available in the Cell Types Database:

```

from allensdk.core.cell_types_cache import CellTypesCache
from allensdk.ephys.extract_cell_features import extract_cell_features
from collections import defaultdict

# initialize the cache
ctc = CellTypesCache(manifest_file='cell_types/manifest.json')

# pick a cell to analyze
specimen_id = 324257146

# download the ephys data and sweep metadata
data_set = ctc.get_ephys_data(specimen_id)
sweeps = ctc.get_ephys_sweeps(specimen_id)

# group the sweeps by stimulus
sweep_numbers = defaultdict(list)
for sweep in sweeps:
    sweep_numbers[sweep['stimulus_name']].append(sweep['sweep_number'])

# calculate features
cell_features = extract_cell_features(data_set,
                                     sweep_numbers['Ramp'],
                                     sweep_numbers['Short Square'],
                                     sweep_numbers['Long Square'])

```

2.2.4 File Formats

This section provides a short description of the file formats used for Allen Cell Types data.

Morphology SWC Files

Morphological neuron reconstructions are available for download as SWC files. The SWC file format is a white-space delimited text file with a standard set of headers. The file lists a set of 3D neuronal compartments, each of which has:

Column	Data Type	Description
id	string	compartment ID
type	integer	compartment type
x	float	3D compartment position (x)
y	float	3D compartment position (y)
z	float	3D compartment position (z)
radius	float	compartment radius
parent	string	parent compartment ID

Comment lines begin with a '#'. Reconstructions in the Allen Cell Types Database can contain the following compartment types:

Type	Description
0	unknown
1	soma
2	axon
3	basal dendrite
4	apical dendrite

The Allen SDK comes with a `swc` Python module that provides helper functions and classes for manipulating SWC files. Consider the following example:

```
import allensdk.core.swc as swc

# if you ran the examples above, you will have a reconstruction here
file_name = 'cell_types/specimen_485909730/reconstruction.swc'
morphology = swc.read_swc(file_name)

# subsample the morphology 3x. root, soma, junctions, and the first child of the root_
↳are preserved.
sparse_morphology = morphology.sparsify(3)

# compartments in the order that they were specified in the file
compartment_list = sparse_morphology.compartment_list

# a dictionary of compartments indexed by compartment id
compartments_by_id = sparse_morphology.compartment_index

# the root soma compartment
soma = morphology.soma

# all compartments are dictionaries of compartment properties
# compartments also keep track of ids of their children
for child in morphology.children_of(soma):
    print(child['x'], child['y'], child['z'], child['radius'])
```

Neurodata Without Borders

The electrophysiology data collected in the Allen Cell Types Database is stored in the [Neurodata Without Borders](#) (NWB) file format. This format, created as part of the [NWB initiative](#), is designed to store a variety of neurophysiology data, including data from intra- and extracellular electrophysiology experiments, optophysiology experiments, as well as tracking and stimulus data. It has a defined schema and metadata labeling system designed so software tools can easily access contained data.

The Allen SDK provides a basic Python class for extracting data from Allen Cell Types Database NWB files. These files store data from intracellular patch-clamp recordings. A stimulus current is presented to the cell and the cell's voltage response is recorded. The file stores both stimulus and response for several experimental trials, here called "sweeps." The following code snippet demonstrates how to extract a sweep's stimulus, response, sampling rate, and estimated spike times:

```
from allensdk.core.nwb_data_set import NwbDataSet

# if you ran the examples above, you will have a NWB file here
file_name = 'cell_types/specimen_485909730/ephys.nwb'
data_set = NwbDataSet(file_name)

sweep_numbers = data_set.get_sweep_numbers()
sweep_number = sweep_numbers[0]
sweep_data = data_set.get_sweep(sweep_number)

# spike times are in seconds relative to the start of the sweep
spike_times = data_set.get_spike_times(sweep_number)

# stimulus is a numpy array in amps
stimulus = sweep_data['stimulus']
```

(continues on next page)

(continued from previous page)

```
# response is a numpy array in volts
reponse = sweep_data['response']

# sampling rate is in Hz
sampling_rate = sweep_data['sampling_rate']

# start/stop indices that exclude the experimental test pulse (if applicable)
index_range = sweep_data['index_range']
```

HDF5 Overview

NWB is implemented in [HDF5](#). HDF5 files provide a hierarchical data storage that mirrors the organization of a file system. Just as a file system has directories and files, and HDF5 file has groups and datasets. The best way to understand an HDF5 (and NWB) file is to open a data file in an HDF5 browser. [HDFView](#) is the recommended browser from the makers of HDF5.

There are HDF5 manipulation libraries for many languages and platforms. MATLAB and Python in particular have strong HDF5 support.

2.3 Mouse Connectivity

The Allen Mouse Brain Connectivity Atlas consists of high-resolution images of axonal projections targeting different anatomic regions or various cell types using Cre-dependent specimens. Each data set is processed through an informatics data analysis pipeline to obtain spatially mapped quantified projection information.

This page describes how to use the SDK to access experimental projection data and metadata. For more information, please visit the Connectivity Atlas [home page](#) and the [API documentation](#)

2.3.1 Structure-Level Projection Data

All AAV projection signal in the Allen Mouse Connectivity Atlas has been registered to the expert-annotated Common Coordinate Framework (CCF) and summarized to structures in the adult mouse structure ontology. Most commonly used for analysis are measures of the density of projection signal in all brain areas for every experiment. This data is available for download and is described in more detail on the [structure unionizes page](#).

2.3.2 Voxel-Level Projection Data

The CCF-registered AAV projection signal is also available for download as a set of 3D volumes for each experiment. The following data volumes are available for download:

- **projection density:** sum of detected projection pixels / sum of all pixels in voxel
- **injection_fraction:** fraction of pixels belonging to manually annotated injection site
- **injection_density:** density of detected projection pixels within the manually annotated injection site
- **data_mask:** binary mask indicating if a voxel contains valid data. Only valid voxels should be used for analysis.

2.3.3 Code Examples

The Mouse Connectivity Jupyter notebook has many code samples to help get started with analysis:

- Download experimental metadata by injection structure and transgenic line
- Download projection signal statistics at a structure level
- Build a structure-to-structure matrix of projection signal values
- Download and visualize gridded projection signal volumes

2.3.4 Mouse Connectivity Cache

The `MouseConnectivityCache` class saves all of the data you can download via the `MouseConnectivityApi` in well known locations so that you don't have to think about file names and directories. It also takes care of knowing if you've already downloaded some files and reads them from disk instead of downloading them again. The following example demonstrates how to download meta data for all experiments with injections in the isocortex and download the projection density volume for one of them:

```
from allensdk.core.mouse_connectivity_cache import MouseConnectivityCache

# tell the cache class what resolution (in microns) of data you want to download
mcc = MouseConnectivityCache(resolution=25)

# use the structure tree class to get information about the isocortex structure
structure_tree = mcc.get_structure_tree()
isocortex_id = structure_tree.get_structures_by_name(['Isocortex'])[0]['id']

# a list of dictionaries containing metadata for non-Cre experiments
experiments = mcc.get_experiments(file_name='non_cre.json',
                                  injection_structure_ids=[isocortex_id])

# download the projection density volume for one of the experiments
pd = mcc.get_projection_density(experiments[0]['id'])
```

2.3.5 File Formats

This section provides a short description of the file formats used for data in the Allen Mouse Connectivity Atlas.

NRRD Files

All of the volumetric data in the connectivity atlas are stored as **NRRD (Nearly Raw Raster Data)** files. A NRRD file consists of a short ASCII header followed by a binary array of data values.

To read these in Python, we recommend the `pynrrd` package. Usage is straightforward:

```
import nrrd

file_name = 'mouse_connectivity/experiment_644250774/projection_density_25.nrrd'
data_array, metadata = nrrd.read(file_name)
```

2.4 Reference Space

Allen Institute atlases and data are registered, when possible, to one of several common reference spaces. Working in such a space allows you to easily compare data across subjects and experimental modalities.

This page documents how to use the Allen SDK to interact with a reference space. For more information and a list of reference spaces, see the [atlas drawings and ontologies API documentation](#) and the [3D reference models API documentation](#). For details about the construction of the Common Coordinate Framework space, see the [CCFv3 whitepaper](#).

2.4.1 Structure Tree

Brain structures in our reference spaces are arranged in trees. The leaf nodes of the tree describe the very fine anatomical divisions of the space, while nodes closer to the root correspond to gross divisions. The `StructureTree` class provides an interface for interacting with a structure tree.

To download a structure tree, use the `allensdk.core.reference_space_cache.ReferenceSpaceCache` class as seen in [this example](#)

2.4.2 Annotation Volumes

An annotation volume is a 3d raster image that segments the reference space into structures. Each voxel in the annotation volume is assigned an integer value that describes the finest structure to which that point in space definitely belongs.

To download a nrrd formatted annotation volume at a specified isometric resolution, use the `allensdk.core.reference_space_cache.ReferenceSpaceCache` class. There is [an example](#) in the notebook.

2.4.3 ReferenceSpaceCache Class

The `allensdk.core.reference_space_cache.ReferenceSpaceCache` class provides a Python interface for downloading structure trees and annotation volumes. It takes care of knowing if you've already downloaded the files and reads them from disk instead of downloading them again.

The class contains methods for working with our reference spaces. Some use cases might include:

- Building an indicator mask for one or more structures
- Viewing the annotation
- Querying the structure graph

Please see the [example notebook](#) for more code samples.

2.5 API Access

The `allensdk.api` package is designed to help retrieve data from the [Allen Brain Atlas API](#). `api` contains methods to help formulate API queries and parse the returned results. There are several pre-made subclasses available that provide pre-made queries specific to certain data sets. Currently there are several subclasses in Allen SDK:

- `CellTypesApi`: data related to the Allen Cell Types Database
- `BiophysicalApi`: data related to biophysical models
- `GlifApi`: data related to GLIF models

- *AnnotatedSectionDataSetsApi*: search for experiments by intensity, density, pattern, and age
- *GridDataApi*: used to download 3-D expression grid data
- *ImageDownloadApi*: download whole or partial two-dimensional images
- *MouseConnectivityApi*: common operations for accessing the Allen Mouse Brain Connectivity Atlas
- *OntologiesApi*: data about neuroanatomical regions of interest
- *ConnectedServices*: schema of Allen Institute Informatics Pipeline services available through the RmaApi
- *RmaApi*: general-purpose HTTP interface to the Allen Institute API data model and services
- *SvgApi*: annotations associated with images as scalable vector graphics (SVG)
- *SynchronizationApi*: data about image alignment
- *TreeSearchApi*: list ancestors or descendants of structure and specimen trees

2.5.1 RMA Database and Service API

One API subclass is the *RmaApi* class. It is intended to simplify constructing an RMA query.

The *RmaApi* is a base class for much of the `allensdk.api.queries` package, but it may be used directly to customize queries or to build queries from scratch.

Often a query will simply request a table of data of one type:

```
from allensdk.api.queries.rma_api import RmaApi

rma = RmaApi()

data = rma.model_query('Atlas',
                       criteria="[name$il'*Mouse*']")
```

This will construct the RMA query url, make the query and parse the resulting JSON into an array of Python dicts with the names, ids and other information about the atlases that can be accessed via the API.

Using the criteria, include and other parameter, specific data can be requested.

```
associations = ''.join(['[id$eq1]',
                       'structure_graph(ontology)', ', ',
                       'graphic_group_labels'])

atlas_data = rma.model_query('Atlas',
                             include=associations,
                             criteria=associations,
                             only=['atlases.id',
                                   'atlases.name',
                                   'atlases.image_type',
                                   'ontologies.id',
                                   'ontologies.name',
                                   'structure_graphs.id',
                                   'structure_graphs.name',
                                   'graphic_group_labels.id',
                                   'graphic_group_labels.name'])
```

Note that a ‘class’ name is used for the first parameter. ‘Association’ names are used to construct the include and criteria parameters nested using parentheses and commas. In the only clause, the ‘table’ form is used, which is

generally a plural lower-case version of the class name. The only clause selects specific ‘fields’ to be returned. The schema that includes the classes, fields, associations and tables can be accessed in JSON form using:

```
# http://api.brain-map.org/api/v2/data.json
schema = rma.get_schema()
for entry in schema:
    data_description = entry['DataDescription']
    clz = list(data_description.keys())[0]
    info = list(data_description.values())[0]
    fields = info['fields']
    associations = info['associations']
    table = info['table']
    print("class: %s" % (clz))
    print("fields: %s" % (','.join(f['name'] for f in fields)))
    print("associations: %s" % (','.join(a['name'] for a in associations)))
    print("table: %s\n" % (table))
```

2.5.2 Using Pandas to Process Query Results

When it is difficult to get data in exactly the required form using only an RMA query, it may be helpful to perform additional operations on the client side. The pandas library can be useful for this.

Data from the API can be read directly into a pandas [Dataframe](#) object.

```
import pandas as pd

structures = pd.DataFrame(
    rma.model_query('Structure',
                    criteria='[graph_id$eq1]',
                    num_rows='all'))
```

Indexing subsets of the data (certain columns, certain rows) is one use of pandas: specifically `.loc`:

```
names_and_acronyms = structures.loc[:, ['name', 'acronym']]
```

and Boolean indexing

```
mea = structures[structures.acronym == 'MEA']
mea_id = mea.iloc[0,:].id
mea_children = structures[structures.parent_structure_id == mea_id]
print(mea_children['name'])
```

`Concatenate`, `merge` and `join` are used to add columns or rows:

When an RMA call contains an include clause, the associated data will be represented as a python dict in a single column. The column may be converted to a proper Dataframe and optionally dropped.

```
criteria_string = "structure_sets[name$eq'Mouse Connectivity - Summary']"
include_string = "ontology"
summary_structures = \
    pd.DataFrame(
        rma.model_query('Structure',
                        criteria=criteria_string,
                        include=include_string,
                        num_rows='all'))

ontologies = \
    pd.DataFrame(
```

(continues on next page)

(continued from previous page)

```
list(summary_structures.ontology)).drop_duplicates()
flat_structures_dataframe = summary_structures.drop(['ontology'], axis=1)
```

Alternatively, it can be accessed using normal python dict and list operations.

```
print(summary_structures.ontology[0]['name'])
```

Pandas Dataframes can be written to a CSV file using `to_csv` and read using `load_csv`.

```
summary_structures[['id',
                    'parent_structure_id',
                    'acronym']].to_csv('summary_structures.csv',
                                      index_label='structure_id')
reread = pd.read_csv('summary_structures.csv')
```

Iteration over a Dataframe of API data can be done in several ways. The `.itertuples` method is one way to do it.

```
for id, name, parent_structure_id in summary_structures[['name',
                                                         'parent_structure_id']].
    .itertuples():
    print("%d %s %d" % (id, name, parent_structure_id))
```

2.5.3 Caching Queries on Disk

`wrap()` has several parameters for querying the API, saving the results as CSV or JSON and reading the results as a pandas dataframe.

```
from allensdk.api.cache import Cache

cache_writer = Cache()
do_cache=True
structures_from_api = \
    cache_writer.wrap(rma.model_query,
                     path='summary.csv',
                     cache=do_cache,
                     model='Structure',
                     criteria='[graph_id$eq1]',
                     num_rows='all')
```

If you change `to_cache` to `False` and run the code again it will read the data from disk rather than executing the query.

2.6 Visual Coding – Neuropixels

The Visual Coding – Neuropixels project uses high-density extracellular electrophysiology (**Ecephys**) probes to record spikes from a wide variety of regions in the mouse brain. Our experiments are designed to study the activity of the visual cortex and thalamus in the context of passive visual stimulation, but these data can be used to address a wide variety of topics.

Spike-sorted data and metadata are available via the AllenSDK as [Neurodata Without Borders](#) files. However, if you're using the AllenSDK to interact with the data, no knowledge of the NWB data format is required.

2.6.1 Getting Started

To jump right in, check out the [quick start guide \(download .ipynb\)](#), which will show you how to download the data, align spikes to a visual stimulus, and decode natural images from neural activity patterns. For a quick summary of experimental design and data access, see the [cheat sheet](#).

If you would like more example code, the [full example notebook \(download .ipynb\)](#) covers all of the ways to access data for each experiment.


Additional tutorials are available on the following topics:

1. [Data access \(download .ipynb\)](#)
2. [Unit quality metrics \(download .ipynb\)](#)
3. [LFP data analysis \(download .ipynb\)](#)
4. [Receptive field mapping \(download .ipynb\)](#)
5. [Optotagging \(download .ipynb\)](#)

For detailed information about the experimental design, data acquisition, and informatics methods, please refer to our [technical whitepaper](#). AllenSDK API documentation is [available here](#).

A note on terminology: Throughout the SDK, we refer to neurons as “units,” because we cannot guarantee that all the spikes assigned to one unit actually originate from a single cell. Unlike in two-photon imaging, where you can visualize each neuron throughout the entire experiment, with electrophysiology we can only “see” a neuron when it fires a spike. If a neuron moves relative to the probe, or if it’s far away from the probe, some of its spikes may get mixed together with those from other neurons. Because of this inherent ambiguity, we provide a variety of quality metrics to allow you to find the right units for your analysis. Even highly contaminated units contain potentially valuable information about brain states, so we didn’t want to leave them out of the dataset. But certain types of analysis require more stringent quality thresholds, to ensure that all of the included units are well isolated from their neighbors.

2.6.2 Data Processing



`_static/neuropixels_data_processing.png`

Neuropixels probes contain 374 or 383 channels that continuously detect voltage fluctuations in the surrounding neural tissue. Each channel is split into two separate data streams, or “bands,” on the probes. The “spike band” is digitized at 30 kHz, and contains information about action potentials fired by neurons directly adjacent to the probe. The “LFP band” is digitized at 2.5 kHz, and records the low-frequency (<1000 Hz) fluctuations that result from synchronized neural activity over a wider area.

To go from the raw spike-band data to NWB files, we perform the following processing steps:

1. Median-subtraction to remove common-mode noise from the continuous traces
2. High-pass filtering (>150 Hz) and whitening across blocks of 32 channels
3. Spike sorting with [Kilosort2](#), to detect spikes and assign them to individual units
4. Computing the mean waveform for each unit
5. Removing units with artifactual waveforms
6. Computing quality metrics for every unit
7. Computing stimulus-specific tuning metrics

For the LFP band, we:


1. Downsample the signals in space and time (every 4th channel and every 2nd sample)
2. High-pass filter at 0.1 Hz to remove the DC offset from each channel
3. Re-reference to channels outside of the brain to remove common-mode noise

The packaged NWB files contain:

1. Spike times, spike amplitudes, mean waveforms, and quality metrics for every unit
2. Information about the visual stimulus
3. Time series of the mouse's running speed, pupil diameter, and pupil position
4. LFP traces for channels in the brain
5. Experiment metadata

All code for data processing and packaging is available in the [ecephys_spike_sorting](#) and the `ecephys` section of the AllenSDK.

2.6.3 Visual Stimulus Sets



`_static/neuropixels_stimulus_sets.png`

A central aim of the Visual Coding – Neuropixels project is to measure the impact of visual stimuli on neurons throughout the mouse visual system. To that end, all mice viewed one of two possible stimulus sets, known as “Brain Observatory 1.1” or “Functional Connectivity”. Both stimulus sets began with a Gabor stimulus flashed at 81 different locations on the screen, used to map receptive fields of visually responsive units. Next, the mice were shown brief flashes of light or dark, to measure the temporal dynamics of the visual response.

The remainder of the visual stimulus set either consisted of the same stimuli shown in the two-photon experiments


(“Brain Observatory 1.1”), or a subset of those stimuli shown with a higher number of repeats. We also added a dot motion stimulus, to allow us to measure the speed tuning of units across the mouse visual system.

2.6.4 AllenSDK 2.0 and Data Compatability

AllenSDK version 2.0 marks a major update to released Visual Coding Neuropixels datasets. Due to newer versions of pynwb/hdmf having issues reading previously released Visual Coding Neuropixels NWB files and due to the significant reorganization of updated NWB file contents, this release contains breaking changes that necessitate a major version revision. NWB files released prior to 6/11/2020 are not guaranteed to work with the 2.0.0 version of AllenSDK. If you cannot or choose not to re-download the updated NWB files, you can continue using a prior version of AllenSDK (< 2.0.0) to access them. However, no further features or bugfixes for AllenSDK (< 2.0.0) are planned. Data released for other projects (Cell Types, Mouse Connectivity, etc.) are *NOT* affected and will *NOT* need to be re-downloaded.

When using the Visual Coding **EcephysProjectCache** from this updated AllenSDK version, if a **ManifestError** is encountered, this indicates that previously downloaded cached data files need to be removed and re-downloaded. The location these files as well as manifest are user defined and are set when instantiating an **EcephysProjectCache**.

2.6.5 Quality Metrics



`_static/neuropixels_quality_metrics.png`

Every NWB file includes a table of quality metrics, which can be used to assess the completeness, contamination, and stability of units in the recording. By default, we won't show you units below a pre-determined quality threshold; we hide any units that are not present for the whole session (`presence_ratio < 0.95`), that include many contaminating spikes (`isi_violations > 0.5`), or are likely missing a large fraction of spikes (`amplitude_cutoff > 0.1`). However, even contaminated or incomplete units contain information about brain states, and may be of interest to analyze. Therefore, the complete units table can be accessed via special flags in the AllenSDK.

In general, we do not make a distinction between ‘single-unit’ and ‘multi-unit’ activity. There is no obvious place to draw a boundary in the overall distributions of quality metrics, and setting a strict cutoff (e.g. `isi_violations = 0`) will remove a lot of potentially valuable data. We prefer to leave it up to the end user to decide what level of contamination is tolerable. But that means you need to be aware that different units will have different levels of cleanliness.

It should also be noted that all of these metrics assume that the spike waveform is stable throughout the experiment. Given that the probe drifts, on average, about 40 microns over the course of the ~3 hour recordings, this assumption is almost never valid. The resulting changes in waveform shape can cause a unit’s quality to fluctuate. If you’re unsure about a unit’s quality, it can be helpful to plot its spike amplitudes over time. This can make it obvious if it’s drifting below threshold, or if it contains spikes from multiple neurons.

Documentation on the various quality metrics can be found in the [ecephys_spike_sorting](#) repository.

For a detailed discussion of the appropriate way to apply each of these metrics, please check out [this tutorial](#) (download [.ipynb](#))

2.6.6 Precomputed Stimulus Metrics

Tables of precomputed metrics are available for download to support population analysis and filtering. The table below describes all of the available metrics. The `get_unit_analysis_metrics()` method will load this table as a [pandas DataFrame](#).

Stimulus	Metric	Field Name
drifting gratings	preferred orientation	pref_ori_dg
	preferred temporal frequency	pref_tf_dg
	global ori. selectivity	g_osi_dg
	global dir. selectivity	g_dsi_dg
	running modulation	run_mod_dg
	running modulation p-value	p_run_mod_dg
	firing rate	firing_rate_dg
	fano factor	fano_dg
	modulation index	mod_idx_dg
	f1/f0	f1_f0_dg
	lifetime sparseness	lifetime_sparseness_dg
	c50 (contrast tuning stimulus)	c50_dg
static gratings	preferred orientation	pref_ori_sg
	preferred spatial frequency	pref_sf_sg
	preferred phase	pref_phase_sg
	global ori. selectivity	g_osi_sg
	running modulation	run_mod_sg
	running modulation p-value	p_run_mod_sg
	firing rate	firing_rate_sg
	fano factor	fano_sg
	lifetime sparseness	lifetime_sparseness_sg
natural scenes	preferred image index	pref_image_ns
	image selectivity	image_selectivity_ns
	running modulation	run_mod_ns
	running modulation p-value	p_run_mod_ns
	firing rate	firing_rate_ns
	fano factor	fano_factor_ns
	lifetime sparseness	lifetime_sparseness_ns
dot motion	preferred speed	pref_speed_dm
	preferred direction	pref_dir_dm

Continued on next page

Table 2 – continued from previous page

Stimulus	Metric	Field Name
	running modulation	run_mod_dm
	running modulation p-value	p_run_mod_dm
	firing rate	firing_rate_dm
	fano factor	fano_factor_dm
	lifetime sparseness	lifetime_sparseness_dm
full-field flashes	on/off ratio	on_off_ratio_fl
	running modulation	run_mod_fl
	running modulation p-value	p_run_mod_fl
	firing rate	firing_rate_fl
	fano factor	fano_factor_fl
	lifetime sparseness	lifetime_sparseness_fl
gabors	RF area	area_rf
	RF elevation	elevation_rf
	RF azimuth	azimuth_rf
	RF p-value	p_value_rf
	running modulation	run_mod_rf
	running modulation p-value	p_run_mod_rf
	firing rate	firing_rate_rf
	fano factor	fano_factor_rf
	lifetime sparseness	lifetime_sparseness_rf

The Allen SDK currently focuses on models generated from electrophysiology data in the Allen Cell Types Database. There are two classes of models available for download: biophysical models and generalize leaky integrate-and-fire models.

3.1 Generalized LIF Models

The Allen Cell Types Database contains Generalized Leaky Integrate and Fire (GLIF) models that simulate the firing behavior of neurons at five levels of complexity. Review the GLIF technical [white paper](#) for details on these models and how their parameters were optimized.

The Allen SDK GLIF simulation module is an explicit time-stepping simulator that evolves a neuron's simulated voltage over the course of an input current stimulus. The module also tracks the neuron's simulated spike threshold and registers action potentials whenever voltage surpasses threshold. Action potentials initiate reset rules that update voltage, threshold, and (optionally) trigger afterspike currents.

The GLIF simulator in this package has a modular architecture that enables users to choose from a number of dynamics and reset rules that update the simulation's voltage, spike threshold, and afterspike currents during the simulation. The GLIF package contains a built-in set of rules, however developers can plug in custom rule implementations provided they follow a simple argument specification scheme.

The Allen SDK GLIF simulator was developed and tested with Python 2.7.9, installed as part of [Anaconda Python](#) distribution version 2.1.0.

The rest of this page provides examples demonstrating how to download models, examples of simulating these models, and general GLIF model documentation.

Note: the GLIF simulator module is still under heavy development and may change significantly in the future.

3.1.1 Downloading GLIF Models

There are two ways to download files necessary to run a GLIF model. The first way is to visit <http://celltypes.brain-map.org> and find cells that have GLIF models available for download. The electrophysiology details page for a cell has a neuronal model download link. Specifically:

1. Click ‘More Options +’ and filter for GLIF models.
2. Click the electrophysiology thumbnail for a cell on the right hand panel.
3. Choose a GLIF model from the ‘Show model responses’ dropdown.
4. Scroll down to the model response click ‘Download model’.

One such link (for a simple LIF neuronal model, ID 566302806), would look like this:

```
http://api.brain-map.org/neuronal_model/download/566302806
```

This link returns .zip archive containing the neuron configuration file and sweep metadata required to simulate the model with stimuli applied to the cell. Specifically, the .zip archive will contain:

- **472423251_neuron_config.json**: JSON config file for the GLIF model
- **ephys_sweeps.json**: JSON with metadata for sweeps presented to the cell
- **neuronal_model.json**: JSON with general metadata for the cell

If you would like to reproduce the model traces seen in the Cell Types Database, you can download an NWB file containing both the stimulus and cell response traces via a ‘Download data’ link on the cell’s electrophysiology page. See the [NWB](#) description section for more details on the NWB file format.

You can also download all of these files, including the cell’s NWB file, using the *GlifApi* class:

```
from allensdk.api.queries.glif_api import GlifApi
from allensdk.core.cell_types_cache import CellTypesCache
import allensdk.core.json_utilities as json_utilities

neuronal_model_id = 566302806

# download model metadata
glif_api = GlifApi()
nm = glif_api.get_neuronal_models_by_id([neuronal_model_id])[0]

# download the model configuration file
nc = glif_api.get_neuron_configs([neuronal_model_id])[neuronal_model_id]
neuron_config = glif_api.get_neuron_configs([neuronal_model_id])
json_utilities.write('neuron_config.json', neuron_config)

# download information about the cell
ctc = CellTypesCache()
ctc.get_ephys_data(nm['specimen_id'], file_name='stimulus.nwb')
ctc.get_ephys_sweeps(nm['specimen_id'], file_name='ephys_sweeps.json')
```

3.1.2 Running a GLIF Simulation

To run a GLIF simulation, the most important file you need is the `neuron_config` JSON file. You can use this file to instantiate a simulator and feed in your own stimulus:

```

import allensdk.core.json_utilities as json_utilities
from allensdk.model.glif.glif_neuron import GlifNeuron

# initialize the neuron
neuron_config = json_utilities.read('neuron_config.json')['566302806']
neuron = GlifNeuron.from_dict(neuron_config)

# make a short square pulse. stimulus units should be in Amps.
stimulus = [ 0.0 ] * 100 + [ 10e-9 ] * 100 + [ 0.0 ] * 100

# important! set the neuron's dt value for your stimulus in seconds
neuron.dt = 5e-6

# simulate the neuron
output = neuron.run(stimulus)

voltage = output['voltage']
threshold = output['threshold']
spike_times = output['interpolated_spike_times']

```

Note: The GLIF simulator does not simulate during action potentials. Instead it inserts NaN values for a fixed number of time steps when voltage surpasses threshold. The simulator skips `neuron.spike_cut_length` time steps after voltage surpasses threshold.

To reproduce the model's traces displayed on the Allen Cell Types Database web page, the Allen SDK provides the `allensdk.core.model.glif.simulate_neuron` module for simulating all sweeps presented to a cell and storing them in the NWB format:

```

import allensdk.core.json_utilities as json_utilities

from allensdk.model.glif.glif_neuron import GlifNeuron
from allensdk.model.glif.simulate_neuron import simulate_neuron

neuron_config = json_utilities.read('neuron_config.json')['566302806']
ephys_sweeps = json_utilities.read('ephys_sweeps.json')
ephys_file_name = 'stimulus.nwb'

neuron = GlifNeuron.from_dict(neuron_config)

sweep_numbers = [ s['sweep_number'] for s in ephys_sweeps if s['stimulus_units'] ==
↳ 'Amps' ]
sweep_numbers = sweep_numbers[:1] # for the sake of a speedy example, just run the_
↳ first one
simulate_neuron(neuron, sweep_numbers, ephys_file_name, ephys_file_name, 0.05)

```

Warning: These stimuli are sampled at a very high resolution (200kHz), and a given cell can have many sweeps. This process can take over an hour.

The `simulate_neuron` function call simulates all sweeps in the NWB file. Because the same NWB file is being used for both input and output, the cell's response traces will be overwritten as stimuli are simulated. `simulate_neuron` optionally accepts a value which will be used to overwrite these NaN values generated during action potentials (in this case 0.05 Volts).

If you would like to run a single sweep instead of all sweeps, try the following:

```
import allensdk.core.json_utilities as json_utilities
from allensdk.model.glif.glif_neuron import GlifNeuron
from allensdk.core.nwb_data_set import NwbDataSet

neuron_config = json_utilities.read('neuron_config.json')['566302806']
ephys_sweeps = json_utilities.read('ephys_sweeps.json')
ephys_file_name = 'stimulus.nwb'

# pull out the stimulus for the current-clamp first sweep
ephys_sweep = next( s for s in ephys_sweeps
                    if s['stimulus_units'] == 'Amps' )
ds = NwbDataSet(ephys_file_name)
data = ds.get_sweep(ephys_sweep['sweep_number'])
stimulus = data['stimulus']

# initialize the neuron
# important! update the neuron's dt for your stimulus
neuron = GlifNeuron.from_dict(neuron_config)
neuron.dt = 1.0 / data['sampling_rate']

# simulate the neuron
output = neuron.run(stimulus)

voltage = output['voltage']
threshold = output['threshold']
spike_times = output['interpolated_spike_times']
```

Note: The dt value provided in the downloadable GLIF neuron configuration files does not correspond to the sampling rate of the original stimulus. Stimuli were subsampled and filtered for parameter optimization. Be sure to overwrite the neuron's dt with the correct sampling rate.

If you would like to plot the outputs of this simulation using numpy and matplotlib, try:

```
import numpy as np
import matplotlib.pyplot as plt

voltage = output['voltage']
threshold = output['threshold']
interpolated_spike_times = output['interpolated_spike_times']
spike_times = output['interpolated_spike_times']
interpolated_spike_voltages = output['interpolated_spike_voltage']
interpolated_spike_thresholds = output['interpolated_spike_threshold']
grid_spike_indices = output['spike_time_steps']
grid_spike_times = output['grid_spike_times']
after_spike_currents = output['AScurrents']

# create a time array for plotting
time = np.arange(len(stimulus))*neuron.dt

plt.figure(figsize=(10, 10))

# plot stimulus
plt.subplot(3,1,1)
plt.plot(time, stimulus)
plt.xlabel('time (s)')
```

(continues on next page)

(continued from previous page)

```

plt.ylabel('current (A)')
plt.title('Stimulus')

# plot model output
plt.subplot(3,1,2)
plt.plot(time, voltage, label='voltage')
plt.plot(time, threshold, label='threshold')

if grid_spike_indices is not None:
    plt.plot(interpolated_spike_times, interpolated_spike_voltages, 'x',
             label='interpolated spike')

    plt.plot((grid_spike_indices-1)*neuron.dt, voltage[grid_spike_indices-1], '.',
             label='last step before spike')

plt.xlabel('time (s)')
plt.ylabel('voltage (V)')
plt.legend(loc=3)
plt.title('Model Response')

# plot after spike currents
plt.subplot(3,1,3)
for ii in range(np.shape(after_spike_currents)[1]):
    plt.plot(time, after_spike_currents[:,ii])
plt.xlabel('time (s)')
plt.ylabel('current (A)')
plt.title('After Spike Currents')

plt.tight_layout()
plt.show()

```

Note: There both interpolated spike times and grid spike times. The grid spike is the first time step where the voltage is higher than the threshold. Note that if you try to plot the voltage at the grid spike indices the output will be NaN. The interpolated spike is the calculated intersection of the threshold and voltage between the time steps.

3.1.3 GLIF Configuration

Instances of the *GlifNeuron* class require many parameters for initialization. Fixed neuron parameters are stored directly as properties on the class instance:

Parameter	Description	Units	Type
EI	resting potential	Volts	float
dt	time duration of each simulation step	seconds	float
R_input	input resistance	Ohms	float
C	capacitance	Farads	float
asc_vector	afterspike current coefficients	Amps	np.array
spike_cut_length	spike duration	time steps	int
th_inf	instantaneous threshold	Volts	float
th_adapt	adapted threshold	Volts	float

Some of these fixed parameters were optimized to fit Allen Cell Types Database electrophysiology data. Optimized coefficients for these parameters are stored by name in the `neuron.coeffs` dictionary. For more details on which

parameters were optimized, please see the technical [white paper](#).

The `GlifNeuron` class has six methods that can be customized: three rules for updating voltage, threshold, and afterspike currents during the simulation; and three rules for updating those values when a spike is detected (voltage surpasses threshold).

Method Type	Description
<code>voltage_dynamics_method</code>	Update simulation voltage for the next time step.
<code>threshold_dynamics_method</code>	Update simulation threshold for the next time step.
<code>AScurrent_dynamics_method</code>	Update afterspike current coefficients for the next time step.
<code>voltage_reset_method</code>	Reset simulation voltage after a spike occurs.
<code>threshold_reset_method</code>	Reset simulation threshold after a spike occurs.
<code>AScurrent_reset_method</code>	Reset afterspike current coefficients after a spike occurs.

The GLIF neuron configuration files available from the Allen Brain Atlas API use built-in methods, however you can supply your own custom method if you like:

```
# define your own custom voltage reset rule
# this one linearly scales the input voltage
def custom_voltage_reset_rule(neuron, voltage_t0, custom_param_a, custom_param_b):
    return custom_param_a * voltage_t0 + custom_param_b

# initialize a neuron from a neuron config file
neuron_config = json_utilities.read('neuron_config.json')['566302806']
neuron = GlifNeuron.from_dict(neuron_config)

# configure a new method and overwrite the neuron's old method
method = neuron.configure_method('custom', custom_voltage_reset_rule,
                                { 'custom_param_a': 0.1, 'custom_param_b': 0.0 })
neuron.voltage_reset_method = method

output = neuron.run(stimulus)
```

Notice that the function is allowed to take custom parameters (here `custom_param_a` and `custom_param_b`), which are configured on method initialization from a dictionary. For more details, see the documentation for the `GlifNeuron` and `GlifNeuronMethod` classes.

3.1.4 Built-in Dynamics Rules

The job of a dynamics rule is to describe how the simulator should update the voltage, spike threshold, and afterspike currents of the simulator at a given simulation time step.

Voltage Dynamics Rules

These methods update the output voltage of the simulation. They all expect a voltage, afterspike current vector, and current injection value to be passed in by the `GlifNeuron`. All other function parameters must be fixed using the `GlifNeuronMethod` class. They all return an updated voltage value.

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_forward_euler()
```

Threshold Dynamics Rules

These methods update the spike threshold of the simulation. They all expect the current threshold and voltage values of the simulation to be passed in by the `GlifNeuron`. All other function parameters must be fixed using the `GlifNeuronMethod` class. They all return an updated threshold value.

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_three_components_exact()
```



```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_spike_component()
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_inf()
```

Afterspike Current Dynamics Rules

These methods expect current afterspike current coefficients, current time step, and time steps of all previous spikes to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated afterspike current array.

```
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_exp()
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_none()
```

3.1.5 Built-in Reset Rules

The job of a reset rule is to describe how the simulator should update the voltage, spike threshold, and afterspike currents of the simulator after the simulator has detected that the simulated voltage has surpassed threshold.

Voltage Reset Rules

These methods update the output voltage of the simulation after voltage has surpassed threshold. They all expect a voltage to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated voltage value.

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_zero()
allensdk.model.glif.glif_neuron_methods.reset_voltage_v_before()
```

Threshold Reset Rules

These methods update the spike threshold of the simulation after a spike has been detected. They all expect the current threshold and the reset voltage value of the simulation to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated threshold value.

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_inf()
allensdk.model.glif.glif_neuron_methods.reset_threshold_three_components()
```

Afterspike Reset Rules

These methods expect current afterspike current coefficients to be passed in by the GlifNeuron. All other function parameters must be fixed using the GlifNeuronMethod class. They all return an updated afterspike current array.

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_none()
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_sum()
```

3.2 Biophysical Models

The Allen Cell Types Database contains biophysical models that characterize the firing behavior of neurons measured in slices through current injection by a somatic whole-cell patch clamp electrode. These models contain a set of 10 active conductances placed at the soma and use the reconstructed 3D morphologies of the modeled neurons. The biophysical modeling [technical white paper](#) contains details on the specific construction of these models and the optimization of the model parameters to match the experimentally-recorded firing behaviors.

The biophysical models are run with the [NEURON](#) simulation environment. The Allen SDK package contains libraries that assist in downloading and setting up the models available on the Allen Institute web site for users to run using NEURON. The examples and scripts provided run on Linux using the bash shell.

3.2.1 Prerequisites

You must have NEURON with the Python interpreter enabled and the Allen SDK installed.

The Allen Institute perisomatic biophysical models were generated using NEURON [version v7.4.rel-1370](#). Instructions for compiling NEURON with the Python interpreter are available from the NEURON team under the heading [Installation with Python as an alternative interpreter](#). The Allen SDK is compatible with Python version 2.7.9, included in the Anaconda 2.1.0 distribution.

Instructions for optional [Docker installation](#) are also available.

Note: Building and installing NEURON with the Python wrapper enabled is not always easy. This page targets users that have a background in NEURON usage and installation.

3.2.2 Downloading Biophysical Models

There are two ways to download files necessary to run a biophysical model. The first way is to visit <http://celltypes.brain-map.org> and find cells that have biophysical models available for download. The electrophysiology details page for a cell has a neuronal model download link. Specifically:

1. Click ‘More Options+’
2. Check ‘Models -> Biophysical - perisomatic’ or ‘Biophysical - all active’
3. Use the Filters, Cell Location and Cell Feature Filters to narrow your results.
4. Click on a Cell Summary to view the Mouse Experiment Electrophysiology.
5. Click the “download data” link to download the NWB stimulus and response file.
6. Click “show model response” and select ‘Biophysical - perisomatic’ or ‘Biophysical - all active’.
7. Scroll down and click the ‘Biophysical - perisomatic’ or ‘Biophysical - all active’ “download model” link.

This may also be done programmatically. The neuronal model id can be found to the left of the corresponding ‘Biophysical - perisomatic’ or ‘Biophysical - all active’ “download model” link.

```
from allensdk.api.queries.biophysical_api import BiophysicalApi

bp = BiophysicalApi()
bp.cache_stimulus = True # change to False to not download the large stimulus NWB file
neuronal_model_id = 472451419 # get this from the web site as above
bp.cache_data(neuronal_model_id, working_directory='neuronal_model')
```

More help can be found in the [online help](#) for the Allen Cell Types Database web application.

3.2.3 Directory Structure

The structure of the directory created looks like this. It includes stimulus files, model parameters, morphology, cellular mechanisms and application configuration.

```
neuronal_model
|-- manifest.json
|-- 472451419_fit.json
```

(continues on next page)

(continued from previous page)

```

|-- Nr5a1-Cre_Ai14_IVSCC_-169248.04.02.01.nwb
|-- Nr5a1-Cre_Ai14_IVSCC_-169248.04.02.01_403165543_m.swc
|-- modfiles
|   |--CaDynamics.mod
|   |--Ca_HVA.mod
|   |--Ca_LVA.mod
|   |--Ih.mod
|   `--...etc.
|
|--x86_64
`---work

```

3.2.4 Running the Simulation (Linux shell prompt)

All of the sweeps available from the web site are included in manifest.json and will be run by default. This can take some time.

```

cd neuronal_model
nrnivmodl ./modfiles # compile the model (only needs to be done once)
python -m allensdk.model.biophysical.runner manifest.json # perisomatic models
python -m allensdk.model.biophysical.runner manifest.json # legacy all-active models
# new all-active models (axon replaced by a 60 micron long 1 micron diameter stub)
python -m allensdk.model.biophysical.runner manifest.json --axon_type stub

```

3.2.5 Selecting a Specific Sweep

The sweeps are listed in manifest.json. You can remove all of the sweep numbers that you do not want run.

3.2.6 Simulation Main Loop

The top level script is in the `run()` method of the `allensdk.model.biophysical.runner` module. The implementation of the method is discussed here step-by-step:

First configure NEURON based on the configuration file, which was read in from the command line at the very bottom of the script.

`run()`:

```

# configure NEURON -- this will infer model type (perisomatic vs. all-active)
utils = Utils.create_utils(description)
h = utils.h

```

The next step is to get the path of the morphology file and pass it to NEURON.

```

# configure model
manifest = description.manifest
morphology_path = description.manifest.get_path('MORPHOLOGY')
utils.generate_morphology(morphology_path.encode('ascii', 'ignore'))
utils.load_cell_parameters()

```

Then read the stimulus and recording configuration and configure NEURON

```
# configure stimulus and recording
stimulus_path = description.manifest.get_path('stimulus_path')
nwb_out_path = manifest.get_path("output")
output = NwbDataSet(nwb_out_path)
run_params = description.data['runs'][0]
sweeps = run_params['sweeps']
junction_potential = description.data['fitting'][0]['junction_potential']
mV = 1.0e-3
```

Loop through the stimulus sweeps and write the output.

```
# run sweeps
for sweep in sweeps:
    utils.setup_iclamp(stimulus_path, sweep=sweep)
    vec = utils.record_values()

    h.finitialize()
    h.run()

    # write to an NWB File
    output_data = (numpy.array(vec['v']) - junction_potential) * mV
    output.set_sweep(sweep, None, output_data)
```

3.2.7 Customization

Much of the code in the perisomatic simulation is not core Allen SDK code. The `runner.py` script largely reads the configuration file and calls into methods in the `Utils` class. `Utils` is a subclass of the `HocUtils` class, which provides access to objects in the NEURON package. The various methods called by the `runner.py` script are implemented here, including: `generate_morphology()`, `load_cell_parameters()`, `setup_iclamp()`, `read_stimulus()` and `record_values()`.

`Utils`:

```
from allensdk.model.biophys_sim.neuron.hoc_utils import HocUtils

.....

class Utils(HocUtils):
    .....

    def __init__(self, description):
        super(Utils, self).__init__(description)
    ....
```

To create a biophysical model using your own software or data, simply model your directory structure on one of the downloaded simulations or one of the examples below. Add your own `runner.py` and `utils.py` module to the simulation directory.

Compile the `.mod` files using NEURON's `nrnivmodl` command (Linux shell):

```
nrnivmodl modfiles
```

Then call your runner script directly, passing in the manifest file to your script:

```
python runner.py manifest.json
```

The output from your simulation and any intermediate files will go in the work directory.

3.2.8 Examples

A minimal example (`simple_example.tgz`) and a multicell example (`multicell_example.tgz`) are available to download as a starting point for your own projects.

Each example provides its own `utils.py` file along with a main script (Linux shell) and supporting configuration files.

`simple_example.tgz`:

```
tar xvzf simple_example.tgz
cd simple
nrnivmodl modfiles
python simple.py
```

`multicell_example.tgz`:

```
tar xvzf multicell_example.tgz
cd multicell
nrnivmodl modfiles
python multi.py
python multicell_diff.py
```

3.2.9 Exporting Output to Text Format or Image

This is an example of using the AllenSDK to save a response voltage to other formats.

```
from allensdk.core.dat_utilities import \
    DatUtilities
from allensdk.core.nwb_data_set import \
    NwbDataSet
import numpy as np
import matplotlib
matplotlib.use("Agg")
import matplotlib.pyplot as plt

nwb_file = '313862020.nwb'
sweep_number = 52
dat_file = '313862020_%d.dat' % (sweep_number)

nwb = NwbDataSet(nwb_file)
sweep = nwb.get_sweep(sweep_number)

# read v and t as numpy arrays
v = sweep['response']
dt = 1.0e3 / sweep['sampling_rate']
num_samples = len(v)
t = np.arange(num_samples) * dt

# save as text file
data = np.transpose(np.vstack((t, v)))
with open (dat_file, "w") as f:
    np.savetxt(f, data)

# save image using matplotlib
fig, ax = plt.subplots(nrows=1, ncols=1)
ax.plot(t, v)
```

(continues on next page)

(continued from previous page)

```
ax.set_title("Sweep %s" % (sweep_number))
fig.savefig('out.png')
```

3.2.10 Model Description Files

Basic Structure

A model description file is simply a JSON object with several sections at the top level and an array of JSON objects within each section.

```
{
  "cell_section": [
    {
      "name": "cell 1",
      "shape": "pyramidal",
      "position": [ 0.1, 0.2, 0.3 ]
    },
    {
      "name": "cell 2",
      "shape": "glial",
      "position": [ 0.1, 0.2, 0.3 ]
    }
  ],
  "extra": [
    { "what": "wood",
      "who": "woodchuck"
    }
  ]
}
```

Even if a section contains no objects or only one object the array brackets must be present.

Objects Within Sections

While no restrictions are enforced on what kinds of objects are stored in a section, some rules of thumb make the file easier to work with.

1. All objects within a section are the same structure. Common operations on a section are to display it as a table, iterate over it, load from or write to a spreadsheet or csv file. These operations are all easier if the section is fairly homogeneous.
2. Objects are not deeply nested. While some shallow nesting is often useful, deep nesting such as a tree structure is not recommended. It makes interoperability with other tools and data formats more difficult.
3. Arrays are allowed, though they should not be deeply nested either.
4. Object member values should be literals. Do not use pickled classes, for example.

Comment Lines

The JSON specification does not allow comments. However, the Allen SDK library applies a preprocessing stage to remove C++-style comments, so they can be used in description files.

Multi-line comments should be surrounded by `/* */` and single-line comments start with `//`. Commented description files will not be recognized by strict json parsers unless the comments are stripped.

commented.json:

```
{
  /*
   * multi-line comment
   */
  "section1": [
    {
      "name": "simon" // single line comment
    }
  ]
}
```

Split Description Files by Section

A model description can be split into multiple files by putting some sections in one file and other sections into another file. This can be useful if you want to put a topology of cells and connections in one file and experimental conditions and stimulus in another file. The resulting structure in memory will behave the same way as if the files were not split. This allows a small experiment to be described in a single file and large experiments to be more modular.

cells.json:

```
{
  "cell_section": [
    {
      "name": "cell 1",
      "shape": "pyramidal",
      "position": [ 0.1, 0.2, 0.3 ]
    },
    {
      "name": "cell 2",
      "shape": "glial",
      "position": [ 0.1, 0.2, 0.3 ]
    }
  ]
}
```

extras.json:

```
{
  "extra": [
    {
      "what": "wood",
      "who": "woodchuck"
    }
  ]
}
```

Split Sections Between Description Files

If two description files containing the same sections are combined, the resulting description will contain objects from both files. This feature allows sub-networks to be described in separate files. The sub-networks can then be composed into a larger network with an additional description of the interconnections.

network1.json:

```
/* A self-contained sub-network */
{
  "cells": [
    { "name": "cell1" },
    { "name": "cell2" }
  ],
  /* intra-network connections */
  "connections": [
    { "source": "cell1", "target" : "cell2" }
  ]
}
```

network2.json:

```
/* Another self-contained sub-network */
{
  "cells": [
    { "name": "cell3" },
    { "name": "cell4" }
  ],
  "connections": [
    { "source": "cell3", "target" : "cell4" }
  ]
}
```

interconnect.json:

```
{
  // the additional connections needed to
  // connect the network1 and network2
  // into a ring topology.
  "connections": [
    { "source": "cell2", "target": "cell3" },
    { "source": "cell4", "target": "cell1" }
  ]
}
```

3.2.11 Resource Manifest

JSON has many advantages. It is widely supported, readable and easy to parse and edit. As data sets get larger or specialized those advantages diminish. Large or complex models and experiments generally need more than a single model description file to completely describe an experiment. A manifest file is a way to describe all of the resources needed within the Allen SDK description format itself.

The manifest section is named “manifest” by default, though it is configurable. The objects in the manifest section each specify a directory, file, or file pattern. Files and directories may be organized in a parent-child relationship.

A Simple Manifest

This is a simple manifest file that specifies the BASEDIR directory using “.”, meaning the current directory:

```
{
  "manifest": [
```

(continues on next page)

(continued from previous page)

```

    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    }
  ] }
}
```

Parent Child Relationships

Adding the optional “parent_key” member to a manifest object creates a parent-child relation. In this case WORKDIR will be found in “./work”:

```

{
  "manifest": [
    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    },
    {
      "key": "WORKDIR",
      "type": "dir",
      "spec": "/work",
      "parent_key": "BASEDIR"
    }
  ] }
}
```

File Spec Patterns

Files can be specified using the type “file” instead of “dir”. If a sequence of many files is needed, the spec may contain patterns to indicate where the sequence number (%d) or string (%s) will be interpolated:

```

{
  "manifest": [
    {
      "key": "BASEDIR",
      "type": "dir",
      "spec": "."
    },
    {
      "key": "voltage_out_cell_path",
      "type": "file",
      "spec": "v_out-cell-%d.dat",
      "parent_key": "BASEDIR"
    }
  ] }
}
```

Split Manifest Files

Manifest files can be split like any description file. This allows the specification of a general directory structure in a shared file and specific files in a separate configuration (i.e. stimulus and working directory)

Extensions

To date, manifest description files have not been used to reference URLs that provide model data, but it is a planned future use case.

3.2.12 Further Reading

- [NEURON](#)
- [Python](#)
- [JSON](#)

CHAPTER 4

Examples

Take a look at the table below for a list of SDK example notebooks and scripts.

Description	Link
Introduction to the Mouse Connectivity Atlas	Jupyter notebook (download .ipynb)
Introduction to the Cell Types Database	Jupyter notebook (download .ipynb)
Introduction to the Brain Observatory	Jupyter notebook (download .ipynb)
Brain Observatory Stimulus Manipulation	Jupyter notebook (download .ipynb)
Brain Observatory Tuning Analysis	Jupyter notebook (download .ipynb)
Brain Observatory Receptive Field Analysis	Jupyter notebook (download .ipynb)
Brain Observatory Cell Specimen ID Mapping	Jupyter notebook (download .ipynb)
Brain Observatory Monitor	Jupyter notebook (download .ipynb)
Dynamic Brain Workshop 2015 experiment detail	Jupyter notebook (download .ipynb)
Stimulating a biophysical model with a square pulse	Jupyter notebook (download .ipynb)
Using a Reference Space	Jupyter notebook (download .ipynb)
Downloading Images	Jupyter notebook (download .ipynb)
Visual Coding Neuropixels Quick Start	Jupyter notebook (download .ipynb)
Visual Coding Neuropixels Reference	Jupyter notebook (download .ipynb)

CHAPTER 5

Authors

Michael Buice @mabuice
Nicholas Cain @nicain
Tom Chartrand @tmchartrand
Kael Dai @kaeldai
Saskia de Vries @saskiad
David Feng @dyf
Tim Fliss @tfliss
Marina Garrett @matchings
Richard Gerkin @rgerkin
Nathan Gouwens @gouwens
Nile Graddis @nilegraddis
Sergey Gratiy @sgratiy
@jennan @jennan
Xiaoxuan Jia @jiaxx
Wesley Jones @wesley-jones
Justin Kiggins @neuromusic
Joe Knox @jknox13
Peter Ledochowitsch @ledochowitsch
Nicholas Mei @njmei
Chris Mochizuki @mochic
Ani Nandi @anirban6908
Gabe Ocker @gocker

Michael Oliver @the-moliver

Doug Ollerenshaw @dougollerenshaw

Jed Perkins @jfperkins

Alex Piet @alexpier

Nick Ponvert @nickponvert

Kat Schelonka @kschelonka

Shu Shi @shus2018

Josh Siegle @jsiegle

Corinne Teeter @corinneteeter

Shreejoy Tripathy @stripathy

Werner Van Geit @wvangeit

Michael Wang @aimichaelwang

Isaak Willett @isaak-willett

Derric Williams @derricw

6.1 Subpackages

6.1.1 allensdk.api package

Subpackages

allensdk.api.queries package

Submodules

allensdk.api.queries.annotated_section_data_sets_api module

class allensdk.api.queries.annotated_section_data_sets_api.**AnnotatedSectionDataSetsApi** (*base*
Bases: *allensdk.api.queries.rma_api.RmaApi*)

See: [Searching Annotated SectionDataSets](#)

get_annotated_section_data_sets (*self*, *structures*, *intensity_values=None*, *density_values=None*, *pattern_values=None*, *age_names=None*)

For a list of target structures, find the SectionDataSet that matches the parameters for intensity_values, density_values, pattern_values, and Age.

Parameters

structure_graph_id [dict of integers] what to retrieve

intensity_values [array of strings, optional] 'High', 'Low', 'Medium' (default)

density_values [array of strings, optional] 'High', 'Low'

pattern_values [array of strings, optional] 'Full'

age_names [array of strings, options] for example 'E11.5', '13.5'

Returns

data [dict] The parsed JSON response message.

Notes

This method uses the non-RMA Annotated SectionDataSet endpoint.

```
get_annotated_section_data_sets_via_rma (self, structures, intensity_values=None,  
                                           density_values=None, pattern_values=None,  
                                           age_names=None)
```

For a list of target structures, find the SectionDataSet that matches the parameters for intensity_values, density_values, pattern_values, and Age.

Parameters

structure_graph_id [dict of integers] what to retrieve

intensity_values [array of strings, optional] intensity values, 'High', 'Low', 'Medium' (default)

density_values [array of strings, optional] density values, 'High', 'Low'

pattern_values [array of strings, optional] pattern values, 'Full'

age_names [array of strings, options] for example 'E11.5', '13.5'

Returns

data [dict] The parsed JSON response message.

Notes

This method uses the RMA endpoint to search annotated SectionDataSet data.

```
get_compound_annotated_section_data_sets (self, queries, fmt='json')
```

Find the SectionDataSet that matches several annotated_section_data_sets queries linked together with a Boolean 'and' or 'or'.

Parameters

queries [array of dicts] dicts with args like build_query

fmt [string, optional] 'json' or 'xml'

Returns

data [dict] The parsed JSON response message.

allensdk.api.queries.biophysical_api module

```
class allensdk.api.queries.biophysical_api.BiophysicalApi (base_uri=None)
```

Bases: *allensdk.api.queries.rma_template.RmaTemplate*

```
BIOPHYSICAL_MODEL_TYPE_IDS = (491455321, 329230710)
```

```
build_rma (self, neuronal_model_id, fmt='json')
```

Construct a query to find all files related to a neuronal model.

Parameters

neuronal_model_id [integer or string representation] key of experiment to retrieve.

fmt [string, optional] json (default) or xml

Returns

string RMA query url.

cache_data (*self*, *neuronal_model_id*, *working_directory=None*)

Take a an experiment id, query the Api RMA to get well-known-files download the files, and store them in the working directory.

Parameters

neuronal_model_id [int or string representation] found in the neuronal_model table in the api

working_directory [string] Absolute path name where the downloaded well-known files will be stored.

create_manifest (*self*, *fit_path=""*, *model_type=""*, *stimulus_filename=""*, *swc_morphology_path=""*, *marker_path=""*, *sweeps=[]*)

Generate a json configuration file with parameters for a a biophysical experiment.

Parameters

fit_path [string] filename of a json configuration file with cell parameters.

stimulus_filename [string] path to an NWB file with input currents.

swc_morphology_path [string] file in SWC format.

sweeps [array of integers] which sweeps in the stimulus file are to be used.

get_neuronal_models (*self*, *specimen_ids*, *num_rows='all'*, *count=False*, *model_type_ids=None*, ***kwargs*)

Fetch all of the biophysically detailed model records associated with a particular specimen_id

Parameters

specimen_ids [list] One or more integer ids identifying specimen records.

num_rows [int, optional] how many records to retrieve. Default is 'all'.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

model_type_ids [list, optional] One or more integer ids identifying categories of neuronal model. Defaults to all-active and perisomatic biophysical_models.

Returns

List of dict Each element is a biophysical model record, containing a unique integer id, the id of the associated specimen, and the id of the model type to which this model belongs.

get_well_known_file_ids (*self*, *neuronal_model_id*)

Query the current RMA endpoint with a neuronal_model id to get the corresponding well known file ids.

Returns

list A list of well known file id strings.

is_well_known_file_type (*self*, *wkf*, *name*)

Check if a structure has the expected name.

Parameters

wkf [dict] A well-known-file structure with nested type information.

name [string] The expected type name

See also:

`read_json` where this helper function is used.

`read_json(self, json_parsed_data)`

Get the list of well_known_file ids from a response body containing nested sample,microarray_slides,well_known_files.

Parameters

json_parsed_data [dict] Response from the Allen Institute Api RMA.

Returns

list of strings Well known file ids.

```
rma_templates = {'model_queries': [{ 'name': 'models_by_specimen', 'description': 's
```

allensdk.api.queries.brain_observatory_api module

```
class allensdk.api.queries.brain_observatory_api.BrainObservatoryApi (base_uri=None,
                                                                    dat-
                                                                    acube_uri=None)
```

Bases: `allensdk.api.queries.rma_template.RmaTemplate`

CELL_MAPPING_ID = 590985414

NWB_FILE_TYPE = 'NWBophys'

OPHYS_ANALYSIS_FILE_TYPE = 'OphysExperimentCellRoiMetricsFile'

OPHYS_EVENTS_FILE_TYPE = 'ObservatoryEventsFile'

`dataframe_query(self, data, filters, primary_key)`

Given a list of dictionary records and a list of filter dictionaries, filter the records using Pandas and return the filtered set of records.

Parameters

data: list of dicts List of dictionaries

filters: list of dicts Each dictionary describes a filtering operation on a field in the dictionary. The general form is { 'field': <field>, 'op': <operation>, 'value': <filter_value(s)> }. For example, you can apply a threshold on the "osi_dg" column with something like this: { 'field': 'osi_dg', 'op': '>', 'value': 1.0 }. See `_QUERY_TEMPLATES` for a full list of operators.

`dataframe_query_string(self, filters)`

Convert a list of cell metric filter dictionaries into a Pandas query string.

`filter_cell_specimens(self, cell_specimens, ids=None, experiment_container_ids=None, include_failed=False, filters=None)`

Filter a list of cell specimen records returned from the `get_cell_metrics` method according some of their properties.

Parameters

cell_specimens: list of dicts List of records returned by the `get_cell_metrics` method.

ids: list of integers Return only records for cells with cell specimen ids in this list

experiment_container_ids: list of integers Return only records for cells that belong to experiment container ids in this list

include_failed: bool Whether to include cells from failed experiment containers

filters: list of dicts Custom query used to reproduce filter sets created in the Allen Brain Observatory web application. The general form is a list of dictionaries each of which describes a filtering operation based on a metric. For more information, see `dataframe_query`.

filter_experiment_containers (*self*, *containers*, *ids=None*, *targeted_structures=None*, *imaging_depths=None*, *cre_lines=None*, *reporter_lines=None*, *transgenic_lines=None*, *include_failed=False*, *simple=False*)

filter_experiments_and_containers (*self*, *objs*, *ids=None*, *targeted_structures=None*, *imaging_depths=None*, *cre_lines=None*, *reporter_lines=None*, *transgenic_lines=None*, *include_failed=False*)

filter_ophys_experiments (*self*, *experiments*, *ids=None*, *experiment_container_ids=None*, *targeted_structures=None*, *imaging_depths=None*, *cre_lines=None*, *reporter_lines=None*, *transgenic_lines=None*, *stimuli=None*, *session_types=None*, *include_failed=False*, *require_eye_tracking=False*, *simple=False*)

get_cell_metrics (*self*, *cell_specimen_ids=None*, **args*, ***kwargs*)
Get cell metrics by id

Parameters

cell_metrics_ids [integer or list of integers, optional] only select specific cell metric records.

Returns

dict [cell metric metadata]

get_cell_specimen_id_mapping (*self*, *file_name*, *mapping_table_id=None*)
Download mapping table from old to new cell specimen IDs.

The mapping table is a CSV file that maps cell specimen ids that have changed between processing runs of the Brain Observatory pipeline.

Parameters

file_name [string] Filename to save locally.

mapping_table_id [integer] ID of the mapping table file. Defaults to the most recent mapping table.

Returns

pandas.DataFrame Mapping table as a DataFrame.

get_column_definitions (*self*, *api_class_name=None*)
Get column definitions

Parameters

api_class_names [string or list of strings, optional] only select specific column definition records.

Returns

dict [column definition metadata]

get_experiment_container_metrics (*self*, *experiment_container_metric_ids=None*)
Get experiment container metrics by id

Parameters

isi_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [isi experiment metadata]

get_experiment_containers (*self*, *experiment_container_ids=None*)

Get experiment container by id

Parameters

experiment_container_ids [integer or list of integers, optional] only select specific experiment containers.

Returns

dict [experiment container metadata]

get_isi_experiments (*self*, *isi_experiment_ids=None*)

Get ISI Experiments by id

Parameters

isi_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [isi experiment metadata]

get_ophys_experiments (*self*, *ophys_experiment_ids=None*)

Get OPhys Experiments by id

Parameters

ophys_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [ophys experiment metadata]

get_stimulus_mappings (*self*, *stimulus_mapping_ids=None*)

Get stimulus mappings by id

Parameters

stimulus_mapping_ids [integer or list of integers, optional] only select specific stimulus mapping records.

Returns

dict [stimulus mapping metadata]

list_column_definition_class_names (*self*)

Get column definitions

Returns

list [api class name strings]

list_isi_experiments (*self*, *isi_ids=None*)

List ISI experiments available through the Allen Institute API

Parameters

neuronal_model_ids [integer or list of integers, optional] only select specific isi experiments.

Returns

```

    dict [neuronal model metadata]

    rma_templates = {'brain_observatory_queries': [{'name': 'list_isi_experiments', 'des
save_ophys_experiment_analysis_data(self, ophys_experiment_id, file_name)
save_ophys_experiment_data(self, ophys_experiment_id, file_name)
save_ophys_experiment_event_data(self, ophys_experiment_id, file_name)
save_ophys_experiment_eye_gaze_data(self, ophys_experiment_id: int, ophys_session_id:
                                int, file_name: str)
simplify_experiment_containers(self, containers)
simplify_ophys_experiments(self, exps)

allensdk.api.queries.brain_observatory_api.find_container_tags(container)
    Custom logic for extracting tags from donor conditions. Filtering out tissueocyte tags.

allensdk.api.queries.brain_observatory_api.find_experiment_acquisition_age(exp)
allensdk.api.queries.brain_observatory_api.find_specimen_cre_line(specimen)
allensdk.api.queries.brain_observatory_api.find_specimen_reporter_line(specimen)
allensdk.api.queries.brain_observatory_api.find_specimen_transgenic_lines(specimen)

```

allensdk.api.queries.cell_types_api module

```

class allensdk.api.queries.cell_types_api.CellTypesApi (base_uri=None)
    Bases: allensdk.api.queries.rma_api.RmaApi

    HUMAN = 'Homo Sapiens'
    MARKER_FILE_TYPE = '3DNeuronMarker'
    MOUSE = 'Mus musculus'
    NWB_FILE_TYPE = 'NWBDownload'
    SWC_FILE_TYPE = '3DNeuronReconstruction'

    filter_cells(self, cells, require_morphology, require_reconstruction, reporter_status, species)
        Filter a list of cell specimens to those that optionally have morphologies or have morphological recon-
        structions.

```

Parameters

```

    cells: list List of cell metadata dictionaries to be filtered

    require_morphology: boolean Filter out cells that have no morphological images.

    require_reconstruction: boolean Filter out cells that have no morphological recon-
        structions.

    reporter_status: list Filter for cells that have a particular cell reporter status

    species: list Filter for cells that belong to one or more species. If None, return all. Must be
        one of [ CellTypesApi.MOUSE, CellTypesApi.HUMAN ].

    filter_cells_api(self, cells, require_morphology=False, require_reconstruction=False, re-
        porter_status=None, species=None, simple=True)

```

get_cell (*self*, *id*)

Query the API for a one cells in the Cell Types Database.

Returns

list Meta data for one cell.

get_ephys_features (*self*)

Query the API for the full table of EphysFeatures for all cells.

get_ephys_sweeps (*self*, *specimen_id*)

Query the API for a list of sweeps for a particular cell in the Cell Types Database.

Parameters

specimen_id: int Specimen ID of a cell.

Returns

list: List of sweep dictionaries belonging to a cell

get_morphology_features (*self*)

Query the API for the full table of morphology features for all cells

Notes

by default the tags column is removed because it isn't useful

list_cells (*self*, *id=None*, *require_morphology=False*, *require_reconstruction=False*, *reporter_status=None*, *species=None*)

Query the API for a list of all cells in the Cell Types Database.

Parameters

id: int ID of a cell. If not provided returns all matching cells.

require_morphology: boolean Only return cells that have morphology images.

require_reconstruction: boolean Only return cells that have morphological reconstructions.

reporter_status: list Return cells that have a particular cell reporter status.

species: list Filter for cells that belong to one or more species. If None, return all. Must be one of [CellTypesApi.MOUSE, CellTypesApi.HUMAN].

Returns

list Meta data for all cells.

list_cells_api (*self*, *id=None*, *require_morphology=False*, *require_reconstruction=False*, *reporter_status=None*, *species=None*)

save_ephys_data (*self*, *specimen_id*, *file_name*)

Save the electrophysiology recordings for a cell as an NWB file.

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the NWB file.

save_reconstruction (*self*, *specimen_id*, *file_name*)

Save the morphological reconstruction of a cell as an SWC file.

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the SWC file.

save_reconstruction_markers (*self, specimen_id, file_name*)

Save the marker file for the morphological reconstruction of a cell. These are comma-delimited files indicating points of interest in a reconstruction (truncation points, early tracing termination, etc).

Parameters

specimen_id: int ID of the specimen, from the Specimens database model in the Allen Institute API.

file_name: str Path to save the marker file.

simplify_cells_api (*self, cells*)

allensdk.api.queries.connected_services module

class allensdk.api.queries.connected_services.**ConnectedServices**

Bases: object

A class representing a schema of informatics web services.

Notes

See [Connected Services and Pipes](#) for a human-readable list of services and parameters.

The URL format is documented at [Service Pipelines](#).

Connected Services only include API services that are accessed via the RMA endpoint using an rma::services stage.

ARRAY = 'array'

BOOLEAN = 'boolean'

FLOAT = 'float'

INTEGER = 'integer'

STRING = 'string'

build_url (*self, service_name, kwargs*)

Create a single stage RMA url from a service name and parameters.

classmethod **schema** ()

Dictionary of service names and parameters.

Notes

See [Connected Services and Pipes](#) for a human-readable list of connected services and their parameters.

allensdk.api.queries.glif_api module

class allensdk.api.queries.glif_api.GlifApi (*base_uri=None*)

Bases: *allensdk.api.queries.rma_template.RmaTemplate*

GLIF_TYPES = [395310498, 395310469, 395310475, 395310479, 471355161]

NWB_FILE_TYPE = None

cache_stimulus_file (*self, output_file_name*)

DEPRECATED Download the NWB file for the current neuronal model and save it to a file.

Parameters

output_file_name: string File name to store the NWB file.

get_ephys_sweeps (*self*)

DEPRECATED Retrieve ephys sweep information out of downloaded metadata for a neuronal model

Returns

list A list of sweeps metadata dictionaries

get_neuron_config (*self, output_file_name=None*)

DEPRECATED Retrieve a model configuration file from the API, optionally save it to disk, and return the contents of that file as a dictionary.

Parameters

output_file_name: string File name to store the neuron configuration (optional).

get_neuron_configs (*self, neuronal_model_ids=None*)

get_neuronal_model (*self, neuronal_model_id*)

DEPRECATED Query the current RMA endpoint with a neuronal_model id to get the corresponding well known files and meta data.

Returns

dict A dictionary containing

get_neuronal_model_templates (*self*)

get_neuronal_models (*self, ephys_experiment_ids=None*)

get_neuronal_models_by_id (*self, neuronal_model_ids=None*)

list_neuronal_models (*self*)

DEPRECATED Query the API for a list of all GLIF neuronal models.

Returns

list Meta data for all GLIF neuronal models.

rma_templates = {'glif_queries': [{'name': 'neuronal_model_templates', 'description':

allensdk.api.queries.grid_data_api module

class allensdk.api.queries.grid_data_api.GridDataApi (*resolution=None,*

base_uri=None)

Bases: *allensdk.api.queries.rma_api.RmaApi*

HTTP Client for the Allen 3-D Expression Grid Data Service.

See: [Downloading 3-D Expression Grid Data](#)


```
DATA_MASK = 'data_mask'
```

```
DENSITY = 'density'
```

```
ENERGY = 'energy'
```

```
INJECTION_DENSITY = 'injection_density'
```

```
INJECTION_ENERGY = 'injection_energy'
```

```
INJECTION_FRACTION = 'injection_fraction'
```

```
INTENSITY = 'intensity'
```

```
PROJECTION_DENSITY = 'projection_density'
```

```
PROJECTION_ENERGY = 'projection_energy'
```

```
download_alignment3d(self, section_data_set_id, num_rows='all', count=False, **kwargs)
```

Download the parameters of the 3D affine transformation mapping this section data set's image-space stack to CCF-space (or vice-versa).

Parameters

section_data_set_id [int] download the parameters for this data set.

Returns

dict : parameters of this section data set's alignment3d

```
download_deformation_field(self, section_data_set_id, header_path=None,
                             voxel_path=None, voxel_type='DeformationFieldVoxels',
                             header_type='DeformationFieldHeader')
```

Download the local alignment parameters for this dataset. This a 3D vector image (3 components) describing a deformable local mapping from CCF voxels to this section data set's affine-aligned image stack.

Parameters

section_data_set_id [int]

Download the deformation field for this data set

header_path [str, optional] If supplied, the deformation field header will be downloaded to this path.

voxel_path [str, optional] If supplied, the deformation field voxels will be downloaded to this path.

voxel_type [str] WellKnownFileType of this dataset's data file

header_type [str] WellKnownFileType of this dataset's header file

```
download_expression_grid_data(self, section_data_set_id, include=None, path=None)
```

Download in zipped metaimage format.

Parameters

section_data_set_id [integer] What to download.

include [list of strings, optional] Image volumes. 'energy' (default), 'density', 'intensity'.

path [string, optional] File name to save as.

Returns

file [3-D expression grid data packaged into a compressed archive file (.zip).]

download_gene_expression_grid_data (*self, section_data_set_id, volume_type, path*)

Download a metainage file containing registered gene expression grid data

Parameters

section_data_set_id [int] Download data from this experiment

volume_type [str] Download this type of data (options are GridDataApi.ENERGY, GridDataApi.DENSITY, GridDataApi.INTENSITY)

path [str] Download to this path

download_projection_grid_data (*self, section_data_set_id, image=None, resolution=None, save_file_path=None*)

Download in NRRD format.

Parameters

section_data_set_id [integer] What to download.

image [list of strings, optional] Image volume. 'projection_density', 'projection_energy', 'injection_fraction', 'injection_density', 'injection_energy', 'data_mask'.

resolution [integer, optional] in microns. 10, 25, 50, or 100 (default).

save_file_path [string, optional] File name to save as.

Notes

See [Downloading 3-D Projection Grid Data](#) for additional documentation.

allensdk.api.queries.image_download_api module

class allensdk.api.queries.image_download_api.**ImageDownloadApi** (*base_uri=None*)

Bases: *allensdk.api.queries.rma_template.RmaTemplate*

HTTP Client to download whole or partial two-dimensional images from the Allen Institute with the Section-Image, AtlasImage and ProjectionImage Download Services.

See [Downloading an Image](#) for more documentation.

COLORMAPS = {'aba': 8, 'aibsmmap_alt': 9, 'blue': 6, 'colormap': 10, 'expression':

atlas_image_query (*self, atlas_id, image_type_name=None*)

List atlas images belonging to a specified atlas

Parameters

atlas_id [integer, optional] Find images from this atlas.

image_type_name [string, optional] Restrict response to images of this type. If not provided, the query will get it from the atlas id.

Returns

list of dict : Each element is an AtlasImage record.

Notes

See [Downloading Atlas Images and Graphics](#) for additional documentation. `allensdk.api.queries.ontologies_api.OntologiesApi.get_atlases()` can also be used to list atlases along with their ids.

download_atlas_image (*self*, *atlas_image_id*, *file_path=None*, ***kwargs*)

download_image (*self*, *image_id*, *file_path=None*, *endpoint=None*, ***kwargs*)

Download whole or partial two-dimensional images from the Allen Institute with the SectionImage or AtlasImage service.

Parameters

image_id [integer] SubImage to download.

file_path [string, optional] where to put it, defaults to image_id.jpg

downsample [int, optional] Number of times to downsample the original image.

quality [int, optional] jpeg quality of the returned image, 0 to 100 (default)

expression [boolean, optional] Request the expression mask for the SectionImage.

view [string, optional] 'expression', 'projection', 'tumor_feature_annotation' or 'tumor_feature_boundary'

top [int, optional] Index of the topmost row of the region of interest.

left :int, optional Index of the leftmost column of the region of interest.

width [int, optional] Number of columns in the output image.

height [int, optional] Number of rows in the output image.

range [list of ints, optional] Filter to specify the RGB channels. low,high,low,high,low,high

colormap [list of floats, optional] Filter to specify the RGB channels. [lower_threshold,colormap] gain 0-1, colormap id is a string from ImageDownload-Api.COLORMAPS

rgb [list of floats, optional] Filter to specify the RGB channels. [red,green,blue] 0-1

contrast [list of floats, optional] Filter to specify contrast parameters. [gain,bias] 0-1

annotation [boolean, optional] Request the annotated AtlasImage

atlas [int, optional] Specify the desired Atlas' annotations.

projection [boolean, optional] Request projection for the specified image.

downsample_dimensions [boolean, optional] Indicates if the width and height should be adjusted to account for downsampling.

Returns

None the file is downloaded and saved to the path.

Notes

By default, an unfiltered full-sized image with the highest quality is returned as a download if no parameters are provided.

'downsample=1' halves the number of pixels of the original image both horizontally and vertically.
range_list = kwargs.get('range', None)

Specifying ‘downsample=2’ quarters the height and width values.

Quality must be an integer from 0, for the lowest quality, up to as high as 100. If it is not specified, it defaults to the highest quality.

Top is specified in full-resolution (largest tier) pixel coordinates. `SectionImage.y` is the default value.

Left is specified in full-resolution (largest tier) pixel coordinates. `SectionImage.x` is the default value.

Width is specified in tier-resolution (desired tier) pixel coordinates. `SectionImage.width` is the default value. It is automatically adjusted when downsampled.

Height is specified in tier-resolution (desired tier) pixel coordinates. `SectionImage.height` is the default value. It is automatically adjusted when downsampled.

The range parameter consists of 6 comma delimited integers that define the lower (0) and upper (4095) bound for each channel in red-green-blue order (i.e. “range=0,1500,0,1000,0,4095”). The default range values can be determined by referring to the following fields on the Equalization model associated with the `SectionDataSet`: `red_lower`, `red_upper`, `green_lower`, `green_upper`, `blue_lower`, `blue_upper`. For more information, see the [Image Controls](#) section of the Allen Mouse Brain Connectivity Atlas: [Projection Dataset](#) help topic. See: ‘Image Download Service’ <<http://help.brain-map.org/display/api/Downloading+an+Image>>_

download_projection_image (*self*, *projection_image_id*, *file_path=None*, ***kwargs*)

download_section_image (*self*, *section_image_id*, *file_path=None*, ***kwargs*)

get_section_data_sets_by_product (*self*, *product_ids*, *include_failed=False*,
num_rows='all', *count=False*, ***kwargs*)

List all of the section data sets produced as part of one or more products

Parameters

product_ids [list of int] Integer specifiers for Allen Institute products. A product is a set of related data.

include_failed [bool, optional] If True, find both failed and passed datasets. Default is False

num_rows [int, optional] how many records to retrieve. Default is ‘all’.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

Returns

list of dict : Each returned element is a section data set record.

Notes

See <http://api.brain-map.org/api/v2/data/query.json?criteria=model::Product> for a list of products.

get_section_image_ranges (*self*, *section_image_ids*, *num_rows='all'*, *count=False*,
as_lists=True, ***kwargs*)

Section images from the Mouse Connectivity Atlas are displayed on connectivity.brain-map.org after having been linearly windowed and leveled. This method obtains parameters defining channelwise upper and lower bounds of the windows used for one or more images.

Parameters

section_image_ids [list of int] Each element is a unique identifier for a section image.

num_rows [int, optional] how many records to retrieve. Default is ‘all’.

count [bool, optional] If True, return a count of the lines found by the query. Default is False.

as_lists [bool, optional] If True, return the window parameters in a list, rather than a dict (this is the format of the range parameter on ImageDownloadApi.download_image). Default is False.

Returns

list of dict or list of list : For each section image id provided, return the window bounds for each channel.

```
rma_templates = {'image_queries': [{'name': 'section_image_ranges', 'description':
```

```
section_image_query (self, section_data_set_id, num_rows='all', count=False, **kwargs)
```

List section images belonging to a specified section data set

Parameters

atlas_id [integer, optional] Find images from this section data set.

num_rows [int] how many records to retrieve. Default is 'all'

count [bool] If True, return a count of the lines found by the query.

Returns

list of dict : Each element is an SectionImage record.

Notes

The SectionDataSet model is used to represent single experiments which produce an array of images. This includes Mouse Connectivity and Mouse Brain Atlas experiments, among other projects. You may see references to the ids of experiments from those projects. These are the same as section data set ids.

allensdk.api.queries.mouse_atlas_api module

```
class allensdk.api.queries.mouse_atlas_api.MouseAtlasApi (base_uri=None)
```

Bases: *allensdk.api.queries.reference_space_api.ReferenceSpaceApi*, *allensdk.api.queries.grid_data_api.GridDataApi*

Downloads Mouse Brain Atlas grid data, reference volumes, and metadata.

```
DEVMOUSE_ATLAS_PRODUCTS = (3,)
```

```
HUMAN_ORGANISM = (1,)
```

```
MOUSE_ATLAS_PRODUCTS = (1,)
```

```
MOUSE_ORGANISM = (2,)
```

```
download_expression_density (self, path, experiment_id)
```

```
download_expression_energy (self, path, experiment_id)
```

```
download_expression_intensity (self, path, experiment_id)
```

```
get_genes (self, organism_ids=None, chromosome_ids=None, **kwargs)
```

Download a list of genes

Parameters

organism_ids [list of int, optional] Filter genes to those appearing in these organisms. Defaults to mouse (2).

chromosome_ids [list of int, optional] Filter genes to those appearing on these chromosomes. Defaults to all.

Returns

list of dict: Each element is a gene record, with a nested chromosome record (also a dict).

get_section_data_sets (*self*, *gene_ids=None*, *product_ids=None*, ***kwargs*)

Download a list of section data sets (experiments) from the Mouse Brain Atlas project.

Parameters

gene_ids [list of int, optional] Filter results based on the genes whose expression was characterized in each experiment. Default is all.

product_ids [list of int, optional] Filter results to a subset of products. Default is the Mouse Brain Atlas.

Returns

list of dict : Each element is a section data set record, with one or more gene records nested in a list.

allensdk.api.queries.mouse_connectivity_api module

class allensdk.api.queries.mouse_connectivity_api.**MouseConnectivityApi** (*base_uri=None*)

Bases: [allensdk.api.queries.reference_space_api.ReferenceSpaceApi](#), [allensdk.api.queries.grid_data_api.GridDataApi](#)

HTTP Client for the Allen Mouse Brain Connectivity Atlas.

See: [Mouse Connectivity API](#)

PRODUCT_IDS = [5, 31]

build_reference_aligned_image_channel_volumes_url (*self*, *data_set_id*)

Construct url to download the red, green, and blue channels aligned to the 25um adult mouse brain reference space volume.

Parameters

data_set_id [integerallensdk.api.queries] aka attachable_id

Notes

See: [Reference-aligned Image Channel Volumes](#) for additional documentation.

calculate_injection_centroid (*self*, *injection_density*, *injection_fraction*, *resolution=25*)

Compute the centroid of an injection site.

Parameters

injection_density: **np.ndarray** The injection density volume of an experiment

injection_fraction: **np.ndarray** The injection fraction volume of an experiment

download_data_mask (*self*, *path*, *experiment_id*, *resolution*)

download_injection_density (*self*, *path*, *experiment_id*, *resolution*)

download_injection_fraction (*self*, *path*, *experiment_id*, *resolution*)

download_projection_density (*self*, *path*, *experiment_id*, *resolution*)

download_reference_aligned_image_channel_volumes (*self*, *data_set_id*,
save_file_path=None)

Returns

The well known file is downloaded

experiment_correlation_search (*self*, ***kwargs*)

Select a seed experiment and a domain over which the similarity comparison is to be made.

Parameters

row [integer] SectionDataSet.id to correlate against.

structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

hemisphere [string, optional] Use 'right' or 'left'. Defaults to both hemispheres.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Correlation Search](#) and [service::mouse_connectivity_correlation](#).

experiment_injection_coordinate_search (*self*, ***kwargs*)

User specifies a seed location within the 3D reference space. The service returns a rank list of experiments by distance of its injection site to the specified seed location.

Parameters

seed_point [list of floats] The coordinates of a point in 3-D SectionDataSet space.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Injection Coordinate Search](#) and `service::mouse_connectivity_injection_coordinate`.

experiment_source_search (*self*, ***kwargs*)

Search over the whole projection signal statistics dataset to find experiments with specific projection profiles.

Parameters

injection_structures [list of integers or strings] Integer Structure.id or String Structure.acronym.

target_domain [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

injection_hemisphere [string, optional] 'right' or 'left', Defaults to both hemispheres.

target_hemisphere [string, optional] 'right' or 'left', Defaults to both hemispheres.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

injection_domain [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Source Search](#), [Target Search](#), and `service::mouse_connectivity_injection_structure`.

experiment_spatial_search (*self*, ***kwargs*)

Displays all SectionDataSets with projection signal density ≥ 0.1 at the seed point. This service also returns the path along the most dense pixels from the seed point to the center of each injection site..

Parameters

seed_point [list of floats] The coordinates of a point in 3-D SectionDataSet space.

transgenic_lines [list of integers or strings, optional] Integer TransgenicLine.id or String TransgenicLine.name. Specify ID 0 to exclude all TransgenicLines.

section_data_sets [list of integers, optional] Ids to filter the results.

injection_structures [list of integers or strings, optional] Integer Structure.id or String Structure.acronym.

primary_structure_only [boolean, optional]

product_ids [list of integers, optional] Integer Product.id

start_row [integer, optional] For paging purposes. Defaults to 0.

num_rows [integer, optional] For paging purposes. Defaults to 2000.

Notes

See [Spatial Search](#) and `service::mouse_connectivity_target_spatial`.

get_experiment_detail (*self*, *experiment_id*)

Retrieve the experiments data.

get_experiments (*self*, *structure_ids*, ***kwargs*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas.

Parameters

structure_ids [integer or list, optional] injection structure

Returns

url [string] The constructed URL

get_experiments_api (*self*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas via the ApiConnectivity table.

Returns

url [string] The constructed URL

get_manual_injection_summary (*self*, *experiment_id*)

Retrieve manual injection summary.

get_projection_image_info (*self*, *experiment_id*, *section_number*)

Fetch meta-information of one projection image.

Parameters

experiment_id [integer]

section_number [integer]

Notes

See: [image examples](#) under [Experimental Overview and Metadata](#) for additional documentation. Download the image using `allensdk.api.queries.image_download_api.ImageDownloadApi.download_section_image()`

get_reference_aligned_image_channel_volumes_url (*self*, *data_set_id*)

Retrieve the download link for a specific data set. Notes — See [Reference-aligned Image Channel Volumes](#) for additional documentation.

get_structure_unionizes (*self*, *experiment_ids*, *is_injection=None*, *structure_name=None*, *structure_ids=None*, *hemisphere_ids=None*, *normalized_projection_volume_limit=None*, *include=None*, *debug=None*, *order=None*)

allensdk.api.queries.ontologies_api module

class `allensdk.api.queries.ontologies_api.OntologiesApi` (*base_uri=None*)

Bases: `allensdk.api.queries.rma_template.RmaTemplate`

See: [Atlas Drawings and Ontologies](#)

get_atlases (*self*)

get_atlases_table (*self*, *atlas_ids=None*, *brief=True*)

List Atlases available through the API with associated ontologies and structure graphs.

Parameters

atlas_ids [integer or list of integers, optional] only select specific atlases

brief [boolean, optional] True (default) requests only name and id fields.

Returns

dict [atlas metadata]

Notes

This query is based on the [table of available Atlases](#). See also: [Class: Atlas](#)

get_structure_graphs (*self*)

get_structure_sets (*self*, *structure_set_ids=None*)

get_structures (*self*, *structure_graph_ids=None*, *structure_graph_names=None*,
structure_set_ids=None, *structure_set_names=None*, *order=*
der=['structures.graph_order'], *num_rows='all'*, *count=False*, ***kwargs*)

Retrieve data about anatomical structures.

Parameters

structure_graph_ids [int or list of ints, optional] database keys to get all structures in particular graphs

structure_graph_names [string or list of strings, optional] list of graph names to narrow the query

structure_set_ids [int or list of ints, optional] database keys to get all structures in a particular set

structure_set_names [string or list of strings, optional] list of set names to narrow the query.

order [list of strings] list of RMA order clauses for sorting

num_rows [int] how many records to retrieve

Returns

dict the parsed json response containing data from the API

Notes

Only one of the methods of limiting the query should be used at a time.

get_structures_with_sets (*self*, *structure_graph_ids*, *order=*
der=['structures.graph_order'],
num_rows='all', *count=False*, ***kwargs*)

Download structures along with the sets to which they belong.

Parameters

structure_graph_ids [int or list of int] Only fetch structure records from these graphs.

order [list of strings] list of RMA order clauses for sorting

num_rows [int] how many records to retrieve

Returns

dict the parsed json response containing data from the API

rma_templates = {'ontology_queries': [{'name': 'structures_by_graph_ids', 'descripti

unpack_structure_set_ancestors (*self*, *structure_dataframe*)

Convert a slash-separated *structure_id_path* field to a list.

Parameters

structure_dataframe [DataFrame] structure data from the API

Returns

None A new column is added to the dataframe containing the ancestor list.

allensdk.api.queries.reference_space_api module

class allensdk.api.queries.reference_space_api.**ReferenceSpaceApi** (*base_uri=None*)

Bases: *allensdk.api.queries.rma_api.RmaApi*

ARA_NISSL = 'ara_nissl'

AVERAGE_TEMPLATE = 'average_template'

CCF_2015 = 'annotation/ccf_2015'

CCF_2016 = 'annotation/ccf_2016'

CCF_2017 = 'annotation/ccf_2017'

CCF_VERSION_DEFAULT = 'annotation/ccf_2017'

DEVMOUSE_2012 = 'annotation/devmouse_2012'

MOUSE_2011 = 'annotation/mouse_2011'

VOXEL_RESOLUTION_100_MICRONS = 100

VOXEL_RESOLUTION_10_MICRONS = 10

VOXEL_RESOLUTION_25_MICRONS = 25

VOXEL_RESOLUTION_50_MICRONS = 50

build_volumetric_data_download_url (*self*, *data_path*, *file_name*, *voxel_resolution=None*,
release=None, *coordinate_framework=None*)

Construct url to download 3D reference model in NRRD format.

Parameters

data_path [string] 'average_template', 'ara_nissl', 'annotation/ccf_{year}', 'annotation/mouse_2011', or 'annotation/devmouse_2012'

voxel_resolution [int] 10, 25, 50 or 100

coordinate_framework [string] 'mouse_ccf' (default) or 'mouse_annotation'

Notes

See: [3-D Reference Models](#) for additional documentation.

download_annotation_volume (*self*, *ccf_version*, *resolution*, *file_name*)

Download the annotation volume at a particular resolution.

Parameters

ccf_version: string Which reference space version to download. Defaults to “annotation/ccf_2017”

resolution: int Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name: string Where to save the annotation volume.

Note: the parameters must be used as positional parameters, not keywords

download_mouse_atlas_volume (*self, age, volume_type, file_name*)

Download a reference volume (annotation, grid annotation, atlas volume) from the mouse brain atlas project

Parameters

age [str] Specify a mouse age for which to download the reference volume

volume_type [str] Specify the type of volume to download

file_name [str] Specify the path to the downloaded volume

download_structure_mask (*self, structure_id, ccf_version, resolution, file_name*)

Download an indicator mask for a specific structure.

Parameters

structure_id [int] Unique identifier for the annotated structure

ccf_version [string] Which reference space version to download. Defaults to “annotation/ccf_2017”

resolution [int] Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name [string] Where to save the downloaded mask.

download_structure_mesh (*self, structure_id, ccf_version, file_name*)

Download a Wavefront obj file containing a triangulated 3d mesh built from an annotated structure.

Parameters

structure_id [int] Unique identifier for the annotated structure

ccf_version [string] Which reference space version to download. Defaults to “annotation/ccf_2017”

file_name [string] Where to save the downloaded mask.

download_template_volume (*self, resolution, file_name*)

Download the registration template volume at a particular resolution.

Parameters

resolution: int Desired resolution to download in microns. Must be 10, 25, 50, or 100.

file_name: string Where to save the registration template volume.

download_volumetric_data (*self, data_path, file_name, voxel_resolution=None, save_file_path=None, release=None, coordinate_framework=None*)

Download 3D reference model in NRRD format.

Parameters

data_path [string] ‘average_template’, ‘ara_nissl’, ‘annotation/ccf_{year}’, ‘annotation/mouse_2011’, or ‘annotation/devmouse_2012’

file_name [string] server-side file name. ‘annotation_10.nrrd’ for example.

voxel_resolution [int] 10, 25, 50 or 100

coordinate_framework [string] 'mouse_ccf' (default) or 'mouse_annotation'

Notes

See: [3-D Reference Models](#) for additional documentation.

allensdk.api.queries.rma_api module

class allensdk.api.queries.rma_api.**RmaApi** (*base_uri=None*)

Bases: [allensdk.api.api.Api](#)

See: [RESTful Model Access \(RMA\)](#)

ALL = 'all'

COUNT = 'count'

CRITERIA = 'rma::criteria'

DEBUG = 'debug'

EQ = '\$eq'

EXCEPT = 'except'

EXCPT = 'excpt'

FALSE = 'false'

INCLUDE = 'rma::include'

IS = '\$is'

MODEL = 'model::'

NUM_ROWS = 'num_rows'

ONLY = 'only'

OPTIONS = 'rma::options'

ORDER = 'order'

PIPE = 'pipe::'

PREVIEW = 'preview'

SERVICE = 'service::'

START_ROW = 'start_row'

TABULAR = 'tabular'

TRUE = 'true'

build_query_url (*self, stage_clauses, fmt='json'*)

Combine one or more RMA query stages into a single RMA query.

Parameters

stage_clauses [list of strings] subqueries

fmt [string, optional] json (default), xml, or csv

Returns

string complete RMA url

build_schema_query (*self*, *clazz=None*, *fmt='json'*)

Build the URL that will fetch the data schema.

Parameters

clazz [string, optional] Name of a specific class or None (default).

fmt [string, optional] json (default) or xml

Returns

url [string] The constructed URL

Notes

If a class is specified, only the schema information for that class will be requested, otherwise the url requests the entire schema.

debug_clause (*self*, *debug_value=None*)

Construct a debug clause for use in an rma::options clause. Parameters ——— debug_value : string or boolean

True, False, None (default) or 'preview'

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

True will request debugging information in the response. False will request no debugging information. None will return an empty clause. 'preview' will request debugging information without the query being run.

filter (*self*, *key*, *value*)

serialize a single RMA query filter clause.

Parameters

key [string] keys for narrowing a query.

value [string] value for narrowing a query.

Returns

string a single filter clause for an RMA query string.

filters (*self*, *filters*)

serialize RMA query filter clauses.

Parameters

filters [dict] keys and values for narrowing a query.

Returns

string filter clause for an RMA query string.

get_schema (*self*, *clazz=None*)

Retrieve schema information.

model_query (*self*, *args, **kwargs)

Construct and execute a model stage of an RMA query string.

Parameters

- model** [string] The top level data type
- filters** [dict] key, value comparisons applied to the top-level model to narrow the results.
- criteria** [string] raw RMA criteria clause to choose what object are returned
- include** [string] raw RMA include clause to return associated objects
- only** [list of strings, optional] to be joined into an rma::options only filter to limit what data is returned
- except** [list of strings, optional] to be joined into an rma::options except filter to limit what data is returned
- except** [list of strings, optional] synonym for except parameter to avoid a reserved word conflict.
- tabular** [list of string, optional] return columns as a tabular data structure rather than a nested tree.
- count** [boolean, optional] False to skip the extra database count query.
- debug** [string, optional] 'true', 'false' or 'preview'
- num_rows** [int or string, optional] how many database rows are returned (may not correspond directly to JSON tree structure)
- start_row** [int or string, optional] which database row is start of returned data (may not correspond directly to JSON tree structure)

Notes

See [RMA Path Syntax](#) for a brief overview of the normalized RMA syntax. Normalized RMA syntax differs from the legacy syntax used in much of the RMA documentation. Using the &debug=true option with an RMA URL will include debugging information in the response, including the normalized query.

model_stage (*self*, *model*, **kwargs)

Construct a model stage of an RMA query string.

Parameters

- model** [string] The top level data type
- filters** [dict] key, value comparisons applied to the top-level model to narrow the results.
- criteria** [string] raw RMA criteria clause to choose what object are returned
- include** [string] raw RMA include clause to return associated objects
- only** [list of strings, optional] to be joined into an rma::options only filter to limit what data is returned
- except** [list of strings, optional] to be joined into an rma::options except filter to limit what data is returned
- tabular** [list of string, optional] return columns as a tabular data structure rather than a nested tree.
- count** [boolean, optional] False to skip the extra database count query.

debug [string, optional] ‘true’, ‘false’ or ‘preview’

num_rows [int or string, optional] how many database rows are returned (may not correspond directly to JSON tree structure)

start_row [int or string, optional] which database row is start of returned data (may not correspond directly to JSON tree structure)

Notes

See [RMA Path Syntax](#) for a brief overview of the normalized RMA syntax. Normalized RMA syntax differs from the legacy syntax used in much of the RMA documentation. Using the `&debug=true` option with an RMA URL will include debugging information in the response, including the normalized query.

only_except_tabular_clause (*self*, *filter_type*, *attribute_list*)

Construct a clause to filter which attributes are returned for use in an `rma::options` clause.

Parameters

filter_type [string] ‘only’, ‘except’, or ‘tabular’

attribute_list [list of strings] for example [‘acronym’, ‘products.name’, ‘structure.id’]

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

The title of tabular columns can be set by adding ‘+as+<title>’ to the attribute. The tabular filter type requests a response that is row-oriented rather than a nested structure. Because of this, the tabular option can mask the lazy query behavior of an `rma::include` clause. The tabular option does not mask the inner-join behavior of an `rma::include` clause. The tabular filter is required for .csv format RMA requests.

options_clause (*self*, ***kwargs*)

build `rma::options` clause.

Parameters

only [list of strings, optional]

except [list of strings, optional]

tabular [list of string, optional]

count [boolean, optional]

debug [string, optional] ‘true’, ‘false’ or ‘preview’

num_rows [int or string, optional]

start_row [int or string, optional]

order_clause (*self*, *order_list=None*)

Construct a debug clause for use in an `rma::options` clause.

Parameters

order_list [list of strings] for example [‘acronym’, ‘products.name+asc’, ‘structure.id+desc’]

Returns

clause [string] The query clause for inclusion in an RMA query URL.

Notes

Optionally adding ‘+asc’ (default) or ‘+desc’ after an attribute will change the sort order.

pipe_stage (*self*, *pipe_name*, *parameters*)

Connect model and service stages via their JSON responses.

Notes

See: [Service Pipelines](#) and [Connected Services and Pipes](#)

quote_string (*self*, *the_string*)

Wrap a clause in single quotes.

Parameters

the_string [string] a clause to be included in an rma query that needs to be quoted

Returns

string input wrapped in single quotes

service_query (*self*, **args*, ***kwargs*)

Construct and Execute a single-stage RMA query to send a request to a connected service.

Parameters

service_name [string] Name of a documented connected service.

parameters [dict] key-value pairs as in the online documentation.

Notes

See: [Service Pipelines](#) and [Connected Services and Pipes](#)

service_stage (*self*, *service_name*, *parameters=None*)

Construct an RMA query fragment to send a request to a connected service.

Parameters

service_name [string] Name of a documented connected service.

parameters [dict] key-value pairs as in the online documentation.

Notes

See: [Service Pipelines](#) and [Connected Services and Pipes](#)

tuple_filters (*self*, *filters*)

Construct an RMA filter clause.

Notes

See [RMA Path Syntax - Square Brackets for Filters](#) for additional documentation.

allensdk.api.queries.rma_pager module

```
class allensdk.api.queries.rma_pager.RmaPager
    Bases: object

    static pager (fn, *args, **kwargs)

allensdk.api.queries.rma_pager.pageable (total_rows=None, num_rows=None)
```

allensdk.api.queries.rma_template module

```
class allensdk.api.queries.rma_template.RmaTemplate (base_uri=None,
                                                    query_manifest=None)
    Bases: allensdk.api.queries.rma_api.RmaApi
    See: Atlas Drawings and Ontologies

    template_query (self, template_name, entry_name, **kwargs)
    to_filter_rhs (self, rhs)
```

allensdk.api.queries.svg_api module

```
class allensdk.api.queries.svg_api.SvgApi (base_uri=None)
    Bases: allensdk.api.api.Api

    build_query (self, section_image_id, groups=None, download=False)
        Build the URL that will fetch meta data for the specified structure.

        Parameters

            section_image_id [integer] Key of the object to be retrieved.

            groups [array of integers] Keys of the group labels to filter the svg types that are returned.

        Returns

            url [string] The constructed URL

    download_svg (self, section_image_id, groups=None, file_path=None)
        Download the svg file

    get_svg (self, section_image_id, groups=None)
        Get the svg document.
```

allensdk.api.queries.synchronization_api module

```
class allensdk.api.queries.synchronization_api.SynchronizationApi (base_uri=None)
    Bases: allensdk.api.api.Api

    HTTP client for image synchronization services uses the image alignment results from the Informatics Data Processing Pipeline. Note: all locations on SectionImages are reported in pixel coordinates and all locations in 3-D ReferenceSpaces are reported in microns.

    See Image to Image Synchronization for additional documentation.

    get_image_to_atlas (self, section_image_id, x, y, atlas_id)
        For a specified Atlas, find the closest annotated SectionImage and (x,y) location as defined by a seed SectionImage and seed (x,y) location.
```

Parameters

section_image_id [integer] Seed for spatial sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

atlas_id [int] Target Atlas for image sync.

Returns

dict The parsed json response

get_image_to_image (*self, section_image_id, x, y, section_data_set_ids*)

For a list of target SectionDataSets, find the closest SectionImage and (x,y) location as defined by a seed SectionImage and seed (x,y) pixel location.

Parameters

section_image_id [integer] Seed for spatial sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

section_data_set_ids [list of integers] Target SectionDataSet IDs for image sync.

Returns

dict The parsed json response

get_image_to_image_2d (*self, section_image_id, x, y, section_image_ids*)

For a list of target SectionImages, find the closest (x,y) location as defined by a seed SectionImage and seed (x,y) location.

Parameters

section_image_id [integer] Seed for image sync.

x [float] Pixel coordinate of the seed location in the seed SectionImage.

y [float] Pixel coordinate of the seed location in the seed SectionImage.

section_image_ids [list of ints] Target SectionImage IDs for image sync.

Returns

dict The parsed json response

get_image_to_reference (*self, section_image_id, x, y*)

For a specified SectionImage and (x,y) location, return the (x,y,z) location in the ReferenceSpace of the associated SectionDataSet.

Parameters

section_image_id [integer] Seed for image sync.

x [float] Pixel coordinate on the specified SectionImage.

y [float] Pixel coordinate on the specified SectionImage.

Returns

dict The parsed json response

get_reference_to_image (*self, reference_space_id, x, y, z, section_data_set_ids*)

For a list of target SectionDataSets, find the closest SectionImage and (x,y) location as defined by a (x,y,z) location in a specified ReferenceSpace.

Parameters

reference_space_id [integer] Seed for spatial sync.

x [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

y [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

z [float] Coordinate (in microns) of the seed location in the seed ReferenceSpace.

section_data_set_ids [list of ints] Target SectionDataSets IDs for image sync.

Returns

dict The parsed json response

get_structure_to_image (*self, section_data_set_id, structure_ids*)

For a list of target structures, find the closest SectionImage and (x,y) location as defined by the centroid of each Structure.

Parameters

section_data_set_id [integer] primary key

structure_ids [list of integers] primary key

Returns

dict The parsed json response

allensdk.api.queries.tree_search_api module

class allensdk.api.queries.tree_search_api.**TreeSearchApi** (*base_uri=None*)

Bases: [allensdk.api.api.Api](#)

See [Searching a Specimen or Structure Tree](#) for additional documentation.

get_tree (*self, kind, db_id, ancestors=None, descendants=None*)

Fetch meta data for the specified structure or specimen.

Parameters

kind [string] ‘Structure’ or ‘Specimen’

db_id [integer] The id of the structure or specimen to search.

ancestors [boolean, optional] whether to include ancestors in the response (defaults to False)

descendants [boolean, optional] whether to include descendants in the response (defaults to False)

Returns

dict parsed json response data

Module contents**Submodules****allensdk.api module**

class allensdk.api.**Api** (*api_base_url_string=None*)

Bases: object

cleanup_truncated_file (*self*, *file_path*)

Helper for removing files.

Parameters

file_path [string] Absolute path including the file name to remove.

construct_well_known_file_download_url (*self*, *well_known_file_id*)

Join data api endpoint and id.

Parameters

well_known_file_id [integer or string representing an integer] well known file id

Returns

string the well-known-file download url for the current api server

See also:

[*retrieve_file_over_http*](#) Can be used to retrieve the file from the url.

default_api_url = 'http://api.brain-map.org'

do_query (*self*, *url_builder_fn*, *json_traversal_fn*, **args*, ***kwargs*)

Bundle an query url construction function with a corresponding response json traversal function.

Parameters

url_builder_fn [function] A function that takes parameters and returns an rma url.

json_traversal_fn [function] A function that takes a json-parsed python data structure and returns data from it.

post [boolean, optional kwarg] True does an HTTP POST, False (default) does a GET

args [arguments] Arguments to be passed to the url builder function.

kwargs [keyword arguments] Keyword arguments to be passed to the rma builder function.

Returns

any type The data extracted from the json response.

Examples

A simple Api subclass example.

do_rma_query (*self*, *rma_builder_fn*, *json_traversal_fn*, **args*, ***kwargs*)

Bundle an RMA query url construction function with a corresponding response json traversal function.

..note:: **Deprecated in AllenSDK 0.9.2** *do_rma_query* will be removed in AllenSDK 1.0, it is replaced by *do_query* because the latter is more general.

Parameters

rma_builder_fn [function] A function that takes parameters and returns an rma url.

json_traversal_fn [function] A function that takes a json-parsed python data structure and returns data from it.

args [arguments] Arguments to be passed to the rma builder function.

kwargs [keyword arguments] Keyword arguments to be passed to the rma builder function.

Returns

any type The data extracted from the json response.

Examples

A simple `Api` subclass example.

```
download_url = 'http://download.alleninstitute.org'
```

```
json_msg_query(self, url, dataframe=False)
```

Common case where the url is fully constructed and the response data is stored in the 'msg' field.

Parameters

url [string] Where to get the data in json form

dataframe [boolean] True converts to a pandas dataframe, False (default) doesn't

Returns

dict or DataFrame returned data; type depends on dataframe option

```
load_api_schema(self)
```

Download the RMA schema from the current RMA endpoint

Returns

dict the parsed json schema message

Notes

This information and other [Allen Brain Atlas Data Portal Data Model](#) documentation is also available as a [Class Hierarchy](#) and [Class List](#).

```
read_data(self, parsed_json)
```

Return the message data from the parsed query.

Parameters

parsed_json [dict] A python structure corresponding to the JSON data returned from the API.

Notes

See [API Response Formats - Response Envelope](#) for additional documentation.

```
retrieve_file_over_http(self, url, file_path, zipped=False)
```

Get a file from the data api and save it.

Parameters

url [string] `Url[R099781a1d33c-1]_` from which to get the file.

file_path [string] Absolute path including the file name to save.

zipped [bool, optional] If true, assume that the response is a zipped directory and attempt to extract contained files into the directory containing `file_path`. Default is False.

See also:

`construct_well_known_file_download_url` Can be used to construct the url.

References

[1]

retrieve_parsed_json_over_http (*self, url, post=False*)

Get the document and put it in a Python data structure

Parameters

url [string] Full API query url.

post [boolean] True does an HTTP POST, False (default) encodes the URL and does a GET

Returns

dict Result document as parsed by the JSON library.

retrieve_xml_over_http (*self, url*)

Get the document and put it in a Python data structure

Parameters

url [string] Full API query url.

Returns

string Unparsed xml string.

set_api_urls (*self, api_base_url_string*)

Set the internal RMA and well known file download endpoint urls based on a api server endpoint.

Parameters

api_base_url_string [string] url of the api to point to

set_default_working_directory (*self, working_directory*)

Set the working directory where files will be saved.

Parameters

working_directory [string] the absolute path string of the working directory.

`allensdk.api.api.stream_file_over_http` (*url, file_path, timeout=(9.05, 31.1)*)

Supply an http get request and stream the response to a file.

Parameters

url [str] Send the request to this url

file_path [str] Stream the response to this path

timeout [float or tuple of float, optional] Specify a timeout for the request. If a tuple, specify separate connect and read timeouts.

`allensdk.api.api.stream_zip_directory_over_http` (*url, directory, members=None, timeout=(9.05, 31.1)*)

Supply an http get request and stream the response to a file.

Parameters

url [str] Send the request to this url

directory [str] Extract the response to this directory

members [list of str, optional] Extract only these files

timeout [float or tuple of float, optional] Specify a timeout for the request. If a tuple, specify separate connect and read timeouts.

allensdk.api.cache module

class allensdk.api.cache.Cache (manifest=None, cache=True, version=None, **kwargs)
Bases: object

add_manifest_paths (self, manifest_builder)

Add cache-class specific paths to the manifest. In derived classes, should call super.

build_manifest (self, file_name)

Creation of default path specifications.

Parameters

file_name [string] where to save it

static cache_csv ()

static cache_csv_dataframe ()

static cache_csv_json ()

static cache_json ()

static cache_json_dataframe ()

static cacher (fn, *args, **kwargs)

make an rma query, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

strategy [string or None, optional] 'create' always generates the data, 'file' loads from disk, 'lazy' queries the server if no file exists, None generates the data and bypasses all caching behavior

pre [function] dfljson->dfljson, takes one data argument and returns filtered version, None for pass-through

post [function] dfljson->?, takes one data argument and returns Object

reader [function, optional] path -> data, default NOP

writer [function, optional] path, data -> None, default NOP

kwargs [objects] passed through to the query function

Returns

Object or None data type depends on fn, reader and/or post methods.

static csv_writer (pth, gen)

get_cache_path (self, file_name, manifest_key, *args)

Helper method for accessing path specs from manifest keys.

Parameters

file_name [string]

manifest_key [string]

args [ordered parameters]

Returns

string or None path

static json_remove_keys (*data, keys*)

static json_rename_columns (*data, new_old_name_tuples=None*)

Convenience method to rename columns in a pandas dataframe.

Parameters

data [dataframe] edited in place.

new_old_name_tuples [list of string tuples (new, old)]

load_csv (*self, path, rename=None, index=None*)

Read a csv file as a pandas dataframe.

Parameters

rename [list of string tuples (new old), optional] columns to rename

index [string, optional] post-rename column to use as the row label.

load_json (*self, path, rename=None, index=None*)

Read a json file as a pandas dataframe.

Parameters

rename [list of string tuples (new old), optional] columns to rename

index [string, optional] post-rename column to use as the row label.

load_manifest (*self, file_name, version=None*)

Read a keyed collection of path specifications.

Parameters

file_name [string] path to the manifest file

Returns

Manifest

manifest_dataframe (*self*)

Convenience method to view manifest as a pandas dataframe.

static nocache_dataframe ()

static nocache_json ()

static pathfinder (*file_name_position,* *secondary_file_name_position=None,*
path_keyword=None)

helper method to find path argument in legacy methods written prior to the @cacheable decorator. Do not use for new @cacheable methods.

Parameters

file_name_position [integer] zero indexed position in the decorated method args where file path may be found.

secondary_file_name_position [integer] zero indexed position in the decorated method args where the file path may be found.

path_keyword [string] kwarg that may have the file path.

Notes

This method is only intended to provide backward-compatibility for some methods that otherwise do not follow the path conventions of the `@cacheable` decorator.

static remove_keys (*data*, *keys=None*)

DataFrame version

static rename_columns (*data*, *new_old_name_tuples=None*)

Convenience method to rename columns in a pandas dataframe.

Parameters

data [dataframe] edited in place.

new_old_name_tuples [list of string tuples (new, old)]

wrap (*self*, *fn*, *path*, *cache*, *save_as_json=True*, *return_dataframe=False*, *index=None*, *rename=None*, ***kwargs*)

make an rma query, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

cache [boolean] True will make the query, False just loads from disk

save_as_json [boolean, optional] True (default) will save data as json, False as csv

return_dataframe [boolean, optional] True will cast the return value to a pandas dataframe, False (default) will not

index [string, optional] column to use as the pandas index

rename [list of string tuples, optional] (new, old) columns to rename

kwargs [objects] passed through to the query function

Returns

dict or DataFrame data type depends on `return_dataframe` option.

Notes

Column renaming happens after the file is reloaded for json

`allensdk.api.cache.cacheable` (*strategy=None*, *pre=None*, *writer=None*, *reader=None*,
post=None, *pathfinder=None*)

decorator for rma queries, save it and return the dataframe.

Parameters

fn [function reference] makes the actual query using kwargs.

path [string] where to save the data

strategy [string or None, optional] 'create' always gets the data from the source (server or generated), 'file' loads from disk, 'lazy' creates the data and saves to file if no file exists, None queries the server and bypasses all caching behavior

pre [function] `dfjson->dfjson`, takes one data argument and returns filtered version, None for pass-through

post [function] `dfjson->?`, takes one data argument and returns Object

reader [function, optional] path -> data, default NOP
writer [function, optional] path, data -> None, default NOP
kwargs [objects] passed through to the query function

Returns

dict or DataFrame data type depends on dataframe option.

Notes

Column renaming happens after the file is reloaded for json

`allensdk.api.cache.get_default_manifest_file(cache_name)`

`allensdk.api.cache.memoize(f)`

Creates an unbound cache of function calls and results. Note that arguments of different types are not cached separately (so `f(3.0)` and `f(3)` are not treated as distinct calls)

Arguments to the cached function must be hashable.

View the cache size with `f.cache_size()`. Clear the cache with `f.cache_clear()`. Access the underlying function with `f.__wrapped__`.

allensdk.api.caching_utilities module

```
allensdk.api.caching_utilities.call_caching(fetch: Callable[[Q], ~Q], write:
                                           Callable[[~Q], NoneType], read:
                                           Union[Callable[[Q], ~P], NoneType] =
                                           None, pre_write: Union[Callable[[~Q],
                                           ~Q], NoneType] = None, cleanup:
                                           Union[Callable[[Q], NoneType], None-
                                           Type] = None, lazy: bool = True, num_tries:
                                           int = 1, failure_message: str = "") →
                                           Union[~P, NoneType]
```

Access data, caching on a local store for future accesses.

Parameters

fetch : Function which pulls data from a remote/expensive source.
write : Function which stores data in a local/inexpensive store.
read : Function which pulls data from a local/inexpensive store.
pre_write : Function applied to obtained data after fetching, but before writing.
cleanup : Function for fixing a failed fetch. e.g. unlinking a partially downloaded file. Exceptions raised by cleanup are not themselves handled
lazy : If True, attempt to read the data from the local/inexpensive store before fetching it. If False, forcibly fetch from the remote/expensive store.
num_tries : How many fetches to attempt before (re)raising an exception. A fetch is failed if reading the result raises an exception.
failure_message : Provides additional context in the event of a failed download. Emitted when retrying, and when a fetch failure occurs after tries are exhausted

Returns

The result of calling read

```
allensdk.api.caching_utilities.one_file_call_caching(path: Union[pathlib.Path, str],
                                                    fetch: Callable[[], ~Q], write:
                                                    Callable[[Union[pathlib.Path,
                                                                    str], ~Q], NoneType], read:
                                                    Union[Callable[[Union[pathlib.Path,
                                                                    str]], ~P], NoneType]
                                                    = None, pre_write:
                                                    Union[Callable[[~Q], ~Q],
                                                    NoneType] = None, cleanup:
                                                    Union[Callable[[], NoneType],
                                                    NoneType] = None, lazy: bool
                                                    = True, num_tries: int = 1,
                                                    failure_message: str = ") →
                                                    Union[~P, NoneType]
```

A call_caching variant where the local store is a single file. See call_caching for complete documentation.

Parameters

path : Path at which the data will be stored

Module contents

Subclasses of allensdk.api.api.Api to implement specific queries to the [Allen Brain Atlas Data Portal](#).

6.1.2 allensdk.brain_observatory package

Subpackages

allensdk.brain_observatory.behavior package

Subpackages

allensdk.brain_observatory.behavior.behavior_ophys_api package

Submodules

allensdk.brain_observatory.behavior.behavior_ophys_api.behavior_ophys_nwb_api module

Module contents

```
class allensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase
    Bases: object

    get_average_projection(self)
    get_cell_specimen_table(self)
    get_corrected_fluorescence_traces(self)
    get_dff_traces(self)
    get_eye_tracking_data(self)
```

```

get_licks(self)
get_max_projection(self)
get_metadata(self)
get_motion_correction(self)
get_ophys_experiment_id(self) → int
get_ophys_timestamps(self)
get_rewards(self)
get_running_data_df(self)
get_running_speed(self)
get_segmentation_mask_image(self)
get_stimulus_presentations(self)
get_stimulus_templates(self)
get_stimulus_timestamps(self)
get_task_parameters(self)
get_trials(self)

```

allensdk.brain_observatory.behavior.internal package

Submodules

allensdk.brain_observatory.behavior.internal.behavior_base module

```

class allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase
    Bases: abc.ABC

```

Abstract base class implementing required methods for interacting with behavior session data.

Child classes should be instantiated with a fetch API that implements these methods.

```
get_behavior_session_id(self) → int
```

Returns the behavior_session_id associated with this experiment, if applicable.

```
get_licks(self) → pandas.core.frame.DataFrame
```

Get lick data from pkl file.

Returns

np.ndarray A dataframe containing lick timestamps.

```
get_rewards(self) → pandas.core.frame.DataFrame
```

Get reward data from pkl file.

Returns

pd.DataFrame A dataframe containing timestamps of delivered rewards.

```
get_running_data_df(self) → pandas.core.frame.DataFrame
```

Get running speed data.

Returns

pd.DataFrame Dataframe containing various signals used to compute running speed.

get_running_speed (*self*) → allensdk.brain_observatory.running_speed.RunningSpeed
Get running speed using timestamps from self.get_stimulus_timestamps.

NOTE: Do not correct for monitor delay.

Returns

RunningSpeed (NamedTuple with two fields)

timestamps [np.ndarray] Timestamps of running speed data samples

values [np.ndarray] Running speed of the experimental subject (in cm / s).

get_stimulus_presentations (*self*) → pandas.core.frame.DataFrame
Get stimulus presentation data.

NOTE: Uses timestamps that do not account for monitor delay.

Returns

pd.DataFrame Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_templates (*self*) → Dict[str, numpy.ndarray]
Get stimulus templates (movies, scenes) for behavior session.

Returns

Dict[str, np.ndarray] A dictionary containing the stimulus images presented during the session. Keys are data set names, and values are 3D numpy arrays.

get_stimulus_timestamps (*self*) → numpy.ndarray
Get stimulus timestamps from pkl file.

NOTE: Located with behavior_session_id

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor that do not account for monitor delay.

get_task_parameters (*self*) → dict
Get task parameters from pkl file.

Returns

dict A dictionary containing parameters used to define the task runtime behavior.

get_trials (*self*) → pandas.core.frame.DataFrame
Get trials from pkl file

Returns

pd.DataFrame A dataframe containing behavioral trial start/stop times, and trial data

allensdk.brain_observatory.behavior.internal.behavior_ophys_base module

class allensdk.brain_observatory.behavior.internal.behavior_ophys_base.**BehaviorOphysBase**
Bases: [allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase](#)

Abstract base class implementing required methods for interacting with behavior+ophys session data.

Child classes should be instantiated with a fetch API that implements these methods.

get_average_projection (*self*) → allensdk.brain_observatory.behavior.image_api.Image

Get an image whose values are the average obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to avg projection image data and metadata.

get_cell_specimen_table (*self*) → pandas.core.frame.DataFrame

Get a cell specimen dataframe containing ROI information about cells identified in an ophys experiment.

Returns

pd.DataFrame Cell ROI information organized into a dataframe. Index is the cell ROI IDs.

get_corrected_fluorescence_traces (*self*) → pandas.core.frame.DataFrame

Get motion-corrected fluorescence traces.

Returns

pd.DataFrame Motion-corrected fluorescence traces organized into a dataframe. Index is the cell ROI IDs.

get_dff_traces (*self*) → pandas.core.frame.DataFrame

Get a table of delta fluorescence over fluorescence traces.

Returns

pd.DataFrame The traces of dff (normalized fluorescence) organized into a dataframe. Index is the cell ROI IDs.

get_max_projection (*self*) → allensdk.brain_observatory.behavior.image_api.Image

Get an image whose values are the maximum obtained values at each pixel of the ophys movie over time.

Returns

allensdk.brain_observatory.behavior.image_api.Image: Array-like interface to max projection image data and metadata.

get_metadata (*self*) → dict

Get behavior+ophys session metadata.

Returns

dict A dictionary of session-specific metadata.

get_motion_correction (*self*) → pandas.core.frame.DataFrame

Get motion correction trace data.

Returns

pd.DataFrame A dataframe containing trace data used during motion correction computation.

get_ophys_session_id (*self*) → Union[int, NoneType]

Returns the ophys_session_id associated with this experiment, if applicable.

get_ophys_timestamps (*self*) → numpy.ndarray

Get optical physiology frame timestamps.

Returns

np.ndarray Timestamps associated with frames captured by the microscope.

get_raw_stimulus_timestamps (*self*) → numpy.ndarray

Get raw stimulus timestamps.

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor without accounting for monitor delay.

get_stimulus_presentations (*self*) → pandas.core.frame.DataFrame

Get stimulus presentation data.

NOTE: Uses monitor delay corrected stimulus timestamps.

Returns

pd.DataFrame Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_timestamps (*self*) → numpy.ndarray

Get stimulus timestamps.

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor after accounting for monitor delay.

allensdk.brain_observatory.behavior.internal.behavior_project_base module

Module contents

allensdk.brain_observatory.behavior.sync package

Submodules

allensdk.brain_observatory.behavior.sync.process_sync module

allensdk.brain_observatory.behavior.sync.process_sync.**calculate_delay** (*sync_data*,
stim_vsync_fall,
sample_frequency)

allensdk.brain_observatory.behavior.sync.process_sync.**filter_digital** (*rising*,
falling,
threshold=0.0001)

Removes short transients from digital signal.

Rising and falling should be same length and units in seconds.

Kwargs: threshold (float): transient width

Module contents

Created on Sunday July 15 2018

@author: marinag

Parameters

Returns

```
allensdk.brain_observatory.behavior.sync.get_behavior_monitoring(dataset: al-
lensdk.brain_observatory.sync_dataset,
permissive:
bool =
False) →
Union[numpy.ndarray,
NoneType]
```

Parameters

permissive [If True, None will be returned if timestamps are not found. If] False, a `KeyError` will be raised

Returns

```
allensdk.brain_observatory.behavior.sync.get_eye_tracking(dataset: al-
lensdk.brain_observatory.sync_dataset.Dataset,
permissive: bool
= False) →
Union[numpy.ndarray,
NoneType]
```

Parameters

permissive [If True, None will be returned if timestamps are not found. If] False, a `KeyError` will be raised

Returns

array of timestamps (floating point; seconds; relative to experiment start) or **None**. If None, no eye tracking timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_lick_times(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool = False)  
→ Union[numpy.ndarray,  
NoneType]
```

Report the timestamps of each detected lick

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no lick timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_ophys_frames(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool  
= False) →  
numpy.ndarray
```

Report the timestamps of each optical physiology video frame

Parameters

dataset [describes experiment timing]

Returns

array of timestamps (floating point; seconds; relative to experiment start).

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Notes

use rising edge for Scientifica, falling edge for Nikon <http://confluence.corp.alleninstitute.org/display/IT/Ophys+Time+Sync> This function uses rising edges

```
allensdk.brain_observatory.behavior.sync.get_raw_stimulus_frames(dataset: al-  
lensdk.brain_observatory.sync_dataset.Dataset,  
permissive: bool  
= False) →  
numpy.ndarray
```

Report the raw timestamps of each stimulus frame. This corresponds to the time at which the psychopy window's flip method returned, but not necessarily to the time at which the stimulus frame was displayed.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

array of timestamps (floating point; seconds; relative to experiment start).

```
allensdk.brain_observatory.behavior.sync.get_stim_photodiode (dataset: al-  

                                                                lensdk.brain_observatory.sync_dataset.Dataset,  

                                                                permissive: bool  

                                                                = False) →  

                                                                Union[List[float],  

                                                                NoneType]
```

Report the timestamps of each detected sync square transition (both black -> white and white -> black) in this experiment.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

array of timestamps (floating point; seconds; relative to experiment start) or None. If None, no photodiode timestamps were found in this sync dataset.

```
allensdk.brain_observatory.behavior.sync.get_stimulus_rebase_function (data,  

                                                                stim-  

                                                                u-  

                                                                lus_timestamps_no_monitor_a
```

Create a rebase function to align times for licks and stimulus timestamps in the “pickle” log with the same events in the event “sync” log.

```
allensdk.brain_observatory.behavior.sync.get_sync_data (sync_path: str, permissive:  

                                                                bool = False) → Dict[str,  

                                                                Union[List, numpy.ndarray,  

                                                                NoneType]]
```

Convenience function for extracting several timestamp arrays from a sync file.

Parameters

sync_path [The hdf5 file here ought to be a Visual Behavior sync output] file. See allensdk.brain_observatory.sync_dataset for more details of this format.

permissive [If True, None will be returned if timestamps are not found. If False, a KeyError will be raised]

Returns

A dictionary with the following keys. All timestamps in seconds: ophys_frames : timestamps of each optical physiology frame lick_times : timestamps of each detected lick ophys_trigger : The time at which ophys acquisition was started eye_tracking : timestamps of each eye tracking video frame behavior_monitoring : timestamps of behavior monitoring video frame stim_photodiode : timestamps of each photodiode transition stimulus_times_no_delay : raw stimulus frame timestamps

Some values may be None. This indicates that the corresponding timestamps were not located in this sync file.

```
allensdk.brain_observatory.behavior.sync.get_trigger (dataset: al-  

                                                                lensdk.brain_observatory.sync_dataset.Dataset,  

                                                                permissive: bool = False) →  

                                                                Union[numpy.ndarray, None-  

                                                                Type]
```

Returns (as a 1-element array) the time at which optical physiology acquisition was started.

Parameters

dataset [describes experiment timing]

permissive [If True, None will be returned if timestamps are not found. If] False, a Key-Error will be raised

Returns

timestamps (floating point; seconds; relative to experiment start) or None. If None, no timestamps were found in this sync dataset.

Notes

Ophys frame timestamps can be recorded before acquisition start when experimenters are setting up the recording session. These do not correspond to acquired ophys frames.

allensdk.brain_observatory.behavior.write_nwb package**Module contents****Submodules****allensdk.brain_observatory.behavior.behavior_data_session module**

class allensdk.brain_observatory.behavior.behavior_data_session.**BehaviorDataSession** (*api: Optional[Type[BehaviorDataSession]] = None*)

Bases: object

behavior_session_id

Unique identifier for this experimental session. :rtype: int

cache_clear (*self*) → None

Convenience method to clear the api cache, if applicable.

classmethod from_lims (*behavior_session_id: int*) → 'BehaviorDataSession'**classmethod from_nwb_path** (*nwb_path: str, **api_kwargs: Any*) → 'BehaviorDataSession'**licks**

Get lick data from pkl file.

Returns

np.ndarray A dataframe containing lick timestamps.

list_api_methods (*self*) → List[Tuple[str, str]]

Convenience method to expose list of API *get* methods. These methods can be accessed by referencing the API used to initialize this BehaviorDataSession via its *api* instance attribute. :rtype: list of tuples, where the first value in the tuple is the method name, and the second value is the method docstring.

metadata

Return metadata about the session. :rtype: dict

ophys_experiment_ids

The unique identifiers for the ophys experiment(s) associated with this behavior session (if one exists)
:rtype: int

ophys_session_id

The unique identifier for the ophys session associated with this behavior session (if one exists) :rtype: int

raw_running_data_df

Get running speed data without the 10Hz low pass filter.

Returns

pd.DataFrame Dataframe containing various signals used to compute running speed, and the unfiltered speed.

raw_running_speed

Get running speed using timestamps from self.get_stimulus_timestamps. Does not apply the 10Hz filter.

NOTE: Do not correct for monitor delay.

Returns**RunningSpeed (NamedTuple with two fields)**

timestamps [np.ndarray] Timestamps of running speed data samples

values [np.ndarray] Running speed of the experimental subject (in cm / s).

rewards

Get reward data from pkl file.

Returns

pd.DataFrame A dataframe containing timestamps of delivered rewards.

running_data_df

Get running speed data. By default applies a 10Hz low pass filter to the data. To get the running speed without the filter, use *raw_running_data_df*.

Returns

pd.DataFrame Dataframe containing various signals used to compute running speed, and the filtered speed.

running_speed

Get running speed using timestamps from self.get_stimulus_timestamps. Applies a 10Hz filter to the speed. To get the unfiltered speed, use *raw_running_speed*.

NOTE: Do not correct for monitor delay.

Returns**RunningSpeed (NamedTuple with two fields)**

timestamps [np.ndarray] Timestamps of running speed data samples

values [np.ndarray] Running speed of the experimental subject (in cm / s).

stimulus_presentations

Get stimulus presentation data.

NOTE: Uses timestamps that do not account for monitor delay.

Returns

pd.DataFrame Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

stimulus_templates

Get stimulus templates (movies, scenes) for behavior session.

Returns

Dict[str, np.ndarray] A dictionary containing the stimulus images presented during the session. Keys are data set names, and values are 3D numpy arrays.

stimulus_timestamps

Get stimulus timestamps from pkl file.

NOTE: Located with behavior_session_id

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor that do not account for monitor delay.

task_parameters

Get task parameters from pkl file.

Returns

dict A dictionary containing parameters used to define the task runtime behavior.

trials

Get trials from pkl file

Returns

pd.DataFrame A dataframe containing behavioral trial start/stop times, and trial data

allensdk.brain_observatory.behavior.behavior_ophys_analysis module

allensdk.brain_observatory.behavior.behavior_ophys_session module

allensdk.brain_observatory.behavior.behavior_project_cache module

allensdk.brain_observatory.behavior.behavior_project_lims_api module

allensdk.brain_observatory.behavior.criteria module

Functions for calculating mtrain state transitions. If criteria are met, return true. Otherwise, return false.

allensdk.brain_observatory.behavior.criteria.consistency_is_key (*session_summary*)

need some way to judge consistency of various parameters

- dprime
- num trials
- hit rate
- fa rate
- lick timing

allensdk.brain_observatory.behavior.criteria.consistent_behavior_within_session (*session_summary*)

need some way to measure consistent performance within a session

- compare peak to overall dprime?

- variance in rolling window dprime?

`allensdk.brain_observatory.behavior.criteria.meets_engagement_criteria(session_summary)`

Returns true if engagement criteria were met for the past 3 days, else false. Args:

`session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak' and 'num_engaged_trials', ordered ascending by training day, for at least 3 days. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function) The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.mostly_useful(trials)`

Returns True if fewer than half the trial time on the last day were aborted trials.

Args: `trials` (pd.DataFrame): Pandas dataframe with columns 'training_day', 'trial_type', and 'trial_length'.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.n_complete(threshold, count)`

For compatibility with original API. If count >= threshold, return True. Otherwise return False. Args:

`threshold` (numeric): Threshold for the count to meet. `count` (numeric): The count to compare to the threshold.

Returns: True if count >= threshold, otherwise False.

`allensdk.brain_observatory.behavior.criteria.no_response_bias(session_summary)`

the mouse meets this criterion if their last session exhibited a response bias between 10% and 90%

Args: `session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'response_bias', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.summer_over(trials)`

Returns true if the maximum value of 'training_day' in the trials dataframe is >= 40, else false.

`allensdk.brain_observatory.behavior.criteria.two_out_of_three_aint_bad(session_summary)`

Returns true if 2 of the last 3 days showed a peak d-prime above 2.

Args: `session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak', ordered ascending by training day, for at least the past 3 days. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by mtrain function). The mtrain implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.whole_lotta_trials(session_summary)`

Mouse meets this criterion if the last session has more than 300 trials. Args:

`session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'num_contingent_trials', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by `mtrain` function). The `mtrain` implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.criteria.yesterday_was_good(session_summary)`

Returns true if the last day showed a peak d-prime above 2 Args:

`session_summary` (pd.DataFrame): Pandas dataframe with daily values for 'dprime_peak', ordered ascending by training day, for at least 1 day. If dataframe is not properly ordered, criterion may not be correctly calculated. This function does not sort the data to preserve prior behavior (sorting column was not required by `mtrain` function). The `mtrain` implementation created the required columns if they didn't exist, so a more informative error is raised here to assist end-users in debugging.

Returns: bool: True if criterion is met, False otherwise

`allensdk.brain_observatory.behavior.dprime` module

`allensdk.brain_observatory.behavior.dprime.get_catch_responses` (*correct_reject=None*,
false_alarm=None,
aborted=None)

`allensdk.brain_observatory.behavior.dprime.get_dprime` (*hit_rate*, *fa_rate*, *sliding_window=100*)

calculates the d-prime for a given hit rate and false alarm rate https://en.wikipedia.org/wiki/Sensitivity_index

Parameters ——— `hit_rate` : float

rate of hits in the True class

fa_rate [float] rate of false alarms in the False class

limits [tuple, optional] limits on extreme values, which distort. default: (0.01,0.99)

d_prime

`allensdk.brain_observatory.behavior.dprime.get_false_alarm_rate` (*correct_reject=None*,
false_alarm=None,
aborted=None,
sliding_window=100)

`allensdk.brain_observatory.behavior.dprime.get_go_responses` (*hit=None*,
miss=None,
aborted=None)

`allensdk.brain_observatory.behavior.dprime.get_hit_rate` (*hit=None*, *miss=None*,
aborted=None, *sliding_window=100*)

`allensdk.brain_observatory.behavior.dprime.get_rolling_dprime` (*rolling_hit_rate*,
rolling_fa_rate,
sliding_window=100)


```
allensdk.brain_observatory.behavior.dprime.get_trial_count_corrected_false_alarm_rate(corrected,
                                                                                   false_alarm_rate,
                                                                                   hit_rate,
                                                                                   miss_rate,
                                                                                   aborted_rate,
                                                                                   sliding_window)

allensdk.brain_observatory.behavior.dprime.get_trial_count_corrected_hit_rate(hit=None,
                                                                              miss=None,
                                                                              aborted=None,
                                                                              sliding_window=100)

allensdk.brain_observatory.behavior.dprime.trial_number_limit(p, N)
```

allensdk.brain_observatory.behavior.eye_tracking_processing module

```
allensdk.brain_observatory.behavior.eye_tracking_processing.compute_circular_area(df_row:
                                                                                  pan-
                                                                                  das.core.series
                                                                                  →
                                                                                  float)
```

Calculate the area of the pupil as a circle using the max of the height/width as radius.

Note: This calculation assumes that the pupil is a perfect circle and any eccentricity is a result of the angle at which the pupil is being viewed.

Parameters

df_row [pd.Series] A row from an eye tracking dataframe containing only “pupil_width” and “pupil_height”.

Returns

float The circular area of the pupil in pixels².

```
allensdk.brain_observatory.behavior.eye_tracking_processing.compute_elliptical_area(df_row:
                                                                                  pan-
                                                                                  das.core.series
                                                                                  →
                                                                                  float)
```

Calculate the area of corneal reflection (cr) or eye ellipse fits using the ellipse formula.

Parameters

df_row [pd.Series] A row from an eye tracking dataframe containing either: “cr_width”, “cr_height” or “eye_width”, “eye_height”

Returns

float The elliptical area of the eye or cr in pixels²

```
allensdk.brain_observatory.behavior.eye_tracking_processing.determine_likely_blinks(eye_areas: pandas.core.series.Series, pupil_areas: pandas.core.series.Series, outliers: pandas.core.series.Series, dilation_frames: int) → pandas.core.series.Series
```

Determine eye tracking frames which contain likely blinks or outliers

Parameters

eye_areas [pd.Series] A pandas series of eye areas.

pupil_areas [pd.Series] A pandas series of pupil areas.

outliers [pd.Series] A pandas series containing bool values of outlier rows.

dilation_frames [int, optional] Determines the number of additional adjacent frames to mark as ‘likely_blink’, by default 2.

Returns

pd.Series A pandas series of bool values that has the same length as the number of eye tracking dataframe rows (frames).

```
allensdk.brain_observatory.behavior.eye_tracking_processing.determine_outliers(data_df: pandas.DataFrame, z_threshold: float) → pandas.core.series.Series
```

Given a dataframe and some z-score threshold return a pandas boolean Series where each entry indicates whether a given row contains at least one outlier (where outliers are calculated along columns).

Parameters

data_df [pd.DataFrame] A dataframe containing only columns where outlier detection is desired. (e.g. “cr_area”, “eye_area”, “pupil_area”)

z_threshold [float] z-score values higher than the z_threshold will be considered outliers.

Returns

pd.Series A pandas boolean Series whose length == len(data_df.index). True denotes that a row in the data_df contains at least one outlier.

```
allensdk.brain_observatory.behavior.eye_tracking_processing.load_eye_tracking_hdf(eye_tracking_file,
path-
lib.Path)
→
pan-
das.core.frame
```

Load a DeepLabCut hdf5 file containing eye tracking data into a dataframe.

Note: The eye tracking hdf5 file contains 3 separate dataframes. One for corneal reflection (cr), eye, and pupil ellipse fits. This function loads and returns this data as a single dataframe.

Parameters

eye_tracking_file [Path] Path to an hdf5 file produced by the DeepLabCut eye tracking pipeline. The hdf5 file will contain the following keys: “cr”, “eye”, “pupil”. Each key has an associated dataframe with the following columns: “center_x”, “center_y”, “height”, “width”, “phi”.

Returns

pd.DataFrame A dataframe containing combined corneal reflection (cr), eyelid (eye), and pupil data. Column names for each field will be renamed by prepending the field name. (e.g. center_x -> eye_center_x)

```
allensdk.brain_observatory.behavior.eye_tracking_processing.process_eye_tracking_data(eye_data,
pan-
das.co
frame_
pan-
das.co
z_thre:
float
=
3.0,
di-
la-
tion_fr
int
=
2)
→
pan-
das.co
```

Processes and refines raw eye tracking data by adding additional computed feature columns.

Parameters

eye_data [pd.DataFrame] A ‘raw’ eye tracking dataframe produced by load_eye_tracking_hdf()

frame_times [pd.Series] A series of frame times acquired from a behavior + ophy session ‘sync file’.

z_threshold [float] z-score values higher than the z_threshold will be considered outliers, by default 3.0.

dilation_frames [int, optional] Determines the number of additional adjacent frames to mark as ‘likely_blink’, by default 2.

Returns

pd.DataFrame A refined eye tracking dataframe that contains additional information about frame times, eye areas, pupil areas, and frames with likely blinks/outliers.

Raises

RuntimeError If the number of sync file frame times does not match the number of eye tracking frames.

allensdk.brain_observatory.behavior.image_api module

```
class allensdk.brain_observatory.behavior.image_api.Image
    Bases: tuple
    Describes a 2D Image
    data [np.ndarray] Image data points
    spacing [tuple] Spacing describes the physical size of each pixel
    unit [str] Physical unit of the spacing (currently constrained to be isotropic)
    data
        Alias for field number 0
    spacing
        Alias for field number 1
    unit
        Alias for field number 2
class allensdk.brain_observatory.behavior.image_api.ImageApi
    Bases: object
    static deserialize (img)
    static serialize (data, spacing, unit)
```

allensdk.brain_observatory.behavior.metadata_processing module

```
allensdk.brain_observatory.behavior.metadata_processing.get_task_parameters (data)
```

allensdk.brain_observatory.behavior.mtrain module

```
class allensdk.brain_observatory.behavior.mtrain.ExtendedTrialSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
    Bases: marshmallow.schema.Schema
    This schema describes the edf core trial structure
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.mtrain.FriendlyDate (format=None,
                                                                **kwargs)
```

```
Bases: marshmallow.fields.Date
```

```
class allensdk.brain_observatory.behavior.mtrain.FriendlyDateTime (format=None,
                                                                    **kwargs)
```

```
Bases: marshmallow.fields.DateTime
```

```
allensdk.brain_observatory.behavior.mtrain.annotate_change_detect (trials)
adds change and detect columns to dataframe
```

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

io.load_trials

```
allensdk.brain_observatory.behavior.mtrain.annotate_trials (trials)
performs multiple annotations:
```

- `annotate_change_detect`
- `fix_change_time`
- `explode_response_window`

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

io.load_trials

```
allensdk.brain_observatory.behavior.mtrain.assign_session_id (trials)
```

adds a column with a unique ID for the session defined as a combination of the mouse ID and startdatetime

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

io.load_trials

```
allensdk.brain_observatory.behavior.mtrain.explode_response_window (trials)
explodes the response_window column in lower & upper columns
```

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

`allensdk.brain_observatory.behavior.mtrain.fix_change_time(trials)`
forces *None* values in the *change_time* column to numpy NaN

Parameters

trials [pandas DataFrame] dataframe of trials

inplace [bool, optional] modify *trials* in place. if False, returns a copy. default: True

See also:

`io.load_trials`

allensdk.brain_observatory.behavior.rewards_processing module

`allensdk.brain_observatory.behavior.rewards_processing.get_rewards(data, stimulus_rebase_function)`

allensdk.brain_observatory.behavior.running_processing module

`allensdk.brain_observatory.behavior.running_processing.calc_deriv(x, time)`
`allensdk.brain_observatory.behavior.running_processing.deg_to_dist(angular_speed: numpy.ndarray) → numpy.ndarray`

Takes the angular speed (radians/s) at each step in radians, and computes the linear speed in cm/s.

Parameters

angular_speed: np.ndarray (1d) 1d array of angular speed in radians/s

Returns

—

np.ndarray (1d) Linear speed in cm/s at each time point.

`allensdk.brain_observatory.behavior.running_processing.get_running_df(data, time: numpy.ndarray, low-pass: bool = True)`

Given the data from the behavior ‘pkl’ file object and a 1d array of timestamps, compute the running speed. Returns a dataframe with the raw voltage data as well as the computed speed at each timestamp. By default, the running speed is filtered with a 10 Hz Butterworth lowpass filter to remove artifacts caused by the rotary encoder.

Parameters

data Deserialized ‘behavior.pkl’ file data

time: np.ndarray (1d) Timestamps for running data measurements

lowpass: bool (default=True) Whether to apply a 10Hz low-pass filter to the running speed data.

Returns

pd.DataFrame Dataframe with an index of timestamps and the following columns:

“speed”: computed running speed “dx”: angular change, computed during data collection “v_sig”: voltage signal from the encoder “v_in”: the theoretical maximum voltage that the encoder

will reach prior to “wrapping”. This should theoretically be 5V (after crossing 5V goes to 0V, or vice versa). In practice the encoder does not always reach this value before wrapping, which can cause transient spikes in speed at the voltage “wraps”.

The raw data are provided so that the user may compute their own speed from source, if desired.

Notes

Though the angular change is available in the raw data (key=”dx”), this method recomputes the angular change from the voltage signal (key=”vsig”) due to very specific, low-level artifacts in the data caused by the encoder. See method docstrings for more detailed information. The raw data is included in the final output in case the end user wants to apply their own corrections and compute running speed from the raw source.

allensdk.brain_observatory.behavior.schemas module

```
class allensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
par-
tial=False,
un-
known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.RaisingSchema`

This schema contains metadata pertaining to behavior.

opts = <marshmallow.schema.SchemaOpts object>

```
class allensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
par-
tial=False,
un-
known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.RaisingSchema`

This schema encompasses task parameters used for behavior or ophys + behavior.

```
neurodata_doc = 'Metadata for behavior or behavior + ophys task parameters'
```

```
neurodata_type = 'BehaviorTaskParameters'
```

```
neurodata_type_inc = 'LabMetaData'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.CompleteOphysBehaviorMetadataSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
partial=False,
unknown=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema`,
`allensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema`

This schema combines fields from behavior, ophys, and subject schemas. Metadata info is passed by the behavior+ophys session in a combined lump containing all the field types.

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema (only=None,
ex-
clude=(),
many=False,
con-
text=None,
load_only=(),
dump_only=(),
partial=False,
unknown=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema`,
`allensdk.brain_observatory.behavior.schemas.OphysMetadataSchema`

This schema contains fields pertaining to ophys+behavior. It is used as a template for generating our custom NWB behavior + ophys extension.

```
neurodata_doc = 'Metadata for behavior + ophys experiments'
```

```
neurodata_type = 'OphysBehaviorMetadata'
```

```
neurodata_type_inc = 'LabMetaData'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.OphysMetadataSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.RaisingSchema`

This schema contains metadata pertaining to optical physiology (ophys).

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.RaisingSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    context=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False, un-
                                                                    known=None)
```

Bases: `marshmallow.schema.Schema`

```
class Meta
```

```
    Bases: object
```

```
    unknown = 'raise'
```

```
opts = <marshmallow.schema.SchemaOpts object>
```

```
class allensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)
```

Bases: `allensdk.brain_observatory.behavior.schemas.RaisingSchema`

This schema contains metadata pertaining to a subject in either a behavior or behavior + ophys experiment.

```
neurodata_doc = 'Metadata for an AIBS behavior or behavior + ophys subject'
neurodata_type = 'BehaviorSubject'
neurodata_type_inc = 'Subject'
opts = <marshmallow.schema.SchemaOpts object>
```

allensdk.brain_observatory.behavior.session_metrics module

`allensdk.brain_observatory.behavior.session_metrics.num_contingent_trials(session_trials)`

Returns the number of “go” and “catch” trials in a training session dataframe. Args:

`session_trials` (pandas.DataFrame): a pandas.DataFrame describing behavior training trials, with the string column “trial_type” describing the type of trial.

Returns (int): Number of “go” and “catch” trials

```
allensdk.brain_observatory.behavior.session_metrics.response_bias(trials, detect_col,
                                                                    trial_types=('go',
                                                                    'catch'))
```

Calculate the response bias for a subset of trial types from a behavioral training dataframe. Args:

trials (pandas.DataFrame): Dataframe containing trial-level information from a behavioral training session. Required columns: “trial_type”, *detect_col*.

detect_col (str): Name of column containing boolean or numeric codings (0/1) for whether or not the mouse had a response.

trial_types (iterable<str>): Iterable containing string trial types to check for the response bias. Trials of types not included in this iterable will be ignored. Default=(“go”, “catch”)

Return: The response bias (or average value of the *detect_col*) for trials in *trial_types*.

allensdk.brain_observatory.behavior.stimulus_processing module

`allensdk.brain_observatory.behavior.stimulus_processing.convert_filepath_caseinsensitive(fil`

```
allensdk.brain_observatory.behavior.stimulus_processing.get_gratings_metadata(stimuli:
                                                                              Dict,
                                                                              start_idx:
                                                                              int
                                                                              =
                                                                              0)
                                                                              →
                                                                              pan-
                                                                              das.core.frame.DataFrame
```

This function returns the metadata for each unique grating that was presented during the experiment. If no gratings were displayed during this experiment it returns an empty dataframe with the expected columns. Parameters
—— stimuli: the stimuli loaded from the experiment pkl file start_idx: the index to start index column

Returns

pd.DataFrame: DataFrame containing the unique stimuli presented during an experiment. The columns contained in this DataFrame are ‘image_category’, ‘image_name’, ‘image_set’,

‘phase’, ‘spatial_frequency’, ‘orientation’, and ‘image_index’. This returns empty if no gratings were presented.

`allensdk.brain_observatory.behavior.stimulus_processing.get_images_dict (pkl)`
→

Dict

Gets the dictionary of images that were presented during an experiment along with image set metadata and the image specific metadata. This function uses the path to the image pkl file to read the images and their metadata from the pkl file and return this dictionary. Parameters ——— pkl: The pkl file containing the data for the stimuli presented during

experiment

Returns

Dict: A dictionary containing keys images, metadata, and image_attributes. These correspond to paths to images to images presented, metadata on the whole set of images, and metadata on specific images, respectively.

`allensdk.brain_observatory.behavior.stimulus_processing.get_stimulus_metadata (pkl)`
→

pan-

das.core.frame.Dat

Gets the stimulus metadata for each type of stimulus presented during the experiment. The metadata is return for gratings, images, and omitted stimuli. Parameters ——— pkl: the pkl file containing the information about what stimuli were

presented during the experiment

Returns

pd.DataFrame: The dataframe containing a row for every stimulus that was presented during the experiment. The row contains the following data, image_category, image_name, image_set, phase, spatial_frequency, orientation, and image index.

`allensdk.brain_observatory.behavior.stimulus_processing.get_stimulus_presentations (data,`
stim-

u-

lus_timesta

→

pan-

das.core.fra

This function retrieves the stimulus presentation dataframe and renames the columns, adds a stop_time column, and set's index to stimulus_presentation_id before sorting and returning the dataframe. :param data: stimulus file associated with experiment id :param stimulus_timestamps: timestamps indicating when stimuli switched

during experiment

Returns stimulus_table: dataframe containing the stimuli metadata as well as what stimuli was presented

`allensdk.brain_observatory.behavior.stimulus_processing.get_stimulus_templates (pkl)`
→

Dict

Gets dictionary of images presented during experimentation Parameters ——— pkl: pkl file containing the data for the presented stimuli

Returns

Dict: Dictionary of images that were presented during the experiment

```
allensdk.brain_observatory.behavior.stimulus_processing.get_visual_stimuli_df(data,  
                                                                           time)  
                                                                           →  
                                                                           pan-  
                                                                           das.core.frame.DataFrame
```

This function loads the stimuli and the omitted stimuli into a dataframe. These stimuli are loaded from the input data, where the `set_log` and `draw_log` contained within are used to calculate the epochs. These epochs are used as `start_frame` and `end_frame` and converted to times by input stimulus timestamps. The omitted stimuli do not have a `end_frame` by design though their duration is always 250ms. :param data: the behavior data file :param time: the stimulus timestamps indicating when each stimuli is

displayed

Returns df: a pandas dataframe containing the stimuli and omitted stimuli that were displayed with their frame, end_frame, start_time, and duration

```
allensdk.brain_observatory.behavior.stimulus_processing.load_pickle(pstream)
```

```
allensdk.brain_observatory.behavior.stimulus_processing.unpack_change_log(change)
```

allensdk.brain_observatory.behavior.trial_masks module

```
allensdk.brain_observatory.behavior.trial_masks.contingent_trials(trials)
```

GO & CATCH trials only

Parameters

trials [pandas DataFrame] dataframe of trials

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

```
allensdk.brain_observatory.behavior.trial_masks.reward_rate(trials, thresh=2.0)
```

masks trials where the reward rate (per minute) is below some threshold.

This de facto omits trials in which the animal was not licking for extended periods or periods when they were licking indiscriminantly.

Parameters

trials [pandas DataFrame] dataframe of trials

thresh [float, optional] threshold under which trials will not be included, default: 2.0

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

```
allensdk.brain_observatory.behavior.trial_masks.trial_types(trials, trial_types)
```

only include trials of certain trial types

Parameters

trials [pandas DataFrame] dataframe of trials

trial_types [list or other iterator]

Returns

mask [pandas Series of booleans, indexed to trials DataFrame]

allensdk.brain_observatory.behavior.trials_processing module

```

allensdk.brain_observatory.behavior.trials_processing.calculate_reward_rate (response_latency=None, start-  
time=None, win-  
dow=0.75, trial_window=25, ini-  
tial_trials=10)

allensdk.brain_observatory.behavior.trials_processing.categorize_one_trial (tr)

allensdk.brain_observatory.behavior.trials_processing.colormap (trial_type,  
response_type)

allensdk.brain_observatory.behavior.trials_processing.create_extended_trials (trials=None, meta-  
data=None, time=None, licks=None)

allensdk.brain_observatory.behavior.trials_processing.data_to_licks (data,  
time)

allensdk.brain_observatory.behavior.trials_processing.data_to_metadata (data,  
time)

allensdk.brain_observatory.behavior.trials_processing.find_licks (reward_times,  
licks, win-  
dow=3.5)

allensdk.brain_observatory.behavior.trials_processing.get_change_time_frame_response_latency (data)

allensdk.brain_observatory.behavior.trials_processing.get_even_sampling (data)
    Get status of even_sampling

```

Parameters

data: Mapping foraging2 experiment output data

Returns

bool: True if even_sampling is enabled

```

allensdk.brain_observatory.behavior.trials_processing.get_extended_trials (data,  
time=None)

allensdk.brain_observatory.behavior.trials_processing.get_image_info_from_trial (trial_log,  
ti)

allensdk.brain_observatory.behavior.trials_processing.get_mouse_id (exp_data)

allensdk.brain_observatory.behavior.trials_processing.get_ori_info_from_trial (trial_log,  
ti)

allensdk.brain_observatory.behavior.trials_processing.get_params (exp_data)

allensdk.brain_observatory.behavior.trials_processing.get_response_latency (change_event,  
trial)

allensdk.brain_observatory.behavior.trials_processing.get_response_type (trials)

```

```
allensdk.brain_observatory.behavior.trials_processing.get_stimulus_attr_changes(stim_dict,  
                                                                              change_frame,  
                                                                              first_frame,  
                                                                              last_frame)
```

Notes

- assumes only two stimuli are ever shown
- converts attr_names to lowercase
- gets the net attr changes from the start of a trial to the end of a trial

```
allensdk.brain_observatory.behavior.trials_processing.get_time(exp_data)
```

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_image_names(trial,  
                                                                              stim-  
                                                                              uli)  
                                                                              →  
                                                                              Dict[str,  
                                                                              str]
```

Gets the name of the stimulus presented at the beginning of the trial and what is it changed to at the end of the trial. Parameters ——— trial: A trial in a behavior ophys session stimuli: The stimuli presentation log for the behavior session

Returns

A dictionary indicating the starting_stimulus and what the stimulus is changed to.

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_lick_times(lick_times,  
                                                                            start_time,  
                                                                            stop_time)
```

extract lick times in time range

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_reward_time(rebased_reward_time,  
                                                                              start_time,  
                                                                              stop_time)
```

extract reward times in time range

```
allensdk.brain_observatory.behavior.trials_processing.get_trial_timing(event_dict:
                                                                    dict,
                                                                    stim-
                                                                    u-
                                                                    lus_presentations_df:
                                                                    pandas.core.frame.DataFrame,
                                                                    licks:
                                                                    List[float],
                                                                    go:
                                                                    bool,
                                                                    catch:
                                                                    bool,
                                                                    auto_rewarded:
                                                                    bool,
                                                                    hit:
                                                                    bool,
                                                                    false_alarm:
                                                                    bool,
                                                                    aborted:
                                                                    bool)
```

Extract a dictionary of trial timing data. See `trial_data_from_log` for a description of the trial types.

Parameters

- event_dict: dict** Dictionary of trial events in the well-known *pkl* file
- stimulus_presentations_df: pd.DataFrame** pandas dataframe of stimulus presentations, from the `get_stimulus_presentations` response for the `BehaviorOphysSession.api`.
- licks: List[float]** list of lick timestamps, from the `get_licks` response for the `BehaviorOphysSession.api`.
- go: bool** True if “go” trial, False otherwise. Mutually exclusive with *catch*.
- catch: bool** True if “catch” trial, False otherwise. Mutually exclusive with *go*.
- auto_rewarded: bool** True if “auto_rewarded” trial, False otherwise.
- hit: bool** True if “hit” trial, False otherwise
- false_alarm: bool** True if “false_alarm” trial, False otherwise
- aborted: bool** True if “aborted” trial, False otherwise

Returns

dict

- start_time: float** The time the trial started (in seconds elapsed from recording start)
- stop_time: float** The time the trial ended (in seconds elapsed from recording start)
- trial_length: float** Duration of the trial in seconds
- response_time: float** The response time, for non-aborted trials. This is equal to the first lick in the trial. For aborted trials or trials without licks, *response_time* is NaN.
- change_frame: int** The frame number that the stimulus changed
- change_time: float** The time in seconds that the stimulus changed

response_latency: float or None The time in seconds between the stimulus change and the animal's lick response, if the trial is a "go", "catch", or "auto_rewarded" type. If the animal did not respond, return *float("inf")*. In all other cases, return None.

Notes

The following parameters are mutually exclusive (exactly one can be true):

hit, miss, false_alarm, aborted, auto_rewarded

```
allensdk.brain_observatory.behavior.trials_processing.get_trials(data,
                                                                licks_df,
                                                                rewards_df,
                                                                stimu-
                                                                lus_presentations_df,
                                                                rebase)
allensdk.brain_observatory.behavior.trials_processing.get_trials_v0(data,
                                                                    time)
allensdk.brain_observatory.behavior.trials_processing.local_time(iso_timestamp,
                                                                time-
                                                                zone=None)
allensdk.brain_observatory.behavior.trials_processing.resolve_initial_image(stimuli,
                                                                            start_frame)
```

Attempts to resolve the initial image for a given start_frame for a trial

Parameters

stimuli: Mapping foraging2 shape stimuli mapping

start_frame: int start frame of the trial

Returns

initial_image_category_name: str stimulus category of initial image

initial_image_group: str group name of the initial image

initial_image_name: str name of the initial image

```
allensdk.brain_observatory.behavior.trials_processing.trial_data_from_log(trial)
Infer trial logic from trial log. Returns a dictionary.
```

- reward volume: volume of water delivered on the trial, in mL

Each of the following values is boolean:

Trial category values are mutually exclusive * go: trial was a go trial (trial with a stimulus change) * catch: trial was a catch trial (trial with a sham stimulus change)

stimulus_change/sham_change are mutually exclusive * stimulus_change: did the stimulus change (True on 'go' trials) * sham_change: stimulus did not change, but response was evaluated (True on 'catch' trials)

Each trial can be one (and only one) of the following: * hit (stimulus changed, animal responded in response window) * miss (stimulus changed, animal did not respond in response window) * false_alarm (stimulus did not change, animal responded in response window) * correct_reject (stimulus did not change, animal did not respond in response window) * aborted (animal responded before change time) * auto_rewarded (reward was automatically delivered following the change. This will bias the animals choice and should not be categorized as hit/miss)

`allensdk.brain_observatory.behavior.trials_processing.validate_trial_condition_exclusivity`

ensure that only one of N possible mutually exclusive trial conditions is True

`allensdk.brain_observatory.behavior.validation` module

Module contents

`allensdk.brain_observatory.ecephys` package

Subpackages

`allensdk.brain_observatory.ecephys.align_timestamps` package

Submodules

`allensdk.brain_observatory.ecephys.align_timestamps.barcode` module

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.extract_barcodes_from_times` (`on_`

`off_`
`in_`
`ter_`
`bar`
`bar`
`cod`
`nbin`

Read barcodes from timestamped rising and falling edges.

Parameters

- on_times** [numpy.ndarray] Timestamps of rising edges on the barcode line
- off_times** [numpy.ndarray] Timestamps of falling edges on the barcode line
- inter_barcode_interval** [numeric, optional] Minimum duration of time between barcodes.
- bar_duration** [numeric, optional] A value slightly shorter than the expected duration of each bar
- barcode_duration_ceiling** [numeric, optional] The maximum duration of a single barcode
- nbits** [int, optional] The bit-depth of each barcode

Returns

- barcode_start_times** [list of numeric] For each detected barcode, the time at which that barcode started
- barcodes** [list of int] For each detected barcode, the value of that barcode as an integer.

Notes

ignores first code in prod (ok, but not intended) ignores first on pulse (intended - this is needed to identify that a barcode is starting)

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.find_matching_index` (*master_barcode*, *probe_barcode*, *alignment_type*=*'start'* or *'end'*)

Given a set of barcodes for the master clock and the probe clock, find the indices of a matching set, either starting from the beginning or the end of the list.

Parameters

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe line. One per barcode

alignment_type [string] *'start'* or *'end'*

Returns

master_barcode_index [int] matching index for master barcodes (None if not found)

probe_barcode_index [int] matching index for probe barcodes (None if not found)

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.get_probe_time_offset` (*master_times*, *master_barcodes*, *probe_times*, *probe_barcode_index*, *acq_start_index*, *local_probe_rate*)

Time offset between master clock and recording probes. For converting probe time to master clock.

Parameters

master_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

acq_start_index [int] sample index of probe acquisition start time

local_probe_rate [float] the probe's apparent sampling rate

Returns

total_time_shift [float] Time at which the probe started acquisition, assessed on the master clock. If < 0, the probe started earlier than the master line.

probe_rate [float] The probe's sampling rate, assessed on the master clock

master_endpoints [iterable] Defines the start and end times of the sync interval on the master clock

`allensdk.brain_observatory.ecephys.align_timestamps.barcode.linear_transform_from_intervals`

Find a scale and translation which aligns two 1d segments

Parameters

master [iterable] Pair of floats defining the master interval. Order is [start, end].

probe [iterable] Pair of floats defining the probe interval. Order is [start, end].

Returns

scale [float] Scale factor. If > 1.0 , the probe clock is running fast compared to the master clock. If < 1.0 , the probe clock is running slow.

translation [float] If > 0 , the probe clock started before the master clock. If < 0 , after.

Notes

solves $(\text{master} + \text{translation}) * \text{scale} = \text{probe}$

for scale and translation

```
allensdk.brain_observatory.ecephys.align_timestamps.barcode.match_barcodes(master_times,
                                                                           mas-
                                                                           ter_barcodes,
                                                                           probe_times,
                                                                           probe_barcodes)
```

Given sequences of barcode values and (local) times on a probe line and a master line, find the time points on each clock corresponding to the first and last shared barcode.

If there's only one probe barcode, only the first matching timepoint is returned.

Parameters

master_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

Returns

probe_interval [np.ndarray] Start and end times of the matched interval according to the probe_clock.

master_interval [np.ndarray] Start and end times of the matched interval according to the master clock

allensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset module

class allensdk.brain_observatory.ecephys.align_timestamps.barcode_sync_dataset.**BarcodeSyncDataset**

Bases: `allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset`

barcode_line

Obtain the index of the barcode line for this dataset.

extract_barcodes (*self*, ***barcode_kwargs*)

Read barcodes and their times from this dataset's barcode line.

Parameters

****barcode_kwargs** : Will be passed to `.barcode.extract_barcodes_from_times`

Returns

times [np.ndarray] The start times of each detected barcode.

codes [np.ndarray] The values of each detected barcode

get_barcode_table (*self*, ***barcode_kwargs*)

A convenience method for getting barcode times and codes in a dictionary.

Notes

This method is deprecated!

allensdk.brain_observatory.ecephys.align_timestamps.channel_states module

allensdk.brain_observatory.ecephys.align_timestamps.channel_states.**extract_barcodes_from_st**

Obtain barcodes from timestamped rising/falling edges.

Parameters

channel_states [numpy.ndarray] Rising and falling edges, denoted 1 and -1

timestamps [numpy.ndarray] Sample index of each event.

sampling_rate [numeric] Samples / second

****barcode_kwargs** : Additional parameters describing the barcodes.

allensdk.brain_observatory.ecephys.align_timestamps.channel_states.**extract_splits_from_stat**

Obtain barcodes from timestamped rising/falling edges.

Parameters

channel_states [numpy.ndarray] Rising and falling edges, denoted 1 and -1

timestamps [numpy.ndarray] Sample index of each event.

sampling_rate [numeric] Samples / second

****barcode_kwargs** : Additional parameters describing the barcodes.

allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer module

class allensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer.**ProbeSynchronizer**

Bases: object

classmethod compute (*master_barcode_times*, *master_barcodes*, *probe_barcode_times*, *probe_barcodes*, *min_time*, *max_time*, *probe_start_index*, *local_probe_sampling_rate*)

Compute a transform from probe samples to master times by aligning barcodes.

Parameters

master_barcode_times [np.ndarray] start times of barcodes (according to the master clock) on the master line. One per barcode.

master_barcodes [np.ndarray] barcode values on the master line. One per barcode

probe_barcode_times [np.ndarray] start times (according to the probe clock) of barcodes on the probe line. One per barcode

probe_barcodes [np.ndarray] barcode values on the probe_line. One per barcode

min_time [Float] time (in seconds) of first barcode to align

max_time [Float] time (in seconds) of last barcode to align

probe_start_index [int] sample index of probe acquisition start time

local_probe_sampling_rate [float] the probe's apparent sampling rate

Returns

ProbeSynchronizer : When called, applies the transform computed here to samples on the probe clock.

sampling_rate_scale

The ratio of the probe's sampling rate assessed on the global clock to the probe's locally assessed sampling rate.

Module contents

allensdk.brain_observatory.ecephys.copy_utility package

Module contents

allensdk.brain_observatory.ecephys.current_source_density package

Module contents

allensdk.brain_observatory.ecephys.ecephys_project_api package

Submodules

`allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api` module

class `allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi`

Bases: `object`

get_channels (*self*, *channel_ids*: `Union[~ArrayLike, NoneType]` = `None`, *probe_ids*: `Union[~ArrayLike, NoneType]` = `None`, *session_ids*: `Union[~ArrayLike, NoneType]` = `None`, *published_at*: `Union[str, NoneType]` = `None`)

get_isi_experiments (*self*, **args*, ***kwargs*)

get_natural_movie_template (*self*, *number*) → `Iterable`

get_natural_scene_template (*self*, *number*) → `Iterable`

get_probe_lfp_data (*self*, *probe_id*: `int`) → `Iterable`

get_probes (*self*, *probe_ids*: `Union[~ArrayLike, NoneType]` = `None`, *session_ids*: `Union[~ArrayLike, NoneType]` = `None`, *published_at*: `Union[str, NoneType]` = `None`)

get_session_data (*self*, *session_id*: `int`) → `Iterable`

get_sessions (*self*, *session_ids*: `Union[~ArrayLike, NoneType]` = `None`, *published_at*: `Union[str, NoneType]` = `None`)

get_unit_analysis_metrics (*self*, *unit_ids*: `Union[~ArrayLike, NoneType]` = `None`, *ecephys_session_ids*: `Union[~ArrayLike, NoneType]` = `None`, *session_types*: `Union[~ArrayLike, NoneType]` = `None`) → `pandas.core.frame.DataFrame`

get_units (*self*, *unit_ids*: `Union[~ArrayLike, NoneType]` = `None`, *channel_ids*: `Union[~ArrayLike, NoneType]` = `None`, *probe_ids*: `Union[~ArrayLike, NoneType]` = `None`, *session_ids*: `Union[~ArrayLike, NoneType]` = `None`, *published_at*: `Union[str, NoneType]` = `None`)

`allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_fixed_api` module

`allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_lims_api` module

`allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_warehouse_api` module

`allensdk.brain_observatory.ecephys.ecephys_project_api.http_engine` module

`allensdk.brain_observatory.ecephys.ecephys_project_api.rma_engine` module

`allensdk.brain_observatory.ecephys.ecephys_project_api.utilities` module

Module contents

`allensdk.brain_observatory.ecephys.ecephys_session_api` package

Submodules

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_nwb1_session_api` module

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_nwb_session_api` module

`allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api` module

class `allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi`

Bases: `object`

get_channels (*self*) → `pandas.core.frame.DataFrame`

get_ecephys_session_id (*self*) → `int`

get_invalid_times (*self*) → `pandas.core.frame.DataFrame`

get_lfp (*self*, *probe_id: int*) → `xarray.core.dataarray.DataArray`

get_mean_waveforms (*self*) → `Dict[int, numpy.ndarray]`

get_metadata (*self*)

get_optogenetic_stimulation (*self*) → `pandas.core.frame.DataFrame`

get_probes (*self*) → `pandas.core.frame.DataFrame`

get_pupil_data (*self*) → `Union[pandas.core.frame.DataFrame, NoneType]`

get_rig_metadata (*self*) → `Union[dict, NoneType]`

get_running_speed (*self*) → `allensdk.brain_observatory.running_speed.RunningSpeed`

get_screen_gaze_data (*self*, *include_filtered_data=False*) → `Union[pandas.core.frame.DataFrame, NoneType]`

get_session_start_time (*self*) → `datetime.datetime`

get_spike_amplitudes (*self*) → `Dict[int, numpy.ndarray]`

get_spike_times (*self*) → `Dict[int, numpy.ndarray]`

get_stimulus_presentations (*self*) → `pandas.core.frame.DataFrame`

get_units (*self*) → `pandas.core.frame.DataFrame`

session_na = `-1`

test (*self*) → `bool`

Module contents

`allensdk.brain_observatory.ecephys.file_io` package

Submodules

allensdk.brain_observatory.ecephys.file_io.continuous_file module

```
class allensdk.brain_observatory.ecephys.file_io.continuous_file.ContinuousFile(data_path,
                                     times-
                                     tamps_path,
                                     to-
                                     tal_num_chann
                                     dtype=<class
                                     'numpy.int16'>
```

Bases: object

Represents a continuous (.dat) file, and its associated timestamps

```
get_lfp_channel_order(self)
    Returns the channel ordering for LFP data extracted from NPX files.

    None
```

```
load(self, memmap=False, memmap_thresh=10000000000.0)
    Reads lfp data and timestamps from the filesystem
```

memmap [bool, optional] If True, the returned data array will be a memory map of the file on disk. Default is True.

memmap_thresh [float, optional] Files above this size in bytes will be memory-mapped, regardless of memmap setting

allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset module

```
class allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset
    Bases: allensdk.brain_observatory.sync_dataset.Dataset
```

```
extract_frame_times(self, strategy, photodiode_cycle=60, frame_keys=('frames', 'stim_vsync'),
                    photodiode_keys=('photodiode', 'stim_photodiode'))
```

```
extract_frame_times_from_photodiode(self, photodiode_cycle=60, frame_keys=('frames',
                                'stim_vsync'), photodiode_keys=('photodiode',
                                'stim_photodiode'))
```

```
extract_frame_times_from_vsyncs(self, photodiode_cycle=60, frame_keys=('frames',
                                'stim_vsync'), photodiode_keys=('photodiode',
                                'stim_photodiode'))
```

```
extract_led_times(self, keys=('LED_sync', 'opto_trial'), fallback_line=18)
```

```
classmethod factory(path)
    Build a new SyncDataset.
```

Parameters

path [str] Filesystem path to the h5 file containing sync information to be loaded.

sample_frequency

allensdk.brain_observatory.ecephys.file_io.stim_file module

```
class allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile(data,
                                                                                       **kwargs)
    Bases: object
```


angular_wheel_rotation

Extract the total rotation of the running wheel on each frame.

angular_wheel_velocity

Extract the mean angular velocity of the running wheel (degrees / s) for each frame.

classmethod factory (*path*, ***kwargs*)**frames_per_second**

Framerate of stimulus presentation

pre_blank_sec

Time (s) before initial stimulus presentation

stimuli

List of dictionaries containing information about individual stimuli

vin**vsig**

Running speed signal voltage

Module contents**allensdk.brain_observatory.ecephys.lfp_subsampling package****Submodules****allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling module**

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.remove_lfp_noise` (*lfp*, *surface_channel*, *channel_numbers*, *channel_max=384*, *channel_limit=380*)

Subtract mean of channels out of brain to remove noise

lfp [numpy.ndarray] 2D array of LFP values (time x channels)

surface_channel [int] Surface channel (relative to original probe)

channel_numbers [numpy.ndarray] Channel numbers in 'lfp' array (relative to original probe)

Returns:

lfp_noise_removed [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.remove_lfp_offset` (*lfp*, *sampling_frequency*, *cut_off_frequency*, *filter_order*)

High-pass filters LFP data to remove offset

lfp [numpy.ndarray] 2D array of LFP values (time x channels)

sampling_frequency [float] Sampling frequency in Hz

cutoff_frequency [float] Cutoff frequency for highpass filter

filter_order [int] Butterworth filter order

Returns:

lfp_filtered [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.select_channels` (*total_channels*,
sur-
face_channel,
sur-
face_padding,
start_channel_offset,
chan-
nel_stride,
chan-
nel_order,
noisy_channels=
dtype=float64),
re-
move_noisy_channels,
ref-
er-
ence_channels=
dtype=float64),
re-
move_references)

Selects a subset of channels for spatial downsampling

total_channels [int] Number of channels in the original data file

surface_channel [int] Index of channel at brain surface

surface_padding [int] Number of channels above surface to save

start_channel_offset [int] First channel to save

channel_stride [int] Number of channels to skip in output

channel_order [np.ndarray] Actual order of LFP channels (needed to account for the bug in NPX extraction)

noisy_channels [numpy.ndarray] Array indicating noisy channels

remove_noisy_channels [bool] Flag to remove noisy channels

reference_channels [numpy.ndarray] Array indicating reference channels

remove_references [bool] Flag to remove reference channels

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.subsample_lfp` (*lfp_raw*,
se-
lected_channels,
sub-
sam-
pling_factor)

Subsamples LFP data

lfp_raw [numpy.ndarray] 2D array of LFP values (time x channels)

selected_channels [numpy.ndarray] Indices of channels to select (spatial subsampling)

downsampling_factor [int] Factor by which to subsample in time

Returns:

lfp_subsampled [numpy.ndarray] New 2D array of LFP values

`allensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.subsample_timestamps` (*timestamp*
sub-
sam-
pling_fac

Subsamples an array of timestamps

timestamps [numpy.ndarray] 1D array of timestamp values

downsampling_factor [int] Factor by which to subsample the timestamps

Returns:

timestamps_sub [numpy.ndarray] New 1D array of timestamps

Module contents

`allensdk.brain_observatory.ecephys.nwb` package

Submodules

`allensdk.brain_observatory.ecephys.nwb.ecephys_nwb_extension_builder` module

Module contents

`allensdk.brain_observatory.ecephys.optotagging_table` package

Module contents

`allensdk.brain_observatory.ecephys.stimulus_analysis` package

Submodules

`allensdk.brain_observatory.ecephys.stimulus_analysis.dot_motion` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.drifting_gratings` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.flashes` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.natural_movies` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.natural_scenes` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.receptive_field_mapping` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.static_gratings` module

`allensdk.brain_observatory.ecephys.stimulus_analysis.stimulus_analysis` module

Module contents

`allensdk.brain_observatory.ecephys.stimulus_table` package

Subpackages

`allensdk.brain_observatory.ecephys.stimulus_table.visualization` package

Submodules

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks` module

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks`.**build_colormap**

`allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks`.**get_blocks** (*table*)

```
allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.main(table_csv_path)
allensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks.plot_blocks(block_id, color, norm)
```

Module contents

Submodules

allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes module

Created on Fri Dec 16 15:11:23 2016

@author: Xiaoxuan Jia

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.apply_display_sequence(sweep_frames_table, frame_display_sequence)

Adjust raw sweep frames for a stimulus based on the display sequence for that stimulus.
```

Parameters

sweep_frames_table [pd.DataFrame] Each row is a sweep. Has two columns, 'start' and 'end', which describe (in frames) when that sweep began and ended.

frame_display_sequence [np.ndarray] 2D array. Rows are display intervals. The 0th column is the start frame of that interval, the 1st the end frame.

Returns

sweep_frames_table [pd.DataFrame] As above, but start and end frames have been adjusted based on the display sequence.

Notes

The frame values in the raw sweep_frames_table are given in 0-indexed offsets from the start of display for this stimulus. This domain only takes into account frames which are part of a display interval for that stimulus, so the frame ids need to be adjusted to lie on the global frame sequence.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.apply_frame_times(stimulus_table, frame_times)

Converts sweep times from frames to seconds.
```

Parameters

stimulus_table [pd.DataFrame] Rows are sweeps. Columns are stimulus parameters as well as start and end frames for each sweep.

frame_times [numpy.ndarray] Gives the time in seconds at which each frame (indices) began.

frames_per_second [numeric, optional] If provided, and `extra_frame_time` is `True`, will be used to calculate the `extra_frame_time`.

extra_frame_time [float, optional] If provided, an additional frame time will be appended. The time will be incremented by `extra_frame_time` from the previous last frame time, to denote the time at which the last frame ended. If `False`, no extra time will be appended. If `None` (default), the increment will be 1.0/fps.

map_columns [tuple of str, optional] Which columns to replace with times. Defaults to 'Start' and 'End'

Returns

stimulus_table [pd.DataFrame] As above, but with `map_columns` values converted to seconds from frames.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.assign_sweep_values(stim_
swee
on='
drop:
tmp_
```

Left joins a stimulus table to a sweep table in order to associate epochs in time with stimulus characteristics.

Parameters

stim_table [pd.DataFrame] Each row is a stimulus epoch, with start and end times and a foreign key onto a particular sweep.

sweep_table [pd.DataFrame] Each row is a sweep. Should have columns in common with the `stim_table` - the resulting table will use values from the `sweep_table`.

on [str, optional] Column on which to join.

drop [bool, optional] If `True` (default), the join column (argument `on`) will be dropped from the output.

tmp_suffix [str, optional] Will be used to identify overlapping columns. Should not appear in the name of any column in either dataframe.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.build_stimuluswise_table
```

Construct a table of sweeps, including their times on the experiment-global clock and the values of each relevant parameter.

Parameters

stimulus [dict] Describes presentation of a stimulus on a particular experiment. Has a number of fields, of which we are using:

stim_path [str] windows file path to the stimulus data

sweep_frames [list of lists] rows are sweeps, columns are start and end frames of that sweep (in the stimulus-specific frame domain). C-order.

sweep_order [list of int] indices are frames, values are the sweep on that frame

display_sequence [list of list]

rows are intervals in which the stimulus was displayed. Columns are start and end times (s, global) of the display. C-order.

dimnames [list of str] Names of parameters for this stimulus (such as “Contrast”)

sweep_table [list of tuple] Each element is a tuple of parameter values (1 per dim-name) describing a single sweep.

seconds_to_frames [function] Converts experiment seconds to frames

start_key [str, optional] key to use for start frame indices. Defaults to ‘Start’

end_key [str, optional] key to use for end frame indices. Defaults to ‘End’

name_key [str, optional] key to use for stimulus name annotations. Defaults to ‘stimulus_name’

block_key [str, optional] key to use for the 0-index position of this stimulus block

get_stimulus_name [function | dict -> str, optional] extracts stimulus name from the stimulus dictionary. Default is read_stimulus_name_from_path

Returns

list of pandas.DataFrame : Each table corresponds to an entry in the display sequence. Rows are sweeps, columns are stimulus parameter values as well as “Start” and “End”.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.create_stim_table (stimuli,
stim-
u-
lus_table
spon-
ta-
neous_a
sort_key
block_k
in-
dex_key
```

Build a full stimulus table

Parameters

stimuli [list of dict] Each element is a stimulus dictionary, as provided by the stim.pkl file.

stimulus_tabler [function] A function which takes a single stimulus dictionary as its argument and returns a stimulus table dataframe.

spontaneous_activity_tabler [function] A function which takes a list of stimulus tables as arguments and returns a list of 0 or more tables describing spontaneous activity sweeps.

sort_key [str, optional] Sort the final stimulus table in ascending order by this key. Defaults to ‘Start’.

Returns

stim_table_full [pandas.DataFrame] Each row is a sweep. Has columns describing (in frames) the start and end times of each sweep. Other columns describe the values of stimulus parameters on those sweeps.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.make_spontaneous_activity_table
```

Fills in frame gaps in a set of stimulus tables. Suitable for use as the `spontaneous_activity_table` in `create_stim_table`.

Parameters

stimulus_tables [list of pd.DataFrame] Input tables - should have `start_key` and `end_key` columns.

start_key [str, optional] Column name for the start of a sweep. Defaults to 'Start'.

end_key [str, optional] Column name for the end of a sweep. Defaults to 'End'.

duration_threshold [numeric or None] If not None (default is 0), remove spontaneous activity sweeps whose duration is less than this threshold.

Returns

list : Either empty, or contains a single pd.DataFrame. The rows of the dataframe are spontaneous activity sweeps.

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.read_stimulus_name_from_filename
```

Obtains a human-readable stimulus name by looking at the filename of the 'stim_path' item.

Parameters

stimulus [dict] must contain a 'stim_path' item.

Returns

str : name of stimulus

```
allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.split_column(table, column, new_columns, drop_old=True)
```

Divides a dataframe column into multiple columns.

Parameters

table [pandas.DataFrame] Columns will be drawn from and assigned to this dataframe. This dataframe will NOT be modified inplace.

column [str] This column will be split.

new_columns [dict, mapping strings to functions] Each key will be the name of a new column, while its value (a function) will be used to build the new column's values. The functions should map from a single value of the original column to a single value of the new column.

drop_old [bool, optional] If True, the original column will be dropped from the table.

Returns

table [pd.DataFrame] The modified table

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities module

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**add_number_to_shuffled_r**

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**collapse_columns**(*table*)
merge, where possible, columns that describe the same parameter. This is pretty conservative - it only matches columns by capitalization and it only overrides nans.

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**drop_empty_columns**(*table*)
Remove from the stimulus table columns whose values are all nan

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**map_column_names**(*table*,
name_ma
ig-
nore_case

allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.**map_stimulus_names**(*table*,
name_
stim_c

Applies a mapping to the stimulus names in a stimulus table

Parameters

table [pd.DataFrame] the input stimulus table

name_map [dict, optional] rename the stimuli according to this mapping

stim_colname: str, optional look in this column for stimulus names

```
allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities.standardize_movie_numbers
```

Natural movie stimuli in visual coding are numbered using words, like “natural_movie_two” rather than “natural_movie_2”. This function ensures that all of the natural movie stimuli in an experiment are named by that convention.

Parameters

table [pd.DataFrame] the incoming stimulus table

movie_re [re.Pattern, optional] regex that matches movie stimulus names

numeral_re [re.Pattern, optional] regex that extracts movie numbers from stimulus names

digit_names [dict, optional] map from numerals to english words

stim_colname [str, optional] the name of the dataframe column that contains stimulus names

Returns

table [pd.DataFrame] the stimulus table with movie numerals having been mapped to english words

allensdk.brain_observatory.ecephys.stimulus_table.output_validation module

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_epoch_duration
```

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_epoch_order (table, threshold, error_message)
```

```
allensdk.brain_observatory.ecephys.stimulus_table.output_validation.validate_max_spontaneous
```

allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction module

```
allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.extract_com
```

Parameters which are not set as sweep_params in the stimulus script (usually because they are not varied during the course of the session) are not output in an easily machine-readable format. This function attempts to recover them by parsing the string repr of the stimulus.

Parameters

stim_repr [str]

The repr of the camstim stimulus object. Served up per-stimulus in the stim pickle.

repr_params_re [re.Pattern] Extracts attributes as “=”-seperated strings

array_re [re.Pattern] Extracts list reprs from numpy array reprs.

Returns

repr_params [dict] dictionary of parameter keys and values extracted from the stim repr.
Where possible, the values are converted to native Python types.

```
allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.extract_st
```

`allensdk.brain_observatory.ecephys.stimulus_table.stimulus_parameter_extraction.parse_stim`

Read the string representation of a psychopy stimulus and extract stimulus parameters.

Parameters

stim_repr [str]

drop_params [tuple]

repr_params_re [re.Pattern]

array_re [re.Pattern]

Returns

dict : maps extracted parameter names to values

Module contents

`allensdk.brain_observatory.ecephys.visualization` package

Module contents

`allensdk.brain_observatory.ecephys.visualization.plot_mean_waveforms` (*mean_waveforms*,
unit_ids,
peak_channels)

Utility for plotting mean waveforms on each unit's peak channel

Parameters

mean_waveforms [dictionary] Maps unit ids to channelwise average spike waveforms for those units

unit_ids [array-like] unique integer identifiers for units to be included

```
allensdk.brain_observatory.ecephys.visualization.plot_spike_counts(data_array,
                                                                    time_coords,
                                                                    cbar_label,
                                                                    title, xlabel=
                                                                    'time
                                                                    relative to
                                                                    stimulus
                                                                    onset
                                                                    (s)', ylabel=
                                                                    'unit',
                                                                    xtick_step=20)
```

Utility for making a simple spike counts plot.

Parameters

data_array [xarray.DataArray] 2D data array unitwise values per time bin. See EcephysSession.sweepwise_spike_counts

```
allensdk.brain_observatory.ecephys.visualization.raster_plot(spike_times,
                                                                figsize=(8, 8),
                                                                cmap=<matplotlib.colors.ListedColormap
                                                                object at
                                                                0x7fc9c3b5f550>,
                                                                title='spike
                                                                raster', cycle_
                                                                colors=False)
```

allensdk.brain_observatory.ecephys.write_nwb package

Module contents

Submodules

allensdk.brain_observatory.ecephys.ecephys_project_cache module

allensdk.brain_observatory.ecephys.ecephys_session module

allensdk.brain_observatory.ecephys.stimulus_sync module

```
allensdk.brain_observatory.ecephys.stimulus_sync.allocate_by_vsync(vs_diff,
                                                                    index,
                                                                    starts,
                                                                    ends,
                                                                    frame_duration,
                                                                    irregularity, cycle)

allensdk.brain_observatory.ecephys.stimulus_sync.assign_to_last(index,
                                                                    starts, ends,
                                                                    frame_duration,
                                                                    irregularity,
                                                                    cycle)
```

```
allensdk.brain_observatory.ecephys.stimulus_sync.compute_frame_times(photodiode_times,  
                                                                    frame_duration,  
                                                                    num_frames,  
                                                                    cycle,  
                                                                    irregu-  
                                                                    lar_interval_policy=<function  
                                                                    as-  
                                                                    sign_to_last  
                                                                    at  
                                                                    0x7fc98a2d9400>)  
  
allensdk.brain_observatory.ecephys.stimulus_sync.correct_on_off_effects(pd_times)
```

Notes

This cannot (without additional info) determine whether an assymmetric offset is odd-long or even-long.

```
allensdk.brain_observatory.ecephys.stimulus_sync.estimate_frame_duration(pd_times,  
                                                                    cy-  
                                                                    cle=60)  
  
allensdk.brain_observatory.ecephys.stimulus_sync.fix_unexpected_edges(pd_times,  
                                                                    ndevs=10,  
                                                                    cy-  
                                                                    cle=60,  
                                                                    max_frame_offset=4)  
  
allensdk.brain_observatory.ecephys.stimulus_sync.flag_unexpected_edges(pd_times,  
                                                                    ndevs=10)  
  
allensdk.brain_observatory.ecephys.stimulus_sync.trim_border_pulses(pd_times,  
                                                                    vs_times,  
                                                                    frame_interval=0.016666666666666666,  
                                                                    num_frames=5)  
  
allensdk.brain_observatory.ecephys.stimulus_sync.trimmed_stats(data, pc-  
                                                                    tiles=(10, 90))
```

Module contents

```
allensdk.brain_observatory.ecephys.get_unit_filter_value(key, pop=True, re-  
                                                         place_none=True,  
                                                         **source)
```

allensdk.brain_observatory.extract_running_speed package

Module contents

allensdk.brain_observatory.gaze_mapping package

Module contents

allensdk.brain_observatory.nwb package

Submodules

allensdk.brain_observatory.nwb.behavior_ophys_nwb_extension_builder module

allensdk.brain_observatory.nwb.metadata module

```
allensdk.brain_observatory.nwb.metadata.create_pynwb_extension_from_schemas (schema_list,
                                                                                   pre-
                                                                                   fix:
                                                                                   str)
```

```
allensdk.brain_observatory.nwb.metadata.extract_from_schema (schema)
```

```
allensdk.brain_observatory.nwb.metadata.load_pynwb_extension (schema, prefix:
                                                                str)
```

allensdk.brain_observatory.nwb.nwb_api module

```
class allensdk.brain_observatory.nwb.nwb_api.NwbApi (path, **kwargs)
    Bases: object
    classmethod from_nwbfile (nwbfile, **kwargs)
    classmethod from_path (path, **kwargs)
    get_image (self, name, module, image_api=None) → SimpleITK.SimpleITK.Image
    get_invalid_times (self) → pandas.core.frame.DataFrame
    get_running_speed (self) → allensdk.brain_observatory.running_speed.RunningSpeed
    get_stimulus_presentations (self) → pandas.core.frame.DataFrame
    nwbfile
    path
```

allensdk.brain_observatory.nwb.nwb_utils module

```
allensdk.brain_observatory.nwb.nwb_utils.get_column_name (table_cols: list, possible_names: set) → str
```

This function returns a column name, given a table with unknown column names and a set of possible column names which are expected. The table column name returned should be the only name contained in the “expected” possible names. :param table_cols: the table columns to search for the possible name within :param possible_names: the names that could exist within the data columns :return: the first entry of the intersection between the possible names

and the names of the columns of the stimulus table

```
allensdk.brain_observatory.nwb.nwb_utils.set_omitted_stop_time (stimulus_table:
                                                                pandas.core.frame.DataFrame,
                                                                omitted_time_duration:
                                                                float = 0.25) →
                                                                None
```

This function sets the stop time for a row that of a stimuli table that is a omitted stimuli. A omitted stimuli is a stimuli where a mouse is shown only a grey screen and these last for 250 milliseconds. These do not include

a stop_time or end_frame as other stimuli in the stimulus table due to design choices. For these stimuli to be added they must have the stop_time calculated and put into the row as data before writing to NWB. :param stimulus_table: pd.DataFrame that contains the stimuli presented to

an experiment subject

Parameters omitted_time_duration – The duration in seconds of the expected length of the omitted stimuli

Returns

stimulus_table_row: returns the same dictionary as inputted but with an additional entry for stop_time.

allensdk.brain_observatory.nwb.schemas module

```
class allensdk.brain_observatory.nwb.schemas.RunningSpeedPathsSchema (only=None,
                                                                    exclude=(),
                                                                    many=False,
                                                                    context=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    partial=False,
                                                                    unknown=None)

Bases: allensdk.brain_observatory.argschema_utilities.RaisingSchema

opts = <marshmallow.schema.SchemaOpts object>
```

Module contents

```
allensdk.brain_observatory.nwb.add_average_image (nwbfile, average_image, image_api=None)
```

```
allensdk.brain_observatory.nwb.add_cell_specimen_table (nwbfile:
                                                         pynwb.file.NWBFile,
                                                         cell_specimen_table: pandas.core.frame.DataFrame)
```

This function takes the cell specimen table and writes the ROIs contained within. It writes these to a new NWB imaging plane based off the previously supplied metadata Parameters ——— nwbfile: NWBFile

this is the in memory NWBFile currently being written to which ROI data is added

cell_specimen_table: pd.DataFrame this is the DataFrame containing the cells segmented from a ophys experiment, stored in json file and loaded. example: /home/nicholasc/projects/allensdk/allensdk/test/

brain_observatory/behavior/cell_specimen_table_789359614.json

Returns

nwbfile: NWBFile The altered in memory NWBFile object that now has a specimen table

```
allensdk.brain_observatory.nwb.add_corrected_fluorescence_traces (nwbfile, corrected_fluorescence_traces)
```



```
allensdk.brain_observatory.nwb.add_dff_traces(nwbfile, dff_traces, ophys_timestamps)
```

```
allensdk.brain_observatory.nwb.add_eye_gaze_data_interfaces(pynwb_container:
    pynwb.core.NWBContainer,
    pupil_areas: pan-
    das.core.series.Series,
    eye_areas: pan-
    das.core.series.Series,
    screen_coordinates:
    pan-
    das.core.frame.DataFrame,
    screen_coordinates_spherical:
    pan-
    das.core.frame.DataFrame,
    synced_timestamps:
    pan-
    das.core.series.Series)
→
pynwb.core.NWBContainer
```

```
allensdk.brain_observatory.nwb.add_eye_gaze_mapping_data_to_nwbfile(nwbfile:
    pynwb.file.NWBFile,
    eye_gaze_data:
    dict) →
    pynwb.file.NWBFile
```

```
allensdk.brain_observatory.nwb.add_eye_tracking_ellipse_fit_data_to_nwbfile(nwbfile:
    pynwb.file.NWBFile,
    eye_dlc_tracking_data:
    dict,
    synced_timestamps:
    pan-
    das.core.series.Series)
→
pynwb.file.NWBFile
```

```
allensdk.brain_observatory.nwb.add_image(nwbfile, image_data, image_name, module_name,
    module_description, image_api=None)
```

```
allensdk.brain_observatory.nwb.add_invalid_times(nwbfile, epochs)
```

Write invalid times to nwbfile if epochs are not empty Parameters ——— nwbfile: pynwb.NWBFile epochs: list of dicts

records of invalid epochs

Returns

pynwb.NWBFile

```
allensdk.brain_observatory.nwb.add_licks(nwbfile, licks)
```

```
allensdk.brain_observatory.nwb.add_max_projection(nwbfile, max_projection, image_api=None)
```

```
allensdk.brain_observatory.nwb.add_metadata(nwbfile, metadata: dict)
```

```
allensdk.brain_observatory.nwb.add_motion_correction(nwbfile, motion_correction)
```

```
allensdk.brain_observatory.nwb.add_rewards(nwbfile, rewards_df)
```

```
allensdk.brain_observatory.nwb.add_running_data_df_to_nwbfile(nwbfile,      run-
                                                                ning_data_df,
                                                                unit_dict,      in-
                                                                dex_key='timestamps')
```

Adds running speed data to an NWBFile as timeseries in acquisition and processing

Parameters

nwbfile [pynwb.NWBFile] File to which running speeds will be written

running_speed [pandas.DataFrame] Contains 'speed' and 'times', 'v_in', 'vsig', 'dx'

unit [str, optional] SI units of running speed values

Returns

nwbfile [pynwb.NWBFile]

```
allensdk.brain_observatory.nwb.add_running_speed_to_nwbfile(nwbfile,      run-
                                                                ning_speed,
                                                                name='speed',
                                                                unit='cm/s')
```

Adds running speed data to an NWBFile as a timeseries in acquisition

Parameters

nwbfile [pynwb.NWBFile] File to which running speeds will be written

running_speed [RunningSpeed] Contains attributes 'values' and 'timestamps'

name [str, optional] used as name of timeseries object

unit [str, optional] SI units of running speed values

Returns

nwbfile [pynwb.NWBFile]

```
allensdk.brain_observatory.nwb.add_segmentation_mask_image(nwbfile,      segmenta-
                                                                tion_mask_image,
                                                                image_api=None)
```

```
allensdk.brain_observatory.nwb.add_stimulus_index(nwbfile,      stimulus_index,
                                                    nwb_template)
```

```
allensdk.brain_observatory.nwb.add_stimulus_presentations(nwbfile, stimulus_table,
                                                            tag='stimulus_time_interval')
```

Adds a stimulus table (defining stimulus characteristics for each time point in a session) to an nwbfile as TimeIntervals.

Parameters

nwbfile [pynwb.NWBFile]

stimulus_table: pd.DataFrame Each row corresponds to an interval of time. Columns define the interval (start and stop time) and its characteristics. Nans in columns with string data will be replaced with the empty strings. Required columns are:

start_time :: the time at which this interval started stop_time :: the time at which this interval ended

tag [str, optional] Each interval in an nwb file has one or more tags. This string will be applied as a tag to all TimeIntervals created here

Returns

nwbfile [pynwb.NWBFile]

```

allensdk.brain_observatory.nwb.add_stimulus_template(nwbfile, image_data, name)
allensdk.brain_observatory.nwb.add_stimulus_timestamps(nwbfile, stimulus_timestamps, module_name='stimulus')
allensdk.brain_observatory.nwb.add_task_parameters(nwbfile, task_parameters)
allensdk.brain_observatory.nwb.add_trials(nwbfile, trials, description_dict={})
allensdk.brain_observatory.nwb.check_nwbfile_version(nwbfile_path: str, desired_minimum_version: str, warning_msg: str)
allensdk.brain_observatory.nwb.create_eye_gaze_mapping_dataframe(eye_gaze_data: dict) → pandas.core.frame.DataFrame
allensdk.brain_observatory.nwb.create_eye_tracking_nwb_processing_module(eye_dlc_tracking_data: dict, synced_timestamps: pandas.core.series.Series) → pynwb.base.ProcessingModule
allensdk.brain_observatory.nwb.create_gaze_mapping_nwb_processing_modules(eye_gaze_data: dict)
allensdk.brain_observatory.nwb.create_stimulus_presentation_time_interval(name: str, description: str, columns_to_add: Iterable) → pynwb.epoch.TimeInterval
allensdk.brain_observatory.nwb.eye_tracking_data_is_valid(eye_dlc_tracking_data: dict, synced_timestamps: pandas.core.series.Series) → bool
allensdk.brain_observatory.nwb.read_eye_dlc_tracking_ellipses(input_path: pathlib.Path) → dict

```

Reads eye tracking ellipse fit data from an h5 file.

Args: input_path (Path): Path to eye tracking ellipse fit h5 file

Returns:

dict: Loaded h5 data. Each 'params' field contains dataframes with] ellipse fit parameters. Dataframes contain 5 columns each consisting of: "center_x", "center_y", "height", "phi", "width"

`allensdk.brain_observatory.nwb.read_eye_gaze_mappings` (*input_path: pathlib.Path*) → dict

Reads eye gaze mapping data from an h5 file.

Args:

input_path (Path): Path to eye gaze mapping h5 data file produced by 'allensdk.brain_observatory.gaze_mapping' module.

Returns:

dict: Loaded h5 data. *_eye_areas: Area of eye (in pixels²) over time *_pupil_areas: Area of pupil (in pixels²) over time *_screen_coordinates: y, x screen coordinates (in cm) over time *_screen_coordinates_spherical: y, x screen coordinates (in deg) over time synced_frame_timestamps: synced timestamps for video frames (in sec)

`allensdk.brain_observatory.nwb.setup_table_for_epochs` (*table, timeseries, tag*)

`allensdk.brain_observatory.nwb.setup_table_for_invalid_times` (*invalid_epochs*)

Create table with invalid times if invalid_epochs are present

Parameters

invalid_epochs: list of dicts of invalid epoch records

Returns

pd.DataFrame of invalid times if epochs are not empty, otherwise return None

allensdk.brain_observatory.ophys package

Subpackages

allensdk.brain_observatory.ophys.trace_extraction package

Module contents

Module contents

allensdk.brain_observatory.receptive_field_analysis package

Submodules

allensdk.brain_observatory.receptive_field_analysis.chisquarperf module

`allensdk.brain_observatory.receptive_field_analysis.chisquarperf.NLL_to_pvalue` (*NLLs, log_base=10.0*)

`allensdk.brain_observatory.receptive_field_analysis.chisquarperf.build_trial_matrix` (*LSN_template, num_trials, on_off_lum, 0*)

Construct indicator arrays for on/off pixels across trials.

Parameters

LSN_template [np.ndarray] Dimensions are (nTrials, nYPixels, nXPixels). Luminance values per pixel and trial. The size of the first dimension may be larger than the num_trials argument (in which case only the first num_trials slices will be used) but may not be smaller.

num_trials [int] The number of trials (left-justified) to build indicators for.

on_off_luminance [array-like, optional] The zeroth element is the luminance value of a pixel when on, the first when off. Defaults are [255, 0].

Returns

trial_mat [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}, nTrials). Boolean values indicate that a pixel was on/off on a particular trial.

```
allensdk.brain_observatory.receptive_field_analysis.chisquareref.chi_square_binary(events,
LSN_template)
allensdk.brain_observatory.receptive_field_analysis.chisquareref.chi_square_within_mask(exclusion_mask,
events_per_pixel,
trials_per_pixel)
```

Determine if cells respond preferentially to on/off pixels in a mask using a chi2 test.

Parameters

exclusion_mask [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer indicator for INCLUSION (!) of a pixel within the testing region.

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Integer values are response counts by cell to on/off luminance at each pixel.

trials_per_pixel [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer values are counts of trials where a pixel is on/off.

Returns

p_vals [np.ndarray] One-dimensional, of length nCells. Float values are p-values for the hypothesis that a given cell has a receptive field within the exclusion mask.

chi [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Values (float) are squared residual event counts divided by expected event counts.

```
allensdk.brain_observatory.receptive_field_analysis.chisquareref.deinterpolate_RF(rf_map,
x_pts,
y_pts,
deg_per_pnt)
```

Downsample an image

Parameters

rf_map [np.ndarray] Input image

x_pts [np.ndarray] Count of sample points along the first (column) axis

y_pts [np.ndarray] Count of sample points along the zeroth (row) axis

deg_per_pnt [numeric] scale factor

Returns

sampled_yx [np.ndarray] Downsampled image

```
allensdk.brain_observatory.receptive_field_analysis.chisquarperf.get_disc_masks(LSN_template,
                                                                                   ra-
                                                                                   dius=3,
                                                                                   on_luminance=255,
                                                                                   off_luminance=0)
```

Obtain an indicator mask surrounding each pixel. The mask is a square, excluding pixels which are coactive on any trial with the main pixel.

Parameters

LSN_template [np.ndarray] Dimensions are (nTrials, nYPixels, nXPixels). Luminance values per pixel and trial.

radius [int] The base mask will be a box whose sides are $2 * \text{radius} + 1$ in length.

on_luminance [int, optional] The value of the luminance for on trials. Default is 255

off_luminance [int, optional] The value of the luminance for off trials. Default is 0

Returns

masks [np.ndarray] Dimensions are (nYPixels, nXPixels, nYPixels, nXPixels). The first 2 dimensions describe the pixel from which the mask was computed. The last 2 serve as the dimensions of the mask images themselves. Masks are binary arrays of type float, with 1 indicating inside, 0 outside.

```
allensdk.brain_observatory.receptive_field_analysis.chisquarperf.get_events_per_pixel(response_np, trial_matrix)
```

Obtain a matrix linking cellular responses to pixel activity.

Parameters

responses_np [np.ndarray] Dimensions are (nTrials, nCells). Boolean values indicate presence/absence of a response on a given trial.

trial_matrix [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}, nTrials). Boolean values indicate that a pixel was on/off on a particular trial.

Returns

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Values for each cell, pixel, and on/off state are the sum of events for that cell across all trials where the pixel was in the on/off state.

```
allensdk.brain_observatory.receptive_field_analysis.chisquarperf.get_expected_events_by_pixel(response_np, trial_matrix, exclusion_mask)
```

Calculate expected number of events per pixel

Parameters

exclusion_mask [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer indicator for INCLUSION (!) of a pixel within the testing region.

events_per_pixel [np.ndarray] Dimensions are (nCells, nYPixels, nXPixels, {on, off}). Integer values are response counts by cell to on/off luminance at each pixel.

trials_per_pixel [np.ndarray] Dimensions are (nYPixels, nXPixels, {on, off}). Integer values are counts of trials where a pixel is on/off.

Returns

np.ndarray : Dimensions (nCells, nYPixels, nXPixels, {on, off}). Float values are pixelwise counts of events expected if events are evenly distributed in mask across trials.

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.get_peak_significance` (*chi_sq*,
LSN_t,
al-,
pha=0)

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.interpolate_RF` (*rf_map*,
deg_per_pnt)

Upsample an image

Parameters

rf_map [np.ndarray] Input image
deg_per_pnt [numeric] scale factor

Returns

interpolated [np.ndarray] Upsampled image

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.locate_median` (*y*,
x)

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.pvalue_to_NLL` (*p_values*,
max_NLL=10.0)

`allensdk.brain_observatory.receptive_field_analysis.chisquarerf.smooth_STA` (*STA*,
gauss_std=0.75,
to-,
tal_degrees=64)

Smooth an image by convolution with a gaussian kernel

Parameters

STA [np.ndarray] Input image
gauss_std [numeric, optional] Standard deviation of the gaussian kernel. Will be applied to the upsampled image, so units are visual degrees. Default is 0.75
total_degrees [int, optional] Size in visual degrees of the input image along its zeroth (row) axis. Used to set the scale factor for up/downsampling.

Returns

STA_smoothed [np.ndarray] Smoothed image

`allensdk.brain_observatory.receptive_field_analysis.eventdetection` module

`allensdk.brain_observatory.receptive_field_analysis.eventdetection.detect_events` (*data*,
cell_index,
stim-,
u-,
lus,
de-,
bug_plots=Fa)

allensdk.brain_observatory.receptive_field_analysis.fit_parameters module

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.add_to_fit_parameters_d`

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.compute_distance` (*center_on*,
center_off)

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.compute_overlap` (*data_fitted*,
data_fitted)

`allensdk.brain_observatory.receptive_field_analysis.fit_parameters.get_gaussian_fit_single`

allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D module

exception `allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.GaussianFitError`
Bases: `RuntimeError`

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.fitgaussian2D` (*data*)
Fit a 2D gaussian to an image

Parameters

data [np.ndarray] input image

Returns

p2 [list] height row mean column mean row standard deviation column standard deviation
rotation

Notes

see `gaussian2D` for details about output values

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.gaussian2D` (*height*,
center_x,
center_y,
width_x,
width_y,
rotation)

Build a function which evaluates a scaled 2d gaussian pdf

Parameters

height [float] scale factor

center_x [float] first coordinate of mean

center_y [float] second coordinate of mean

width_x [float] standard deviation along x axis

width_y [float] standard deviation along y axis

rotation [float] degrees clockwise by which to rotate the gaussian

Returns

rotgauss: fn parameters are x and y positions (row/column semantics are set by your inputs to this function). Return value is the scaled gaussian pdf evaluated at the argued point.

`allensdk.brain_observatory.receptive_field_analysis.fitgaussian2D.moments2(data)`
Treating input image data as an independent multivariate gaussian, estimate mean and standard deviations

Parameters

data [np.ndarray] 2d numpy array.

Returns

height [float] The maximum observed value in the data

y [float] Mean row index

x [float] Mean column index

width_y [float] The standard deviation along the mean row

width_x [float] The standard deviation along the mean column

None : This function returns an instance of None.

Notes

uses original method from website for finding center

allensdk.brain_observatory.receptive_field_analysis.postprocessing module

`allensdk.brain_observatory.receptive_field_analysis.postprocessing.get_gaussian_fit(rf)`

`allensdk.brain_observatory.receptive_field_analysis.postprocessing.run_postprocessing(data, rf)`

allensdk.brain_observatory.receptive_field_analysis.receptive_field module

`allensdk.brain_observatory.receptive_field_analysis.receptive_field.compute_receptive_field`

`allensdk.brain_observatory.receptive_field_analysis.receptive_field.compute_receptive_field`

```
allensdk.brain_observatory.receptive_field_analysis.receptive_field.events_to_pvalues_no_f
```

```
allensdk.brain_observatory.receptive_field_analysis.receptive_field.get_attribute_dict(rf)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.print_summary(rf)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.read_h5_group(g)
allensdk.brain_observatory.receptive_field_analysis.receptive_field.read_receptive_field_f
allensdk.brain_observatory.receptive_field_analysis.receptive_field.write_receptive_field_t
```

allensdk.brain_observatory.receptive_field_analysis.tools module

```
allensdk.brain_observatory.receptive_field_analysis.tools.dict_generator(indict,
                                                                           pre=None)
allensdk.brain_observatory.receptive_field_analysis.tools.list_of_dicts_to_dict_of_lists(lis
allensdk.brain_observatory.receptive_field_analysis.tools.read_h5_group(g)
```

allensdk.brain_observatory.receptive_field_analysis.utilities module

```
allensdk.brain_observatory.receptive_field_analysis.utilities.convolve(img,
                                                                           sigma=4)
    2D Gaussian convolution
allensdk.brain_observatory.receptive_field_analysis.utilities.get_A(data,
                                                                           stimu-
                                                                           lus)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_A_blur(data,
                                                                           stim-
                                                                           u-
                                                                           lus)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_attribute_dict(rf)
allensdk.brain_observatory.receptive_field_analysis.utilities.get_components(receptive_field_data
allensdk.brain_observatory.receptive_field_analysis.utilities.get_shuffle_matrix(data,
                                                                           event_vector,
                                                                           A,
                                                                           num-
                                                                           ber_of_shuffle
                                                                           re-
                                                                           sponse_detect
```

```
allensdk.brain_observatory.receptive_field_analysis.utilities.get_sparse_noise_epoch_mask_
```

```
allensdk.brain_observatory.receptive_field_analysis.utilities.smooth(x, win-
                                                                    dow_len=11,
                                                                    win-
                                                                    dow='hanning',
                                                                    mode='valid')
```

smooth the data using a window with requested size.

This method is based on the convolution of a scaled window with the signal. The signal is prepared by introducing reflected copies of the signal (with the window size) in both ends so that transient parts are minimized in the beginning and end part of the output signal.

input: x: the input signal window_len: the dimension of the smoothing window; should be an odd integer
window: the type of window from 'flat', 'hanning', 'hamming', 'bartlett', 'blackman'

flat window will produce a moving average smoothing.

output: the smoothed signal

example:

```
t=linspace(-2,2,0.1) x=sin(t)+randn(len(t))*0.1 y=smooth(x)
```

see also:

numpy.hanning, numpy.hamming, numpy.bartlett, numpy.blackman, numpy.convolve scipy.signal.lfilter

TODO: the window parameter could be the window itself if an array instead of a string NOTE: length(output) != length(input), to correct this: return y[(window_len/2-1):-(window_len/2)] instead of just y.

```
allensdk.brain_observatory.receptive_field_analysis.utilities.upsample_image_to_degrees (img
```

allensdk.brain_observatory.receptive_field_analysis.visualization module

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_chi_square_summary (n
```

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_ellipses (gaussian_fit_di
                                                                    ax=None,
                                                                    show=True,
                                                                    close=True,
                                                                    save_file_name
                                                                    color='b')
```

Example Usage: oeid, cell_index, stimulus = 512176430, 12, 'locally_sparse_noise' brain_observatory_cache = BrainObservatoryCache() data_set = brain_observatory_cache.get_ophys_experiment_data(oeid) lsn = LocallySparseNoise(data_set, stimulus) result = compute_receptive_field_with_postprocessing(data_set, cell_index, stimulus, alpha=.05, number_of_shuffles=5000) plot_ellipses(result['off']['gaussian_fit'], color='r')

```
allensdk.brain_observatory.receptive_field_analysis.visualization.plot_fields(on_data,
                                                                              off_data,
                                                                              on_axes,
                                                                              off_axes,
                                                                              cbar_axes=None,
                                                                              clim=None,
                                                                              cmap='magma')

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_gaussian_fit(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=
                                                                              cmap='m

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_mask(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=None,
                                                                              cmap='magma')

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_msr_summary(lsn,
                                                                              cell_index,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=N
                                                                              cmap=Non

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_p_values(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=None,
                                                                              cmap='magma

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_receptive_field_data

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_rts_blur_summary(rf_data,
                                                                              ax_
                                                                              ax_
                                                                              ax_
                                                                              cma

allensdk.brain_observatory.receptive_field_analysis.visualization.plot_rts_summary(rf_data,
                                                                              ax_on,
                                                                              ax_off,
                                                                              ax_cbar=N
                                                                              cmap='ma

allensdk.brain_observatory.receptive_field_analysis.visualization.pvalue_to_NLL(p_values,
                                                                              max_NLL=10.0
```

Module contents

allensdk.brain_observatory.sync_utilities package

Module contents

```
allensdk.brain_observatory.sync_utilities.get_synchronized_frame_times(session_sync_file:
                                                                    path-
                                                                    lib.Path,
                                                                    sync_line_label_keys:
                                                                    Tu-
                                                                    ple[str,
                                                                    ...],
                                                                    trim_after_spike:
                                                                    bool
                                                                    =
                                                                    True)
                                                                    →
                                                                    pan-
                                                                    das.core.series.Series
```

Get experimental frame times from an experiment session sync file.

Parameters

session_sync_file [Path] Path to an ephys session sync file. The sync file contains rising/falling edges from a daq system which indicates when certain events occur (so they can be related to each other).

sync_line_label_keys [Tuple[str, ...]] Line label keys to get times for. See class attributes of `allensdk.brain_observatory.sync_dataset.Dataset` for a listing of possible keys.

trim_after_spike [bool = True] If True, will call `trim_discontiguous_times` on the frame times before returning them, which will detect any spikes in the data and remove all elements for the list which come after the spike.

Returns

pd.Series An array of times when frames for the eye tracking camera were acquired.

```
allensdk.brain_observatory.sync_utilities.trim_discontiguous_times(times,
                                                                    thresh-
                                                                    old=100)
```

allensdk.brain_observatory.visualization package

Module contents

```
allensdk.brain_observatory.visualization.plot_running_speed(timestamps,    val-
                                                                    ues,    start_index=0,
                                                                    stop_index=None,
                                                                    step=1,    yla-
                                                                    bel='running speed
                                                                    (cm/s)', xlabel='time
                                                                    (s)', title=None)
```

Make a simple plot of a running speed trace

Parameters

timestamps [numpy.ndarray] Times at which running speed samples were collected

values [numpy.ndarray] Running speed values (by default: linear cm / s with negative values indicating backwards movement)

Submodules

allensdk.brain_observatory.argschema_utilities module

```
class allensdk.brain_observatory.argschema_utilities.ArgSchemaParserPlus (*args,  
                                                                    **kwargs)  
    Bases: argschema.argschema_parser.ArgSchemaParser
```

```
class allensdk.brain_observatory.argschema_utilities.InputFile (default=<marshmallow.missing>,  
                                                                    attribute=None,  
                                                                    data_key=None,  
                                                                    error=None,  
                                                                    validate=None,  
                                                                    required=False,  
                                                                    al-  
                                                                    low_none=None,  
                                                                    load_only=False,  
                                                                    dump_only=False,  
                                                                    miss-  
                                                                    ing=<marshmallow.missing>,  
                                                                    er-  
                                                                    ror_messages=None,  
                                                                    **metadata)  
    Bases: marshmallow.fields.String
```

A marshmallow String field subclass which deserializes json str fields that represent a desired input path to pathlib.Path. Also performs read access checking.

```
class allensdk.brain_observatory.argschema_utilities.OutputFile (default=<marshmallow.missing>,  
                                                                    at-  
                                                                    tribute=None,  
                                                                    data_key=None,  
                                                                    error=None,  
                                                                    vali-  
                                                                    date=None,  
                                                                    re-  
                                                                    quired=False,  
                                                                    al-  
                                                                    low_none=None,  
                                                                    load_only=False,  
                                                                    dump_only=False,  
                                                                    miss-  
                                                                    ing=<marshmallow.missing>,  
                                                                    er-  
                                                                    ror_messages=None,  
                                                                    **metadata)  
    Bases: marshmallow.fields.String
```

A marshmallow String field subclass which deserializes json str fields that represent a desired output file path to a pathlib.Path. Also performs write access checking.

```

class allensdk.brain_observatory.argschema_utilities.RaisingSchema (only=None,
                                                                    ex-
                                                                    clude=(),
                                                                    many=False,
                                                                    con-
                                                                    text=None,
                                                                    load_only=(),
                                                                    dump_only=(),
                                                                    par-
                                                                    tial=False,
                                                                    un-
                                                                    known=None)

    Bases: argschema.schemas.DefaultSchema

    class Meta
        Bases: object

        unknown = 'raise'

        opts = <marshmallow.schema.SchemaOpts object>

allensdk.brain_observatory.argschema_utilities.check_read_access (path)
allensdk.brain_observatory.argschema_utilities.check_write_access (filepath, al-
                                                                    low_exists=False)
allensdk.brain_observatory.argschema_utilities.check_write_access_dir (dirpath)
allensdk.brain_observatory.argschema_utilities.check_write_access_overwrite (path)
allensdk.brain_observatory.argschema_utilities.optional_lims_inputs (argv, in-
                                                                    put_schema,
                                                                    out-
                                                                    put_schema,
                                                                    lims_input_getter)
allensdk.brain_observatory.argschema_utilities.write_or_print_outputs (data,
                                                                    parser)

```

allensdk.brain_observatory.brain_observatory_exceptions module

```

exception allensdk.brain_observatory.brain_observatory_exceptions.BrainObservatoryAnalysisException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.EpochSeparationException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.MissingStimulusException
    Bases: Exception

exception allensdk.brain_observatory.brain_observatory_exceptions.NoEyeTrackingException
    Bases: Exception

```

allensdk.brain_observatory.brain_observatory_plotting module

```

allensdk.brain_observatory.brain_observatory_plotting.plot_drifting_grating_traces (dg,
                                                                    save_dir)
    saves figures with a Ori X TF grid of mean responses

```

```
allensdk.brain_observatory.brain_observatory_plotting.plot_lsn_traces (lsn,  
                                                                    save_dir,  
                                                                    suf-  
                                                                    fix="")  
allensdk.brain_observatory.brain_observatory_plotting.plot_ns_traces (nsa,  
                                                                    save_dir)  
allensdk.brain_observatory.brain_observatory_plotting.plot_running_a (dg,  
                                                                    nm1,  
                                                                    nm3,  
                                                                    save_dir)  
allensdk.brain_observatory.brain_observatory_plotting.plot_sg_traces (sg,  
                                                                    save_dir)
```

allensdk.brain_observatory.chisquare_categorical module

Created on Wed Jun 5 15:52:22 2019

@author: dan

```
allensdk.brain_observatory.chisquare_categorical.advance_combination (curr_combination,  
                                                                    op-  
                                                                    tions_per_column)  
allensdk.brain_observatory.chisquare_categorical.chisq_from_stim_table (stim_table,  
                                                                    columns,  
                                                                    mean_sweep_events,  
                                                                    num_shuffles=1000,  
                                                                    ver-  
                                                                    bose=False)  
allensdk.brain_observatory.chisquare_categorical.compute_chi (observed,      ex-  
                                                                    pected)  
allensdk.brain_observatory.chisquare_categorical.compute_chi_shuffle (mean_sweep_events,  
                                                                    sweep_categories,  
                                                                    num_shuffles=1000)  
allensdk.brain_observatory.chisquare_categorical.compute_expected (mean_sweep_events,  
                                                                    sweep_conditions)  
allensdk.brain_observatory.chisquare_categorical.compute_observed (mean_sweep_events,  
                                                                    sweep_conditions)  
allensdk.brain_observatory.chisquare_categorical.make_category_dummy (sweep_categories)  
allensdk.brain_observatory.chisquare_categorical.stim_table_to_categories (stim_table,  
                                                                    columns,  
                                                                    ver-  
                                                                    bose=False)
```

allensdk.brain_observatory.circle_plots module

```
class allensdk.brain_observatory.circle_plots.CoronaPlotter (angle_start=270,  
                                                                    plot_scale=1.2,  
                                                                    inner_radius=0.3,  
                                                                    *args, **kwargs)
```

Bases: `allensdk.brain_observatory.circle_plots.PolarPlotter`


```

infer_dims (self, category_data)

plot (self, category_data, data=None, clim=None, cmap=<matplotlib.colors.LinearSegmentedColormap
      object at 0x7fc9a4777208>)

set_dims (self, categories)

show_arrow (self, color=None)

show_circle (self, color=None)

class allensdk.brain_observatory.circle_plots.FanPlotter (group_scale=0.9, *args,
                                                         **kwargs)
    Bases: allensdk.brain_observatory.circle_plots.PolarPlotter

    static for_drifting_gratings ()

    static for_static_gratings ()

    infer_dims (self, r_data, angle_data, group_data)

    plot (self, r_data, angle_data, group_data=None, data=None,
          cmap=<matplotlib.colors.LinearSegmentedColormap object at 0x7fc9a4777208>, clim=None,
          rmap=None, rlim=None, axis_color=None, label_color=None)

    set_dims (self, rs, angles, groups)

    show_angle_labels (self, angles=None, labels=None, color=None, offset=0.05, fontdict=None)

    show_axes (self, angles=None, radii=None, closed=False, color=None)

    show_group_labels (self, groups=None, color=None, fontdict=None)

    show_r_labels (self, radii=None, labels=None, color=None, offset=0.1, fontdict=None)

class allensdk.brain_observatory.circle_plots.PolarPlotter (direction=-1, an-
                                                            gle_start=0, cir-
                                                            cle_scale=1.1, in-
                                                            ner_radius=None,
                                                            plot_center=(0.0,
                                                            0.0), plot_scale=0.9)

    Bases: object

    DIR_CCW = 1

    DIR_CW = -1

    finalize (self)

class allensdk.brain_observatory.circle_plots.TrackPlotter (direction=-1, an-
                                                            gle_start=270.0,
                                                            inner_radius=0.45,
                                                            ring_length=None,
                                                            *args, **kwargs)

    Bases: allensdk.brain_observatory.circle_plots.PolarPlotter

    plot (self, data, clim=None, cmap=<matplotlib.colors.LinearSegmentedColormap object at
          0x7fc9a4777208>, mean_cmap=<matplotlib.colors.LinearSegmentedColormap object at
          0x7fc98a9db048>, norm=None)

    show_arrow (self, color=None)

allensdk.brain_observatory.circle_plots.add_angle_labels (ax, angles, labels, ra-
                                                          dius, color=None, font-
                                                          dict=None, offset=0.05)

```

```
allensdk.brain_observatory.circle_plots.add_arrow(ax, radius, start_angle, end_angle,
                                                  color=None, width=18.0)
allensdk.brain_observatory.circle_plots.angle_lines(angles, inner_radius,
                                                    outer_radius)
allensdk.brain_observatory.circle_plots.build_hex_pack(n)
allensdk.brain_observatory.circle_plots.hex_pack(radius, n)
allensdk.brain_observatory.circle_plots.make_pincushion_plot(data, trials, on,
                                                            nrows, ncols,
                                                            clim=None,
                                                            color_map=None,
                                                            radius=None)
allensdk.brain_observatory.circle_plots.polar_line_circles(radii, theta,
                                                           start_r=0)
allensdk.brain_observatory.circle_plots.polar_linspace(radius, start_angle,
                                                       stop_angle, num, end-
                                                       point=False, degrees=True)
    Evenly distributed list of x,y coordinates from an input range of angles and a radius in polar coordinates.
allensdk.brain_observatory.circle_plots.polar_to_xy(angles, radius)
    Convert an array of angles (in radians) and a radius in polar coordinates to an array of x,y coordinates.
allensdk.brain_observatory.circle_plots.radial_arcs(rs, start_theta, end_theta)
allensdk.brain_observatory.circle_plots.radial_circles(rs)
allensdk.brain_observatory.circle_plots.reset_hex_pack()
allensdk.brain_observatory.circle_plots.rings_in_hex_pack(ct)
allensdk.brain_observatory.circle_plots.spiral_trials(radii, x=0.0, y=0.0)
allensdk.brain_observatory.circle_plots.spiral_trials_polar(r, theta, radii, off-
                                                            set=None)
allensdk.brain_observatory.circle_plots.wedge_ring(N, inner_radius, outer_radius,
                                                    start=0, stop=360)
```

allensdk.brain_observatory.demixer module

```
allensdk.brain_observatory.demixer.demix_time_dep_masks(raw_traces:
                                                         numpy.ndarray,
                                                         stack: numpy.ndarray,
                                                         masks: numpy.ndarray,
                                                         max_block_size: int
                                                         = 1000) → Tuple[
                                                         numpy.ndarray, list]
```

Demix traces of potentially overlapping masks extracted from a single 2p recording.

Parameters

- **raw_traces** – 2d array of traces for each mask, of dimensions (n, t), where *t* is the number of time points and *n* is the number of masks.
- **stack** – 3d array representing a 1p recording movie, of dimensions (t, H, W) or corresponding hdf5 dataset.
- **masks** – 3d array of binary roi masks, of shape (n, H, W), where *n* is the number of masks, and HW are the dimensions of an individual frame in the movie *stack*.

Max_block_size int representing maximum number of movie frames to read at a time (-1 for full length *t* of *stack*) (the default is 1000)

Returns Tuple of demixed traces and whether each frame was skipped in the demixing calculation.

```
allensdk.brain_observatory.demixer.find_negative_baselines(trace)
allensdk.brain_observatory.demixer.find_negative_transients_threshold(trace,
                                                                    win-
                                                                    dow=500,
                                                                    length=10,
                                                                    std_devs=3)
allensdk.brain_observatory.demixer.find_zero_baselines(traces)
allensdk.brain_observatory.demixer.identify_valid_masks(mask_array)
allensdk.brain_observatory.demixer.plot_negative_baselines(raw_traces,
                                                            demix_traces,
                                                            mask_array,
                                                            roi_ids_mask,
                                                            plot_dir, ext='png')
allensdk.brain_observatory.demixer.plot_negative_transients(raw_traces,
                                                            demix_traces,
                                                            valid_roi,
                                                            mask_array,
                                                            roi_ids_mask,
                                                            plot_dir, ext='png')
allensdk.brain_observatory.demixer.plot_overlap_masks_lengthOne(roi_ind, masks,
                                                                savefile=None,
                                                                weighted=False)
allensdk.brain_observatory.demixer.plot_traces(raw_trace, demix_trace, roi_id, roi_ind,
                                                save_file)
allensdk.brain_observatory.demixer.plot_transients(roi_ind, t_trans, masks, traces,
                                                    demix_traces, savefile)
allensdk.brain_observatory.demixer.rolling_window(trace, window=500)
```

Parameters

- **trace** –
- **window** –

Returns

allensdk.brain_observatory.dff module

```
allensdk.brain_observatory.dff.calculate_dff(traces, dff_computation_cb=None,
                                              save_plot_dir=None)
```

Apply dF/F computation to a set of traces.

The default computation method is `compute_dff_windowed_median()` using default window parameters.

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

dff_computation_cb [function] Function that takes traces as an argument and returns an array of the same shape that is the calculated dF/F.

save_plot_dir [str] Directory to save dF/F plots to. By default no plots are saved.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.compute_dff_windowed_median(traces, median_kernel_long=5401, median_kernel_short=101, noise_stds=None, n_small_baseline_frames=None, **kwargs)
```

Compute dF/F of a set of traces with median filter detrending.

The operation is basically:

```
T_long = windowed_median(T) # long timescale kernel
T_dff1 = (T - T_long) / elementwise_max(T_long, noise_std(T))
T_short = windowed_median(T_dff1) # short timescale kernel
T_dff = T_dff1 - elementwise_min(T_short, 2.5*noise_std(T_dff1))
```

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

median_kernel_long [int] Window size to use for long timescale median detrending.

median_kernel_short [int] Window size to use for short timescale median detrending.

noise_stds [list] List that will contain noise_std(T_dff1) for each trace. The value for each trace will be appended to the list if provided.

n_small_baseline_frames [list] List that will contain the number of frames for each trace where the long-timescale median window is less than noise_std(T). The value for each trace will be appended to the list if provided.

kwargs: Additional keyword arguments are passed to `noise_std()`.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.compute_dff_windowed_mode(traces, mode_kernelsize=5400, mean_kernelsize=3000)
```

Compute dF/F of a set of traces using a low-pass windowed-mode operator.

The operation is basically:

```
T_mm = windowed_mean(windowed_mode(T))
T_dff = (T - T_mm) / T_mm
```

Parameters

traces [np.ndarray] 2D array of traces to be analyzed.

mode_kernelsize [int] Window size to use for windowed_mode.

mean_kernelsize [int] Window size to use for windowed_mean.

Returns

dff [np.ndarray] 2D array of dF/F traces.

```
allensdk.brain_observatory.dff.main()
```

```
allensdk.brain_observatory.dff.movingaverage(x, kernelsize, y)
```

Compute the windowed average of an array.

Parameters

x [np.ndarray] Array to be analyzed

kernelsize [int] Size of the moving window

y [np.ndarray] Output array to store the results

```
allensdk.brain_observatory.dff.movingmode_fast(x, kernelsize, y)
```

Compute the windowed mode of an array. A running mode is initialized with a histogram of values over the initial kernelsize/2 values. The mode is then updated as the kernel moves by adding and subtracting values from the histogram.

Parameters

x [np.ndarray] Array to be analyzed

kernelsize [int] Size of the moving window

y [np.ndarray] Output array to store the results

```
allensdk.brain_observatory.dff.noise_std(x, noise_kernel_length=31, positive_peak_scale=1.5, outlier_std_scale=2.5)
```

Robust estimate of the standard deviation of the trace noise.

```
allensdk.brain_observatory.dff.plot_onetrace(dff, fc)
```

Debug plotting function

```
allensdk.brain_observatory.dff.robust_std(x)
```

Robust estimate of standard deviation.

Estimate of the standard deviation using the median absolute deviation of x.

allensdk.brain_observatory.drifting_gratings module

```
class allensdk.brain_observatory.drifting_gratings.DriftingGratings(data_set, **kwargs)
```

Bases: `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`

Perform tuning analysis specific to drifting gratings stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

static from_analysis_file (data_set, analysis_file)

get_noise_correlation (self, corr='spearman')

get_peak (self)

Computes metrics related to each cell's peak response condition.

Returns

Pandas data frame containing the following columns (_dg suffix is for drifting grating):

- ori_dg (orientation)
- tf_dg (temporal frequency)
- reliability_dg
- osi_dg (orientation selectivity index)
- dsi_dg (direction selectivity index)
- peak_dff_dg (peak dF/F)
- ptest_dg
- p_run_dg
- run_modulation_dg
- cv_dg (circular variance)

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# orientations, # temporal frequencies, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant response ($p < 0.05$) to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

number_ori

number_tf

open_star_plot (*self*, *cell_specimen_id*=None, *include_labels*=False, *cell_index*=None)

orivals

plot_direction_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_orientation_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_direction (*self*, *include_labels*=False, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_temporal_frequency (*self*, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim x repetition for noise correlations

tfvals

allensdk.brain_observatory.findlevel module

`allensdk.brain_observatory.findlevel.findlevel` (*inwave, threshold, direction='both'*)

allensdk.brain_observatory.locally_sparse_noise module

class `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` (*data_set, stimulus=None, **kwargs*)

Bases: `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`

Perform tuning analysis specific to the locally sparse noise stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

stimulus: string Name of locally sparse noise stimulus. See `brain_observatory.stimulus_info`.

nrows: int Number of rows in the stimulus template

ncol: int Number of columns in the stimulus template

LSN

LSN_GREY = 127

LSN_OFF = 0

LSN_OFF_SCREEN = 64

LSN_ON = 255

LSN_mask

cell_index_receptive_field_analysis_data

extralength

static from_analysis_file (*data_set, analysis_file, stimulus*)

get_mean_response (*self*)

get_peak (*self*)

Implemented by subclasses.

get_receptive_field (*self*)

Calculates receptive fields for each cell

get_receptive_field_analysis_data (*self*)

Calculates receptive fields for each cell

get_receptive_field_attribute_df (*self*)

interlength

mean_response

static merge_mean_response (*rc1, rc2*)

Move out of this class, to session analysis

open_pincushion_plot (*self, on, cell_specimen_id=None, color_map=None, cell_index=None*)

```
plot_cell_receptive_field(self, on, cell_specimen_id=None, color_map=None, clim=None,
                           mask=None, cell_index=None, scalebar=True)
plot_population_receptive_field(self, color_map='RdPu', clim=None, mask=None,
                                scalebar=True)
plot_receptive_field_analysis_data(self, cell_index, **kwargs)
populate_stimulus_table(self)
    Implemented by subclasses.
static read_cell_index_receptive_field_analysis(file_handle, prefix, path=None)
receptive_field
static save_cell_index_receptive_field_analysis(cell_index_receptive_field_analysis_data,
                                                new_nwb, prefix)

sort_trials(self)
sweeplength
```

allensdk.brain_observatory.natural_movie module

```
class allensdk.brain_observatory.natural_movie.NaturalMovie(data_set,
                                                             movie_name,
                                                             **kwargs)
Bases: allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
```

Perform tuning analysis specific to natural movie stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

movie_name: string

one of [stimulus_info.NATURAL_MOVIE_ONE, stimulus_info.NATURAL_MOVIE_TWO,
stimulus_info.NATURAL_MOVIE_THREE]

static from_analysis_file (data_set, analysis_file, movie_name)

get_peak (self)

Computes properties of the peak response condition for each cell.

Returns

Pandas data frame with the below fields. A suffix of “nm1”, “nm2” or “nm3” is appended to the field name

on which of three movie clips was presented.

- peak_nm1 (frame with peak response)
- response_variability_nm1

get_sweep_response (self)

Returns the dF/F response for each cell

Returns

Numpy array

open_track_plot (self, cell_specimen_id=None, cell_index=None)

populate_stimulus_table (*self*)

Implemented by subclasses.

sweep_response

sweeplength

allensdk.brain_observatory.natural_scenes module

class allensdk.brain_observatory.natural_scenes.**NaturalScenes** (*data_set*,
***kwargs*)

Bases: *allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis*

Perform tuning analysis specific to natural scenes stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

extralength

static from_analysis_file (*data_set*, *analysis_file*)

get_noise_correlation (*self*, *corr*='spearman')

get_peak (*self*)

Computes metrics about peak response condition for each cell.

Returns

Pandas data frame with the following fields ('_ns' suffix is for natural scene):

- scene_ns (scene number)
- reliability_ns
- peak_dff_ns (peak dF/F)
- ptest_ns
- p_run_ns
- run_modulation_ns
- time_to_peak_ns

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# scenes, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant ($p < 0.05$) response to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

interlength

number_scenes

open_corona_plot (*self*, *cell_specimen_id*=None, *cell_index*=None)

plot_time_to_peak (*self*, *p_value_max*=0.05, *color_map*=<matplotlib.colors.LinearSegmentedColormap object at 0x7fc9b307bcf8>)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim x repetition for noise correlations

sweeplength

allensdk.brain_observatory.observatory_plots module

class allensdk.brain_observatory.observatory_plots.**DimensionPatchHandler** (*vals*,
start_color,
end_color,
**args*,
***kwargs*)

Bases: object

dim_color (*self*, *index*)

legend_artist (*self*, *legend*, *orig_handle*, *fontsize*, *handlebox*)

allensdk.brain_observatory.observatory_plots.**figure_in_px** (*w*, *h*, *file_name*,
dpi=96.0, *transparent*=False)

allensdk.brain_observatory.observatory_plots.**finalize_no_axes** (*pad*=0.0)

allensdk.brain_observatory.observatory_plots.**finalize_no_labels** (*pad*=0.3, *legend*=False)

allensdk.brain_observatory.observatory_plots.**finalize_with_axes** (*pad*=0.3)

allensdk.brain_observatory.observatory_plots.**float_label** (*n*)

allensdk.brain_observatory.observatory_plots.**plot_cell_correlation** (*sig_corrs*,
labels,
colors,
scale=15)

allensdk.brain_observatory.observatory_plots.**plot_combined_speed** (*binmed_resp_vis*,
binmed_dx_vis,
binmed_resp_sp,
binmed_dx_sp,
evoked_color,
spont_color)

allensdk.brain_observatory.observatory_plots.**plot_condition_histogram** (*vals*,
bins,
color='#ccccdd')

allensdk.brain_observatory.observatory_plots.**plot_mask_outline** (*mask*, *ax*,
color='k')

```

allensdk.brain_observatory.observatory_plots.plot_pupil_location(xy_deg, s=1,
                                                                    c=None,
                                                                    cmap=<matplotlib.colors.LinearSegmentedColormap: 0x7fc98a0a84a8>,
                                                                    object at
                                                                    0x7fc98a0a84a8>,
                                                                    edge-
                                                                    color=", in-
                                                                    clude_labels=True)

allensdk.brain_observatory.observatory_plots.plot_radial_histogram(angles,
                                                                    counts,
                                                                    all_angles=None,
                                                                    in-
                                                                    clude_labels=False,
                                                                    off-
                                                                    set=180.0,
                                                                    direction=-
                                                                    1,
                                                                    closed=False,
                                                                    color='#ccccdd')

allensdk.brain_observatory.observatory_plots.plot_receptive_field(rf,
                                                                    color_map=None,
                                                                    clim=None,
                                                                    mask=None,
                                                                    out-
                                                                    line_color='#cccccc',
                                                                    scale-
                                                                    bar=True)

allensdk.brain_observatory.observatory_plots.plot_representational_similarity(rs,
                                                                    dims=None,
                                                                    dim_labels=None,
                                                                    col-
                                                                    ors=None,
                                                                    dim_order=None,
                                                                    la-
                                                                    bels=True)

allensdk.brain_observatory.observatory_plots.plot_selectivity_cumulative_histogram(sis,
                                                                    xla-
                                                                    bel,
                                                                    si_range=[
                                                                    1.5],
                                                                    n_hist_bins,
                                                                    color='#cccccc')

allensdk.brain_observatory.observatory_plots.plot_speed(binned_resp, binned_dx,
                                                                    num_bins, color)

allensdk.brain_observatory.observatory_plots.plot_time_to_peak(msrs, ttps,
                                                                    t_start, t_end,
                                                                    stim_start,
                                                                    stim_end,
                                                                    cmap)

```

```
allensdk.brain_observatory.observatory_plots.population_correlation_scatter(sig_corrs,  
                                                                           noise_corrs,  
                                                                           la-  
                                                                           bels,  
                                                                           col-  
                                                                           ors,  
                                                                           scale=15)
```

allensdk.brain_observatory.r_neuropil module

```
class allensdk.brain_observatory.r_neuropil.NeuropilSubtract (lam=0.05, dt=1.0,  
                                                             folds=4)
```

Bases: object

TODO: docs

```
estimate_error (self, r)
```

Estimate error values for a given r for each fold and return the mean.

```
fit (self, r_range=[0.0, 2.0], iterations=3, dr=0.1, dr_factor=0.1)
```

Estimate error values for a range of r values. Identify a new r range around the minimum error values and repeat multiple times. TODO: docs

```
fit_block_coordinate_desc (self, r_init=5.0, min_delta_r=1e-08)
```

```
set_F (self, F_M, F_N)
```

Break the F_M and F_N traces into the number of folds specified in the class constructor and normalize each fold of F_M and R_N relative to F_N.

```
allensdk.brain_observatory.r_neuropil.ab_from_T (T, lam, dt)
```

```
allensdk.brain_observatory.r_neuropil.ab_from_diagonals (mat_dict)
```

Constructs value for `scipy.linalg.solve_banded`

Parameters

mat_dict: dictionary of diagonals keyed by offsets

Returns

ab: value for `scipy.linalg.solve_banded`

```
allensdk.brain_observatory.r_neuropil.alpha_filter (A=1.0, alpha=0.05, beta=0.25,  
                                                    T=100)
```

```
allensdk.brain_observatory.r_neuropil.error_calc (F_M, F_N, F_C, r)
```

```
allensdk.brain_observatory.r_neuropil.error_calc_outlier (F_M, F_N, F_C, r)
```

```
allensdk.brain_observatory.r_neuropil.estimate_contamination_ratios (F_M,  
                                                                           F_N,  
                                                                           lam=0.05,  
                                                                           folds=4,  
                                                                           itera-  
                                                                           tions=3,  
                                                                           r_range=[0.0,  
                                                                           2.0],  
                                                                           dr=0.1,  
                                                                           dr_factor=0.1)
```

Calculates neuropil contamination of ROI

Parameters

F_M: ROI trace F_N: Neuropil trace

Returns

dictionary: key-value pairs

- 'r': the contamination ratio – corrected trace = $M - r * N$
- 'err': RMS error
- 'min_error': minimum error
- 'bounds_error': boolean. True if error or R are outside tolerance

`allensdk.brain_observatory.r_neuropil.get_diagonals_from_sparse(mat)`

Returns a dictionary of diagonals keyed by offsets

Parameters

mat: scipy.sparse matrix

Returns

dictionary: diagonals keyed by offsets

`allensdk.brain_observatory.r_neuropil.normalize_F(F_M, F_N)`

`allensdk.brain_observatory.r_neuropil.synthesize_F(T, af1, af2, p1=0.05, p2=0.1)`

Build a synthetic F_C, F_M, F_N, and r of length T TODO: docs

`allensdk.brain_observatory.r_neuropil.validate_with_synthetic_F(T, N)`

Compute N synthetic traces of length T with known values of r, then estimate r. TODO: docs

allensdk.brain_observatory.roi_masks module

class `allensdk.brain_observatory.roi_masks.Mask` (*image_w*, *image_h*, *label*,
mask_group)

Bases: object

Abstract class to represent image segmentation mask. Its two main subclasses are RoiMask and NeuropilMask. The former represents the mask of a region of interest (ROI), such as a cell observed in 2-photon imaging. The latter represents the neuropil around that cell, and is useful when subtracting the neuropil signal from the measured ROI signal.

This class should not be instantiated directly.

Parameters

image_w: integer Width of image that ROI resides in

image_h: integer Height of image that ROI resides in

label: text User-defined text label to identify mask

mask_group: integer User-defined number to help put masks into different categories

get_mask_plane (*self*)

Returns mask content on full-size image plane

Returns

numpy 2D array [*img_rows*][*img_cols*]

init_by_pixels (*self*, *border*, *pix_list*)

Initialize mask using a list of mask pixels

Parameters

border: float[4] Coordinates defining useable area of image. See `create_roi_mask()`

pix_list: integer[][2] List of pixel coordinates (x,y) that define the mask

overlaps_motion_border

class `allensdk.brain_observatory.roi_masks.NeuropilMask` (*w, h, label, mask_group*)

Bases: `allensdk.brain_observatory.roi_masks.Mask`

init_by_mask (*self, border, array*)

Initialize mask using spatial mask

Parameters

border: float[4] Border widths on the [right, left, down, up] sides. The resulting neuropil mask will not include pixels falling into a border.

array: integer[image height][image width] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

class `allensdk.brain_observatory.roi_masks.RoiMask` (*image_w, image_h, label, mask_group*)

Bases: `allensdk.brain_observatory.roi_masks.Mask`

init_by_mask (*self, border, array*)

Initialize mask using spatial mask

Parameters

border: float[4] Coordinates defining useable area of image. See `create_roi_mask()`.

roi_mask: integer[image height][image width] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

`allensdk.brain_observatory.roi_masks.calculate_roi_and_neuropil_traces` (*movie_h5, roi_mask_list, motion_border*)

get roi and neuropil masks

`allensdk.brain_observatory.roi_masks.calculate_traces` (*stack, mask_list, block_size=1000*)

Calculates the average response of the specified masks in the image stack

Parameters

stack: float[image height][image width] Image stack that masks are applied to

mask_list: list<Mask> List of masks

Returns

float[number masks][number frames] This is the average response for each Mask in each image frame

`allensdk.brain_observatory.roi_masks.create_neuropil_mask` (*roi, border, combined_binary_mask, label=None*)

Convenience function to create and initialize a Neuropil mask. Neuropil masks are defined as the region around an ROI, up to 13 pixels out, that does not include other ROIs

Parameters

roi: RoiMask object The ROI that the neuropil masks will be based on

border: float[4] Border widths on the [right, left, down, up] sides. The resulting neuropil mask will not include pixels falling into a border.

combined_binary_mask List of pixel coordinates (x,y) that define the mask

combined_binary_mask: integer[image_h][image_w] Image-sized array that shows the position of all ROIs in the image. ROI masks should have a value of one. Background pixels must be zero. In other words, the combined_binary_mask is a bitmap union of all ROI masks

label: text User-defined text label to identify the mask

Returns

NeuropilMask object

```
allensdk.brain_observatory.roi_masks.create_roi_mask(image_w, image_h, border,
                                                    pix_list=None,
                                                    roi_mask=None, label=None,
                                                    mask_group=-1)
```

Convenience function to create and initialize an RoiMask

Parameters

image_w: integer Width of image that ROI resides in

image_h: integer Height of image that ROI resides in

border: float[4] Coordinates defining useable area of image. If the entire image is usable, and masks are valid anywhere in the image, this should be [0, 0, 0, 0]. The following constants help describe the array order:

RIGHT_SHIFT = 0

LEFT_SHIFT = 1

DOWN_SHIFT = 2

UP_SHIFT = 3

When parts of the image are unusable, for example due to motion correction shifting of different image frames, the border array should store the usable image area

pix_list: integer[][2] List of pixel coordinates (x,y) that define the mask

roi_mask: integer[image_h][image_w] Image-sized array that describes the mask. Active parts of the mask should have values >0. Background pixels must be zero

label: text User-defined text label to identify mask

mask_group: integer User-defined number to help put masks into different categories

Returns

RoiMask object

```
allensdk.brain_observatory.roi_masks.create_roi_mask_array(rois)
```

Create full image mask array from list of RoiMasks.

Parameters

rois: list<RoiMask> List of roi masks.

Returns

np.ndarray: NxWxH array Boolean array of len(rois) image masks.

`allensdk.brain_observatory.roi_masks.validate_mask(mask)`
Check a given roi or neuropil mask for (a subset of) disqualifying problems.

`allensdk.brain_observatory.running_speed` module

class `allensdk.brain_observatory.running_speed.RunningSpeed`
Bases: `tuple`
Describes the rate at which an experimental subject ran during a session.
values `[np.ndarray]` running speed (cm/s) at each sample point
timestamps `[np.ndarray]` The time at which each sample was collected (s).
timestamps
Alias for field number 0
values
Alias for field number 1

`allensdk.brain_observatory.session_analysis` module

class `allensdk.brain_observatory.session_analysis.SessionAnalysis` (*nwb_path*,
save_path)
Bases: `object`

Run all of the stimulus-specific analyses associated with a single experiment session.

Parameters

nwb_path: string, path to NWB file

save_path: string, path to HDF5 file to store outputs. Recommended NOT to modify the NWB file.

append_experiment_metrics (*self*, *metrics*)
Extract stimulus-agnostic metrics from an experiment into a dictionary

append_metadata (*self*, *df*)
Append the metadata fields from the NWB file as columns to a `pd.DataFrame`

append_metrics_drifting_grating (*self*, *metrics*, *dg*)
Extract metrics from the DriftingGratings peak response table into a dictionary.

append_metrics_locally_sparse_noise (*self*, *metrics*, *lsn*)
Extract metrics from the LocallySparseNoise peak response table into a dictionary.

append_metrics_natural_movie_one (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_ONE) peak response table into a dictionary.

append_metrics_natural_movie_three (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_THREE) peak response table into a dictionary.

append_metrics_natural_movie_two (*self*, *metrics*, *nma*)
Extract metrics from the NaturalMovie(stimulus_info.NATURAL_MOVIE_TWO) peak response table into a dictionary.

append_metrics_natural_scene (*self*, *metrics*, *ns*)
Extract metrics from the NaturalScenes peak response table into a dictionary.

append_metrics_static_grating (*self, metrics, sg*)

Extract metrics from the StaticGratings peak response table into a dictionary.

save_session_a (*self, dg, nm1, nm3, peak*)

Save the output of session A analysis to self.save_path.

Parameters

dg: DriftingGratings instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm3: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_THREE

peak: pd.DataFrame The combined peak response property table created in self.session_a().

save_session_b (*self, sg, nm1, ns, peak*)

Save the output of session B analysis to self.save_path.

Parameters

sg: StaticGratings instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

ns: NaturalScenes instance

peak: pd.DataFrame The combined peak response property table created in self.session_b().

save_session_c (*self, lsn, nm1, nm2, peak*)

Save the output of session C analysis to self.save_path.

Parameters

lsn: LocallySparseNoise instance

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm2: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_TWO

peak: pd.DataFrame The combined peak response property table created in self.session_c().

save_session_c2 (*self, lsn4, lsn8, nm1, nm2, peak*)

Save the output of session C2 analysis to self.save_path.

Parameters

lsn4: LocallySparseNoise instance This LocallySparseNoise instance should have been created with self.stimulus = stimulus_info.LOCALLY_SPARSE_NOISE_4DEG.

lsn8: LocallySparseNoise instance This LocallySparseNoise instance should have been created with self.stimulus = stimulus_info.LOCALLY_SPARSE_NOISE_8DEG.

nm1: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_ONE

nm2: NaturalMovie instance This NaturalMovie instance should have been created with movie_name = stimulus_info.NATURAL_MOVIE_TWO

peak: pd.DataFrame The combined peak response property table created in `self.session_c2()`.

session_a (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural movie one, natural movie three, and drifting gratings. The input NWB be for a `stimulus_info.THREE_SESSION_A` experiment.

Parameters

plot_flag: bool Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: bool Whether to save the output of analysis to `self.save_path` upon completion.

session_b (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural scenes, static gratings, and natural movie one. The input NWB be for a `stimulus_info.THREE_SESSION_B` experiment.

Parameters

plot_flag: bool Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: bool Whether to save the output of analysis to `self.save_path` upon completion.

session_c (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for natural movie one, natural movie two, and locally sparse noise. The input NWB be for a `stimulus_info.THREE_SESSION_C` experiment.

Parameters

plot_flag: bool Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: bool Whether to save the output of analysis to `self.save_path` upon completion.

session_c2 (*self*, *plot_flag=False*, *save_flag=True*)

Run stimulus-specific analysis for locally sparse noise (4 deg.), locally sparse noise (8 deg.), natural movie one, and natural movie two. The input NWB be for a `stimulus_info.THREE_SESSION_C2` experiment.

Parameters

plot_flag: bool Whether to generate `brain_observatory_plotting` work plots after running analysis.

save_flag: bool Whether to save the output of analysis to `self.save_path` upon completion.

verify_roi_lists_equal (*self*, *roi1*, *roi2*)

TODO: replace this with simpler numpy comparisons

`allensdk.brain_observatory.session_analysis.main()`

`allensdk.brain_observatory.session_analysis.multi_dataframe_merge(dfs)`

merge a number of `pd.DataFrames` into a single dataframe on their index columns. If any columns are duplicated, prefer the first occuring instance of the column

```
allensdk.brain_observatory.session_analysis.run_session_analysis(nwb_path,
                                                                save_path,
                                                                plot_flag=False,
                                                                save_flag=True)
```

Inspect an NWB file to determine which experiment session was run and compute all stimulus-specific analyses.

Parameters

nwb_path: `string` Path to NWB file.

save_path: `string` path to save results. Recommended NOT to use NWB file.

plot_flag: `bool` Whether to save brain_observatory_plotting work plots.

save_flag: `bool` Whether to save results to save_path.

allensdk.brain_observatory.session_api_utils module

```
class allensdk.brain_observatory.session_api_utils.ParamsMixin(ignore: set =
                                                                {'api'})
```

Bases: `object`

This mixin adds parameter management functionality to the class it is mixed into.

This mixin expects that the class it is mixed into will have an `__init__` with type annotated parameters. It also expects for the class to have semi-private attributes of the `__init__` type annotated parameters.

Example:

SomeClassWhereParamManagementIsDesired(ParamsMixin):

```
# Managed params should be typed (with simple types if possible)! def __init__(self,
param_to_ignore, a_param_1: int, a_param_2: float,
    b_param_1: list):
    # Parameters can be ignored by the mixin super().__init__(ignore={'param_to_ignore'})
    # Pay attention to the naming scheme! self._a_param_1 = a_param_1 self._a_param_2
    = a_param_2 self._b_param_1 = b_param_1
    ...
```

After being mixed in, methods like `'get_params'`, `'set_params'`, `'needs_data_refresh'`, and `'clear_updated_params'` will be available.

clear_updated_params (*self*, *data_params: set*)
This method clears 'updated params' whose data have been updated

get_params (*self*) → `Dict[str, Any]`
Get managed params and their values

needs_data_refresh (*self*, *data_params: set*) → `bool`
Check if specific params have been updated via *set_params()*

set_params (*self*, ***params*)
Set managed params

allensdk.brain_observatory.session_api_utils.is_equal (*a: Any*, *b: Any*) → `bool`
Function to deal with checking if two variables of possibly mixed types have the same value.

allensdk.brain_observatory.static_gratings module

class allensdk.brain_observatory.static_gratings.**StaticGratings** (*data_set*,
***kwargs*)
Bases: *allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis*

Perform tuning analysis specific to static gratings stimulus.

Parameters

data_set: BrainObservatoryNwbDataSet object

extralength

static_from_analysis_file (*data_set*, *analysis_file*)

get_noise_correlation (*self*, *corr*='spearman')

get_peak (*self*)

Computes metrics related to each cell's peak response condition.

Returns

Panda data frame with the following fields (**_sg** suffix is

for static grating):

- **ori_sg** (orientation)
- **sf_sg** (spatial frequency)
- **phase_sg**
- **response_variability_sg**
- **osi_sg** (orientation selectivity index)
- **peak_dff_sg** (peak dF/F)
- **ptest_sg**
- **time_to_peak_sg**

get_representational_similarity (*self*, *corr*='spearman')

get_response (*self*)

Computes the mean response for each cell to each stimulus condition. Return is a (# orientations, # spatial frequencies, # phasees, # cells, 3) np.ndarray. The final dimension contains the mean response to the condition (index 0), standard error of the mean of the response to the condition (index 1), and the number of trials with a significant response ($p < 0.05$) to that condition (index 2).

Returns

Numpy array storing mean responses.

get_signal_correlation (*self*, *corr*='spearman')

interlength

number_ori

number_phase

number_sf

open_fan_plot (*self*, *cell_specimen_id*=None, *include_labels*=False, *cell_index*=None)

orivals

phasevals

plot_orientation_selectivity (*self*, *si_range*=[0, 1.5], *n_hist_bins*=50, *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_orientation (*self*, *include_labels*=False, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_preferred_spatial_frequency (*self*, *si_range*=[0, 1.5], *color*='#ccccdd',
p_value_max=0.05, *peak_dff_min*=3)

plot_time_to_peak (*self*, *p_value_max*=0.05, *color_map*=<matplotlib.colors.LinearSegmentedColormap
object at 0x7fc9b307bcf8>)

populate_stimulus_table (*self*)

Implemented by subclasses.

reshape_response_array (*self*)

Returns response array in cells x stim conditions x repetition for noise correlations

this is a re-organization of the mean sweep response table

sfvals

sweeplength

allensdk.brain_observatory.stimulus_analysis module

class allensdk.brain_observatory.stimulus_analysis.**StimulusAnalysis** (*data_set*)

Bases: object

Base class for all response analysis code. Subclasses are responsible for computing metrics and traces relevant to a particular stimulus. The base class contains methods for organizing sweep responses row of a stimulus stable (get_sweep_response). Subclasses implement the get_response method, computes the mean sweep response to all sweeps for a each stimulus condition.

Parameters

data_set: BrainObservatoryNwbDataSet instance

speed_tuning: boolean, deprecated Whether or not to compute speed tuning histograms

acquisition_rate

binned_cells_sp

binned_cells_vis

binned_dx_sp

binned_dx_vis

cell_id

celltraces

dfftraces

dxcm

dxtime

get_fluorescence (*self*)

get_peak (*self*)

Implemented by subclasses.

get_response (*self*)

Implemented by subclasses.

get_speed_tuning (*self*, *binsize*)

Calculates speed tuning, spontaneous versus visually driven. The return is a 5-tuple of speed and dF/F histograms.

binned_dx_sp: (bins,2) np.ndarray of running speeds binned during spontaneous activity stimulus. The first bin contains all speeds below 1 cm/s. Dimension 0 is mean running speed in the bin. Dimension 1 is the standard error of the mean.

binned_cells_sp: (bins,2) np.ndarray of fluorescence during spontaneous activity stimulus. First bin contains all data for speeds below 1 cm/s. Dimension 0 is mean fluorescence in the bin. Dimension 1 is the standard error of the mean.

binned_dx_vis: (bins,2) np.ndarray of running speeds outside of spontaneous activity stimulus. The first bin contains all speeds below 1 cm/s. Dimension 0 is mean running speed in the bin. Dimension 1 is the standard error of the mean.

binned_cells_vis: np.ndarray of fluorescence outside of spontaneous activity stimulus. First bin contains all data for speeds below 1 cm/s. Dimension 0 is mean fluorescence in the bin. Dimension 1 is the standard error of the mean.

peak_run: pd.DataFrame of speed-related properties of a cell.

Returns

tuple: binned_dx_sp, binned_cells_sp, binned_dx_vis, binned_cells_vis, peak_run

get_sweep_response (*self*)

Calculates the response to each sweep in the stimulus table for each cell and the mean response. The return is a 3-tuple of:

- **sweep_response**: pd.DataFrame of response dF/F traces organized by cell (column) and sweep (row)
- **mean_sweep_response**: mean values of the traces returned in **sweep_response**
- **pval**: p value from 1-way ANOVA comparing response during sweep to response prior to sweep

Returns

3-tuple: sweep_response, mean_sweep_response, pval

mean_sweep_response

numbercells

peak

peak_run

plot_representational_similarity (*self*, *repsim*, *stimulus=False*)

plot_running_speed_histogram (*self*, *xlim=None*, *nbins=None*)

plot_speed_tuning (*self*, *cell_specimen_id=None*, *cell_index=None*, *evoked_color='#b30000'*, *spontaneous_color='#0000b3'*)

populate_stimulus_table (*self*)

Implemented by subclasses.

pval

response

```

roi_id
row_from_cell_id (self, csid=None, idx=None)
stim_table
sweep_response
timestamps

```

`allensdk.brain_observatory.stimulus_analysis.nonraising_ks_2samp` (*data1*, *data2*,
***kwargs*
scipy.stats.ks_2samp now raises a `ValueError` if one of the input arrays is of length 0. Previously it signaled this case by returning nans. This function restores the prior behavior.

`allensdk.brain_observatory.stimulus_info` module

```

class allensdk.brain_observatory.stimulus_info.BinaryIntervalSearchTree (search_list)
    Bases: object
    add (self, input_list, tmp=None)
    static from_df (input_df)
    search (self, fi, tmp=None)

```

```

class allensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor (experiment_geometry=None)
    Bases: allensdk.brain_observatory.stimulus_info.Monitor
    http://help.brain-map.org/display/observatory/Documentation?preview=/10616846/10813485/VisualCoding\_VisualStimuli.pdf https://www.cnet.com/products/asus-pa248q/specs/
    grating_to_screen (self, phase, spatial_frequency, orientation, **kwargs)
    lsn_image_to_screen (self, img, **kwargs)
    pixels_to_visual_degrees (self, n, **kwargs)
    visual_degrees_to_pixels (self, vd, **kwargs)
    warp_image (self, img, **kwargs)

```

```

class allensdk.brain_observatory.stimulus_info.ExperimentGeometry (distance,
                                                                    mon_height_cm,
                                                                    mon_width_cm,
                                                                    mon_res,
                                                                    eyepoint)
    Bases: object
    generate_warp_coordinates (self)
    warp_coordinates

```

```

class allensdk.brain_observatory.stimulus_info.Monitor (n_pixels_r,      n_pixels_c,
                                                         panel_size, spatial_unit)
    Bases: object
    aspect_ratio
    get_mask (self)
    grating_to_screen (self, phase, spatial_frequency, orientation, distance_from_monitor,
                       p2p_amp=256, baseline=127, translation=(0, 0))
    height

```

```
lsn_image_to_screen (self, img, stimulus_type, origin='lower', background_color=127, translation=(0, 0))
map_stimulus (self, source_stimulus_coordinate, source_stimulus_type, target_stimulus_type)
mask
natural_movie_image_to_screen (self, img, origin='lower', translation=(0, 0))
natural_scene_image_to_screen (self, img, origin='lower', translation=(0, 0))
panel_size
pixel_size
pixels_to_visual_degrees (self, n, distance_from_monitor, small_angle_approximation=True)
set_spatial_unit (self, new_unit)
show_image (self, img, ax=None, show=True, mask=False, warp=False, origin='lower')
spatial_frequency_to_pix_per_cycle (self, spatial_frequency, distance_from_monitor)
visual_degrees_to_pixels (self, vd, distance_from_monitor, small_angle_approximation=True)
width
class allensdk.brain_observatory.stimulus_info.StimulusSearch (nwb_dataset)
    Bases: object
    search (self, fi)
allensdk.brain_observatory.stimulus_info.all_stimuli ()
    Return a list of all stimuli in the data set
allensdk.brain_observatory.stimulus_info.get_spatial_grating (height=None, aspect_ratio=None, ori=None, pix_per_cycle=None, phase=None, p2p_amp=2, baseline=0)
allensdk.brain_observatory.stimulus_info.get_spatio_temporal_grating (t, tempo-ral_frequency=None, **kwargs)
allensdk.brain_observatory.stimulus_info.lsn_coordinate_to_monitor_coordinate (lsn_coordinate, monitor_shape, stimulus_type)
allensdk.brain_observatory.stimulus_info.make_display_mask (display_shape=(1920, 1200))
    Build a display-shaped mask that indicates which pixels are on screen after warping the stimulus.
```



```

allensdk.brain_observatory.stimulus_info.map_monitor_coordinate_to_stimulus_coordinate(moni-
                                                                                       i-
                                                                                       tor_s
                                                                                       stim-
                                                                                       u-
                                                                                       lus_t

allensdk.brain_observatory.stimulus_info.map_monitor_coordinate_to_template_coordinate(moni-
                                                                                       i-
                                                                                       tor_s
                                                                                       tem-
                                                                                       plate

allensdk.brain_observatory.stimulus_info.map_stimulus(source_stimulus_coordinate,
                                                       source_stimulus_type, tar-
                                                       get_stimulus_type, moni-
                                                       tor_shape)

allensdk.brain_observatory.stimulus_info.map_stimulus_coordinate_to_monitor_coordinate(temp
                                                                                       mon-
                                                                                       i-
                                                                                       tor_s
                                                                                       stim-
                                                                                       u-
                                                                                       lus_t

allensdk.brain_observatory.stimulus_info.map_template_coordinate_to_monitor_coordinate(temp
                                                                                       mon-
                                                                                       i-
                                                                                       tor_s
                                                                                       tem-
                                                                                       plate

allensdk.brain_observatory.stimulus_info.mask_stimulus_template(template_display_coords,
                                                                tem-
                                                                plate_shape,
                                                                dis-
                                                                play_mask=None,
                                                                thresh-
                                                                old=1.0)

```

Build a mask for a stimulus template of a given shape and display coordinates that indicates which part of the template is on screen after warping.

Parameters

- template_display_coords:** **list** list of (x,y) display coordinates
- template_shape:** **tuple** (width,height) of the display template
- display_mask:** **np.ndarray** boolean 2D mask indicating which display coordinates are on screen after warping.
- threshold:** **float** Fraction of pixels associated with a template display coordinate that should remain on screen to count as belonging to the mask.

Returns

- tuple:** (template mask, pixel fraction)

```
allensdk.brain_observatory.stimulus_info.monitor_coordinate_to_lsn_coordinate(monitor_coordinate,
                                     monitor_shape,
                                     monitor_type,
                                     lsn_shape,
                                     lsn_type)
```

```
allensdk.brain_observatory.stimulus_info.monitor_coordinate_to_natural_movie_coordinate(monitor, monitor_i-, monitor)
```

```
allensdk.brain_observatory.stimulus_info.natural_movie_coordinate_to_monitor_coordinate(natural_movie_coordinate, monitor_coordinate)
```

```
allensdk.brain_observatory.stimulus_info.natural_scene_coordinate_to_monitor_coordinate(natural_scene_coordinate, monitor_coordinate)
```

```
allensdk.brain_observatory.stimulus_info.rotate(X, Y, theta)
```

```
allensdk.brain_observatory.stimulus_info.sessions_with_stimulus(stimulus)
```

Return the names of the sessions that contain a given stimulus.

```
allensdk.brain_observatory.stimulus_info.stimuli_in_session(session, al-  
low_unknown=True)
```

Return a list what stimuli are available in a given session.

Parameters

session: string Must be one of: [stimulus_info.THREE_SESSION_A, stimulus_info.THREE_SESSION_B, stimulus_info.THREE_SESSION_C, stimulus_info.THREE_SESSION_C2]

```
allensdk.brain_observatory.stimulus_info.translate_image_and_fill(img, trans-  
lation=(0,  
0))
```

```
allensdk.brain_observatory.stimulus_info.warp_stimulus_coords(vertices, distance=15.0,
mon_height_cm=32.5,
mon_width_cm=51.0,
mon_res=(1920,
1200), eye-
point=(0.5, 0.5))
```

For a list of screen vertices, provides a corresponding list of texture coordinates.

Parameters

vertices: `numpy.ndarray` `[[x0,y0], [x1,y1], ...]` A set of vertices to convert to texture positions.

distance: float distance from the monitor in cm.

mon_height cm: float monitor height in cm

mon_width_cm: float monitor width in cm

mon_res: tuple monitor resolution (x,y)

eyepoint: tuple

Returns

np.ndarray x,y coordinates shaped like the input that describe what pixel coordinates are displayed on the input coordinates after warping the stimulus.

allensdk.brain_observatory.sync_dataset module

dataset.py

Dataset object for loading and unpacking an HDF5 dataset generated by sync.py

@author: derriew

Allen Institute for Brain Science

Dependencies

numpy <http://www.numpy.org/> h5py <http://www.h5py.org/>

class allensdk.brain_observatory.sync_dataset.**Dataset** (*path*)
Bases: object

A sync dataset. Contains methods for loading and parsing the binary data.

Parameters

path [str] Path to HDF5 file.

Examples

```
>>> dset = Dataset('my_h5_file.h5')
>>> logger.info(dset.meta_data)
>>> dset.stats()
>>> dset.close()
```

```
>>> with Dataset('my_h5_file.h5') as d:
...     logger.info(dset.meta_data)
...     dset.stats()
```

BEHAVIOR_TRACKING_KEYS = ('beh_frame_received', 'cam1_exposure', 'behavior_monitoring')

DEPRECATED_KEYS = {'behavior_monitoring', 'cam1_exposure', 'cam2_exposure', 'eye_tracking'}

EYE_TRACKING_KEYS = ('eye_frame_received', 'cam2_exposure', 'eyetracking', 'eye_tracking')

FRAME_KEYS = ('frames', 'stim_vsync')

OPTOGENETIC_STIMULATION_KEYS = ('LED_sync', 'opto_trial')

PHOTODIODE_KEYS = ('photodiode', 'stim_photodiode')

analog_meta_data

close (*self*)

Closes the dataset.

duty_cycle (*self*, *line*)

Doesn't work right now. Freezes python for some reason.

Returns the duty cycle of a line.

frequency (*self*, *line*, *edge*='rising')

Returns the average frequency of a line.

get_all_bits (*self*)

Returns the data for all bits.

get_all_events (*self*)

Returns all counter values and their cooresponding IO state.

get_all_times (*self*, *units*='samples')

Returns all counter values.

Parameters

units [str] Return times in 'samples' or 'seconds'

get_analog_channel (*self*, *channel*, *start_time*=0.0, *stop_time*=None, *downsample*=1)

Returns the data from the specified analog channel between the timepoints.

Args: *channel* (int, str): desired channel index or label *start_time* (Optional[float]): start time in seconds
stop_time (Optional[float]): stop time in seconds *downsample* (Optional[int]): downsample factor

Returns: ndarray: slice of data for specified channel

Raises: KeyError: no analog data present

get_analog_meta (*self*)

Returns the metadata for the analog data.

get_bit (*self*, *bit*)

Returns the values for a specific bit.

Parameters

bit [int] Bit to return.

get_bit_changes (*self*, *bit*)

Returns the first derivative of a specific bit. Data points are 1 on rising edges and 255 on falling edges.

Parameters

bit [int] Bit for which to return changes.

get_edges (*self*, *kind*: str, *keys*: Union[str, Sequence[str]], *units*: str = 'seconds', *permissive*: bool = False) → Union[numpy.ndarray, NoneType]

Utility function for extracting edge times from a line

Parameters

kind [One of "rising", "falling", or "all". Should this method return] timestamps for rising, falling or both edges on the appropriate line

keys [These will be checked in sequence. Timestamps will be returned] for the first which is present in the line labels

units [one of "seconds", "samples", or "indices". The returned] "time"stamps will be given in these units.

raise_missing [If True and no matching line is found, a KeyError will] be raised

Returns

An array of edge times. If `raise_missing` is `False` and none of the keys were found, returns `None`.

Raises

KeyError [none of the provided keys were found among this dataset's] line labels

get_events_by_bit (*self*, *bit*, *units*='samples')

Returns all counter values for transitions (both rising and falling) for a specific bit.

Parameters

bit [int] Bit for which to return events.

get_events_by_line (*self*, *line*, *units*='samples')

Returns all counter values for transitions (both rising and falling) for a specific line.

Parameters

line [str] Line for which to return events.

get_falling_edges (*self*, *line*, *units*='samples')

Returns the counter values for the falling edges for a specific bit or line.

Parameters

line [str] Line for which to return edges.

get_line (*self*, *line*)

Returns the values for a specific line.

Parameters

line [str] Line to return.

get_line_changes (*self*, *line*)

Returns the first derivative of a specific line. Data points are 1 on rising edges and 255 on falling edges.

Parameters

line [(str)] Line name for which to return changes.

get_nearest (*self*, *source*, *target*, *source_edge*='rising', *target_edge*='rising', *direction*='previous', *units*='indices')

For all values of the source line, finds the nearest edge from the target line.

By default, returns the indices of the target edges.

Args: *source* (str, int): desired source line *target* (str, int): desired target line *source_edge* [Optional(str)]: "rising" or "falling" source edges *target_edge* [Optional(str)]: "rising" or "falling" target edges *direction* (str): "previous" or "next". Whether to prefer the

previous edge or the following edge.

units (str): "indices"

get_rising_edges (*self*, *line*, *units*='samples')

Returns the counter values for the rizing edges for a specific bit or line.

Parameters

line [str] Line for which to return edges.

line_stats (*self*, *line*, *print_results=True*)

Quick-and-dirty analysis of a bit.

##TODO: Split this up into smaller functions.

load (*self*, *path*)

Loads an hdf5 sync dataset.

Parameters

path [str] Path to hdf5 file.

period (*self*, *line*, *edge='rising'*)

Returns a dictionary with avg, min, max, and st of period for a line.

plot_all (*self*, *start_time*, *stop_time*, *auto_show=True*)

Plot all active bits.

Yikes. Come up with a better way to show this.

plot_bit (*self*, *bit*, *start_time=0.0*, *end_time=None*, *auto_show=True*, *axes=None*, *name=""*)

Plots a specific bit at a specific time period.

plot_bits (*self*, *bits*, *start_time=0.0*, *end_time=None*, *auto_show=True*)

Plots a list of bits.

plot_line (*self*, *line*, *start_time=0.0*, *end_time=None*, *auto_show=True*)

Plots a specific line at a specific time period.

plot_lines (*self*, *lines*, *start_time=0.0*, *end_time=None*, *auto_show=True*)

Plots specific lines at a specific time period.

sample_freq

stats (*self*)

Quick-and-dirty analysis of all bits. Prints a few things about each bit where events are found.

`allensdk.brain_observatory.sync_dataset.get_bit(uint_array, bit)`

Returns a bool array for a specific bit in a uint ndarray.

Parameters

uint_array [(numpy.ndarray)] The array to extract bits from.

bit [(int)] The bit to extract.

`allensdk.brain_observatory.sync_dataset.unpack_uint32(uint32_array, endian='L')`

Unpacks an array of 32-bit unsigned integers into bits.

Default is least significant bit first.

***Not currently used by sync dataset because get_bit is better and does** basically the same thing. I'm just leaving it in because it could potentially account for endianness and possibly have other uses in the future.

Module contents

```
class allensdk.brain_observatory.JSONEncoder(*, skipkeys=False, ensure_ascii=True,
                                             check_circular=True, allow_nan=True,
                                             sort_keys=False, indent=None, separa-
                                             tors=None, default=None)
```

Bases: `json.encoder.JSONEncoder`

default (*self*, *o*)

Implement this method in a subclass such that it returns a serializable object for *o*, or calls the base implementation (to raise a `TypeError`).

For example, to support arbitrary iterators, you could implement default like this:

```
def default(self, o):
    try:
        iterable = iter(o)
    except TypeError:
        pass
    else:
        return list(iterable)
    # Let the base class default method raise the TypeError
    return JSONEncoder.default(self, o)
```

`allensdk.brain_observatory.dict_to_indexed_array` (*dc*, *order=None*)

Given a dictionary and an ordered arr, build a concatenation of the dictionary's values and an index describing how that concatenation can be unpacked

`allensdk.brain_observatory.hook` (*json_dict*)

6.1.3 allensdk.config package

Subpackages

allensdk.config.app package

Submodules

allensdk.config.app.application_config module

```
class allensdk.config.app.application_config.ApplicationConfig(defaults,
                                                                name='app',
                                                                help='Run ap-
                                                                plication.', de-
                                                                fault_log_config=None)
```

Bases: `object`

Convenience class that handles of application configuration from environment variables, .conf files and the command line using Python standard libraries and formats.

apply_configuration_from_command_line (*self*, *parsed_args*)

Read application configuration variables from the command line.

Unassigned variables are left unchanged if previously assigned, set to their default values, or None if no default is specified at init time. Assigned variables will overwrite the previous value.

see: <https://docs.python.org/2/howto/argparse.html>

Parameters

parsed_args [dict] the arguments as parsed from the command line.

apply_configuration_from_environment (*self*)

Read application configuration variables from the environment.

The variable names are upper case and have a prefix defined by the application.

See: <https://docs.python.org/2/library/os.html>

apply_configuration_from_file (*self*, *config_file_path*)

Read application configuration variables from a .conf file.

Unassigned variables are set to their default values or None if no default is specified at init time. The variables are found in a section named by the application.

Parameters

config_file_path [string] path to to an INI (.conf) or JSON format application config file.

Returns

see: <https://docs.python.org/2/library/configparser.html>

create_argparser (*self*)

Initialization for the command-line parsing stage.

An application specific prefix is applied to argument names.

Parameters

prog [string] Application specific prefix for argument names.

description [string] A brief 'help' description of the application.

Returns

argParse.ArgumentParser The initialized argument parser object.

Notes

Defaults are set at the first environment reading. Command line args only override them when present

from_json_file (*self*, *json_path*)

Read an application configuration from a JSON format file.

Parameters

json_path [string] Path to the JSON file.

Returns

string An application configuration in INI format

from_json_string (*self*, *json_string*)

Read a configuration from a JSON format string.

Parameters

json_string [string] A JSON-formatted string containing an application configuration.

Returns

string An application configuration in INI format

load (*self*, *command_line_args*, *disable_existing_loggers=True*)

Load application configuration options, first from the environment, then from the configuration file, then from the command line.

Each stage of loading can override the previous stage.

Parameters

command_line_args [dict] Parameters passed to the application.

disable_existing_loggers [boolean] Reset the logging system or not.

Returns

fileConfig Configuration object with all levels applied

parse_command_line_args (*self*, *args*)

Simply call the internal argparser object.

Parameters

args [array] Parameters passed to the application.

Returns

Namespace Parsed parameters.

to_config_string (*self*, *description*)

Create a configuration string from a dict.

Parameters

description [dict] Configuration options for an application.

Returns

string Equivalent configuration as an INI format string

Notes

The Python configparser library natively supports this functionality in Python 3.

Module contents

allensdk.config.app is a package that assists in configuring application software, as opposed to domain-specific configuration.

allensdk.config.model package

Subpackages

allensdk.config.model.formats package

Submodules

allensdk.config.model.formats.hdf5_util module

```
class allensdk.config.model.formats.hdf5_util.Hdf5Util
    Bases: object

    read (self, file_path)

    write (self, file_path, m)
```

allensdk.config.model.formats.json_description_parser module

```
class allensdk.config.model.formats.json_description_parser.JsonDescriptionParser
    Bases: allensdk.config.model.description_parser.DescriptionParser

    log = <Logger allensdk.config.model.formats.json_description_parser (WARNING)>

    read (self, file_path, description=None, section=None, **kwargs)
        Parse a complete or partial configuration.
```

Parameters

json_string [string] Input to parse.

description [Description, optional] Where to put the parsed configuration. If None a new one is created.

section [string, optional] Where to put the parsed configuration within the description.

Returns

Description The input description with parsed configuration added.

Section is only specified for “bare” objects that are to be added to a section array.

```
read_string (self, json_string, description=None, section=None, **kwargs)
    Parse a complete or partial configuration.
```

Parameters

json_string [string] Input to parse.

description [Description, optional] Where to put the parsed configuration. If None a new one is created.

section [string, optional] Where to put the parsed configuration within the description.

Returns

Description The input description with parsed configuration added.

Section is only specified for “bare” objects that are to be added to a section array.

```
write (self, filename, description)
    Write the description to a JSON file.
```

Parameters

description [Description] Object to write.

```
write_string (self, description)
    Write the description to a JSON string.
```

Parameters

description [Description] Object to write.

Returns

string JSON serialization of the input.

allensdk.config.model.formats.pycfg_description_parser module

class allensdk.config.model.formats.pycfg_description_parser.**PycfgDescriptionParser**
 Bases: *allensdk.config.model.description_parser.DescriptionParser*

log = <Logger allensdk.config.model.formats.pycfg_description_parser (WARNING)>

read (*self*, *pycfg_file_path*, *description=None*, *section=None*, ***kwargs*)
 Read a serialized description from a Python (.pycfg) file.

Parameters

filename [string] Name of the .pycfg file.

Returns

Description Configuration object.

read_string (*self*, *python_string*, *description=None*, *section=None*, ***kwargs*)
 Read a serialized description from a Python (.pycfg) string.

Parameters

python_string [string] Python string with a serialized description.

Returns

Description Configuration object.

write (*self*, *filename*, *description*)
 Write the description to a Python (.pycfg) file.

Parameters

filename [string] Name of the file to write.

write_string (*self*, *description*)
 Write the description to a pretty-printed Python string.

Parameters

description [Description] Configuration object to write.

Module contents**Submodules****allensdk.config.model.description module**

class allensdk.config.model.description.**Description**
 Bases: object

fix_unary_sections (*self*, *section_names=None*)
 Wrap section contents that don't have the proper array surrounding them in an array.

Parameters

section_names [list of strings, optional] Keys of sections that might not be in array form.

is_empty (*self*)

Check if anything is in the object.

Returns

boolean true if self.data is missing or empty

unpack (*self*, *data*, *section=None*)

Read the manifest and other stand-alone configuration structure, or insert a configuration object into a section of an existing configuration.

Parameters

data [dict] A configuration object including top level sections, or an configuration object to be placed within a section.

section [string, optional.] If this is present, place data within an existing section array.

unpack_manifest (*self*, *data*)

Pull the manifest configuration section into a separate place.

Parameters

data [dict] A configuration structure that still has a manifest section.

update_data (*self*, *data*, *section=None*)

Merge configuration data possibly from multiple files.

Parameters

data [dict] Configuration structure to add.

section [string, optional] What configuration section to read it into if the file does not specify.

allensdk.config.model.description_parser module

class allensdk.config.model.description_parser.**DescriptionParser**

Bases: object

log = <Logger allensdk.config.model.description_parser (WARNING)>

parser_for_extension (*self*, *filename*)

Choose a subclass that can read the format.

Parameters

filename [string] For the extension.

Returns

DescriptionParser Appropriate subclass.

read (*self*, *file_path*, *description=None*, *section=None*, ***kwargs*)

Parse data needed for a simulation.

Parameters

description [dict] Configuration from parsing previous files.

section [string, optional] What configuration section to read it into if the file does not specify.

read_string (*self*, *data_string*, *description=None*, *section=None*, *header=None*)

Parse data needed for a simulation from a string.

write (*self*, *filename*, *description*)
Save the configuration.

Parameters

filename [string] Name of the file to write.

Module contents

Submodules

allensdk.config.manifest module

class allensdk.config.manifest.**Manifest** (*config=None*, *relative_base_dir='.'*, *version=None*)

Bases: object

Manages the location of external files referenced in an Allen SDK configuration

DIR = 'dir'

DIRNAME = 'dir_name'

FILE = 'file'

VERSION = 'manifest_version'

add_file (*self*, *file_key*, *file_name*, *dir_key=None*, *path_format=None*)
Insert a new file entry.

Parameters

file_key [string] Reference to the entry.

file_name [string] Substitutions of the %s, %d style allowed.

dir_key [string] Reference to the parent directory entry.

path_format [string, optional] File type for further parsing.

add_path (*self*, *key*, *path*, *path_type='dir'*, *absolute=True*, *path_format=None*, *parent_key=None*)
Insert a new entry.

Parameters

key [string] Identifier for referencing the entry.

path [string] Specification for a path using %s, %d style substitution.

path_type [string enumeration] 'dir' (default) or 'file'

absolute [boolean] Is the spec relative to the process current directory.

path_format [string, optional] Indicate a known file type for further parsing.

parent_key [string] Refer to another entry.

add_paths (*self*, *path_info*)
add information about paths stored in the manifest.

Parameters

path_info [dict] Information about the new paths

as_dataframe (*self*)

check_dir (*self*, *path_key*, *do_exit=False*)

Verify a directories existence or optionally exit.

Parameters

path_key [string] Reference to the entry.

do_exit [boolean] What to do if the directory is not present.

create_dir (*self*, *path_key*)

Make a directory for an entry.

Parameters

path_key [string] Reference to the entry.

get_format (*self*, *path_key*)

Retrieve the type of a path entry.

Parameters

path_key [string] reference to the entry

Returns

string File type.

get_path (*self*, *path_key*, **args*)

Retrieve an entry with substitutions.

Parameters

path_key [string] Refer to the entry to retrieve.

args [any types, optional] arguments to be substituted into the path spec for %s, %d, etc.

Returns

string Path with parent structure and substitutions applied.

load_config (*self*, *config*, *version=None*)

Load paths into the manifest from an Allen SDK config section.

Parameters

config [Config] Manifest section of an Allen SDK config.

log = <Logger allensdk.config.manifest (WARNING)>

resolve_paths (*self*, *description_dict*, *suffix='_key'*)

Walk input items and expand those that refer to a manifest entry.

Parameters

description_dict [dict] Any entries with key names ending in suffix will be expanded.

suffix [string] Indicates the entries to be expanded.

classmethod safe_make_parent_dirs (*file_name*)

Create a parent directories for file.

Parameters

file_name [string]

Returns

leftmost [string] most rootward directory created

classmethod `safe_mkdir(directory)`

Create path if not already there.

Parameters

directory [string] create it if it doesn't exist

Returns

leftmost [string] most rootward directory created

exception `allensdk.config.manifest.ManifestVersionError(message, version, found_version)`

Bases: Exception

outdated

allensdk.config.manifest_builder module

class `allensdk.config.manifest_builder.ManifestBuilder`

Bases: object

add_path (*self*, *key*, *spec*, *typename='dir'*, *parent_key=None*, *format=None*)

add_section (*self*, *name*, *contents*)

as_dataframe (*self*)

df_columns = ['key', 'parent_key', 'spec', 'type', 'format']

from_dataframe (*self*, *df*)

get_config (*self*)

get_manifest (*self*)

set_version (*self*, *value*)

write_json_file (*self*, *path*, *overwrite=False*)

write_json_string (*self*)

Module contents

`allensdk.config.enable_console_log(level=None)`

configure allensdk logging to output to the console.

Parameters

level [int] logging level 0-50 (logging.INFO, logging.DEBUG, etc.)

Notes

See: [Logging Cookbook](#)

6.1.4 allensdk.core package

Subpackages

allensdk.core.lazy_property package

Submodules

allensdk.core.lazy_property.lazy_property module

```
class allensdk.core.lazy_property.lazy_property.LazyProperty (api_method, wrap-  
                                                         pers=(), *args,  
                                                         **kwargs)  
  
    Bases: object  
    calculate (self)
```

allensdk.core.lazy_property.lazy_property_mixin module

```
class allensdk.core.lazy_property.lazy_property_mixin.LazyPropertyMixin  
    Bases: object  
    LazyProperty
```

Module contents

Submodules

allensdk.core.auth_config module

allensdk.core.authentication module

```
class allensdk.core.authentication.CredentialProvider  
    Bases: abc.ABC  
    METHOD = 'custom'  
    provide (self, credential)  
  
class allensdk.core.authentication.DbCredentials (dbname, user, host, port, password)  
    Bases: tuple  
    dbname  
        Alias for field number 0  
    host  
        Alias for field number 2  
    password  
        Alias for field number 4  
    port  
        Alias for field number 3
```


user

Alias for field number 1

```
class allensdk.core.authentication.EnvCredentialProvider(environ: Optional[Dict[str, Any]] = None)
```

Bases: `allensdk.core.authentication.CredentialProvider`

Provides credentials from environment variables for variables listed in CREDENTIAL_KEYS.

METHOD = 'env'

provide (*self, credential*)

```
allensdk.core.authentication.credential_injector(credential_map: Dict[str, Any], provider: Union[allensdk.core.authentication.CredentialProvider, NoneType] = None)
```

Decorator used to inject credentials from another source if not explicitly provided in the function call. This function will only supply values for keyword arguments. All keys defined in *credential_map* must correspond to keyword arguments in the function signature.

Parameters

credential_map: Dict[Str: Any] Dictionary where the keys are the keyword of a credential kwarg passed to the decorated function, and the values are the name of the credential in the credential provider (see CREDENTIAL_KEYS).

Example of *credential_map* for PostgresQueryMixin connecting to LIMS database:

```
{ "dbname": "LIMS_DBNAME", "user": "LIMS_USER", "host": "LIMS_HOST", "password": "LIMS_PASSWORD", "port": "LIMS_PORT" }
```

provider: Optional[CredentialProvider] Subclass of CredentialProvider to provide credentials to the wrapped function. If left unspecified, will default to EnvCredentialProvider, which provides credentials from environment variables.

```
allensdk.core.authentication.get_credential_provider()
```

```
allensdk.core.authentication.set_credential_provider(provider)
```

allensdk.core.brain_observatory_cache module

```
class allensdk.core.brain_observatory_cache.BrainObservatoryCache(cache=True, manifest_file=None, base_uri=None, api=None)
```

Bases: `allensdk.api.cache.Cache`

Cache class for storing and accessing data from the Brain Observatory. By default, this class will cache any downloaded metadata or files in well known locations defined in a manifest file. This behavior can be disabled.

Parameters

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. `get_ophys_experiment_data(file_name='file.nwb')`).

manifest_file: string File name of the manifest to be read. Default is "brain_observatory_manifest.json".

Attributes

api: BrainObservatoryApi instance The object used for making API queries related to the Brain Observatory.

```
ANALYSIS_DATA_KEY = 'ANALYSIS_DATA'
CELL_SPECIMENS_KEY = 'CELL_SPECIMENS'
EVENTS_DATA_KEY = 'EVENTS_DATA'
EXPERIMENTS_KEY = 'EXPERIMENTS'
EXPERIMENT_CONTAINERS_KEY = 'EXPERIMENT_CONTAINERS'
EXPERIMENT_DATA_KEY = 'EXPERIMENT_DATA'
EYE_GAZE_DATA_KEY = 'EYE_GAZE_DATA'
MANIFEST_VERSION = '1.3'
STIMULUS_MAPPINGS_KEY = 'STIMULUS_MAPPINGS'
```

build_manifest (*self*, *file_name*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

get_all_cre_lines (*self*)

Return a list of all cre driver lines in the data set.

get_all_imaging_depths (*self*)

Return a list of all imaging depths in the data set.

get_all_reporter_lines (*self*)

Return a list of all reporter lines in the data set.

get_all_session_types (*self*)

Return a list of all stimulus sessions in the data set.

get_all_stimuli (*self*)

Return a list of all stimuli in the data set.

get_all_targeted_structures (*self*)

Return a list of all targeted structures in the data set.

get_cell_specimens (*self*, *file_name=None*, *ids=None*, *experiment_container_ids=None*, *include_failed=False*, *simple=True*, *filters=None*)

Return cell specimens that have certain properties.

Parameters

file_name: string File name to save/read the cell specimens. If *file_name* is *None*, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is *None*.

ids: list List of cell specimen ids.

experiment_container_ids: list List of experiment container ids.

include_failed: bool Whether to include cells from failed experiment containers

simple: boolean Whether or not to simplify the dictionary properties returned by this method to a more concise subset.

filters: list of dicts List of filter dictionaries. The Allen Brain Observatory web site can generate filters in this format to reproduce a filtered set of cells found there. To see what these look like, visit <http://observatory.brain-map.org/visualcoding>, perform a cell search and apply some filters (e.g. find cells in a particular area), then click the “view these cells in the AllenSDK” link on the bottom-left of the search results page. This will take you to a page that contains a code sample you can use to apply those same filters via this argument. For more detail on the filter syntax, see `BrainObservatoryApi.dataframe_query`.

Returns

list of dictionaries

get_experiment_containers (*self*, *file_name=None*, *ids=None*, *targeted_structures=None*, *imaging_depths=None*, *cre_lines=None*, *reporter_lines=None*, *transgenic_lines=None*, *include_failed=False*, *simple=True*)

Get a list of experiment containers matching certain criteria.

Parameters

file_name: string File name to save/read the experiment containers. If *file_name* is *None*, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is *None*.

ids: list List of experiment container ids.

targeted_structures: list List of structure acronyms. Must be in the list returned by `BrainObservatoryCache.get_all_targeted_structures()`.

imaging_depths: list List of imaging depths. Must be in the list returned by `BrainObservatoryCache.get_all_imaging_depths()`.

cre_lines: list List of cre lines. Must be in the list returned by `BrainObservatoryCache.get_all_cre_lines()`.

reporter_lines: list List of reporter lines. Must be in the list returned by `BrainObservatoryCache.get_all_reporter_lines()`.

transgenic_lines: list List of transgenic lines. Must be in the list returned by `BrainObservatoryCache.get_all_cre_lines()` or `BrainObservatoryCache.get_all_reporter_lines()`.

include_failed: boolean Whether or not to include failed experiment containers.

simple: boolean Whether or not to simplify the dictionary properties returned by this method to a more concise subset.

Returns

list of dictionaries

get_nwb_filepath (*self*, *ophys_experiment_id=None*)

get_ophys_experiment_analysis (*self*, *ophys_experiment_id*, *stimulus_type*, *file_name=None*)

Download the h5 analysis file for a stimulus set, for a particular *ophys_experiment* (if it hasn’t already been downloaded) and return a data accessor object.

Parameters

file_name: string File name to save/read the data set. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

ophys_experiment_id: int id of the ophys_experiment to retrieve

stimulus_name: str stimulus type; should be an element of self.list_stimuli()

Returns

BrainObservatoryNwbDataSet

get_ophys_experiment_data (*self, ophys_experiment_id, file_name=None*)

Download the NWB file for an ophys_experiment (if it hasn't already been downloaded) and return a data accessor object.

Parameters

file_name: string File name to save/read the data set. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

ophys_experiment_id: integer id of the ophys_experiment to retrieve

Returns

BrainObservatoryNwbDataSet

get_ophys_experiment_events (*self, ophys_experiment_id, file_name=None*)

Download the npz events file for an ophys_experiment if it hasn't already been downloaded and return the events array.

Parameters

file_name: string File name to save/read the data set. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

ophys_experiment_id: int id of the ophys_experiment to retrieve events for

Returns

events: numpy.ndarray [N_cells,N_times] array of events.

get_ophys_experiment_stimuli (*self, experiment_id*)

For a single experiment, return the list of stimuli present in that experiment.

get_ophys_experiments (*self, file_name=None, ids=None, experiment_container_ids=None, targeted_structures=None, imaging_depths=None, cre_lines=None, reporter_lines=None, transgenic_lines=None, stimuli=None, session_types=None, cell_specimen_ids=None, include_failed=False, require_eye_tracking=False, simple=True*)

Get a list of ophys experiments matching certain criteria.

Parameters

file_name: string File name to save/read the ophys experiments. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

ids: list List of ophys experiment ids.

experiment_container_ids: list List of experiment container ids.

targeted_structures: list List of structure acronyms. Must be in the list returned by `BrainObservatoryCache.get_all_targeted_structures()`.

imaging_depths: list List of imaging depths. Must be in the list returned by `BrainObservatoryCache.get_all_imaging_depths()`.

cre_lines: list List of cre lines. Must be in the list returned by `BrainObservatoryCache.get_all_cre_lines()`.

reporter_lines: list List of reporter lines. Must be in the list returned by `BrainObservatoryCache.get_all_reporter_lines()`.

transgenic_lines: list List of transgenic lines. Must be in the list returned by `BrainObservatoryCache.get_all_cre_lines()` or `BrainObservatoryCache.get_all_reporter_lines()`.

stimuli: list List of stimulus names. Must be in the list returned by `BrainObservatoryCache.get_all_stimuli()`.

session_types: list List of stimulus session type names. Must be in the list returned by `BrainObservatoryCache.get_all_session_types()`.

cell_specimen_ids: list Only include experiments that contain cells with these ids.

include_failed: boolean Whether or not to include experiments from failed experiment containers.

simple: boolean Whether or not to simplify the dictionary properties returned by this method to a more concise subset.

require_eye_tracking: boolean If True, only return experiments that have eye tracking results. Default: False.

Returns

list of dictionaries

get_ophys_pupil_data (*self*, *ophys_experiment_id: int*, *file_name: str = None*, *suppress_pupil_data: bool = True*) → `pandas.core.frame.DataFrame`

Download the h5 eye gaze mapping file for an ophys_experiment if it hasn't already been downloaded and return it as a `pandas.DataFrame`.

Parameters

file_name: string File name to save/read the data set. If *file_name* is None, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

ophys_experiment_id: int id of the ophys_experiment to retrieve pupil data for.

suppress_pupil_data: bool Whether or not to suppress pupil data from dataset. Default is True.

Returns

`pd.DataFrame`

If 'suppress_eye_gaze_data' is set to 'False':

Contains raw/filtered columns for gaze mapping: **_eye_area*

**_pupil_area* **_screen_coordinates_x_cm* **_screen_coordinates_y_cm*

**_screen_coordinates_spherical_x_deg* **_screen_coordinates_spherical_y_deg*

Otherwise: An empty pandas DataFrame

allensdk.core.brain_observatory_nwb_data_set module

```
class allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet (nwb_file)
    Bases: object

    FILE_METADATA_MAPPING = {'age': 'general/subject/age', 'device_string': 'general/dev
    MOTION_CORRECTION_DATASETS = ['MotionCorrection/2p_image_series/xy_translations', 'Mot
    PIPELINE_DATASET = 'brain_observatory_pipeline'
    STIMULUS_TABLE_TYPES = {'abstract_feature_series': ['drifting_gratings', 'static_grat
    SUPPORTED_PIPELINE_VERSION = '3.0'

    get_cell_specimen_ids (self)
        Returns an array of cell IDs for all cells in the file

        Returns

            cell specimen IDs: list

    get_cell_specimen_indices (self, cell_specimen_ids)
        Given a list of cell specimen ids, return their index based on their order in this file.

        Parameters

            cell_specimen_ids: list of cell specimen ids

    get_corrected_fluorescence_traces (self, cell_specimen_ids=None)
        Returns an array of demixed and neuropil-corrected fluorescence traces for all ROIs and the timestamps
        for each datapoint

        Parameters

            cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is
                None (default) then all are returned

        Returns

            timestamps: 2D numpy array Timestamp for each fluorescence sample
            traces: 2D numpy array Corrected fluorescence traces for each cell

    get_demixed_traces (self, cell_specimen_ids=None)
        Returns an array of demixed fluorescence traces for all ROIs and the timestamps for each datapoint

        Parameters

            cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is
                None (default) then all are returned

        Returns

            timestamps: 2D numpy array Timestamp for each fluorescence sample
            traces: 2D numpy array Demixed fluorescence traces for each cell

    get_dff_traces (self, cell_specimen_ids=None)
        Returns an array of dF/F traces for all ROIs and the timestamps for each datapoint

        Parameters

            cell_specimen_ids: list or array (optional) List of cell IDs to return data for. If this is
                None (default) then all are returned

        Returns
```

timestamps: 2D numpy array Timestamp for each fluorescence sample

dF/F: 2D numpy array dF/F values for each cell

get_fluorescence_timestamps (*self*)

Returns an array of timestamps in seconds for the fluorescence traces

get_fluorescence_traces (*self*, *cell_specimen_ids=None*)

Returns an array of fluorescence traces for all ROI and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Fluorescence traces for each cell

get_locally_sparse_noise_stimulus_template (*self*, *stimulus*, *mask_off_screen=True*)

Return an array of the stimulus template for the specified stimulus.

Parameters

stimulus: string

Which locally sparse noise stimulus to retrieve. Must be one of:

stimulus_info.LOCALLY_SPARSE_NOISE	stimu-
lus_info.LOCALLY_SPARSE_NOISE_4DEG	stimu-
lus_info.LOCALLY_SPARSE_NOISE_8DEG	

mask_off_screen: boolean Set off-screen regions of the stimulus to LocallySparseNoise.LSN_OFF_SCREEN.

Returns

tuple: (template, off-screen mask)

get_max_projection (*self*)

Returns the maximum projection image for the 2P movie.

Returns

max projection: np.ndarray

get_metadata (*self*)

Returns a dictionary of meta data associated with each experiment, including Cre line, specimen number, visual area imaged, imaging depth

Returns

metadata: dictionary

get_motion_correction (*self*)

Returns a Panda DataFrame containing the x- and y- translation of each image used for image alignment

get_neuropil_r (*self*, *cell_specimen_ids=None*)

Returns a scalar value of r for neuropil correction of fluorescence traces

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then results for all are returned

Returns

r: 1D numpy array, `len(r)=len(cell_specimen_ids)` Scalar for neuropil subtraction for each cell

get_neuropil_traces (*self*, *cell_specimen_ids=None*)

Returns an array of neuropil fluorescence traces for all ROIs and the timestamps for each datapoint

Parameters

cell_specimen_ids: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Neuropil fluorescence traces for each cell

get_pupil_location (*self*, *as_spherical=True*)

Returns the x, y pupil location.

Parameters

as_spherical [bool] Whether to return the location as spherical (default) or not. If true, the result is altitude and azimuth in degrees, otherwise it is x, y in centimeters. (0,0) is the center of the monitor.

Returns

(timestamps, location) Timestamps is an (Nx1) array of timestamps in seconds. Location is an (Nx2) array of spatial location.

get_pupil_size (*self*)

Returns the pupil area in pixels.

Returns

(timestamps, areas) Timestamps is an (Nx1) array of timestamps in seconds. Areas is an (Nx1) array of pupil areas in pixels.

get_roi_ids (*self*)

Returns an array of IDs for all ROIs in the file

Returns

ROI IDs: list

get_roi_mask (*self*, *cell_specimen_ids=None*)

Returns an array of all the ROI masks

Parameters

cell specimen IDs: list or array (optional) List of cell IDs to return traces for. If this is None (default) then all are returned

Returns

List of ROI_Mask objects

get_roi_mask_array (*self*, *cell_specimen_ids=None*)

Return a numpy array containing all of the ROI masks for requested cells. If *cell_specimen_ids* is omitted, return all masks.

Parameters

cell_specimen_ids: list List of cell specimen ids. Default None.

Returns

np.ndarray: NxWxH array, where N is number of cells

get_running_speed (*self*)

Returns the mouse running speed in cm/s

get_session_type (*self*)

Returns the type of experimental session, presently one of the following: three_session_A, three_session_B, three_session_C

Returns

session type: string

get_stimulus (*self*, *frame_ind*)

get_stimulus_epoch_table (*self*)

Returns a pandas dataframe that summarizes the stimulus epoch duration for each acquisition time index in the experiment

Parameters

None

Returns

timestamps: 2D numpy array Timestamp for each fluorescence sample

traces: 2D numpy array Fluorescence traces for each cell

get_stimulus_table (*self*, *stimulus_name*)

Return a stimulus table given a stimulus name

Notes

For more information, see: http://help.brain-map.org/display/observatory/Documentation?preview=/10616846/10813485/VisualCoding_VisualStimuli.pdf

get_stimulus_template (*self*, *stimulus_name*)

Return an array of the stimulus template for the specified stimulus.

Parameters

stimulus_name: string Must be one of the strings returned by list_stimuli().

Returns

stimulus table: pd.DataFrame

list_stimuli (*self*)

Return a list of the stimuli presented in the experiment.

Returns

stimuli: list of strings

number_of_cells

Number of cells in the experiment

save_analysis_arrays (*self*, **datasets*)

save_analysis_dataframes (*self*, **tables*)

stimulus_search

```
allensdk.core.brain_observatory_nwb_data_set.align_running_speed(dxcm, dx-  
time, times-  
tamps)
```

If running speed timestamps differ from fluorescence timestamps, adjust by inserting NaNs to running speed.

Returns

tuple: dxcm, **dx**time

```
allensdk.core.brain_observatory_nwb_data_set.get_epoch_mask_list(st, threshold,  
max_cuts=2)
```

Convenience function to cut a stim table into multiple epochs

Parameters

- **st** – input stimtable
- **threshold** – threshold on the max duration of a subepoch
- **max_cuts** – maximum number of allowed epochs to cut into

Returns epoch_mask_list, a list of indices that define the start and end of sub-epochs

allensdk.core.cache_method_utilities module

```
class allensdk.core.cache_method_utilities.CachedInstanceMethodMixin
```

Bases: object

```
cache_clear(self)
```

Calls *cache_clear* method on all bound methods in this instance (where valid). Intended to clear calls cached with the *memoize* decorator. Note that this will also clear functions decorated with *lru_cache* and *lfu_cache* in this class (or any other function with *cache_clear* attribute).

allensdk.core.cell_types_cache module

```
class allensdk.core.cell_types_cache.CellTypesCache(cache=True, mani-  
fest_file=None, base_uri=None)
```

Bases: *allensdk.api.cache.Cache*

Cache class for storing and accessing data from the Cell Types Database. By default, this class will cache any downloaded metadata or files in well known locations defined in a manifest file. This behavior can be disabled.

Parameters

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. *get_ephys_data*(file_name='file.nwb')).

manifest_file: string File name of the manifest to be read. Default is "cell_types_manifest.json".

Attributes

api: CellTypesApi instance The object used for making API queries related to the Cell Types Database

```
CELLS_KEY = 'CELLS'
```

```
EPHYS_DATA_KEY = 'EPHYS_DATA'
```

```
EPHYS_FEATURES_KEY = 'EPHYS_FEATURES'
```

```
EPHYS_SWEEPS_KEY = 'EPHYS_SWEEPS'
```

```
MANIFEST_VERSION = '1.1'
```

```
MARKER_KEY = 'MARKER'
```

```
MORPHOLOGY_FEATURES_KEY = 'MORPHOLOGY_FEATURES'
```

```
RECONSTRUCTION_KEY = 'RECONSTRUCTION'
```

```
build_manifest(self, file_name)
```

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

```
get_all_features(self, dataframe=False, require_reconstruction=True)
```

Download morphology and electrophysiology features for all cells and merge them into a single table.

Parameters

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

require_reconstruction: boolean Only return ephys and morphology features for cells that have reconstructions. Default True.

```
get_cells(self, file_name=None, require_morphology=False, require_reconstruction=False, reporter_status=None, species=None, simple=True)
```

Download metadata for all cells in the database and optionally return a subset filtered by whether or not they have a morphology or reconstruction.

Parameters

file_name: string File name to save/read the cell metadata as JSON. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

require_morphology: boolean Filter out cells that have no morphological images.

require_reconstruction: boolean Filter out cells that have no morphological reconstructions.

reporter_status: list Filter for cells that have one or more cell reporter statuses.

species: list Filter for cells that belong to one or more species. If None, return all. Must be one of [CellTypesApi.MOUSE, CellTypesApi.HUMAN].

```
get_ephys_data(self, specimen_id, file_name=None)
```

Download electrophysiology traces for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

NwbDataSet A class instance with helper methods for retrieving stimulus and response traces out of an NWB file.

```
get_ephys_features(self, dataframe=False, file_name=None)
```

Download electrophysiology features for all cells in the database.

Parameters

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

get_ephys_sweeps (*self, specimen_id, file_name=None*)

Download sweep metadata for a single cell specimen.

Parameters

specimen_id: int ID of a cell.

get_morphology_features (*self, dataframe=False, file_name=None*)

Download morphology features for all cells with reconstructions in the database.

Parameters

file_name: string File name to save/read the ephys features metadata as CSV. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

dataframe: boolean Return the output as a Pandas DataFrame. If False, return a list of dictionaries.

get_reconstruction (*self, specimen_id, file_name=None*)

Download and open a reconstruction for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the reconstruction SWC. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

Morphology A class instance with methods for accessing morphology compartments.

get_reconstruction_markers (*self, specimen_id, file_name=None*)

Download and open a reconstruction marker file for a single cell in the database.

Parameters

specimen_id: int The ID of a cell specimen to download.

file_name: string File name to save/read the reconstruction marker. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

Returns

Morphology A class instance with methods for accessing morphology compartments.

class allensdk.core.cell_types_cache.**ReporterStatus**

Bases: object

Valid strings for filtering by cell reporter status.

INDETERMINATE = None

NA = None

```
NEGATIVE = 'negative'
```

```
POSITIVE = 'positive'
```

allensdk.core.dat_utilities module

```
class allensdk.core.dat_utilities.DatUtilities
```

```
Bases: object
```

```
classmethod save_voltage(output_path, v, t)
```

Save a single voltage output result into a simple text format.

The output file is one t v pair per line.

Parameters

output_path [string] file name for output

v [numpy array] voltage

t [numpy array] time

allensdk.core.exceptions module

```
exception allensdk.core.exceptions.DataFrameIndexError(msg,
                                                         caught_exception=None)
```

```
Bases: LookupError
```

More verbose method for accessing invalid rows or columns in a dataframe. Should be used when an index error is thrown on a dataframe.

```
exception allensdk.core.exceptions.DataFrameKeyError(msg, caught_exception=None)
```

```
Bases: LookupError
```

More verbose method for accessing invalid rows or columns in a dataframe. Should be used when a keyerror is thrown on a dataframe.

```
exception allensdk.core.exceptions.MissingDataError
```

```
Bases: ValueError
```

allensdk.core.h5_utilities module

```
allensdk.core.h5_utilities.decode_bytes(bytes_dataset, encoding='UTF-8')
```

Convert the elements of a dataset of bytes to str

```
allensdk.core.h5_utilities.h5_object_matcher_relname_in(relnames,
                                                         h5_object_name,
                                                         h5_object)
```

Asks if an h5 object's relative name (the final section of its absolute name) is contained within a provided array

Parameters

relnames [array-like] Relative names against which to match

h5_object_name [str] Full name (path from origin) of h5 object

h5_object [h5py.Group, h5py.Dataset] Check this object's relative name

Returns

bool : whether the match succeeded

h5_object [h5py.group, h5py.Dataset] the argued object

`allensdk.core.h5_utilities.keyed_locate_h5_objects` (*matcher_cbs*, *h5_file*,
start_node=None)

Traverse an h5 file and build up a dictionary mapping supplied keys to located objects

`allensdk.core.h5_utilities.load_datasets_by_relnames` (*relnames*, *h5_file*, *start_node*)

A convenience function for finding and loading into memory one or more datasets from an h5 file

`allensdk.core.h5_utilities.locate_h5_objects` (*matcher_cb*, *h5_file*, *start_node=None*)

Traverse an h5 file and return objects matching supplied criteria

`allensdk.core.h5_utilities.traverse_h5_file` (*callback*, *h5_file*, *start_node=None*)

Traverse an h5 file and apply a callback to each node

allensdk.core.json_utilities module

class `allensdk.core.json_utilities.JsonComments`

Bases: object

classmethod `read_file` (*file_name*)

classmethod `read_string` (*json_string*)

classmethod `remove_comments` (*json_string*)

Strip single and multiline javascript-style comments.

Parameters

json [string] Json string with javascript-style comments.

Returns

string Copy of the input with comments removed.

Note: A JSON decoder MAY accept and ignore comments.

classmethod `remove_multiline_comments` (*json_string*)

Rebuild input without substrings matching `/.../`.

Parameters

json_string [string] may or may not contain multiline comments.

Returns

string Copy of the input without the comments.

`allensdk.core.json_utilities.json_handler` (*obj*)

Used by `write_json` convert a few non-standard types to things that the json package can handle.

`allensdk.core.json_utilities.read` (*file_name*)

Shortcut reading JSON from a file.

`allensdk.core.json_utilities.read_url` (*url*, *method='POST'*)

`allensdk.core.json_utilities.read_url_get` (*url*)

Transform a JSON contained in a file into an equivalent nested python dict.

Parameters

url [string] where to get the json.

Returns

dict Python version of the input

**Note: if the input is a bare array or literal, for example,
the output will be of the corresponding type.**

```
allensdk.core.json_utilities.read_url_post(url)
```

Transform a JSON contained in a file into an equivalent nested python dict.

Parameters

url [string] where to get the json.

Returns

dict Python version of the input

**Note: if the input is a bare array or literal, for example,
the output will be of the corresponding type.**

```
allensdk.core.json_utilities.write(file_name, obj)
```

Shortcut for writing JSON to a file. This also takes care of serializing numpy and data types.

```
allensdk.core.json_utilities.write_string(obj)
```

Shortcut for writing JSON to a string. This also takes care of serializing numpy and data types.

allensdk.core.mouse_connectivity_cache module

```
class allensdk.core.mouse_connectivity_cache.MouseConnectivityCache(resolution=None,
                                                                    cache=True,
                                                                    manifest_file=None,
                                                                    ccf_version=None,
                                                                    base_uri=None,
                                                                    version=None)
```

Bases: `allensdk.core.reference_space_cache.ReferenceSpaceCache`

Cache class for storing and accessing data related to the adult mouse Connectivity Atlas. By default, this class will cache any downloaded metadata or files in well known locations defined in a manifest file. This behavior can be disabled.

Parameters

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

ccf_version: string Desired version of the Common Coordinate Framework. This affects the annotation volume (`get_annotation_volume`) and structure masks (`get_structure_mask`). Must be one of (`MouseConnectivityApi.CCF_2015`, `MouseConnectivityApi.CCF_2016`). Default: `MouseConnectivityApi.CCF_2016`

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. `get_projection_density(file_name='file.nrrd')`).

manifest_file: string File name of the manifest to be read. Default is "mouse_connectivity_manifest.json".

Attributes

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

api: MouseConnectivityApi instance Used internally to make API queries.

ALIGNMENT3D_KEY = 'ALIGNMENT3D'

DATA_MASK_KEY = 'DATA_MASK'

DEFAULT_STRUCTURE_SET_IDS = (167587189,)

DEFORMATION_FIELD_HEADER_KEY = 'DEFORMATION_FIELD_HEADER'

DEFORMATION_FIELD_VOXEL_KEY = 'DEFORMATION_FIELD_VOXELS'

DFMFLD_RESOLUTIONS = (25,)

EXPERIMENTS_KEY = 'EXPERIMENTS'

INJECTION_DENSITY_KEY = 'INJECTION_DENSITY'

INJECTION_FRACTION_KEY = 'INJECTION_FRACTION'

MANIFEST_VERSION = 1.3

PROJECTION_DENSITY_KEY = 'PROJECTION_DENSITY'

STRUCTURE_UNIONIZES_KEY = 'STRUCTURE_UNIONIZES'

SUMMARY_STRUCTURE_SET_ID = 167587189

add_manifest_paths (*self*, *manifest_builder*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

default_structure_ids

filter_experiments (*self*, *experiments*, *cre=None*, *injection_structure_ids=None*)

Take a list of experiments and filter them by cre status and injection structure.

Parameters

cre: boolean or list If True, return only cre-positive experiments. If False, return only cre-negative experiments. If None, return all experiments. If list, return all experiments with cre line names in the supplied list. Default None.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If None, return all experiments. Default None.

filter_structure_unionizes (*self*, *unionizes*, *is_injection=None*, *structure_ids=None*, *include_descendants=False*, *hemisphere_ids=None*)

Take a list of unionizes and return a subset of records filtered by injection status, structure, and hemisphere.

Parameters

is_injection: boolean If True, only return unionize records that disregard non-injection pixels. If False, only return unionize records that disregard injection pixels. If None, return all records. Default None.

structure_ids: list Only return unionize records for a set of structures. If None, return all records. Default None.

include_descendants: **boolean** Include all descendant records for specified structures. Default False.

hemisphere_ids: **list** Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If None, include all records [1, 2, 3]. Default None.

get_affine_parameters (*self*, *section_data_set_id*, *direction='trv'*, *file_name=None*)

Extract the parameters of the 3D affine tranformation mapping this section data set's image-space stack to CCF-space (or vice-versa).

Parameters

section_data_set_id [int] download the parameters for this data set.

direction [str, optional]

Valid options are:

trv ["transform from reference to volume". Maps CCF points to image space points. If you are] resampling data into CCF, this is the direction you want.

tvr : "transform from volume to reference". Maps image space points to CCF points.

file_name [str] If provided, store the downloaded file here.

Returns

alignment [numpy.ndarray]

4 X 3 matrix. In order to transform a point [X_1, X_2, X_3] run `np.dot([X_1, X_2, X_3, 1], alignment)`. In

to build a SimpleITK affine transform run: `transform = sitk.AffineTransform(3)`
`transform.SetParameters(alignment.flatten())`

get_data_mask (*self*, *experiment_id*, *file_name=None*)

Read a data mask volume for a single experiment. Download it first if it doesn't exist. Data mask is a binary mask of voxels that have valid data. Only use valid data in analysis!

Parameters

experiment_id: **int** ID of the experiment to download/read. This corresponds to `section_data_set_id` in the API.

file_name: **string** File name to store the template volume. If it already exists, it will be read from this file. If `file_name` is None, the `file_name` will be pulled out of the manifest. Default is None.

get_deformation_field (*self*, *section_data_set_id*, *header_path=None*, *voxel_path=None*)

Extract the local alignment parameters for this dataset. This a 3D vector image (3 components) describing a deformable local mapping from CCF voxels to this section data set's affine-aligned image stack.

Parameters

section_data_set_id [int]

Download the deformation field for this data set

header_path [str, optional] If supplied, the deformation field header will be downloaded to this path.

voxel_path [str, optiona] If supplied, the deformation field voxels will be downloaded to this path.

Returns

numpy.ndarray: 3D X 3 component vector array (origin 0, 0, 0; 25-micron isometric resolution) defining a deformable transformation from CCF-space to affine-transformed image space.

get_experiment_structure_unionizes (*self*, *experiment_id*, *file_name=None*, *is_injection=None*, *structure_ids=None*, *include_descendants=False*, *hemisphere_ids=None*)

Retrieve the structure unionize data for a specific experiment. Filter by structure, injection status, and hemisphere.

Parameters

experiment_id: int ID of the experiment of interest. Corresponds to `section_data_set_id` in the API.

file_name: string File name to save/read the experiments list. If `file_name` is `None`, the `file_name` will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is `None`.

is_injection: boolean If `True`, only return unionize records that disregard non-injection pixels. If `False`, only return unionize records that disregard injection pixels. If `None`, return all records. Default `None`.

structure_ids: list Only return unionize records for a specific set of structures. If `None`, return all records. Default `None`.

include_descendants: boolean Include all descendant records for specified structures. Default `False`.

hemisphere_ids: list Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If `None`, include all records [1, 2, 3]. Default `None`.

get_experiments (*self*, *dataframe=False*, *file_name=None*, *cre=None*, *injection_structure_ids=None*)

Read a list of experiments that match certain criteria. If caching is enabled, this will save the whole (unfiltered) list of experiments to a file.

Parameters

dataframe: boolean Return the list of experiments as a Pandas DataFrame. If `False`, return a list of dictionaries. Default `False`.

file_name: string File name to save/read the structures table. If `file_name` is `None`, the `file_name` will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is `None`.

cre: boolean or list If `True`, return only cre-positive experiments. If `False`, return only cre-negative experiments. If `None`, return all experiments. If list, return all experiments with cre line names in the supplied list. Default `None`.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If `None`, return all experiments. Default `None`.

get_injection_density (*self*, *experiment_id*, *file_name=None*)

Read an injection density volume for a single experiment. Download it first if it doesn't exist. Injection density is the proportion of projecting pixels in a grid voxel only including pixels that are part of the injection site in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_injection_fraction (*self, experiment_id, file_name=None*)

Read an injection fraction volume for a single experiment. Download it first if it doesn't exist. Injection fraction is the proportion of pixels in the injection site in a grid voxel in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_projection_density (*self, experiment_id, file_name=None*)

Read a projection density volume for a single experiment. Download it first if it doesn't exist. Projection density is the proportion of projecting pixels in a grid voxel in [0,1].

Parameters

experiment_id: int ID of the experiment to download/read. This corresponds to section_data_set_id in the API.

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_projection_matrix (*self, experiment_ids, projection_structure_ids=None, hemisphere_ids=None, parameter='projection_volume', dataframe=False*)

get_structure_unionizes (*self, experiment_ids, is_injection=None, structure_ids=None, include_descendants=False, hemisphere_ids=None*)

Get structure unionizes for a set of experiment IDs. Filter the results by injection status, structure, and hemisphere.

Parameters

experiment_ids: list List of experiment IDs. Corresponds to section_data_set_id in the API.

is_injection: boolean If True, only return unionize records that disregard non-injection pixels. If False, only return unionize records that disregard injection pixels. If None, return all records. Default None.

structure_ids: list Only return unionize records for a specific set of structures. If None, return all records. Default None.

include_descendants: boolean Include all descendant records for specified structures. Default False.

hemisphere_ids: list Only return unionize records that disregard pixels outside of a hemisphere. or set of hemispheres. Left = 1, Right = 2, Both = 3. If None, include all records [1, 2, 3]. Default None.

rank_structures (*self, experiment_ids, is_injection, structure_ids=None, hemisphere_ids=None, rank_on='normalized_projection_volume', n=5, threshold=0.01*)

Produces one or more (per experiment) ranked lists of brain structures, using a specified data field.

Parameters

- experiment_ids** [list of int] Obtain injection_structures for these experiments.
- is_injection** [boolean] Use data from only injection (or non-injection) unionizes.
- structure_ids** [list of int, optional] Consider only these structures. It is a good idea to make sure that these structures are not spatially overlapping; otherwise your results will contain redundant information. Defaults to the summary structures - a brain-wide list of nonoverlapping mid-level structures.
- hemisphere_ids** [list of int, optional] Consider only these hemispheres (1: left, 2: right, 3: both). Like with structures, you might get redundant results if you select overlapping options. Defaults to [1, 2].
- rank_on** [str, optional] Rank unionize data using this field (descending). Defaults to `normalized_projection_volume`.
- n** [int, optional] Return only the top n structures.
- threshold** [float, optional] Consider only records whose data value - specified by the `rank_on` parameter - exceeds this value.

Returns

- list** : Each element (1 for each input experiment) is a list of dictionaries. The dictionaries describe the top injection structures in descending order. They are specified by their `structure` and `hemisphere_id` fields and additionally report the value specified by the `rank_on` parameter.

allensdk.core.nwb_data_set module

```
class allensdk.core.nwb_data_set.NwbDataSet (file_name, spike_time_key=None)
```

Bases: object

A very simple interface for extracting electrophysiology data from an NWB file.

```
DEPRECATED_SPIKE_TIMES = 'aibs_spike_times'
```

```
SPIKE_TIMES = 'spike_times'
```

```
fill_sweep_responses (self, fill_value=0.0, sweep_numbers=None, extend_experiment=False)
```

Fill sweep response arrays with a single value.

Parameters

- fill_value: float** Value used to fill sweep response array
- sweep_numbers: list** List of integer sweep numbers to be filled (default all sweeps)
- extend_experiment: bool** If True, extend experiment epoch length to the end of the sweep (undo any truncation)

```
get_experiment_sweep_numbers (self)
```

Get all of the sweep numbers for experiment epochs in the file, not including test sweeps.

```
get_pipeline_version (self)
```

Returns the AI pipeline version number, stored in the metadata field 'generated_by'. If that field is missing, version 0.0 is returned.

Returns

- int tuple: (major, minor)**

get_spike_times (*self*, *sweep_number*, *key=None*)

Return any spike times stored in the NWB file for a sweep.

Parameters

sweep_number: `int` index to access

key [`string`] label where the spike times are stored (default `NwbDataSet.SPIKE_TIMES`)

Returns

list list of spike times in seconds relative to the start of the sweep

get_sweep (*self*, *sweep_number*)

Retrieve the stimulus, response, index_range, and sampling rate for a particular sweep. This method hides the NWB file's distinction between a "Sweep" and an "Experiment". An experiment is a subset of a sweep that excludes the initial test pulse. It also excludes any erroneous response data at the end of the sweep (usually for ramp sweeps, where recording was terminated mid-stimulus).

Some sweeps do not have an experiment, so full data arrays are returned. Sweeps that have an experiment return full data arrays (include the test pulse) with any erroneous data trimmed from the back of the sweep.

Parameters

sweep_number: `int`

Returns

dict A dictionary with 'stimulus', 'response', 'index_range', and 'sampling_rate' elements. The index range is a 2-tuple where the first element indicates the end of the test pulse and the second index is the end of valid response data.

get_sweep_metadata (*self*, *sweep_number*)

Retrieve the sweep level metadata associated with each sweep. Includes information on stimulus parameters like its name and amplitude as well as recording quality metadata, like access resistance and seal quality.

Parameters

sweep_number: `int`

Returns

dict A dictionary with 'aibs_stimulus_amplitude_pa', 'aibs_stimulus_name', 'gain', 'initial_access_resistance', 'seal' elements. These specific fields are ones encoded in the original AIBS in vitro .nwb files.

get_sweep_numbers (*self*)

Get all of the sweep numbers in the file, including test sweeps.

set_spike_times (*self*, *sweep_number*, *spike_times*, *key=None*)

Set or overwrite the spikes times for a sweep.

Parameters

sweep_number [`int`] index to access

key [`string`] where the times are stored (default `NwbDataSet.SPIKE_TIME`)

spike_times: `np.array` array of spike times in seconds

set_sweep (*self*, *sweep_number*, *stimulus*, *response*)

Overwrite the stimulus or response of an NWB file. If the supplied arrays are shorter than stored arrays, they are padded with zeros to match the original data size.

Parameters

sweep_number: int

stimulus: np.array Overwrite the stimulus with this array. If None, stimulus is unchanged.

response: np.array Overwrite the response with this array. If None, response is unchanged.

allensdk.core.obj_utilities module

`allensdk.core.obj_utilities.parse_obj(lines)`

Parse a wavefront obj file into a triplet of vertices, normals, and faces. This parser is specific to obj files generated from our annotation volumes

Parameters

lines [list of str] Lines of input obj file

Returns

vertices [np.ndarray] Dimensions are (nSamples, nCoordinates=3). Locations in the reference space of vertices

vertex_normals [np.ndarray] Dimensions are (nSample, nElements=3). Vectors normal to vertices.

face_vertices [np.ndarray] Dimensions are (sample, nVertices=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertices that make up each face.

face_normals [np.ndarray] Dimensions are (sample, nNormals=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertex normals that make up each face.

Notes

This parser is specialized to the obj files that the Allen Institute for Brain Science generates from our own structure annotations.

`allensdk.core.obj_utilities.read_obj(path)`

allensdk.core.ontology module

class `allensdk.core.ontology.Ontology(df)`

Bases: object

Note: Deprecated from 0.12.5 *Ontology* has been replaced by *StructureTree*.

get_child_ids (*self*, *structure_ids*)

Find the set of ids that are immediate children of one or more structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

set Set of child structure ids

get_children (*self*, *structure_ids*)

Find the set of structures that are immediate children of one or more structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

pandas.DataFrame Set of child structures

get_descendant_ids (*self*, *structure_ids*)

Find the set of the ids of structures that are descendants of one or more structures. The returned set will include the input structure ids.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

set Set of descendant structure ids.

get_descendants (*self*, *structure_ids*)

Find the set of structures that are descendants of one or more structures. The returned set will include the input structures.

Parameters

structure_ids: iterable Any iterable type that contains structure ids that can be cast to integers.

Returns

pandas.DataFrame Set of descendant structures.

structure_descends_from (*self*, *child_id*, *parent_id*)

Return whether one structure id is a descendant of another structure id.

allensdk.core.ophys_experiment_session_id_mapping module

allensdk.core.reference_space module

class allensdk.core.reference_space.**ReferenceSpace** (*structure_tree*, *annotation*, *resolution*)

Bases: object

static check_and_write (*base_dir*, *structure_id*, *fn*)

A many_structure_masks callback that writes the mask to a nrrd file if the file does not already exist.

check_coverage (*self*, *structure_ids*, *domain_mask*)

Determines whether a spatial domain is completely covered by structures in a set.

Parameters

structure_ids [list of int] Specifies the set of structures to check.

domain_mask [numpy ndarray] Same shape as annotation. 1 inside the mask, 0 out. Specifies spatial domain.

Returns

numpy ndarray : 1 where voxels are missing from the candidate, 0 where the candidate exceeds the domain

direct_voxel_counts (*self*)

Determines the number of voxels directly assigned to one or more structures.

Returns

dict : Keys are structure ids, values are the number of voxels directly assigned to those structures.

direct_voxel_map

downsample (*self*, *target_resolution*)

Obtain a smaller reference space by downsampling

Parameters

target_resolution [tuple of numeric] Resolution in microns of the output space.

interpolator [string] Method used to interpolate the volume. Currently only 'nearest' is supported

Returns

ReferenceSpace : A new ReferenceSpace with the same structure tree and a downsampled annotation.

export_itksnap_labels (*self*, *id_type*=<class 'numpy.uint16'>, *label_description_kwargs*=None)

Produces itksnap labels, remapping large ids if needed.

Parameters

id_type [np.integer, optional] Used to determine the type of the output annotation and whether ids need to be remapped to smaller values.

label_description_kwargs [dict, optional] Keyword arguments passed to StructureTree.export_label_description

Returns

np.ndarray : Annotation volume, remapped if needed

pd.DataFrame label_description dataframe

get_slice_image (*self*, *axis*, *position*, *cmap*=None)

Produce a AxBx3 RGB image from a slice in the annotation

Parameters

axis [int] Along which to slice the annotation volume. 0 is coronal, 1 is horizontal, and 2 is sagittal.

position [int] In microns. Take the slice from this far along the specified axis.

cmap [dict, optional] Keys are structure ids, values are rgb triplets. Defaults to structure rgb_triplets.

Returns

np.ndarray : RGB image array.

Notes

If you assign a custom colormap, make sure that you take care of the background in addition to the structures.

make_structure_mask (*self*, *structure_ids*, *direct_only=False*)

Return an indicator array for one or more structures

Parameters

structure_ids [list of int] Make a mask that indicates the union of these structures' voxels

direct_only [bool, optional] If True, only include voxels directly assigned to a structure in the mask. Otherwise include voxels assigned to descendants.

Returns

numpy ndarray : Same shape as annotation. 1 inside mask, 0 outside.

many_structure_masks (*self*, *structure_ids*, *output_cb=None*, *direct_only=False*)

Build one or more structure masks and do something with them

Parameters

structure_ids [list of int] Specify structures to be masked

output_cb [function, optional] Must have the following signature: `output_cb(structure_id, fn)`. On each requested id, `fn` will be curried to make a mask for that id. Defaults to returning the structure id and mask.

direct_only [bool, optional] If True, only include voxels directly assigned to a structure in the mask. Otherwise include voxels assigned to descendants.

Yields

Return values of **output_cb** called on each **structure_id**, **structure_mask** pair.

Notes

`output_cb` is called on every yield, so any side-effects (such as writing to a file) will be carried out regardless of what you do with the return values. You do actually have to iterate through the output, though.

remove_unassigned (*self*, *update_self=True*)

Obtains a structure tree consisting only of structures that have at least one voxel in the annotation.

Parameters

update_self [bool, optional] If True, the contained structure tree will be replaced,

Returns

list of dict : elements are filtered structures

static return_mask_cb (*structure_id*, *fn*)

A basic callback for `many_structure_masks`

total_voxel_counts (*self*)

Determines the number of voxels assigned to a structure or its descendants

Returns

dict : Keys are structure ids, values are the number of voxels assigned to structures' descendants.

total_voxel_map

validate_structures (*self*, *structure_ids*, *domain_mask*)

Determines whether a set of structures produces an exact and nonoverlapping tiling of a spatial domain

Parameters

structure_ids [list of int] Specifies the set of structures to check.

domain_mask [numpy ndarray] Same shape as annotation. 1 inside the mask, 0 out.
Specifies spatial domain.

Returns

set : Ids of structures that are the ancestors of other structures in the supplied set.

numpy ndarray : Indicator for missing voxels.

write_itksnap_labels (*self*, *annotation_path*, *label_path*, ***kwargs*)

Generate a label file (nrrd) and a label_description file (csv) for use with ITKSnap

Parameters

annotation_path [str] write generated label file here

label_path [str] write generated label_description file here

****kwargs** : will be passed to self.export_itksnap_labels

allensdk.core.reference_space_cache module

class allensdk.core.reference_space_cache.**ReferenceSpaceCache** (*resolution*, *reference_space_key*, ***kwargs*)

Bases: *allensdk.api.cache.Cache*

ANNOTATION_KEY = 'ANNOTATION'

MANIFEST_VERSION = 1.2

REFERENCE_SPACE_VERSION_KEY = 'REFERENCE_SPACE_VERSION'

STRUCTURES_KEY = 'STRUCTURES'

STRUCTURE_MASK_KEY = 'STRUCTURE_MASK'

STRUCTURE_MESH_KEY = 'STRUCTURE_MESH'

STRUCTURE_TREE_KEY = 'STRUCTURE_TREE'

TEMPLATE_KEY = 'TEMPLATE'

add_manifest_paths (*self*, *manifest_builder*)

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

get_annotation_volume (*self*, *file_name=None*)

Read the annotation volume. Download it first if it doesn't exist.

Parameters

file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_reference_space (*self, structure_file_name=None, annotation_file_name=None*)

Build a ReferenceSpace from this cache's annotation volume and structure tree. The ReferenceSpace does operations that relate brain structures to spatial domains.

Parameters

structure_file_name: string File name to save/read the structures table. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

annotation_file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

get_structure_mask (*self, structure_id, file_name=None, annotation_file_name=None*)

Read a 3D numpy array shaped like the annotation volume that has non-zero values where voxels belong to a particular structure. This will take care of identifying substructures.

Parameters

structure_id: int ID of a structure.

file_name: string File name to store the structure mask. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

annotation_file_name: string File name to store the annotation volume. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

Notes

This method downloads structure masks from the Allen Institute. To make your own locally, see `ReferenceSpace.many_structure_masks`.

get_structure_mesh (*self, structure_id, file_name=None*)

Obtain a 3D mesh specifying the surface of an annotated structure.

Parameters

structure_id: int ID of a structure.

file_name: string File name to store the structure mesh. If it already exists, it will be read from this file. If file_name is None, the file_name will be pulled out of the manifest. Default is None.

Returns

vertices [np.ndarray] Dimensions are (nSamples, nCoordinates=3). Locations in the reference space of vertices

vertex_normals [np.ndarray] Dimensions are (nSample, nElements=3). Vectors normal to vertices.

face_vertices [np.ndarray] Dimensions are (sample, nVertices=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertices that make up each face.

face_normals [np.ndarray] Dimensions are (sample, nNormals=3). References are given in indices (0-indexed here, but 1-indexed in the file) of vertex normals that make up each face.

Notes

These meshes are meant for 3D visualization and as such have been smoothed. If you are interested in performing quantitative analyses, we recommend that you use the structure masks instead.

get_structure_tree (*self*, *file_name=None*, *structure_graph_id=1*)

Read the list of adult mouse structures and return an StructureTree instance.

Parameters

file_name: string File name to save/read the structures table. If *file_name* is None, the *file_name* will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

structure_graph_id: int Build a tree using structure only from the identified structure graph.

get_template_volume (*self*, *file_name=None*)

Read the template volume. Download it first if it doesn't exist.

Parameters

file_name: string File name to store the template volume. If it already exists, it will be read from this file. If *file_name* is None, the *file_name* will be pulled out of the manifest. Default is None.

classmethod validate_structure_id (*structure_id*)

classmethod validate_structure_ids (*structure_ids*)

allensdk.core.simple_tree module

class allensdk.core.simple_tree.SimpleTree (*nodes*, *node_id_cb*, *parent_id_cb*)

Bases: object

ancestor_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' ancestors

Parameters

node_ids [list of hashable] Items are ids of nodes whose ancestors you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' ancestors' ids.

Notes

Given the tree: A -> B -> C

 '-> D

The ancestors of C are [C, B, A]. The ancestors of A are [A]. The ancestors of D are [D, A]

ancestors (*self*, *node_ids*)

Get one or more nodes' ancestor nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose ancestors will be found.

Returns

list of list of dict : Items are lists of ancestor nodes corresponding to argued ids.

child_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' children

Parameters

node_ids [list of hashable] Items are ids of nodes whose children you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' children's ids.

children (*self*, *node_ids*)

Get one or more nodes' child nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose children will be found.

Returns

list of list of dict : Items are lists of child nodes corresponding to argued ids.

descendant_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' descendants

Parameters

node_ids [list of hashable] Items are ids of nodes whose descendants you wish to find.

Returns

list of list of hashable : Items are lists of input nodes' descendants' ids.

Notes

Given the tree: A -> B -> C

‘-> D

The descendants of A are [B, C, D]. The descendants of C are [].

descendants (*self*, *node_ids*)

Get one or more nodes' descendant nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose descendants will be found.

Returns

list of list of dict : Items are lists of descendant nodes corresponding to argued ids.

filter_nodes (*self*, *criterion*)

Obtain a list of nodes filtered by some criterion

Parameters

criterion [function | node dict => bool] Only nodes for which criterion returns true will be returned.

Returns

list of dict : Items are node dictionaries that passed the filter.

node (*self*, *node_ids=None*)

node_ids (*self*)

Obtain the node ids of each node in the tree

Returns

list : elements are node ids

nodes (*self*, *node_ids=None*)

Get one or more nodes' full dictionaries from their ids.

Parameters

node_ids [list of hashable] Items are ids of nodes to be returned. Default is all.

Returns

list of dict : Items are nodes corresponding to argued ids.

nodes_by_property (*self*, *key*, *values*, *to_fn=None*)

Get nodes by a specified property

Parameters

key [hashable or function] The property used for lookup. Should be unique. If a function, will be invoked on each node.

values [list] Select matching elements from the lookup.

to_fn [function, optional] Defines the outputs, on a per-node basis. Defaults to returning the whole node.

Returns

list : outputs, 1 for each input value.

parent (*self*, *node_ids*)

parent_id (*self*, *node_ids*)

parent_ids (*self*, *node_ids*)

Obtain the ids of one or more nodes' parents

Parameters

node_ids [list of hashable] Items are ids of nodes whose parents you wish to find.

Returns

list of hashable : Items are ids of input nodes' parents in order.

parents (*self*, *node_ids*)

Get one or mode nodes' parent nodes

Parameters

node_ids [list of hashable] Items are ids of nodes whose parents will be found.

Returns

list of dict : Items are parents of nodes corresponding to argued ids.

value_map (*self*, *from_fn*, *to_fn*)

Obtain a look-up table relating a pair of node properties across nodes

Parameters

from_fn [function | node dict => hashable value] The keys of the output dictionary will be obtained by calling from_fn on each node. Should be unique.

to_fn [function | node_dict => value] The values of the output function will be obtained by calling to_fn on each node.

Returns

dict : Maps the node property defined by from_fn to the node property defined by to_fn across nodes.

allensdk.core.sitk_utilities module

`allensdk.core.sitk_utilities.fix_array_dimensions(array, ncomponents=1)`

Convenience function that reorders ndarray dimensions for io with SimpleITK

Parameters

array [np.ndarray] The array to be reordered

ncomponents [int, optional] Number of components per pixel, default 1.

Returns

np.ndarray : Reordered array

`allensdk.core.sitk_utilities.get_sitk_image_information(image)`

Extract information about a SimpleITK image

Parameters

image [sitk.Image] Extract information about this image.

Returns

dict : Extracted information. Includes spacing, origin, size, direction, and number of components per pixel

`allensdk.core.sitk_utilities.read_ndarray_with_sitk(path)`

Read a numpy array from a file using SimpleITK

Parameters

path [str] Read from this path

Returns

image [np.ndarray] Obtained array

information [dict] Additional information about the array

`allensdk.core.sitk_utilities.set_sitk_image_information(image, information)`

Set information on a SimpleITK image

Parameters

image [sitk.Image] Set information on this image.

information [dict] Stores information to be set. Supports spacing, origin, direction. Also checks (but cannot set) size and number of components per pixel

`allensdk.core.sitk_utilities.write_ndarray_with_sitk(array, path, **information)`

Write a numpy array to a file using SimpleITK

Parameters

array [np.ndarray] Array to be written.

path [str] Write to here

****information** [dict] Contains additional information to be stored in the image file. See `set_sitk_image_information` for more information.

allensdk.core.structure_tree module

class `allensdk.core.structure_tree.StructureTree(nodes)`

Bases: `allensdk.core.simple_tree.SimpleTree`

static `clean_structures(structures, whitelist=None, data_transforms=None, re-names=None)`

Convert structures_with_sets query results into a form that can be used to construct a StructureTree

Parameters

structures [list of dict] Each element describes a structure. Should have a structure id path field (str values) and a structure_sets field (list of dict).

whitelist [list of str, optional] Only these fields will be included in the final structure record. Default is the output of `StructureTree.whitelist`.

data_transforms [dict, optional] Keys are str field names. Values are functions which will be applied to the data associated with those fields. Default is to map colors from hex to rgb and convert the structure id path to a list of int.

renames [dict, optional] Controls the field names that appear in the output structure records. Default is to map 'color_hex_triplet' to 'rgb_triplet'.

Returns

list of dict : structures, after conversion of structure_id_path and structure_sets

static `collect_sets(structure)`

Structure sets may be specified by full records or id. This method collects all of the structure set records/ids in a structure record and replaces them with a single list of id records.

static `data_transforms()`

export_label_description (*self*, *alphas=None*, *exclude_label_vis=None*, *exclude_mesh_vis=None*, *label_key='acronym'*)

Produces an itksnap label_description table from this structure tree

Parameters

alphas [dict, optional] Maps structure ids to alpha levels. Optional - will only use provided ids.

exclude_label_vis [list, optional] The structures denoted by these ids will not be visible in ITKSnap.

exclude_mesh_vis [list, optional] The structures denoted by these ids will not have visible meshes in ITKSnap.

label_key: str, optional Use this column for display labels.

Returns

pd.DataFrame : Contains data needed for loading as an ITKSnap label description file.

get_ancestor_id_map (*self*)

Get a dictionary mapping structure ids to ancestor ids across all nodes.

Returns

dict : Keys are structure ids. Values are lists of ancestor ids.

get_colormap (*self*)

Get a dictionary mapping structure ids to colors across all nodes.

Returns

dict : Keys are structure ids. Values are RGB lists of integers.

get_id_acronym_map (*self*)

Get a dictionary mapping structure acronyms to ids across all nodes.

Returns

dict : Keys are structure acronyms. Values are structure ids.

get_name_map (*self*)

Get a dictionary mapping structure ids to names across all nodes.

Returns

dict : Keys are structure ids. Values are structure name strings.

get_structure_sets (*self*)

Lists all unique structure sets that are assigned to at least one structure in the tree.

Returns

list of int : Elements are ids of structure sets.

get_structures_by_acronym (*self*, *acronyms*)

Obtain a list of brain structures from their acronyms

Parameters

names [list of str] Get structures corresponding to these acronyms.

Returns

list of dict : Each item describes a structure.

get_structures_by_id (*self*, *structure_ids*)

Obtain a list of brain structures from their structure ids

Parameters

structure_ids [list of int] Get structures corresponding to these ids.

Returns

list of dict : Each item describes a structure.

get_structures_by_name (*self*, *names*)

Obtain a list of brain structures from their names,

Parameters

names [list of str] Get structures corresponding to these names.

Returns

list of dict : Each item describes a structure.

get_structures_by_set_id (*self*, *structure_set_ids*)

Obtain a list of brain structures from by the sets that contain them.

Parameters

structure_set_ids [list of int] Get structures belonging to these structure sets.

Returns

list of dict : Each item describes a structure.

has_overlaps (*self*, *structure_ids*)

Determine if a list of structures contains structures along with their ancestors

Parameters

structure_ids [list of int] Check this set of structures for overlaps

Returns

set : Ids of structures that are the ancestors of other structures in the supplied set.

static hex_to_rgb (*hex_color*)

Convert a hexadecimal color string to a uint8 triplet

Parameters

hex_color [string] Must be 6 characters long, unless it is 7 long and the first character is #. If hex_color is a triplet of int, it will be returned unchanged.

Returns

list of int : 3 characters long - 1 per two characters in the input string.

static path_to_list (*path*)

Structure id paths are sometimes formatted as “/”-seperated strings. This method converts them to a list of integers, if needed.

static renames ()

structure_descends_from (*self*, *child_id*, *parent_id*)

Tests whether one structure descends from another.

Parameters

child_id [int] Id of the putative child structure.

parent_id [int] Id of the putative parent structure.

Returns

bool : True if the structure specified by child_id is a descendant of the one specified by parent_id. Otherwise False.

static whitelist ()

allensdk.core.swc module

class allensdk.core.swc.**Compartment** (*args, **kwargs)

Bases: dict

A dictionary class storing information about a single morphology node

print_node (*self*)

print out compartment information with field names

```
class allensdk.core.swc.Marker (*args, **kwargs)
```

Bases: dict

Simple dictionary class for handling reconstruction marker objects.

CUT_DENDRITE = 10

NO_RECONSTRUCTION = 20

SPACING = [0.1144, 0.1144, 0.28]

```
class allensdk.core.swc.Morphology (compartment_list=None, compartment_index=None)
```

Bases: object

Keep track of the list of compartments in a morphology and provide a few helper methods (soma, tree information, pruning, etc).

APICAL_DENDRITE = 4

AXON = 2

BASAL_DENDRITE = 3

DENDRITE = 3

NODE_TYPES = [1, 2, 3, 3, 4]

SOMA = 1

```
append (self, node_list)
```

Add additional nodes to this Morphology. Those nodes must originate from another morphology object.

Parameters

node_list: list of Morphology nodes

```
apply_affine (self, aff, scale=None)
```

Apply an affine transform to all compartments in this morphology. Node radius is adjusted as well.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

```
change_parent (self, child, parent)
```

Change the parent of a node. The child node is adjusted to point to the new parent, the child is taken off of the previous parent's child list, and it is added to the new parent's child list.

Parameters

child: integer or Morphology Object The ID of the child node, or the child node itself

parent: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

Nothing**children_of** (*self*, *seg*)

Returns a list of the children of the specified node

Parameters**seg: integer or Morphology Object** The ID of the parent node, or the parent node itself**Returns****A list of the child morphology objects. If the ID of the parent node is invalid, None is returned.****compartment_index**

Return the compartment index. This is a property to ensure that the compartment list and compartment index are in sync.

compartment_index_by_type (*self*, *compartment_type*)

Return an dictionary of compartments indexed by id that all have a particular compartment type.

Parameters**compartment_type: int** Desired compartment type**Returns****A dictionary of Morphology Objects, indexed by ID****compartment_list**

Return the compartment list. This is a property to ensure that the compartment list and compartment index are in sync.

compartment_list_by_type (*self*, *compartment_type*)

Return an list of all compartments having the specified compartment type.

Parameters**compartment_type: int** Desired compartment type**Returns****A list of of Morphology Objects****convert_type** (*self*, *old_type*, *new_type*)

Converts all compartments from one type to another. Nodes of the original type are not affected so this procedure can also be used as a merge procedure.

Parameters**old_type: enum** The compartment type to be changed. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE**new_type: enum** The target compartment type. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE**delete_tree** (*self*, *n*)

Delete tree, and all of its compartments, from the morphology.

Parameters**n: Integer** The tree number to delete

find (*self*, *x*, *y*, *z*, *dist*, *node_type=None*)

Returns a list of Morphology Objects located within ‘dist’ of coordinate (x,y,z). If node_type is specified, the search will be constrained to return only nodes of that type.

Parameters

x, y, z: float The x,y,z coordinates from which to search around

dist: float The search radius

node_type: enum (optional) One of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE or APICAL_DENDRITE

Returns

A list of all Morphology Objects matching the search criteria

node (*self*, *n*)

Returns the morphology node having the specified ID.

Parameters

n: integer ID of desired node

Returns

A morphology object having the specified ID, or None if such a node doesn’t exist

num_nodes

Return the number of compartments in the morphology.

num_trees

Return the number of trees in the morphology. A tree is defined as everything following from a single root compartment.

parent_of (*self*, *seg*)

Returns parent of the specified node.

Parameters

seg: integer or Morphology Object The ID of the child node, or the child node itself

Returns

A morphology object, or None if no parent exists or if the specified node ID doesn’t exist

root

[deprecated] Returns root node of soma, if present. Use ‘soma’ instead of ‘root’

save (*self*, *file_name*)

Write this morphology out to an SWC file

Parameters

file_name: string desired name of your SWC file

soma

Returns root node of soma, if present

sparsify (*self*, *modulo*, *compress_ids=False*)

Return a new Morphology object that has a given number of non-leaf, non-root nodes removed. IDs can be reassigned so as to be continuous.

Parameters

modulo: int keep 1 out of every modulo nodes.

compress_ids: boolean Reassign ids so that ids are continuous (no missing id numbers).

Returns

Morphology A new morphology instance

strip_all_other_types (*self*, *node_type*, *keep_soma=True*)

Strips everything from the morphology except for the specified type. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The compartment type to keep in the morphology. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

keep_soma: Boolean (optional) True (default) if soma nodes should remain in the morphology, and False if the soma should also be stripped

strip_type (*self*, *node_type*)

Strips all compartments of the specified type from the morphology. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The compartment type to strip from the morphology. Use one of the following constants: SOMA, AXON, DENDRITE, BASAL_DENDRITE, or APICAL_DENDRITE

stumpify_axon (*self*, *count=10*)

Remove all axon compartments except the first ‘count’ nodes, as counted from the connected axon root.

Parameters

count: Integer The length of the axon ‘stump’, in number of compartments

tree (*self*, *n*)

Returns a list of all Morphology Nodes within the specified tree. A tree is defined as a fully connected graph of nodes. Each tree has exactly one root.

Parameters

n: integer ID of desired tree

Returns

A list of all morphology objects in the specified tree, or None

if the tree doesn’t exist

write (*self*, *file_name*)

`allensdk.core.swc.read_marker_file` (*file_name*)

read in a marker file and return a list of dictionaries

`allensdk.core.swc.read_swc` (*file_name*, *columns='NOT_USED'*, *numeric_columns='NOT_USED'*) *nu-*

Read in an SWC file and return a Morphology object.

Parameters

file_name: string SWC file name.

Returns

Morphology A Morphology instance.

allensdk.core.typing module

class allensdk.core.typing.SupportsStr
Bases: typing._Protocol
Classes that support the __str__ method

Module contents**6.1.5 allensdk.ephys package****Submodules****allensdk.ephys.ephys_extractor module**

class allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor(*ramps_ext*,
short_squares_ext,
long_squares_ext,
subthresh_min_amp=-100)

Bases: object

SAG_TARGET = -100.0

SUBTHRESH_MAX_AMP = 0

as_dict (*self*)
Create dict of cell features.

cell_features (*self*)

long_squares_features (*self*, *option=None*)

long_squares_stim_amps (*self*, *option=None*)

process (*self*, *keys=None*)
Processes features. Can take a specific key (or set of keys) to do a subset of processing.

ramps_features (*self*, *all=False*)

short_squares_features (*self*)

```
class allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor (t=None,  
                                                             v=None,  
                                                             i=None,  
                                                             start=None,  
                                                             end=None,  
                                                             filter=10.0,  
                                                             dv_cutoff=20.0,  
                                                             max_interval=0.005,  
                                                             min_height=2.0,  
                                                             min_peak=-  
                                                             30.0,  
                                                             thresh_frac=0.05,  
                                                             base-  
                                                             line_interval=0.1,  
                                                             base-  
                                                             line_detect_thresh=0.3,  
                                                             id=None)
```

Bases: object

Feature calculation for a sweep (voltage and/or current time series).

as_dict (*self*)

Create dict of features and spikes.

burst_metrics (*self*)

Find bursts and return max “burstiness” index (normalized max rate in burst vs out).

Returns

max_burstiness_index [max “burstiness” index across detected bursts]

num_bursts [number of bursts detected]

delay_metrics (*self*)

Calculates ratio of latency to dominant time constant of rise before spike

Returns

delay_ratio [ratio of latency to tau (higher means more delay)]

tau [dominant time constant of rise before spike]

estimate_sag (*self*, *peak_width=0.005*)

Calculate the sag in a hyperpolarizing voltage response.

Parameters

peak_width [window width to get more robust peak estimate in sec (default 0.005)]

Returns

sag [fraction that membrane potential relaxes back to baseline]

estimate_time_constant (*self*)

Calculate the membrane time constant by fitting the voltage response with a single exponential.

Returns

tau [membrane time constant in seconds]

is_spike_feature_affected_by_clipping (*self*, *key*)

pause_metrics (*self*)

Estimate average number of pauses and average fraction of time spent in a pause

Attempts to detect pauses with a variety of conditions and averages results together.

Pauses that are consistently detected contribute more to estimates.

Returns

avg_n_pauses [average number of pauses detected across conditions]

avg_pause_frac [average fraction of interval (between start and end) spent in a pause]

max_reliability [max fraction of times most reliable pause was detected given weights tested]

n_max_rel_pauses [number of pauses detected with *max_reliability*]

process_new_spike_feature (*self*, *feature_name*, *feature_func*, *affected_by_clipping=False*)

Add new spike-level feature calculation function

The function should take this sweep extractor as its argument. Its results can be accessed by calling the method `spike_feature(<feature_name>)`.

process_new_sweep_feature (*self*, *feature_name*, *feature_func*)

Add new sweep-level feature calculation function

The function should take this sweep extractor as its argument. Its results can be accessed by calling the method `sweep_feature(<feature_name>)`.

process_spikes (*self*)

Perform spike-related feature analysis

set_stimulus_amplitude_calculator (*self*, *function*)

spike_feature (*self*, *key*, *include_clipped=False*, *force_exclude_clipped=False*)

Get specified feature for every spike.

Parameters

key [feature name]

include_clipped: return values for every identified spike, even when clipping means they will be incorrect/

Returns

spike_feature_values [ndarray of features for each spike]

spike_feature_keys (*self*)

Get list of every available spike feature.

spikes (*self*)

Get all features for each spike as a list of records.

stimulus_amplitude (*self*)

sweep_feature (*self*, *key*, *allow_missing=False*)

Get sweep-level feature (*key*).

Parameters

key [name of sweep-level feature]

allow_missing [return np.nan if key is missing for sweep (default False)]

Returns

sweep_feature [sweep-level feature value]

sweep_feature_keys (*self*)

Get list of every available sweep-level feature.

voltage_deflection (*self*, *deflect_type=None*)

Measure deflection (min or max, between start and end if specified).

Parameters

deflect_type [measure minimal ('min') or maximal ('max') voltage deflection] If not specified, it will check to see if the current (i) is positive or negative between start and end, then choose 'max' or 'min', respectively If the current is not defined, it will default to 'min'.

Returns

deflect_v [peak]

deflect_index [index of peak deflection]

```
class allensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor (t_set=None,  
                                                                v_set=None,  
                                                                i_set=None,  
                                                                start=None,  
                                                                end=None,  
                                                                filter=10.0,  
                                                                dv_cutoff=20.0,  
                                                                max_interval=0.005,  
                                                                min_height=2.0,  
                                                                min_peak=-30.0,  
                                                                thresh_frac=0.05,  
                                                                base_line_interval=0.1,  
                                                                base_line_detect_thresh=0.3,  
                                                                id_set=None)
```

Bases: object

classmethod from_sweeps (*sweep_list*)

Initialize EphysSweepSetFeatureExtractor object with a list of pre-existing sweep feature extractor objects.

process_spikes (*self*)

Analyze spike features for all sweeps.

spike_feature_averages (*self*, *key*)

Get nparray of average spike-level feature (*key*) for all sweeps

sweep_features (*self*, *key*, *allow_missing=False*)

Get nparray of sweep-level feature (*key*) for all sweeps

Parameters

key [name of sweep-level feature]

allow_missing [return np.nan if key is missing for sweep (default False)]

Returns

sweep_feature [nparray of sweep-level feature values]

sweeps (*self*)

Get list of EphysSweepFeatureExtractor objects.

```
allensdk.ephys.ephys_extractor.cell_extractor_for_nwb(dataset, ramps,
                                                    short_squares, long_squares,
                                                    subthresh_min_amp=-100)
```

Initialize EphysCellFeatureExtractor object from NWB data set

Parameters

dataset [NwbDataSet]

ramps [list of sweep numbers of ramp sweeps]

short_squares [list of sweep numbers of short square sweeps]

long_squares [list of sweep numbers of long square sweeps]

```
allensdk.ephys.ephys_extractor.extractor_for_nwb_sweeps(dataset, sweep_numbers,
                                                         fixed_start=None,
                                                         fixed_end=None,
                                                         dv_cutoff=20.0,
                                                         thresh_frac=0.05)
```

```
allensdk.ephys.ephys_extractor.fit_fit_slope(ext)
```

Fit the rate and stimulus amplitude to a line and return the slope of the fit.

```
allensdk.ephys.ephys_extractor.input_resistance(ext)
```

Estimate input resistance in MOhms, assuming all sweeps in passed extractor are hyperpolarizing responses.

```
allensdk.ephys.ephys_extractor.membrane_time_constant(ext)
```

Average the membrane time constant values estimated from each sweep in passed extractor.

```
allensdk.ephys.ephys_extractor.reset_long_squares_start(when)
```

allensdk.ephys.ephys_features module

exception allensdk.ephys.ephys_features.FeatureError

Bases: Exception

Generic Python-exception-derived object raised by feature detection functions.

```
allensdk.ephys.ephys_features.adaptation_index(isis)
```

Calculate adaptation index of *isis*.

```
allensdk.ephys.ephys_features.analyze_trough_details(v, t, spike_indexes,
                                                       peak_indexes, clipped=None,
                                                       end=None, filter=10.0,
                                                       heavy_filter=1.0,
                                                       term_frac=0.01,
                                                       adp_thresh=0.5, tol=0.5,
                                                       flat_interval=0.002,
                                                       adp_max_delta_t=0.005,
                                                       adp_max_delta_v=10.0,
                                                       dvd_t=None)
```

Analyze trough to determine if an ADP exists and whether the reset is a ‘detour’ or ‘direct’

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]
peak_indexes [numpy array of spike peak indexes]
end [end of time window (optional)]
filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (default 1)]
heavy_filter [lower cutoff frequency for 4-pole low-pass Bessel filter in kHz (default 1)]
thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]
adp_thresh: minimum dV/dt in V/s to exceed to be considered to have an ADP (optional, default 1.5)

tol [tolerance for evaluating whether Vm drops appreciably further after end of spike (default 1.0 mV)]

flat_interval: if the trace is flat for this duration, stop looking for an ADP (default 0.002 s)

adp_max_delta_t: max possible ADP delta t (default 0.005 s)

adp_max_delta_v: max possible ADP delta v (default 10 mV)

dvdvdt [pre-calculated time-derivative of voltage (optional)]

Returns

isi_types [numpy array of isi reset types (direct or detour)]

fast_trough_indexes [numpy array of indexes at the start of the trough (i.e. end of the spike)]

adp_indexes [numpy array of adp indexes (np.nan if there was no ADP in that ISI)]

slow_trough_indexes [numpy array of indexes at the minimum of the slow phase of the trough] (if there wasn't just a fast phase)

`allensdk.ephys.ephys_features.average_rate(t, spikes, start, end)`
Calculate average firing rate during interval between *start* and *end*.

Parameters

t [numpy array of times in seconds]

spikes [numpy array of spike indexes]

start [start of time window for spike detection]

end [end of time window for spike detection]

Returns

avg_rate [average firing rate in spikes/sec]

`allensdk.ephys.ephys_features.average_voltage(v, t, start=None, end=None)`
Calculate average voltage between start and end.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

start [start of time window for spike detection (optional, default None)]

end [end of time window for spike detection (optional, default None)]

Returns

v_avg [average voltage]

`allensdk.ephys.ephys_features.calculate_dvdt` (*v*, *t*, *filter=None*)

Low-pass filters (if requested) and differentiates voltage by time.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default None)]

Returns

dvdt [numpy array of time-derivative of voltage (V/s = mV/ms)]

`allensdk.ephys.ephys_features.check_thresholds_and_peaks` (*v*, *t*, *spike_indexes*,
peak_indexes, *upstroke_indexes*,
end=None,
max_interval=0.005,
thresh_frac=0.05, *filter=10.0*, *dvdt=None*,
tol=1.0)

Validate thresholds and peaks for set of spikes

Check that peaks and thresholds for consecutive spikes do not overlap Spikes with overlapping thresholds and peaks will be merged.

Check that peaks and thresholds for a given spike are not too far apart.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]

peak_indexes [numpy array of indexes of spike peaks]

upstroke_indexes [numpy array of indexes of spike upstrokes]

max_interval [maximum allowed time between start of spike and time of peak in sec (default 0.005)]

thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

tol [tolerance for returning to threshold in mV (optional, default 1)]

Returns

spike_indexes [numpy array of modified spike indexes]

peak_indexes [numpy array of modified spike peak indexes]

upstroke_indexes [numpy array of modified spike upstroke indexes]

clipped [numpy array of clipped status of spikes]

`allensdk.ephys.ephys_features.detect_bursts` (*isis*, *isi_types*, *fast_tr_v*, *fast_tr_t*, *slow_tr_v*,
slow_tr_t, *thr_v*, *tol=0.5*, *pause_cost=1.0*)

Detect bursts in spike train.

Parameters

isis [numpy array of n interspike intervals]
isi_types [numpy array of n interspike interval types]
fast_tr_v [numpy array of fast trough voltages for the n + 1 spikes of the train]
fast_tr_t [numpy array of fast trough times for the n + 1 spikes of the train]
slow_tr_v [numpy array of slow trough voltages for the n + 1 spikes of the train]
slow_tr_t [numpy array of slow trough times for the n + 1 spikes of the train]
thr_v [numpy array of threshold voltages for the n + 1 spikes of the train]
tol [tolerance for the difference in slow trough voltages and thresholds (default 0.5 mV)] Used to identify “delay” interspike intervals that occur within a burst

Returns

bursts [list of bursts] Each item in list is a tuple of the form (burst_index, start, end) where *burst_index* is a comparison index between the highest instantaneous rate within the burst vs the highest instantaneous rate outside the burst. *start* is the index of the first ISI of the burst, and *end* is the ISI index immediately following the burst.

`allensdk.ephys.ephys_features.detect_pauses(isis, isi_types, cost_weight=1.0)`

Determine which ISIs are “pauses” in ongoing firing.

Pauses are unusually long ISIs with a “detour reset” among “direct resets”.

Parameters

isis [numpy array of interspike intervals]
isi_types [numpy array of interspike interval types (‘direct’ or ‘detour’)]
cost_weight [weight for cost function for calling an ISI a pause] Higher cost weights lead to fewer ISIs identified as pauses. The cost function also depends on the difference between the duration of the “pause” ISIs and the average duration and standard deviation of “non-pause” ISIs.

Returns

pauses [numpy array of indices corresponding to pauses in *isis*]

`allensdk.ephys.ephys_features.detect_putative_spikes(v, t, start=None, end=None, filter=10.0, dv_cutoff=20.0)`

Perform initial detection of spikes and return their indexes.

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
start [start of time window for spike detection (optional)]
end [end of time window for spike detection (optional)]
filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]
dv_cutoff [minimum dV/dt to qualify as a spike in V/s (optional, default 20)]
dvdtdt [pre-calculated time-derivative of voltage (optional)]

Returns

putative_spikes [numpy array of preliminary spike indexes]

```
allensdk.ephys.ephys_features.estimate_adjusted_detection_parameters(v_set,
                                                                    t_set,
                                                                    inter-
                                                                    val_start,
                                                                    inter-
                                                                    val_end,
                                                                    fil-
                                                                    ter=10)
```

Estimate adjusted values for spike detection by analyzing a period when the voltage changes quickly but passively (due to strong current stimulation), which can result in spurious spike detection results.

Parameters

v_set [list of numpy arrays of voltage time series in mV]
t_set [list of numpy arrays of times in seconds]
interval_start [start of analysis interval (sec)]
interval_end [end of analysis interval (sec)]

Returns

new_dv_cutoff [adjusted dv/dt cutoff (V/s)]
new_thresh_frac [adjusted fraction of avg upstroke to find threshold]

```
allensdk.ephys.ephys_features.filter_putative_spikes(v, t, spike_indexes,
                                                    peak_indexes, min_height=2.0,
                                                    min_peak=-30.0, filter=10.0,
                                                    dvdt=None)
```

Filter out events that are unlikely to be spikes based on:

- Voltage failing to go down between peak and the next spike's threshold
- Height (threshold to peak)
- Absolute peak level

Parameters

v [numpy array of voltage time series in mV]
t [numpy array of times in seconds]
spike_indexes [numpy array of preliminary spike indexes]
peak_indexes [numpy array of indexes of spike peaks]
min_height [minimum acceptable height from threshold to peak in mV (optional, default 2)]
min_peak [minimum acceptable absolute peak level in mV (optional, default -30)]
filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]
dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

spike_indexes [numpy array of threshold indexes]
peak_indexes [numpy array of peak indexes]

```
allensdk.ephys.ephys_features.find_downstroke_indexes(v, t, peak_indexes,  
trough_indexes,  
clipped=None, filter=10.0,  
dvdt=None)
```

Find indexes of minimum voltage (troughs) between spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

peak_indexes [numpy array of spike peak indexes]

trough_indexes [numpy array of threshold indexes]

clipped: boolean array - False if spike not clipped by edge of window

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

downstroke_indexes [numpy array of downstroke indexes]

```
allensdk.ephys.ephys_features.find_peak_indexes(v, t, spike_indexes, end=None)
```

Find indexes of spike peaks.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of preliminary spike indexes]

end [end of time window for spike detection (optional)]

```
allensdk.ephys.ephys_features.find_time_index(t, t_0)
```

Find the index value of a given time (t_0) in a time series (t).

```
allensdk.ephys.ephys_features.find_trough_indexes(v, t, spike_indexes, peak_indexes,  
clipped=None, end=None)
```

Find indexes of minimum voltage (trough) between spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]

peak_indexes [numpy array of spike peak indexes]

end [end of time window (optional)]

Returns

trough_indexes [numpy array of threshold indexes]

```
allensdk.ephys.ephys_features.find_upstroke_indexes(v, t, spike_indexes, peak_indexes,  
filter=10.0, dvdt=None)
```

Find indexes of maximum upstroke of spike.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of preliminary spike indexes]

peak_indexes [numpy array of indexes of spike peaks]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdv [pre-calculated time-derivative of voltage (optional)]

Returns

upstroke_indexes [numpy array of upstroke indexes]

`allensdk.ephys.ephys_features.find_widths(v, t, spike_indexes, peak_indexes, trough_indexes, clipped=None)`

Find widths at half-height for spikes.

Widths are only returned when heights are defined

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

spike_indexes [numpy array of spike indexes]

peak_indexes [numpy array of spike peak indexes]

trough_indexes [numpy array of trough indexes]

Returns

widths [numpy array of spike widths in sec]

`allensdk.ephys.ephys_features.fit_membrane_time_constant(v, t, start, end, min_rsme=0.0001)`

Fit an exponential to estimate membrane time constant between start and end

Parameters

v [numpy array of voltages in mV]

t [numpy array of times in seconds]

start [start of time window for exponential fit]

end [end of time window for exponential fit]

min_rsme: minimal acceptable root mean square error (default 1e-4)

Returns

a, inv_tau, y0 [Coefficients of equation $y_0 + a * \exp(-\text{inv_tau} * x)$]

returns np.nan for values if fit fails

`allensdk.ephys.ephys_features.fit_prespike_time_constant(v, t, start, spike_time, dv_limit=-0.001, tau_limit=0.3)`

Finds the dominant time constant of the pre-spike rise in voltage

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

start [start of voltage rise (seconds)]

spike_time [time of first spike (seconds)]

dv_limit [dV/dt cutoff (default -0.001)] Shortens fit window if rate of voltage drop exceeds this limit

tau_limit [upper bound for slow time constant (seconds, default 0.3)] If the slower time constant of a double-exponential fit is twice that of the faster and exceeds this limit, the faster one will be considered the dominant one

Returns

tau [dominant time constant (seconds)]

`allensdk.ephys.ephys_features.get_isis(t, spikes)`

Find interspike intervals in sec between spikes (as indexes).

`allensdk.ephys.ephys_features.has_fixed_dt(t)`

Check that all time intervals are identical.

`allensdk.ephys.ephys_features.latency(t, spikes, start)`

Calculate time to the first spike.

`allensdk.ephys.ephys_features.norm_diff(a)`

Calculate average of $(a[i] - a[i+1]) / (a[i] + a[i+1])$.

`allensdk.ephys.ephys_features.norm_sq_diff(a)`

Calculate average of $(a[i] - a[i+1])^2 / (a[i] + a[i+1])^2$.

`allensdk.ephys.ephys_features.refine_threshold_indexes(v, t, upstroke_indexes, thresh_frac=0.05, filter=10.0, dvdt=None)`

Refine threshold detection of previously-found spikes.

Parameters

v [numpy array of voltage time series in mV]

t [numpy array of times in seconds]

upstroke_indexes [numpy array of indexes of spike upstrokes (for threshold target calculation)]

thresh_frac [fraction of average upstroke for threshold calculation (optional, default 0.05)]

filter [cutoff frequency for 4-pole low-pass Bessel filter in kHz (optional, default 10)]

dvdt [pre-calculated time-derivative of voltage (optional)]

Returns

threshold_indexes [numpy array of threshold indexes]

`allensdk.ephys.extract_cell_features` module

`allensdk.ephys.extract_cell_features.extract_cell_features(data_set, ramp_sweep_numbers, short_square_sweep_numbers, long_square_sweep_numbers, sub-thresh_min_amp=None)`

`allensdk.ephys.extract_cell_features.extract_sweep_features(data_set, sweeps_by_type)`

```
allensdk.ephys.extract_cell_features.get_ramp_stim_characteristics(i, t)
    Identify the start time and start index of a ramp sweep.

allensdk.ephys.extract_cell_features.get_square_stim_characteristics(i, t,
                                                                    no_test_pulse=False)
    Identify the start time, duration, amplitude, start index, and end index of a square stimulus. This assumes that
    there is a test pulse followed by the stimulus square.

allensdk.ephys.extract_cell_features.get_stim_characteristics(i, t,
                                                             no_test_pulse=False)
    Identify the start time, duration, amplitude, start index, and end index of a general stimulus. This assumes that
    there is a test pulse followed by the stimulus square.

allensdk.ephys.extract_cell_features.mean_features_spike_zero(sweeps)
    Compute mean feature values for the first spike in list of extractors
```

allensdk.ephys.feature_extractor module

```
class allensdk.ephys.feature_extractor.EphysFeatureExtractor
    Bases: object

    adaptation_index(self, spikes, stim_end)
    calculate_trough(self, spike, v, curr, t, next_idx)
    isicv(self, spikes)
    process_instance(self, name, v, curr, t, onset, dur, stim_name)
    push_summary(self, new_summary)
    score_feature_set(self, set_num)
    summarize(self, summary)

class allensdk.ephys.feature_extractor.EphysFeatures(name)
    Bases: object

    clone(self, param_dict)
    print_out(self)
```

Module contents

6.1.6 allensdk.internal package

Subpackages

allensdk.internal.api package

Subpackages

allensdk.internal.api.queries package

Submodules

allensdk.internal.api.queries.biophysical_module_api module

class allensdk.internal.api.queries.biophysical_module_api.**BiophysicalModuleApi** (*base_uri=None*)
Bases: *allensdk.api.queries.rma_template.RmaTemplate*

get_neuronal_model_runs (*self, neuronal_model_run_ids=None*)
List Neuronal Model Rusn available through LIMS with associated info needed to run in NEURON.

Parameters

neuronal_model_run_ids [integer or list of integers, optional] only select specific neuronal_model_runs.

Returns

dict [neuronal model run metadata]

get_neuronal_models (*self, neuronal_model_ids=None*)
List Neuronal Models available through LIMS with associated info needed to run in NEURON.

Parameters

neuronal_model_ids [integer or list of integers, optional] only select specific neuronal_models.

Returns

dict [neuronal model metadata]

rma_templates = {'biophysical_lims_queries': [{'name': 'neuronal_model_runs_by_ids',

allensdk.internal.api.queries.biophysical_module_reader module

class allensdk.internal.api.queries.biophysical_module_reader.**BiophysicalModuleReader**
Bases: *object*

MOD_FILE_TYPE_ID = 292178729

MORPHOLOGY_TYPE_ID = 303941301

STIMULUS_CONTENT_TYPE = None

fit_parameters_file_entries (*self*)
read the fit_parameter file path from the lims result corresponding to the stimulus file :return: well_known_file entries :rtype: array of dicts

fit_parameters_path (*self*)
Get the path to the fit parameters file from the lims result. :return: path to file :rtype: string

lims_working_directory (*self*)
While this is the same directory as the neuronal_model_run directory, it can be mocked out for testing if the other directory is read only.

mod_file_entries (*self*)
read the NERUON .mod file entries from the lims result corresponding to the NeuronModel :return: well known file entries :rtype: array of dicts

mod_file_paths (*self*)
Get the paths to the mod files from the lims result. :return: paths to mod files :rtype: array of strings

model_type (*self*)
TODO: comment

morphology_file_entries (*self*)

read the well known file paths from the lims result corresponding to the morphology

Returns

array of dicts: well known file entries

morphology_path (*self*)

Get the path to the morphology file from the lims result. :return: path to morphology file :rtype: string

neuronal_model_run_dir (*self*)

read the directory path where output goes from the lims optimization config json

Returns

string: directory path

read_json (*self*, *path*)

read_json_string (*self*, *json_string*)

read_lims_file (*self*, *lims_path*)

read_lims_message (*self*, *message*, *lims_path*)

set_workflow_state (*self*, *state*)

stimulus_file_entries (*self*)

read the well known file path from the lims result corresponding to the stimulus file :return: well_known_file entries :rtype: array of dicts

stimulus_path (*self*)

Get the path to the stimulus file from the lims result. :return: path to stimulus file :rtype: string

sweep_entries (*self*)

read the sweep entries from the lims result corresponding to the stimulus :return: stimulus sweep entries :rtype: array of dicts

sweep_numbers (*self*)

Get the stimulus sweep numbers from the lims result :return: list of sweep numbers :rtype: array of ints

sweep_numbers_by_type (*self*)

to_manifest (*self*, *manifest_path*=None)

update_well_known_file (*self*, *path*, *well_known_file_type_id*=None)

write_file (*self*, *path*)

allensdk.internal.api.queries.grid_data_api_prerelease module

class allensdk.internal.api.queries.grid_data_api_prerelease.**GridDataApiPrerelease** (*storage_dir*, *res-*

o-
lu-
tion=None
base_uri=l

Bases: *allensdk.api.queries.grid_data_api.GridDataApi*

Client for retrieving prereleased mouse connectivity data from lims.

Parameters

base_uri [string, optional] Does not affect pulling from lims.

file_name [string, optional] File name to save/read storage_directories dict. Passed to Grid-DataApiPrerelease constructor.

GRID_DATA_DIRECTORY = 'grid'

download_projection_grid_data (*self, path, experiment_id, file_name*)
Copy data from path to file_name.

Parameters

path [string] path to file in shared directory (copy source)

experiment_id [int] image series id.

file_name [string] path to file destination (copy target)

classmethod from_file_name (*file_name, cache=True, **kwargs*)
Alternative constructor using cache path file_name.

Parameters

file_name [string] Path where storage_directories will be saved.

****kwargs** Keyword arguments to be supplied to `__init__`

Returns

cls [instance of GridDataApiPrerelease]

allensdk.internal.api.queries.mouse_connectivity_api_prerelease module

class allensdk.internal.api.queries.mouse_connectivity_api_prerelease.**MouseConnectivityApi**

Bases: *allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi*

Client for retrieving prereleased mouse connectivity data from lims.

Parameters

base_uri [string, optional] Does not affect pulling from lims.

file_name [string, optional] File name to save/read storage_directories dict. Passed to Grid-DataApiPrerelease constructor.

download_data_mask (*self, path, experiment_id, resolution*)

download_injection_density (*self, path, experiment_id, resolution*)

download_injection_fraction (*self, path, experiment_id, resolution*)

download_projection_density (*self, path, experiment_id, resolution*)

get_experiments (*self*)

Fetch experiment metadata from the Mouse Brain Connectivity Atlas.

Parameters

structure_ids [integer or list, optional] injection structure

Returns

url [string] The constructed URL

get_structure_unionizes (*self*)

allensdk.internal.api.queries.optimize_config_reader module

```
class allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader
    Bases: object
```

```
MOD_FILE_TYPE_ID = 292178729
```

```
MORPHOLOGY_TYPE_ID = 303941301
```

```
NEURONAL_MODEL_PARAMETERS = 329230374
```

```
STIMULUS_CONTENT_TYPE = None
```

```
build_manifest (self, manifest_path=None)
```

```
lims_working_directory (self)
```

While this is the same directory as the optimize directory, it can be mocked out for testing if the optimize directory is write only.

```
mod_file_entries (self)
```

read the NERUON .mod file entries from the lims result corresponding to the NeuronModel :return: well known file entries :rtype: array of dicts

```
mod_file_paths (self)
```

Get the paths to the mod files from the lims result. :return: paths to mod files :rtype: array of strings

```
morphology_file_entries (self)
```

read the well known file paths from the lims result corresponding to the morphology

Returns

array of dicts: well known file entries

```
morphology_path (self)
```

Get the path to the morphology file from the lims result. :return: path to morphology file :rtype: string

```
neuronal_model_optimize_dir (self)
```

read the directory path where output goes from the lims optimization config json

Returns

string: directory path

```
output_directory (self)
```

```
read_json (self, path)
```

```
read_json_string (self, json_string)
```

```
read_lims_file (self, lims_path)
```

```
read_lims_message (self, message, lims_path)
```

```
stimulus_file_entries (self)
```

read the well known file path from the lims result corresponding to the stimulus file :return: well_known_file entries :rtype: array of dicts

```
stimulus_path (self)
```

Get the path to the stimulus file from the lims result. :return: path to stimulus file :rtype: string

```
sweep_entries (self)
```

read the sweep entries from the lims result corresponding to the stimulus :return: stimulus sweep entries :rtype: array of dicts

```
sweep_numbers (self)
```

Get the stimulus sweep numbers from the lims result :return: list of sweep numbers :rtype: array of ints

```
to_manifest (self, manifest_path=None)
update_well_known_file (self, path, well_known_file_type_id=None)
write_file (self, path)
```

allensdk.internal.api.queries.pre_release module

```
class allensdk.internal.api.queries.pre_release.BrainObservatoryApiPreRelease (base_uri=None,
                                                                              dat-
                                                                              acube_uri=None)
```

Bases: *allensdk.api.queries.brain_observatory_api.BrainObservatoryApi*

```
get_cell_metrics (self)
    Get cell metrics by id
```

Parameters

cell_metrics_ids [integer or list of integers, optional] only select specific cell metric records.

Returns

dict [cell metric metadata]

```
get_experiment_containers (self)
    Get experiment container by id
```

Parameters

experiment_container_ids [integer or list of integers, optional] only select specific experiment containers.

Returns

dict [experiment container metadata]

```
get_ophys_experiments (self)
    Get OPhys Experiments by id
```

Parameters

ophys_experiment_ids [integer or list of integers, optional] only select specific experiments.

Returns

dict [ophys experiment metadata]

Module contents

Submodules

allensdk.internal.api.api_prerelease module

```
class allensdk.internal.api.api_prerelease.ApiPrerelease (api_base_url_string=None)
```

Bases: *allensdk.api.api.Api*

Extends allensdk.api.api to copy files ‘locally’ from shared storage.

retrieve_file_from_storage (*self*, *storage_path*, *save_file_path*)
Copy data from path to file_name.

Parameters

storage_path [string] path to file in shared directory (copy source)

save_file_name [string] path to file destination (copy target)

allensdk.internal.api.behavior_data_lims_api module

```
class allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi (behavior_session_id:
    int,
    lims_credentials:
    Optional[allensdk.core.authentication.LimsCredentials] =
    None,
    mtrain_credentials:
    Optional[allensdk.core.authentication.MtrainCredentials] =
    None)

Bases:
    allensdk.core.cache_method_utilities.CachedInstanceMethodMixin,
    allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase
```

get_age (*self*) → str

Returns age code of the subject. :rtype: str

get_behavior_session_id (*self*) → int

Getter to be consistent with BehaviorOphysLimsApi.

get_behavior_session_uuid (*self*) → Union[int, NoneType]

get_behavior_stimulus_file (*self*) → str

Return the path to the StimulusPickle file for a session. :rtype: str

get_birth_date (*self*) → <method 'date' of 'datetime.datetime' objects>

Returns the birth date of the animal. :rtype: datetime.date

get_driver_line (*self*) → List[str]

Returns the genotype name(s) of the driver line(s). :rtype: list

get_experiment_date (*self*) → datetime.datetime

Return timestamp the behavior stimulus file began recording in UTC :rtype: datetime

get_external_specimen_name (*self*) → int

Returns the LabTracks ID :rtype: int

get_full_genotype (*self*) → str

Return the name of the subject's genotype :rtype: str

get_licks (*self*) → pandas.core.frame.DataFrame

Get lick data from pkl file. This function assumes that the first sensor in the list of lick_sensors is the desired lick sensor. If this changes we need to update to get the proper line.

Since licks can occur outside of a trial context, the lick times are extracted from the vsyncs and the frame number in *lick_events*. Since we don't have a timestamp for when in "experiment time" the vsync stream

starts (from `self.get_stimulus_timestamps`), we compute it by fitting a linear regression (frame number x time) for the *start_trial* and *end_trial* events in the *trial_log*, to true up these time streams.

Returns `pd.DataFrame` – A dataframe containing lick timestamps

get_metadata (*self*) → `Dict[str, Any]`

Return metadata about the session. :rtype: dict

get_reporter_line (*self*) → `List[str]`

Returns the genotype name(s) of the reporter line(s). :rtype: list

get_rewards (*self*) → `pandas.core.frame.DataFrame`

Get reward data from pkl file, based on pkl file timestamps (not sync file).

Returns `pd.DataFrame` – A dataframe containing timestamps of delivered rewards.

get_rig_name (*self*) → `str`

Returns the name of the experimental rig. :rtype: str

get_running_data_df (*self*, *lowpass=True*) → `pandas.core.frame.DataFrame`

Get running speed data.

Returns `pd.DataFrame` – dataframe containing various signals used to compute running speed.

get_running_speed (*self*, *lowpass=True*) → `allensdk.brain_observatory.running_speed.RunningSpeed`

Get running speed using timestamps from `self.get_stimulus_timestamps`.

NOTE: Do not correct for monitor delay.

Returns `RunningSpeed` – a `NamedTuple` containing the subject's timestamps and running speeds (in cm/s)

get_sex (*self*) → `str`

Returns sex of the animal (M/F) :rtype: str

get_stimulus_frame_rate (*self*) → `float`

get_stimulus_name (*self*) → `str`

Returns the name of the stimulus set used for the session. :rtype: str

get_stimulus_presentations (*self*) → `pandas.core.frame.DataFrame`

Get stimulus presentation data.

NOTE: Uses timestamps that do not account for monitor delay.

Returns `pd.DataFrame` – Table whose rows are stimulus presentations (i.e. a given image, for a given duration, typically 250 ms) and whose columns are presentation characteristics.

get_stimulus_templates (*self*) → `Dict[str, numpy.ndarray]`

Get stimulus templates (movies, scenes) for behavior session.

Returns

Dict[str, np.ndarray] A dictionary containing the stimulus images presented during the session. Keys are data set names, and values are 3D numpy arrays.

get_stimulus_timestamps (*self*) → `numpy.ndarray`

Get stimulus timestamps (vsyncs) from pkl file. Align to the (frame, time) points in the trial events.

NOTE: Located with `behavior_session_id`. Does not use the `sync_file` which requires `ophys_session_id`.

Returns

np.ndarray Timestamps associated with stimulus presentations on the monitor that do no account for monitor delay.

get_task_parameters (*self*) → dict
Get task parameters from pkl file.

Returns

dict A dictionary containing parameters used to define the task runtime behavior.

get_trials (*self*) → pandas.core.frame.DataFrame
Get trials from pkl file

Returns

pd.DataFrame A dataframe containing behavioral trial start/stop times, and trial data

allensdk.internal.api.behavior_lims_api module

```
class allensdk.internal.api.behavior_lims_api.BehaviorLimsApi (behavior_experiment_id:
                                                                int,
                                                                lims_credentials:
                                                                Optional[allensdk.core.authentication.DbCredentials] = None)

Bases: object

static behavior_session_id_to_foraging_id (behavior_session_id)
    maps behavior_session_id to foraging_id

static foraging_id_to_behavior_session_id (foraging_id)
    maps foraging_id to behavior_session_id

classmethod from_foraging_id (foraging_id: str, lims_credentials:
                                Union[allensdk.core.authentication.DbCredentials, None] = None)

get_behavior_experiment_id (self)

get_behavior_stimulus_file (self)

get_extended_trials (self)
```

allensdk.internal.api.behavior_ophys_api module

```
class allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi (ophys_experiment_id:
                                                                int,
                                                                lims_credentials:
                                                                Optional[allensdk.core.authentication.DbCredentials] = None)

Bases: allensdk.internal.api.ophys_lims_api.OphysLimsApi, allensdk.
        brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase

classmethod get_api_list_by_container_id (container_id)

get_average_projection (self, image_api=None)

get_behavior_session_uuid (self)

get_behavior_stimulus_file (self)

static get_containers_df (only_passed=True)
```

```
get_corrected_fluorescence_traces (self)
get_dff_traces (self)
get_experiment_container_id (self)
get_extended_trials (self)
get_eye_tracking (self, z_threshold: float = 3.0, dilation_frames: int = 2)
get_eye_tracking_filepath (self)
get_licks (self)
get_metadata (self)
get_motion_correction (self)
get_nwb_filepath (self)
static get_ophys_experiment_df ()
get_ophys_frame_rate (self)
get_ophys_timestamps (self)
get_rewards (self)
get_running_data_df (self, lowpass=True)
get_running_speed (self, lowpass=True)
get_stimulus_frame_rate (self)
get_stimulus_presentations (self)
get_stimulus_rebase_function (self)
get_stimulus_templates (self)
get_stimulus_timestamps (self)
get_sync_data (self)
get_sync_licks (self)
get_task_parameters (self)
get_trials (self)
```

allensdk.internal.api.lims_api module

```
class allensdk.internal.api.lims_api.LimsApi (lims_credentials: Optional[allensdk.core.authentication.DbCredentials]
                                              = None)
    Bases: object
    get_behavior_tracking_video_filepath_df (self)
    get_experiment_id (self)
    get_eye_tracking_video_filepath_df (self)
```

allensdk.internal.api.mtrain_api module

```

class allensdk.internal.api.mtrain_api.MtrainApi (api_base='http://mtrain:5000')
    Bases: object

    get_behavior_training_df (self, LabTracks_ID=None)

    get_current_stage (self, LabTracks_ID)

    get_df (self, table_name, get_obj=None, **kwargs)

    get_page (self, table_name, get_obj=None, filters=[], **kwargs)

    get_session (self, behavior_session_uuid=None, behavior_session_id=None)

    get_subjects (self)

class allensdk.internal.api.mtrain_api.MtrainSqlApi (dbname=None, user=None,
                                                    host=None, password=None,
                                                    port=None)

    Bases: object

    get_behavior_training_df (self, LabTracks_ID)

    get_subjects (self)

```

allensdk.internal.api.ophys_lims_api module

```

class allensdk.internal.api.ophys_lims_api.OphysLimsApi (ophys_experiment_id:
                                                         int, lims_credentials: Op-
                                                         tional[allensdk.core.authentication.DbCredentials]
                                                         = None)
    Bases: allensdk.core.cache_method_utilities.CachedInstanceMethodMixin

    get_age (self)

    get_average_intensity_projection_image_file (self)

    get_behavior_session_id (self) → Union[int, NoneType]
        Returns the behavior_session_id associated with this experiment, if applicable.

    get_cell_roi_ids (self)

    get_cell_specimen_table (self)

    get_demix_file (self)

    get_dff_file (self)

    get_driver_line (self)

    get_equipment_id (self)

    get_experiment_date (self)

    get_external_specimen_name (self, ophys_experiment_id=None)

    get_field_of_view_shape (self)

    get_foraging_id (self)

    get_full_genotype (self)

    get_imaging_depth (self)

```

get_imaging_plane_group (*self*) → Union[int, NoneType]
Get the imaging plane group number. This is a numeric index that indicates the order that the frames were acquired when there is more than one frame acquired concurrently. Relevant for mesoscope data timestamps, as the laser jumps between plane groups during the scan. Will be None for non-mesoscope data.

get_max_projection (*self*, *image_api=None*)

get_max_projection_file (*self*)

get_metadata (*self*)

get_motion_corrected_image_stack_file (*self*)

get_nwb_filepath (*self*)

get_objectlist_file (*self*)

get_ophys_cell_segmentation_run_id (*self*)

get_ophys_experiment_dir (*self*)

get_ophys_experiment_id (*self*)

get_plane_group_count (*self*) → int
Gets the total number of plane groups in the session. This is required for resampling ophys timestamps for mesoscope data. Will be 0 if the scope did not capture multiple concurrent frames. See *get_imaging_plane_group* for more info.

get_raw_cell_specimen_table_dict (*self*)

get_raw_dff_data (*self*)

get_reporter_line (*self*)

get_rig_name (*self*)

get_rigid_motion_transform_file (*self*)

get_sex (*self*)

get_stimulus_name (*self*)

get_surface_2p_pixel_size_um (*self*)

get_sync_file (*self*, *ophys_experiment_id=None*)

get_targeted_structure (*self*)

get_workflow_state (*self*)

Module contents

exception allensdk.internal.api.OneOrMoreResultExpectedError
Bases: RuntimeError

class allensdk.internal.api.PostgresQueryMixin (*, *dbname, user, host, password, port*)
Bases: object

fetchall (*self*, *query*, *strict=True*)

fetchone (*self*, *query*, *strict=True*)

get_connection (*self*)

get_cursor (*self*)

select (*self*, *query*)

select_one (*self*, *query*)

`allensdk.internal.api.psycpg2_select` (*query*, *database*, *host*, *port*, *username*, *password*)

allensdk.internal.brain_observatory package

Subpackages

allensdk.internal.brain_observatory.resources package

Module contents

Submodules

allensdk.internal.brain_observatory.annotated_region_metrics module

Module for calculating annotated region metrics from ISI data

`allensdk.internal.brain_observatory.annotated_region_metrics.create_region_mask` (*image_shape*,
x,
y,
width,
height,
mask)

Create mask for region on retinotopic map

Parameters

image_shape [tuple] (height, width) of retinotopic map

x [int] x offset of region mask within retinotopic map

y [int] y offset of region mask within retinotopic map

width [int] width of region mask

height [int] height of region mask

mask [list] region mask as a list of lists

Returns

numpy.ndarray Region mask

`allensdk.internal.brain_observatory.annotated_region_metrics.eccentricity` (*az*,
alt,
az_center,
alt_center)

Compute eccentricity

Parameters

az [numpy.ndarray] Azimuth retinotopic map

alt [numpy.ndarray] Altitude retinotopic map

az_center [float] Azimuth value to use as center of eccentricity map

alt_center [float] Altitude value to use as center of eccentricity map

Returns**numpy.ndarray** Eccentricity map

```
allensdk.internal.brain_observatory.annotated_region_metrics.get_metrics(altitude_phase,
                                                                           az-
                                                                           imuth_phase,
                                                                           x=None,
                                                                           y=None,
                                                                           width=None,
                                                                           height=None,
                                                                           mask=None,
                                                                           al-
                                                                           ti-
                                                                           tude_scale=0.322,
                                                                           az-
                                                                           imuth_scale=0.383)
```

Calculate annotated region metrics

```
allensdk.internal.brain_observatory.annotated_region_metrics.retinotopy_metric(mask,
                                                                           isi_map)
```

Compute retinotopic metrics for a responding area

Parameters**mask** [numpy.ndarray] Mask representing the area over which to calculate metrics**isi_map** [numpy.ndarray] Retinotopic map**Returns****(float, float, float, float) tuple** min, max, range, bias of retinotopic map over masked region**allensdk.internal.brain_observatory.demix_report module**

```
allensdk.internal.brain_observatory.demix_report.background_trace(trace,
                                                                    save_dir,
                                                                    data_set=None)

allensdk.internal.brain_observatory.demix_report.compute_correlations(dm,
                                                                    movie_path,
                                                                    movie_dataset)

allensdk.internal.brain_observatory.demix_report.compute_correlations_without_masks(dm)

allensdk.internal.brain_observatory.demix_report.compute_non_overlap_masks(dm)

allensdk.internal.brain_observatory.demix_report.compute_non_overlap_traces(dm,
                                                                    movie_path,
                                                                    movie_dataset)

allensdk.internal.brain_observatory.demix_report.correlation_report(dm,
                                                                    save_dir,
                                                                    with-
                                                                    out_masks=True)
```

parameters: dm: [DeMix object] without_masks: boolean


```
allensdk.internal.brain_observatory.demix_report.plot_masks(dm, save_dir,
                                                           movie_file,
                                                           movie_dataset,
                                                           window=150,
                                                           add_background=True)
```

allensdk.internal.brain_observatory.demixer module

```
allensdk.internal.brain_observatory.demixer.demix_time_dep_masks(raw_traces,
                                                                    stack,
                                                                    masks)
```

Parameters

- **raw_traces** – extracted traces
- **stack** – movie (same length as traces)
- **masks** – binary roi masks

Returns

 demixed traces

```
allensdk.internal.brain_observatory.demixer.find_negative_baselines(trace)
allensdk.internal.brain_observatory.demixer.find_negative_transients_threshold(trace,
                                                                                win-
                                                                                dow=500,
                                                                                length=10,
                                                                                std_devs=3)
allensdk.internal.brain_observatory.demixer.find_zero_baselines(traces)
allensdk.internal.brain_observatory.demixer.identify_valid_masks(mask_array)
allensdk.internal.brain_observatory.demixer.plot_negative_baselines(raw_traces,
                                                                    demix_traces,
                                                                    mask_array,
                                                                    roi_ids_mask,
                                                                    plot_dir,
                                                                    ext='png')
allensdk.internal.brain_observatory.demixer.plot_negative_transients(raw_traces,
                                                                    demix_traces,
                                                                    valid_roi,
                                                                    mask_array,
                                                                    roi_ids_mask,
                                                                    plot_dir,
                                                                    ext='png')
allensdk.internal.brain_observatory.demixer.plot_overlap_masks_lengthOne(roi_ind,
                                                                    masks,
                                                                    save-
                                                                    file=None,
                                                                    weighted=False)
allensdk.internal.brain_observatory.demixer.plot_traces(raw_trace, demix_trace,
                                                         roi_id, roi_ind, save_file)
```

```
allensdk.internal.brain_observatory.demixer.plot_transients (roi_ind,      t_trans,
                                                            masks,      traces,
                                                            demix_traces, save-
                                                            file)
```

```
allensdk.internal.brain_observatory.demixer.rolling_window (trace, window=500)
```

Parameters

- **trace** –
- **window** –

Returns

allensdk.internal.brain_observatory.eye_calibration module

```
class allensdk.internal.brain_observatory.eye_calibration.EyeCalibration (monitor_position=array([
    8.62,
    3.16]),
    mon-
    i-
    tor_rotations=array([0.,
    0.,
    0.]),
    led_position=array([25.8,
    -
    6.12,
    3.21]),
    cam-
    era_position=array([13.,
    0.,
    0.]),
    cam-
    era_rotations=array([0.,
    0.,
    0.22863813]),
    eye_radius=0.1682,
    cm_per_pixel=0.0010199)
```

Bases: object

Class for performing eye-tracking calibration.

Provides methods for estimating the position of the pupil in 3D space and projecting the gaze onto the monitor in both 3D space and monitor space given the experimental geometry.

Parameters

- monitor_position** [numpy.ndarray] [x,y,z] position of monitor in cm.
- monitor_rotations** [numpy.ndarray] [x,y,z] rotations of monitor in radians.
- led_position** [numpy.ndarray] [x,y,z] position of LED in cm.
- camera_position** [numpy.ndarray] [x,y,z] position of camera in cm.
- camera_rotations** [numpy.ndarray] [x,y,z] rotations for camera in radians. X and Y must be 0.
- eye_radius** [float] Radius of the eye in cm.

cm_per_pixel [float] Pixel size of eye-tracking camera.

compute_area (*self*, *pupil_parameters*)

Compute the area of the pupil.

Assume the pupil is a circle, and that as it moves off-axis with the camera the observed ellipse major axis remains the diameter of the circle.

Parameters

pupil_parameters [numpy.ndarray] [nx5] array of pupil parameters.

Returns

numpy.ndarray [nx1] array of pupil areas in estimated pixels.

static cr_position_in_mouse_eye_coordinates (*led_position*, *eye_radius*)

Determine the 3D position of the corneal reflection.

The eye is modeled as a spherical mirror, so the reflection appears to be half the radius of the eye from the origin along the eye-LED axis.

Parameters

led_position [numpy.ndarray] [x,y,z] position of the LED in eye coordinates.

eye_radius [float] Radius of the eye in centimeters.

Returns

numpy.ndarray [x,y,z] location of the corneal reflection in eye coordinates.

pupil_position_in_mouse_eye_coordinates (*self*, *pupil_parameters*, *cr_parameters*)

Compute the 3D pupil position in mouse eye coordinates.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_parameters [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimates in eye coordinates.

pupil_position_on_monitor_in_cm (*self*, *pupil_parameters*, *cr_parameters*)

Compute the pupil position on the monitor in cm.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_parameters [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimates in eye coordinates.

pupil_position_on_monitor_in_degrees (*self*, *pupil_parameters*, *cr_parameters*)

Get pupil position on monitor measured in visual degrees.

Parameters

pupil_parameters [numpy.ndarray] Array of pupil parameters for each eye tracking frame.

cr_params [numpy.ndarray] Array of corneal reflection parameters for each eye tracking frame.

Returns

numpy.ndarray Pupil position estimate in visual degrees.

`allensdk.internal.brain_observatory.eye_calibration.base_object_to_eye_rotation_matrix(object_position)`
Rotation matrix to rotate base object frame to eye coordinates.

By convention, any other object's coordinate frame before rotations is set with positive Z pointing from the object's position back to the origin of the eye coordinate system, with X parallel to the eye X-Y plane.

Parameters

object_position [np.ndarray] [x, y, z] position of object in eye coordinates.

Returns

numpy.ndarray [3x3] rotation matrix.

`allensdk.internal.brain_observatory.eye_calibration.object_norm_eye_coordinates(object_position, x_rotation, y_rotation, z_rotation)`
Get the normal vector for the object plane in eye coordinates.

Parameters

object_position [numpy.ndarray] [x, y, z] location of the object in eye coordinates.

x_rotation [float] Rotation about the x-axis in radians.

y_rotation [float] Rotation about the y-axis in radians.

z_rotation [float] Rotation about the z-axis in radians.

Returns

numpy.ndarray Endpoint of the object plane vector in eye coordinates.

`allensdk.internal.brain_observatory.eye_calibration.object_rotation_matrix(x_rotation, y_rotation, z_rotation)`
Rotation matrix in object coordinate frame.

The rotation matrix for rotating the object coordinate frame from the initial position. This is done by rotating around x, then around y', then around z''.

Parameters

x_rotation [float] Rotation about x axis in radians.

y_rotation [float] Rotation about y axis in radians.

z_rotation [float] Rotation about z axis in radians.

Returns

numpy.ndarray [3x3] rotation matrix.

`allensdk.internal.brain_observatory.eye_calibration.project_to_plane(plane_normal, plane_point, points)`
Project from the origin through points onto a plane.

Parameters

plane_normal [numpy.ndarray] [x, y, z] normal unit vector to the plane.

plane_point [numpy.ndarray] [x, y, z] point on the plane.

points [numpy.ndarray] [nx3] points in space through which to project.

Returns

numpy.ndarray [nx3] points projected on the plane.

allensdk.internal.brain_observatory.fit_ellipse module

```
class allensdk.internal.brain_observatory.fit_ellipse.FitEllipse (min_points,  
max_iter,  
threshold,  
num_close)
```

Bases: object

choose_inliers (*self, candidate_points*)

fit_ellipse (*self, inlier_points*)

outlier_cost (*self, outlier_points, params*)

ransac_fit (*self, candidate_points*)

allensdk.internal.brain_observatory.fit_ellipse.**ellipse_angle_of_rotation** (*a*)

allensdk.internal.brain_observatory.fit_ellipse.**ellipse_angle_of_rotation2** (*a*)

allensdk.internal.brain_observatory.fit_ellipse.**ellipse_axis_length** (*a*)

allensdk.internal.brain_observatory.fit_ellipse.**ellipse_center** (*a*)

allensdk.internal.brain_observatory.fit_ellipse.**fit_ellipse** (*candidate_points*)

allensdk.internal.brain_observatory.fit_ellipse.**rotate_vector** (*y, x, theta*)

allensdk.internal.brain_observatory.fit_ellipse.**test_fit** ()

allensdk.internal.brain_observatory.frame_stream module

```
class allensdk.internal.brain_observatory.frame_stream.CvInputStream (movie_path,  
num_frames=None,  
block_size=1,  
cache_frames=False)
```

Bases: object

close (*self*)

open (*self*)

```
class allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream(movie_path,
                                                                    frame_shape,
                                                                    ffm-
                                                                    peg_bin='ffmpeg',
                                                                    num_frames=None,
                                                                    block_size=1,
                                                                    cache_frames=False,
                                                                    pro-
                                                                    cess_frame_cb=None)

    Bases: allensdk.internal.brain_observatory.frame_stream.FrameInputStream

    close (self)

    create_images (self, output_directory, image_type)

    open (self)

class allensdk.internal.brain_observatory.frame_stream.FfmpegOutputStream(frame_shape,
                                                                    ffm-
                                                                    peg_bin='ffmpeg',
                                                                    block_size=1)

    Bases: allensdk.internal.brain_observatory.frame_stream.FrameOutputStream

    close (self)

    open (self, movie_path)

class allensdk.internal.brain_observatory.frame_stream.FrameInputStream(movie_path,
                                                                    num_frames=None,
                                                                    block_size=1,
                                                                    cache_frames=False,
                                                                    pro-
                                                                    cess_frame_cb=None)

    Bases: object

    close (self)

    create_images (self, output_directory, image_type)

    open (self)

class allensdk.internal.brain_observatory.frame_stream.FrameOutputStream(block_size=1)

    Bases: object

    close (self)

    open (self, movie_path)

    write (self, frame)

class allensdk.internal.brain_observatory.frame_stream.ImageOutputStream(block_size=1)

    Bases: allensdk.internal.brain_observatory.frame_stream.FrameOutputStream
```

allensdk.internal.brain_observatory.itracker module

allensdk.internal.brain_observatory.itracker_utils module

```
allensdk.internal.brain_observatory.itracker_utils.default_ray(n)

allensdk.internal.brain_observatory.itracker_utils.eccentricity(a1, a2)
```

```
allensdk.internal.brain_observatory.itracker_utils.filter_bad_params(params,
                                                                    frame_width,
                                                                    frame_height)
```

Replace positions outside image with nan

```
allensdk.internal.brain_observatory.itracker_utils.generate_rays(image_array,
                                                                    seed_pixel)
```

```
allensdk.internal.brain_observatory.itracker_utils.initial_cr_point(image_array,
                                                                    bbox=None)
```

bbox is a tuple of (xmin, xmax, ymin, ymax)

```
allensdk.internal.brain_observatory.itracker_utils.initial_pupil_point(image_array,
                                                                    bbox=None)
```

bbox is a tuple of (xmin, xmax, ymin, ymax)

```
allensdk.internal.brain_observatory.itracker_utils.medfilt_custom(x,          ker-
                                                                    nel_size=3)
```

This median filter returns 'nan' whenever any value in the kernel width is 'nan' and the median otherwise

```
allensdk.internal.brain_observatory.itracker_utils.median_absolute_deviation(a,
                                                                    con-
                                                                    sis-
                                                                    tency_constant=1.4)
```

Calculate the median absolute deviation of a univariate dataset.

Parameters

a [numpy.ndarray] Sample data.

consistency_constant [float] Constant to make the MAD a consistent estimator of the population standard deviation (1.4826 for a normal distribution).

Returns

float Median absolute deviation of the data.

```
allensdk.internal.brain_observatory.itracker_utils.post_process_cr(cr_params)
```

This will replace questionable values of the CR x and y position with 'nan'

- 1) threshold ellipse area by 99th percentile area distribution
- 2) median filter using custom median filter
- 3) remove deviations from discontinuous jumps

The 'nan' values likely represent obscured CRs, secondary reflections, merges with the secondary reflection, or visual distortions due to the whisker or deformations of the eye

```
allensdk.internal.brain_observatory.itracker_utils.post_process_pupil(pupil_params)
```

Filter pupil parameters to replace outliers with nan

Parameters

pupil_params [numpy.ndarray] (Nx5) array of pupil parameters [x, y, angle, axis1, axis2].

Returns

numpy.ndarray Pupil parameters with outliers replaced with nan

```
allensdk.internal.brain_observatory.itracker_utils.rotate_ray(ray, theta)
```

```
allensdk.internal.brain_observatory.itracker_utils.sobel_grad(image_array)
```

allensdk.internal.brain_observatory.mask_set module

```
class allensdk.internal.brain_observatory.mask_set.MaskSet (masks)
```

```
    Bases: object
```

```
    close (self, mask_idxs, max_dist)
```

```
    close_sets (self, set_size, max_dist)
```

```
    count
```

```
    detect_duplicates (self, overlap_threshold)
```

```
    detect_unions (self, set_size=2, max_dist=10, threshold=0.7)
```

```
    distance (self, mask_idxs)
```

```
    intersection (self, mask_idxs)
```

```
    intersection_size (self, mask_idxs)
```

```
    mask (self, mask_idx)
```

```
    mask_is_union_of_set (self, mask_idx, set_idxs, threshold)
```

```
    overlap_fraction (self, idx0, idx1)
```

```
    size (self, mask_idx)
```

```
    union (self, mask_idxs)
```

```
    union_size (self, mask_idxs)
```

```
allensdk.internal.brain_observatory.mask_set.bb_dist (bbs)
```

```
allensdk.internal.brain_observatory.mask_set.make_bbs (masks)
```

allensdk.internal.brain_observatory.ophys_session_decomposition module

```
allensdk.internal.brain_observatory.ophys_session_decomposition.export_frame_to_hdf5 (raw_filename,  
                                                                                       data_hdf5_filename,  
                                                                                       auxiliary_hdf5_filename,  
                                                                                       frame_number,  
                                                                                       compression,  
                                                                                       compression_level,  
                                                                                       compression_type,  
                                                                                       compression_options)
```

Export a frame from raw to hdf5.

Data with the channel_description *data* is stored in the *data_hdf5_filename*, while any other data is stored in the *auxiliary_hdf5_filename*

```
allensdk.internal.brain_observatory.ophys_session_decomposition.load_frame (raw_filename,  
                                                                                   json_meta,  
                                                                                   use_memmap=False)
```

Load a frame of a multi-frame raw file.


```
allensdk.internal.brain_observatory.ophys_session_decomposition.open_view_on_binary(file_like,
dtype=<dtype='<float64'>',
'numpy.uint32',
mode='r',
offset=0,
set=0,
shape=None,
or-
der='C',
strides=None)
```

Open a view into a memory-mapped binary file.

Parameters

file_like [{string, file object}] File to open.

dtype [numpy.dtype] Numpy dtype to open the memory-mapped array as.

mode [string] Mode to open the file in.

offset [integer] Offset (in bytes) into the file at which to start the memory map.

shape [(tuple, list)] Shape of the array.

order [{"C", "F"}] C or Fortran ordering.

strides [(tuple, list)] Strides along each axis for reading the array.

Returns

numpy.memmap Strided view into memory-mapped array.

```
allensdk.internal.brain_observatory.ophys_session_decomposition.read_strided(filename,
dtype,
offset,
set,
shape,
strides)
```

Load a frame without memory-mapping.

allensdk.internal.brain_observatory.roi_filter module

allensdk.internal.brain_observatory.roi_filter_utils module

```
allensdk.internal.brain_observatory.roi_filter_utils.CRITERIA()
```

```
class allensdk.internal.brain_observatory.roi_filter_utils.TrainingLabelClassifier(criteria)
Bases: object
```

Very basic threshold_based classifier.

Has a decision function that is just the number of distinct criteria met by the classifier. Criteria are defined as a list of strings used with pandas.DataFrame.eval.

Parameters

criteria [list] List of evaluation strings.

decision_function (*self*, *X*)

Get the distance from the decision boundary.

Parameters

X [array-like] Features for each ROI.

Returns

T [array-like] Distance for each sample from the decision boundary.

class `allensdk.internal.brain_observatory.roi_filter_utils.TrainingMultiLabelClassifier` (*crit*

Bases: `object`

Multilabel classifier using groups of `TrainingLabelClassifiers`.

This was used to generate labeling for training the original SVM for classification.

Parameters

criteria [dictionary] Label names and criteria for each label.

get_eXcluded (*self, X*)

Get the calculated value of the eXcluded column.

This is useful for comparison with the original classifier implementation.

Parameters

X [pandas.DataFrame] Object features from the object list file.

Returns

numpy.ndarray Calculated eXcluded score from the classifier.

label_data (*self, X, as_columns=True*)

Generate labels for each row in X.

Parameters

X [pandas.DataFrame] Object features from the object list file.

Returns

numpy.ndarray Array of label codes representing the combination of labels found for each row.

`allensdk.internal.brain_observatory.roi_filter_utils.calculate_max_border` (*motion_df,*
max_shift)

Calculate motion boundary from frame offsets.

When the motion correction algorithm fails to find sufficient matches, it generates very large frame offsets. The use of *max_shift* avoids filtering too many cells due to the large offsets, with the tradeoff that those frames will be noise.

Parameters

motion_df [pandas.DataFrame] Dataframe containing the x, y offsets from motion correction.

max_shift [float] Maximum shift to allow when considering motion correction. Any larger shifts are considered outliers.

Returns

list [right_shift, left_shift, down_shift, up_shift]

`allensdk.internal.brain_observatory.roi_filter_utils.get_indices_by_distance` (*object_list_points,*
mask_points)

Find indices of nearest neighbor matches.

Require a distance of 0 (perfect match) and a unique match between masks and object_list entries.

`allensdk.internal.brain_observatory.roi_filter_utils.get_rois` (*segmentation_stack*,
border=None)

Extract a list of rois from the segmentation data array.

Parameters

segmentation_stack [numpy.ndarray] The array from the maxInt_masks file showing the object masks.

border [list] [right_shift, left_shift, down_shift, up_shift] bounding box determined from motion correction.

Returns

list List of RoiMask objects.

`allensdk.internal.brain_observatory.roi_filter_utils.order_rois_by_object_list` (*object_data*,
rois)

Reorder rois by matching bounding boxes to object list.

Parameters

object_data [pandas.DataFrame] Object list data.

rois [list] List of RoiMasks.

Returns

list The list of rois reordered to index the same as *object_data*.

`allensdk.internal.brain_observatory.run_itracker` module

`allensdk.internal.brain_observatory.time_sync` module

class `allensdk.internal.brain_observatory.time_sync.OphysTimeAligner` (*sync_file*,
*scan-
ner=None*,
dff_file=None,
*stimu-
lus_pkl=None*,
eye_video=None,
*behav-
ior_video=None*,
long_stim_threshold=0.2)

Bases: `object`

behavior_video_timestamps

corrected_behavior_video_timestamps

corrected_eye_video_timestamps

corrected_ophys_timestamps

corrected_stim_timestamps

dataset

eye_video_timestamps

ophys_timestamps

Get the timestamps for the ophys data.

stim_timestamps

```
allensdk.internal.brain_observatory.time_sync.corrected_video_timestamps(video_name,  
                                                                           times-  
                                                                           tamps,  
                                                                           data_length)
```

```
allensdk.internal.brain_observatory.time_sync.get_alignment_array(ref, other,  
                                                                    int_method=<ufunc  
                                                                    'floor'>)
```

Generate an alignment array

```
allensdk.internal.brain_observatory.time_sync.get_keys(sync_dset:          al-  
                                                         lensdk.brain_observatory.sync_dataset.Dataset)  
                                                         → dict
```

Gets the correct keys for the sync file by searching the sync file line labels. Removes key from the dictionary if it is not in the sync dataset line labels. Args:

sync_dset: The sync dataset to search for keys within

Returns:

key_dict: dictionary of key value pairs for finding data in the sync file

```
allensdk.internal.brain_observatory.time_sync.get_ophys_data_length(filename)
```

```
allensdk.internal.brain_observatory.time_sync.get_photodiode_events(sync_dset,  
                                                                      photodi-  
                                                                      ode_key)
```

Returns the photodiode events with the start/stop indicators and the window init flash stripped off. These transitions occur roughly ~1.0s apart, since the sync square changes state every N frames (where N = 60, and frame rate is 60 Hz). Because there are no markers for when the first transition of this type started, we estimate based on the event intervals. For the first valid event, find the first two events that both meet the following criteria:

The next event occurs ~1.0s later

First the last valid event, find the first two events that both meet the following criteria:

The last valid event occurred ~1.0s before

```
allensdk.internal.brain_observatory.time_sync.get_real_photodiode_events(sync_dset,  
                                                                           pho-  
                                                                           to-  
                                                                           di-  
                                                                           ode_key,  
                                                                           anomaly_threshold=0.5)
```

Gets the photodiode events with the anomalies removed.

```
allensdk.internal.brain_observatory.time_sync.get_stim_data_length(filename:  
                                                                    str) → int
```

Get stimulus data length from .pkl file.

Parameters

filename [str] Path of stimulus data .pkl file.

Returns

int Stimulus data length.

```
allensdk.internal.brain_observatory.time_sync.get_video_length(filename)
```

```
allensdk.internal.brain_observatory.time_sync.monitor_delay(sync_dset,
                                                           stim_times,    pho-
                                                           todiod_key,  transi-
                                                           tion_frame_interval=60,
                                                           max_monitor_delay=0.07)
```

Calculate monitor delay.

Module contents

allensdk.internal.core package

Submodules

allensdk.internal.core.lims_pipeline_module module

```
class allensdk.internal.core.lims_pipeline_module.PipelineModule(description="",
                                                                parser=None)
```

Bases: object

args

input_data (*self*)

write_output_data (*self*, *data*)

```
allensdk.internal.core.lims_pipeline_module.default_argument_parser(description="")
```

```
allensdk.internal.core.lims_pipeline_module.run_module(module,          input_data,
                                                         storage_directory,    op-
                                                         tional_args=None,
                                                         python='/shared/utis.x86_64/python-
                                                         2.7/bin/python',
                                                         sdk_path='/shared/bioapps/infoapps/lims2_modules
                                                         local=False, pbs=None)
```

allensdk.internal.core.lims_utilities module

```
allensdk.internal.core.lims_utilities.append_well_known_file(wkfs,          path,
                                                             wkf_type_id=None,
                                                             con-
                                                             tent_type=None)
```

```
allensdk.internal.core.lims_utilities.connect(user='limsreader',          host='limsdb2',
                                              database='lims2',    password='limsro',
                                              port=5432)
```

```
allensdk.internal.core.lims_utilities.convert_from_titan_linux(file_name)
```

```
allensdk.internal.core.lims_utilities.get_input_json(object_id, object_class, strat-
                                                         egy_class,          host='lims2',
                                                         **kwargs)
```

```
allensdk.internal.core.lims_utilities.get_well_known_file_by_name(wkfs,    file-
                                                                    name)
```

```
allensdk.internal.core.lims_utilities.get_well_known_file_by_type(wkfs,
                                                                    wkf_type_id)
```

```
allensdk.internal.core.lims_utilities.get_well_known_files_by_name(wkfs, file-
                                                                    name)
allensdk.internal.core.lims_utilities.get_well_known_files_by_type(wkfs,
                                                                    wkf_type_id)
allensdk.internal.core.lims_utilities.linux_to_windows(file_name)
allensdk.internal.core.lims_utilities.query(query, user='limsreader', host='limsdb2',
                                             database='lims2', password='limsro',
                                             port=5432)
allensdk.internal.core.lims_utilities.safe_system_path(file_name)
allensdk.internal.core.lims_utilities.select(cursor, query)
```

allensdk.internal.core.mouse_connectivity_cache_prerelease module

```
class allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePre
```

Bases: *allensdk.core.mouse_connectivity_cache.MouseConnectivityCache*

Extends MouseConnectivityCache to use prereleased data from lims.

Parameters

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

ccf_version: string Desired version of the Common Coordinate Framework. This affects the annotation volume (*get_annotation_volume*) and structure masks (*get_structure_mask*). Must be one of (*MouseConnectivityApi.CCF_2015*, *MouseConnectivityApi.CCF_2016*). Default: *MouseConnectivityApi.CCF_2016*

cache: boolean Whether the class should save results of API queries to locations specified in the manifest file. Queries for files (as opposed to metadata) must have a file location. If caching is disabled, those locations must be specified in the function call (e.g. *get_projection_density(file_name='file.nrrd')*).

manifest_file: string File name of the manifest to be read. Default is "mouse_connectivity_manifest.json".

Attributes

resolution: int Resolution of grid data to be downloaded when accessing projection volume, the annotation volume, and the annotation volume. Must be one of (10, 25, 50, 100). Default is 25.

api: MouseConnectivityApiPrerelease instance Used internally to make API queries.

```
EXPERIMENTS_PRERELEASE_KEY = 'EXPERIMENTS_PRERELEASE'
```

```
STORAGE_DIRECTORIES_PRERELEASE_KEY = 'STORAGE_DIRECTORIES_PRERELEASE'
```

```
add_manifest_paths (self, manifest_builder)
```

Construct a manifest for this Cache class and save it in a file.

Parameters

file_name: string File location to save the manifest.

```
filter_experiments (self, experiments, cre=None, injection_structure_ids=None, age=None, gender=None, workflow_state=None, workflows=None, project_code=None)
```

Take a list of experiments and filter them by cre status and injection structure.

Parameters

cre: boolean or list If True, return only cre-positive experiments. If False, return only cre-negative experiments. If None, return all experiments. If list, return all experiments with cre line names in the supplied list. Default None.

injection_structure_ids: list Only return experiments that were injected in the structures provided here. If None, return all experiments. Default None.

age [list] Only return experiments with specimens with ages provided here. If None, return all experiments. Default None.

```
get_experiments (self, dataframe=False, file_name=None, cre=None, injection_structure_ids=None, age=None, gender=None, workflow_state=None, workflows=None, project_code=None)
```

Read a list of experiments.

If caching is enabled, this will save the whole (unfiltered) list of experiments to a file.

Parameters

dataframe: boolean Return the list of experiments as a Pandas DataFrame. If False, return a list of dictionaries. Default False.

file_name: string File name to save/read the structures table. If file_name is None, the file_name will be pulled out of the manifest. If caching is disabled, no file will be saved. Default is None.

allensdk.internal.core.simpletree module

```
class allensdk.internal.core.simpletree.SimpleTree (nodes, node_id_cb, parent_id_cb)
```

Bases: object

```
ancestor_ids (self, nid)
```

```
ancestors (self, nid)
```

```
child_ids (self, nid)
```

```
children (self, nid)
```

```
descendant_ids (self, nid)
```

```
descendants (self, nid)
```

```
node (self, nid)
```

```
node_ids (self)
```

```
nodes (self, nids=None)
```

```
parent (self, nid)
```

```
parent_id(self, nid)
```

allensdk.internal.core.swc module

```
class allensdk.internal.core.swc.Marker(*args, **kwargs)
```

Bases: dict

Simple dictionary class for handling reconstruction marker objects.

```
CUT_DENDRITE = 10
```

```
NO_RECONSTRUCTION = 20
```

```
SPACING = [0.1144, 0.1144, 0.28]
```

```
allensdk.internal.core.swc.read_marker_file(file_name)
```

read in a marker file and return a list of dictionaries

```
allensdk.internal.core.swc.read_swc(file_name)
```

Read in an SWC file and return a Morphology object.

Parameters

file_name: string SWC file name.

Returns

Morphology A Morphology instance.

Module contents

allensdk.internal.ephys package

Submodules

allensdk.internal.ephys.core_feature_extract module

```
allensdk.internal.ephys.core_feature_extract.extract_data(data, nwb_file)
```

```
allensdk.internal.ephys.core_feature_extract.filter_sweeps(sweeps, types=None,
                                                            passed_only=True,
                                                            iclamp_only=True)
```

```
allensdk.internal.ephys.core_feature_extract.filtered_sweep_numbers(sweeps,
                                                                      types=None,
                                                                      passed_only=True,
                                                                      iclamp_only=True)
```

```
allensdk.internal.ephys.core_feature_extract.find_coarse_long_square_amp_delta(sweeps,
                                                                                dec-
                                                                                i-
                                                                                mals=0)
```

Find the delta between amplitudes of coarse long square sweeps. Includes failed sweeps.

```
allensdk.internal.ephys.core_feature_extract.find_stim_start(stim, idx0=0)
```

Find the index of the first nonzero positive or negative jump in an array.

Parameters

stim: np.ndarray Array to be searched

idx0: int Start searching with this index (default: 0).

Returns

int

```
allensdk.internal.ephys.core_feature_extract.find_sweep_stim_start (data_set,
                                                                    sweep_number)
allensdk.internal.ephys.core_feature_extract.generate_output_cell_features (cell_features,
                                                                    sweep_features,
                                                                    sweep_index)
allensdk.internal.ephys.core_feature_extract.nan_get (obj, key)
    Return a value from a dictionary. If it does not exist, return None. If it is NaN, return None
allensdk.internal.ephys.core_feature_extract.save_qc_figures (qc_fig_dir,
                                                                    nwb_file,      out-
                                                                    put_data,
                                                                    plot_cell_figures)
allensdk.internal.ephys.core_feature_extract.update_output_sweep_features (cell_features,
                                                                    sweep_features,
                                                                    sweep_index)
```

allensdk.internal.ephys.plot_qc_figures module

```
allensdk.internal.ephys.plot_qc_figures.exp_curve (x, a, inv_tau, y0)
    Function used for tau curve fitting
allensdk.internal.ephys.plot_qc_figures.get_features (sweep_features,
                                                                    sweep_number)
allensdk.internal.ephys.plot_qc_figures.get_spikes (sweep_features, sweep_number)
allensdk.internal.ephys.plot_qc_figures.get_time_string ()
allensdk.internal.ephys.plot_qc_figures.load_experiment (file_name,
                                                                    sweep_number)
allensdk.internal.ephys.plot_qc_figures.main ()
allensdk.internal.ephys.plot_qc_figures.make_cell_html (image_files,
                                                                    ephys_roi_result, file_name,
                                                                    relative_sweep_link)
allensdk.internal.ephys.plot_qc_figures.make_cell_page (nwb_file, ephys_roi_result,
                                                                    working_dir,
                                                                    save_cell_plots=True)
allensdk.internal.ephys.plot_qc_figures.make_sweep_html (sweep_files, file_name)
allensdk.internal.ephys.plot_qc_figures.make_sweep_page (nwb_file, ephys_roi_result,
                                                                    working_dir)
allensdk.internal.ephys.plot_qc_figures.mask_nulls (data)
allensdk.internal.ephys.plot_qc_figures.plot_cell_figures (nwb_file,
                                                                    ephys_roi_result,
                                                                    image_dir, sizes)
```

```
allensdk.internal.ephys.plot_qc_figures.plot_fi_curve_figures(nwb_file,  
                                                                cell_features,  
                                                                lims_features,  
                                                                sweep_features,  
                                                                image_dir, sizes,  
                                                                cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_hero_figures(nwb_file, cell_features,  
                                                            lims_features,  
                                                            sweep_features,  
                                                            image_dir, sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_images(ephys_roi_result, image_dir,  
                                                    sizes, image_sets)  
  
allensdk.internal.ephys.plot_qc_figures.plot_instantaneous_threshold_thumbnail(nwb_file,  
                                                                                sweep_numbers,  
                                                                                cell_features,  
                                                                                lims_features,  
                                                                                sweep_features,  
                                                                                color='red')  
  
allensdk.internal.ephys.plot_qc_figures.plot_long_square_summary(nwb_file,  
                                                                    cell_features,  
                                                                    lims_features,  
                                                                    sweep_features)  
  
allensdk.internal.ephys.plot_qc_figures.plot_ramp_figures(nwb_file,  
                                                            cell_specimen,  
                                                            cell_features,  
                                                            lims_features,  
                                                            sweep_features,  
                                                            image_dir, sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_rheo_figures(nwb_file, cell_features,  
                                                            lims_features,  
                                                            sweep_features,  
                                                            image_dir, sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_sag_figures(nwb_file, cell_features,  
                                                            lims_features,  
                                                            sweep_features,  
                                                            image_dir, sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures.plot_short_square_figures(nwb_file,  
                                                                    cell_features,  
                                                                    lims_features,  
                                                                    sweep_features,  
                                                                    image_dir,  
                                                                    sizes,  
                                                                    cell_image_files)
```

```

allensdk.internal.ephys.plot_qc_figures.plot_single_ap_values (nwb_file,
                                                             sweep_numbers,
                                                             lims_features,
                                                             sweep_features,
                                                             cell_features,
                                                             type_name)

allensdk.internal.ephys.plot_qc_figures.plot_subthreshold_long_square_figures (nwb_file,
                                                                                  cell_features,
                                                                                  lims_features,
                                                                                  sweep_features,
                                                                                  im-
                                                                                  age_dir,
                                                                                  sizes,
                                                                                  cell_image_files)

allensdk.internal.ephys.plot_qc_figures.plot_sweep_figures (nwb_file,
                                                             ephys_roi_result,
                                                             image_dir, sizes)

allensdk.internal.ephys.plot_qc_figures.plot_sweep_set_summary (nwb_file, high-
                                                                light_sweep_number,
                                                                sweep_numbers,
                                                                high-
                                                                light_color='#0779BE',
                                                                back-
                                                                ground_color='#dddddd')

allensdk.internal.ephys.plot_qc_figures.plot_sweep_value_figures (cell_specimen,
                                                                    image_dir,
                                                                    sizes,
                                                                    cell_image_files)

allensdk.internal.ephys.plot_qc_figures.save_figure (fig, image_name, im-
                                                         age_set_name, image_dir,
                                                         sizes, image_sets, scalew=1,
                                                         scaleh=1, ext='jpg')

```

allensdk.internal.ephys.plot_qc_figures3 module

```

allensdk.internal.ephys.plot_qc_figures3.exp_curve (x, a, inv_tau, y0)
    Function used for tau curve fitting

allensdk.internal.ephys.plot_qc_figures3.get_features (sweep_features,
                                                         sweep_number)

allensdk.internal.ephys.plot_qc_figures3.get_spikes (sweep_features, sweep_number)

allensdk.internal.ephys.plot_qc_figures3.get_time_string ()

allensdk.internal.ephys.plot_qc_figures3.load_experiment (file_name,
                                                            sweep_number)

allensdk.internal.ephys.plot_qc_figures3.make_cell_html (image_files, file_name,
                                                            relative_sweep_link,
                                                            specimen_info, fields)

```

```
allensdk.internal.ephys.plot_qc_figures3.make_cell_page(nwb_file,    cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         sweep_info,
                                                         well_known_files,
                                                         specimen_info,    work-
ing_dir,    fields_to_show,
                                                         save_cell_plots=True)
```

nwb_file: name of nwb file (string)

cell_features:

rheo_features: dict containing extracted features from rheobase sweep

sweep_features:

sweep_info:

well_known_files: LIMS-output information containing graphics file names

working_dir:

save_cell_plots:

```
allensdk.internal.ephys.plot_qc_figures3.make_sweep_html(sweep_files, file_name)
```

```
allensdk.internal.ephys.plot_qc_figures3.make_sweep_page(nwb_file,    working_dir,
                                                         sweep_data)
```

```
allensdk.internal.ephys.plot_qc_figures3.mask_nulls(data)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_cell_figures(nwb_file,
                                                         cell_features,
                                                         sweep_features,
                                                         rheo_features,    im-
age_dir,    sweep_info,
                                                         sizes)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_fi_curve_figures(nwb_file,
                                                         cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         image_dir,
                                                         sizes,
                                                         cell_image_files)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_hero_figures(nwb_file,
                                                         cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         image_dir,    sizes,
                                                         cell_image_files)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_images(well_known_files,    image_dir,
                                                         sizes, image_sets)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_instantaneous_threshold_thumbnail(nwb_file,
                                                         sweep_numbers,
                                                         cell_features,
                                                         rheo_features,
                                                         sweep_features,
                                                         color='red')
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_long_square_summary(nwb_file,  
                                                                    cell_features,  
                                                                    rheo_features,  
                                                                    sweep_features)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_ramp_figures(nwb_file, sweep_info,  
                                                            cell_features,  
                                                            rheo_features,  
                                                            sweep_features,  
                                                            image_dir,      sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_rheo_figures(nwb_file,  
                                                            cell_features,  
                                                            rheo_features,  
                                                            sweep_features,  
                                                            image_dir,      sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_sag_figures(nwb_file, cell_features,  
                                                            rheo_features,  
                                                            sweep_features,  
                                                            image_dir,      sizes,  
                                                            cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_short_square_figures(nwb_file,  
                                                                    cell_features,  
                                                                    rheo_features,  
                                                                    sweep_features,  
                                                                    im-  
                                                                    age_dir,  
                                                                    sizes,  
                                                                    cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_single_ap_values(nwb_file,  
                                                                    sweep_numbers,  
                                                                    rheo_features,  
                                                                    sweep_features,  
                                                                    cell_features,  
                                                                    type_name)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_subthreshold_long_square_figures(nwb_file,  
                                                                    cell_features,  
                                                                    rheo_features,  
                                                                    sweep_features,  
                                                                    im-  
                                                                    age_dir,  
                                                                    sizes,  
                                                                    cell_image_files)  
  
allensdk.internal.ephys.plot_qc_figures3.plot_sweep_figures(nwb_file,  
                                                                    sweep_data,      im-  
                                                                    age_dir, sizes)
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_sweep_set_summary(nwb_file, high-  
                                                                light_sweep_number,  
                                                                sweep_numbers,  
                                                                high-  
                                                                light_color='#0779BE',  
                                                                back-  
                                                                ground_color='#dddddd')
```

```
allensdk.internal.ephys.plot_qc_figures3.plot_sweep_value_figures(sweep_info,  
                                                                    image_dir,  
                                                                    sizes,  
                                                                    cell_image_files)
```

```
allensdk.internal.ephys.plot_qc_figures3.save_figure(fig, image_name, im-  
                                                       age_set_name, image_dir,  
                                                       sizes, image_sets, scalew=1,  
                                                       scaleh=1, ext='jpg')
```

Module contents

allensdk.internal.model package

Subpackages

allensdk.internal.model.biophysical package

Subpackages

allensdk.internal.model.biophysical.fits package

Subpackages

allensdk.internal.model.biophysical.fits.fit_styles package

Module contents

Module contents

allensdk.internal.model.biophysical.passive_fitting package

Subpackages

allensdk.internal.model.biophysical.passive_fitting.passive package

Module contents

Submodules

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.arg_parser()
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.main()
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit.process_inputs(parser)
```

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit2 module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit2.main()
```

allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit_elec module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit_elec.main()
```

allensdk.internal.model.biophysical.passive_fitting.neuron_utils module

```
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.get_h()
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.load_morphology(filename)
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.parse_neuron_output(output_s)
allensdk.internal.model.biophysical.passive_fitting.neuron_utils.read_neuron_fit_stdout(functi
```

allensdk.internal.model.biophysical.passive_fitting.output_grabber module

```
class allensdk.internal.model.biophysical.passive_fitting.output_grabber.OutputGrabber(stream, thread)
```

Bases: object

Class used to grab standard output or another stream.

escape_char = '\x08'

readOutput(self)

Read the stream data (one byte at a time) and save the text in *capturedtext*.

start(self)

Start capturing the stream data.

stop(self)

Stop capturing the stream data and save the text in *capturedtext*.

allensdk.internal.model.biophysical.passive_fitting.preprocess module

```
allensdk.internal.model.biophysical.passive_fitting.preprocess.get_cap_check_indices(i)
allensdk.internal.model.biophysical.passive_fitting.preprocess.get_passive_fit_data(cap_checks, data_set)
allensdk.internal.model.biophysical.passive_fitting.preprocess.main()
```

Module contents

Submodules

`allensdk.internal.model.biophysical.biophysical_archiver` module

```
class allensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver (archive_dir)
    Bases: object
    archive_cell (self, ephys_result_id, specimen_id, template, neuronal_model_id)
    get_cells (self)
    get_neuronal_models (self, specimen_ids)
    get_stimulus_file (self, neuronal_model_id)
    get_template_names (self)
```

`allensdk.internal.model.biophysical.check_fi_shift` module

```
allensdk.internal.model.biophysical.check_fi_shift.calculate_fi_curves (data_set,
                                                                    sweeps)
allensdk.internal.model.biophysical.check_fi_shift.estimate_fi_shift (data_set,
                                                                    sweeps)
```

`allensdk.internal.model.biophysical.deap_utils` module

```
class allensdk.internal.model.biophysical.deap_utils.Utils (description)
    Bases: allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils
    actual_parameters_from_normalized (self, params)
    calculate_feature_errors (self, t_ms, v, i)
    generate_morphology (self, morph_filename)
    insert_iclamp (self)
    load_cell_parameters (self)
    normalize_actual_parameters (self, params)
    record_values (self)
    set_actual_parameters (self, params)
    set_iclamp_params (self, amp, delay, dur)
    set_normalized_parameters (self, params)
```

`allensdk.internal.model.biophysical.ephys_utils` module

```
allensdk.internal.model.biophysical.ephys_utils.get_step_stim_characteristics (i,
                                                                    t)
allensdk.internal.model.biophysical.ephys_utils.get_sweep_v_i_t_from_set (data_set,
                                                                    sweep_number)
```


`allensdk.internal.model.biophysical.ephys_utils.get_sweeps_of_type` (*sweep_type*,
sweeps)

allensdk.internal.model.biophysical.fit_stage_1 module

allensdk.internal.model.biophysical.fit_stage_2 module

allensdk.internal.model.biophysical.make_deap_fit_json module

class `allensdk.internal.model.biophysical.make_deap_fit_json.Report` (*top_level_description*,
fit_type)

Bases: object

best_fit_value (*self*)

check_org_selections_for_noise_block (*self*)

gather_from_seeds (*self*)

generate_fit_file (*self*)

make_fit_json_file (*self*)

setup_model (*self*)

allensdk.internal.model.biophysical.neuron_parallel module

allensdk.internal.model.biophysical.optimize module

allensdk.internal.model.biophysical.run_optimize module

class `allensdk.internal.model.biophysical.run_optimize.RunOptimize` (*input_json*,
out-put_json)

Bases: object

copy_local (*self*)

generate_manifest_lims (*self*, *lims_json_path*, *manifest_path*)

generate_manifest_rma (*self*, *neuronal_model_id*, *manifest_path*, *api_url=None*)

info (*self*, *lims_json_path*)

return a string that a bash script can use to find the working directory, etc. to clean up.

load_manifest (*self*)

make_fit (*self*)

nrnivmodl (*self*)

start_specimen (*self*)

`allensdk.internal.model.biophysical.run_optimize.main` (*command*, *input_json*, *out-put_json*)

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.internal.model.biophysical.run_optimize_workflow module**allensdk.internal.model.biophysical.run_passive_fit module**

`allensdk.internal.model.biophysical.run_passive_fit.main(limit, manifest_path)`

`allensdk.internal.model.biophysical.run_passive_fit.run_passive_fit(description)`

allensdk.internal.model.biophysical.run_simulate_lims module

class `allensdk.internal.model.biophysical.run_simulate_lims.RunSimulateLims` (*input_json, out-put_json*)

Bases: `allensdk.model.biophysical.run_simulate.RunSimulate`

copy_local (*self*)

generate_manifest_lims (*self, lims_data_path, manifest_path*)

generate_manifest_rma (*self, neuronal_model_run_id, manifest_path, api_url=None*)

`allensdk.internal.model.biophysical.run_simulate_lims.main` (*command, lims_strategy_json, lims_response_json*)

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.internal.model.biophysical.run_simulate_workflow module**Module contents****allensdk.internal.model.glif package****Submodules****allensdk.internal.model.glif.ASGLM module****allensdk.internal.model.glif.MLIN module**

`allensdk.internal.model.glif.MLIN.MLIN` (*voltage, current, res, cap, dt, MAKE_PLOT=False, SHOW_PLOT=False, BLOCK=False, PUBLICATION_PLOT=False*)

voltage, current input:

voltage: numpy array of voltage with test pulse cut out
current: numpy array of stimulus with test pulse cut out

`allensdk.internal.model.glif.MLIN.autocorr` (*x*)

`allensdk.internal.model.glif.MLIN.exp_decay` (*time, amp, tau*)

`allensdk.internal.model.glif.MLIN.expsymm_cdf` (*v, dv*)

`allensdk.internal.model.glif.MLIN.expsymm_pdf` (*v, dv*)

`allensdk.internal.model.glif.MLIN.find_bin_center(edges)`

`allensdk.internal.model.glif.are_two_lists_of_arrays_the_same` module

`allensdk.internal.model.glif.are_two_lists_of_arrays_the_same.are_two_lists_of_arrays_the_same`

returns False if to lists of arrays are different. otherwise the function returns True.

`allensdk.internal.model.glif.configure_model` module

`allensdk.internal.model.glif.error_functions` module

`allensdk.internal.model.glif.error_functions.MLIN_list_error(param_guess, experiment, input_data)`

`allensdk.internal.model.glif.find_spikes` module

`allensdk.internal.model.glif.find_spikes.align_and_cut_spikes(voltage_list, current_list, dt, spike_window=None)`

This function aligns the spikes to some criteria and returns a current and voltage trace of of the spike over a time window. Also returns zero crossing, and threshold in reference to the aligned spikes.

`allensdk.internal.model.glif.find_spikes.find_spikes_list(voltage_list, dt)`

`allensdk.internal.model.glif.find_spikes.find_spikes_list_old(voltage_list, dt)`

`allensdk.internal.model.glif.find_spikes.find_spikes_old(v, dt)`

`allensdk.internal.model.glif.find_spikes.find_spikes_ssq_list(voltage_list, dt, dv_cutoff, thresh_frac)`

`allensdk.internal.model.glif.find_sweeps` module

exception `allensdk.internal.model.glif.find_sweeps.MissingSweepException`

Bases: Exception

`allensdk.internal.model.glif.find_sweeps.find_long_square_sweeps(sweeps)`

`allensdk.internal.model.glif.find_sweeps.find_noise_sweeps(sweeps)`

Find 1) the noise1 sweeps 2) the noise2 sweeps 4) all noise sweeps

`allensdk.internal.model.glif.find_sweeps.find_ramp_sweeps(sweeps)`

Find 1) all ramp sweeps

2) all subthreshold ramps

3) all superthreshold ramps

`allensdk.internal.model.glif.find_sweeps.find_ramp_to_rheo_sweeps(sweeps)`

`allensdk.internal.model.glif.find_sweeps.find_ranked_sweep(sweep_list, key, reverse=False)`

```
allensdk.internal.model.glif.find_sweeps.find_short_square_sweeps(sweeps)
```

Find 1) all of the subthreshold short square sweeps

2) all of the superthreshold short square sweeps

3) the subthresholds short square sweep with maximum stimulus amplitude

```
allensdk.internal.model.glif.find_sweeps.find_sweeps(sweep_list)
```

```
allensdk.internal.model.glif.find_sweeps.get_sweep_numbers(sweep_list)
```

```
allensdk.internal.model.glif.find_sweeps.get_sweeps_by_name(sweeps,
                                                             sweep_type)
```

```
allensdk.internal.model.glif.find_sweeps.main()
```

```
allensdk.internal.model.glif.find_sweeps.organize_sweeps_by_name(sweeps,
                                                                  name)
```

```
allensdk.internal.model.glif.find_sweeps.parse_arguments()
```

allensdk.internal.model.glif.glif_experiment module

```
class allensdk.internal.model.glif.glif_experiment.GlifExperiment(neuron, dt,
                                                                  stim_list,
                                                                  resp_list,
                                                                  spike_time_steps,
                                                                  grid_spike_times,
                                                                  grid_spike_voltages,
                                                                  param_fit_names,
                                                                  **kwargs)
```

Bases: object

neuron_parameter_count (*self*)

run (*self*, *param_guess*)

This code will run the loaded neuron model in reference to the target neuron spikes. inputs:

self: is the instance of the neuron model and parameters alone with the values of the target spikes.

NOTE the values in each array of the self.gridSpikeIndexTarge_list and the self.interpolated_spike_times are in reference to the time start of of the stim in each induvidual array (not the universal time)

param_guess: array of scalars of the values that will be inserted into the mapping function below.

returns:

voltage_list: list of array of voltage values. **NOTE: IF THE MODEL NEURON SPIKES BEFORE THE TARG**
NOT BE CALCULATED THEREFORE THE RESULTING VECTOR WILL NOT BE AS
LONG AS THE TARGET AND ALSO WILL NOT MAKE SENSE WITH THE STIMULUS
UNLESS YOU CUT IT AND OUTPUT IT TOO.

grid_spike_times_list: interpolated_spike_time_list: an array of the actual times of the spikes.
NOTE: THESE TIMES ARE CALCULATED BY ADDING THE

TIME OF THE INDIVIDUAL SPIKE TO THE TIME OF THE LAST SPIKE.

gridISIFromLastTargSpike_list: list of arrays of spike times of the model in reference to the last target (biologi
spike (not in reference to sweep start)

interpolatedISIFromLastTargSpike_list: list of arrays of spike times of the model in reference to the last target spike (not in reference to sweep start)

voltageOfModelAtGridBioSpike_list: list of arrays of scalars that contain the voltage of the model neuron when the target or bio neuron spikes. **theshOfModelAtGridBioSpike_list:** list of arrays of scalars that contain the threshold of the model neuron when the target or bio neuron spikes.

run_base_model (*self*, *param_guess*)

This code will run the loaded neuron model. inputs:

self: is the instance of the neuron model and parameters alone with the values of the target spikes.

NOTE the values in each array of the `self.gridSpikeIndexTarge_list` and the `self.interpolated_spike_times` are in reference to the time start of of the stim in each induvidual array (not the universal time)

`param_guess`: array of scalars of the values that will be inserted into the mapping function below.

returns:

voltage_list: list of array of voltage values. **NOTE: IF THE MODEL NEURON SPIKES BEFORE THE TARG**
NOT BE CALCULATED THEREFORE THE RESULTING VECTOR WILL NOT BE AS
LONG AS THE TARGET AND ALSO WILL NOT MAKE SENSE WITH THE STIMULUS
UNLESS YOU CUT IT AND OUTPUT IT TOO.

`gridTime_list`: `interpolatedTime_list`: an array of the actual times of the spikes. NOTE: THESE
TIMES ARE CALCULATED BY ADDING THE

TIME OF THE INDIVIDUAL SPIKE TO THE TIME OF THE LAST SPIKE.

grid_ISI_list: list of arrays of spike times of the model in reference to the last target (biological)
spike (not in reference to sweep start)

interpolated_ISI_list: list of arrays of spike times of the model in reference to the last target (biological)
spike (not in reference to sweep start)

`grid_spike_voltage_list`: list of arrays of scalars that contain the voltage of the model neuron when
the target or bio neuron spikes. `grid_spike_threshold_list`: list of arrays of scalars that contain the
threshold of the model neuron when the target or bio neuron spikes.

set_neuron_parameters (*self*, *param_guess*)

Maps the parameter guesses to the coefficients of the model. input:

`param_guess` is vector of values. It is assumed that the length will be

allensdk.internal.model.glif.glif_optimizer module

```
class allensdk.internal.model.glif.glif_optimizer.GlifOptimizer (experiment, dt,  
outer_iterations,  
inner_iterations,  
sigma_outer,  
sigma_inner,  
param_fit_names,  
stim, xtol,  
ftol, internal_iterations,  
bessel, error_function=None,  
error_function_data=None,  
init_params=None)
```

Bases: `object`

evaluate (*self, x, dt_multiplier=100*)

initiate_unique_seed (*self, seed=None*)

randomize_parameter_values (*self, values, sigma*)

run_many (*self, iteration_finished_callback=None, seed=None*)

run_once (*self, param0*)

@param *param0*: a list of the initial guesses for the optimizer @return: tuple including parameters that optimize function and value - see `fmin` docs

run_once_bound (*self, low_bound, high_bound*)

@param *low_bound*: a scalar initial guess for the optimizer @param *high_bound*: a scalar high bound for the optimizer @return: tuple including parameters that optimize function and value - see `fmin` docs

to_dict (*self*)

allensdk.internal.model.glif.glif_optimizer_neuron module

```
exception allensdk.internal.model.glif.glif_optimizer_neuron.GlifBadInitializationException
```

Bases: `Exception`

Exception raised when voltage is above threshold at the beginning of a sweep. i.e. probably caused by the optimizer.

```
exception allensdk.internal.model.glif.glif_optimizer_neuron.GlifNeuronException (message,  
data)
```

Bases: `Exception`

Exception for catching simulation errors and reporting intermediate data.

```
class allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron (*args,  
**kwargs)
```

Bases: `allensdk.model.glif.glif_neuron.GlifNeuron`

Contains methods for running the neuron model in a “forced-spike” paradigm used during optimization.

TYPE = `'GLIF'`

classmethod `from_dict(d)`

classmethod `from_dict_legacy(d)`

run_until_biological_spike (*self, voltage_t0, threshold_t0, AScurrents_t0, stimulus, response, start_index, after_end_index, bio_spike_time_steps*)

Run the neuron simulation over a segment of a stimulus given initial conditions for use in the “forced spike” optimization paradigm. [Note: the section of stimulus is meant to be between two biological neuron spikes. Thus the stimulus is during the interspike interval (ISI)]. The model is simulated until either the model spikes or the end of the segment is reached. If the model does not spike, a spike time is extrapolated past the end of the simulation segment.

This function also returns the initial conditions for the subsequent stimulus segment. In the forced spike paradigm there are several ways

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

stimulus [np.ndarray] the full stimulus array (not just the segment of data being simulated)

response [np.ndarray] the full response array (not just the segment of data being simulated)

start_index [int] index of global stimulus at which to start simulation

after_end_index [int] index of global stimulus *after* the last index to be simulated

bio_spike_time_steps [list] time steps of input spikes

Returns

dict

a dictionary containing: ‘voltage’: simulated voltage value ‘threshold’: simulated threshold values ‘AScurrent_matrix’: afterspike current values during the simulation ‘grid_model_spike_time’: model spike time (in units of dt) ‘interpolated_model_spike_time’: model spike time (in units of dt) interpolated between time steps ‘voltage_t0’: reset voltage value to be used in subsequent simulation interval ‘threshold_t0’: reset threshold value to be used in subsequent simulation interval ‘AScurrents_t0’: reset afterspike current value to be used in subsequent simulation interval ‘grid_bio_spike_model_voltage’: model voltage at the time of the input spike ‘grid_bio_spike_model_threshold’: model threshold at the time of the input spike

run_with_biological_spikes (*self, stimulus, response, bio_spike_time_steps*)

Run the neuron simulation over a stimulus, but do not allow the model to spike on its own. Rather, force the simulation to spike and reset at a given set of spike indices. Dynamics rules are applied between spikes regardless of the simulated voltage and threshold values. Reset rules are applied only at input spike times. This is used during optimization to force the model to follow the spikes of biological data. The model is optimized in this way so that history effects due to spiking can be adequately modeled. For example, every time the model spikes a new set of afterspike currents will be initiated. To ensure that afterspike currents can be optimized, we force them to be initiated at the time of the biological spike.

Parameters

stimulus [np.ndarray] vector of scalar current values

responses [np.ndarray] vector of scalar voltage values

bio_spike_time_steps [list] spike time step indices

Returns

dict

a dictionary containing: 'voltage': simulated voltage values, 'threshold': simulated threshold values, 'AScurrent_matrix': afterspike currents during the simulation, 'grid_model_spike_times': spike times of the model aligned to the simulation grid (when it would have spiked), 'interpolated_model_spike_times': spike times of the model linearly interpolated between time steps, 'grid_ISI': interspike interval between grid model spike times, 'interpolated_ISI': interspike interval between interpolated model spike times, 'grid_bio_spike_model_voltage': voltage of the model at biological/input spike times, 'grid_bio_spike_model_threshold': voltage of the model at biological/input spike times interpolated between time steps

to_dict (*self*)

Convert the neuron to a serializable dictionary.

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_model_spike_from_endpoints (n
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_model_spike_from_endpoints_s
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_spike_time (dt,  
                                                                           num_time_steps,  
                                                                           thresh-  
                                                                           old_t0,  
                                                                           thresh-  
                                                                           old_t1,  
                                                                           volt-  
                                                                           age_t0,  
                                                                           volt-  
                                                                           age_t1)
```

Given two voltage and threshold values and an interval between them, extrapolate a spike time by intersecting lines the thresholds and voltages.


```
allensdk.internal.model.glif.glif_optimizer_neuron.extrapolate_spike_voltage(dt,
                                                                    num_time_steps,
                                                                    thresh-
                                                                    old_t0,
                                                                    thresh-
                                                                    old_t1,
                                                                    volt-
                                                                    age_t0,
                                                                    volt-
                                                                    age_t1)
```

Given two voltage and threshold values and an interval between them, extrapolate a spike time by intersecting lines the thresholds and voltages.

```
allensdk.internal.model.glif.glif_optimizer_neuron.find_first_model_spike(voltage,
                                                                    thresh-
                                                                    old,
                                                                    volt-
                                                                    age_t1,
                                                                    thresh-
                                                                    old_t1,
                                                                    dt)
```

```
allensdk.internal.model.glif.glif_optimizer_neuron.interpolate_spike_voltage(dt,
                                                                    time_step,
                                                                    thresh-
                                                                    old_t0,
                                                                    thresh-
                                                                    old_t1,
                                                                    volt-
                                                                    age_t0,
                                                                    volt-
                                                                    age_t1)
```

Given two voltage and threshold values, the dt between them and the initial time step, interpolate a spike time within the dt interval by intersecting the two lines.

allensdk.internal.model.glif.optimize_neuron module

allensdk.internal.model.glif.plotting module

Written by Corinne Teeter 3-31-14

```
allensdk.internal.model.glif.plotting.checkPreprocess(originalStim_list,    pro-
                                                                    cessedStim_list,    orig-
                                                                    inalVoltage_list,    pro-
                                                                    cessedVoltage_list,    config,
                                                                    blockME=False)
```

```
allensdk.internal.model.glif.plotting.checkSpikeCutting(originalStim_list,    cut-
                                                                    Stim_list,    originalVolt-
                                                                    age_list,    cutVoltage_list,
                                                                    allindOfNonSpiking_list,
                                                                    config, blockME=False)
```

```
allensdk.internal.model.glif.plotting.plotLineRegress1(slope, intercept, r, xlim)
```

```
allensdk.internal.model.glif.plotting.plotLineRegressRed(slope, intercept, r, xlim)
```

```
allensdk.internal.model.glif.plotting.plotSpikes(voltage_list, spike_ind_list, dt,  
                                                  blockME=False, method=False)
```

allensdk.internal.model.glif.preprocess_neuron module

allensdk.internal.model.glif.rc module

```
allensdk.internal.model.glif.rc.least_squares_RCEl_calc_tested(voltage_list,  
                                                                current_list, dt)
```

Calculate resistance, capacitance and resting potential by performing least squares on current and voltage.

Parameters

voltage_list: list of arrays voltage responses for several sweep repeats

current_list: list of arrays current injections for several sweep repeats

dt: float time step size in voltage and current traces

Returns

r_list: list of floats each value corresponds to the resistance of a sweep

c_list: list of floats each value corresponds to the capacitance of a sweep

el_list: list of floats each value corresponds to the resting potential of a sweep

allensdk.internal.model.glif.spike_cutting module

```
allensdk.internal.model.glif.spike_cutting.calc_spike_cut_and_v_reset_via_expvar_residuals
```

This function calculates where the spike should be cut based on explained variance. The goal is to find a model where the voltage after a spike maximally explains the voltage before a spike. This will also specify the voltage reset rule inputs:

spike_determination_method: string specifying the method used to find threshold
all_current_list: list of current (list of current traces injected into neuron)
all_voltage_list: list of voltages (list of voltage trace)

The change is that if the slope is greater than one or intercept is greater than zero it forces it. Regardless of required force the residuals are used.

```
allensdk.internal.model.glif.spike_cutting.plotLineRegress1(slope, intercept, r,  
                                                            xlim)
```

```
allensdk.internal.model.glif.spike_cutting.plotLineRegressRed(slope, intercept, r,  
                                                                xlim)
```

allensdk.internal.model.glif.threshold_adaptation module

```
allensdk.internal.model.glif.threshold_adaptation.calc_spike_component_of_threshold_from_m
```

Calculate the spike components of the threshold by fitting a decaying exponential function to data to threshold versus time since last spike in the multiblip data. The exponential is forced to decay to the local `th_inf` (calculated as the mean all of the threshold values of the first spikes in each individual triblip stimulus). For each multiblip stimulus in a stimulus set if there is more than one spike the difference in voltages from the first and second spike are plotted versus the separation in time. Note that this algorithm should only be implemented on multiblips sweeps where the neuron spike on the first and second blip. Since there is no easy way to do this, this erroneous data should not be provided to this algorithm (i.e. is should be visually checked and eliminated the preprocessor should hold back this data manually for now.)

#TODO: check to see if this is still true. Notes: The standard SDK spike detection algorithm does not work with the multiblip stimulus due to artifacts when the stimulus turns on and off. Please see the `find_multiblip_spikes` module for more information.

Input:

multi_SS: dictionary contains multiblip information such as current and stimulus

dt: float time step in seconds

Returns:

const_to_add_to_thresh_for_reset: float amplitude of the exponential fit otherwise known as `a_spike`. Note that this is without any spike cutting

decay_const: float decay constant of exponential. Note the function fit is a negative exponential which will mean this value will either have to be negated when it is used or the functions used will have to have to include the negative.

`thresh_inf`: float

```
allensdk.internal.model.glif.threshold_adaptation.exp_fit_c(t, a1, k1, const)
```

```
allensdk.internal.model.glif.threshold_adaptation.exp_force_c(t_const, a1, k1)
```

```
allensdk.internal.model.glif.threshold_adaptation.fit_avoltage_bvoltage(x,
                                                                    v_trace_list,
                                                                    El_list,
                                                                    spike_cut_length,
                                                                    all_spikeInd_list,
                                                                    th_inf,
                                                                    dt,
                                                                    a_spike,
                                                                    b_spike,
                                                                    fake=False)
                                                                    th_trace,
```

This is a version of `fit_avoltage_bvoltage_debug` that does not require the

`v_component_of_thresh_trace`, and `spike_component_of_thresh_trace` needed for debugging. A test should be run to make sure the same output comes out from this and the debug function

This function returns the squared error for the difference between the ‘known’ voltage component of the threshold obtained from the biological neuron and the voltage component of the threshold of the model obtained with the input parameters (so that the minimum can be searched for via `fmin`). The overall threshold is the sum of threshold infinity the spike component of the threshold and the voltage component of the threshold. Therefore threshold infinity and the spike component of the threshold must be subtracted from the threshold of the neuron in order to isolate the voltage component of the threshold. In the evaluation of the model the actual voltage of the neuron is used so that any errors in the other components of the model will not influence the fits here (for example, if a afterspike current was estimated incorrectly)

Notes: * The spike component of the threshold is subtracted from the

voltage which means that the voltage component of the threshold should only be added to rules.

- **b_spike was fit using a negative value in the function therefore the negative is placed in the equation.**
- **values in this function are in ‘real’ voltage as opposed to voltage** relative to resting potential.
- **current injection during the spike is not taken into account. This seems reasonable as the** ion channels are open during this time and injected current may not greatly influence the neuron.

x: numpy array `x[0]`=a_voltage input, `x[1]` is b_voltage_input, `x[2]` is th_inf

v_trace_list: list of numpy arrays voltage traces (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

El_list: list of floats reversal potential (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

spike_cut_length: int number of indices removed after initiation of a spike

all_spikeInd_list: list of numpy arrays indices of spike trains

th_inf: float threshold infinity (`v_trace`, `El`, and `th_inf` must be in the same frame of reference)

dt: float size of time step (SI units)

a_spike: float amplitude of spike component of threshold.

b_spike: float decay constant in spike component of the threshold

fake: Boolean if True makes uses the voltage value of spike step-1 because there is not a voltage value at the spike step because it is set to nan in the simulator.

```
allensdk.internal.model.glif.threshold_adaptation.fit_avoltage_bvoltage_th(x,
                                                                           v_trace_list,
                                                                           El_list,
                                                                           spike_cut_length,
                                                                           all_spikeInd_list,
                                                                           dt,
                                                                           a_spike,
                                                                           b_spike,
                                                                           fake=False)
```

This is a version of `fit_avoltage_bvoltage_th_debug` that does not require the `th_trace`, `v_component_of_thresh_trace`, and `spike_component_of_thresh_trace` needed for debugging. A test should be run to make sure the same output comes out from this and the debug function

This function returns the squared error for the difference between the ‘known’ voltage component of the threshold obtained from the biological neuron and the voltage component of the threshold of the model obtained with the input parameters (so that the minimum can be searched for via `fmin`). The overall threshold is the sum of

threshold infinity the spike component of the threshold and the voltage component of the threshold. Therefore threshold infinity and the spike component of the threshold must be subtracted from the threshold of the neuron in order to isolate the voltage component of the threshold. In the evaluation of the model the actual voltage of the neuron is used so that any errors in the other components of the model will not influence the fits here (for example, if a afterspike current was estimated incorrectly)

Notes: * The spike component of the threshold is subtracted from the

voltage which means that the voltage component of the threshold should only be added to rules.

- **b_spike was fit using a negative value in the function therefore the negative is placed in the equation.**
- **values in this function are in ‘real’ voltage as opposed to voltage** relative to resting potential.
- **current injection during the spike is not taken into account. This seems reasonable as the** ion channels are open during this time and injected current may not greatly influence the neuron.

x: numpy array x[0]=a_voltage input, x[1] is b_voltage_input, x[2] is th_inf

v_trace_list: list of numpy arrays voltage traces (v_trace, El, and th_inf must be in the same frame of reference)

El_list: list of floats reversal potential (v_trace, El, and th_inf must be in the same frame of reference)

spike_cut_length: int number of indicies removed after initiation of a spike

all_spikeInd_list: list of numpy arrays indicies of spike trains

dt: float size of time step (SI units)

a_spike: float amplitude of spike component of threshold.

b_spike: float decay constant in spike component of the threshold

fake: Boolean if True makes uses the voltage value of spike step-1 because there is not a voltage value at the spike step because it is set to nan in the simulator.

`allensdk.internal.model.glif.threshold_adaptation.get_peaks(voltage, above-Value=0)`

This function was written by Corinne Teeter and calculates the action potential peaks of a voltage equation” inputs

voltage: numpy array of voltages aboveValue: scalar voltage value over which voltage is considered a spike.

outputs: peakInd: array of indicies of peaks

Module contents

Submodules

`allensdk.internal.model.AIC` module

`allensdk.internal.model.AIC.AIC(RSS, k, n)`

Computes the Akaike Information Criterion.

RSS-residual sum of squares of the fitting errors. k - number of fitted parameters. n - number of observations.

```
allensdk.internal.model.AIC.AICc (RSS, k, n)
```

Corrected AIC. formula from Wikipedia.

```
allensdk.internal.model.AIC.BIC (RSS, k, n)
```

Bayesian information criterion or Schwartz information criterion. Formula from wikipedia.

allensdk.internal.model.GLM module

```
allensdk.internal.model.GLM.create_basis_IPSP (neye, ncos, kpeaks, ks, DTsim, t0, I_stim,  
nkt, flag_exp, npcut)
```

```
allensdk.internal.model.GLM.ff (x, c, dc)
```

```
allensdk.internal.model.GLM.invnl (x)
```

```
allensdk.internal.model.GLM.makeBasis_StimKernel (kbasprs, nkt)
```

```
allensdk.internal.model.GLM.makeBasis_StimKernel_exp (kbasprs, nkt)
```

```
allensdk.internal.model.GLM.makeFitStruct_GLM (dtsim, kbasprs, nkt, flag_exp)
```

```
allensdk.internal.model.GLM.nlin (x)
```

```
allensdk.internal.model.GLM.normalizecols (A)
```

```
allensdk.internal.model.GLM.sameconv (A, B)
```

allensdk.internal.model.data_access module

```
allensdk.internal.model.data_access.load_sweep (file_name, sweep_number, de-  
sired_dt=None, cut=0, bessel=False)
```

load a data sweep and do specified data processing. Inputs:

file_name: **string** name of .nwb data file

sweep_number: number specifying the sweep to be loaded

desired_dt: the size of the time step the data should be subsampled to

cut: indicie of which to start reporting data (i.e. cut off data before this indicie)

bessel: **dictionary** contains parameters 'N' and 'Wn' to implement standard python bessel filtering

Returns:

dictionary containing voltage: array current: array dt: time step of the returned data start_idx: the index at which the first stimulus starts (excluding the test pulse)

```
allensdk.internal.model.data_access.load_sweeps (file_name, sweep_numbers, dt=None,  
cut=0, bessel=False)
```

load sweeps and do specified data processing. Inputs:

file_name: **string** name of .nwb data file

sweep_numbers: sweep numbers to be loaded

desired_dt: the size of the time step the data should be subsampled to

cut: indicie of which to start reporting data (i.e. cut off data before this indicie)

bessel: **dictionary** contains parameters 'N' and 'Wn' to implement standard python bessel filtering

Returns:

dictionary containing voltage: list of voltage trace arrays current: list of current trace arrays dt: list of time step corresponding to each array of the returned data start_idx: list of the indices at which the first stimulus starts (excluding the test pulse) in each returned sweep

```
allensdk.internal.model.data_access.subsample_data(data, method, present_time_step,
                                                    desired_time_step)
```

Module contents

allensdk.internal.morphology package

Submodules

allensdk.internal.morphology.compartment module

```
class allensdk.internal.morphology.compartment.Compartment (node1, node2)
    Bases: object
```

allensdk.internal.morphology.morphology module

```
class allensdk.internal.morphology.morphology.Morphology (node_list=None)
    Bases: object
```

Keep track of the list of nodes in a morphology and provide a few helper methods (soma, tree information, pruning, etc).

APICAL_DENDRITE = 4

AXON = 2

BASAL_DENDRITE = 3

NODE_TYPES = [1, 2, 3, 4]

SOMA = 1

append (*self*, *nodes*)

Add additional nodes to this Morphology. Those nodes must originate from another morphology object.

Parameters

nodes: list of Morphology nodes

apply_affine (*self*, *aff*, *scale=None*)

Apply an affine transform to all nodes in this morphology. Compartment radius is adjusted as well.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

apply_affine_only_rotation (*self*, *aff*)

Apply an affine transform to all nodes in this morphology. Only the rotation element of the transform is performed (i.e., although the entire transformation and translation matrix is supplied, only the rotation element is used). The morphology is translated to the point where the soma root is at 0,0,0.

Format of the affine matrix is:

[x0 y0 z0] [tx] [x1 y1 z1] [ty] [x2 y2 z2] [tz]

where the left 3x3 the matrix defines the affine rotation and scaling, and the right column is the translation vector.

The matrix must be collapsed and stored in a list as follows:

[x0 y0, z0, x1, y1, z1, x2, y2, z2, tx, ty, tz]

Parameters

aff: 3x4 array of floats (python 2D list, or numpy 2D array) the transformation matrix

change_parent (*self*, *child*, *parent*)

Change the parent of a node. The child node is adjusted to point to the new parent, the child is taken off of the previous parent's child list, and it is added to the new parent's child list.

Parameters

child: integer or Morphology Object The ID of the child node, or the child node itself

parent: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

Nothing

children_of (*self*, *seg*)

Returns a list of the children of the specified node

Parameters

seg: integer or Morphology Object The ID of the parent node, or the parent node itself

Returns

A list of the child morphology objects. If the ID of the parent node is invalid, None is returned.

clone (*self*)

Create a clone (deep copy) of this morphology

compartment (*self*, *n*)

Returns the morphology Compartment having the specified ID.

Parameters

n: integer ID of desired compartment

Returns

A morphology object having the specified ID, or None if such a

node doesn't exist

compartment_list

convert_type (*self*, *from_type*, *to_type*)

Convert all nodes in morphology from one type to another

Parameters

from_type: enum The node type that will be eliminated and replaced. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

to_type: enum The new type that will replace it. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

delete_tree (*self*, *n*)

Delete tree, and all of its nodes, from the morphology.

Parameters

n: Integer The tree number to delete

find (*self*, *x*, *y*, *z*, *dist*, *node_type=None*)

Returns a list of Morphology Objects located within 'dist' of coordinate (x,y,z). If node_type is specified, the search will be constrained to return only nodes of that type.

Parameters

x, y, z: float The x,y,z coordinates from which to search around

dist: float The search radius

node_type: enum (optional) One of the following constants: SOMA, AXON, BASAL_DENDRITE or APICAL_DENDRITE

Returns

A list of all Morphology Objects matching the search criteria

get_dimensions (*self*)

Returns tuple of overall width, height and depth of morphology. WARNING: if locations of nodes in morphology are manipulated then this value can become incorrect. It can be reset and recalculated by programmitically setting self.dims to None.

Returns

3 real arrays: [width, height, depth], [min_x, min_y, min_z],

[max_x, max_y, max_z]

node (*self*, *n*)

Returns the morphology node having the specified ID.

Parameters

n: integer ID of desired node

Returns

A morphology node having the specified ID, or None if such a node doesn't exist

node_list

Return the node list. This is a property to ensure that the node list and node index are in sync.

node_list_by_type (*self*, *node_type*)

Return an list of all nodes having the specified node type.

Parameters

node_type: int Desired node type

Returns

A list of of Morphology Objects

num_nodes

Return the number of nodes in the morphology.

num_trees

Return the number of trees in the morphology. A tree is defined as everything following from a single root node.

parent_of (*self*, *seg*)

Returns parent of the specified node.

Parameters

seg: integer or Morphology Object The ID of the child node, or the child node itself

Returns

A morphology object, or None if no parent exists or if the specified node ID doesn't exist

save (*self*, *file_name*)

Write this morphology out to an SWC file

Parameters

file_name: string desired name of your SWC file

soma_root (*self*)

Returns root node of soma, if present

sparsify (*self*, *modulo*)

Return a new Morphology object that has a given number of non-leaf, non-root nodes removed.

Parameters

modulo: int keep 1 out of every modulo nodes.

Returns

Morphology A new morphology instance

strip_all_other_types (*self*, *node_type*, *keep_soma=True*)

Strips everything from the morphology except for the specified type. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The node type to keep in the morphology. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

keep_soma: Boolean (optional) True (default) if soma nodes should remain in the morphology, and False if the soma should also be stripped

strip_type (*self*, *node_type*)

Strips all nodes of the specified type from the morphology. Parent and child relationships are updated accordingly, creating new roots when necessary.

Parameters

node_type: enum The node type to strip from the morphology. Use one of the following constants: SOMA, AXON, BASAL_DENDRITE, or APICAL_DENDRITE

stumpify_axon (*self*, *count=10*)

Remove all axon nodes except the first ‘count’ nodes, as counted from the connected axon root.

Parameters

count: Integer The length of the axon ‘stump’, in number of nodes

to_dict (*self*)

Returns a dictionary of Node objects. These Nodes are a copy of the Morphology. Modifying them will not modify anything in the Morphology itself.

tree (*self*, *n*)

Returns a list of all Morphology nodes within the specified tree. A tree is defined as a fully connected graph of nodes. Each tree has exactly one root.

Parameters

n: integer ID of desired tree

Returns

**A list of all morphology objects in the specified tree, or None
if the tree doesn’t exist**

write (*self*, *file_name*)

allensdk.internal.morphology.morphvis module

class allensdk.internal.morphology.morphvis.**MorphologyColors**

Bases: object

set_apical_color (*self*, *r*, *g*, *b*)

set_axon_color (*self*, *r*, *g*, *b*)

set_basal_color (*self*, *r*, *g*, *b*)

set_soma_color (*self*, *r*, *g*, *b*)

allensdk.internal.morphology.morphvis.**calculate_scale** (*morph*, *pix_width*,
pix_height)

Calculates scaling factor and x,y insets required to auto-scale and center morphology into box with specified numbers of pixels

Parameters

morph: AISDK Morphology object

pix_width: int

Number of image pixels on X axis

pix_height: int

Number of image pixels on Y axis

Returns

real, real, real

First return value is the scaling factor. Second is the number of pixels needed to adjust x-coordinates so that the morphology is horizontally centered. Third is the number of pixels needed to adjust the y-coordinates so that the morphology is vertically centered.

```
allensdk.internal.morphology.morphvis.create_image(w, h, color=None, alpha=False)
allensdk.internal.morphology.morphvis.draw_density_hist(img, morph, vert_scale,
                                                         inset_left=0, inset_right=0, inset_top=0,
                                                         inset_bottom=0,
                                                         num_bins=None, colors=None)
```

Draws density histogram onto image. When no scaling is applied, and no insets are provided, the coordinates of the morphology are used directly – i.e., 100 in morphology coordinates is equal to 100 pixels.

The scale factor is multiplied to morphology coordinates before being drawn. If `scale_factor=2` then 50 in morphology coordinates is 100 pixels. Left and top insets shift the coordinate axes for drawing. E.g., if `left=10` and `top=5` then 0,0 in morphology coordinates is 10,5 in pixel space. Bottom and right insets are ignored.

If `scale_to_fit` is set then scale factor is ignored. The morphology is scaled to be the maximum size that fits in the image, taking into account insets. In a 100x100 image, if all insets=10, then the image is scaled to fit into the center 80x80 pixel area, and nothing is drawn in the inset border areas.

Axons are drawn before soma and dendrite compartments.

Parameters

img: PIL image object

morph: AISDK Morphology object

vert_scale: real

This is the amount required to multiply to a morphology y-coordinate to convert it to relative cortical depth (on [0,1]).

This is the inverse of the cortical thickness.

inset_*: real

This is the number of pixels to use as border on top/bottom/right/left. If `scale_to_fit` is false then only the top/left values are used, as the `scale_factor` will determine how large the morphology is (it can be drawn beyond insets and even beyond image boundaries)

num_bins: int

The number of bins in the histogram

colors: MorphologyColors object

This is the color scheme used to draw the morphology. If `colors=None` then default coloring is used

Returns

Histogram arrays: [hist, hist2, hist3, hist4]

**where hist is the histogram of all neurites, and hist[234] are
the histograms of SWC types 2,3,4**

```
allensdk.internal.morphology.morphvis.draw_morphology(img, morph, inset_left=0,
                                                         inset_right=0, inset_top=0, inset_bottom=0,
                                                         scale_to_fit=False,
                                                         scale_factor=1.0, colors=None)
```

Draws morphology onto image When no scaling is applied, and no insets are provided, the coordinates of the morphology are used directly – i.e., 100 in morphology coordinates is equal to 100 pixels.

The scale factor is multiplied to morphology coordinates before being drawn. If `scale_factor=2` then 50 in morphology coordinates is 100 pixels. Left and top insets shift the coordinate axes for drawing. E.g., if `left=10` and `top=5` then 0,0 in morphology coordinates is 10,5 in pixel space. Bottom and right insets are ignored.

If `scale_to_fit` is set then scale factor is ignored. The morphology is scaled to be the maximum size that fits in the image, taking into account insets. In a 100x100 image, if all insets=10, then the image is scaled to fit into the center 80x80 pixel area, and nothing is drawn in the inset border areas.

Axons are drawn before soma and dendrite compartments.

Parameters

img: PIL image object

morph: AISDK Morphology object

inset_*: real

This is the number of pixels to use as border on top/bottom/
right/left. If `scale_to_fit` is false then only the top/left
values are used, as the `scale_factor` will determine how
large the morphology is (it can be drawn beyond insets and even
beyond image boundaries)

scale_to_fit: boolean

If true then morphology is scaled to the inset area of the
image and `scale_factor` is ignored. Morphology is centered
in the image in the sense that the top/bottom and left/right
edges of the morphology are equidistant from image borders.

scale_factor: real

A scalar amount that is multiplied to morphology coordinates
before drawing

colors: MorphologyColors object

This is the color scheme used to draw the morphology. If
`colors=None` then default coloring is used

Returns

2-dimensional array, the pixel coordinates of the soma root [x,y]

allensdk.internal.morphology.node module

class allensdk.internal.morphology.node.**Node** (*n, t, x, y, z, r, pn, **kwargs*)

Bases: object

Represents node in SWC morphology file

classmethod **from_dict** (*d*)

short_string (*self*)

create string with node information in succinct, single-line form

to_dict (*self*)

Convert the node into a serializable dictionary

allensdk.internal.morphology.node.**euclidean_distance** (*node1, node2*)

allensdk.internal.morphology.node.**midpoint** (*node1, node2*)

allensdk.internal.morphology.validate_swc module

class allensdk.internal.morphology.validate_swc.**TestNode** (*n, t, x, y, z, r, pn*)

Bases: object

allensdk.internal.morphology.validate_swc.**main** ()

allensdk.internal.morphology.validate_swc.**resave_swc** (*orig_swc, new_file*)

Reads SWC file into AllenSDK Morphology object and resaves it. This can fix some problems in an SWC file that may disrupt other software tools reading the file (e.g., NEURON)

Parameters

orig_swc: string Name of SWC file to read

new_file: string Name of output SWC file

allensdk.internal.morphology.validate_swc.**validate_swc** (*swc_file*)

Tests SWC files for compatibility with AllenSDK

To be compatible with NEURON, SWC files must have the following properties:

- 1) a single root node with parent ID '-1'
- 2) sequentially increasing ID numbers
- 3) immediate children of the soma cannot branch

To be compatible with feature analysis, SWC files can only have node types in the range 1-4:

1 = soma 2 = axon 3 = [basal] dendrite 4 = apical dendrite

Module contents

allensdk.internal.mouse_connectivity package

Subpackages

allensdk.internal.mouse_connectivity.interval_unionize package

Submodules

allensdk.internal.mouse_connectivity.interval_unionize.cav_unionize module

allensdk.internal.mouse_connectivity.interval_unionize.cav_unionizer module

allensdk.internal.mouse_connectivity.interval_unionize.data_utilities module

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_cav_density(cav_data)`

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_injection_data(injection_data)`
 Read nrrd files containing injection signal data

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_projection_data(projection_data)`
 Read nrrd files containing global signal data

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.get_sum_pixel_intensities(sum_pixel_intensities)`
`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.load_annotation(annotation_data)`
 Read data files segmenting the reference space into regions of valid and invalid data, then further among brain structures

`allensdk.internal.mouse_connectivity.interval_unionize.data_utilities.read(path)`

allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer module

class `allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer.IntervalUnionizer`
 Bases: object

direct_unionize (*self*, *data_arrays*, *pre_sorted=False*, ***kwargs*)
 Obtain unionize records from directly annotated regions.

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened arrays.

sorted [bool, optional] If False, data arrays will be sorted.

extract_data (*self*, *data_arrays*, *low*, *high*, ***kwargs*)
 Given flattened data arrays and a specified interval, generate summary data

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened, sorted arrays.

low [int] Index at which interval of interest begins. Inclusive.

high [int] Index at which interval of interest ends. Exclusive.

postprocess_unionizes (*self*, *raw_unionizes*, ***kwargs*)

Carry out additional calculations/formatting derivative of core unionization.

Parameters

raw_unionizes [list of unionizes] Each entry is a unionize record.

classmethod propagate_record (*child_record*, *ancestor_record*, *copy_all=False*)

Updates one unionize corresponding to a rootward structure with information from a unionize corresponding to a leafward structure

Parameters

child_record [unionize] Data will be drawn from this record

ancestor_record [unionize] This record will be updated

classmethod propagate_to_bilateral (*lateral_unionizes*)

classmethod propagate_unionizes (*direct_unionizes*, *ancestor_id_map*)

Structures are arranged in a tree, whose leafward-oriented edges indicate physical containment. This method updates rootward unionize records with information from leafward ones.

Parameters

direct_unionizes [list of unionizes] Each entry is a unionize record produced from a collection of directly labeled voxels in the segmentation volume.

ancestor_id_map [dict] Keys are structure ids. Values are ids of all structures rootward in

the tree, including the key node

Returns

output_unionizes [list of unionizes] Contains completed unionize records at all depths in the structure tree

classmethod record_cb ()

setup_interval_map (*self*, *annotation*)

Build a map from structure ids to intervals in the sorted flattened reference space.

Parameters

annotation [np.ndarray] Segmentation label array.

sort_data_arrays (*self*, *data_arrays*)

Apply the precomputed sort to flattened data arrays

Parameters

data_arrays [dict] Keys identify types of data volume. Values are flattened, unsorted arrays.

Returns

dict : As input, but values are sorted

allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_cav module

allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic module

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.get`

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.get`

`allensdk.internal.mouse_connectivity.interval_unionize.run_tissuecyte_unionize_classic.run`

allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record module

class `allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.Ti`

Bases: `allensdk.internal.mouse_connectivity.interval_unionize.unionize_record.Unionize`

direct_sum_projection_pixels

max_voxel_density

max_voxel_index

output (*self*, *output_spacing_iso*, *volume_scale*, *target_shape*, *sort*)
Generate derived data for this unionize

Parameters

output_spacing_iso [numeric] Isometric spacing of reference space in microns

volume_scale [numeric] Scale factor mapping pixels to microns³

target_shape [array-like of numeric] Shape of reference space

projection_density

projection_energy

projection_intensity

propagate (*self*, *ancestor*, *copy_all=False*)
Update a rootward unionize with data from this unionize record

Parameters

ancestor [TissuecyteBaseUnionize] will be updated

Returns

ancestor [TissuecyteBaseUnionize]

set_max_voxel (*self*, *density_array*, *low*)
Find the voxel of greatest density in this unionizes spatial domain

Parameters

density_array [ndarray] Float values are densities per voxel

low [int] index in full flattened, sorted array of starting voxel

sum_pixel_intensity

sum_pixels

sum_projection_pixel_intensity

sum_projection_pixels**class** allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.**Ti**

Bases: *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteBaseUnionize*

calculate (*self, low, high, data_arrays*)**class** allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.**Ti**

Bases: *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionize_record.TissuecyteBaseUnionize*

calculate (*self, low, high, data_arrays, ij_record*)**allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer module****class** allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer.**Tissuecyt**

Bases: *allensdk.internal.mouse_connectivity.interval_unionize.interval_unionizer.IntervalUnionizer*

A specialization of the IntervalUnionizer set up for unionizing Tissuecyte-derived projection data.

extract_data (*self, data_arrays, low, high*)

As parent

postprocess_unionizes (*self, raw_unionizes, image_series_id, output_spacing_iso, volume_scale, target_shape, sort*)

As parent

classmethod propagate_record (*child_record, ancestor_record, copy_all=False*)

As parent

classmethod record_cb ()**allensdk.internal.mouse_connectivity.interval_unionize.unionize_record module****class** allensdk.internal.mouse_connectivity.interval_unionize.unionize_record.**Unionize** (*args, **kwargs)

Bases: object

Abstract base class for unionize records.

calculate (*self, *args, **kwargs*)**output** (*self, *args, **kwargs*)**propagate** (*self, ancestor, copy_all, *args, **kwargs*)**slice_arrays** (*self, low, high, data_arrays*)

Extract a slice from several aligned arrays

Parameters

low [int] start of slice, inclusive

high [int] end of slice, exclusive

data_arrays [dict] keys are varieties of data. values are sorted, flattened data arrays

Module contents

allensdk.internal.mouse_connectivity.projection_thumbnail package

Submodules

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip module

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**apply_c**

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**blend_v**

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**do_blur**

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**handle**

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**max_cb**

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**run** (vol
imi
ima
ro-
ta-
tion
col-
orn

allensdk.internal.mouse_connectivity.projection_thumbnail.generate_projection_strip.**simple**

allensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet module

```
class allensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet.ImageSheet
    Bases: object
    append (self, new_cell)
    apply (self, fn, *args, **kwargs)
```

```
static build_from_image (image, n, axis)  
copy (self)  
get_output (self, axis)
```

allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions module

```
allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.convert_axis  
allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.max_project
```

```
allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions.template_pro
```

allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities module

```
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.blend (ima  
                                                                                               wei
```

Parameters

image_stack :: list of np.ndarray The images to be blended. Shapes cannot differ

weight_stack :: list of np.ndarray The weight of each image at each pixel. Will be normal-
ized.

```
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.convert_c
```

Generates a matplotlib continuous colormap on [0, 1] from a discrete colormap at N evenly spaced points.

Parameters

data [list of list] Sublists are [r, g, b].

Returns

matplotlib.colors.LinearSegmentedColormap Gamma is 1. Output space is 3 X [0, 1]

```
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.minmax_n
```

```
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.normalize
```

```
allensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities.sitk_safe
```

allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector module

```

class allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector.VolumePro
    Bases: object

    build_rotation_transform (self, from_axis, to_axis, angle)

    extract (self, cb, volume=None)

    classmethod fixed_factory (volume, size)

    rotate (self, from_axis, to_axis, angle)

    rotate_and_extract (self, from_axes, to_axes, angles, cb)

    classmethod safe_factory (volume)

```

allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities module

```

allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_center
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_diagonal
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_image_p
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_get_size_p
allensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities.sitk_paste_into

```

Module contents**allensdk.internal.mouse_connectivity.tissuecyte_stitching package****Submodules****allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher module**

```

class allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.Stitcher (image_dimensions,
                                         tiles,
                                         average_image_tiles,
                                         channels)

    Bases: object

    run (self, cb=<built-in function array>)

    stitch (self, slice_image, stitched_indicator, tile, cb=<built-in function array>)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.blend_component_from_pos

```

Obtains a normalized component of the blend, which describes depth of overlap along a specified axis in a specified direction

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_blend(indicator_region,
                                                                              stup,
                                                                              cb=<built-in
                                                                              function
                                                                              ar-
                                                                              ray>)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_blend_component(indicator_region,
                                                                                          lg,
                                                                                          axis,
                                                                                          mesh)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_indicator_bound_point(indicator_region,
                                                                                              lg,
                                                                                              axis,
                                                                                              mesh)
```

Finds the index of first change in a binary mask along a specified axis in a specified direction

```
allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.get_overall_blend(indicator_region,
                                                                                          lg,
                                                                                          axis,
                                                                                          mesh)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.initialize_image(dimension,
                                                                                      nchan-
                                                                                      nels,
                                                                                      dtype,
                                                                                      or-
                                                                                      der='C')

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.initialize_images(dimension,
                                                                                      nchan-
                                                                                      nels)

allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.make_blended_tile(blend,
                                                                                      tile,
                                                                                      cur-
                                                                                      rent_region)
```

allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile module

```
class allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile(index,
                                                                              im-
                                                                              age,
                                                                              is_missing,
                                                                              bounds,
                                                                              chan-
                                                                              nel,
                                                                              size,
                                                                              mar-
                                                                              gins,
                                                                              *args,
                                                                              **kwargs)
```

Bases: object

```
apply_average_tile(self, average_tile)

apply_average_tile_to_self(self, average_tile)

average_tile_is_untrimmed(self, average_tile)
```

```

get_image_region (self)
get_missing_path (self)
initialize_image (self)
trim (self, image)
trim_self (self)

```

Module contents

Module contents

allensdk.internal.pipeline_modules package

Subpackages

allensdk.internal.pipeline_modules.gbm package

Submodules

allensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_run_records module

```

allensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_run_records.main (analysis_records_  

                                                                 db_host,  

                                                                 db_port,  

                                                                 db_name,  

                                                                 db_user,  

                                                                 db_passwd)

```

allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap module

```

allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_gene_fpkm_table (analysis_ru  

    Creates a a matrix ("rows x columns = genes x samples") of fpkm gene expression values for each particular  

    (gene, sample) pair. Rows are sorted by entrez_id and columns are by rna_well_id
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_genes_for_transcripts (an  

    Creates a list that contains the associated gene for each transcript sorted alphabetically
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_sample_metadata (sample_me  

    Creates a table of sample metadata sorted by rna_well_id
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_transcript_fpkm_table (an  

    Creates a a matrix ("rows x columns = transcripts x samples") of fpkm gene expression values for each particular  

    (transcript, sample) pair. Rows are sorted by transcript id and columns are by rna_well_id
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.create_transcripts_for_genes (an  

    Creates a list that contains the associated transcript for each gene sorted by entrez_id
allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap.main ()

```

allensdk.internal.pipeline_modules.gbm.generate_gbm_sample_metadata module

```
allensdk.internal.pipeline_modules.gbm.generate_gbm_sample_metadata.main(sample_metadata_json_lo-  
db_host,  
db_port,  
db_name,  
db_user,  
db_passwd)
```

Module contents

Submodules

allensdk.internal.pipeline_modules.run_annotated_region_metrics module

Run annotated region metrics calculations

```
allensdk.internal.pipeline_modules.run_annotated_region_metrics.debug(region_id,  
stor-  
age_directory='./',  
lo-  
cal=True,  
sdk_path='/data/informatics/C  
script_path='/data/informatics  
lims_host='lims2')  
  
allensdk.internal.pipeline_modules.run_annotated_region_metrics.load_arrays(h5_file)  
allensdk.internal.pipeline_modules.run_annotated_region_metrics.main()
```

allensdk.internal.pipeline_modules.run_demixing module

```
allensdk.internal.pipeline_modules.run_demixing.assert_exists(file_name)  
allensdk.internal.pipeline_modules.run_demixing.debug(experiment_id, local=False)  
allensdk.internal.pipeline_modules.run_demixing.get_path(obj, key, check_exists)  
allensdk.internal.pipeline_modules.run_demixing.main()  
allensdk.internal.pipeline_modules.run_demixing.parse_input(data, ex-  
clude_labels)
```

allensdk.internal.pipeline_modules.run_dff_computation module

```
allensdk.internal.pipeline_modules.run_dff_computation.main()  
allensdk.internal.pipeline_modules.run_dff_computation.parse_input(data)
```

allensdk.internal.pipeline_modules.run_eye_tracking module

allensdk.internal.pipeline_modules.run_neuropil_correction module

```
allensdk.internal.pipeline_modules.run_neuropil_correction.adjust_r_for_negativity(r,
                                                                                    F_C,
                                                                                    F_M,
                                                                                    F_N)

allensdk.internal.pipeline_modules.run_neuropil_correction.debug(experiment_id,
                                                                    local=False)

allensdk.internal.pipeline_modules.run_neuropil_correction.debug_plot(file_name,
                                                                        roi_trace,
                                                                        neu-
                                                                        ropil_trace,
                                                                        cor-
                                                                        rected_trace,
                                                                        r,
                                                                        r_vals=None,
                                                                        err_vals=None)

allensdk.internal.pipeline_modules.run_neuropil_correction.main()
```

allensdk.internal.pipeline_modules.run_observatory_analysis module

```
allensdk.internal.pipeline_modules.run_observatory_analysis.debug(experiment_ids,
                                                                    lo-
                                                                    cal=False,
                                                                    OUT-
                                                                    PUT_DIR='/data/informatics/CAM/',
                                                                    SDK_PATH='/data/informatics/CAM',
                                                                    wall-
                                                                    time='10:00:00',
                                                                    python='/shared/utils.x86_64/python
                                                                    2.7/bin/python',
                                                                    queue='braintv')

allensdk.internal.pipeline_modules.run_observatory_analysis.get_experiment_nwb_file(experiment)

allensdk.internal.pipeline_modules.run_observatory_analysis.get_experiment_session(experiment)

allensdk.internal.pipeline_modules.run_observatory_analysis.main()
```

allensdk.internal.pipeline_modules.run_observatory_container_thumbnails module

allensdk.internal.pipeline_modules.run_observatory_thumbnails module

```
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_cell_plots(cell_specimen_id,
                                         pre-
                                         fix,
                                         as-
                                         spect,
                                         con-
                                         figs,
                                         out-
                                         put_dir,
                                         axes=None,
                                         trans-
                                         par-
                                         ent=False)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_correlation_plots(data_source,
                                              anal-
                                              y-
                                              sis_file,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_drifting_gratings(dga,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_experiment_thumbnails(name,
                                              a,
                                              y,
                                              s,
                                              a,
                                              h,
                                              t,
                                              t)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_eye_tracking_plots(data_source,
                                              con-
                                              figs,
                                              out-
                                              put_dir)

allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_locally_sparse_noise(location,
                                              co-
                                              figs,
                                              ou-
                                              pu-
                                              on)
```

```
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_natural_movie(nma,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir,  
                                          name)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_natural_scenes(nsa,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_plots(prefix,  
                                     as-  
                                     pect,  
                                     con-  
                                     figs,  
                                     out-  
                                     put_dir,  
                                     axes=None,  
                                     trans-  
                                     par-  
                                     ent=False)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_receptive_field(lsna,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_speed_tuning(analysis,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_static_gratings(sga,  
                                          con-  
                                          figs,  
                                          out-  
                                          put_dir)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.build_type(nwb_file,  
                                   data_file,  
                                   con-  
                                   figs,  
                                   out-  
                                   put_dir,  
                                   type_name)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.debug(experiment_id,  
                             plots=None,  
                             lo-  
                             cal=False)  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_analysis_file  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_files(experiment  
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_experiment_nwb_file(experim
```

```
allensdk.internal.pipeline_modules.run_observatory_thumbnails.get_input_data(experiment_id)
allensdk.internal.pipeline_modules.run_observatory_thumbnails.lsna_check_hvas(data_set,
                                                                               data_file)
allensdk.internal.pipeline_modules.run_observatory_thumbnails.main()
allensdk.internal.pipeline_modules.run_observatory_thumbnails.parse_input(data)
```

allensdk.internal.pipeline_modules.run_ophys_eye_calibration module

```
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.debug(experiment_id,
                                                                    lo-
                                                                    cal=False)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.get_wkf(wkf_type,
                                                                    experi-
                                                                    ment_id)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.main()
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.parse_input_data(data)
allensdk.internal.pipeline_modules.run_ophys_eye_calibration.write_output(filename,
                                                                    po-
                                                                    si-
                                                                    tion_degrees,
                                                                    po-
                                                                    si-
                                                                    tion_cm,
                                                                    ar-
                                                                    eas)
```

allensdk.internal.pipeline_modules.run_ophys_session_decomposition module

```
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.convert_frame(conversion_de
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.create_fake_metadata(exp
raw
cha
nel
wid
heig
item
size
n_p

allensdk.internal.pipeline_modules.run_ophys_session_decomposition.debug(experiment_id,
                                                                    lo-
                                                                    cal=False,
                                                                    raw_path=None)
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.main()
allensdk.internal.pipeline_modules.run_ophys_session_decomposition.parse_input(data)
    Load all input data from the input json.
```

allensdk.internal.pipeline_modules.run_ophys_time_sync module

class allensdk.internal.pipeline_modules.run_ophys_time_sync.**TimeSyncOutputs**

Bases: tuple

Schema for synchronization outputs

behavior_alignment

Alias for field number 12

behavior_delta

Alias for field number 5

behavior_times

Alias for field number 9

experiment_id

Alias for field number 0

eye_alignment

Alias for field number 11

eye_delta

Alias for field number 4

eye_times

Alias for field number 8

ophys_delta

Alias for field number 2

ophys_times

Alias for field number 6

stimulus_alignment

Alias for field number 10

stimulus_delay

Alias for field number 1

stimulus_delta

Alias for field number 3

stimulus_times

Alias for field number 7

class allensdk.internal.pipeline_modules.run_ophys_time_sync.**TimeSyncWriter**(*output_h5_path:*

str,
output_json_path:
Optional[str]
=
None)

Bases: object

validate_paths(*self*)

Determines whether we can actually write to the specified paths, allowing for creation of intermediate directories. It is a good idea to run this before doing any heavy calculations!

write(*self*, *outputs:* allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs)

Convenience for writing both an output h5 and (if applicable) an output json.

Parameters

outputs [the data to be written]

write_output_h5 (*self*, *outputs*)

Write (mainly) heaviweight data to an h5 file.

Parameters

outputs [the data to be written]

write_output_json (*self*, *outputs*)

Write lightweight data to a json

Parameters

outputs [the data to be written]

`allensdk.internal.pipeline_modules.run_ophys_time_sync.check_stimulus_delay` (*obt_delay*:
float,
min_delay:
float,
max_delay:
float)

Raise an exception if the monitor delay is not within specified bounds

Parameters

obt_delay [obtained monitor delay (s)]

min_delay [lower threshold (s)]

max_delay [upper threshold (s)]

`allensdk.internal.pipeline_modules.run_ophys_time_sync.main()`

`allensdk.internal.pipeline_modules.run_ophys_time_sync.run_ophys_time_sync` (*aligner*:
al-
lensdk.internal.brain_d
ex-
per-
i-
ment_id:
int,
min_stimulus_delay:
float,
max_stimulus_delay:
float)
→
al-
lensdk.internal.pipeline

Carry out synchronization of timestamps across the data streams of an ophys experiment.

Parameters

aligner [drives alignment. See OphysTimeAligner for details of the] attributes and properties that must be implemented.

experiment_id [unique identifier for the experiment being aligned]

min_stimulus_delay [reject alignment run (raise a ValueError) if the] calculated monitor delay is below this value (s).

max_stimulus_delay [reject alignment run (raise a ValueError) if the] calculated monitor delay is above this value (s).

Returns

A **TimeSyncOutputs** (see definition for more information) of output parameters and arrays of aligned timestamps.

`allensdk.internal.pipeline_modules.run_roi_filter` module

`allensdk.internal.pipeline_modules.run_tissuecyte_projection_thumbnail_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_stitching_classic` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_cav_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_counts_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_from_json` module

`allensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic_from_json.main()`

Module contents

Module contents

6.1.7 allensdk.model package

Subpackages

`allensdk.model.biophys_sim` package

Subpackages

`allensdk.model.biophys_sim.neuron` package

Submodules

`allensdk.model.biophys_sim.neuron.hoc_utils` module

class `allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils` (*description*)

Bases: `object`

A helper class for containing references to NEUORN.

Attributes

h [object] The NEURON hoc object.

nrn [object] The NEURON python object.

neuron [module] The NEURON module.

```
h = None
initialize_hoc(self)
    Basic setup for NEURON.
neuron = None
nrn = None
```

Module contents

allensdk.model.biophys_sim.scripts package

Module contents

Submodules

allensdk.model.biophys_sim.bps_command module

```
allensdk.model.biophys_sim.bps_command.choose_bps_command(command='bps_simple',
                                                            conf_file=None)
allensdk.model.biophys_sim.bps_command.run_module(description, module_name, func-
                                                    tion_name)
```

allensdk.model.biophys_sim.config module

```
class allensdk.model.biophys_sim.config.Config
```

Bases: *allensdk.config.app.application_config.ApplicationConfig*

```
load(self, config_path, disable_existing_logs=False)
```

Parse the application configuration then immediately load the model configuration files.

Parameters

disable_existing_logs [boolean, optional] If false (default) leave existing logs after configuration.

```
read_model_description(self)
```

parse the model_file field of the application configuration and read the files.

The model_file field of the application configuration is first split at commas, since it may list more than one file.

The files may be uris of the form `file:filename?section=name`, in which case a bare configuration object is read from filename into the configuration section with key 'name'.

A simple filename without a section option is treated as a standard multi-section configuration file.

Returns

description [Description] Configuration object.

Module contents

allensdk.model.biophysical package

Submodules

allensdk.model.biophysical.run_simulate module

class allensdk.model.biophysical.run_simulate.**RunSimulate** (*input_json*, *output_json*)

Bases: object

load_manifest (*self*)

nrnivmodl (*self*)

simulate (*self*)

allensdk.model.biophysical.run_simulate.**main** (*command*, *lims_strategy_json*, *lims_response_json*)

Entry point for module. :param command: select behavior, nrnivmodl or simulate :type command: string
:param lims_strategy_json: path to json file output from lims. :type lims_strategy_json: string :param
lims_response_json: path to json file returned to lims. :type lims_response_json: string

allensdk.model.biophysical.runner module

allensdk.model.biophysical.runner.**load_description** (*args_dict*)

Read configurations.

Parameters

args_dict [dict] Parsed arguments dictionary with following keys.

manifest_file [string] .json file with containing the experiment configuration

axon_type [string] Axon handling for the all-active models

Returns

Config Object with all information needed to run the experiment.

allensdk.model.biophysical.runner.**prepare_nwb_output** (*nwb_stimulus_path*, *nwb_result_path*)

Copy the stimulus file, zero out the recorded voltages and spike times.

Parameters

nwb_stimulus_path [string] NWB file name

nwb_result_path [string] NWB file name

allensdk.model.biophysical.runner.**run** (*args*, *sweeps=None*, *procs=6*)

Main function for simulating sweeps in a biophysical experiment.

Parameters

args [dict] Parsed arguments to run the experiment.

procs [int] number of sweeps to simulate simultaneously.

sweeps [list] list of experiment sweep numbers to simulate. If None, simulate all sweeps.

`allensdk.model.biophysical.runner.run_sync(description, sweeps=None)`

Single-process main function for simulating sweeps in a biophysical experiment.

Parameters

description [Config] All information needed to run the experiment.

sweeps [list] list of experiment sweep numbers to simulate. If None, simulate all sweeps.

`allensdk.model.biophysical.runner.save_nwb(output_path, v, sweep, sweeps_by_type)`

Save a single voltage output result into an existing sweep in a NWB file. This is intended to overwrite a recorded trace with a simulated voltage.

Parameters

output_path [string] file name of a pre-existing NWB file.

v [numpy array] voltage

sweep [integer] which entry to overwrite in the file.

allensdk.model.biophysical.utils module

class `allensdk.model.biophysical.utils.AllActiveUtils(description, axon_type)`

Bases: `allensdk.model.biophysical.utils.Utils`

generate_morphology (*self*, *morph_filename*)

Load a neuroLucida or swc-format cell morphology file.

Parameters

morph_filename [string] Path to morphology.

load_cell_parameters (*self*)

Configure a neuron after the cell morphology has been loaded.

class `allensdk.model.biophysical.utils.Utils(description)`

Bases: `allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils`

A helper class for NEURON functionality needed for biophysical simulations.

Attributes

h [object] The NEURON hoc object.

nrn [object] The NEURON python object.

neuron [module] The NEURON module.

generate_morphology (*self*, *morph_filename*)

Load a swc-format cell morphology file.

Parameters

morph_filename [string] Path to swc.

get_recorded_data (*self*, *vec*)

Extract recorded voltages and timestamps given the recorded Vector instance. If `self.stimulus_sampling_rate` is smaller than `self.simulation_sampling_rate`, resample to `self.stimulus_sampling_rate`.

Parameters

vec [neuron.Vector] constructed by `self.record_values`

Returns

dict with two keys: 'v' = numpy.ndarray with voltages, 't' = numpy.ndarray with timestamps

load_cell_parameters (*self*)

Configure a neuron after the cell morphology has been loaded.

static nearest_neuron_sampling_rate (*hz*, *target_hz=40000*)

read_stimulus (*self*, *stimulus_path*, *sweep=0*)

Load current values for a specific experiment sweep and setup simulation and stimulus sampling rates.

NOTE: NEURON only allows simulation timestamps of multiples of 40KHz. To avoid aliasing, we set the simulation sampling rate to the least common multiple of the stimulus sampling rate and 40KHz.

Parameters

stimulus_path [string] NWB file name

sweep [integer, optional] sweep index

record_values (*self*)

Set up output voltage recording.

setup_iclamp (*self*, *stimulus_path*, *sweep=0*)

Assign a current waveform as input stimulus.

Parameters

stimulus_path [string] NWB file name

update_default_cell_hoc (*self*, *description*, *default_cell_hoc='cell.hoc'*)

replace the default 'cell.hoc' path in the manifest with 'cell.hoc' packaged within AllenSDK if it does not exist

`allensdk.model.biophysical.utils.create_utils` (*description*, *model_type=None*)

Factory method to create a Utils subclass.

Parameters

description [Config instance] used to initialize Utils subclass

model_type [string] Must be one of [PERISOMATIC_TYPE, ALL_ACTIVE_TYPE]. If none, defaults to PERISOMATIC_TYPE

Returns

Utils instance

Module contents**allensdk.model.glif package****Submodules****allensdk.model.glif.glif_neuron module**

exception `allensdk.model.glif.glif_neuron.GlifBadResetException` (*message*, *dv*)

Bases: `Exception`

Exception raised when voltage is still above threshold after a reset rule is applied.

```
class allensdk.model.glif.glif_neuron.GlifNeuron (El, dt, asc_tau_array, R_input, C,  
                                                asc_amp_array, spike_cut_length,  
                                                th_inf, th_adapt, coeffs, AScurrent_dynamics_method,  
                                                voltage_dynamics_method,  
                                                threshold_dynamics_method,  
                                                AScurrent_reset_method,  
                                                voltage_reset_method,  
                                                threshold_reset_method,  
                                                init_voltage,  
                                                init_threshold,  
                                                init_AScurrents,  
                                                **kwargs)
```

Bases: object

Implements the current-based Mihalas Neiber GLIF neuron. Simulations model the voltage, threshold, and afterspike currents of a neuron given an input stimulus. A set of modular dynamics rules are applied until voltage crosses threshold, at which point a set of modular reset rules are applied. See `glif_neuron_methods.py` for a list of what options there are for voltage, threshold, and afterspike current dynamics and reset rules.

Parameters

El [float]
resting potential

dt [float] duration between time steps

asc_tau_array: np.ndarray TODO

R_input [float] input resistance

C [float] capacitance

asc_amp_array [np.ndarray] afterspike current vector. one element per element of `asc_tau_array`.

spike_cut_length [int] how many time steps to replace with NaNs when a spike occurs.

th_inf [float] instantaneous threshold

coeffs [dict] dictionary coefficients premultiplied to neuron properties during simulation. used for optimization.

AScurrent_dynamics_method [dict] dictionary containing the ‘name’ of the afterspike current dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

voltage_dynamics_method [dict] dictionary containing the ‘name’ of the voltage dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

threshold_dynamics_method [dict] dictionary containing the ‘name’ of the threshold dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

AScurrent_reset_method [dict] dictionary containing the ‘name’ of the afterspike current dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

voltage_reset_method [dict] dictionary containing the ‘name’ of the voltage dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

threshold_reset_method [dict] dictionary containing the ‘name’ of the threshold dynamics method to use and a ‘params’ dictionary parameters to pass to that function.

init_voltage [float] initial voltage value

init_threshold [float] initial spike threshold value

init_AScurrents [np.ndarray] initial afterspike current vector. one element per element of asc_tau_array.

TYPE = 'GLIF'

append_threshold_components (*self, spike, voltage*)

static configure_library_method (*method_type, params*)

Create a GlifNeuronMethod instance out of a library of functions organized by type name. This refers to the METHOD_LIBRARY in glif_neuron_methods.py, which lays out the available functions that can be used for dynamics and reset rules.

Parameters

method_type [string] the name of a function category (e.g. 'AScurrent_dynamics_method' for the afterspike current dynamics methods)

params [dict] a dictionary with two members. 'name': the string name of function you want, and 'params': parameters you want to pass to that function

Returns

GlifNeuronMethod a GlifNeuronMethod instance

static configure_method (*method_name, method, method_params*)

Create a GlifNeuronMethod instance given a name, a function, and function parameters. This is just a shortcut to the GlifNeuronMethod constructor.

Parameters

method_name [string] name for referring to this method later

method [function] a python function

method_parameters [dict] function arguments whose values should be fixed

Returns

GlifNeuronMethod a GlifNeuronMethod instance

dynamics (*self, voltage_t0, threshold_t0, AScurrents_t0, inj, time_step, spike_time_steps*)

Update the voltage, threshold, and afterspike currents of the neuron for a single time step.

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

inj [float] the current value of the current injection into the neuron

time_step [int] the current time step of the neuron simulation

spike_time_steps [list] a list of all of the time steps of spikes in the neuron

Returns

tuple voltage_t1 (voltage at next time step), threshold_t1 (threshold at next time step), AScurrents_t1 (afterspike currents at next time step)

classmethod from_dict (*d*)

reset (*self*, *voltage_t0*, *threshold_t0*, *AScurrents_t0*)

Apply reset rules to the neuron's voltage, threshold, and afterspike currents assuming a spike has occurred (voltage is above threshold).

Parameters

voltage_t0 [float] the current voltage of the neuron

threshold_t0 [float] the current spike threshold level of the neuron

AScurrents_t0 [np.ndarray] the current state of the afterspike currents in the neuron

Returns

tuple *voltage_t1* (voltage at next time step), *threshold_t1* (threshold at next time step), *AScurrents_t1* (afterspike currents at next time step)

run (*self*, *stim*)

Run neuron simulation over a given stimulus. This steps through the stimulus applying dynamics equations. After each step it checks if voltage is above threshold. If so, *self.spike_cut_length* NaNs are inserted into the output voltages, reset rules are applied to the voltage, threshold, and afterspike currents, and the simulation resumes.

Parameters

stim [np.ndarray] vector of scalar current values

Returns

dict

a dictionary containing: 'voltage': simulated voltage values, 'threshold': threshold values during the simulation, 'AScurrents': afterspike current values during the simulation, 'grid_spike_times': spike times (in units of *self.dt*) aligned to simulation time steps, 'interpolated_spike_times': spike times (in units of *self.dt*) linearly interpolated between time steps, 'spike_time_steps': the indices of grid spike times, 'interpolated_spike_voltage': voltage of the simulation at interpolated spike times, 'interpolated_spike_threshold': threshold of the simulation at interpolated spike times

tau_m

to_dict (*self*)

Convert the neuron to a serializable dictionary.

`allensdk.model.glif.glif_neuron.interpolate_spike_time` (*dt*, *time_step*, *threshold_t0*,
threshold_t1, *voltage_t0*,
voltage_t1)

Given two voltage and threshold values, the *dt* between them and the initial time step, interpolate a spike time within the *dt* interval by intersecting the two lines.

`allensdk.model.glif.glif_neuron.interpolate_spike_value` (*dt*, *interpo-*
lated_spike_time_offset,
v0, *v1*)

Take a value at two adjacent time steps and linearly interpolate what the value would be at an offset between the two time steps.

`allensdk.model.glif.glif_neuron.line_crossing_x` (*dx*, *a0*, *a1*, *b0*, *b1*)

Find the *x* value of the intersection of two lines.

`allensdk.model.glif.glif_neuron.line_crossing_y` (*dx*, *a0*, *a1*, *b0*, *b1*)

Find the *y* value of the intersection of two lines.

allensdk.model.glif.glif_neuron_methods module

The methods in this module are used for configuring dynamics and reset rules for the GlifNeuron. For more details on how to use these methods, see [Generalized LIF Models](#).

```
class allensdk.model.glif.glif_neuron_methods.GlifNeuronMethod(method_name,
                                                                method,
                                                                method_params)
```

Bases: object

A simple class to keep track of the name and parameters associated with a neuron method. This class is initialized with a name, function, and parameters to pass to the function. The function then has those passed parameters fixed to a partial function using `functools.partial`. This class then mimics a function itself using the `__call__` convention. Parameters that are not fixed in this way are assumed to be passed into the method when it is called. If the passed parameters contain an argument that is not part of the function signature, an exception will be raised.

Parameters

method_name [string] A shorthand name that will be used to reference this method in the *GlifNeuron*.

method [function] A python function to be called when this instance is called.

method_params [dict] A dictionary mapping function arguments to values for values that should be fixed.

```
modify_parameter(self, param, operator)
```

Modify a function parameter needs to be modified after initialization.

Parameters

param [string] the name of the parameter to modify

operator [callable] a function or lambda that returns the desired modified value

Returns

type the new value of the variable that was just modified.

```
to_dict(self)
```

```
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_exp(neuron, AS-
                                                                currents_t0,
                                                                time_step,
                                                                spike_time_steps)
```

Exponential afterspike current dynamics method takes a current at t0 and returns the current at a time step later.

```
allensdk.model.glif.glif_neuron_methods.dynamics_AScurrent_none(neuron, AS-
                                                                currents_t0,
                                                                time_step,
                                                                spike_time_steps)
```

This method always returns zeros for the afterspike currents, regardless of input.

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_inf(neuron, thresh-
                                                                old_t0,    volt-
                                                                age_t0,    AS-
                                                                currents_t0,
                                                                inj)
```

Set threshold to the neuron's instantaneous threshold.

Parameters

neuron [class]
threshold_t0 [not used here]
voltage_t0 [not used here]
AScurrents_t0 [not used here]
inj [not used here]
AScurrents_t0 [not used here]
inj [not used here]

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_spike_component(neuron,  
                                                                           thresh-  
                                                                           old_t0,  
                                                                           volt-  
                                                                           age_t0,  
                                                                           AS-  
                                                                           cur-  
                                                                           rents_t0,  
                                                                           inj,  
                                                                           a_spike,  
                                                                           b_spike,  
                                                                           a_voltage,  
                                                                           b_voltage)
```

Analytical solution for spike component of threshold. The threshold will adapt via a component initiated by a spike which decays as an exponential. The component is in reference to threshold infinity and are recorded in the neuron's threshold components. The voltage component of the threshold is set to zero in the threshold components because it is zero here. The third component refers to `th_inf` which is added separately as opposed to being included in the voltage component of the threshold as is done in equation 2.1 of Mihalas and Nieber 2009. Threshold infinity is removed for simple optimization.

Parameters

neuron [class]
threshold_t0 [float] threshold input to function
voltage_t0 [float] voltage input to function
AScurrents_t0 [vector] values of after spike currents
inj [float] current injected into the neuron

```
allensdk.model.glif.glif_neuron_methods.dynamics_threshold_three_components_exact(neuron,  
                                                                           thresh-  
                                                                           old_t0,  
                                                                           volt-  
                                                                           age_t0,  
                                                                           AS-  
                                                                           cur-  
                                                                           rents_t0,  
                                                                           inj,  
                                                                           a_spike,  
                                                                           b_spike,  
                                                                           a_voltage,  
                                                                           b_voltage)
```

Analytical solution for threshold dynamics. The threshold will adapt via two mechanisms: 1. a voltage dependent adaptation. 2. a component initiated by a spike which decays as an exponential. These two component are

in reference to threshold infinity and are recorded in the neuron's threshold components. The third component refers to `th_inf` which is added separately as opposed to being included in the voltage component of the threshold as is done in equation 2.1 of Mihalas and Nieber 2009. Threshold infinity is removed for simple optimization.

Parameters

neuron [class]
threshold_t0 [float] threshold input to function
voltage_t0 [float] voltage input to function
AScurrents_t0 [vector] values of after spike currents
inj [float] current injected into the neuron

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_exact(neuron,
                                                                      volt-
                                                                      age_t0,
                                                                      AS-
                                                                      cur-
                                                                      rents_t0,
                                                                      inj)
```

(TODO) Linear voltage dynamics.

```
allensdk.model.glif.glif_neuron_methods.dynamics_voltage_linear_forward_euler(neuron,
                                                                                volt-
                                                                                age_t0,
                                                                                AS-
                                                                                cur-
                                                                                rents_t0,
                                                                                inj)
```

(TODO) Linear voltage dynamics.

```
allensdk.model.glif.glif_neuron_methods.max_of_line_and_const(x, b, c, d)
Find the maximum of a value and a position on a line
```

Parameters

x: float x position on line 1
c: float slope of line 1
d: float y-intercept of line 1
b: float y-intercept of line 2

Returns

float the max of a line value and a constant

```
allensdk.model.glif.glif_neuron_methods.min_of_line_and_zero(x, c, d)
Find the minimum of a value and a position on a line
```

Parameters

x: float x position on line 1
c: float slope of line 1
d: float y-intercept of line 1
b: float y-intercept of line 2

Returns

float the max of a line value and a constant

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_none(neuron, AScurrents_t0)
```

Reset afterspike currents to zero.

```
allensdk.model.glif.glif_neuron_methods.reset_AScurrent_sum(neuron, AScurrents_t0, r)
```

Reset afterspike currents by adding summed exponentials. Left over currents from last spikes as well as newly initiated currents from current spike. Currents amplitudes in `neuron.asc_amp_array` need to be the amplitudes advanced though the spike cutting. I.e. In the preprocessor if the after spike currents are calculated via the GLM from spike initiation the amplitude at the time after the spike cutting needs to be calculated and `neuron.asc_amp_array` needs to be set to this value.

Parameters

r [np.ndarray] a coefficient vector applied to the afterspike currents

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_inf(neuron, threshold_t0, voltage_v1)
```

Reset the threshold to instantaneous threshold.

```
allensdk.model.glif.glif_neuron_methods.reset_threshold_three_components(neuron, threshold_t0, voltage_v1, a_spike, b_spike)
```

This method calculates the two components of the threshold: a spike (fast) component and a voltage (slow) component. The `threshold_components` vectors are then updated so that the traces match the voltage, current, and total threshold traces. The spike component of the threshold decays via an exponential fit specified by the amplitude `a_spike` and the time constant `b_spike` fit via the multiblip data. The voltage component does not change during the duration of the spike. The spike component are threshold component are summed along with threshold infinity to return the total threshold. Note that in the current implementation `a_spike` is added to the last value of the `threshold_components` which means that `a_spike` is the amplitude after spike cutting (if there is any).

Inputs:

neuron: class contains attributes of the neuron

threshold_t0, voltage_t0: float are not used but are here for consistency with other methods

a_spike: float amplitude of the exponential decay of spike component of threshold after spike cutting has been implemented.

b_spike: float amplitude of the exponential decay of spike component of threshold

Outputs:

Returns: float the total threshold which is the sum of the spike component of threshold, the voltage component of threshold and threshold infinity (with it's corresponding coefficient)

neuron.threshold_components: dictionary containing

a_spike: list vector of spiking component of threshold that corresponds to the voltage, current, and total threshold traces

b_spike: list

vector of voltage component of threshold that corresponds to the voltage, current, and total threshold traces.

Note that this function can be changed to use `a_spike` at the time of the spike and then have the spike component plus the residual decay thought the spike. There are benefits and drawbacks to this. This potential change would be beneficial as it perhaps makes more biological sense for the threshold to go up at the time of spike if the traces are ever used. Also this would mean that `a_spike` would not have to be adjusted thought the spike cutting after the multiblip fit. However the current implementation makes sense in that it is similar to how afterspike currents are implemented.

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_v_before(neuron, voltage_t0, a, b)
```

Reset voltage to the previous value with a scale and offset applied.

Parameters

a [float] voltage scale constant

b [float] voltage offset constant

```
allensdk.model.glif.glif_neuron_methods.reset_voltage_zero(neuron, voltage_t0)
```

Reset voltage to zero.

```
allensdk.model.glif.glif_neuron_methods.spike_component_of_threshold_exact(th0, b_spike, t)
```

Spike component of threshold modeled as an exponential decay. Implemented here as exact analytical solution.

Parameters

th0 [float] threshold input to function

b_spike [float] decay constant of exponential

t [float or array] time step if used in an Euler setup time if used analytically

```
allensdk.model.glif.glif_neuron_methods.spike_component_of_threshold_forward_euler(th_t0, b_spike, dt)
```

Spike component of threshold modeled as an exponential decay. Implemented here for forward Euler

Parameters

th_t0 [float] threshold input to function

b_spike [float] decay constant of exponential

dt [float] time step

```
allensdk.model.glif.glif_neuron_methods.voltage_component_of_threshold_exact(th0, v0, I, t, a_voltage, b_voltage, C, g, EI)
```

Note this function is the exact formulation; however, `dt` is used because `t0` is the initial time and `dt` is the time the function is exactly evaluated at. Note: that here, this equation is in reference to `th_inf`. Therefore `th0` is the total threshold-`thr_inf` (`threshold_inf` replaced with 0 in the equation to be verbose). This is done so that `th_inf` can be optimized without affecting this function.

Parameters

th0 [float] threshold input to function

v0 [float] voltage input to function

I [float] total current entering neuron (note if there are after spike currents these must be included in this value)

t [float or array] time step if used in an Euler setup time if used analytically

a_voltage [float] constant a

b_voltage [float] constant b

C [float] capacitance

g [float] conductance (1/resistance)

El [float] reversal potential

```
allensdk.model.glif.glif_neuron_methods.voltage_component_of_threshold_forward_euler(th_t0,
                                                                                       v_t0,
                                                                                       dt,
                                                                                       a_voltage,
                                                                                       b_voltage,
                                                                                       El)
```

Equation 2.1 of Mihalas and Nieber, 2009 implemented for use in forward Euler. Note here all variables are in reference to threshold infinity. Therefore thr_inf is zero here (replaced threshold_inf with 0 in the equation to be verbose). This is done so that th_inf can be optimized without affecting this function.

Parameters

th_t0 [float] threshold input to function

v_t0 [float] voltage input to function

dt [float] time step

a_voltage [float] constant a

b_voltage [float] constant b

El [float] reversal potential

allensdk.model.glif.simulate_neuron module

```
allensdk.model.glif.simulate_neuron.load_sweep(file_name, sweep_number)
```

Load the stimulus for a sweep from file.

```
allensdk.model.glif.simulate_neuron.main()
```

```
allensdk.model.glif.simulate_neuron.parse_arguments()
```

Use argparse to get required arguments from the command line

```
allensdk.model.glif.simulate_neuron.simulate_neuron(neuron, sweep_numbers,
                                                    input_file_name, out-
                                                    put_file_name, spike_cut_value)
```

```
allensdk.model.glif.simulate_neuron.simulate_sweep(neuron, stimulus,
                                                    spike_cut_value)
```

Simulate a neuron given a stimulus and initial conditions.

```
allensdk.model.glif.simulate_neuron.simulate_sweep_from_file(neuron,
                                                             sweep_number,
                                                             input_file_name,
                                                             output_file_name,
                                                             spike_cut_value)
```

Load a sweep stimulus, simulate the response, and write it out.

```
allensdk.model.glif.simulate_neuron.write_sweep_response(file_name,
                                                         sweep_number, response,
                                                         spike_times)
```

Overwrite the response in a file.

Module contents

A Generalized Linear Integrate and Fire (GLIF) neuron modeling package. Use this code to run the GLIF models available in the Allen Cell Types Atlas. See [Generalized LIF Models](#) for more details.

Module contents

6.1.8 allensdk.morphology package

Submodules

allensdk.morphology.validate_swc module

```
allensdk.morphology.validate_swc.main()
allensdk.morphology.validate_swc.validate_swc(swc_file)
```

To be compatible with NEURON, SWC files must have the following properties:

- 1) a single root node with parent ID '-1'
- 2) sequentially increasing ID numbers
- 3) immediate children of the soma cannot branch

Module contents

6.1.9 allensdk.mouse_connectivity package

Subpackages

allensdk.mouse_connectivity.grid package

Subpackages

allensdk.mouse_connectivity.grid.subimage package

Submodules

allensdk.mouse_connectivity.grid.subimage.base_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.base_subimage.IntensitySubImage(reduce_level,
                                         in_dims,
                                         in_spacing,
                                         coarse_spacing,
                                         in-
                                         ten-
                                         sity_paths,
                                         *args,
                                         **kwargs)
```

Bases: *allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage*

```
get_intensity(self)
```

```
required_intensities = []
```

```
setup_images(self)
```

```
class allensdk.mouse_connectivity.grid.subimage.base_subimage.PolygonSubImage(reduce_level,
                                         in_dims,
                                         in_spacing,
                                         coarse_spacing,
                                         poly-
                                         gon_info,
                                         *args,
                                         **kwargs)
```

Bases: *allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage*

```
get_polygons(self)
```

```
optional_polys = []
```

```
required_polys = []
```

```
setup_images(self)
```

```
class allensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentationSubImage(reduce_level,
                                         in_dims,
                                         in_spacing,
                                         coarse_spac-
                                         seg-
                                         men-
                                         ta-
                                         tion_paths,
                                         *args,
                                         **kwargs)
```

Bases: *allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage*

```
extract_injection_from_segmentation(self, segmentation_name='segmentation', injec-
                                         tion_name='injection')
```

Notes

Currently, the segmentation uses a series of codes to map 8-bit values onto meaningful classifications. The code for signal pixels is a 1 in at least one of the 5 rightmost bits.

```
extract_signal_from_segmentation(self, segmentation_name='segmentation', sig-
                                         nal_name='signal')
```

Notes

Currently, the segmentation uses a series of codes to map 8-bit values onto meaningful classifications. The code for signal pixels is a 1 in the leftmost bit.

In some cases, bit 5 indicates that the pixel was not removed in a posfiltering process. Optionally, this postfilter can be applied in gridding.

get_segmentation (*self*)

process_segmentation (*self*)

read_segmentation_image (*self*, *segmentation_name*='segmentation')

Notes

We downsample in memory rather than using the jp2 pyramid because the segmentation is a label image.

required_segmentations = []

setup_images (*self*)

```
class allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage (reduce_level,
                                                                    in_dims,
                                                                    in_spacing,
                                                                    coarse_spacing,
                                                                    *args,
                                                                    **kwargs)
```

Bases: object

apply_mask (*self*, *image_name*, *mask_name*, *positive*=True)

apply_pixel_counter (*self*, *accumulator_name*, *image*)

binarize (*self*, *image_name*)

compute_coarse_planes (*self*)

make_pixel_counter (*self*)

pixel_counter

setup_images (*self*)

```
allensdk.mouse_connectivity.grid.subimage.base_subimage.run_subimage (input_data)
```

allensdk.mouse_connectivity.grid.subimage.cav_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage (reduce_level,
                                                                    in_dims,
                                                                    in_spacing,
                                                                    coarse_spacing,
                                                                    polygon_info,
                                                                    *args,
                                                                    **kwargs)
```

Bases: `allensdk.mouse_connectivity.grid.subimage.base_subimage.`

`PolygonSubImage`

compute_coarse_planes (*self*)

```
required_polys = ['missing_tile', 'cav_tracer']
```

`allensdk.mouse_connectivity.grid.subimage.classic_subimage` module

```
class allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage(reduce_level,  
in_dims,  
in_spacing,  
coarse_spacing,  
poly-  
gon_info,  
seg-  
men-  
ta-  
tion_paths,  
in-  
ten-  
sity_paths,  
in-  
jec-  
tion_polygon,  
*args,  
**kwargs)
```

```
Bases:      allensdk.mouse_connectivity.grid.subimage.base_subimage.  
IntensitySubImage, allensdk.mouse_connectivity.grid.subimage.base_subimage.  
SegmentationSubImage,      allensdk.mouse_connectivity.grid.subimage.  
base_subimage.PolygonSubImage
```

```
compute_coarse_planes(self)
```

```
compute_injection(self)
```

```
compute_intensity(self)
```

```
compute_projection(self)
```

```
compute_sum_pixels(self)
```

```
optional_polys = ['aav_tracer']
```

```
process_segmentation(self)
```

```
required_intensities = ['green']
```

```
required_polys = ['missing_tile', 'no_signal', 'aav_exclusion']
```

```
required_segmentations = ['segmentation']
```


allensdk.mouse_connectivity.grid.subimage.count_subimage module

```
class allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage (reduce_level,
                                                                    in_dims,
                                                                    in_spacing,
                                                                    coarse_spacing,
                                                                    poly-
                                                                    gon_info,
                                                                    seg-
                                                                    men-
                                                                    ta-
                                                                    tion_paths,
                                                                    in-
                                                                    jec-
                                                                    tion_polygon_key=''
                                                                    *args,
                                                                    **kwargs)
```

```
Bases:      allensdk.mouse_connectivity.grid.subimage.base_subimage.
             SegmentationSubImage,      allensdk.mouse_connectivity.grid.subimage.
             base_subimage.PolygonSubImage
```

```
compute_coarse_planes (self)
compute_injection (self)
compute_projection (self)
compute_sum_pixels (self)
process_segmentation (self)
required_polys = ['missing_tile', 'no_signal', 'aav_exclusion']
required_segmentations = ['segmentation']
```

Module contents

```
allensdk.mouse_connectivity.grid.subimage.run_subimage (input_data)
```

allensdk.mouse_connectivity.grid.utilities package**Submodules****allensdk.mouse_connectivity.grid.utilities.downsampling_utilities module**

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.apply_divisions (image,
                                                                    win-
                                                                    dow_size)
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.block_average (volume,
                                                                    fac-
                                                                    tor)
```

```
allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.conv(image,
                                                                    fac-
                                                                    tor,
                                                                    win-
                                                                    dow_size)

allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.downsample_average(volume,
                                                                    cur-
                                                                    rent_spa-
                                                                    ce,
                                                                    tar-
                                                                    get_spa-
                                                                    ce)

allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.extract(image,
                                                                    fac-
                                                                    tor,
                                                                    win-
                                                                    dow_size,
                                                                    win-
                                                                    dow_step,
                                                                    out-
                                                                    put_shape)

allensdk.mouse_connectivity.grid.utilities.downsampling_utilities.window_average(volume,
                                                                    fac-
                                                                    tor)
```

allensdk.mouse_connectivity.grid.utilities.image_utilities module

```
allensdk.mouse_connectivity.grid.utilities.image_utilities.block_apply(in_image,
                                                                    out_shape,
                                                                    dtype,
                                                                    blocks,
                                                                    fn)

allensdk.mouse_connectivity.grid.utilities.image_utilities.build_affine_transform(aff_params)

allensdk.mouse_connectivity.grid.utilities.image_utilities.build_composite_transform(dfmfield,
                                                                    aff_params)

allensdk.mouse_connectivity.grid.utilities.image_utilities.compute_coarse_parameters(in_dims,
                                                                    in_spacing,
                                                                    out_spacing,
                                                                    reduce_level)

allensdk.mouse_connectivity.grid.utilities.image_utilities.grid_image_blocks(in_shape,
                                                                    in_spacing,
                                                                    out_spacing)

allensdk.mouse_connectivity.grid.utilities.image_utilities.image_from_array(array,
                                                                    spacing,
                                                                    origin,
                                                                    gin=True)
```

```

allensdk.mouse_connectivity.grid.utilities.image_utilities.new_image(dims,
                                                                    spacing,
                                                                    dtype,
                                                                    origin=True)

allensdk.mouse_connectivity.grid.utilities.image_utilities.np_sitk_convert(np_type)
allensdk.mouse_connectivity.grid.utilities.image_utilities.rasterize_polygons(shape,
                                                                    scale,
                                                                    polys)

allensdk.mouse_connectivity.grid.utilities.image_utilities.read_intensity_image()
allensdk.mouse_connectivity.grid.utilities.image_utilities.read_segmentation_image(path)
allensdk.mouse_connectivity.grid.utilities.image_utilities.resample_into_volume(image,
                                                                    transform,
                                                                    z,
                                                                    vol,
                                                                    dtype=8)

allensdk.mouse_connectivity.grid.utilities.image_utilities.resample_volume(volume,
                                                                    dims,
                                                                    spacing,
                                                                    interpolator=None,
                                                                    transform=None)

allensdk.mouse_connectivity.grid.utilities.image_utilities.set_image_spacing(image,
                                                                    spacing,
                                                                    origin=True)

allensdk.mouse_connectivity.grid.utilities.image_utilities.sitk_np_convert(sitk_type)
allensdk.mouse_connectivity.grid.utilities.image_utilities.write_volume(volume,
                                                                    name,
                                                                    prefix=None,
                                                                    specify_resolution=None,
                                                                    extension='.nrrd',
                                                                    paths=None)

```

Module contents

allensdk.mouse_connectivity.grid.writers package

Module contents

```
allensdk.mouse_connectivity.grid.writers.cav_writer(gridder, grid_prefix, accumula-
                                                    tor_prefix, **kwargs)
allensdk.mouse_connectivity.grid.writers.classic_writer(gridder,      grid_prefix,
                                                         accumulator_prefix, tar-
                                                         get_spacings, **kwargs)
allensdk.mouse_connectivity.grid.writers.count_writer(gridder, grid_prefix, accumu-
                                                         lator_prefix, target_spacings,
                                                         **kwargs)
allensdk.mouse_connectivity.grid.writers.handle_pyramid(isg, key, target_spacings,
                                                         prefix, paths)
allensdk.mouse_connectivity.grid.writers.ratio_and_pyramid(isg, num, den, out,
                                                         accumulator_prefix,
                                                         grid_prefix,      tar-
                                                         get_spacings, paths)
```

Submodules

allensdk.mouse_connectivity.grid.image_series_gridder module

```
class allensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder(in_dims,
                                                                                   in_spacing,
                                                                                   out_dims,
                                                                                   out_spacing,
                                                                                   re-
                                                                                   duce_level,
                                                                                   subim-
                                                                                   ages,
                                                                                   subim-
                                                                                   age_kwargs,
                                                                                   npro-
                                                                                   cesses,
                                                                                   affine_params,
                                                                                   dfm-
                                                                                   fld_path)
```

Bases: object

```
accumulator_to_numpy(self, key, cb)
build_coarse_grids(self)
consume_volume(self, key, cb)
initialize_coarse_volume(self, key, dtype)
make_ratio_volume(self, num_key, den_key, ratio_key)
    assume parents numpified
paste_slice(self, key, index, slice_array)
```

```

paste_subimage (self, index, output)
    Inserts planar accumulators into coarse grid volumes

resample_volume (self, key)

set_coarse_grid_parameters (self)

setup_subimages (self)

transform

```

Module contents

Module contents

6.1.10 allensdk.test_utilities package

Submodules

allensdk.test_utilities.custom_comparators module

```

class allensdk.test_utilities.custom_comparators.WhitespaceStrippedString (string:
                                                                    str,
                                                                    whites-
                                                                    pace_chars:
                                                                    str
                                                                    =
                                                                    '\s',
                                                                    ASCII:
                                                                    bool
                                                                    =
                                                                    False)

```

Bases: object

Comparator class to compare strings that have been stripped of whitespace. By default removes any unicode whitespace character that matches the regex `s`, (which includes `[` `]`, and other unicode whitespace characters).

allensdk.test_utilities.regression_fixture module

```
allensdk.test_utilities.regression_fixture.get_list_of_path_dict()
```

allensdk.test_utilities.temp_dir module

```
allensdk.test_utilities.temp_dir.temp_dir(request)
```

Module contents

6.2 Submodules

6.2.1 allensdk.deprecated module

`allensdk.deprecated.class_deprecated(message=None)`

`allensdk.deprecated.deprecated(message=None)`

`allensdk.deprecated.legacy(message=None)`

6.2.2 allensdk.tmp module

6.3 Module contents


exception `allensdk.OneResultExpectedError`

Bases: `RuntimeError`

`allensdk.one(x)`

The Allen Software Development Kit houses source code for reading and processing Allen Brain Atlas data. The Allen SDK focuses on the Allen Brain Observatory, Cell Types Database, and Mouse Brain Connectivity Atlas.

Attention: As of October 2019, we have dropped Python 2 support and any files with a py2 dependency (for example analysis files) have been updated.




`_static/sdk_cam.png`

Allen Brain Observatory

The [Allen Brain Observatory](#) is a data resource for understanding sensory processing in the mouse visual cortex. This study systematically measures visual responses in multiple cortical areas and layers using two-photon calcium imaging of GCaMP6-labeled neurons targeted using Cre driver lines. Response characterizations include orientation tuning, spatial and temporal frequency tuning, temporal dynamics, and spatial receptive field structure.

The mean fluorescence traces for all segmented cells are available in the Neurodata Without Borders file format ([NWB files](#)). These files contain standardized descriptions of visual stimuli to support stimulus-specific tuning analysis. The Allen SDK provides code to:

- download and organize experiment data according to cortical area, imaging depth, and Cre line
- remove the contribution of neuropil signal from fluorescence traces
- access (or compute) dF/F traces based on the neuropil-corrected traces
- perform stimulus-specific tuning analysis (e.g. drifting grating direction tuning)

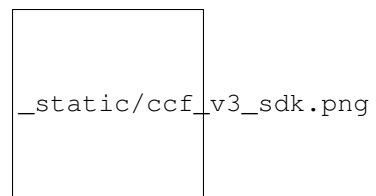


Allen Cell Types Database

The [Allen Cell Types Database](#) contains electrophysiological and morphological characterizations of individual neurons in the mouse primary visual cortex. The Allen SDK provides Python code for accessing electrophysiology measurements ([NWB files](#)) for all neurons and morphological reconstructions ([SWC files](#)) for a subset of neurons.

The Database also contains two classes of models fit to this data set: biophysical models produced using the NEURON simulator and generalized leaky integrate and fire models (GLIFs) produced using custom Python code provided with this toolkit.

The Allen SDK provides sample code demonstrating how to download neuronal model parameters from the Allen Brain Atlas API and run your own simulations using stimuli from the Allen Cell Types Database or custom current injections:

- *Biophysical Models*
- *Generalized LIF Models*

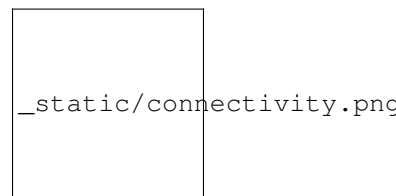

Allen Mouse Brain Connectivity Atlas

The [Allen Mouse Brain Connectivity Atlas](#) is a high-resolution map of neural connections in the mouse brain. Built on an array of transgenic mice genetically engineered to target specific cell types, the Atlas comprises a unique compendium of projections from selected neuronal populations throughout the brain. The primary data of the Atlas consists of high-resolution images of axonal projections targeting different anatomic regions or various cell types using Cre-dependent specimens. Each data set is processed through an informatics data analysis pipeline to obtain spatially mapped quantified projection information.

The Allen SDK provides Python code for accessing experimental metadata along with projection signal volumes registered to a common coordinate framework. This framework has structural annotations, which allows users to compute structure-level signal statistics.

See the [mouse connectivity section](#) for more details.

CHAPTER 10

What's New - 2.3.1 (October 13, 2020)

As of the 2.3.1 release:

- (Internal) Fixed a write_nwb bug for behavior ophys experiments involving the BehaviorOphysJsonApi expecting a mesoscope-specific method.

What's New - 2.3.0 (October 9, 2020)

As of the 2.3.0 release:

- Visual behavior running speed is now low-pass filtered at 10Hz. The raw running speed data is still available. The running speed is corrected for encoder threshold crossing artifacts.
- Support for stimulus gratings for visual behavior.
- Fixed an eye-tracking sync problem.
- Updates to some visual behavior pynwb implementations.
- Adds load sync data for individual plane sets to relate accurate event timings to mesoscope data.
- Adds public API method to access the `behavior_session_id` from an instance of `BehaviorOphysSession`.

What's New - 2.2.0 (September 3, 2020)

As of the 2.2.0 release:

- AllenSDK HTTP engine streaming requests now include a progress bar
- *ImportError: cannot import name 'MultiContainerInterface' from 'hdmf.container'* errors should now be resolved (by removing explicit version bounds on the *hdmf* package).
- The optical physiology 2-photon trace demixer has been modified to be more memory friendly and should no longer result in out of memory errors when trying to demix very large movie stacks.

CHAPTER 13

What's New - 2.1.0 (July 16, 2020)

As of the 2.1.0 release:

- behavior ophys nwb files can now be written using updated pynwb and hdmf
- A warning has been added if you are using AllenSDK with outdated NWB files
- A new documentation file has been added which will contain Visual Behavior specific terms for quick lookup

What's New - 2.0.0 (June 11, 2020)

As of the 2.0.0 release:

- pynwb and hdmf version requirements have been made less strict
- The organization of data for ecephys neuropixels Neurodata Without Borders (NWB) files has been significantly changed to conform with NWB specifications and best practices
- CCF locations for ecephys neuropixels electrodes are now written to NWB files
- Examples for accessing eye tracking ellipse fit and screen gaze location data have been added to ecephys example notebooks

Important Note: Due to newer versions of pynwb/hdmf having issues reading previously released Visual Coding Neuropixels NWB files and due to the significant reorganization of their NWB file contents, this release contains breaking changes that necessitate a major version revision. NWB files released prior to 6/11/2020 are not guaranteed to work with the 2.0.0 version of AllenSDK. If you cannot or choose not to re-download the updated NWB files, you can continue using a prior version of AllenSDK (< 2.0.0) to access them. However, no further features or bugfixes for AllenSDK (< 2.0.0) are planned. Data released for other projects (Cell Types, Mouse Connectivity, etc.) are *NOT* affected and will *NOT* need to be re-downloaded

CHAPTER 15

Previous Release Notes

- 1.8.0
- 1.7.1
- 1.7.0
- 1.6.0
- 1.5.0
- 1.4.0
- 1.3.0
- 1.2.0
- 1.1.1
- 1.1.0
- 1.0.2
- 0.16.3
- 0.16.2
- 0.16.1
- 0.16.0
- 0.14.5
- 0.14.4
- 0.14.3
- 0.14.2
- 0.13.2
- 0.13.1
- 0.13.0

- 0.12.4

Bibliography

[1] Allen Brain Atlas Data Portal: [Downloading a WellKnownFile](#).

Python Module Index

a

allensdk, 338
allensdk.api, 80
allensdk.api.api, 72
allensdk.api.cache, 76
allensdk.api.caching_utilities, 79
allensdk.api.queries, 72
allensdk.api.queries.annotated_section_data_sets_api, 43
allensdk.api.queries.biophysical_api, 44
allensdk.api.queries.brain_observatory_api, 46
allensdk.api.queries.cell_types_api, 49
allensdk.api.queries.connected_services, 51
allensdk.api.queries.glif_api, 52
allensdk.api.queries.grid_data_api, 52
allensdk.api.queries.image_download_api, 54
allensdk.api.queries.mouse_atlas_api, 57
allensdk.api.queries.mouse_connectivity_api, 58
allensdk.api.queries.ontologies_api, 61
allensdk.api.queries.reference_space_api, 63
allensdk.api.queries.rma_api, 65
allensdk.api.queries.rma_pager, 70
allensdk.api.queries.rma_template, 70
allensdk.api.queries.svg_api, 70
allensdk.api.queries.synchronization_api, 70
allensdk.api.queries.tree_search_api, 72
allensdk.brain_observatory, 179
allensdk.brain_observatory.argschema_utilities, 146
allensdk.brain_observatory.behavior, 109
allensdk.brain_observatory.behavior.behavior_data_sets_api, 88
allensdk.brain_observatory.behavior.behavior_ophys_api, 80
allensdk.brain_observatory.behavior.criteria, 90
allensdk.brain_observatory.behavior.dprime, 92
allensdk.brain_observatory.behavior.eye_tracking_page, 93
allensdk.brain_observatory.behavior.image_api, 96
allensdk.brain_observatory.behavior.internal, 84
allensdk.brain_observatory.behavior.internal.behavior_data_sets_api, 81
allensdk.brain_observatory.behavior.internal.behavior_data_sets_api, 82
allensdk.brain_observatory.behavior.metadata_processing, 96
allensdk.brain_observatory.behavior.mtrain, 96
allensdk.brain_observatory.behavior.rewards_processing, 98
allensdk.brain_observatory.behavior.running_processes, 98
allensdk.brain_observatory.behavior.schemas, 99
allensdk.brain_observatory.behavior.session_metrics, 102
allensdk.brain_observatory.behavior.stimulus_processing, 102
allensdk.brain_observatory.behavior.sync, 84
allensdk.brain_observatory.behavior.sync.process_sync, 84
allensdk.brain_observatory.behavior.trial_masks, 84
allensdk.brain_observatory.behavior.trials_processing, 84

b

105 allensdk.brain_observatory.ecephys.stimulus_table.
 allensdk.brain_observatory.behavior.write_nwb, 125
 88 allensdk.brain_observatory.ecephys.stimulus_table.
 allensdk.brain_observatory.brain_observatory_exceptions, 126
 147 allensdk.brain_observatory.ecephys.stimulus_table.
 allensdk.brain_observatory.brain_observatory_plotting, 127
 147 allensdk.brain_observatory.ecephys.stimulus_table.v
 allensdk.brain_observatory.chisquare_categorical, 121
 148 allensdk.brain_observatory.ecephys.stimulus_table.v
 allensdk.brain_observatory.circle_plots, 120
 148 allensdk.brain_observatory.ecephys.visualization,
 allensdk.brain_observatory.demixer, 150 128
 allensdk.brain_observatory.dff, 151 allensdk.brain_observatory.ecephys.write_nwb,
 allensdk.brain_observatory.drifting_gratings, 129
 153 allensdk.brain_observatory.extract_running_speed,
 allensdk.brain_observatory.ecephys, 130 130
 allensdk.brain_observatory.ecephys.alignment_utils, 155
 113 allensdk.brain_observatory.findlevel,
 allensdk.brain_observatory.ecephys.alignment_utils.brain_observatory.gaze_mapping,
 109 130
 allensdk.brain_observatory.ecephys.alignment_utils.brain_observatory.gaze_mapping,
 111 155
 allensdk.brain_observatory.ecephys.alignment_utils.brain_observatory.natural_movie,
 112 156
 allensdk.brain_observatory.ecephys.alignment_utils.brain_observatory.natural_scenes,
 113 157
 allensdk.brain_observatory.ecephys.copy_allensdk.brain_observatory.nwb, 132
 113 allensdk.brain_observatory.nwb.behavior_ophys_nwb_e
 allensdk.brain_observatory.ecephys.current_source_density, 134
 113 allensdk.brain_observatory.nwb.metadata,
 allensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api,
 114 allensdk.brain_observatory.nwb.nwb_api,
 allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api,
 115 allensdk.brain_observatory.nwb.nwb_utils,
 allensdk.brain_observatory.ecephys.file_io, 131
 117 allensdk.brain_observatory.nwb.schemas,
 allensdk.brain_observatory.ecephys.file_io.continuous_file, 132
 116 allensdk.brain_observatory.observatory_plots,
 allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset, 158
 116 allensdk.brain_observatory.ophys, 136
 allensdk.brain_observatory.ecephys.file_io.allensdk.brain_observatory.ophys.trace_extraction,
 116 136
 allensdk.brain_observatory.ecephys.lfp_subsampling, 160
 119 allensdk.brain_observatory.r_neuropil,
 allensdk.brain_observatory.ecephys.lfp_subsampling.brain_observatory.receptive_field_analysis,
 117 144
 allensdk.brain_observatory.ecephys.optotagging, 136
 119 allensdk.brain_observatory.receptive_field_analysis,
 allensdk.brain_observatory.ecephys.stimulus_sync, 139
 129 allensdk.brain_observatory.receptive_field_analysis,
 allensdk.brain_observatory.ecephys.stimulus_table, 140
 128 allensdk.brain_observatory.receptive_field_analysis,
 allensdk.brain_observatory.ecephys.stimulus_table.ecephys_prespikes, 140
 121 allensdk.brain_observatory.receptive_field_analysis

[allensdk.brain_observatory.receptive_field_utilities](#), 141
[allensdk.brain_observatory.receptive_field_utilities](#), 141
[allensdk.brain_observatory.receptive_field_utilities](#), 142
[allensdk.brain_observatory.receptive_field_utilities](#), 142
[allensdk.brain_observatory.receptive_field_utilities](#), 143
[allensdk.brain_observatory.roi_masks](#), 161
[allensdk.brain_observatory.running_speed](#), 164
[allensdk.brain_observatory.session_analysis](#), 164
[allensdk.brain_observatory.session_api_utilities](#), 167
[allensdk.brain_observatory.static_gratings](#), 168
[allensdk.brain_observatory.stimulus_analysis](#), 169
[allensdk.brain_observatory.stimulus_info](#), 171
[allensdk.brain_observatory.sync_dataset](#), 175
[allensdk.brain_observatory.sync_utilities](#), 145
[allensdk.brain_observatory.visualization](#), 145
C
[allensdk.config](#), 187
[allensdk.config.app](#), 181
[allensdk.config.app.application_config](#), 179
[allensdk.config.manifest](#), 185
[allensdk.config.manifest_builder](#), 187
[allensdk.config.model](#), 185
[allensdk.config.model.description](#), 183
[allensdk.config.model.description_parser](#), 184
[allensdk.config.model.formats](#), 183
[allensdk.config.model.formats.hdf5_util](#), 182
[allensdk.config.model.formats.json_description_parser](#), 182
[allensdk.config.model.formats.pycfg_description_parser](#), 183
[allensdk.core](#), 227
[allensdk.core.auth_config](#), 188
[allensdk.core.authentication](#), 188
[allensdk.core.brain_observatory_cache](#), 189
[allensdk.core.brain_observatory_nwb_data_set](#), 194
[allensdk.core.cache_utilities](#), 198
[allensdk.core.cell_types_cache](#), 198
[allensdk.core.dat_utilities](#), 201
[allensdk.core.exceptions](#), 201
[allensdk.core.h5_utilities](#), 201
[allensdk.core.visualization_utilities](#), 202
[allensdk.core.lazy_property](#), 188
[allensdk.core.lazy_property.lazy_property](#), 188
[allensdk.core.lazy_property.lazy_property_mixin](#), 188
[allensdk.core.mouse_connectivity_cache](#), 203
[allensdk.core.nwb_data_set](#), 208
[allensdk.core.obj_utilities](#), 210
[allensdk.core.ontology](#), 210
[allensdk.core.ophys_experiment_session_id_mapping](#), 211
[allensdk.core.reference_space](#), 211
[allensdk.core.reference_space_cache](#), 214
[allensdk.core.simple_tree](#), 216
[allensdk.core.sitk_utilities](#), 219
[allensdk.core.structure_tree](#), 220
[allensdk.core.swc](#), 222
[allensdk.core.typing](#), 227
d
[allensdk.deprecated](#), 338
e
[allensdk.ephys](#), 239
[allensdk.ephys.ephys_extractor](#), 227
[allensdk.ephys.ephys_features](#), 231
[allensdk.ephys.extract_cell_features](#), 238
[allensdk.ephys.feature_extractor](#), 239
i
[allensdk.internal](#), 315
[allensdk.internal.api](#), 250
[allensdk.internal.api.api_prerelease](#), 244
[allensdk.internal.api.behavior_data_lims_api](#), 245
[allensdk.internal.api.behavior_lims_api](#), 247
[allensdk.internal.api.behavior_ophys_api](#), 247
[allensdk.internal.api.lims_api](#), 248
[allensdk.internal.api.mtrain_api](#), 249

[allensdk.internal.api.ophys_lims_api](#), 249
[allensdk.internal.api.queries](#), 244
[allensdk.internal.api.queries.biophysical_model_reader](#), 240
[allensdk.internal.api.queries.biophysical_model_reader](#), 240
[allensdk.internal.api.queries.grid_data_api_pre_release](#), 241
[allensdk.internal.api.queries.mouse_connectivity_pre_release](#), 242
[allensdk.internal.api.queries.optimize_configuration_pre_release](#), 243
[allensdk.internal.api.queries.pre_release](#), 244
[allensdk.internal.brain_observatory](#), 265
[allensdk.internal.brain_observatory.annotated_region_metrics](#), 251
[allensdk.internal.brain_observatory.demix_report](#), 252
[allensdk.internal.brain_observatory.demixer](#), 253
[allensdk.internal.brain_observatory.eye_calibration](#), 254
[allensdk.internal.brain_observatory.fit_ellipse](#), 257
[allensdk.internal.brain_observatory.frame_stream](#), 257
[allensdk.internal.brain_observatory.itracker_utils](#), 258
[allensdk.internal.brain_observatory.mask_set](#), 260
[allensdk.internal.brain_observatory.ophys_session_decomposition](#), 260
[allensdk.internal.brain_observatory.resources](#), 251
[allensdk.internal.brain_observatory.roi_filter](#), 261
[allensdk.internal.brain_observatory.time_sync](#), 263
[allensdk.internal.core](#), 268
[allensdk.internal.core.lims_pipeline_modules](#), 265
[allensdk.internal.core.lims_utilities](#), 265
[allensdk.internal.core.mouse_connectivity_checker_pre_release](#), 266
[allensdk.internal.core.simpletree](#), 267
[allensdk.internal.core.swc](#), 268
[allensdk.internal.ephys](#), 274
[allensdk.internal.ephys.core_feature_extract](#), 268
[allensdk.internal.ephys.plot_qc_figures](#), 269
[allensdk.internal.ephys.plot_qc_figures3](#), 271
[allensdk.internal.model](#), 291
[allensdk.internal.model.AIC](#), 289
[allensdk.internal.model.biophysical](#), 278
[allensdk.internal.model.biophysical.biophysical_array_reader](#), 276
[allensdk.internal.model.biophysical.check_fit_shift](#), 276
[allensdk.internal.model.biophysical.deap_utils](#), 276
[allensdk.internal.model.biophysical.ephys_utils](#), 276
[allensdk.internal.model.biophysical.fits](#), 274
[allensdk.internal.model.biophysical.fits.fit_styles](#), 274
[allensdk.internal.model.biophysical.make_deap_fit](#), 277
[allensdk.internal.model.biophysical.passive_fitting](#), 276
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.passive_fitting](#), 275
[allensdk.internal.model.biophysical.run_optimize](#), 277
[allensdk.internal.model.biophysical.run_optimize_w](#), 277
[allensdk.internal.model.biophysical.run_passive_fit](#), 278
[allensdk.internal.model.biophysical.run_simulate_l](#), 278
[allensdk.internal.model.biophysical.run_simulate_w](#), 278
[allensdk.internal.model.data_access](#), 290
[allensdk.internal.model.glif](#), 289
[allensdk.internal.model.glif.are_two_lists_of_arrays](#), 279
[allensdk.internal.model.glif.error_functions](#), 279
[allensdk.internal.model.glif.find_spikes](#), 279
[allensdk.internal.model.glif.find_sweeps](#), 279

279
 allensdk.internal.model.glif.glif_experiment, 305
 280
 allensdk.internal.model.glif.glif_optimizer, 307
 282
 allensdk.internal.model.glif.glif_optimizer_neuron, 305
 282
 allensdk.internal.model.glif.MLIN, 278 306
 allensdk.internal.model.glif.plotting, allensdk.internal.pipeline_modules, 315
 285
 allensdk.internal.model.glif.rc, 286 308
 allensdk.internal.model.glif.spike_cutting, allensdk.internal.pipeline_modules.gbm.generate_gbr
 286 307
 allensdk.internal.model.glif.threshold_adapter, allensdk.internal.pipeline_modules.gbm.generate_gbr
 287 307
 allensdk.internal.model.GLM, 290 allensdk.internal.pipeline_modules.gbm.generate_gbr
 allensdk.internal.morphology, 298 308
 allensdk.internal.morphology.compartment, allensdk.internal.pipeline_modules.run_annotated_re
 291 308
 allensdk.internal.morphology.morphology, allensdk.internal.pipeline_modules.run_demixing,
 291 308
 allensdk.internal.morphology.morphvis, allensdk.internal.pipeline_modules.run_dff_computat
 295 308
 allensdk.internal.morphology.node, 298 allensdk.internal.pipeline_modules.run_neuropil_co
 allensdk.internal.morphology.validate_sw, 309
 298
 allensdk.internal.mouse_connectivity, 309
 307
 allensdk.internal.mouse_connectivity.interval_uni
 303 310
 allensdk.internal.mouse_connectivity.interval_uni
 299 312
 allensdk.internal.mouse_connectivity.interval_uni
 299 312
 allensdk.internal.mouse_connectivity.interval_uni
 301 313
 allensdk.internal.mouse_connectivity.interval_uni
 301 315
 allensdk.internal.mouse_connectivity.interval_uni
 302 316
 allensdk.internal.mouse_connectivity.interval_uni
 302 317
 allensdk.internal.mouse_connectivity.projection
 305 316
 allensdk.internal.mouse_connectivity.projection
 303 316
 allensdk.internal.mouse_connectivity.projection
 303 315
 allensdk.internal.mouse_connectivity.projection
 304 319
 allensdk.internal.mouse_connectivity.projection
 304 317
 allensdk.internal.mouse_connectivity.projection
 305 318
 allensdk.internal.mouse_connectivity.projection
 305 319

m

`allensdk.model.glif.glif_neuron`, [319](#)
`allensdk.model.glif.glif_neuron_methods`,
 [323](#)
`allensdk.model.glif.simulate_neuron`, [328](#)
`allensdk.morphology`, [329](#)
`allensdk.morphology.validate_swc`, [329](#)
`allensdk.mouse_connectivity`, [337](#)
`allensdk.mouse_connectivity.grid`, [337](#)
`allensdk.mouse_connectivity.grid.image_series_gridder`,
 [336](#)
`allensdk.mouse_connectivity.grid.subimage`,
 [333](#)
`allensdk.mouse_connectivity.grid.subimage.base_subimage`,
 [330](#)
`allensdk.mouse_connectivity.grid.subimage.cav_subimage`,
 [331](#)
`allensdk.mouse_connectivity.grid.subimage.classic_subimage`,
 [332](#)
`allensdk.mouse_connectivity.grid.subimage.count_subimage`,
 [333](#)
`allensdk.mouse_connectivity.grid.utilities`,
 [335](#)
`allensdk.mouse_connectivity.grid.utilities.downsampling_utilities`,
 [333](#)
`allensdk.mouse_connectivity.grid.utilities.image_utilities`,
 [334](#)
`allensdk.mouse_connectivity.grid.writers`,
 [336](#)

t

`allensdk.test_utilities`, [338](#)
`allensdk.test_utilities.custom_comparators`,
 [337](#)
`allensdk.test_utilities.regression_fixture`,
 [337](#)
`allensdk.test_utilities.temp_dir`, [337](#)

A

`ab_from_diagonals()` (in module `allensdk.brain_observatory.r_neuropil`), 160
`ab_from_T()` (in module `allensdk.brain_observatory.r_neuropil`), 160
`accumulator_to_numpy()` (in module `allensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`), 336
`acquisition_rate` (in module `allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`), 169
`actual_parameters_from_normalized()` (in module `allensdk.internal.model.biophysical.deap_utils.Utils`), 276
`adaptation_index()` (in module `allensdk.ephys.feature_extractor.EphysFeatureExtractor`), 239
`adaptation_index()` (in module `allensdk.ephys.ephys_features`), 231
`add()` (`allensdk.brain_observatory.stimulus_info.BinaryIntervalSearchTree`), 171
`add_angle_labels()` (in module `allensdk.brain_observatory.circle_plots`), 149
`add_arrow()` (in module `allensdk.brain_observatory.circle_plots`), 149
`add_average_image()` (in module `allensdk.brain_observatory.nwb`), 132
`add_cell_specimen_table()` (in module `allensdk.brain_observatory.nwb`), 132
`add_corrected_fluorescence_traces()` (in module `allensdk.brain_observatory.nwb`), 132
`add_dff_traces()` (in module `allensdk.brain_observatory.nwb`), 133
`add_eye_gaze_data_interfaces()` (in module `allensdk.brain_observatory.nwb`), 133
`add_eye_gaze_mapping_data_to_nwbfile()` (in module `allensdk.brain_observatory.nwb`), 133
`add_eye_tracking_ellipse_fit_data_to_nwbfile()` (in module `allensdk.brain_observatory.nwb`), 133
`add_file()` (`allensdk.config.manifest.Manifest`), 185
`add_image()` (in module `allensdk.brain_observatory.nwb`), 133
`add_invalid_times()` (in module `allensdk.brain_observatory.nwb`), 133
`add_licks()` (in module `allensdk.brain_observatory.nwb`), 133
`add_manifest_paths()` (`allensdk.api.cache.Cache`), 76
`add_manifest_paths()` (`allensdk.core.mouse_connectivity_cache.MouseConnectivityCache`), 204
`add_manifest_paths()` (`allensdk.core.reference_space_cache.ReferenceSpaceCache`), 214
`add_manifest_paths()` (`allensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease`), 267
`add_max_projection()` (in module `allensdk.brain_observatory.nwb`), 133
`add_metadata()` (in module `allensdk.brain_observatory.nwb`), 133
`add_motion_correction()` (in module `allensdk.brain_observatory.nwb`), 133
`add_number_to_shuffled_movie()` (in module `allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities`), 125
`add_path()` (`allensdk.config.manifest.Manifest`), 185
`add_path()` (`allensdk.config.manifest_builder.ManifestBuilder`), 187
`add_paths()` (`allensdk.config.manifest.Manifest`), 185
`add_rewards()` (in module `allensdk.brain_observatory.nwb`), 133
`add_running_data_df_to_nwbfile()` (in module `allensdk.brain_observatory.nwb`), 133

ule `allensdk.brain_observatory.nwb`), 133
`add_running_speed_to_nwbfile()` (in module `allensdk.brain_observatory.nwb`), 134
`add_section()` (in module `allensdk.config.manifest_builder.ManifestBuilder` method), 187
`add_segmentation_mask_image()` (in module `allensdk.brain_observatory.nwb`), 134
`add_stimulus_index()` (in module `allensdk.brain_observatory.nwb`), 134
`add_stimulus_presentations()` (in module `allensdk.brain_observatory.nwb`), 134
`add_stimulus_template()` (in module `allensdk.brain_observatory.nwb`), 134
`add_stimulus_timestamps()` (in module `allensdk.brain_observatory.nwb`), 135
`add_task_parameters()` (in module `allensdk.brain_observatory.nwb`), 135
`add_to_fit_parameters_dict_single()` (in module `allensdk.brain_observatory.receptive_field_analysis_fit_parameters`), 140
`add_trials()` (in module `allensdk.brain_observatory.nwb`), 135
`adjust_r_for_negativity()` (in module `allensdk.internal.pipeline_modules.run_neuropil_corrections`), 309
`advance_combination()` (in module `allensdk.brain_observatory.chisquare_categorical`), 148
`AIC()` (in module `allensdk.internal.model.AIC`), 289
`AICc()` (in module `allensdk.internal.model.AIC`), 289
`align_and_cut_spikes()` (in module `allensdk.internal.model.glif.find_spikes`), 279
`align_running_speed()` (in module `allensdk.core.brain_observatory_nwb_data_set`), 197
`ALIGNMENT3D_KEY` (in module `allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), 204
`ALL` (`allensdk.api.queries.rma_api.RmaApi` attribute), 65
`all_stimuli()` (in module `allensdk.brain_observatory.stimulus_info`), 172
`AllActiveUtils` (class in module `allensdk.model.biophysical.utils`), 318
`allensdk` (module), 338
`allensdk.api` (module), 80
`allensdk.api.api` (module), 72
`allensdk.api.cache` (module), 76
`allensdk.api.caching_utilities` (module), 79
`allensdk.api.queries` (module), 72
`allensdk.api.queries.annotated_section_data_sets_api` (module), 43
`allensdk.api.queries.biophysical_api` (module), 44
`allensdk.api.queries.brain_observatory_api` (module), 46
`allensdk.api.queries.cell_types_api` (module), 49
`allensdk.api.queries.connected_services` (module), 51
`allensdk.api.queries.glif_api` (module), 52
`allensdk.api.queries.grid_data_api` (module), 52
`allensdk.api.queries.image_download_api` (module), 54
`allensdk.api.queries.mouse_atlas_api` (module), 57
`allensdk.api.queries.mouse_connectivity_api` (module), 58
`allensdk.api.queries.ontologies_api` (module), 61
`allensdk.api.queries.reference_space_api` (module), 63
`allensdk.api.queries.rma_api` (module), 65
`allensdk.api.queries.rma_pager` (module), 70
`allensdk.api.queries.rma_template` (module), 70
`allensdk.api.queries.svg_api` (module), 70
`allensdk.api.queries.synchronization_api` (module), 70
`allensdk.api.queries.tree_search_api` (module), 72
`allensdk.brain_observatory` (module), 179
`allensdk.brain_observatory.argschema_utilities` (module), 146
`allensdk.brain_observatory.behavior` (module), 109
`allensdk.brain_observatory.behavior.behavior_data_cache` (module), 88
`allensdk.brain_observatory.behavior.behavior_ophys` (module), 80
`allensdk.brain_observatory.behavior.criteria` (module), 90
`allensdk.brain_observatory.behavior.dprime` (module), 92
`allensdk.brain_observatory.behavior.eye_tracking_plots` (module), 93
`allensdk.brain_observatory.behavior.image_api` (module), 96
`allensdk.brain_observatory.behavior.internal` (module), 84
`allensdk.brain_observatory.behavior.internal.behavior_data_cache` (module), 81
`allensdk.brain_observatory.behavior.internal.behavior_ophys` (module), 81

(module), 82	(module), 113
allensdk.brain_observatory.behavior.metadata	allensdk.brain_observatory.ecephys.ecephys_project
(module), 96	(module), 114
allensdk.brain_observatory.behavior.mtrac	allensdk.brain_observatory.ecephys.ecephys_session
(module), 96	(module), 115
allensdk.brain_observatory.behavior.reward	allensdk.brain_observatory.ecephys.file_io
(module), 98	(module), 117
allensdk.brain_observatory.behavior.running	allensdk.brain_observatory.ecephys.file_io.continuous
(module), 98	(module), 116
allensdk.brain_observatory.behavior.schemas	allensdk.brain_observatory.ecephys.file_io.ecephys
(module), 99	(module), 116
allensdk.brain_observatory.behavior.sessions	allensdk.brain_observatory.ecephys.file_io.stim_files
(module), 102	(module), 116
allensdk.brain_observatory.behavior.stimulus	allensdk.brain_observatory.ecephys.lfp_subsampling
(module), 102	(module), 119
allensdk.brain_observatory.behavior.syncallensdk	allensdk.brain_observatory.ecephys.lfp_subsampling
(module), 84	(module), 117
allensdk.brain_observatory.behavior.syncapi	allensdk.brain_observatory.ecephys.optotagging_table
(module), 84	(module), 119
allensdk.brain_observatory.behavior.trials	allensdk.brain_observatory.ecephys.stimulus_sync
(module), 104	(module), 129
allensdk.brain_observatory.behavior.trials	allensdk.brain_observatory.ecephys.stimulus_table
(module), 105	(module), 128
allensdk.brain_observatory.behavior.write	allensdk.brain_observatory.ecephys.stimulus_table.
(module), 88	(module), 121
allensdk.brain_observatory.brain_observatory	allensdk.brain_observatory.ecephys.stimulus_table.m
(module), 147	(module), 125
allensdk.brain_observatory.brain_observatory	allensdk.brain_observatory.ecephys.stimulus_table.
(module), 147	(module), 126
allensdk.brain_observatory.chisquare_categorical	allensdk.brain_observatory.ecephys.stimulus_table.
(module), 148	(module), 127
allensdk.brain_observatory.circle_plots	allensdk.brain_observatory.ecephys.stimulus_table.v
(module), 148	(module), 121
allensdk.brain_observatory.demixer (module),	allensdk.brain_observatory.ecephys.stimulus_table.v
151	(module), 120
allensdk.brain_observatory.dff	allensdk.brain_observatory.ecephys.visualization
(module), 151	(module), 128
allensdk.brain_observatory.drifting_gratings	allensdk.brain_observatory.ecephys.write_nwb
(module), 153	(module), 129
allensdk.brain_observatory.ecephys (module),	allensdk.brain_observatory.extract_running_speed
130	(module), 130
allensdk.brain_observatory.ecephys.alignment	allensdk.brain_observatory.findlevel
(module), 113	(module), 155
allensdk.brain_observatory.ecephys.alignment	allensdk.brain_observatory.gaze_mapping
(module), 109	(module), 130
allensdk.brain_observatory.ecephys.alignment	allensdk.brain_observatory.gaze_mapping
(module), 111	(module), 130
allensdk.brain_observatory.ecephys.alignment	allensdk.brain_observatory.gaze_mapping
(module), 112	(module), 130
allensdk.brain_observatory.ecephys.alignment	allensdk.brain_observatory.gaze_mapping
(module), 113	(module), 130
allensdk.brain_observatory.ecephys.copy_allensdk	allensdk.brain_observatory.nwb
(module), 113	(module), 132
allensdk.brain_observatory.ecephys.current	allensdk.brain_observatory.nwb.behavior_ophys_nwb
(module), 113	(module), 132

(module), 131
 allensdk.brain_observatory.nwb.metadata (module), 131
 allensdk.brain_observatory.nwb.nwb_api (module), 131
 allensdk.brain_observatory.nwb.nwb_utils (module), 131
 allensdk.brain_observatory.nwb.schemas (module), 132
 allensdk.brain_observatory.observatory_plots (module), 158
 allensdk.brain_observatory.ophys (module), 136
 allensdk.brain_observatory.ophys.trace_extraction (module), 136
 allensdk.brain_observatory.r_neuropil (module), 160
 allensdk.brain_observatory.receptive_field_analysis (module), 144
 allensdk.brain_observatory.receptive_field_analysis (module), 136
 allensdk.brain_observatory.receptive_field_analysis (module), 139
 allensdk.brain_observatory.receptive_field_analysis (module), 140
 allensdk.brain_observatory.receptive_field_analysis (module), 140
 allensdk.brain_observatory.receptive_field_analysis (module), 141
 allensdk.brain_observatory.receptive_field_analysis (module), 141
 allensdk.brain_observatory.receptive_field_analysis (module), 142
 allensdk.brain_observatory.receptive_field_analysis (module), 142
 allensdk.brain_observatory.receptive_field_analysis (module), 143
 allensdk.brain_observatory.roi_masks (module), 161
 allensdk.brain_observatory.running_speed (module), 164
 allensdk.brain_observatory.session_analysis (module), 164
 allensdk.brain_observatory.session_api_utils (module), 167
 allensdk.brain_observatory.static_gratings (module), 168
 allensdk.brain_observatory.stimulus_analysis (module), 169
 allensdk.brain_observatory.stimulus_info (module), 171
 allensdk.brain_observatory.sync_dataset (module), 175
 allensdk.brain_observatory.sync_utilities (module), 145
 allensdk.brain_observatory.visualization (module), 145
 allensdk.config (module), 187
 allensdk.config.app (module), 181
 allensdk.config.app.application_config (module), 179
 allensdk.config.manifest (module), 185
 allensdk.config.manifest_builder (module), 187
 allensdk.config.model (module), 185
 allensdk.config.model.description (module), 183
 allensdk.config.model.description_parser (module), 184
 allensdk.config.model.formats (module), 183
 allensdk.config.model.formats.hdf5_util (module), 182
 allensdk.config.model.formats.json_description_parser (module), 182
 allensdk.config.model.formats.pycfg_description_parser (module), 183
 allensdk.config.model.parameters (module), 187
 allensdk.core.auth_config (module), 188
 allensdk.core.authentication_and_session (module), 188
 allensdk.core.brain_observatory_cache (module), 189
 allensdk.core.data_processing (module), 189
 allensdk.core.brain_observatory_nwb_data_set (module), 194
 allensdk.core.receptive_field (module), 194
 allensdk.core.cache_method_utilities (module), 198
 allensdk.core.cell_types_cache (module), 198
 allensdk.core.is_utilities (module), 201
 allensdk.core.dat_utilities (module), 201
 allensdk.core.exceptions (module), 201
 allensdk.core.h5_utilities (module), 201
 allensdk.core.json_utilities (module), 202
 allensdk.core.lazy_property (module), 188
 allensdk.core.lazy_property.lazy_property (module), 188
 allensdk.core.lazy_property.lazy_property_mixin (module), 188
 allensdk.core.mouse_connectivity_cache (module), 203
 allensdk.core.nwb_data_set (module), 208
 allensdk.core.obj_utilities (module), 210
 allensdk.core.ontology (module), 210
 allensdk.core.ophys_experiment_session_id_mapping (module), 211
 allensdk.core.reference_space (module), 211
 allensdk.core.reference_space_cache (module), 214

[allensdk.core.simple_tree \(module\), 216](#)
[allensdk.core.sitk_utilities \(module\), 219](#)
[allensdk.core.structure_tree \(module\), 220](#)
[allensdk.core.swc \(module\), 222](#)
[allensdk.core.typing \(module\), 227](#)
[allensdk.deprecated \(module\), 338](#)
[allensdk.ephys \(module\), 239](#)
[allensdk.ephys.ephys_extractor \(module\), 227](#)
[allensdk.ephys.ephys_features \(module\), 231](#)
[allensdk.ephys.extract_cell_features \(module\), 238](#)
[allensdk.ephys.feature_extractor \(module\), 239](#)
[allensdk.internal \(module\), 315](#)
[allensdk.internal.api \(module\), 250](#)
[allensdk.internal.api.api_prerelease \(module\), 244](#)
[allensdk.internal.api.behavior_data_lims_api \(module\), 245](#)
[allensdk.internal.api.behavior_lims_api \(module\), 247](#)
[allensdk.internal.api.behavior_ophys_api \(module\), 247](#)
[allensdk.internal.api.lims_api \(module\), 248](#)
[allensdk.internal.api.mtrain_api \(module\), 249](#)
[allensdk.internal.api.ophys_lims_api \(module\), 249](#)
[allensdk.internal.api.queries \(module\), 244](#)
[allensdk.internal.api.queries.biophysical_module \(module\), 240](#)
[allensdk.internal.api.queries.biophysical_module_reader \(module\), 240](#)
[allensdk.internal.api.queries.grid_data_api_pre_release \(module\), 241](#)
[allensdk.internal.api.queries.mouse_connectivity_pre_release \(module\), 242](#)
[allensdk.internal.api.queries.optimize_config_pre_release \(module\), 243](#)
[allensdk.internal.api.queries.pre_release \(module\), 244](#)
[allensdk.internal.brain_observatory \(module\), 265](#)
[allensdk.internal.brain_observatory.annotated_metrics \(module\), 251](#)
[allensdk.internal.brain_observatory.demix_report \(module\), 252](#)
[allensdk.internal.brain_observatory.demixer \(module\), 253](#)
[allensdk.internal.brain_observatory.eye_calibration \(module\), 254](#)
[allensdk.internal.brain_observatory.fit_ellipse \(module\), 257](#)
[allensdk.internal.brain_observatory.frame_stream \(module\), 257](#)
[allensdk.internal.brain_observatory.itracker_utils \(module\), 258](#)
[allensdk.internal.brain_observatory.mask_set \(module\), 260](#)
[allensdk.internal.brain_observatory.ophys_session \(module\), 260](#)
[allensdk.internal.brain_observatory.resources \(module\), 251](#)
[allensdk.internal.brain_observatory.roi_filter_utils \(module\), 261](#)
[allensdk.internal.brain_observatory.time_sync \(module\), 263](#)
[allensdk.internal.core \(module\), 268](#)
[allensdk.internal.core.lims_pipeline_module \(module\), 265](#)
[allensdk.internal.core.lims_utilities \(module\), 265](#)
[allensdk.internal.core.mouse_connectivity_cache_pre_release \(module\), 266](#)
[allensdk.internal.core.simpletree \(module\), 267](#)
[allensdk.internal.core.swc \(module\), 268](#)
[allensdk.internal.ephys \(module\), 274](#)
[allensdk.internal.ephys.core_feature_extractor \(module\), 268](#)
[allensdk.internal.ephys.plot_qc_figures \(module\), 269](#)
[allensdk.internal.ephys.plot_qc_figures3 \(module\), 271](#)
[allensdk.internal.model \(module\), 291](#)
[allensdk.internal.model.AIC \(module\), 289](#)
[allensdk.internal.model.biophysical \(module\), 278](#)
[allensdk.internal.model.biophysical.biophysical_array \(module\), 276](#)
[allensdk.internal.model.biophysical.biophysical_pre_release \(module\), 276](#)
[allensdk.internal.model.biophysical.check_fit_shift \(module\), 276](#)
[allensdk.internal.model.biophysical.deap_utils \(module\), 276](#)
[allensdk.internal.model.biophysical.ephys_utils \(module\), 276](#)
[allensdk.internal.model.biophysical.fits \(module\), 274](#)
[allensdk.internal.model.biophysical.fits_metrics \(module\), 274](#)
[allensdk.internal.model.biophysical.fits.fit_styles \(module\), 274](#)
[allensdk.internal.model.biophysical.make_deap_fit \(module\), 277](#)
[allensdk.internal.model.biophysical.passive_fitting \(module\), 276](#)

allensdk.internal.model.biophysical.passive_fit	(module), 291
(module), 275	allensdk.internal.morphology.morphology
allensdk.internal.model.biophysical.passive_fit	(module), 291
(module), 275	allensdk.internal.morphology.morphvis
allensdk.internal.model.biophysical.passive_fit	(module), 295
(module), 275	allensdk.internal.morphology.node (mod-
allensdk.internal.model.biophysical.passive_fit	ule), 298
(module), 275	allensdk.internal.morphology.validate_swc
allensdk.internal.model.biophysical.passive_fit	(module), 298
(module), 275	allensdk.internal.mouse_connectivity
allensdk.internal.model.biophysical.passive_fit	(module), 298
(module), 274	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.passive_fit	(module), 298
(module), 275	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.run_optimize	(module), 299
(module), 277	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.run_optimize	(module), 299w
(module), 278	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.run_passive	(module), 301
(module), 278	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.run_simulation	(module), 301
(module), 278	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.biophysical.run_simulation	(module), 303w
(module), 278	allensdk.internal.mouse_connectivity.interval_union
allensdk.internal.model.data_access	(module), 302
(module), 290	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif (module), 289	(module), 305
allensdk.internal.model.glif.are_two_lists_of_neurons_the_same	(module), 303
(module), 279	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif.error_function	(module), 303
(module), 279	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif.find_spikes	(module), 304
(module), 279	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif.find_sweeps	(module), 304
(module), 279	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif.glif_experiment_name	(module), 305
(module), 280	allensdk.internal.mouse_connectivity.projection_thru
allensdk.internal.model.glif.glif_optimization	(module), 305
(module), 282	allensdk.internal.mouse_connectivity.tissuecyte_stimuli
allensdk.internal.model.glif.glif_optimization	(module), 307
(module), 282	allensdk.internal.mouse_connectivity.tissuecyte_stimuli
allensdk.internal.model.glif.MLIN (module), 278	(module), 305
allensdk.internal.model.glif.plotting	(module), 306
(module), 285	allensdk.internal.pipeline_modules (module), 315
allensdk.internal.model.glif.rc (module), 286	allensdk.internal.pipeline_modules.gbm
(module), 286	(module), 308
allensdk.internal.model.glif.spike_cutting	(module), 307
(module), 287	allensdk.internal.pipeline_modules.gbm.generate_gbm
allensdk.internal.model.glif.threshold_adapter	(module), 307
(module), 287	allensdk.internal.pipeline_modules.gbm.generate_gbm
allensdk.internal.model.GLM (module), 290	(module), 307
allensdk.internal.morphology (module), 298	allensdk.internal.pipeline_modules.gbm.generate_gbm
allensdk.internal.morphology.compartmental	allensdk.internal.pipeline_modules.gbm.generate_gbm

(module), 308
 allensdk.internal.pipeline_modules.run_anahd_sdk_region_connectivity.grid.subimage
 (module), 308
 allensdk.internal.pipeline_modules.run_dahlsdk.mouse_connectivity.grid.subimage.base_sub
 (module), 308
 allensdk.internal.pipeline_modules.run_diff_sdk.mouse_connectivity.grid.subimage.cav_subin
 (module), 308
 allensdk.internal.pipeline_modules.run_nahp_sdk.mouse_connectivity.grid.subimage.classic_
 (module), 309
 allensdk.internal.pipeline_modules.run_obs_sdk.mouse_connectivity.grid.subimage.count_sub
 (module), 309
 allensdk.internal.pipeline_modules.run_obs_sdk.mouse_connectivity.grid.utilities
 (module), 310
 allensdk.internal.pipeline_modules.run_obs_sdk.mouse_connectivity.grid.utilities.downsamp
 (module), 312
 allensdk.internal.pipeline_modules.run_obs_sdk.mouse_connectivity.grid.utilities.image_ut
 (module), 312
 allensdk.internal.pipeline_modules.run_obs_sdk.mouse_connectivity.grid.writers
 (module), 313
 allensdk.internal.pipeline_modules.run_tasks_sdk.task_initialize_tasks (module), 315
 (module), 315
 allensdk.model (module), 329
 allensdk.model.biophys_sim (module), 317
 allensdk.model.biophys_sim.bps_command
 (module), 316
 allensdk.model.biophys_sim.config (mod-
 ule), 316
 allensdk.model.biophys_sim.neuron (mod-
 ule), 316
 allensdk.model.biophys_sim.neuron.hoc_util_alpha_filter() (in module al-
 lensdk.brain_observatory.r_neuropil), 160
 allensdk.model.biophys_sim.scripts (mod-
 ule), 316
 allensdk.model.biophysical (module), 319
 allensdk.model.biophysical.run_simulate
 (module), 317
 allensdk.model.biophysical.runner (mod-
 ule), 317
 allensdk.model.biophysical.utils (mod-
 ule), 318
 allensdk.model.glif (module), 329
 allensdk.model.glif.glif_neuron (module),
 319
 allensdk.model.glif.glif_neuron_methods
 (module), 323
 allensdk.model.glif.simulate_neuron
 (module), 328
 allensdk.morphology (module), 329
 allensdk.morphology.validate_swc (mod-
 ule), 329
 allensdk.mouse_connectivity (module), 337
 allensdk.mouse_connectivity.grid (mod-
 ule), 337
 allensdk.mouse_connectivity.grid.image_series_attribute (module), 336
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.subimage
 (module), 333
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.subimage.cav_subin
 (module), 331
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.subimage.classic_
 (module), 332
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.subimage.count_sub
 (module), 333
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.utilities
 (module), 335
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.utilities.downsamp
 (module), 333
 allensdk.mouse_connectivity.grid.image_series_attribute.grid.utilities.image_ut
 (module), 334
 allensdk.mouse_connectivity.grid.writers
 (module), 336
 allensdk.test_utilities (module), 338
 allensdk.test_utilities.custom_comparators
 (module), 337
 allensdk.test_utilities.regression_fixture
 (module), 337
 allensdk.test_utilities.temp_dir (mod-
 ule), 337
 allocate_by_vsync() (in module al-
 lensdk.brain_observatory.ecephys.stimulus_sync),
 129
 analog_meta_data (al-
 lensdk.brain_observatory.sync_dataset.Dataset
 attribute), 175
 ANALYSIS_DATA_KEY (al-
 lensdk.core.brain_observatory_cache.BrainObservatoryCache
 attribute), 190
 analyze_trough_details() (in module al-
 lensdk.ephys.ephys_features), 231
 ancestor_ids() (al-
 lensdk.core.simple_tree.SimpleTree method),
 216
 ancestor_ids() (al-
 lensdk.internal.core.simpletree.SimpleTree
 method), 267
 ancestors() (allensdk.core.simple_tree.SimpleTree
 method), 216
 ancestors() (allensdk.internal.core.simpletree.SimpleTree
 method), 267
 angle_lines() (in module al-
 lensdk.brain_observatory.circle_plots), 150
 angular_wheel_rotation (al-
 lensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOneP
 attribute), 116

angular_wheel_velocity (al-
lensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile
attribute), 117

annotate_change_detect() (in module al-
lensdk.brain_observatory.behavior.mtrain),
97

annotate_trials() (in module al-
lensdk.brain_observatory.behavior.mtrain),
97

AnnotatedSectionDataSetsApi (class in al-
lensdk.api.queries.annotated_section_data_sets_api),
43

ANNOTATION_KEY (al-
lensdk.core.reference_space_cache.ReferenceSpaceCache
attribute), 214

Api (class in allensdk.api.api), 72

APICAL_DENDRITE (allensdk.core.swc.Morphology
attribute), 223

APICAL_DENDRITE (al-
lensdk.internal.morphology.morphology.Morphology
attribute), 291

ApiPrerelease (class in al-
lensdk.internal.api.api_prerelease), 244

append() (allensdk.core.swc.Morphology method),
223

append() (allensdk.internal.morphology.morphology.Morphology
method), 291

append() (allensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet
method), 303

append_experiment_metrics() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metadata() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_drifting_grating() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_locally_sparse_noise() (allensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_natural_movie_one() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_natural_movie_three() (allensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_natural_movie_two() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_natural_scene() (al-
lensdk.brain_observatory.session_analysis.SessionAnalysis
method), 164

append_metrics_static_grating() (al-

lensdk.brain_observatory.session_analysis.SessionAnalysis
OnePickleStimFile
append_threshold_components() (al-
lensdk.model.glif.glif_neuron.GlifNeuron
method), 321

append_well_known_file() (in module al-
lensdk.internal.core.lims_utilities), 265

ApplicationConfig (class in al-
lensdk.config.app.application_config), 179

apply() (allensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet
method), 303

apply_affine() (allensdk.core.swc.Morphology
method), 223

apply_affine() (al-
lensdk.internal.morphology.morphology.Morphology
method), 291

apply_affine_only_rotation() (al-
lensdk.internal.morphology.morphology.Morphology
method), 292

apply_average_tile() (al-
lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile
method), 306

apply_average_tile_to_self() (al-
lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile
method), 306

apply_colormap() (in module al-
lensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet
method), 303

apply_configuration_from_command_line() (allensdk.config.app.application_config.ApplicationConfig
method), 179

apply_configuration_from_environment() (allensdk.config.app.application_config.ApplicationConfig
method), 180

apply_configuration_from_file() (al-
lensdk.config.app.application_config.ApplicationConfig
method), 180

apply_display_sequence() (in module al-
lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.PreSpikesTable
method), 331

apply_divisions() (in module al-
lensdk.mouse_connectivity.grid.utilities.downsampling_utilities),
333

apply_frame_times() (in module al-
lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes.PreSpikesTable
method), 331

apply_mask() (allensdk.mouse_connectivity.grid.subimage.base_subimage.BaseSubImage
method), 331

apply_pixel_counter() (al-
lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage
method), 331

APPLY_SSL (allensdk.api.queries.reference_space_api.ReferenceSpaceApi
attribute), 63

archive_cell() (al-

[lensdk.internal.model.biophysical.biophysical_archive.BiophysicalArchiveMorphology attribute](#)), 223
[method](#)), 276
[AXON \(allensdk.internal.morphology.morphology.Morphology attribute\)](#), 291
[are_two_lists_of_arrays_the_same\(\)](#) (in module [al-lensdk.internal.model.glif.are_two_lists_of_arrays_the_same](#)), 279
[arg_parser\(\)](#) (in module [al-lensdk.internal.brain_observatory.demix_report](#)), [lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit](#)), 275
[args \(allensdk.internal.core.lims_pipeline_module.PipelineModule attribute\)](#), 111
[ArgSchemaParserPlus](#) (class in [al-lensdk.brain_observatory.argschema_utilities](#)), 146
[ARRAY \(allensdk.api.queries.connected_services.ConnectedServices attribute\)](#), 223
[as_dataframe\(\)](#) ([allensdk.config.manifest.Manifest method](#)), 185
[as_dataframe\(\)](#) ([al-lensdk.config.manifest_builder.ManifestBuilder method](#)), 187
[as_dict\(\)](#) ([allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor method](#)), 227
[as_dict\(\)](#) ([allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor method](#)), 228
[aspect_ratio \(allensdk.brain_observatory.stimulus_info.Monitor attribute\)](#), 171
[assert_exists\(\)](#) (in module [al-lensdk.internal.pipeline_modules.run_demixing](#)), 308
[assign_session_id\(\)](#) (in module [al-lensdk.brain_observatory.behavior.mtrain](#)), 97
[assign_sweep_values\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spikes](#)), 122
[assign_to_last\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.stimulus_sync](#)), 129
[atlas_image_query\(\)](#) ([al-lensdk.api.queries.image_download_api.ImageDownloadApi method](#)), 54
[autocorr\(\)](#) (in module [al-lensdk.internal.model.glif.MLIN](#)), 278
[average_rate\(\)](#) (in module [al-lensdk.ephys.ephys_features](#)), 232
[AVERAGE_TEMPLATE](#) ([al-lensdk.api.queries.reference_space_api.ReferenceSpaceApi attribute](#)), 63
[average_tile_is_untrimmed\(\)](#) ([al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile_tile method](#)), 306
[average_voltage\(\)](#) (in module [al-lensdk.ephys.ephys_features](#)), 232
[barcode_line \(allensdk.brain_observatory.ecephys.align_timestamps.barcode_line\)](#), 111
[BarcodeSyncDataset](#) (class in [al-lensdk.brain_observatory.ecephys.align_timestamps.barcode_sync](#)), 111
[BASAL_DENDRITE \(allensdk.core.swc.Morphology attribute\)](#), 223
[BASAL_DENDRITE](#) ([al-lensdk.internal.morphology.morphology.Morphology attribute](#)), 291
[base_object_to_eye_rotation_matrix\(\)](#) (in module [al-lensdk.internal.brain_observatory.eye_calibration](#)), 260
[bb_dist\(\)](#) (in module [al-lensdk.internal.brain_observatory.mask_set](#)), 260
[behavior_alignment](#) ([al-lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync attribute](#)), 313
[behavior_delta](#) ([al-lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync attribute](#)), 313
[behavior_session_id](#) ([al-lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession attribute](#)), 88
[behavior_session_id_to_foraging_id\(\)](#) ([allensdk.internal.api.behavior_lims_api.BehaviorLimsApi static method](#)), 247
[behavior_times](#) ([al-lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync attribute](#)), 313
[behavior_tracking_keys](#) ([al-lensdk.brain_observatory.sync_dataset.Dataset attribute](#)), 175
[behavior_video_timestamps](#) ([al-lensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute](#)), 263
[BehaviorBase](#) (class in [al-lensdk.brain_observatory.behavior.internal.behavior_base](#)), 81
[BehaviorDataLimsApi](#) (class in [al-lensdk.internal.api.behavior_data_lims_api](#)), 245
[BehaviorDataSession](#) (class in [al-lensdk.brain_observatory.behavior.behavior_data_session](#)), 129

88
 BehaviorLimsApi (class in al-
 lensdk.internal.api.behavior_lims_api), 247
 BehaviorMetadataSchema (class in al-
 lensdk.brain_observatory.behavior.schemas),
 99
 BehaviorOphysApiBase (class in al-
 lensdk.brain_observatory.behavior.behavior_ophys_api),
 80
 BehaviorOphysBase (class in al-
 lensdk.brain_observatory.behavior.internal.behavior_ophys_base),
 82
 BehaviorOphysLimsApi (class in al-
 lensdk.internal.api.behavior_ophys_api),
 247
 BehaviorTaskParametersSchema (class in al-
 lensdk.brain_observatory.behavior.schemas),
 99
 best_fit_value() (al-
 lensdk.internal.model.biophysical.make_deap_fit_pre_release
 method), 277
 BIC() (in module allensdk.internal.model.AIC), 290
 binarize() (allensdk.mouse_connectivity.grid.subimage.base_subimage
 method), 331
 BinaryIntervalSearchTree (class in al-
 lensdk.brain_observatory.stimulus_info),
 171
 binned_cells_sp (al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 attribute), 169
 binned_cells_vis (al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 attribute), 169
 binned_dx_sp (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 attribute), 169
 binned_dx_vis (al-
 lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 attribute), 169
 BIOPHYSICAL_MODEL_TYPE_IDS (al-
 lensdk.api.queries.biophysical_api.BiophysicalApi
 attribute), 44
 BiophysicalApi (class in al-
 lensdk.api.queries.biophysical_api), 44
 BiophysicalArchiver (class in al-
 lensdk.internal.model.biophysical.biophysical_archiver),
 276
 BiophysicalModuleApi (class in al-
 lensdk.internal.api.queries.biophysical_module_api),
 240
 BiophysicalModuleReader (class in al-
 lensdk.internal.api.queries.biophysical_module_reader),
 240
 blend() (in module al-
 lensdk.internal.mouse_connectivity.projection_thumbnail.visualization_utilities),
 304
 blend_component_from_point() (in module al-
 lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher),
 305
 blend_with_background() (in module al-
 lensdk.internal.mouse_connectivity.projection_thumbnail.generator),
 303
 block_apply() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 334
 block_image() (in module al-
 lensdk.mouse_connectivity.grid.utilities.downsampling_utilities),
 333
 BOOLEAN (allensdk.api.queries.connected_services.ConnectedServices
 attribute), 51
 BrainObservatoryAnalysisException, 147
 BrainObservatoryApi (class in al-
 lensdk.api.queries.brain_observatory_api),
 46
 BrainObservatoryApiPreRelease (class in al-
 lensdk.internal.api.queries.pre_release), 244
 BrainObservatoryCache (class in al-
 lensdk.brain_observatory.brain_observatory_cache), 189
 BrainObservatoryMonitor (class in al-
 lensdk.brain_observatory.stimulus_info),
 171
 BrainObservatoryNwbDataSet (class in al-
 lensdk.core.brain_observatory_nwb_data_set),
 194
 build_affine_transform() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 334
 build_cell_plots() (in module al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 310
 build_coarse_grids() (al-
 lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder
 method), 336
 build_colormap() (in module al-
 lensdk.brain_observatory.ecephys.stimulus_table.visualization_utilities),
 120
 build_composite_transform() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities),
 334
 build_correlation_plots() (in module al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 310
 build_drifting_gratings() (in module al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 310
 build_experiment_thumbnails() (in module al-
 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 310

`build_eye_tracking_plots()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 310
`build_from_image()` (in module `allensdk.internal.mouse_connectivity.projection_thumbnail_image_sheet.ImageSheet` static method), 303
`build_hex_pack()` (in module `allensdk.brain_observatory.circle_plots`), 150
`build_locally_sparse_noise()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 310
`build_manifest()` (`allensdk.api.cache.Cache` method), 76
`build_manifest()` (`allensdk.core.brain_observatory_cache.BrainObservatoryCache` method), 190
`build_manifest()` (`allensdk.core.cell_types_cache.CellTypesCache` method), 199
`build_manifest()` (`allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader` method), 243
`build_natural_movie()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 310
`build_natural_scenes()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 311
`build_plots()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 311
`build_query()` (`allensdk.api.queries.svg_api.SvgApi` method), 70
`build_query_url()` (`allensdk.api.queries.rma_api.RmaApi` method), 65
`build_receptive_field()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 311
`build_reference_aligned_image_channel_volumes_url()` (`allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 58
`build_rma()` (`allensdk.api.queries.biophysical_api.BiophysicalApi` method), 44
`build_rotation_transform()` (`allensdk.internal.mouse_connectivity.projection_thumbnail_volume_projection.VolumeProjection` method), 305
`build_schema_query()` (`allensdk.api.queries.rma_api.RmaApi` method), 66
`build_speed_tuning()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 311
`build_static_gratings()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 287
`build_stimuluswise_table()` (in module `allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spike_table.PreSpikeTable` static method), 122
`build_trial_matrix()` (in module `allensdk.brain_observatory.receptive_field_analysis.chisquarperf`), 136
`build_type()` (in module `allensdk.internal.pipeline_modules.run_observatory_thumbnails`), 311
`build_url()` (`allensdk.api.queries.connected_services.ConnectedServices` method), 51
`build_volumetric_data_download_url()` (`allensdk.api.queries.reference_space_api.ReferenceSpaceApi` method), 63
`burst_metrics()` (`allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor` method), 228
`Cache` (class in `allensdk.api.cache`), 76
`cache_clear()` (`allensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession` method), 88
`cache_clear()` (`allensdk.core.cache_method_utilities.CachedInstanceMethodMixin` method), 198
`cache_csv()` (`allensdk.api.cache.Cache` static method), 76
`cache_csv_dataframe()` (`allensdk.api.cache.Cache` static method), 76
`cache_csv_json()` (`allensdk.api.cache.Cache` static method), 76
`cache_data()` (`allensdk.api.queries.biophysical_api.BiophysicalApi` method), 45
`cache_json()` (`allensdk.api.cache.Cache` static method), 76
`cache_json_dataframe()` (`allensdk.api.cache.Cache` static method), 76
`cache_stimulus_file()` (`allensdk.api.queries.glif_api.GlifApi` method), 52
`cacheable()` (in module `allensdk.api.cache`), 78
`CachedInstanceMethodMixin` (class in `allensdk.core.cache_method_utilities`), 198
`cacher()` (`allensdk.api.cache.Cache` static method), 76
`calc_deriv()` (in module `allensdk.brain_observatory.behavior.running_processing`), 98
`calc_spike_component_of_threshold_from_multiblip()` (in module `allensdk.internal.model.glif.threshold_adaptation`), 287

calc_spike_cut_and_v_reset_via_expvar_residuals()	(in module al-lensdk.internal.model.glif.spike_cutting), 286	calc_subimage()	(class in al-lensdk.mouse_connectivity.grid.subimage.cav_subimage), 331
calculate()	(allensdk.core.lazy_property.lazy_property.LazyProperty method), 188	calculate()	(allensdk.api.queries.reference_space_api.ReferenceSpaceApi attribute), 63
calculate()	(allensdk.internal.mouse_connectivity.internal_unionize_reference_space_api.ReferenceSpaceApi method), 302	calculate()	(allensdk.api.queries.reference_space_api.ReferenceSpaceApi attribute), 63
calculate()	(allensdk.internal.mouse_connectivity.internal_unionize_reference_space_api.ReferenceSpaceApi method), 302	calculate()	(allensdk.api.queries.reference_space_api.ReferenceSpaceApi attribute), 63
calculate()	(allensdk.internal.mouse_connectivity.internal_unionize_record.Unionize method), 302	calculate()	(allensdk.api.queries.reference_space_api.ReferenceSpaceApi attribute), 63
calculate_delay()	(in module al-lensdk.brain_observatory.behavior.sync.process_sync), 84	cell_extractor_for_nwb()	(in module al-lensdk.ephys.ephys_extractor), 231
calculate_dff()	(in module al-lensdk.brain_observatory.dff), 151	cell_features()	(al-lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor method), 227
calculate_dvdt()	(in module al-lensdk.ephys.ephys_features), 233	cell_id()	(allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis attribute), 169
calculate_feature_errors()	(al-lensdk.internal.model.biophysical.deap_utils.Utils method), 276	cell_index_receptive_field_analysis_data()	(allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise attribute), 155
calculate_fi_curves()	(in module al-lensdk.internal.model.biophysical.check_fi_shift), 276	CELL_MAPPING_ID()	(al-lensdk.api.queries.brain_observatory_api.BrainObservatoryApi attribute), 46
calculate_injection_centroid()	(al-lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi method), 58	CELLS_KEY()	(al-lensdk.core.brain_observatory_cache.BrainObservatoryCache attribute), 190
calculate_max_border()	(in module al-lensdk.internal.brain_observatory.roi_filter_utils), 262	CELLS_KEY()	(allensdk.core.cell_types_cache.CellTypesCache attribute), 198
calculate_reward_rate()	(in module al-lensdk.brain_observatory.behavior.trials_processing), 105	celltraces()	(allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis attribute), 169
calculate_roi_and_neuropil_traces()	(in module al-lensdk.brain_observatory.roi_masks), 162	CellTypesApi()	(class in al-lensdk.api.queries.cell_types_api), 49
calculate_scale()	(in module al-lensdk.internal.morphology.morphvis), 295	CellTypesCache()	(class in al-lensdk.core.cell_types_cache), 198
calculate_traces()	(in module al-lensdk.brain_observatory.roi_masks), 162	change_parent()	(allensdk.core.swc.Morphology method), 223
calculate_trough()	(al-lensdk.ephys.feature_extractor.EphysFeatureExtractor method), 239	change_parent()	(al-lensdk.internal.morphology.morphology.Morphology method), 292
call_caching()	(in module al-lensdk.api.caching_utilities), 79	check_and_write()	(al-lensdk.core.reference_space.ReferenceSpace static method), 211
CamStimOnePickleStimFile()	(class in al-lensdk.brain_observatory.ecephys.file_io.stim_file), 116	check_coverage()	(al-lensdk.core.reference_space.ReferenceSpace method), 211
categorize_one_trial()	(in module al-lensdk.brain_observatory.behavior.trials_processing), 105	check_dir()	(allensdk.config.manifest.Manifest method), 185
cav_writer()	(in module al-lensdk.mouse_connectivity.grid.writers), 386	check_nwbfile_version()	(in module al-lensdk.brain_observatory.nwb), 135
		check_org_selections_for_noise_block()	(allensdk.internal.model.biophysical.make_deap_fit_json.Report method), 211

method), 277

check_read_access() (in module *allensdk.brain_observatory.argschema_utilities*), 147

check_stimulus_delay() (in module *allensdk.internal.pipeline_modules.run_ophys_time_sync*), 314

check_thresholds_and_peaks() (in module *allensdk.ephys.ephys_features*), 233

check_write_access() (in module *allensdk.brain_observatory.argschema_utilities*), 147

check_write_access_dir() (in module *allensdk.brain_observatory.argschema_utilities*), 147

check_write_access_overwrite() (in module *allensdk.brain_observatory.argschema_utilities*), 147

checkPreprocess() (in module *allensdk.internal.model.glif.plotting*), 285

checkSpikeCutting() (in module *allensdk.internal.model.glif.plotting*), 285

chi_square_binary() (in module *allensdk.brain_observatory.receptive_field_analysis.chisquare_categorical*), 137

chi_square_within_mask() (in module *allensdk.brain_observatory.receptive_field_analysis.chisquare_categorical*), 137

child_ids() (*allensdk.core.simple_tree.SimpleTree* method), 217

child_ids() (*allensdk.internal.core.simpletree.SimpleTree* method), 267

children() (*allensdk.core.simple_tree.SimpleTree* method), 217

children() (*allensdk.internal.core.simpletree.SimpleTree* method), 267

children_of() (*allensdk.core.swc.Morphology* method), 224

children_of() (*allensdk.internal.morphology.morphology.Morphology* method), 292

chisq_from_stim_table() (in module *allensdk.brain_observatory.chisquare_categorical*), 148

choose_bps_command() (in module *allensdk.model.biophys_sim.bps_command*), 316

choose_inliers() (*allensdk.internal.brain_observatory.fit_ellipse.FitEllipse* method), 257

class_deprecated() (in module *allensdk.deprecated*), 338

classic_writer() (in module *allensdk.mouse_connectivity.grid.writers*), 336

ClassicSubImage (class in *allensdk.mouse_connectivity.grid.subimage.classic_subimage*), 332

clone_structures() (*allensdk.core.structure_tree.StructureTree* static method), 220

cleanup_truncated_file() (*allensdk.api.api.Api* method), 72

clear_updated_params() (*allensdk.brain_observatory.session_api_utils.ParamsMixin* method), 167

clone() (*allensdk.ephys.feature_extractor.EphysFeatures* method), 239

clone() (*allensdk.internal.morphology.morphology.Morphology* method), 292

close() (*allensdk.brain_observatory.sync_dataset.Dataset* method), 175

close() (*allensdk.internal.brain_observatory.frame_stream.CvInputStream* method), 257

close() (*allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream* method), 258

close() (*allensdk.internal.brain_observatory.frame_stream.FfmpegOutputStream* method), 258

close() (*allensdk.internal.brain_observatory.frame_stream.FrameInputStream* method), 258

close() (*allensdk.internal.brain_observatory.frame_stream.FrameOutputStream* method), 258

close() (*allensdk.internal.brain_observatory.mask_set.MaskSet* method), 260

close_sets() (*allensdk.internal.brain_observatory.mask_set.MaskSet* method), 260

collapse_columns() (in module *allensdk.brain_observatory.ecephys.stimulus_table.naming_utilities*), 125

collect_sets() (*allensdk.core.structure_tree.StructureTree* static method), 220

colormap() (in module *allensdk.brain_observatory.behavior.trials_processing*), 105

COLORMAPS (*allensdk.api.queries.image_download_api.ImageDownloadAPI* attribute), 54

Compartment (class in *allensdk.core.swc*), 222

Compartment (class in *allensdk.internal.morphology.compartment*), 291

compartment() (*allensdk.internal.morphology.morphology.Morphology* method), 292

compartment_index (*allensdk.core.swc.Morphology* attribute), 224

compartment_index_by_type() (*allensdk.internal.morphology.morphology.Morphology* method), 292

`lensdk.core.swc.Morphology` method), 224
`compartment_list` (`allensdk.core.swc.Morphology` attribute), 224
`compartment_list` (`allensdk.internal.morphology.morphology.Morphology` attribute), 293
`compartment_list_by_type()` (`allensdk.core.swc.Morphology` method), 224
`CompleteOphysBehaviorMetadataSchema` (class in `allensdk.brain_observatory.behavior.schemas`), 100
`compute()` (`allensdk.brain_observatory.ecephys.align_timestamps.probe_sync_probe_sync` class method), 113
`compute_area()` (`allensdk.internal.brain_observatory.eye_calibration.EyeCalibration` method), 255
`compute_chi()` (in module `allensdk.brain_observatory.chisquare_categorical`), 148
`compute_chi_shuffle()` (in module `allensdk.brain_observatory.chisquare_categorical`), 148
`compute_circular_area()` (in module `allensdk.brain_observatory.behavior.eye_tracking_processing`), 93
`compute_coarse_parameters()` (in module `allensdk.mouse_connectivity.grid.utilities.image_utilities`), 334
`compute_coarse_planes()` (`allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage` method), 331
`compute_coarse_planes()` (`allensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage` method), 331
`compute_coarse_planes()` (`allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage` method), 332
`compute_coarse_planes()` (`allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage` method), 333
`compute_correlations()` (in module `allensdk.internal.brain_observatory.demix_report`), 252
`compute_correlations_without_masks()` (in module `allensdk.internal.brain_observatory.demix_report`), 252
`compute_dff_windowed_median()` (in module `allensdk.brain_observatory.dff`), 152
`compute_dff_windowed_mode()` (in module `allensdk.brain_observatory.dff`), 152
`compute_distance()` (in module `allensdk.brain_observatory.receptive_field_analysis.fit_parameters`), 140
`compute_elliptical_area()` (in module `allensdk.brain_observatory.behavior.eye_tracking_processing`), 93
`compute_expected()` (in module `allensdk.brain_observatory.chisquare_categorical`), 148
`compute_frame_times()` (in module `allensdk.brain_observatory.ecephys.stimulus_sync`), 129
`compute_injection()` (`allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage` method), 332
`compute_injection()` (`allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage` method), 332
`compute_intensity()` (`allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage` method), 332
`compute_non_overlap_masks()` (in module `allensdk.internal.brain_observatory.demix_report`), 252
`compute_non_overlap_traces()` (in module `allensdk.internal.brain_observatory.demix_report`), 252
`compute_observed()` (in module `allensdk.brain_observatory.chisquare_categorical`), 148
`compute_overlap()` (in module `allensdk.brain_observatory.receptive_field_analysis.fit_parameters`), 140
`compute_projection()` (`allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage` method), 332
`compute_projection()` (`allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage` method), 332
`compute_receptive_field()` (in module `allensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis`), 141
`compute_receptive_field_with_postprocessing()` (in module `allensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis`), 141
`compute_sum_pixels()` (`allensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage` method), 332
`compute_sum_pixels()` (`allensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage` method), 333
`Config` (class in `allensdk.model.biophys_sim.config`), 316
`configure_library_method()` (`allensdk.model.glif.glif_neuron.GlifNeuron` method), 316

<code>static method), 321</code>	<code>copy_local()</code> (allensdk.internal.model.biophysical.run_optimize.RunOptimize method), 277
<code>configure_method()</code> (allensdk.model.glif.glif_neuron.GlifNeuron static method), 321	<code>copy_local()</code> (allensdk.internal.model.biophysical.run_simulate_lims.RunSimulateLims method), 278
<code>connect()</code> (in module allensdk.internal.core.lims_utilities), 265	<code>CoronaPlotter</code> (class in allensdk.brain_observatory.circle_plots), 148
<code>ConnectedServices</code> (class in allensdk.api.queries.connected_services), 51	<code>correct_on_off_effects()</code> (in module allensdk.brain_observatory.ecephys.stimulus_sync), 130
<code>consistency_is_key()</code> (in module allensdk.brain_observatory.behavior.criteria), 90	<code>corrected_behavior_video_timestamps</code> (allensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute), 263
<code>consistent_behavior_within_session()</code> (in module allensdk.brain_observatory.behavior.criteria), 90	<code>corrected_eye_video_timestamps</code> (allensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute), 263
<code>construct_well_known_file_download_url()</code> (allensdk.api.api.Api method), 73	<code>corrected_ophys_timestamps</code> (allensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute), 263
<code>consume_volume()</code> (allensdk.mouse_connectivity.grid.image_series_grinder.ImageSeriesGrinder method), 336	<code>corrected_video_timestamps</code> (allensdk.internal.brain_observatory.time_sync.OphysTimeAligner attribute), 263
<code>contingent_trials()</code> (in module allensdk.brain_observatory.behavior.trial_masks), 104	<code>corrected_video_timestamps()</code> (in module allensdk.internal.brain_observatory.time_sync), 264
<code>ContinuousFile</code> (class in allensdk.brain_observatory.ecephys.file_io.continuous_file), 116	<code>create_correlation_report()</code> (in module allensdk.internal.brain_observatory.demix_report), 252
<code>conv()</code> (in module allensdk.mouse_connectivity.grid.utilities.downsampling_utilities), 333	<code>create_manifest()</code> (allensdk.api.queries.rma_api.RmaApi attribute), 65
<code>convert_axis()</code> (in module allensdk.internal.mouse_connectivity.projection_thumbnail_projection_utilities), 304	<code>count</code> (allensdk.internal.brain_observatory.mask_set.MaskSet attribute), 260
<code>convert_discrete_colormap()</code> (in module allensdk.internal.mouse_connectivity.projection_thumbnail_visualization_utilities), 304	<code>count_writer()</code> (in module allensdk.mouse_connectivity.grid.writers), 333
<code>convert_filepath_caseinsensitive()</code> (in module allensdk.brain_observatory.behavior.stimulus_processing), 102	<code>CountSubImage</code> (class in allensdk.mouse_connectivity.grid.subimage.count_subimage), 333
<code>convert_frame()</code> (in module allensdk.internal.pipeline_modules.run_ophys_session_decomposition), 312	<code>create_position_in_mouse_eye_coordinates()</code> (allensdk.internal.brain_observatory.eye_calibration.EyeCalibration static method), 255
<code>convert_from_titan_linux()</code> (in module allensdk.internal.core.lims_utilities), 265	<code>create_session_decomposition_user()</code> (allensdk.config.app.application_config.ApplicationConfig method), 180
<code>convert_type()</code> (allensdk.core.swc.Morphology method), 224	<code>create_basis_IPSP()</code> (in module allensdk.internal.model.GLM), 290
<code>convert_type()</code> (allensdk.internal.morphology.morphology.Morphology method), 293	<code>create_dir()</code> (allensdk.config.manifest.Manifest method), 186
<code>convolve()</code> (in module allensdk.brain_observatory.receptive_field_analysis_utilities), 142	<code>create_extended_trials()</code> (in module allensdk.brain_observatory.behavior.trials_processing), 105
<code>copy()</code> (allensdk.internal.mouse_connectivity.projection_thumbnail_image_to_image_conversion_utilities), 304	<code>create_eye_gaze_mapping_dataframe()</code> (in module allensdk.brain_observatory.nwb), 135
	<code>create_image_sheet_processing_module()</code> (in module allensdk.brain_observatory.nwb), 135

[135](#)
`create_fake_metadata()` (in module `allensdk.internal.pipeline_modules.run_ophys_session_deconvolution`), 312
[312](#)
`create_gaze_mapping_nwb_processing_modules()` (in module `allensdk.brain_observatory.nwb`), [135](#)
`create_gene_fpkm_table()` (in module `allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap`), [307](#)
`create_genes_for_transcripts()` (in module `allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap`), [307](#)
`create_image()` (in module `allensdk.internal.morphology.morphvis`), [296](#)
`create_images()` (`allensdk.internal.brain_observatory.frame_stream.FfmpegInputStream` static method), [258](#)
`create_images()` (`allensdk.internal.brain_observatory.frame_stream.FrameInputStream` static method), [258](#)
`create_manifest()` (`allensdk.api.queries.biophysical_api.BiophysicalApi` static method), [45](#)
`create_neuropil_mask()` (in module `allensdk.brain_observatory.roi_masks`), [162](#)
`create_pynwb_extension_from_schemas()` (in module `allensdk.brain_observatory.nwb.metadata`), [131](#)
`create_region_mask()` (in module `allensdk.internal.brain_observatory.annotated_region_metrics`), [251](#)
`create_roi_mask()` (in module `allensdk.brain_observatory.roi_masks`), [163](#)
`create_roi_mask_array()` (in module `allensdk.brain_observatory.roi_masks`), [163](#)
`create_sample_metadata()` (in module `allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap`), [307](#)
`create_stim_table()` (in module `allensdk.brain_observatory.ecephys.stimulus_table.ecephys_pre_spikes`), [123](#)
`create_stimulus_presentation_time_interval()` (in module `allensdk.brain_observatory.nwb`), [135](#)
`create_transcript_fpkm_table()` (in module `allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap`), [307](#)
`create_transcripts_for_genes()` (in module `allensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap`), [200](#)

[307](#)
`create_utils()` (in module `allensdk.internal.biophysical_utils`), [319](#)
`credential_injector()` (in module `allensdk.core.authentication`), [189](#)
`CredentialProvider` (class in `allensdk.core.authentication`), [188](#)
`CRITERIA` (`allensdk.api.queries.rma_api.RmaApi` attribute), [65](#)
`CRITERIA()` (in module `allensdk.internal.brain_observatory.roi_filter_utils`), [261](#)
`create_cache()` (`allensdk.api.cache.Cache` static method), [76](#)
`CUT_DENDRITE` (`allensdk.core.swc.Marker` attribute), [223](#)
`CUT_DENDRITE` (`allensdk.internal.core.swc.Marker` attribute), [268](#)
`CvInputStream` (class in `allensdk.internal.brain_observatory.frame_stream`), [257](#)

D

`Data` (`allensdk.brain_observatory.behavior.image_api.Image` attribute), [96](#)
`DATA_MASK` (`allensdk.api.queries.grid_data_api.GridDataApi` attribute), [52](#)
`DATA_MASK_KEY` (`allensdk.core.mouse_connectivity_cache.MouseConnectivityCache` attribute), [204](#)
`data_to_ticks()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [105](#)
`data_to_metadata()` (in module `allensdk.brain_observatory.behavior.trials_processing`), [105](#)
`data_transforms()` (`allensdk.core.structure_tree.StructureTree` static method), [220](#)
`data_query()` (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` static method), [46](#)
`data_query_string()` (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` static method), [46](#)
`DataFrameIndexError`, [201](#)
`DataFrameKeyError`, [201](#)
`dataset` (`allensdk.internal.brain_observatory.time_sync.OphysTimeAlign` attribute), [263](#)
`DatUtilities` (class in `allensdk.core.dat_utilities`), [175](#)

DbCredentials (class in *allensdk.core.authentication*), 188

dbname (*allensdk.core.authentication.DbCredentials* attribute), 188

DEBUG (*allensdk.api.queries.rma_api.RmaApi* attribute), 65

debug () (in module *allensdk.internal.pipeline_modules.run_annotated_region_metrics*), 308

debug () (in module *allensdk.internal.pipeline_modules.run_demixing*), 308

debug () (in module *allensdk.internal.pipeline_modules.run_neuropil_correction*), 309

debug () (in module *allensdk.internal.pipeline_modules.run_observatory_analysis*), 309

debug () (in module *allensdk.internal.pipeline_modules.run_observatory_demixing*), 311

debug () (in module *allensdk.internal.pipeline_modules.run_ophys_eye_calibration*), 312

debug () (in module *allensdk.internal.pipeline_modules.run_ophys_session_decomposition*), 312

debug_clause () (*allensdk.api.queries.rma_api.RmaApi* method), 66

debug_plot () (in module *allensdk.internal.pipeline_modules.run_neuropil_correction*), 309

decision_function () (*allensdk.internal.brain_observatory.roi_filter_utils.TrainingLabelClassifier* method), 261

decode_bytes () (in module *allensdk.core.h5_utilities*), 201

default () (*allensdk.brain_observatory.JSONEncoder* method), 179

default_api_url (*allensdk.api.api.Api* attribute), 73

default_argument_parser () (in module *allensdk.internal.core.lims_pipeline_module*), 265

default_ray () (in module *allensdk.internal.brain_observatory.itracker_utils*), 258

default_structure_ids (*allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* attribute), 204

DEFAULT_STRUCTURE_SET_IDS (*allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* attribute), 204

DEFORMATION_FIELD_HEADER_KEY (*allensdk.brain_observatory.behavior.image_api.ImageApi*

lensdk.core.mouse_connectivity_cache.MouseConnectivityCache attribute), 204

DEFORMATION_FIELD_VOXEL_KEY (*allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* attribute), 204

deg_to_dist () (in module *allensdk.brain_observatory.behavior.running_processing*), 283

deinterpolate_RF () (in module *allensdk.brain_observatory.receptive_field_analysis.chisquarperf*), 137

delay_metrics () (*allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor* method), 228

delete_tree () (*allensdk.core.swc.Morphology* method), 224

deserialize () (*allensdk.internal.morphology.morphology.Morphology* method), 293

demix_dep_masks () (in module *allensdk.brain_observatory.demixer*), 150

demix_time_dep_masks () (in module *allensdk.internal.brain_observatory.demixer*), 253

DENDRITE (*allensdk.core.swc.Morphology* attribute), 223

DENSITY (*allensdk.api.queries.grid_data_api.GridDataApi* attribute), 53

deprecated () (in module *allensdk.deprecated*), 338

DEPRECATED_KEYS (*allensdk.brain_observatory.sync_dataset.Dataset* attribute), 175

DEPRECATED_SPIKE_TIMES (*allensdk.core.nwb_data_set.NwbDataSet* attribute), 338

descendant_ids () (*allensdk.core.simple_tree.SimpleTree* method), 217

descendant_ids () (*allensdk.internal.core.simpletree.SimpleTree* method), 267

descendants () (*allensdk.core.simple_tree.SimpleTree* method), 217

descendants () (*allensdk.internal.core.simpletree.SimpleTree* method), 267

Description (class in *allensdk.config.model.description*), 183

DescriptionParser (class in *allensdk.config.model.description_parser*), 184

deserialize () (*allensdk.brain_observatory.behavior.image_api.ImageApi*

static method), 96

detect_bursts() (in module *allensdk.ephys.ephys_features*), 233

detect_duplicates() (in module *allensdk.internal.brain_observatory.mask_set.MaskSet* method), 260

detect_events() (in module *allensdk.brain_observatory.receptive_field_analysis.event_detection*), 139

detect_pauses() (in module *allensdk.ephys.ephys_features*), 234

detect_putative_spikes() (in module *allensdk.ephys.ephys_features*), 234

detect_unions() (in module *allensdk.internal.brain_observatory.mask_set.MaskSet* method), 260

determine_likely_blinks() (in module *allensdk.brain_observatory.behavior.eye_tracking_processing*), 93

determine_outliers() (in module *allensdk.brain_observatory.behavior.eye_tracking_processing*), 94

DEVMOUSE_2012 (in module *allensdk.api.queries.reference_space_api.ReferenceSpaceApi* attribute), 63

DEVMOUSE_ATLAS_PRODUCTS (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi* attribute), 57

df_columns(*allensdk.config.manifest_builder.ManifestBuilder* attribute), 187

dfftraces(*allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis* attribute), 169

DFMFLD_RESOLUTIONS (in module *allensdk.core.mouse_connectivity_cache.MouseConnectivityCache* attribute), 204

dict_generator() (in module *allensdk.brain_observatory.receptive_field_analysis.tools*), 142

dict_to_indexed_array() (in module *allensdk.brain_observatory*), 179

dim_color() (*allensdk.brain_observatory.observatory_plots.DimensionPatchHandler* method), 158

DimensionPatchHandler (class in *allensdk.brain_observatory.observatory_plots*), 158

DIR (*allensdk.config.manifest.Manifest* attribute), 185

DIR_CCW (*allensdk.brain_observatory.circle_plots.PolarPlotter* attribute), 149

DIR_CW (*allensdk.brain_observatory.circle_plots.PolarPlotter* attribute), 149

direct_sum_projection_pixels (in module *allensdk.internal.mouse_connectivity.interval_unionize.tissue_detection* attribute), 301

direct_unionize() (in module *allensdk.internal.mouse_connectivity.interval_unionize.interval_unionize* method), 299

direct_voxel_counts() (in module *allensdk.core.reference_space.ReferenceSpace* method), 212

direct_voxel_map (in module *allensdk.core.reference_space.ReferenceSpace* attribute), 212

DIRNAME (*allensdk.config.manifest.Manifest* attribute), 185

distance() (*allensdk.internal.brain_observatory.mask_set.MaskSet* method), 260

do_blur() (in module *allensdk.internal.mouse_connectivity.projection_thumbnail.generator*), 303

do_query() (*allensdk.api.api.Api* method), 73

do_rma_query() (*allensdk.api.api.Api* method), 73

download_alignment3d() (in module *allensdk.api.queries.grid_data_api.GridDataApi* method), 53

download_annotation_volume() (in module *allensdk.api.queries.reference_space_api.ReferenceSpaceApi* method), 63

download_atlas_image() (in module *allensdk.api.queries.image_download_api.ImageDownloadApi* method), 55

download_data_mask() (in module *allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi* method), 58

download_data_mask() (in module *allensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease* method), 242

download_deformation_field() (in module *allensdk.api.queries.grid_data_api.GridDataApi* method), 53

download_expression_density() (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi* method), 57

download_expression_energy() (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi* method), 57

download_expression_grid_data() (in module *allensdk.api.queries.grid_data_api.GridDataApi* method), 53

download_expression_intensity() (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi* method), 57

download_gene_expression_grid_data() (in module *allensdk.api.queries.grid_data_api.GridDataApi* method), 53

download_image() (in module *allensdk.api.queries.image_download_api.ImageDownloadApi* method), 55

download_injection_density() (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi* method), 57

`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` (class in `al-`
`method`), 58 `lensdk.api.mouse_connectivity.grid.utilities.downsampling_utilities`,
334
`download_injection_density()` (alias of `draw_density_hist()` (in module `al-`
`lensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease`), 296
`method`), 242 `draw_morphology()` (in module `al-`
`lensdk.internal.morphology.morphvis`), 297
`download_injection_fraction()` (alias of `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method`), 58 `Gratings` (class in `al-`
`lensdk.brain_observatory.drifting_gratings`),
153
`download_injection_fraction()` (alias of `lensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease`
`method`), 242 `lensdk.brain_observatory.ecephys.stimulus_table.naming_utilities`
`download_mouse_atlas_volume()` (alias of `lensdk.api.queries.reference_space_api.ReferenceSpaceApi`
`method`), 64 `sync_dataset.Dataset`
`method`), 175
`download_projection_density()` (alias of `dxcn()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` attribute), 169
`method`), 59 `dxtime()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`
attribute), 169
`download_projection_density()` (alias of `lensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease`
`method`), 242 `GlifNeuron`
`method`), 321
`download_projection_grid_data()` (alias of `dynamics_ACurrent_exp()` (in module `al-`
`lensdk.api.queries.grid_data_api.GridDataApi` `lensdk.model.glif.glif_neuron_methods`), 323
`method`), 54 `dynamics_ACurrent_none()` (in module `al-`
`lensdk.model.glif.glif_neuron_methods`), 323
`download_projection_grid_data()` (alias of `lensdk.internal.api.queries.grid_data_api_prerelease.GridDataApiPrerelease`
`method`), 242 `inf()` (in module `al-`
`lensdk.model.glif.glif_neuron_methods`), 323
`download_projection_image()` (alias of `dynamics_threshold_spike_component()`
`lensdk.api.queries.image_download_api.ImageDownloadApi` in module `al-`
`method`), 56 `lensdk.model.glif.glif_neuron_methods`),
324
`download_reference_aligned_image_channel_volumes()` (alias of `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method`), 59 `threshold_three_components_exact()`
(in module `al-`
`lensdk.model.glif.glif_neuron_methods`),
324
`download_section_image()` (alias of `lensdk.api.queries.image_download_api.ImageDownloadApi`
`method`), 56 `dynamics_voltage_linear_exact()` (in mod-
ule `allensdk.model.glif.glif_neuron_methods`),
325
`download_structure_mask()` (alias of `lensdk.api.queries.reference_space_api.ReferenceSpaceApi`
`method`), 64 `dynamics_voltage_linear_forward_euler()`
(in module `al-`
`lensdk.api.queries.reference_space_api.ReferenceSpaceApi` `lensdk.model.glif.glif_neuron_methods`),
325
`download_structure_mesh()` (alias of `lensdk.api.queries.reference_space_api.ReferenceSpaceApi` `lensdk.model.glif.glif_neuron_methods`),
`method`), 64
`download_svg()` (alias of `lensdk.api.queries.svg_api.SvgApi` `method`), 70
`download_template_volume()` (alias of `lensdk.internal.brain_observatory.annotated_region_metrics`),
251
`lensdk.api.queries.reference_space_api.ReferenceSpaceApi` `method`), 64
`eccentricity()` (in module `al-`
`lensdk.internal.brain_observatory.itracker_utils`),
258
`download_url()` (`allensdk.api.api.Api` attribute), 74
`download_volumetric_data()` (alias of `lensdk.api.queries.reference_space_api.ReferenceSpaceApi`
`method`), 64 `ProjectApi` (class in `al-`
`lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_`
114
`downsample()` (`allensdk.core.reference_space.ReferenceSpace` `method`), 212
`EcephysSessionApi` (class in `al-`
`lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_`
114
`downsample_average()` (in module `al-`

- 115
 EcephysSyncDataset (class in *allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset*), 160
- 116
 ellipse_angle_of_rotation() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 257
- ellipse_angle_of_rotation2() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 257
- ellipse_axis_length() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 257
- ellipse_center() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 257
- enable_console_log() (in module *allensdk.config*), 187
- ENERGY (*allensdk.api.queries.grid_data_api.GridDataApi* attribute), 53
- EnvCredentialProvider (class in *allensdk.core.authentication*), 189
- EPHYS_DATA_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 198
- EPHYS_FEATURES_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 198
- EPHYS_SWEEPS_KEY (*allensdk.core.cell_types_cache.CellTypesCache* attribute), 199
- EphysCellFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 227
- EphysFeatureExtractor (class in *allensdk.ephys.feature_extractor*), 239
- EphysFeatures (class in *allensdk.ephys.feature_extractor*), 239
- EphysSweepFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 227
- EphysSweepSetFeatureExtractor (class in *allensdk.ephys.ephys_extractor*), 230
- EpochSeparationException, 147
- EQ (*allensdk.api.queries.rma_api.RmaApi* attribute), 65
- error_calc() (in module *allensdk.brain_observatory.r_neuropil*), 160
- error_calc_outlier() (in module *allensdk.brain_observatory.r_neuropil*), 160
- escape_char (*allensdk.internal.model.biophysical.passive_fitting.dynapark_graph_builder.GraphBuilder* attribute), 275
- estimate_adjusted_detection_parameters() (in module *allensdk.ephys.ephys_features*), 235
- estimate_contamination_ratios() (in module *allensdk.brain_observatory.r_neuropil*), 160
- estimate_error() (*allensdk.brain_observatory.r_neuropil.NeuropilSubtract* method), 160
- estimate_fit_shift() (in module *allensdk.internal.model.biophysical.check_fit_shift*), 276
- estimate_frame_duration() (in module *allensdk.brain_observatory.ecephys.stimulus_sync*), 130
- estimate_sag() (*allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor* method), 228
- estimate_time_constant() (*allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor* method), 228
- euclidean_distance() (in module *allensdk.internal.morphology.node*), 298
- evaluate() (*allensdk.internal.model.glif.glif_optimizer.GlifOptimizer* method), 282
- EVENTS_DATA_KEY (*allensdk.core.brain_observatory_cache.BrainObservatoryCache* attribute), 190
- events_to_pvalues_no_fdr_correction() (in module *allensdk.brain_observatory.receptive_field_analysis.receptive_field*), 141
- EXCEPT (*allensdk.api.queries.rma_api.RmaApi* attribute), 65
- EXCPT (*allensdk.api.queries.rma_api.RmaApi* attribute), 65
- exp_curve() (in module *allensdk.internal.ephys.plot_qc_figures*), 269
- exp_curve() (in module *allensdk.internal.ephys.plot_qc_figures3*), 271
- exp_decay() (in module *allensdk.internal.model.glif.MLIN*), 278
- exp_fit_c() (in module *allensdk.internal.model.glif.threshold_adaptation*), 287
- exp_force_c() (in module *allensdk.internal.model.glif.threshold_adaptation*), 287
- EXPERIMENT_CONTAINERS_KEY (*allensdk.core.brain_observatory_cache.BrainObservatoryCache* attribute), 190
- experiment_correlation_search() (*allensdk.brain_observatory_connectivity_api.MouseConnectivityApi* method), 59
- EXPERIMENT_DATA_KEY (*allensdk.core.brain_observatory_cache.BrainObservatoryCache* attribute), 190
- experiment_id (*allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync* attribute), 160

`attribute`), 313
`experiment_injection_coordinate_search()` (in module `al-`
`lensdk.brain_observatory.ecephys.align_timestamps.barcode`),
`(allensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method)`, 59
`experiment_source_search()` (al- `lensdk.ephys.extract_cell_features`), 238
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method)`, 60
`experiment_spatial_search()` (al- `lensdk.brain_observatory.ecephys.stimulus_table.stimulus_param`
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method)`, 60
`ExperimentGeometry` (class in al- `lensdk.internal.mouse_connectivity.interval_unionize.interval_un`
`lensdk.brain_observatory.stimulus_info`),
171
`EXPERIMENTS_KEY` (al- `lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_t`
`lensdk.core.brain_observatory_cache.BrainObservatoryCache`
`attribute`), 190
`EXPERIMENTS_KEY` (al- `lensdk.internal.ephys.core_feature_extract`),
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache`
`attribute`), 204
`EXPERIMENTS_PRERELEASE_KEY` (al- `lensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E`
`lensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease`
`attribute`), 266
`explode_response_window()` (in module al- `(allensdk.brain_observatory.ecephys.file_io.ecephys_sync_data`
`lensdk.brain_observatory.behavior.mtrain`),
97
`export_frame_to_hdf5()` (in module al- `lensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E`
`lensdk.internal.brain_observatory.ophys_session_decomposition`
`method)`, 116
`export_itksnap_labels()` (al- `lensdk.brain_observatory.nwb.metadata`),
`lensdk.core.reference_space.ReferenceSpace`
`method)`, 212
`export_label_description()` (al- `(allensdk.mouse_connectivity.grid.subimage.base_subimage.Segm`
`lensdk.core.structure_tree.StructureTree`
`method)`, 220
`expsymm_cdf()` (in module al- `lensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E`
`lensdk.internal.model.glif.MLIN`), 278
`expsymm_pdf()` (in module al- `lensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E`
`lensdk.internal.model.glif.MLIN`), 278
`ExtendedTrialSchema` (class in al- `lensdk.brain_observatory.ecephys.align_timestamps.channel_stat`
`lensdk.brain_observatory.behavior.mtrain`),
96
`extract()` (`allensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector.VolumeProjector`
`method)`, 305
`extract()` (in module al- (in module al-
`lensdk.mouse_connectivity.grid.utilities.downsampling_utilities`
334
`extract_barcodes()` (al- `extract_sweep_features()` (in module al-
`lensdk.brain_observatory.ecephys.align_timestamps.barcode`
`method)`, 111
`extract_barcodes_from_states()` (in module al- `lensdk.ephys.extract_barcode_effects`
`(in module al-`
`lensdk.brain_observatory.ecephys.align_timestamps.channel_statistics`), 155
`extract_barcodes_from_times()` `extralength(allensdk.brain_observatory.natural_scenes.NaturalScenes`
`attribute)`, 157

extralength (allensdk.brain_observatory.static_gratings.StaticGratings attribute), 168

extrapolate_model_spike_from_endpoints (FfmpegOutputStream (class in allensdk.internal.brain_observatory.frame_stream), (in module al- 257
lensdk.internal.model.glif.glif_optimizer_neuron), 258
284 figure_in_px (in module al-
extrapolate_model_spike_from_endpoints_single (allensdk.brain_observatory.observatory_plots), (in module al- 158
lensdk.internal.model.glif.glif_optimizer_neuron), FILE (allensdk.config.manifest.Manifest attribute), 185
284 FILE_METADATA_MAPPING (al-
extrapolate_spike_time (in module al- lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet attribute), 194
284 lensdk.internal.model.glif.glif_optimizer_neuron), fill_sweep_responses (al-
extrapolate_spike_voltage (in module al- lensdk.core.nwb_data_set.NwbDataSet
284 lensdk.internal.model.glif.glif_optimizer_neuron), method), 208
filter (allensdk.api.queries.rma_api.RmaApi
eye_alignment (al- method), 66
lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs (in module al-
attribute), 313 lensdk.internal.brain_observatory.itracker_utils,
eye_delta (allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs 258
attribute), 313 filter_cell_specimens (al-
EYE_GAZE_DATA_KEY (al- lensdk.api.queries.brain_observatory_api.BrainObservatoryApi
lensdk.core.brain_observatory_cache.BrainObservatoryCache method), 46
attribute), 190 filter_cells (al-
eye_times (allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutputs 258
attribute), 313 filter_cell_types (allensdk.api.queries.cell_types_api.CellTypesApi
method), 49
eye_tracking_data_is_valid (in module al- filter_cells_api (al-
lensdk.brain_observatory.nwb), 135 lensdk.api.queries.cell_types_api.CellTypesApi
method), 49
EYE_TRACKING_KEYS (al- filter_digital (in module al-
lensdk.brain_observatory.sync_dataset.Dataset 84
attribute), 175 lensdk.brain_observatory.behavior.sync.process_sync),
eye_video_timestamps (al- 84
lensdk.internal.brain_observatory.time_sync.OphysTimeAligner experiment_containers (al-
attribute), 263 lensdk.api.queries.brain_observatory_api.BrainObservatoryApi
method), 47
EyeCalibration (class in al- filter_experiments (al-
lensdk.internal.brain_observatory.eye_calibration), 254 lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
method), 204
filter_experiments (al-
F filter_experiments (al-
factory (allensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.EcephysSyncDataset class method), 116
lensdk.core.mouse_connectivity_cache.prerelease.MouseConnectivityCache
method), 267
factory (allensdk.brain_observatory.ecephys.file_io.stim_file_cache.StimFileCache class method), 117
lensdk.api.queries.brain_observatory_api.BrainObservatoryApi
method), 47
FALSE (allensdk.api.queries.rma_api.RmaApi attribute), 65
filter_nodes (al-
FanPlotter (class in al- lensdk.core.simple_tree.SimpleTree method), 217
lensdk.brain_observatory.circle_plots), 149
FeatureError, 231
fetchall (allensdk.internal.api.PostgresQueryMixin filter_ophys_experiments (al-
method), 250 lensdk.api.queries.brain_observatory_api.BrainObservatoryApi
method), 47
fetchone (allensdk.internal.api.PostgresQueryMixin filter_putative_spikes (in module al-
method), 250 lensdk.ephys.ephys_features), 235
ff (in module allensdk.internal.model.GLM), 290
filter_structure_unionizes (al-
FfmpegInputStream (class in al- lensdk.core.mouse_connectivity_cache.MouseConnectivityCache

method), 204

filter_sweeps() (in module *allensdk.internal.ephys.core_feature_extract*), 268

filtered_sweep_numbers() (in module *allensdk.internal.ephys.core_feature_extract*), 268

filters() (*allensdk.api.queries.rma_api.RmaApi* method), 66

finalize() (*allensdk.brain_observatory.circle_plots.PolarPlotter* method), 149

finalize_no_axes() (in module *allensdk.brain_observatory.observatory_plots*), 158

finalize_no_labels() (in module *allensdk.brain_observatory.observatory_plots*), 158

finalize_with_axes() (in module *allensdk.brain_observatory.observatory_plots*), 158

find() (*allensdk.core.swc.Morphology* method), 224

find() (*allensdk.internal.morphology.morphology.Morphology* method), 293

find_bin_center() (in module *allensdk.internal.model.glif.MLIN*), 278

find_coarse_long_square_amp_delta() (in module *allensdk.internal.ephys.core_feature_extract*), 268

find_container_tags() (in module *allensdk.api.queries.brain_observatory_api*), 49

find_downstroke_indexes() (in module *allensdk.ephys.ephys_features*), 235

find_experiment_acquisition_age() (in module *allensdk.api.queries.brain_observatory_api*), 49

find_first_model_spike() (in module *allensdk.internal.model.glif.glif_optimizer_neuron*), 285

find_licks() (in module *allensdk.brain_observatory.behavior.trials_processing*), 105

find_long_square_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_matching_index() (in module *allensdk.brain_observatory.ecephys.align_timestamps_barcode*), 109

find_negative_baselines() (in module *allensdk.brain_observatory.demixer*), 151

find_negative_baselines() (in module *allensdk.internal.brain_observatory.demixer*), 253

find_negative_transients_threshold() (in module *allensdk.brain_observatory.demixer*), 151

find_negative_transients_threshold() (in module *allensdk.internal.brain_observatory.demixer*), 253

find_noise_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_peak_indexes() (in module *allensdk.ephys.ephys_features*), 236

find_ramp_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_ramp_to_rheo_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_ranked_sweep() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_short_square_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 279

find_specimen_cre_line() (in module *allensdk.api.queries.brain_observatory_api*), 49

find_specimen_reporter_line() (in module *allensdk.api.queries.brain_observatory_api*), 49

find_specimen_transgenic_lines() (in module *allensdk.api.queries.brain_observatory_api*), 49

find_spikes_list() (in module *allensdk.internal.model.glif.find_spikes*), 279

find_spikes_list_old() (in module *allensdk.internal.model.glif.find_spikes*), 279

find_spikes_old() (in module *allensdk.internal.model.glif.find_spikes*), 279

find_spikes_ssq_list() (in module *allensdk.internal.model.glif.find_spikes*), 279

find_stim_start() (in module *allensdk.internal.ephys.core_feature_extract*), 268

find_sweep_stim_start() (in module *allensdk.internal.ephys.core_feature_extract*), 269

find_sweeps() (in module *allensdk.internal.model.glif.find_sweeps*), 280

find_time_index() (in module *allensdk.ephys.ephys_features*), 236

find_upstroke_indexes() (in module *allensdk.ephys.ephys_features*), 236

find_widths() (in module *allensdk.ephys.ephys_features*), 237

find_zero_baselines() (in module *allensdk.ephys.ephys_features*), 237

<code>lensdk.brain_observatory.demixer()</code> , 151	<code>class method)</code> , 305
<code>find_zero_baselines()</code> (in module <code>al-</code> <code>lensdk.internal.brain_observatory.demixer</code>), 253	<code>flag_unexpected_edges()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.ecephys.stimulus_sync</code>), 130
<code>findlevel()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.findlevel</code>), 155	<code>FLOAT (allensdk.api.queries.connected_services.ConnectedServices</code> <code>attribute)</code> , 51
<code>fit()</code> (<code>allensdk.brain_observatory.r_neuropil.NeuropilSubtract</code> <code>method</code>), 160	<code>float_label()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.observatory_plots</code>), 158
<code>fit_avoltage_bvoltage()</code> (in module <code>al-</code> <code>lensdk.internal.model.glif.threshold_adaptation</code>), 287	<code>for_drifting_gratings()</code> (<code>al-</code> <code>lensdk.brain_observatory.circle_plots.FanPlotter</code> <code>static method</code>), 149
<code>fit_avoltage_bvoltage_th()</code> (in module <code>al-</code> <code>lensdk.internal.model.glif.threshold_adaptation</code>), 288	<code>for_static_gratings()</code> (<code>al-</code> <code>lensdk.brain_observatory.circle_plots.FanPlotter</code> <code>static method</code>), 149
<code>fit_block_coordinate_desc()</code> (<code>al-</code> <code>lensdk.brain_observatory.r_neuropil.NeuropilSubtract</code> <code>method</code>), 160	<code>fractaging_id_to_behavior_session_id()</code> (<code>allensdk.internal.api.behavior_lims_api.BehaviorLimsApi</code> <code>static method</code>), 247
<code>fit_ellipse()</code> (<code>al-</code> <code>lensdk.internal.brain_observatory.fit_ellipse.FitEllipse</code> <code>method</code>), 257	<code>FRAME_KEYS (allensdk.brain_observatory.sync_dataset.Dataset</code> <code>attribute)</code> , 175
<code>fit_ellipse()</code> (in module <code>al-</code> <code>lensdk.internal.brain_observatory.fit_ellipse</code>), 257	<code>frame_time_offset()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.sync</code>), 84
<code>fit_fi_slope()</code> (in module <code>al-</code> <code>lensdk.ephys.ephys_extractor</code>), 231	<code>FrameInputStream</code> (class in <code>al-</code> <code>lensdk.internal.brain_observatory.frame_stream</code>), 258
<code>fit_membrane_time_constant()</code> (in module <code>al-</code> <code>lensdk.ephys.ephys_features</code>), 237	<code>FrameOutputStream</code> (class in <code>al-</code> <code>lensdk.internal.brain_observatory.frame_stream</code>), 258
<code>fit_parameters_file_entries()</code> (<code>al-</code> <code>lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader</code> <code>method</code>), 240	<code>frames_per_second</code> (<code>al-</code> <code>lensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOneP</code> <code>distributedModuleReader</code> <code>attribute</code>), 176
<code>fit_parameters_path()</code> (<code>al-</code> <code>lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader</code> <code>method</code>), 240	<code>frequency()</code> (<code>allensdk.brain_observatory.sync_dataset.Dataset</code> <code>method</code>), 176
<code>fit_prespike_time_constant()</code> (in module <code>al-</code> <code>lensdk.ephys.ephys_features</code>), 237	<code>FriendlyDate</code> (class in <code>al-</code> <code>lensdk.brain_observatory.behavior.mtrain</code>), 97
<code>FitEllipse</code> (class in <code>al-</code> <code>lensdk.internal.brain_observatory.fit_ellipse</code>), 257	<code>FriendlyDateTime</code> (class in <code>al-</code> <code>lensdk.brain_observatory.behavior.mtrain</code>), 97
<code>fitgaussian2D()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.receptive_field_analysis.fitgaussian2D</code>), 140	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.drifting_gratings.DriftingGratings</code> <code>static method</code>), 153
<code>fix_array_dimensions()</code> (in module <code>al-</code> <code>lensdk.core.sitk_utilities</code>), 219	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise</code> <code>static method</code>), 155
<code>fix_change_time()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.behavior.mtrain</code>), 98	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.natural_movie.NaturalMovie</code> <code>static method</code>), 156
<code>fix_unary_sections()</code> (<code>al-</code> <code>lensdk.config.model.description.Description</code> <code>method</code>), 183	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.natural_scenes.NaturalScenes</code> <code>static method</code>), 157
<code>fix_unexpected_edges()</code> (in module <code>al-</code> <code>lensdk.brain_observatory.ecephys.stimulus_sync</code>), 130	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.natural_scenes.NaturalScenes</code> <code>static method</code>), 157
<code>fixed_factory()</code> (<code>al-</code> <code>lensdk.internal.mouse_connectivity.projection_thumbnail.volume_projection</code> <code>method</code>), 157	<code>from_analysis_file()</code> (<code>al-</code> <code>lensdk.brain_observatory.natural_scenes.NaturalScenes</code> <code>static method</code>), 157

static method), 168

from_dataframe() (allensdk.config.manifest_builder.ManifestBuilder method), 187

from_df() (allensdk.brain_observatory.stimulus_info.BehaviorDataSessionStimulusInfo static method), 171

from_dict() (allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron class method), 282

from_dict() (allensdk.internal.morphology.node.Node class method), 298

from_dict() (allensdk.model.glif.glif_neuron.GlifNeuron class method), 321

from_dict_legacy() (allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron class method), 283

from_file_name() (allensdk.internal.api.queries.grid_data_api_prerelease.GridDataApiPrerelease class method), 242

from_foraging_id() (allensdk.internal.api.behavior_lims_api.BehaviorLimsApi class method), 247

from_json_file() (allensdk.config.app.application_config.ApplicationConfig method), 180

from_json_string() (allensdk.config.app.application_config.ApplicationConfig method), 180

from_lims() (allensdk.brain_observatory.behavior.behavior_data_session_stimulus_info.ExperimentGeometry class method), 88

from_nwb_path() (allensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession class method), 88

from_nwbfile() (allensdk.brain_observatory.nwb.nwb_api.NwbApi class method), 131

from_path() (allensdk.brain_observatory.nwb.nwb_api.NwbApi class method), 131

from_sweeps() (allensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor class method), 230

static method), 277

generate_manifest_lims() (allensdk.internal.model.biophysical.run_simulate_lims.RunSimulate method), 278

generate_manifest_rma() (allensdk.internal.model.biophysical.run_optimize.RunOptimize method), 277

generate_manifest_rma() (allensdk.internal.model.biophysical.run_simulate_lims.RunSimulate method), 278

generate_morphology() (allensdk.internal.model.biophysical.deap_utils.Utils method), 276

generate_morphology() (allensdk.model.biophysical.utils.AllActiveUtils method), 318

generate_grid_data_api_prerelease() (allensdk.model.biophysical.utils.Utils method), 318

generate_output_cell_features() (in module allensdk.internal.ephys.core.feature_extract), 269

generate_rays() (in module allensdk.internal.brain_observatory.itracker_utils), 259

generate_warp_coordinates() (allensdk.brain_observatory.receptive_field_analysis.utilities), 171

get_A() (in module allensdk.brain_observatory.receptive_field_analysis.utilities), 142

get_A_blur() (in module allensdk.brain_observatory.receptive_field_analysis.utilities), 142

get_affine_parameters() (allensdk.core.mouse_connectivity_cache.MouseConnectivityCache method), 205

get_age() (allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi method), 245

get_age() (allensdk.internal.api.ophys_lims_api.OphysLimsApi method), 249

get_alignment_array() (in module allensdk.internal.brain_observatory.time_sync), 264

get_all_bits() (allensdk.brain_observatory.sync_dataset.Dataset method), 176

get_all_cre_lines() (allensdk.core.brain_observatory_cache.BrainObservatoryCache method), 190

get_all_events() (allensdk.brain_observatory.sync_dataset.Dataset method), 176

gather_from_seeds() (allensdk.internal.model.biophysical.make_deap_fit_json.Report method), 277

gaussian2D() (in module allensdk.brain_observatory.receptive_field_analysis.fit_gaussian2d), 140

GaussianFitError, 140

generate_fit_file() (allensdk.internal.model.biophysical.make_deap_fit_json.Report method), 277

generate_manifest_lims() (allensdk.internal.model.biophysical.run_optimize.RunOptimize method), 176

<code>get_all_features()</code> <i>lensdk.core.cell_types_cache.CellTypesCache</i> method), 199	(al- <code>get_attribute_dict()</code> (in module <i>al- lensdk.brain_observatory.receptive_field_analysis.utilities</i>), 142
<code>get_all_imaging_depths()</code> <i>lensdk.core.brain_observatory_cache.BrainObservatoryCache</i> method), 190	(al- <code>get_average_intensity_projection_image_file()</code> <i>allensdk.internal.api.ophys_lims_api.OphysLimsApi</i> method), 249
<code>get_all_reporter_lines()</code> <i>lensdk.core.brain_observatory_cache.BrainObservatoryCache</i> method), 190	(al- <code>get_average_projection()</code> (al- <i>lensdk.brain_observatory.behavior.behavior_ophys_api.Behavior</i> method), 80
<code>get_all_session_types()</code> <i>lensdk.core.brain_observatory_cache.BrainObservatoryCache</i> method), 190	(al- <code>get_average_projection()</code> (al- <i>lensdk.brain_observatory.behavior.internal.behavior_ophys_base</i> method), 83
<code>get_all_stimuli()</code> <i>lensdk.core.brain_observatory_cache.BrainObservatoryCache</i> method), 190	(al- <code>get_average_projection()</code> (al- <i>lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi</i> method), 247
<code>get_all_targeted_structures()</code> <i>lensdk.core.brain_observatory_cache.BrainObservatoryCache</i> method), 190	(al- <code>get_barcode_table()</code> (al- <i>lensdk.brain_observatory.ecephys.align_timestamps.barcode_syn</i> method), 112
<code>get_all_times()</code> <i>lensdk.brain_observatory.sync_dataset.Dataset</i> method), 176	(al- <code>get_behavior_experiment_id()</code> (al- <i>lensdk.internal.api.behavior_lims_api.BehaviorLimsApi</i> method), 247
<code>get_analog_channel()</code> <i>lensdk.brain_observatory.sync_dataset.Dataset</i> method), 176	(al- <code>get_behavior_monitoring()</code> (in module <i>al- lensdk.brain_observatory.behavior.sync</i>), 85
<code>get_analog_meta()</code> <i>lensdk.brain_observatory.sync_dataset.Dataset</i> method), 176	(al- <code>get_behavior_session_id()</code> (al- <i>lensdk.brain_observatory.behavior.internal.behavior_base.Behav</i> method), 81
<code>get_ancestor_id_map()</code> <i>lensdk.core.structure_tree.StructureTree</i> method), 220	(al- <code>get_behavior_session_id()</code> (al- <i>lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp</i> method), 245
<code>get_ancestor_id_map()</code> (in module <i>al- lensdk.internal.mouse_connectivity.interval_unionize.run_tissu</i> 301	(al- <code>get_behavior_session_id()</code> (al- <i>lensdk.internal.api.ophys_lims_api.OphysLimsApi</i> method), 249
<code>get_annotated_section_data_sets()</code> <i>lensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi</i> method), 43	(al- <code>get_behavior_session_uuid()</code> (al- <i>lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp</i> method), 245
<code>get_annotated_section_data_sets_via_rma()</code> <i>lensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi</i> method), 44	(al- <code>get_behavior_session_uuid()</code> (al- <i>lensdk.brain_observatory.behavior.ophys_api.BehaviorOphysLimsApi</i> method), 247
<code>get_annotation_volume()</code> <i>lensdk.core.reference_space_cache.ReferenceSpaceCache</i> method), 214	(al- <code>get_behavior_stimulus_file()</code> (al- <i>lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp</i> method), 245
<code>get_api_list_by_container_id()</code> <i>lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi</i> class method), 247	(al- <code>get_behavior_stimulus_file()</code> (al- <i>lensdk.internal.api.behavior_lims_api.BehaviorLimsApi</i> method), 247
<code>get_atlases()</code> <i>lensdk.api.queries.ontologies_api.OntologiesApi</i> method), 61	(al- <code>get_behavior_stimulus_file()</code> (al- <i>lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi</i> method), 247
<code>get_atlases_table()</code> <i>lensdk.api.queries.ontologies_api.OntologiesApi</i> method), 61	(al- <code>get_behavior_tracking_video_filepath_df()</code> <i>(allensdk.internal.api.lims_api.LimsApi</i> method), 248
<code>get_attribute_dict()</code> (in module <i>al- lensdk.brain_observatory.receptive_field_analysis.receptive_fie</i> 142	(al- <code>get_behavior_training_df()</code> (al- <i>lensdk.internal.api.mtrain_api.MtrainApi</i> method), 249

`get_behavior_training_df()` (al- `lensdk.internal.api.mtrain_api.MtrainSqlApi` method), 249
`get_birth_date()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` method), 245
`get_bit()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 176
`get_bit()` (in module `al- lensdk.brain_observatory.sync_dataset`), 178
`get_bit_changes()` (al- `lensdk.brain_observatory.sync_dataset.Dataset` method), 176
`get_blend()` (in module `al- lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitching`), 305
`get_blend_component()` (in module `al- lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitching`), 306
`get_blocks()` (in module `al- lensdk.brain_observatory.ecephys.stimulus_table.visualization.blocks`), 120
`get_cache_path()` (`allensdk.api.cache.Cache` method), 76
`get_cap_check_indices()` (in module `al- lensdk.internal.model.biophysical.passive_fitting.preprocess`), 275
`get_catch_responses()` (in module `al- lensdk.brain_observatory.behavior.dprime`), 92
`get_cav_density()` (in module `al- lensdk.internal.mouse_connectivity.interval_unionize.data_utilities`), 299
`get_cell()` (`allensdk.api.queries.cell_types_api.CellTypesApi` method), 49
`get_cell_metrics()` (al- `lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 47
`get_cell_metrics()` (al- `lensdk.internal.api.queries.pre_release.BrainObservatoryApiPreRelease` method), 244
`get_cell_roi_ids()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 249
`get_cell_specimen_id_mapping()` (al- `lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 47
`get_cell_specimen_ids()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 194
`get_cell_specimen_indices()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 194
`get_cell_specimen_table()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi` method), 80
`get_cell_specimen_table()` (al- `lensdk.brain_observatory.behavior.internal.behavior_ophys_base_api.BehaviorOphysBaseApi` method), 83
`get_cell_specimen_table()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 249
`get_cell_specimens()` (al- `lensdk.core.brain_observatory_cache.BrainObservatoryCache` method), 190
`get_cells()` (`allensdk.core.cell_types_cache.CellTypesCache` method), 199
`get_cells()` (`allensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver` method), 276
`get_change_time_frame_response_latency()` (`lensdk.brain_observatory.behavior.trials_processing`), 105
`get_channels()` (al- `lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api` method), 114
`get_channels()` (al- `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api` method), 115
`get_child_ids()` (`allensdk.core.ontology.Ontology` method), 210
`get_children()` (`allensdk.core.ontology.Ontology` method), 211
`get_colormap()` (al- `lensdk.core.structure_tree.StructureTree` method), 221
`get_column_definitions()` (al- `lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 47
`get_column_name()` (in module `al- lensdk.brain_observatory.nwb.nwb_utils`), 131
`get_components()` (in module `al- lensdk.brain_observatory.receptive_field_analysis.utilities`), 142
`get_compound_annotated_section_data_sets()` (`allensdk.api.queries.annotated_section_data_sets_api.AnnotatedSectionDataSetsApi` method), 44
`get_config()` (`allensdk.config.manifest_builder.ManifestBuilder` method), 187
`get_connection()` (al- `lensdk.internal.api.PostgresQueryMixin` method), 247
`get_containers_df()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` method), 247
`get_corrected_fluorescence_traces()` (al-

`lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase`
`method`), 80

`get_corrected_fluorescence_traces()` (al- `lensdk.brain_observatory.r_neuropil`), 161
`lensdk.brain_observatory.behavior.internal.behavior_ophys_base.BehaviorOphysBase` (al-
`method`), 83

`get_corrected_fluorescence_traces()` (al- `lensdk.internal.morphology.morphology.Morphology`
`method`), 293

`get_corrected_fluorescence_traces()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (in module al-
`method`), 194

`get_corrected_fluorescence_traces()` (al- `lensdk.brain_observatory.receptive_field_analysis.chisquarperf`),
137

`get_credential_provider()` (in module al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysApiBase` (in module al-
`method`), 247

`get_driver_line()` (in module al- `lensdk.brain_observatory.behavior.dprime`),
92

`get_current_stage()` (al- `lensdk.core.authentication`), 189

`get_driver_line()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`
`method`), 245

`get_cursor()` (al- `lensdk.internal.api.mtrain_api.MtrainApi` `method`), 249

`get_driver_line()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi`
`method`), 249

`get_data_mask()` (al- `lensdk.internal.api.PostgresQueryMixin` `method`), 250

`get_ecephys_session_id()` (al- `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` (al-
`method`), 205

`get_default_manifest_file()` (in module al- `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_`
`method`), 115

`get_edges()` (al- `lensdk.brain_observatory.sync_dataset.Dataset`
`method`), 176

`get_deformation_field()` (al- `lensdk.api.cache`), 79

`get_ecephys_data()` (al- `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` (al-
`method`), 205

`get_ecephys_features()` (al- `lensdk.core.cell_types_cache.CellTypesCache` (al-
`method`), 199

`get_demix_file()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` `method`), 249

`get_ecephys_features()` (al- `lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 50

`get_demixed_traces()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` `method`), 249

`get_ecephys_features()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (al-
`method`), 194

`get_ecephys_features()` (al- `lensdk.core.cell_types_cache.CellTypesCache` (al-
`method`), 199

`get_descendant_ids()` (al- `lensdk.core.ontology.Ontology` `method`), 211

`get_ecephys_sweeps()` (al- `lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 50

`get_descendants()` (al- `lensdk.core.ontology.Ontology` `method`), 211

`get_ecephys_sweeps()` (al- `lensdk.api.queries.glif_api.GlifApi` `method`),
52

`get_df()` (al- `lensdk.internal.api.mtrain_api.MtrainApi` `method`), 249

`get_ecephys_sweeps()` (al- `lensdk.core.cell_types_cache.CellTypesCache`
`method`), 200

`get_dff_file()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` `method`), 249

`get_epoch_mask_list()` (in module al- `lensdk.core.brain_observatory_nwb_data_set`),
198

`get_dff_traces()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase` (al-
`method`), 80

`get_ophys_data()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` `method`), 249

`get_dff_traces()` (al- `lensdk.brain_observatory.behavior.internal.behavior_ophys_base.BehaviorOphysBase` (in module al-
`method`), 83

`get_dff_traces()` (al- `lensdk.brain_observatory.behavior.trials_processing`),
105

`get_dff_traces()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` (al-
`method`), 194

`get_dff_traces()` (al- `lensdk.brain_observatory.sync_dataset.Dataset`
`method`), 177

`get_dff_traces()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysApiBase` (al-
`method`), 80

`get_ophys_data_by_line()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` `method`), 249

`lensdk.brain_observatory.sync_dataset.Dataset` 309
`method`), 177 `get_experiment_structure_unionizes()`
`get_events_per_pixel()` (in module `al-` `(allensdk.core.mouse_connectivity_cache.MouseConnectivityCache)`), 206
`lensdk.brain_observatory.receptive_field_analysis.chisquare_method`), 206
138 `get_experiment_sweep_numbers()` (al-
`get_excluded()` (al- `lensdk.core.nwb_data_set.NwbDataSet` (al-
`lensdk.internal.brain_observatory.roi_filter_utils.TrainingModelClassifier` (al-
`method`), 262 `get_experiments()` (al-
`get_expected_events_by_pixel()` `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
(in module `al-` `method`), 61
`lensdk.brain_observatory.receptive_field_analysis.get_experiments()` (al-
138 `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache`
`get_experiment_analysis_file()` `method`), 206
(in module `al-` `get_experiments()` (al-
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` `lensdk.internal.api.queries.mouse_connectivity_api_prerelease.MouseConnectivityApiPrerelease`
311 `method`), 242
`get_experiment_container_id()` (al- `get_experiments()` (al-
`lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` `lensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease`
`method`), 248 `method`), 267
`get_experiment_container_metrics()` (al- `get_experiments_api()` (al-
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`
`method`), 47 `method`), 61
`get_experiment_containers()` (al- `get_extended_trials()` (al-
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.internal.api.behavior_lims_api.BehaviorLimsApi`
`method`), 48 `method`), 247
`get_experiment_containers()` (al- `get_extended_trials()` (al-
`lensdk.core.brain_observatory_cache.BrainObservatoryCache` `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`
`method`), 191 `method`), 248
`get_experiment_containers()` (al- `get_extended_trials()` (in module `al-`
`lensdk.internal.api.queries.pre_release.BrainObservatoryApiPrerelease` `lensdk.brain_observatory.behavior.trials_processing`),
`method`), 244 105
`get_experiment_date()` (al- `get_external_specimen_name()` (al-
`lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`
`method`), 245 `method`), 245
`get_experiment_date()` (al- `get_external_specimen_name()` (al-
`lensdk.internal.api.ophys_lims_api.OphysLimsApi` `lensdk.internal.api.ophys_lims_api.OphysLimsApi`
`method`), 249 `method`), 249
`get_experiment_detail()` (al- `get_eye_tracking()` (al-
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`
`method`), 61 `method`), 248
`get_experiment_files()` (in module `al-` `get_eye_tracking()` (in module `al-`
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` `lensdk.brain_observatory.behavior.sync`),
311 85
`get_experiment_id()` (al- `get_eye_tracking_data()` (al-
`lensdk.internal.api.lims_api.LimsApi` `method`), `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysLimsApi`
248 `method`), 80
`get_experiment_nwb_file()` (in module `al-` `get_eye_tracking_filepath()` (al-
`lensdk.internal.pipeline_modules.run_observatory_analysis` `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`
309 `method`), 248
`get_experiment_nwb_file()` (in module `al-` `get_eye_tracking_video_filepath_df()`
`lensdk.internal.pipeline_modules.run_observatory_thumbnail` `lensdk.internal.api.lims_api.LimsApi`
311 `method`), 248
`get_experiment_session()` (in module `al-` `get_falling_edges()` (al-
`lensdk.internal.pipeline_modules.run_observatory_analysis` `lensdk.brain_observatory.sync_dataset.Dataset`

method), 177

get_false_alarm_rate() (in module *allensdk.brain_observatory.behavior.dprime*), 92

get_features() (in module *allensdk.internal.ephys.plot_qc_figures*), 269

get_features() (in module *allensdk.internal.ephys.plot_qc_figures3*), 271

get_field_of_view_shape() (in module *allensdk.internal.api.ophys_lims_api.OphysLimsApi*), 249

get_fluorescence() (in module *allensdk.brain_observatory.stimulus_analysis.stimulus_analysis*), 169

get_fluorescence_timestamps() (in module *allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet*), 195

get_fluorescence_traces() (in module *allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet*), 195

get_foraging_id() (in module *allensdk.internal.api.ophys_lims_api.OphysLimsApi*), 249

get_format() (in module *allensdk.config.manifest.Manifest*), 186

get_full_genotype() (in module *allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi*), 245

get_full_genotype() (in module *allensdk.internal.api.ophys_lims_api.OphysLimsApi*), 249

get_gaussian_fit() (in module *allensdk.brain_observatory.receptive_field_analysis.postprocessing*), 141

get_gaussian_fit_single_channel() (in module *allensdk.brain_observatory.receptive_field_analysis.postprocessing*), 140

get_genes() (in module *allensdk.api.queries.mouse_atlas_api.MouseAtlasApi*), 57

get_go_responses() (in module *allensdk.brain_observatory.behavior.dprime*), 92

get_gratings_metadata() (in module *allensdk.brain_observatory.behavior.stimulus_processing*), 102

get_h() (in module *allensdk.internal.model.biophysical.passive_fitting.neuron_utils*), 275

get_hit_rate() (in module *allensdk.brain_observatory.behavior.dprime*), 92

get_id_acronym_map() (in module *allensdk.core.structure_tree.StructureTree*), 221

get_image() (in module *allensdk.brain_observatory.nwb.nwb_api.NwbApi*), 131

get_image_info_from_trial() (in module *allensdk.brain_observatory.behavior.trials_processing*), 105

get_image_region() (in module *allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile*), 306

get_image_to_atlas() (in module *allensdk.api.queries.synchronization_api.SynchronizationApi*), 70

get_image_to_image() (in module *allensdk.api.queries.synchronization_api.SynchronizationApi*), 71

get_images_dict() (in module *allensdk.brain_observatory.behavior.stimulus_processing*), 103

get_imaging_depth() (in module *allensdk.internal.api.ophys_lims_api.OphysLimsApi*), 249

get_imaging_plane_group() (in module *allensdk.internal.api.ophys_lims_api.OphysLimsApi*), 249

get_indicator_bound_point() (in module *allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher*), 306

get_indices_by_distance() (in module *allensdk.internal.brain_observatory.roi_filter_utils*), 262

get_parameters_on_data() (in module *allensdk.internal.mouse_connectivity.interval_unionize.data_utilities*), 299

get_injection_density() (in module *allensdk.core.mouse_connectivity_cache.MouseConnectivityCache*), 206

get_injection_fraction() (in module *allensdk.core.mouse_connectivity_cache.MouseConnectivityCache*), 207

get_input_data() (in module *allensdk.internal.pipeline_modules.run_observatory_thumbnails*), 152

get_input_json() (in module *allensdk.internal.core.lims_utilities*), 265

get_intensity() (in module *allensdk.mouse_connectivity.grid.subimage.base_subimage.Intensity*), 330

`get_invalid_times()` (al- `get_max_projection()` (al-
`lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi` (al-
`method`), 115 `method`), 81
`get_invalid_times()` (al- `get_max_projection()` (al-
`lensdk.brain_observatory.nwb.nwb_api.NwbApi` `lensdk.brain_observatory.behavior.internal.behavior_ophys_base`
`method`), 131 `method`), 83
`get_isi_experiments()` (al- `get_max_projection()` (al-
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryN`
`method`), 48 `method`), 195
`get_isi_experiments()` (al- `get_max_projection()` (al-
`lensdk.brain_observatory.ecephys.ecephys_project_api.EcephysProjectApi` `lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi`
`method`), 114 `method`), 250
`get_isis()` (in `module` al- `get_max_projection_file()` (al-
`lensdk.ephys.ephys_features`), 238 `lensdk.internal.api.ophys_lims_api.OphysLimsApi`
`method`), 250
`get_keys()` (in `module` al-
`lensdk.internal.brain_observatory.time_sync`), `get_mean_response()` (al-
264 `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
`get_lfp()` (allensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi) `get_mean_waveforms()` (al-
`method`), 115 `lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi`
`get_lfp_channel_order()` (al- `lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi`
`lensdk.brain_observatory.ecephys.file_io.continuous_file.ContinuousFile` `method`), 116
`method`), 116
`get_click_times()` (in `module` al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi`
`lensdk.brain_observatory.behavior.sync`), `method`), 81
85 `get_metadata()` (al-
`get_clicks()` (allensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi) `lensdk.brain_observatory.behavior.internal.behavior_ophys_base`
`method`), 80 `method`), 83
`get_clicks()` (allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase) (al-
`method`), 81 `lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi`
`get_clicks()` (allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi) `method`), 195
`method`), 245 `get_metadata()` (al-
`get_clicks()` (allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi) `lensdk.brain_observatory.nwb_data_set.BrainObservatoryN`
`method`), 248 `method`), 195
`get_line()` (allensdk.brain_observatory.sync_dataset.Dataset) `get_metadata()` (al-
`method`), 177 `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`
`get_line_changes()` (al- `method`), 246
`lensdk.brain_observatory.sync_dataset.Dataset` `get_metadata()` (al-
`method`), 177 `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`
`method`), 248
`get_list_of_path_dict()` (in `module` al- `get_metadata()` (al-
`lensdk.test_utilities.regression_fixture`), 337 `lensdk.internal.api.ophys_lims_api.OphysLimsApi`
`get_locally_sparse_noise_stimulus_template()` (allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet) `method`), 195
`method`), 195
`get_metrics()` (in `module` al-
`get_manifest()` (al- `lensdk.internal.brain_observatory.annotated_region_metrics`),
`lensdk.config.manifest_builder.ManifestBuilder` 252
`method`), 187
`get_missing_path()` (al-
`get_manual_injection_summary()` (al- `lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile`
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi`), 307
`method`), 61
`get_morphology_features()` (al-
`get_mask()` (allensdk.brain_observatory.stimulus_info.Monitor) `lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 171 `method`), 50
`get_mask_plane()` (al- `get_morphology_features()` (al-
`lensdk.brain_observatory.roi_masks.Mask` `lensdk.core.cell_types_cache.CellTypesCache`
`method`), 161 `method`), 200

`get_motion_corrected_image_stack_file()` `get_neuronal_models()` (al-
lensdk.internal.api.ophys_lims_api.OphysLimsApi *lensdk.internal.model.biophysical.biophysical_archiver.Biophysical*
method), 250 *method*), 276

`get_motion_correction()` (al- `get_neuronal_models_by_id()` (al-
lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi *lensdk.api.glif_api.GlifApi* *method*),
method), 81 52

`get_motion_correction()` (al- `get_neuropil_r()` (al-
lensdk.brain_observatory.behavior.internal.behavior_ophys_lims_api.BehaviorOphysLimsApi *lensdk.brain_observatory.nwb_data_set.BrainObservatoryNwbData*
method), 83 *method*), 195

`get_motion_correction()` (al- `get_neuropil_traces()` (al-
lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbData *lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbData*
method), 195 *method*), 196

`get_motion_correction()` (al- `get_noise_correlation()` (al-
lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi *lensdk.brain_observatory.drifting_gratings.DriftingGratings*
method), 248 *method*), 153

`get_mouse_id()` (in *module* al- `get_noise_correlation()` (al-
lensdk.brain_observatory.behavior.trials_processing), *lensdk.brain_observatory.natural_scenes.NaturalScenes*
105 *method*), 157

`get_name_map()` (al- `get_noise_correlation()` (al-
lensdk.core.structure_tree.StructureTree *lensdk.brain_observatory.static_gratings.StaticGratings*
method), 221 *method*), 168

`get_natural_movie_template()` (al- `get_nwb_filepath()` (al-
lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi *lensdk.brain_observatory_cache.BrainObservatoryCache*
method), 114 *method*), 191

`get_natural_scene_template()` (al- `get_nwb_filepath()` (al-
lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.EcephysProjectApi *lensdk.brain_observatory_cache.BrainObservatoryCache*
method), 114 *method*), 248

`get_nearest()` (al- `get_nwb_filepath()` (al-
lensdk.brain_observatory.sync_dataset.Dataset *lensdk.internal.api.ophys_lims_api.OphysLimsApi*
method), 177 *method*), 250

`get_neuron_config()` (al- `get_objectlist_file()` (al-
lensdk.api.queries.glif_api.GlifApi *method*), *lensdk.internal.api.ophys_lims_api.OphysLimsApi*
52 *method*), 250

`get_neuron_configs()` (al- `get_ophys_cell_segmentation_run_id()` (al-
lensdk.api.queries.glif_api.GlifApi *method*), (al-
52 *method*), 250

`get_neuronal_model()` (al- `get_ophys_data_length()` (in *module* al-
lensdk.api.queries.glif_api.GlifApi *method*), *lensdk.internal.brain_observatory.time_sync*),
52 264

`get_neuronal_model_runs()` (al- `get_ophys_experiment_analysis()` (al-
lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi *lensdk.brain_observatory_cache.BrainObservatoryCache*
method), 240 *method*), 191

`get_neuronal_model_templates()` (al- `get_ophys_experiment_data()` (al-
lensdk.api.queries.glif_api.GlifApi *method*), *lensdk.core.brain_observatory_cache.BrainObservatoryCache*
52 *method*), 192

`get_neuronal_models()` (al- `get_ophys_experiment_df()` (al-
lensdk.api.queries.biophysical_api.BiophysicalApi *lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi*
method), 45 *static method*), 248

`get_neuronal_models()` (al- `get_ophys_experiment_dir()` (al-
lensdk.api.queries.glif_api.GlifApi *method*), *lensdk.internal.api.ophys_lims_api.OphysLimsApi*
52 *method*), 250

`get_neuronal_models()` (al- `get_ophys_experiment_events()` (al-
lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi *lensdk.brain_observatory_cache.BrainObservatoryCache*
method), 240 *method*), 192

[get_ophys_experiment_id\(\)](#) (al- [get_params\(\)](#) (in module [al-lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase](#) method), 81 [lensdk.brain_observatory.behavior.trials_processing](#), 105)

[get_ophys_experiment_id\(\)](#) (al- [get_passive_fit_data\(\)](#) (in module [al-lensdk.internal.api.ophys_lims_api.OphysLimsApi](#) method), 250 [lensdk.internal.model.biophysical.passive_fitting.preprocess](#), 275)

[get_ophys_experiment_stimuli\(\)](#) (al- [get_path\(\)](#) ([allensdk.config.manifest.Manifest](#) method), 186 [lensdk.core.brain_observatory_cache.BrainObservatoryCache](#) method), 192

[get_ophys_experiments\(\)](#) (al- [get_path\(\)](#) (in module [al-lensdk.internal.pipeline_modules.run_demixing](#)), [lensdk.api.queries.brain_observatory_api.BrainObservatoryApi](#) method), 48

[get_ophys_experiments\(\)](#) (al- [get_peak\(\)](#) ([allensdk.brain_observatory.drifting_gratings.DriftingGratings](#) method), 153 [lensdk.core.brain_observatory_cache.BrainObservatoryCache](#) method), 155

[get_ophys_experiments\(\)](#) (al- [get_peak\(\)](#) ([allensdk.brain_observatory.natural_movie.NaturalMovie](#) method), 156 [lensdk.internal.api.queries.pre_release.BrainObservatoryApiPreRelease](#) method), 244

[get_ophys_frame_rate\(\)](#) (al- [get_peak\(\)](#) ([allensdk.brain_observatory.natural_scenes.NaturalScenes](#) method), 157 [lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi](#) method), 168

[get_ophys_frames\(\)](#) (in module [al-lensdk.brain_observatory.behavior.sync](#), 86 [get_peak\(\)](#) ([allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis](#) method), 169)

[get_ophys_pupil_data\(\)](#) (al- [get_peak_significance\(\)](#) (in module [al-lensdk.brain_observatory.receptive_field_analysis.chisquarerf](#)), [lensdk.core.brain_observatory_cache.BrainObservatoryCache](#) method), 193

[get_ophys_session_id\(\)](#) (al- [get_peaks\(\)](#) (in module [al-lensdk.internal.model.glif.threshold_adaptation](#)), [lensdk.brain_observatory.behavior.internal.behavior_ophys_api.BehaviorOphysBase](#) method), 83

[get_ophys_timestamps\(\)](#) (al- [get_phodiode_events\(\)](#) (in module [al-lensdk.internal.brain_observatory.time_sync](#)), [lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase](#) method), 81

[get_ophys_timestamps\(\)](#) (al- [get_pipeline_version\(\)](#) (al- [lensdk.core.nwb_data_set.NwbDataSet](#) method), 208 [lensdk.brain_observatory.behavior.internal.behavior_ophys_api.BehaviorOphysBase](#) method), 83

[get_ophys_timestamps\(\)](#) (al- [get_plane_group_count\(\)](#) (al- [lensdk.internal.api.ophys_lims_api.OphysLimsApi](#) method), 250 [lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi](#) method), 248

[get_optogenetic_stimulation\(\)](#) (al- [get_polygons\(\)](#) (al- [lensdk.mouse_connectivity.grid.subimage.base_subimage.Polygons](#) method), 300 [lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi](#) method), 115

[get_ori_info_from_trial\(\)](#) (in module [al-lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api](#) method), 114 [lensdk.brain_observatory.behavior.trials_processing](#), 105

[get_output\(\)](#) ([allensdk.internal.mouse_connectivity.projection_threshold_image_shape_image_shape](#) align_timestamps.barcode), 304 [lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api](#), 110

[get_overall_blend\(\)](#) (in module [al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitching](#) method), 114 [lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitching](#), 306

[get_page\(\)](#) ([allensdk.internal.api.mtrain_api.MtrainApi](#) method), 115 [lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api](#) method), 249

[get_params\(\)](#) ([allensdk.brain_observatory.session_api_utils.ParallelMklm](#) internal.mouse_connectivity.interval_unionize.data_utility method), 167 [lensdk.internal.mouse_connectivity.interval_unionize.data_utility](#), 299)

`get_projection_density()` (al- `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` method), 61
`get_projection_image_info()` (al- `lensdk.core.reference_space_cache.ReferenceSpaceCache` method), 215
`get_projection_matrix()` (al- `lensdk.core.mouse_connectivity_cache.MouseConnectivityCache` method), 71
`get_pupil_data()` (al- `lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi` method), 115
`get_pupil_location()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 196
`get_pupil_size()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 196
`get_ramp_stim_characteristics()` (in module `al-
lensdk.ephys.extract_cell_features`), 238
`get_raw_cell_specimen_table_dict()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 250
`get_raw_dff_data()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 250
`get_raw_stimulus_frames()` (in module `al-
lensdk.brain_observatory.behavior.sync`), 86
`get_raw_stimulus_timestamps()` (al- `lensdk.brain_observatory.behavior.internal.behavior_ophys_base.BehaviorOphysBase` method), 83
`get_real_photodiode_events()` (in module `al-
lensdk.internal.brain_observatory.time_sync`), 264
`get_receptive_field()` (al- `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 155
`get_receptive_field_analysis_data()` (al- `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 155
`get_receptive_field_attribute_df()` (al- `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 155
`get_reconstruction()` (al- `lensdk.core.cell_types_cache.CellTypesCache` method), 200
`get_reconstruction_markers()` (al- `lensdk.core.cell_types_cache.CellTypesCache` method), 200
`get_recorded_data()` (al- `lensdk.model.biophysical.utils.Utils` method), 318
`get_reference_aligned_image_channel_volumes()` (al- `lensdk.api.queries.mouse_connectivity_api.MouseConnectivityCache` method), 61
`get_reference_space()` (al- `lensdk.core.reference_space_cache.ReferenceSpaceCache` method), 215
`get_reference_to_image()` (al- `lensdk.api.queries.synchronization_api.SynchronizationApi` method), 71
`get_reporter_line()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` method), 246
`get_reporter_line()` (al- `lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 250
`get_representational_similarity()` (al- `lensdk.brain_observatory.drifting_gratings.DriftingGratings` method), 157
`get_representational_similarity()` (al- `lensdk.brain_observatory.natural_scenes.NaturalScenes` method), 157
`get_representational_similarity()` (al- `lensdk.brain_observatory.static_gratings.StaticGratings` method), 168
`get_response()` (al- `lensdk.brain_observatory.drifting_gratings.DriftingGratings` method), 154
`get_response()` (al- `lensdk.brain_observatory.natural_scenes.NaturalScenes` method), 157
`get_response()` (al- `lensdk.brain_observatory.static_gratings.StaticGratings` method), 168
`get_response()` (al- `lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` method), 170
`get_response_latency()` (in module `al-
lensdk.brain_observatory.behavior.trials_processing`), 105
`get_response_type()` (in module `al-
lensdk.brain_observatory.behavior.trials_processing`), 105
`get_rewards()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi` method), 81
`get_rewards()` (al- `lensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase` method), 81
`get_rewards()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` method), 246
`get_rewards()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` method), 248
`get_rewards()` (in module `al-`

`lensdk.brain_observatory.behavior.rewards_processing`), `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.NwbApi`
 98 `method`), 115
`get_rig_metadata()` (al- `get_running_speed()` (al-
`lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.NwbApi`
`method`), 115 `method`), 131
`get_rig_name()` (al- `get_running_speed()` (al-
`lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet`
`method`), 246 `method`), 197
`get_rig_name()` (al- `get_running_speed()` (al-
`lensdk.internal.api.ophys_lims_api.OphysLimsApi` `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`
`method`), 250 `method`), 246
`get_rigid_motion_transform_file()` (al- `get_running_speed()` (al-
`lensdk.internal.api.ophys_lims_api.OphysLimsApi` `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`
`method`), 250 `method`), 248
`get_rising_edges()` (al- `get_schema()` (`allensdk.api.queries.rma_api.RmaApi`
`lensdk.brain_observatory.sync_dataset.Dataset` `method`), 66
`method`), 177 `get_screen_gaze_data()` (al-
`get_roi_ids()` (al- `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.NwbApi`
`lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` `method`), 196
`method`), 196 `get_section_data_sets()` (al-
`get_roi_mask()` (al- `lensdk.api.queries.mouse_atlas_api.MouseAtlasApi`
`lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` `method`), 196
`method`), 196 `get_section_data_sets_by_product()` (al-
`get_roi_mask_array()` (al- `lensdk.api.queries.image_download_api.ImageDownloadApi`
`lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` `method`), 196
`method`), 196 `get_section_image_ranges()` (al-
`get_rois()` (in module al- `lensdk.api.queries.image_download_api.ImageDownloadApi`
`lensdk.internal.brain_observatory.roi_filter_utils`), `method`), 56
 262 `get_segmentation()` (al-
`get_rolling_dprime()` (in module al- `lensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentImageBase`
`lensdk.brain_observatory.behavior.dprime`), `method`), 331
 92 `get_segmentation_mask_image()` (al-
`get_running_data_df()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi`
`lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase` `method`), 81
`method`), 81 `get_session()` (al-
`get_running_data_df()` (al- `lensdk.internal.api.mtrain_api.MtrainApi`
`lensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase` `method`), 249
`method`), 81 `get_session_data()` (al-
`get_running_data_df()` (al- `lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.NwbApi`
`lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`), 114
`method`), 246 `get_session_start_time()` (al-
`get_running_data_df()` (al- `lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.NwbApi`
`lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi`), 115
`method`), 248 `get_session_type()` (al-
`get_running_df()` (in module al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet`
`lensdk.brain_observatory.behavior.running_processing`), `method`), 197
 98 `get_sessions()` (al-
`get_running_speed()` (al- `lensdk.brain_observatory.ecephys.ecephys_project_api.ecephys_project_api.NwbApi`
`lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase` `method`), 81
`method`), 81 `get_sex()` (`allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi`
`lensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase` `method`), 246
`method`), 82 `method`), 250
`get_running_speed()` (al- `get_shuffle_matrix()` (in module al-

[lensdk.brain_observatory.receptive_field_analysis_utilitiesstim_photodiode\(\)](#) (in module [al-](#)
[142](#) [lensdk.brain_observatory.behavior.sync](#)),
[get_signal_correlation\(\)](#) (al- [86](#)
[lensdk.brain_observatory.drifting_gratings.DriftingGratingsstimulus\(\)](#) (al-
[method](#)), [154](#) [lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryN](#)
[get_signal_correlation\(\)](#) (al- [method](#)), [197](#)
[lensdk.brain_observatory.natural_scenes.NaturalScenesstimulus_attr_changes\(\)](#) (in module [al-](#)
[method](#)), [157](#) [lensdk.brain_observatory.behavior.trials_processing](#)),
[get_signal_correlation\(\)](#) (al- [105](#)
[lensdk.brain_observatory.static_gratings.StaticGratingsstimulus_epoch_table\(\)](#) (al-
[method](#)), [168](#) [lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryN](#)
[get_sitk_image_information\(\)](#) (in module [al-](#) [method](#)), [197](#)
[lensdk.core.sitk_utilities](#)), [219](#) [get_stimulus_file\(\)](#) (al-
[get_slice_image\(\)](#) (al- [lensdk.internal.model.biophysical.biophysical_archiver.Biophysic](#)
[lensdk.core.reference_space.ReferenceSpace](#) [method](#)), [276](#)
[method](#)), [212](#) [get_stimulus_frame_rate\(\)](#) (al-
[get_sparse_noise_epoch_mask_list\(\)](#) [lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp](#)
(a- [method](#)), [246](#)
(in module [al-](#) [lensdk.brain_observatory.receptive_field_analysis_utilitiesstimulus_frame_rate\(\)](#) (al-
[142](#) [lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi](#)
[method](#)), [248](#)
[get_spatial_grating\(\)](#) (in module [al-](#) [get_stimulus_mappings\(\)](#) (al-
[lensdk.brain_observatory.stimulus_info](#)), [lensdk.api.queries.brain_observatory_api.BrainObservatoryApi](#)
[172](#) [method](#)), [48](#)
[get_spatio_temporal_grating\(\)](#) (in module [allensdk.brain_observatory.stimulus_info](#)), [172](#) [get_stimulus_metadata\(\)](#) (in module [al-](#)
[get_speed_tuning\(\)](#) (al- [lensdk.brain_observatory.behavior.stimulus_processing](#)),
[lensdk.brain_observatory.stimulus_analysis.StimulusAnalysisstimulus_metadata\(\)](#) (in module [al-](#)
[method](#)), [170](#) [103](#)
[get_spike_amplitudes\(\)](#) (al- [lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp](#)
[lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi](#)
[method](#)), [115](#) [get_stimulus_name\(\)](#) (al-
[get_spike_times\(\)](#) (al- [lensdk.internal.api.ophys_lims_api.OphysLimsApi](#)
[lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi](#)
[method](#)), [115](#) [get_stimulus_name\(\)](#) (al-
[method](#)), [250](#)
[get_spike_times\(\)](#) (al- [lensdk.brain_observatory.behavior.behavior_ophys_api.Behavior](#)
[lensdk.core.nwb_data_set.NwbDataSet](#) [method](#)), [81](#)
[method](#)), [208](#) [get_stimulus_presentations\(\)](#) (al-
[get_spikes\(\)](#) (in module [al-](#) [lensdk.brain_observatory.behavior.internal.behavior_base.Behav](#)
[lensdk.internal.ephys.plot_qc_figures](#)), [269](#) [method](#)), [82](#)
[get_spikes\(\)](#) (in module [al-](#) [get_stimulus_presentations\(\)](#) (al-
[lensdk.internal.ephys.plot_qc_figures3](#)), [lensdk.brain_observatory.behavior.internal.behavior_ophys_base](#)
[271](#) [method](#)), [84](#)
[get_square_stim_characteristics\(\)](#) (in module [allensdk.ephys.extract_cell_features](#)), [239](#) [get_stimulus_presentations\(\)](#) (al-
[239](#) [lensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_](#)
[method](#)), [115](#)
[get_step_stim_characteristics\(\)](#) (in module [al-](#) [get_stimulus_presentations\(\)](#) (al-
[lensdk.internal.model.biophysical.ephys_utils](#)), [276](#) [lensdk.brain_observatory.nwb.nwb_api.NwbApi](#)
[method](#)), [131](#)
[get_stim_characteristics\(\)](#) (in module [al-](#) [get_stimulus_presentations\(\)](#) (al-
[lensdk.ephys.extract_cell_features](#)), [239](#) [lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsAp](#)
[method](#)), [246](#)
[get_stim_data_length\(\)](#) (in module [al-](#) [get_stimulus_presentations\(\)](#) (al-
[lensdk.internal.brain_observatory.time_sync](#)), [264](#) [lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi](#)
[method](#)), [248](#)

`get_stimulus_presentations()` (in module `allensdk.brain_observatory.behavior.stimulus_processing`), 103
`lensdk.api.queries.ontologies_api.OntologiesApi` method), 62
`get_structure_sets()` (al- `lensdk.core.structure_tree.StructureTree` method), 221
`get_stimulus_rebase_function()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` method), 248
`get_stimulus_rebase_function()` (in module `allensdk.brain_observatory.behavior.sync`), 87
`lensdk.api.queries.synchronization_api.SynchronizationApi` method), 72
`get_stimulus_table()` (al- `get_structure_tree()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 197
`lensdk.core.reference_space_cache.ReferenceSpaceCache` method), 216
`get_stimulus_template()` (al- `get_structure_unionizes()` (al- `lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` method), 197
`lensdk.api.queries.mouse_connectivity_api.MouseConnectivityApi` method), 61
`get_stimulus_templates()` (al- `get_structure_unionizes()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApi` method), 81
`lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase` method), 207
`get_stimulus_templates()` (al- `get_structure_unionizes()` (al- `lensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase` method), 82
`lensdk.api.queries.mouse_connectivity_api_prerelease.M` method), 242
`get_stimulus_templates()` (al- `get_structures()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` method), 246
`lensdk.api.queries.ontologies_api.OntologiesApi` method), 62
`get_stimulus_templates()` (al- `get_structures_by_acronym()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` method), 248
`lensdk.core.structure_tree.StructureTree` method), 221
`get_stimulus_templates()` (in module `al- get_structures_by_id()` (al- `lensdk.brain_observatory.behavior.stimulus_processing`), 103
`lensdk.core.structure_tree.StructureTree` method), 221
`get_stimulus_timestamps()` (al- `get_structures_by_name()` (al- `lensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysApiBase` method), 81
`lensdk.core.structure_tree.StructureTree` method), 221
`get_stimulus_timestamps()` (al- `get_structures_by_set_id()` (al- `lensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase` method), 82
`lensdk.core.structure_tree.StructureTree` method), 221
`get_stimulus_timestamps()` (al- `get_structures_with_sets()` (al- `lensdk.brain_observatory.behavior.internal.behavior_ophys_api.BehaviorOphysLimsApi` method), 84
`lensdk.api.queries.ontologies_api.OntologiesApi` method), 62
`get_stimulus_timestamps()` (al- `get_subjects()` (al- `lensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi` method), 246
`lensdk.internal.api.mtrain_api.MtrainApi` method), 249
`get_stimulus_timestamps()` (al- `get_subjects()` (al- `lensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi` method), 248
`lensdk.internal.api.mtrain_api.MtrainSqlApi` method), 249
`get_structure_graphs()` (al- `get_sum_pixel_intensities()` (in module `al- lensdk.api.queries.ontologies_api.OntologiesApi` method), 62
`lensdk.internal.mouse_connectivity.interval_unionize.data_utilitie` 299
`get_structure_mask()` (al- `get_sum_pixels()` (in module `al- lensdk.core.reference_space_cache.ReferenceSpaceCache` method), 215
`lensdk.internal.mouse_connectivity.interval_unionize.data_utilitie` 299
`get_structure_mesh()` (al- `get_surface_2p_pixel_size_um()` (al- `lensdk.core.reference_space_cache.ReferenceSpaceCache` method), 215
`lensdk.internal.api.ophys_lims_api.OphysLimsApi` method), 250
`get_structure_sets()` (al- `get_svg()` (`allensdk.api.queries.svg_api.SvgApi`

method), 70

get_sweep() (allensdk.core.nwb_data_set.NwbDataSet method), 209

get_sweep_metadata() (allensdk.core.nwb_data_set.NwbDataSet method), 209

get_sweep_numbers() (allensdk.core.nwb_data_set.NwbDataSet method), 209

get_sweep_numbers() (in module allensdk.internal.model.glif.find_sweeps), 280

get_sweep_response() (allensdk.brain_observatory.natural_movie.NaturalMovie method), 156

get_sweep_response() (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis method), 170

get_sweep_v_i_t_from_set() (in module allensdk.internal.model.biophysical.ephys_utils), 276

get_sweeps_by_name() (in module allensdk.internal.model.glif.find_sweeps), 280

get_sweeps_of_type() (in module allensdk.internal.model.biophysical.ephys_utils), 276

get_sync_data() (allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi method), 248

get_sync_data() (in module allensdk.brain_observatory.behavior.sync), 87

get_sync_file() (allensdk.internal.api.ophys_lims_api.OphysLimsApi method), 250

get_sync_licks() (allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi method), 248

get_synchronized_frame_times() (in module allensdk.brain_observatory.sync_utilities), 145

get_targeted_structure() (allensdk.internal.api.ophys_lims_api.OphysLimsApi method), 250

get_task_parameters() (allensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysLimsApi method), 81

get_task_parameters() (allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase method), 82

get_task_parameters() (allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi method), 246

get_task_parameters() (allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi method), 248

get_task_parameters() (in module allensdk.brain_observatory.behavior.metadata_processing), 96

get_template_names() (allensdk.internal.model.biophysical.biophysical_archiver.BiophysicalArchiver method), 276

get_template_volume() (allensdk.core.reference_space_cache.ReferenceSpaceCache method), 216

get_time() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_time_string() (in module allensdk.internal.ephys.plot_qc_figures), 269

get_time_string() (in module allensdk.internal.ephys.plot_qc_figures3), 271

get_tree() (allensdk.api.queries.tree_search_api.TreeSearchApi method), 72

get_trial_count_corrected_false_alarm_rate() (in module allensdk.brain_observatory.behavior.dprime), 92

get_trial_count_corrected_hit_rate() (in module allensdk.brain_observatory.behavior.dprime), 92

get_trial_image_names() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trial_lick_times() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trial_reward_time() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trial_timing() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trials() (allensdk.brain_observatory.behavior.behavior_ophys_api.BehaviorOphysLimsApi method), 81

get_trials() (allensdk.brain_observatory.behavior.internal.behavior_base.BehaviorBase method), 82

get_trials() (allensdk.internal.api.behavior_data_lims_api.BehaviorDataLimsApi method), 247

get_trials() (allensdk.internal.api.behavior_ophys_api.BehaviorOphysLimsApi method), 248

get_trials() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trials_v0() (in module allensdk.brain_observatory.behavior.trials_processing), 106

get_trigger() (in module allensdk.brain_observatory.behavior.trials_processing), 106

`lensdk.brain_observatory.behavior.sync`, 87
`get_unit_analysis_metrics()` (in module `lensdk.brain_observatory.ecephys.ecephys_project_api`), 114
`get_unit_filter_value()` (in module `lensdk.brain_observatory.ecephys`), 130
`get_units()` (in module `lensdk.brain_observatory.ecephys.ecephys_project_api`), 114
`get_units()` (in module `lensdk.brain_observatory.ecephys.ecephys_session_api`), 115
`get_video_length()` (in module `lensdk.internal.brain_observatory.time_sync`), 264
`get_visual_stimuli_df()` (in module `lensdk.brain_observatory.behavior.stimulus_processing`), 104
`get_volume_scale()` (in module `lensdk.internal.mouse_connectivity.interval_unionize.run_tibsstudy_tutorial`), 301
`get_well_known_file_by_name()` (in module `allensdk.internal.core.lims_utilities`), 265
`get_well_known_file_by_type()` (in module `allensdk.internal.core.lims_utilities`), 265
`get_well_known_file_ids()` (in module `allensdk.api.queries.biophysical_api`), 45
`get_well_known_files_by_name()` (in module `allensdk.internal.core.lims_utilities`), 265
`get_well_known_files_by_type()` (in module `allensdk.internal.core.lims_utilities`), 266
`get_wkf()` (in module `lensdk.internal.pipeline_modules.run_ophys_eye_calibration`), 312
`get_workflow_state()` (in module `lensdk.internal.api.ophys_lims_api`), 250
`GLIF_TYPES` (in module `allensdk.api.queries.glif_api`), 52
`GlifApi` (class in `allensdk.api.queries.glif_api`), 52
`GlifBadInitializationException`, 282
`GlifBadResetException`, 319
`GlifExperiment` (class in `lensdk.internal.model.glif.glif_experiment`), 280
`GlifNeuron` (class in `allensdk.model.glif.glif_neuron`), 319
`GlifNeuronException`, 282
`GlifNeuronMethod` (class in `lensdk.model.glif.glif_neuron_methods`), 323
`GlifOptimizer` (class in `lensdk.internal.model.glif.glif_optimizer`), 282
`GlifOptimizerNeuron` (class in `lensdk.internal.model.glif.glif_optimizer_neuron`), 282
`GlifProjectApi` (class in `lensdk.brain_observatory.ecephys.ecephys_project_api`), 114
`GratingToScreen` (class in `lensdk.brain_observatory.stimulus_info`), 171
`grating_to_screen()` (in module `lensdk.brain_observatory.stimulus_info`), 171
`GridDataApi` (class in `lensdk.api.queries.grid_data_api`), 242
`GridDataApiPrerelease` (class in `lensdk.api.queries.grid_data_api_prerelease`), 241
`grid_image_blocks()` (in module `lensdk.mouse_connectivity.grid.utilities.image_utilities`), 334
`GridDataApi` (class in `lensdk.api.queries.grid_data_api`), 52
`GridDataApiPrerelease` (class in `lensdk.api.queries.grid_data_api_prerelease`), 241
`h` (in module `allensdk.model.biophys_sim.neuron.hoc_utils`), 315
`h5_object_matcher_relname_in()` (in module `allensdk.core.h5_utilities`), 201
`handle_output_image()` (in module `lensdk.internal.mouse_connectivity.projection_thumbnail.generator`), 303
`handle_pyramid()` (in module `lensdk.mouse_connectivity.grid.writers`), 336
`hdf5_util` (class in `lensdk.config.model.formats.hdf5_util`), 182
`height` (in module `allensdk.brain_observatory.stimulus_info.monitor`), 171
`hex_pack()` (in module `lensdk.brain_observatory.circle_plots`), 150
`hex_to_rgb()` (in module `allensdk.core.structure_tree.structure_tree`), 222
`HocUtils` (class in `allensdk.model.biophys_sim.neuron.hoc_utils`), 315
`hook()` (in module `allensdk.brain_observatory`), 179
`host` (in module `allensdk.core.authentication.db_credentials`), 188
`HUMAN` (in module `allensdk.api.queries.cell_types_api`), 49
`HUMAN_ORGANISM` (in module `allensdk.api.queries.cell_types_api`), 49

[lensdk.api.queries.mouse_atlas_api.MouseAtlasApi](#) [259](#)
[attribute](#)), 57

I

[identify_valid_masks\(\)](#) (in module [al-lensdk.brain_observatory.demixer](#)), 151
[identify_valid_masks\(\)](#) (in module [al-lensdk.internal.brain_observatory.demixer](#)), 253
[Image](#) (class in [allensdk.brain_observatory.behavior.image_api](#)), 96
[image_from_array\(\)](#) (in module [al-lensdk.mouse_connectivity.grid.utilities.image_utilities](#)), 334
[ImageApi](#) (class in [al-lensdk.brain_observatory.behavior.image_api](#)), 96
[ImageDownloadApi](#) (class in [al-lensdk.api.queries.image_download_api](#)), 54
[ImageOutputStream](#) (class in [al-lensdk.internal.brain_observatory.frame_stream](#)), 258
[ImageSeriesGridder](#) (class in [al-lensdk.mouse_connectivity.grid.image_series_gridder](#)), 336
[ImageSheet](#) (class in [al-lensdk.internal.mouse_connectivity.projection_thumbnail.image_sheet](#)), 303
[INCLUDE](#) ([allensdk.api.queries.rma_api.RmaApi](#) [attribute](#)), 65
[INDETERMINATE](#) ([al-lensdk.core.cell_types_cache.ReporterStatus](#) [attribute](#)), 200
[infer_dims\(\)](#) ([allensdk.brain_observatory.circle_plots.CirclePlots](#) [method](#)), 148
[infer_dims\(\)](#) ([allensdk.brain_observatory.circle_plots.CirclePlots](#) [method](#)), 149
[info\(\)](#) ([allensdk.internal.model.biophysical.run_optimizer.RunOptimizer](#) [method](#)), 277
[init_by_mask\(\)](#) ([al-lensdk.brain_observatory.roi_masks.NeuropilMask](#) [method](#)), 162
[init_by_mask\(\)](#) ([al-lensdk.brain_observatory.roi_masks.RoiMask](#) [method](#)), 162
[init_by_pixels\(\)](#) ([al-lensdk.brain_observatory.roi_masks.Mask](#) [method](#)), 161
[initial_cr_point\(\)](#) (in module [al-lensdk.internal.brain_observatory.itracker_utils](#)), 259
[initial_pupil_point\(\)](#) (in module [al-lensdk.internal.brain_observatory.itracker_utils](#)),

[initialize_coarse_volume\(\)](#) ([al-lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder](#) [method](#)), 336
[initialize_hoc\(\)](#) ([al-lensdk.model.biophys_sim.neuron.hoc_utils.HocUtils](#) [method](#)), 316
[initialize_image\(\)](#) ([al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile](#) [method](#)), 307
[initialize_image\(\)](#) (in module [al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)), 306
[initialize_images\(\)](#) (in module [al-lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)), 306
[initiate_unique_seed\(\)](#) ([al-lensdk.internal.model.glif.glif_optimizer.GlifOptimizer](#) [method](#)), 282
[INJECTION_DENSITY](#) ([al-lensdk.api.queries.grid_data_api.GridDataApi](#) [attribute](#)), 53
[INJECTION_DENSITY_KEY](#) ([al-lensdk.core.mouse_connectivity_cache.MouseConnectivityCache](#) [attribute](#)), 204
[INJECTION_ENERGY](#) ([al-lensdk.api.queries.grid_data_api.GridDataApi](#) [attribute](#)), 53
[INJECTION_FRACTION](#) ([al-lensdk.api.queries.grid_data_api.GridDataApi](#) [attribute](#)), 53
[INJECTION_FRACTION_KEY](#) ([al-lensdk.core.mouse_connectivity_cache.MouseConnectivityCache](#) [attribute](#)), 204
[InjectData\(\)](#) ([allensdk.internal.core.lims_pipeline_module.PipelineModule](#) [method](#)), 265
[input_resistance\(\)](#) (in module [al-lensdk.ephys.ephys_extractor](#)), 231
[InputFile](#) (class in [al-lensdk.brain_observatory.argschema_utilities](#)), 146
[insert_iclamp\(\)](#) ([al-lensdk.internal.model.biophysical.deap_utils.Utils](#) [method](#)), 276
[INTEGER](#) ([allensdk.api.queries.connected_services.ConnectedServices](#) [attribute](#)), 51
[INTENSITY](#) ([allensdk.api.queries.grid_data_api.GridDataApi](#) [attribute](#)), 53
[IntensitySubImage](#) (class in [al-lensdk.mouse_connectivity.grid.subimage.base_subimage](#)), 330
[interlength](#) ([allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise](#) [attribute](#)), 155
[interlength](#) ([allensdk.brain_observatory.natural_scenes.NaturalScenes](#)

[attribute](#)), 157
[interlength](#) ([allensdk.brain_observatory.static_gratings.StaticGratings](#)
[attribute](#)), 168
[interpolate_RF](#) () (in module [al-](#)
[lensdk.brain_observatory.receptive_field_analysis](#)[Lehisquarerf](#)),
139
[interpolate_spike_time](#) () (in module [al-](#)
[lensdk.model.glif.glif_neuron](#)), 322
[interpolate_spike_value](#) () (in module [al-](#)
[lensdk.model.glif.glif_neuron](#)), 322
[interpolate_spike_voltage](#) () (in module [al-](#)
[lensdk.internal.model.glif.glif_optimizer_neuron](#)),
285
[intersection](#) () (al-
[lensdk.internal.brain_observatory.mask_set.MaskSet](#)
[method](#)), 260
[intersection_size](#) () (al-
[lensdk.internal.brain_observatory.mask_set.MaskSet](#)
[method](#)), 260
[IntervalUnionizer](#) (class in al-
[lensdk.internal.mouse_connectivity.interval_unionize_interval_unionizer](#)),
299
[invnl](#) () (in module [allensdk.internal.model.GLM](#)), 290
[IS](#) ([allensdk.api.queries.rma_api.RmaApi](#) attribute), 65
[is_empty](#) () ([allensdk.config.model.description.Description](#)
[method](#)), 183
[is_equal](#) () (in module [al-](#)
[lensdk.brain_observatory.session_api_utils](#)),
167
[is_spike_feature_affected_by_clipping](#) ()
([allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor](#)
[method](#)), 228
[is_well_known_file_type](#) () (al-
[lensdk.api.queries.biophysical_api.BiophysicalApi](#)
[method](#)), 45
[isicv](#) () ([allensdk.ephys.feature_extractor.EphysFeatureExtractor](#)
[method](#)), 239

J

[json_handler](#) () (in module [al-](#)
[lensdk.core.json_utilities](#)), 202
[json_msg_query](#) () ([allensdk.api.api.Api](#) method),
74
[json_remove_keys](#) () ([allensdk.api.cache.Cache](#)
[static method](#)), 77
[json_rename_columns](#) () (al-
[lensdk.api.cache.Cache](#) static method), 77
[JsonComments](#) (class in [allensdk.core.json_utilities](#)),
202
[JsonDescriptionParser](#) (class in al-
[lensdk.config.model.formats.json_description_parser](#)),
182
[JSONEncoder](#) (class in [allensdk.brain_observatory](#)),
179

K

[keyed_locate_h5_objects](#) () (in module [al-](#)
[lensdk.core.h5_utilities](#)), 202
[label_data](#) () ([allensdk.internal.brain_observatory.roi_filter_utils.Train](#)
[method](#)), 262
[latency](#) () (in module [allensdk.ephys.ephys_features](#)),
238
[LazyProperty](#) ([allensdk.core.lazy_property.lazy_property_mixin.LazyPr](#)
[attribute](#)), 188
[LazyProperty](#) (class in al-
[lensdk.core.lazy_property.lazy_property](#)),
188
[LazyPropertyMixin](#) (class in al-
[lensdk.core.lazy_property.lazy_property_mixin](#)),
188
[least_squares_RCEl_calc_tested](#) () (in mod-
ule [allensdk.internal.model.glif.rc](#)), 286
[legacy](#) () (in module [allensdk.deprecated](#)), 338
[legend_artist](#) () (al-
[lensdk.brain_observatory.observatory_plots.DimensionPatchHan](#)
[method](#)), 158
[licks](#) ([allensdk.brain_observatory.behavior.behavior_data_session.Behav](#)
[attribute](#)), 88
[lims_working_directory](#) () (al-
[lensdk.internal.api.queries.biophysical_module_reader.Biophysic](#)
[method](#)), 240
[lims_working_directory](#) () (al-
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConf](#)
[method](#)), 243
[LimsApi](#) (class in [allensdk.internal.api.lims_api](#)), 248
[line_crossing_x](#) () (in module [al-](#)
[lensdk.model.glif.glif_neuron](#)), 322
[line_crossing_y](#) () (in module [al-](#)
[lensdk.model.glif.glif_neuron](#)), 322
[line_stats](#) () ([allensdk.brain_observatory.sync_dataset.Dataset](#)
[method](#)), 178
[linear_transform_from_intervals](#) ()
(in module [al-](#)
[lensdk.brain_observatory.ecephys.align_timestamps.barcode](#)),
110
[linux_to_windows](#) () (in module [al-](#)
[lensdk.internal.core.lims_utilities](#)), 266
[list_api_methods](#) () (al-
[lensdk.brain_observatory.behavior.behavior_data_session.Behav](#)
[method](#)), 88
[list_cells](#) () ([allensdk.api.queries.cell_types_api.CellTypesApi](#)
[method](#)), 50
[list_cells_api](#) () (al-
[lensdk.api.queries.cell_types_api.CellTypesApi](#)
[method](#)), 50
[list_column_definition_class_names](#) ()
([allensdk.api.queries.brain_observatory_api.BrainObservatoryAp](#)

method), 48

list_isi_experiments() (al-
lensdk.api.queries.brain_observatory_api.BrainObservatoryApi
method), 48

list_neuronal_models() (al-
lensdk.api.queries.glif_api.GlifApi method),
52

list_of_dicts_to_dict_of_lists() (in module al-
lensdk.brain_observatory.receptive_field_analysis_tools),
142

list_stimuli() (al-
lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet
method), 197

load() (allensdk.brain_observatory.ecephys.file_io.continuous_file.ContinuousFile
method), 116

load() (allensdk.brain_observatory.sync_dataset.Dataset
method), 178

load() (allensdk.config.app.application_config.ApplicationConfig
method), 180

load() (allensdk.model.biophys_sim.config.Config
method), 316

load_annotation() (in module al-
lensdk.internal.mouse_connectivity.interval_unionization_utilities),
299

load_api_schema() (allensdk.api.api.Api method),
74

load_arrays() (in module al-
lensdk.internal.pipeline_modules.run_annotated_region_metrics),
308

load_cell_parameters() (al-
lensdk.internal.model.biophysical.deap_utils.Utils
method), 276

load_cell_parameters() (al-
lensdk.model.biophysical.utils.AllActiveUtils
method), 318

load_cell_parameters() (al-
lensdk.model.biophysical.utils.Utils method),
319

load_config() (allensdk.config.manifest.Manifest
method), 186

load_csv() (allensdk.api.cache.Cache method), 77

load_datasets_by_relnames() (in module al-
lensdk.core.h5_utilities), 202

load_description() (in module al-
lensdk.model.biophysical.runner), 317

load_experiment() (in module al-
lensdk.internal.ephys.plot_qc_figures), 269

load_experiment() (in module al-
lensdk.internal.ephys.plot_qc_figures3),
271

load_eye_tracking_hdf() (in module al-
lensdk.brain_observatory.behavior.eye_tracking_processing),
94

load_frame() (in module al-
lensdk.internal.brain_observatory.ophys_session_decomposition),
246

load_json() (allensdk.api.cache.Cache method), 77

load_manifest() (allensdk.api.cache.Cache
method), 77

load_manifest() (al-
lensdk.internal.model.biophysical.run_optimize.RunOptimize
method), 277

load_manifest() (al-
lensdk.model.biophysical.run_simulate.RunSimulate
method), 317

load_nwb_extension() (in module al-
lensdk.internal.model.biophysical.passive_fitting.neuron_utils),
104

load_pickle() (in module al-
lensdk.brain_observatory.behavior.stimulus_processing),
104

load_synwb_extension() (in module al-
lensdk.brain_observatory.nwb.metadata),
131

load_sweep() (in module al-
lensdk.internal.model.data_access), 290

load_utilities() (in module al-
lensdk.model.glif.simulate_neuron), 328

load_sweeps() (in module al-
lensdk.internal.model.data_access), 290

local_time() (in module al-
lensdk.brain_observatory.behavior.trials_processing),
108

LocallySparseNoise (class in al-
lensdk.brain_observatory.locally_sparse_noise),
155

locate_h5_objects() (in module al-
lensdk.core.h5_utilities), 202

locate_median() (in module al-
lensdk.brain_observatory.receptive_field_analysis.chisquarerrf),
139

log (allensdk.config.manifest.Manifest attribute), 186

log (allensdk.config.model.description_parser.DescriptionParser
attribute), 184

log (allensdk.config.model.formats.json_description_parser.JsonDescriptionParser
attribute), 182

log (allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser
attribute), 183

long_squares_features() (al-
lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor
method), 227

long_squares_stim_amps() (al-
lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor
method), 227

LSN (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise
attribute), 155

lsn_coordinate_to_monitor_coordinate()

(in module al- 278
 lensdk.brain_observatory.stimulus_info), main() (in module al-
 172 lensdk.internal.model.glif.find_sweeps), 280
 LSN_GREY (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise (in module al-
 attribute), 155 lensdk.internal.morphology.validate_swc),
 lsn_image_to_screen() (al- 298
 lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor (in module al-
 method), 171 lensdk.internal.pipeline_modules.gbm.generate_gbm_analysis_ru
 lsn_image_to_screen() (al- 307
 lensdk.brain_observatory.stimulus_info.Monitor main() (in module al-
 method), 171 lensdk.internal.pipeline_modules.gbm.generate_gbm_heatmap),
 LSN_mask (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise (in module al-
 attribute), 155 main() (in module al-
 LSN_OFF (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise (in module al-
 attribute), 155 lensdk.internal.pipeline_modules.gbm.generate_gbm_sample_me
 LSN_OFF_SCREEN (al- main() (in module al-
 lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise (in module al-
 attribute), 155 lensdk.internal.pipeline_modules.run_annotated_region_metrics),
 LSN_ON (allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise (in module al-
 attribute), 155 lensdk.internal.pipeline_modules.run_demixing),
 lsna_check_hvas() (in module al- 308
 lensdk.internal.pipeline_modules.run_observatory_thumbnail(nails), (in module al-
 312 lensdk.internal.pipeline_modules.run_dff_computation),
 308
M main() (in module al-
 main() (in module allensdk.brain_observatory.dff), 153 lensdk.internal.pipeline_modules.run_neuropil_correction),
 main() (in module al- 309
 lensdk.brain_observatory.ecephys.stimulus_table.visualization.view_blocks), module al-
 120 lensdk.internal.pipeline_modules.run_observatory_analysis),
 main() (in module al- 309
 lensdk.brain_observatory.session_analysis), main() (in module al-
 166 lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 main() (in module al- 312
 lensdk.internal.ephys.plot_qc_figures), 269 main() (in module al-
 main() (in module al- lensdk.internal.pipeline_modules.run_ophys_eye_calibration),
 lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit),
 275 main() (in module al-
 main() (in module al- lensdk.internal.pipeline_modules.run_ophys_session_decomposit
 lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit2),
 275 main() (in module al-
 main() (in module al- lensdk.internal.pipeline_modules.run_ophys_time_sync),
 lensdk.internal.model.biophysical.passive_fitting.neuron_passive_fit_elec),
 275 main() (in module al-
 main() (in module al- lensdk.internal.pipeline_modules.run_tissuecyte_unionize_classic
 lensdk.internal.model.biophysical.passive_fitting.preprocess), 315
 275 main() (in module al-
 main() (in module al- lensdk.model.biophysical.run_simulate),
 lensdk.internal.model.biophysical.run_optimize), 317
 277 main() (in module al-
 main() (in module al- lensdk.model.glif.simulate_neuron), 328
 lensdk.internal.model.biophysical.run_passive_fit), main() (in module allensdk.morphology.validate_swc),
 278 329
 main() (in module al- make_bbs() (in module al-
 lensdk.internal.model.biophysical.run_simulate_lims), lensdk.internal.brain_observatory.mask_set),

260
 make_blended_tile() (in module al-
 lensdk.internal.mouse_connectivity.tissuecyte_stimulus_table), 290
 306
 make_category_dummy() (in module al-
 lensdk.brain_observatory.chisquare_categorical), 185
 148
 make_cell_html() (in module al-
 lensdk.internal.ephys.plot_qc_figures), 269
 make_cell_html() (in module al-
 lensdk.internal.ephys.plot_qc_figures3),
 271
 make_cell_page() (in module al-
 lensdk.internal.ephys.plot_qc_figures), 269
 make_cell_page() (in module al-
 lensdk.internal.ephys.plot_qc_figures3),
 271
 make_display_mask() (in module al-
 lensdk.brain_observatory.stimulus_info),
 172
 make_fit() (allensdk.internal.model.biophysical.run_optimize.RunOptimize
 method), 277
 make_fit_json_file() (al-
 lensdk.internal.model.biophysical.make_deep_fit_json.Report
 method), 277
 make_pincushion_plot() (in module al-
 lensdk.brain_observatory.circle_plots), 150
 make_pixel_counter() (al-
 lensdk.mouse_connectivity.grid.subimage.base_subimage.Subimage
 method), 331
 make_ratio_volume() (al-
 lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder
 method), 336
 make_spontaneous_activity_tables() (in
 module al-
 lensdk.brain_observatory.ecephys.stimulus_table.ecephys_preprocessing
 124
 make_structure_mask() (al-
 lensdk.core.reference_space.ReferenceSpace
 method), 213
 make_sweep_html() (in module al-
 lensdk.internal.ephys.plot_qc_figures), 269
 make_sweep_html() (in module al-
 lensdk.internal.ephys.plot_qc_figures3),
 272
 make_sweep_page() (in module al-
 lensdk.internal.ephys.plot_qc_figures), 269
 make_sweep_page() (in module al-
 lensdk.internal.ephys.plot_qc_figures3),
 272
 makeBasis_StimKernel() (in module al-
 lensdk.internal.model.GLM), 290
 makeBasis_StimKernel_exp() (in module al-
 lensdk.internal.model.GLM), 290
 makeFitStruct_GLM() (in module al-
 lensdk.internal.model.GLM), 290
 ManifestBuilder (class in al-
 lensdk.config.manifest_builder), 187
 ManifestVersionError, 187
 manifest_dataframe() (allensdk.api.cache.Cache
 method), 77
 MANIFEST_VERSION (al-
 lensdk.core.brain_observatory_cache.BrainObservatoryCache
 attribute), 190
 MANIFEST_VERSION (al-
 lensdk.core.cell_types_cache.CellTypesCache
 attribute), 199
 MANIFEST_VERSION (al-
 lensdk.core.mouse_connectivity_cache.MouseConnectivityCache
 attribute), 204
 MANIFEST_VERSION (al-
 lensdk.core.reference_space_cache.ReferenceSpaceCache
 attribute), 214
 ManifestBuilder (class in al-
 lensdk.config.manifest_builder), 187
 ManifestVersionError, 187
 make_masks() (al-
 lensdk.core.reference_space.ReferenceSpace
 method), 213
 make_report() (in module al-
 lensdk.brain_observatory.ecephys.stimulus_table.naming_utilities
 125
 map_monitor_coordinate_to_stimulus_coordinate() (in
 module al-
 lensdk.brain_observatory.stimulus_info),
 172
 map_monitor_coordinate_to_template_coordinate() (al-
 lensdk.brain_observatory.stimulus_info),
 173
 map_stimulus() (al-
 lensdk.brain_observatory.stimulus_info.Monitor
 method), 172
 map_stimulus() (in module al-
 lensdk.brain_observatory.stimulus_info),
 173
 map_stimulus_coordinate_to_monitor_coordinate() (in
 module al-
 lensdk.brain_observatory.stimulus_info),
 173
 map_stimulus_names() (in module al-
 lensdk.brain_observatory.ecephys.stimulus_table.naming_utilities
 125
 map_template_coordinate_to_monitor_coordinate() (in
 module al-
 lensdk.brain_observatory.stimulus_info),
 173
 Marker (class in allensdk.core.swc), 222
 Marker (class in allensdk.internal.core.swc), 268
 MARKER_FILE_TYPE (al-

[lensdk.api.queries.cell_types_api.CellTypesApi](#) (class in [lensdk.api.queries.cell_types_api](#)), 49
[MARKER_KEY](#) ([allensdk.core.cell_types_cache.CellTypesCache](#) attribute), 199
[mask](#) ([allensdk.brain_observatory.stimulus_info.Monitor](#) attribute), 172
[Mask](#) (class in [allensdk.brain_observatory.roi_masks](#)), 161
[mask\(\)](#) ([allensdk.internal.brain_observatory.mask_set.MaskSet](#) method), 260
[mask_is_union_of_set\(\)](#) (method in [allensdk.internal.brain_observatory.mask_set.MaskSet](#)), 260
[mask_nulls\(\)](#) (in module [allensdk.internal.ephys.plot_qc_figures](#)), 269
[mask_nulls3\(\)](#) (in module [allensdk.internal.ephys.plot_qc_figures3](#)), 272
[mask_stimulus_template\(\)](#) (in module [allensdk.brain_observatory.stimulus_info](#)), 173
[MaskSet](#) (class in [allensdk.internal.brain_observatory.mask_set](#)), 260
[match_barcodes\(\)](#) (in module [allensdk.brain_observatory.ecephys.align_timestamps_barcode](#)), 111
[max_cb\(\)](#) (in module [allensdk.internal.mouse_connectivity.projection_thumbnail_generator.projection_str\(p\)](#)), 303
[max_of_line_and_const\(\)](#) (in module [allensdk.model.glif.glif_neuron_methods](#)), 325
[max_projection\(\)](#) (in module [allensdk.internal.mouse_connectivity.projection_thumbnail.projection_functions](#)), 304
[max_voxel_density](#) (in module [allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_base_unionize.record.TissuecyteBaseUnionize](#) attribute), 301
[max_voxel_index](#) (in module [allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_base_unionize.record.TissuecyteBaseUnionize](#) attribute), 301
[mean_features_spike_zero\(\)](#) (in module [allensdk.ephys.extract_cell_features](#)), 239
[mean_response](#) (in module [allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise](#) attribute), 155
[mean_sweep_response](#) (in module [allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis](#) attribute), 170
[medfilt_custom\(\)](#) (in module [allensdk.internal.brain_observatory.itracker_utils](#)), 259
[median_absolute_deviation\(\)](#) (in module [allensdk.internal.brain_observatory.itracker_utils](#)), 259
[meets_engagement_criteria\(\)](#) (in module [allensdk.brain_observatory.behavior.criteria](#)), 91
[membrane_time_constant\(\)](#) (in module [allensdk.ephys.ephys_extractor](#)), 231
[memoize\(\)](#) (in module [allensdk.api.cache](#)), 79
[merge_mean_response\(\)](#) (static method in [allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise](#)), 155
[Metadata](#) ([allensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession](#) attribute), 88
[METHOD](#) ([allensdk.core.authentication.CredentialProvider](#) attribute), 188
[METHOD](#) ([allensdk.core.authentication.EnvCredentialProvider](#) attribute), 189
[midpoint\(\)](#) (in module [allensdk.internal.morphology.node](#)), 298
[min_of_line_and_zero\(\)](#) (in module [allensdk.model.glif.glif_neuron_methods](#)), 325
[minmax_norm\(\)](#) (in module [allensdk.internal.mouse_connectivity.projection_thumbnail.visualization](#)), 304
[MissingDataError](#), 201
[MissingStimulusException](#), 147
[MissingSweepException](#), 279
[MLIN\(\)](#) (in module [allensdk.internal.model.glif.MLIN](#)), 278
[model_generate_projection_str\(p\)](#) (in module [allensdk.internal.model.glif.error_functions](#)), 279
[mod_file_entries\(\)](#) (in module [allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) attribute), 240
[mod_file_entries\(\)](#) (in module [allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) attribute), 243
[MOD_FILE_TYPE_ID](#) (in module [allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) attribute), 240
[MOD_FILE_TYPE_ID](#) (in module [allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) attribute), 243
[MODEL](#) ([allensdk.api.queries.rma_api.RmaApi](#) attribute), 65
[model_query\(\)](#) (in module [allensdk.api.queries.rma_api.RmaApi](#) method), 67
[model_stage\(\)](#) (in module [allensdk.api.queries.rma_api.RmaApi](#) method), 67

[lensdk.api.queries.rma_api.RmaApi](#) (method), [MOTION_CORRECTION_DATASETS](#) (al-
[67](#) [lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet](#))
[model_type\(\)](#) ([allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#)
method), [240](#) [MOUSE](#) ([allensdk.api.queries.cell_types_api.CellTypesApi](#)
method), [240](#)
[modify_parameter\(\)](#) (al- [attribute](#)), [49](#)
[lensdk.model.glif.glif_neuron_methods.GlifNeuronMethod](#) ([2011](#) ([allensdk.api.queries.reference_space_api.ReferenceSpaceApi](#)
method), [323](#) [attribute](#)), [63](#)
[moments2\(\)](#) (in module al- [MOUSE_ATLAS_PRODUCTS](#) (al-
[lensdk.brain_observatory.receptive_field_analysis.fitgaussian_2d](#) ([allensdk.api.queries.mouse_atlas_api.MouseAtlasApi](#)
[141](#) [attribute](#)), [57](#)
[Monitor](#) (class in al- [MOUSE_ORGANISM](#) (al-
[lensdk.brain_observatory.stimulus_info](#)), [lensdk.api.queries.mouse_atlas_api.MouseAtlasApi](#)
[171](#) [attribute](#)), [57](#)
[monitor_coordinate_to_lsn_coordinate\(\)](#) [MouseAtlasApi](#) (class in al-
(in module al- [lensdk.api.queries.mouse_atlas_api](#)), [57](#)
[lensdk.brain_observatory.stimulus_info](#)), [MouseConnectivityApi](#) (class in al-
[173](#) [lensdk.api.queries.mouse_connectivity_api](#)),
[monitor_coordinate_to_natural_movie_coordinate\(\)](#) ([58](#)
(in module al- [MouseConnectivityApiPrerelease](#) (class in al-
[lensdk.brain_observatory.stimulus_info](#)), [lensdk.internal.api.queries.mouse_connectivity_api_prerelease](#)),
[174](#) [242](#)
[monitor_delay\(\)](#) (in module al- [MouseConnectivityCache](#) (class in al-
[lensdk.internal.brain_observatory.time_sync](#)), [lensdk.core.mouse_connectivity_cache](#)),
[264](#) [203](#)
[Morphology](#) (class in [allensdk.core.swc](#)), [223](#) [MouseConnectivityCachePrerelease](#)
[Morphology](#) (class in al- (class in al-
[lensdk.internal.morphology.morphology](#)), [lensdk.internal.core.mouse_connectivity_cache_prerelease](#)),
[291](#) [266](#)
[MORPHOLOGY_FEATURES_KEY](#) (al- [movingaverage\(\)](#) (in module al-
[lensdk.core.cell_types_cache.CellTypesCache](#) [lensdk.brain_observatory.dff](#)), [153](#)
[attribute](#)), [199](#) [movingmode_fast\(\)](#) (in module al-
[morphology_file_entries\(\)](#) (al- [lensdk.brain_observatory.dff](#)), [153](#)
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) ([allensdk.internal.api.mtrain_api](#)),
method), [240](#) [249](#)
[morphology_file_entries\(\)](#) (al- [MtrainSqlApi](#) (class in al-
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) ([allensdk.internal.api.mtrain_api](#)), [249](#)
method), [243](#) [multi_dataframe_merge\(\)](#) (in module al-
[morphology_path\(\)](#) (al- [lensdk.brain_observatory.session_analysis](#)),
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#)
method), [241](#) [166](#)
[morphology_path\(\)](#) (al- **N**
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) (in module al-
method), [243](#) [lensdk.brain_observatory.behavior.criteria](#)),
[MORPHOLOGY_TYPE_ID](#) (al- [91](#)
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) [ReporterStatus](#) at-
[attribute](#)), [240](#) [tribute](#)), [200](#)
[MORPHOLOGY_TYPE_ID](#) (al- [nan_get\(\)](#) (in module al-
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) ([allensdk.internal.ephys.core.feature_extract](#)),
[attribute](#)), [243](#) [269](#)
[MorphologyColors](#) (class in al- [natural_movie_coordinate_to_monitor_coordinate\(\)](#)
[lensdk.internal.morphology.morphvis](#)), [295](#) (in module al-
[mostly_useful\(\)](#) (in module al- [lensdk.brain_observatory.stimulus_info](#)),
[lensdk.brain_observatory.behavior.criteria](#)), [174](#)
[91](#) [natural_movie_image_to_screen\(\)](#) (al-

`lensdk.brain_observatory.stimulus_info.Monitor` (method), 172
`lensdk.brain_observatory.stimulus_info.Monitor` (method), 280
`natural_scene_coordinate_to_monitor_coordinate` (`lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader` (in module `lensdk.brain_observatory.stimulus_info`), 174), 243
`natural_scene_image_to_screen()` (`lensdk.brain_observatory.stimulus_info.Monitor` (method), 172), 243
`NaturalMovie` (class in `lensdk.brain_observatory.natural_movie`), 156
`NaturalScenes` (class in `lensdk.brain_observatory.natural_scenes`), 157
`nearest_neuron_sampling_rate()` (`lensdk.model.biophysical.utils.Utils` (static method), 319), 334
`needs_data_refresh()` (`lensdk.brain_observatory.session_api_utils.ParamsMixin` (method), 167), 334
`NEGATIVE` (`allensdk.core.cell_types_cache.ReporterStatusNO_RECONSTRUCTION` (attribute), 200), 290
`neurodata_doc` (`lensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema` (attribute), 100), 290
`neurodata_doc` (`lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema` (attribute), 101), 290
`neurodata_doc` (`lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema` (attribute), 102), 290
`neurodata_type` (`lensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema` (attribute), 100), 290
`neurodata_type` (`lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema` (attribute), 101), 290
`neurodata_type` (`lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema` (attribute), 102), 290
`neurodata_type_inc` (`lensdk.brain_observatory.behavior.schemas.BehaviorTaskParametersSchema` (attribute), 100), 290
`neurodata_type_inc` (`lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema` (attribute), 101), 290
`neurodata_type_inc` (`lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema` (attribute), 102), 290
`neuron` (`allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils` (attribute), 316), 223
`neuron_parameter_count()` (`lensdk.internal.model.glif.glif_experiment.GlifExperiment` (attribute), 291), 223

`nodes()` (`allensdk.core.simple_tree.SimpleTree` method), 218
`nodes()` (`allensdk.internal.core.simpletree.SimpleTree` method), 267
`nodes_by_property()` (`allensdk.core.simple_tree.SimpleTree` method), 218
`NoEyeTrackingException`, 147
`noise_std()` (in module `allensdk.brain_observatory.dff`), 153
`nonraising_ks_2samp()` (in module `allensdk.brain_observatory.stimulus_analysis`), 171
`norm_diff()` (in module `allensdk.ephys.ephys_features`), 238
`norm_sq_diff()` (in module `allensdk.ephys.ephys_features`), 238
`normalize_actual_parameters()` (`allensdk.internal.model.biophysical.deap_utils.Utils` method), 276
`normalize_F()` (in module `allensdk.brain_observatory.r_neuropil`), 161
`normalize_intensity()` (in module `allensdk.internal.mouse_connectivity.projection_thumbnail_visualization_utilities`), 304
`normalizecols()` (in module `allensdk.internal.model.GLM`), 290
`np_sitk_convert()` (in module `allensdk.mouse_connectivity.grid.utilities.image_utilities`), 335
O
`nrn` (`allensdk.model.biophys_sim.neuron.hoc_utils.HocUtils` attribute), 316
`nrnivmodl()` (`allensdk.internal.model.biophysical.run_optimize.RunOptimize` method), 277
`nrnivmodl()` (`allensdk.model.biophysical.run_simulate.RunSimulate` method), 317
`num_contingent_trials()` (in module `allensdk.brain_observatory.behavior.session_metrics`), 102
`num_nodes` (`allensdk.core.swc.Morphology` attribute), 225
`num_nodes` (`allensdk.internal.morphology.morphology.Morphology` attribute), 294
`NUM_ROWS` (`allensdk.api.queries.rma_api.RmaApi` attribute), 65
`num_trees` (`allensdk.core.swc.Morphology` attribute), 225
`num_trees` (`allensdk.internal.morphology.morphology.Morphology` attribute), 294
`number_of_cells` (`allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` attribute), 197
`number_ori` (`allensdk.brain_observatory.drifting_gratings.DriftingGratings` attribute), 154
`number_ori` (`allensdk.brain_observatory.static_gratings.StaticGratings` attribute), 168
`number_phase` (`allensdk.brain_observatory.static_gratings.StaticGratings` attribute), 168
`number_scenes` (`allensdk.brain_observatory.natural_scenes.NaturalScenes` attribute), 157
`number_sf` (`allensdk.brain_observatory.static_gratings.StaticGratings` attribute), 168
`number_tf` (`allensdk.brain_observatory.drifting_gratings.DriftingGratings` attribute), 154
`numbercells` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` attribute), 170
`NWB_FILE_TYPE` (`allensdk.api.queries.brain_observatory_api.BrainObservatoryApi` attribute), 46
`NWB_FILE_TYPE` (`allensdk.api.queries.cell_types_api.CellTypesApi` attribute), 49
`NWB_FILE_TYPE` (`allensdk.api.queries.glif_api.GlifApi` attribute), 52
`NwbApi` (class in `allensdk.brain_observatory.nwb.nwb_api`), 211
`NwbDataSet` (class in `allensdk.core.nwb_data_set`), 208
`nwbfile` (`allensdk.brain_observatory.nwb.nwb_api.NwbApi` attribute), 131
`object_norm_eye_coordinates()` (in module `allensdk.internal.brain_observatory.eye_calibration`), 256
`object_norm_eye_coordinates()` (in module `allensdk.internal.brain_observatory.eye_calibration`), 256
`object_norm_eye_coordinates()` (in module `allensdk.internal.brain_observatory.eye_calibration`), 256
`one()` (in module `allensdk`), 338
`one_file_call_caching()` (in module `allensdk.api.caching_utilities`), 80
`OneOrMoreResultExpectedError`, 250
`OneOrMoreResultExpectedError`, 338
`ONLY` (`allensdk.api.queries.rma_api.RmaApi` attribute), 65
`only_except_tabular_clause()` (`allensdk.api.queries.rma_api.RmaApi` method), 68
`OntologiesApi` (class in `allensdk.api.queries.ontologies_api`), 61
`Ontology` (class in `allensdk.core.ontology`), 210
`open()` (`allensdk.brain_observatory.frame_stream.CvInputStream` method), 257
`open()` (`allensdk.brain_observatory.frame_stream.FfmpegInputStream` method), 258

`open()` (`allensdk.internal.brain_observatory.frame_stream.FrameStream` `FrameOutputStream` `Reader` (class in `al-`
`method`), 258 `lensdk.internal.api.queries.optimize_config_reader`),
`open()` (`allensdk.internal.brain_observatory.frame_stream.FrameInputStream`
`method`), 258 `optional_lims_inputs()` (in module `al-`
`open()` (`allensdk.internal.brain_observatory.frame_stream.FrameOutputStream` `lensdk.brain_observatory.argschema_utilities`),
`method`), 258 147
`open_corona_plot()` (al- `optional_polys` (al-
`lensdk.brain_observatory.natural_scenes.NaturalScenes` `lensdk.mouse_connectivity.grid.subimage.base_subimage.Polygon`
`method`), 157 `attribute`), 330
`open_fan_plot()` (al- `optional_polys` (al-
`lensdk.brain_observatory.static_gratings.StaticGratings` `lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubimage`
`method`), 168 `attribute`), 332
`open_pincushion_plot()` (al- `OPTIONS` (`allensdk.api.queries.rma_api.RmaApi` `at-`
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`, 65
`method`), 155 `options_clause()` (al-
`open_star_plot()` (al- `lensdk.api.queries.rma_api.RmaApi` `method`),
`lensdk.brain_observatory.drifting_gratings.DriftingGratings`, 68
`method`), 154 `OPTOGENETIC_STIMULATION_KEYS` (al-
`open_track_plot()` (al- `lensdk.brain_observatory.sync_dataset.Dataset`
`lensdk.brain_observatory.natural_movie.NaturalMovie` `attribute`), 175
`method`), 156 `opts` (`allensdk.brain_observatory.argschema_utilities.RaisingSchema`
`open_view_on_binary()` (in module `al-` `attribute`), 147
`lensdk.internal.brain_observatory.ophys_session_data_session` `lensdk.brain_observatory.behavior.mtrain.ExtendedTrialSchema`
260 `attribute`), 96
`OPHYS_ANALYSIS_FILE_TYPE` (al- `opts` (`allensdk.brain_observatory.behavior.schemas.BehaviorMetadataSchema`
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `attribute`), 99
`attribute`), 46 `opts` (`allensdk.brain_observatory.behavior.schemas.BehaviorTaskParameters`
`ophys_delta` (`allensdk.internal.pipeline_modules.run_ophys_time_sync` `lensdk.internal.pipeline_modules.sync_outputs`
`attribute`), 313 `opts` (`allensdk.brain_observatory.behavior.schemas.CompleteOphysBehaviorSchema`
`OPHYS_EVENTS_FILE_TYPE` (al- `attribute`), 100
`lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` `lensdk.brain_observatory.behavior.schemas.OphysBehaviorMetadataSchema`
`attribute`), 46 `attribute`), 101
`ophys_experiment_ids` (al- `opts` (`allensdk.brain_observatory.behavior.schemas.OphysMetadataSchema`
`lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession` `attribute`), 88
`attribute`), 88 `opts` (`allensdk.brain_observatory.behavior.schemas.RaisingSchema`
`ophys_session_id` (al- `attribute`), 101
`lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession` `lensdk.brain_observatory.behavior.schemas.SubjectMetadataSchema`
`attribute`), 89 `attribute`), 102
`ophys_times` (`allensdk.internal.pipeline_modules.run_ophys_time_sync` `lensdk.internal.pipeline_modules.sync_outputs`
`attribute`), 313 `attribute`), 132
`ophys_timestamps` (al- `ORDER` (`allensdk.api.queries.rma_api.RmaApi` `attribute`),
`lensdk.internal.brain_observatory.time_sync.OphysTimeAligner`
`attribute`), 263 `order_clause()` (al-
`OphysBehaviorMetadataSchema` (class in `al-` `lensdk.api.queries.rma_api.RmaApi` `method`),
`lensdk.brain_observatory.behavior.schemas`), 68
100 `order_rois_by_object_list()` (in module `al-`
`OphysLimsApi` (class in `al-` `lensdk.internal.brain_observatory.roi_filter_utils`),
`lensdk.internal.api.ophys_lims_api`), 249 263
`OphysMetadataSchema` (class in `al-` `organize_sweeps_by_name()` (in module `al-`
`lensdk.brain_observatory.behavior.schemas`), `lensdk.internal.model.glif.find_sweeps`), 280
101 `orivals` (`allensdk.brain_observatory.drifting_gratings.DriftingGratings`
`OphysTimeAligner` (class in `al-` `attribute`), 154
`lensdk.internal.brain_observatory.time_sync`), `orivals` (`allensdk.brain_observatory.static_gratings.StaticGratings`
263 `attribute`), 168

outdated (allensdk.config.manifest.ManifestVersionError.parse_arguments() (in module al-
 attribute), 187
 outlier_cost() (al- parse_command_line_args() (al-
 lensdk.internal.brain_observatory.fit_ellipse.FitEllipse lensdk.config.app.application_config.ApplicationConfig
 method), 257
 method), 181
 output() (allensdk.internal.mouse_connectivity.interval_unionize_unionize_record.TissotBaseUnionize
 method), 301
 output() (allensdk.internal.mouse_connectivity.interval_unionize_unionize_record.Unionize
 method), 302
 output_directory() (al- parse_input() (in module al-
 lensdk.internal.pipeline_modules.run_demixing),
 lensdk.internal.pipeline_modules.run_dff_computation),
 method), 208
 parse_input() (in module al-
 OutputFile (class in al- lensdk.internal.pipeline_modules.run_observatory_thumbnails),
 lensdk.brain_observatory.argschema_utilities), 312
 146
 OutputGrabber (class in al- lensdk.internal.pipeline_modules.run_ophys_session_decomposit
 lensdk.internal.model.biophysical.passive_fitting.output_grabber),
 275
 overlap_fraction() (al- parse_input_data() (in module al-
 lensdk.internal.brain_observatory.mask_set.MaskSet 312
 method), 260
 parse_neuron_output() (in module al-
 overlaps_motion_border (al- lensdk.internal.model.biophysical.passive_fitting.neuron_utils),
 lensdk.brain_observatory.roi_masks.Mask 275
 attribute), 162
 parse_obj() (in module allensdk.core.obj_utilities),
 210
 parse_stim_repr() (in module al-
 P
 lensdk.brain_observatory.ecephys.stimulus_table.stimulus_param
 pageable() (in module al- 127
 lensdk.api.queries.rma_pager), 70
 pager() (allensdk.api.queries.rma_pager.RmaPager parser_for_extension() (al-
 static method), 70
 lensdk.config.model.description_parser.DescriptionParser
 panel_size (allensdk.brain_observatory.stimulus_info.Monitor method), 184
 attribute), 172
 password (allensdk.core.authentication.DbCredentials
 ParamsMixin (class in al- attribute), 188
 lensdk.brain_observatory.session_api_utils),
 167
 paste_slice() (al-
 parent() (allensdk.core.simple_tree.SimpleTree lensdk.mouse_connectivity.grid.image_series_gridder.ImageSerie
 method), 218
 method), 336
 parent() (allensdk.internal.core.simpletree.SimpleTree paste_subimage() (al-
 method), 267
 lensdk.mouse_connectivity.grid.image_series_gridder.ImageSerie
 method), 336
 parent_id() (allensdk.core.simple_tree.SimpleTree path (allensdk.brain_observatory.nwb.nwb_api.NwbApi
 method), 218
 attribute), 131
 parent_id() (allensdk.internal.core.simpletree.SimpleTree path_to_list() (al-
 method), 267
 lensdk.core.structure_tree.StructureTree static
 method), 222
 parent_ids() (allensdk.core.simple_tree.SimpleTree pathfinder() (allensdk.api.cache.Cache static
 method), 218
 method), 77
 parent_of() (allensdk.core.swc.Morphology pause_metrics() (al-
 method), 225
 lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor
 parent_of() (allensdk.internal.morphology.morphology.Morphology method), 228
 method), 294
 parents() (allensdk.core.simple_tree.SimpleTree peak (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 method), 218
 attribute), 170
 parse_arguments() (in module al- peak_run (allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis
 lensdk.internal.model.glif.find_sweeps), 280
 attribute), 170
 period() (allensdk.brain_observatory.sync_dataset.Dataset

`method`), 178
`phasevals` (`allensdk.brain_observatory.static_gratings.StaticGratings` attribute), 168
`PHOTODIODE_KEYS` (`allensdk.brain_observatory.sync_dataset.Dataset` attribute), 175
`PIPE` (`allensdk.api.queries.rma_api.RmaApi` attribute), 65
`pipe_stage()` (`allensdk.api.queries.rma_api.RmaApi` method), 69
`PIPELINE_DATASET` (`allensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet` attribute), 194
`PipelineModule` (class in `allensdk.internal.core.lims_pipeline_module`), 265
`pixel_counter` (`allensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage` attribute), 331
`pixel_size` (`allensdk.brain_observatory.stimulus_info.Monitor` attribute), 172
`pixels_to_visual_degrees()` (`allensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor` method), 171
`pixels_to_visual_degrees()` (`allensdk.brain_observatory.stimulus_info.Monitor` method), 172
`plot()` (`allensdk.brain_observatory.circle_plots.CoronaPlotter` method), 149
`plot()` (`allensdk.brain_observatory.circle_plots.FanPlotter` method), 149
`plot()` (`allensdk.brain_observatory.circle_plots.TrackPlotter` method), 149
`plot_all()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 178
`plot_bit()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 178
`plot_bits()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 178
`plot_blocks()` (in module `allensdk.brain_observatory.ecephys.stimulus_table.visualization`), 121
`plot_cell_correlation()` (in module `allensdk.brain_observatory.observatory_plots`), 158
`plot_cell_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 269
`plot_cell_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_cell_receptive_field()` (`allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 155
`plot_chi_square_summary()` (in module `allensdk.brain_observatory.receptive_field_analysis.visualization`), 143
`plot_combined_speed()` (in module `allensdk.brain_observatory.observatory_plots`), 158
`plot_condition_histogram()` (in module `allensdk.brain_observatory.observatory_plots`), 158
`plot_direction_selectivity()` (`allensdk.brain_observatory.drifting_gratings.DriftingGratings` method), 154
`plot_drifting_traces()` (in module `allensdk.brain_observatory.brain_observatory_plotting`), 147
`plot_ellipses()` (in module `allensdk.brain_observatory.receptive_field_analysis.visualization`), 143
`plot_fi_curve_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 269
`plot_fi_curve_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_gaussian_fit()` (in module `allensdk.brain_observatory.receptive_field_analysis.visualization`), 143
`plot_hero_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_hero_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_images()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_images()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_instantaneous_threshold_thumbnail()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_instantaneous_threshold_thumbnail()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_line()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 178
`plot_lines()` (`allensdk.brain_observatory.sync_dataset.Dataset` method), 178
`plot_long_square_summary()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_long_square_summary()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 272
`plot_lsn_traces()` (in module `allensdk.brain_observatory.receptive_field_analysis.visualization`), 143

`lensdk.brain_observatory.brain_observatory_plotting`), `method`), 169
`plot_mask()` (in module `al-` `lensdk.brain_observatory.static_gratings.StaticGratings` `method`), 169
`plot_mask_outline()` (in module `al-` `lensdk.brain_observatory.drifting_gratings.DriftingGratings` `method`), 154
`plot_masks()` (in module `al-` `lensdk.internal.brain_observatory.demix_report`), 158
`plot_mean_waveforms()` (in module `al-` `lensdk.brain_observatory.ecephys.visualization`), 128
`plot_msr_summary()` (in module `al-` `lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
`plot_negative_baselines()` (in module `al-` `lensdk.brain_observatory.demixer`), 151
`plot_negative_baselines()` (in module `al-` `lensdk.internal.brain_observatory.demixer`), 253
`plot_negative_transients()` (in module `al-` `lensdk.brain_observatory.demixer`), 151
`plot_negative_transients()` (in module `al-` `lensdk.internal.brain_observatory.demixer`), 253
`plot_ns_traces()` (in module `al-` `lensdk.brain_observatory.brain_observatory_plotting`), 148
`plot_onetrace()` (in module `al-` `lensdk.brain_observatory.dff`), 153
`plot_orientation_selectivity()` (in module `al-` `lensdk.brain_observatory.drifting_gratings.DriftingGratings` `method`), 154
`plot_orientation_selectivity()` (in module `al-` `lensdk.brain_observatory.static_gratings.StaticGratings` `method`), 169
`plot_overlap_masks_lengthOne()` (in module `al-` `lensdk.brain_observatory.demixer`), 151
`plot_overlap_masks_lengthOne()` (in module `al-` `lensdk.internal.brain_observatory.demixer`), 253
`plot_p_values()` (in module `al-` `lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
`plot_population_receptive_field()` (in module `al-` `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` `method`), 156
`plot_preferred_direction()` (in module `al-` `lensdk.brain_observatory.drifting_gratings.DriftingGratings` `method`), 154
`plot_preferred_orientation()` (in module `al-` `lensdk.brain_observatory.static_gratings.StaticGratings` `method`), 170
`plot_preferred_spatial_frequency()` (in module `al-` `lensdk.brain_observatory.static_gratings.StaticGratings` `method`), 169
`plot_preferred_temporal_frequency()` (in module `al-` `lensdk.brain_observatory.drifting_gratings.DriftingGratings` `method`), 154
`plot_pupil_location()` (in module `al-` `lensdk.brain_observatory.observatory_plots`), 158
`plot_radial_histogram()` (in module `al-` `lensdk.brain_observatory.observatory_plots`), 159
`plot_ramp_figures()` (in module `al-` `lensdk.internal.ephys.plot_qc_figures`), 270
`plot_ramp_figures()` (in module `al-` `lensdk.internal.ephys.plot_qc_figures3`), 273
`plot_receptive_field()` (in module `al-` `lensdk.brain_observatory.observatory_plots`), 159
`plot_receptive_field_analysis_data()` (in module `al-` `lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` `method`), 156
`plot_receptive_field_data()` (in module `al-` `lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
`plot_representational_similarity()` (in module `al-` `lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` `method`), 170
`plot_representational_similarity()` (in module `al-` `lensdk.brain_observatory.observatory_plots`), 159
`plot_rheo_figures()` (in module `al-` `lensdk.internal.ephys.plot_qc_figures`), 270
`plot_rheo_figures()` (in module `al-` `lensdk.internal.ephys.plot_qc_figures3`), 273
`plot_rts_blur_summary()` (in module `al-` `lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
`plot_rts_summary()` (in module `al-` `lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
`plot_running_a()` (in module `al-` `lensdk.brain_observatory.brain_observatory_plotting`), 144
`plot_running_speed()` (in module `al-` `lensdk.brain_observatory.visualization`), 145
`plot_running_speed_histogram()` (in module `al-` `lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` `method`), 170

`plot_sag_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_sag_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_selectivity_cumulative_histogram()` (in module `allensdk.brain_observatory.observatory_plots`), 159
`plot_sg_traces()` (in module `allensdk.brain_observatory.brain_observatory_plotting`), 148
`plot_short_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_short_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_single_ap_values()` (in module `allensdk.internal.ephys.plot_qc_figures`), 270
`plot_single_ap_values()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_speed()` (in module `allensdk.brain_observatory.observatory_plots`), 159
`plot_speed_tuning()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` method), 170
`plot_spike_counts()` (in module `allensdk.brain_observatory.ecephys.visualization`), 128
`plot_subthreshold_long_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 271
`plot_subthreshold_long_square_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_sweep_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 271
`plot_sweep_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_sweep_set_summary()` (in module `allensdk.internal.ephys.plot_qc_figures`), 271
`plot_sweep_set_summary()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 273
`plot_sweep_value_figures()` (in module `allensdk.internal.ephys.plot_qc_figures`), 271
`plot_sweep_value_figures()` (in module `allensdk.internal.ephys.plot_qc_figures3`), 274
`plot_time_to_peak()` (`allensdk.brain_observatory.natural_scenes.NaturalScenes` method), 158
`plot_time_to_peak()` (`allensdk.brain_observatory.static_gratings.StaticGratings` method), 169
`plot_time_to_peak()` (in module `allensdk.brain_observatory.observatory_plots`), 159
`plot_traces()` (in module `allensdk.brain_observatory.demixer`), 151
`plot_traces()` (in module `allensdk.internal.brain_observatory.demixer`), 253
`plot_transients()` (in module `allensdk.brain_observatory.demixer`), 151
`plot_transients()` (in module `allensdk.internal.brain_observatory.demixer`), 253
`plotLineRegress1()` (in module `allensdk.internal.model.glif.plotting`), 285
`plotLineRegress1()` (in module `allensdk.internal.model.glif.spike_cutting`), 286
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.plotting`), 285
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.spike_cutting`), 286
`plotLineRegressRed()` (in module `allensdk.internal.model.glif.plotting`), 285
`polar_line_circles()` (in module `allensdk.brain_observatory.circle_plots`), 150
`polar_linspace()` (in module `allensdk.brain_observatory.circle_plots`), 150
`polar_to_xy()` (in module `allensdk.brain_observatory.circle_plots`), 150
`PolarPlotter` (class in `allensdk.brain_observatory.circle_plots`), 149
`PolygonSubImage` (class in `allensdk.mouse_connectivity.grid.subimage.base_subimage`), 330
`populate_stimulus_table()` (`allensdk.brain_observatory.drifting_gratings.DriftingGratings` method), 154
`populate_stimulus_table()` (`allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` method), 156
`populate_stimulus_table()` (`allensdk.brain_observatory.natural_movie.NaturalMovie` method), 156
`populate_stimulus_table()` (`allensdk.brain_observatory.natural_scenes.NaturalScenes` method), 158
`populate_stimulus_table()` (`allensdk.brain_observatory.static_gratings.StaticGratings` method), 169
`populate_stimulus_table()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis` method), 169

`method`), 170
`population_correlation_scatter()` (in module `al-`
`lensdk.brain_observatory.observatory_plots`),
159
`port` (`allensdk.core.authentication.DbCredentials` at-
tribute), 188
`POSITIVE` (`allensdk.core.cell_types_cache.ReporterStatus`
attribute), 201
`post_process_cr()` (in module `al-`
`lensdk.internal.brain_observatory.itracker_utils`),
259
`post_process_pupil()` (in module `al-`
`lensdk.internal.brain_observatory.itracker_utils`),
259
`PostgresQueryMixin` (class in `al-`
`lensdk.internal.api`), 250
`postprocess_unionizes()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize`
method), 300
`postprocess_unionizes()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize`
method), 302
`pre_blank_sec` (al-
`lensdk.brain_observatory.ecephys.file_io.stim_file.CamStim`
attribute), 117
`prepare_nwb_output()` (in module `al-`
`lensdk.model.biophysical.runner`), 317
`PREVIEW` (`allensdk.api.queries.rma_api.RmaApi` at-
tribute), 65
`print_node()` (`allensdk.core.swc.Compartment`
method), 222
`print_out()` (`allensdk.ephys.feature_extractor.EphysFeature`
method), 239
`print_summary()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.receptive_field`),
142
`ProbeSynchronizer` (class in `al-`
`lensdk.brain_observatory.ecephys.align_timestamps.probe_synchronizer`),
113
`process()` (`allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor`
method), 227
`process_eye_tracking_data()` (in module `al-`
`lensdk.brain_observatory.behavior.eye_tracking_processing`),
95
`process_inputs()` (in module `al-`
`lensdk.internal.model.biophysical.passive_fitting.neuron_packing`),
275
`process_instance()` (al-
`lensdk.ephys.feature_extractor.EphysFeatureExtractor`
method), 239
`process_new_spike_feature()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`process_new_sweep_feature()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`process_segmentation()` (al-
`lensdk.mouse_connectivity.grid.subimage.base_subimage.Segment`
method), 331
`process_segmentation()` (al-
`lensdk.mouse_connectivity.grid.subimage.classic_subimage.Classic`
method), 332
`process_segmentation()` (al-
`lensdk.mouse_connectivity.grid.subimage.count_subimage.Count`
method), 333
`process_spikes()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`process_spikes()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor`
method), 230
`PROJECTION_DENSITY` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
attribute), 58
`project_to_plane()` (in module `al-`
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
method), 256
`PROJECTION_DENSITY` (al-
`lensdk.brain_observatory.ecephys.file_io.stim_file.CamStim`
attribute), 53
`projection_density` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
attribute), 301
`PROJECTION_DENSITY_KEY` (al-
`lensdk.core.mouse_connectivity_cache.MouseConnectivityCache`
attribute), 204
`PROJECTION_ENERGY` (al-
`lensdk.api.queries.grid_data_api.GridDataApi`
attribute), 53
`projection_energy` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
attribute), 301
`propagate()` (`allensdk.internal.mouse_connectivity.interval_unionize.tis`
method), 301
`propagate_record()` (`allensdk.internal.mouse_connectivity.interval_unionize.ur`
method), 302
`propagate_record()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_un`
class method), 300
`propagate_record()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_u`
class method), 302
`propagate_to_bilateral()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_un`
class method), 300

`propagate_unionizes()` (al- `RaisingSchema.Meta` (class in al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize`), 101
 class method), 300
`provide()` (`allensdk.core.authentication.CredentialProvider`), 188
 method), 188
`provide()` (`allensdk.core.authentication.EnvCredentialProvider`), 227
 method), 189
`psycpg2_select()` (in module al- `lensdk.internal.model.glif.glif_optimizer.GlifOptimizer`
`lensdk.internal.api`), 251 method), 282
`pupil_position_in_mouse_eye_coordinates()` (`allensdk.internal.brain_observatory.eye_calibration.EyeCalibration`), 255
 method), 255
`pupil_position_on_monitor_in_cm()` (al- `ransac_fit()` (`allensdk.internal.brain_observatory.fit_ellipse.FitEllipse`), 257
`lensdk.internal.brain_observatory.eye_calibration.EyeCalibration`), 255
 method), 255
`pupil_position_on_monitor_in_degrees()` (`allensdk.internal.brain_observatory.eye_calibration.EyeCalibration`), 255
 method), 255
`push_summary()` (al- `lensdk.mouse_connectivity.grid.utilities.image_utilities`), 335
`lensdk.ephys.feature_extractor.EphysFeatureExtractor` method), 239
`pval()` (`allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis`), 170
 attribute), 170
`pvalue_to_NLL()` (in module al- `raw_running_data_df` (al-
`lensdk.brain_observatory.receptive_field_analysis.chisquare_df`), 139
 attribute), 89
`pvalue_to_NLL()` (in module al- `raw_running_speed` (al-
`lensdk.brain_observatory.receptive_field_analysis.visualization`), 144
 attribute), 89
`PycfgDescriptionParser` (class in al- `read()` (`allensdk.config.model.description_parser.DescriptionParser`
`lensdk.config.model.formats.pycfg_description_parser`), method), 184
 183
`read()` (`allensdk.config.model.formats.hdf5_util.Hdf5Util` method), 182
`read()` (`allensdk.config.model.formats.json_description_parser.JsonDescriptionParser` method), 182
`read()` (`allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser` method), 183
`read()` (in module `allensdk.core.json_utilities`), 202
`read()` (in module al- `lensdk.internal.mouse_connectivity.interval_unionize.data_utilities`), 299
`read_cell_index_receptive_field_analysis()` (`allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise` static method), 156
`read_data()` (`allensdk.api.api.Api` method), 74
`read_eye_dlc_tracking_ellipses()` (in module `allensdk.brain_observatory.nwb`), 135
`read_eye_gaze_mappings()` (in module al- `lensdk.brain_observatory.nwb`), 135
`read_file()` (`allensdk.core.json_utilities.JsonComments` class method), 202
`read_h5_group()` (in module al- `lensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis`), 142
 142

`read_h5_group()` (in module `allensdk.brain_observatory.receptive_field_analysis.tools`), 142
`read_stimulus()` (in module `allensdk.model.biophysical.utils.Utils`), 319
`read_intensity_image()` (in module `allensdk.mouse_connectivity.grid.utilities.image_utilities`), 335
`read_stimulus_name_from_path()` (in module `allensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spike_times`), 319
`read_json()` (`allensdk.api.queries.biophysical_api.BiophysicalAPI`), 24
`read_strided()` (in module `allensdk.model.biophysical.module_reader.BiophysicalModuleReader`), 261
`read_json()` (`allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader`), 241
`read_json()` (`allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader`), 243
`read_strided()` (in module `allensdk.config.model.description_parser.DescriptionParser`), 184
`read_json_string()` (in module `allensdk.config.model.formats.json_description_parser.JsonDescriptionParser`), 182
`read_json_string()` (in module `allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser`), 183
`read_lims_file()` (in module `allensdk.core.json_utilities.JsonComments`), 202
`read_lims_file()` (in module `allensdk.core.swc`), 226
`read_lims_message()` (in module `allensdk.core.json_utilities`), 202
`read_lims_message()` (in module `allensdk.core.swc`), 226
`read_url_get()` (in module `allensdk.core.json_utilities`), 202
`read_url_post()` (in module `allensdk.core.json_utilities`), 203
`read_output()` (`allensdk.internal.model.biophysical.passive_fitting.output_utils`), 300
`receptive_field` (in module `allensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`), 156
`RECONSTRUCTION_KEY` (in module `allensdk.core.cell_types_cache.CellTypesCache`), 199
`read_ndarray_with_sitk()` (in module `allensdk.core.sitk_utilities`), 219
`record_cb()` (`allensdk.internal.mouse_connectivity.interval_unionize.interval_unionize_utils`), 300
`record_cb()` (`allensdk.internal.mouse_connectivity.interval_unionize.interval_unionize_utils`), 302
`record_values()` (in module `allensdk.internal.model.biophysical.deap_utils.Utils`), 276
`read_receptive_field_from_h5()` (in module `allensdk.brain_observatory.receptive_field_analysis.receptive_field_analysis`), 142
`read_segmentation_image()` (in module `allensdk.mouse_connectivity.grid.subimage.base_subimage_segmentation_subimage`), 331
`read_segmentation_image()` (in module `allensdk.mouse_connectivity.grid.utilities.image_utilities`), 335
`REFERENCE_SPACE_KEY` (in module `allensdk.core.reference_space_cache.ReferenceSpaceCache`), 214
`REFERENCE_SPACE` (in module `allensdk.core.reference_space`), 211

ReferenceSpaceApi (class in al-
 lensdk.api.queries.reference_space_api),
 63 required_segmentations (al-
 lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage attribute), 332

ReferenceSpaceCache (class in al-
 lensdk.core.reference_space_cache), 214 required_segmentations (al-
 lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage attribute), 333

refine_threshold_indexes() (in module al-
 lensdk.ephys.ephys_features), 238 resample_into_volume() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities), 335

remove_comments() (al-
 lensdk.core.json_utilities.JsonComments class method), 202 resample_volume() (al-
 lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder method), 337

remove_keys() (allensdk.api.cache.Cache static
 method), 78 resample_volume() (in module al-
 lensdk.mouse_connectivity.grid.utilities.image_utilities), 335

remove_lfp_noise() (in module al-
 lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.Subsampling), 117 resave_swc() (in module al-
 lensdk.internal.morphology.validate_swc), 298

remove_lfp_offset() (in module al-
 lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling.Subsampling), 117 reset() (allensdk.model.glif.glif_neuron.GlifNeuron
 method), 321

remove_multiline_comments() (al-
 lensdk.core.json_utilities.JsonComments class method), 202 reset_AScurrent_none() (in module al-
 lensdk.model.glif.glif_neuron_methods), 326

remove_unassigned() (al-
 lensdk.core.reference_space.ReferenceSpace method), 213 reset_AScurrent_sum() (in module al-
 lensdk.model.glif.glif_neuron_methods), 326

rename_columns() (allensdk.api.cache.Cache static
 method), 78 reset_hex_pack() (in module al-
 lensdk.brain_observatory.circle_plots), 150

renames() (allensdk.core.structure_tree.StructureTree
 static method), 222 reset_long_squares_start() (in module al-
 lensdk.ephys.ephys_extractor), 231

Report (class in al-
 lensdk.internal.model.biophysical.make_deap_fit_json), 277 reset_threshold_inf() (in module al-
 lensdk.model.glif.glif_neuron_methods), 326

ReporterStatus (class in al-
 lensdk.core.cell_types_cache), 200 reset_threshold_three_components() (in module al-
 lensdk.model.glif.glif_neuron_methods), 326

required_intensities (al-
 lensdk.mouse_connectivity.grid.subimage.base_subimage.BaseSubImage attribute), 330 reset_voltage_zero() (in module al-
 lensdk.model.glif.glif_neuron_methods), 327

required_intensities (al-
 lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage attribute), 332 reshape_response_array() (al-
 lensdk.brain_observatory.drifting_gratings.DriftingGratings method), 154

required_polys (al-
 lensdk.mouse_connectivity.grid.subimage.base_subimage.BaseSubImage attribute), 330 reshape_response_array() (al-
 lensdk.brain_observatory.natural_scenes.NaturalScenes method), 158

required_polys (al-
 lensdk.mouse_connectivity.grid.subimage.cav_subimage.CavSubImage attribute), 331 reshape_response_array() (al-
 lensdk.brain_observatory.static_gratings.StaticGratings method), 169

required_polys (al-
 lensdk.mouse_connectivity.grid.subimage.classic_subimage.ClassicSubImage attribute), 332 resolve_initial_image() (in module al-
 lensdk.brain_observatory.behavior.trials_processing), 108

required_polys (al-
 lensdk.mouse_connectivity.grid.subimage.count_subimage.CountSubImage attribute), 333

required_segmentations (al-
 lensdk.mouse_connectivity.grid.subimage.base_subimage.BaseSubImage attribute), 331

<code>resolve_paths()</code>	(<code>allensdk.config.manifest.Manifest</code>	<code>attribute</code>), 170
<code>method</code>), 186	<code>RoiMask</code>	(<code>class</code> in <code>al-</code>
<code>response</code> (<code>allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis</code>)	<code>attribute</code>), 162	
<code>rolling_window()</code>	(in <code>module</code> <code>al-</code>	
<code>response_bias()</code>	(in <code>module</code> <code>al-</code>	<code>lensdk.brain_observatory.demixer</code>), 151
<code>lensdk.brain_observatory.behavior.session_metrics</code>)	<code>rolling_window()</code>	(in <code>module</code> <code>al-</code>
102		<code>lensdk.internal.brain_observatory.demixer</code>),
<code>retinotopy_metric()</code>	(in <code>module</code> <code>al-</code>	254
<code>lensdk.internal.brain_observatory.annotated_region_on_stimulus</code>)	<code>attribute</code>), 225	
252	<code>rotate()</code>	(<code>allensdk.internal.mouse_connectivity.projection_thumbnail.vo</code>
<code>retrieve_file_from_storage()</code>	(<code>al-</code>	<code>method</code>), 305
<code>lensdk.internal.api.api_prerelease.ApiPrerelease</code>	<code>rotate()</code>	(in <code>module</code> <code>al-</code>
<code>method</code>), 244		<code>lensdk.brain_observatory.stimulus_info</code>),
<code>retrieve_file_over_http()</code>	(<code>al-</code>	174
<code>lensdk.api.api.Api</code> <code>method</code>), 74	<code>rotate_and_extract()</code>	(<code>al-</code>
<code>retrieve_parsed_json_over_http()</code>	(<code>al-</code>	<code>lensdk.internal.mouse_connectivity.projection_thumbnail.volume</code>
<code>lensdk.api.api.Api</code> <code>method</code>), 75	<code>method</code>), 305	
<code>retrieve_xml_over_http()</code>	(<code>allensdk.api.api.Api</code> <code>method</code>), 75	<code>rotate_ray()</code>
		(in <code>module</code> <code>al-</code>
<code>return_mask_cb()</code>	(<code>al-</code>	<code>lensdk.internal.brain_observatory.itracker_utils</code>),
<code>lensdk.core.reference_space.ReferenceSpace</code>	<code>rotate_vector()</code>	(in <code>module</code> <code>al-</code>
<code>static method</code>), 213		<code>lensdk.internal.brain_observatory.fit_ellipse</code>),
<code>reward_rate()</code>	(in <code>module</code> <code>al-</code>	257
<code>lensdk.brain_observatory.behavior.trial_masks</code>),	<code>row_from_cell_id()</code>	(<code>al-</code>
104		<code>lensdk.brain_observatory.stimulus_analysis.StimulusAnalysis</code>
<code>rewards</code> (<code>allensdk.brain_observatory.behavior.behavior_data_session_behavior_data_session</code>)	<code>attribute</code>), 89	<code>run()</code>
<code>rings_in_hex_pack()</code>	(in <code>module</code> <code>al-</code>	<code>method</code>), 280
<code>lensdk.brain_observatory.circle_plots</code>), 150	<code>run()</code>	(<code>allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher.S</code>
<code>rma_templates</code>	(<code>al-</code>	<code>method</code>), 305
<code>lensdk.api.queries.biophysical_api.BiophysicalApi</code>	<code>run()</code>	(<code>allensdk.model.glif.glif_neuron.GlifNeuron</code>
<code>attribute</code>), 46		<code>method</code>), 322
<code>rma_templates</code>	(<code>al-</code>	<code>run()</code>
<code>lensdk.api.queries.brain_observatory_api.BrainObservatoryApi</code>	<code>run()</code>	(in <code>module</code> <code>al-</code>
<code>attribute</code>), 49		<code>lensdk.internal.mouse_connectivity.interval_unionize.run_tissuec</code>
<code>rma_templates</code> (<code>allensdk.api.queries.glif_api.GlifApi</code>	<code>run()</code>	(in <code>module</code> <code>al-</code>
<code>attribute</code>), 52		<code>lensdk.internal.mouse_connectivity.projection_thumbnail.generat</code>
<code>rma_templates</code>	(<code>al-</code>	303
<code>lensdk.api.queries.image_download_api.ImageDownloadApi</code>	<code>run()</code>	(in <code>module</code> <code>allensdk.model.biophysical.runner</code>),
<code>attribute</code>), 57		317
<code>rma_templates</code>	(<code>al-</code>	<code>run_base_model()</code>
<code>lensdk.api.queries.ontologies_api.OntologiesApi</code>	<code>run_base_model()</code>	(<code>al-</code>
<code>attribute</code>), 63		<code>lensdk.internal.model.glif.glif_experiment.GlifExperiment</code>
<code>rma_templates</code>	(<code>al-</code>	<code>method</code>), 281
<code>lensdk.internal.api.queries.biophysical_module_api.BiophysicalModuleApi</code>	<code>run_many()</code>	(<code>allensdk.internal.model.glif.glif_optimizer.GlifOptimizer</code>
<code>attribute</code>), 240	<code>run_module()</code>	(in <code>module</code> <code>al-</code>
<code>RmaApi</code> (<code>class</code> in <code>allensdk.api.queries.rma_api</code>), 65		<code>lensdk.internal.core.lims_pipeline_module</code>),
<code>RmaPager</code> (<code>class</code> in <code>allensdk.api.queries.rma_pager</code>),	<code>run_module()</code>	265
70		(in <code>module</code> <code>al-</code>
<code>RmaTemplate</code>	(<code>class</code> in <code>al-</code>	<code>lensdk.model.biophys_sim.bps_command</code>),
<code>lensdk.api.queries.rma_template</code>), 70		316
<code>robust_std()</code>	(in <code>module</code> <code>al-</code>	<code>run_once()</code>
<code>lensdk.brain_observatory.dff</code>), 153		(<code>allensdk.internal.model.glif.glif_optimizer.GlifOptimizer</code>
<code>roi_id</code> (<code>allensdk.brain_observatory.stimulus_analysis.StimulusAnalysis</code>)	<code>bound()</code>	(<code>al-</code>

[lensdk.internal.model.glif.glif_optimizer.GlifOptimizer](#) (al-
 method), 282
[lensdk.config.manifest.Manifest](#) class method),
 run_ophys_time_sync() (in module al- 186
[lensdk.internal.pipeline_modules.run_ophys_time_sync](#)),
 314 (allensdk.config.manifest.Manifest
 class method), 186
 run_passive_fit() (in module al- safe_system_path() (in module al-
[lensdk.internal.model.biophysical.run_passive_fit](#)), [lensdk.internal.core.lims_utilities](#)), 266
 278
 SAG_TARGET ([allensdk.ephys.ephys_extractor.EphysCellFeatureExtractor](#)
 run_postprocessing() (in module al- attribute), 227
[lensdk.brain_observatory.receptive_field_analysis.postprocessing](#)), (in module al-
 141 [lensdk.internal.model.GLM](#)), 290
 run_session_analysis() (in module al- sample_freq ([allensdk.brain_observatory.sync_dataset.Dataset](#)
[lensdk.brain_observatory.session_analysis](#)), attribute), 178
 166
 run_subimage() (in module al- sample_frequency (al-
[lensdk.mouse_connectivity.grid.subimage](#)), [lensdk.brain_observatory.ecephys.file_io.ecephys_sync_dataset.E](#)
 333 attribute), 116
 run_subimage() (in module al- sampling_rate_scale (al-
[lensdk.mouse_connectivity.grid.subimage.base_subimage](#)), [lensdk.brain_observatory.ecephys.align_timestamps.probe_sync](#)
 331 attribute), 113
 run_sync() (in module al- save() ([allensdk.core.swc.Morphology](#) method), 225
[lensdk.model.biophysical.runner](#)), 317 save() ([allensdk.internal.morphology.morphology.Morphology](#)
 method), 294
 run_until_biological_spike() (al- save_analysis_arrays() (al-
[lensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron](#)
 method), 283 [allensdk.brain_observatory_nwb_data_set.BrainObservatoryN](#)
 method), 197
 run_with_biological_spikes() (al- save_analysis_dataframes() (al-
[lensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron](#)
 method), 283 [allensdk.brain_observatory_nwb_data_set.BrainObservatoryN](#)
 method), 197
 running_data_df (al- save_cell_index_receptive_field_analysis() (al-
[lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession](#)
 attribute), 89 [allensdk.brain_observatory.locally_sparse_noise.LocallySparseN](#)
 static method), 156
 running_speed (al- save_ephys_data() (al-
[lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession](#)
 attribute), 89 [allensdk.brain_observatory.cell_types_api.CellTypesApi](#)
 method), 50
 RunningSpeed (class in al- save_figure() (in module al-
[lensdk.brain_observatory.running_speed](#)), [lensdk.internal.ephys.plot_qc_figures](#)), 271
 164
 RunningSpeedPathsSchema (class in al- save_figure() (in module al-
[lensdk.brain_observatory.nwb.schemas](#)), [lensdk.internal.ephys.plot_qc_figures3](#)),
 132 274
 RunOptimize (class in al- save_nwb() (in module al-
[lensdk.internal.model.biophysical.run_optimize](#)), [lensdk.model.biophysical.runner](#)), 318
 277
 RunSimulate (class in al- save_ophys_experiment_analysis_data() (al-
[lensdk.model.biophysical.run_simulate](#)), [allensdk.api.queries.brain_observatory_api.BrainObservatoryAp](#)
 317 method), 49
 RunSimulateLims (class in al- save_ophys_experiment_data() (al-
[lensdk.internal.model.biophysical.run_simulate_lims](#)), [lensdk.api.queries.brain_observatory_api.BrainObservatoryApi](#)
 278 method), 49
 save_ophys_experiment_event_data() (al-
[lensdk.api.queries.brain_observatory_api.BrainObservatoryApi](#)
 method), 49
 save_ophys_experiment_eye_gaze_data() (al-
[allensdk.api.queries.brain_observatory_api.BrainObservatoryAp](#)
 safe_factory() (al- [lensdk.internal.mouse_connectivity.projection_thumbnail.volume_projector.VolumeProjector](#)
 class method), 305 save_qc_figures() (in module al-

`lensdk.internal.ephys.core_feature_extract()`, 269
`save_reconstruction()` (al-
`lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 50
`save_reconstruction_markers()` (al-
`lensdk.api.queries.cell_types_api.CellTypesApi`
`method`), 51
`save_session_a()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 165
`save_session_b()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 165
`save_session_c()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 165
`save_session_c2()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 165
`save_voltage()` (al-
`lensdk.core.dat_utilities.DatUtilities`
`method`), 201
`schema()` (al-
`lensdk.api.queries.connected_services.ConnectedServices`
`class method`), 51
`score_feature_set()` (al-
`lensdk.ephys.feature_extractor.EphysFeatureExtractor`
`method`), 239
`search()` (al-
`lensdk.brain_observatory.stimulus_info.BinarySearchTree`
`method`), 171
`search()` (al-
`lensdk.brain_observatory.stimulus_info.StimulusSearchTree`
`method`), 172
`section_image_query()` (al-
`lensdk.api.queries.image_download_api.ImageDownloadApi`
`method`), 57
`SegmentationSubImage` (class in al-
`lensdk.mouse_connectivity.grid.subimage.base_subimage`), 330
`select()` (al-
`lensdk.internal.api.PostgresQueryMixin`
`method`), 250
`select()` (in module al-
`lensdk.internal.core.lims_utilities`), 266
`select_channels()` (in module al-
`lensdk.brain_observatory.ecephys.lfp_subsampling.subsampling`
`method`), 118
`select_one()` (al-
`lensdk.internal.api.PostgresQueryMixin`
`method`), 251
`serialize()` (al-
`lensdk.brain_observatory.behavior.image_api.ImageApi`
`static method`), 96
`SERVICE` (al-
`lensdk.api.queries.rma_api.RmaApi`
`attribute`), 65
`service_query()` (al-
`lensdk.api.queries.rma_api.RmaApi`
`method`), 69
`service_stage()` (al-
`lensdk.api.queries.rma_api.RmaApi`
`method`), 69
`session_a()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 166
`session_b()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 166
`session_c()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 166
`session_c2()` (al-
`lensdk.brain_observatory.session_analysis.SessionAnalysis`
`method`), 166
`session_na()` (al-
`lensdk.brain_observatory.ecephys.ecephys_session_api.EcephysSessionApi`
`attribute`), 115
`SessionAnalysis` (class in al-
`lensdk.brain_observatory.session_analysis`), 164
`sessions_with_stimulus()` (in module al-
`lensdk.brain_observatory.stimulus_info`), 174
`set_actual_parameters()` (al-
`lensdk.internal.model.biophysical.deap_utils.Utils`
`method`), 276
`set_api_urls()` (al-
`lensdk.api.api.Api`
`method`), 75
`set_apical_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 295
`set_axon_color()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 295
`set_base_search_tree()` (al-
`lensdk.internal.morphology.morphvis.MorphologyColors`
`method`), 295
`set_coarse_grid_parameters()` (al-
`lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeriesGridder`
`method`), 337
`set_credential_provider()` (in module al-
`lensdk.core.authentication`), 189
`set_default_working_directory()` (al-
`lensdk.api.api.Api`
`method`), 75
`set_dims()` (al-
`lensdk.brain_observatory.circle_plots.CoronaPlotter`
`method`), 149
`set_dims()` (al-
`lensdk.brain_observatory.circle_plots.FanPlotter`
`method`), 149
`set_F()` (al-
`lensdk.brain_observatory.r_neuropil.NeuropilSubtract`
`method`), 160
`set_iclamp_params()` (al-
`lensdk.internal.model.biophysical.deap_utils.Utils`
`method`), 276
`set_image_spacing()` (in module al-
`lensdk.mouse_connectivity.grid.utilities.image_utilities`), 335
`set_max_voxel()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_utilities`
`method`), 301
`set_neuron_parameters()` (al-

`lensdk.internal.model.glif.glif_experiment.GlifExperiment` method), 281
`set_normalized_parameters()` (alias `lensdk.internal.model.biophysical.deap_utils.Utils` method), 276
`set_omitted_stop_time()` (in module `lensdk.brain_observatory.nwb.nwb_utils`), 131
`set_params()` (`allensdk.brain_observatory.session_api_utils.ParamsMixin` method), 167
`set_sitk_image_information()` (in module `lensdk.core.sitk_utilities`), 219
`set_soma_color()` (alias `lensdk.internal.morphology.morphvis.MorphologyColors` method), 295
`set_spatial_unit()` (alias `lensdk.brain_observatory.stimulus_info.Monitor` method), 172
`set_spike_times()` (alias `lensdk.core.nwb_data_set.NwbDataSet` method), 209
`set_stimulus_amplitude_calculator()` (alias `lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor` method), 229
`set_sweep()` (`allensdk.core.nwb_data_set.NwbDataSet` method), 209
`set_version()` (alias `lensdk.config.manifest_builder.ManifestBuilder` method), 187
`set_workflow_state()` (alias `lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader` method), 241
`setup_iclamp()` (alias `lensdk.model.biophysical.utils.Utils` method), 319
`setup_images()` (alias `lensdk.mouse_connectivity.grid.subimage.base_subimage.IntensitySubImage` method), 330
`setup_images()` (alias `lensdk.mouse_connectivity.grid.subimage.base_subimage.PolygonsSubImage` method), 330
`setup_images()` (alias `lensdk.mouse_connectivity.grid.subimage.base_subimage.SegmentationSubImage` method), 331
`setup_images()` (alias `lensdk.mouse_connectivity.grid.subimage.base_subimage.SubImage` method), 331
`setup_interval_map()` (alias `lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize.IntervalUnionizer` method), 300
`setup_model()` (alias `lensdk.internal.model.biophysical.make_deap_fit_report.make_report_neuron()` method), 277
`setup_subimages()` (alias `lensdk.mouse_connectivity.grid.image_series_gridder.ImageSeries` method), 337
`setup_table_for_epochs()` (in module `lensdk.brain_observatory.nwb`), 136
`setup_table_for_invalid_times()` (in module `allensdk.brain_observatory.nwb`), 136
`sfvals` (`allensdk.brain_observatory.static_gratings.StaticGratings` attribute), 169
`split_params_mixer_features()` (alias `lensdk.ephys.ephys_extractor.EphysCellFeatureExtractor` method), 227
`short_string()` (alias `lensdk.internal.morphology.node.Node` method), 298
`show_angle_labels()` (alias `lensdk.brain_observatory.circle_plots.FanPlotter` method), 149
`show_arrow()` (`allensdk.brain_observatory.circle_plots.CoronaPlotter` method), 149
`show_arrow()` (`allensdk.brain_observatory.circle_plots.TrackPlotter` method), 149
`show_axes()` (`allensdk.brain_observatory.circle_plots.FanPlotter` method), 149
`show_circle()` (alias `lensdk.brain_observatory.circle_plots.CoronaPlotter` method), 149
`show_group_labels()` (alias `lensdk.brain_observatory.circle_plots.FanPlotter` method), 149
`show_image()` (`allensdk.brain_observatory.stimulus_info.Monitor` method), 172
`show_r_labels()` (alias `lensdk.brain_observatory.circle_plots.FanPlotter` method), 149
`simple_rotation()` (in module `lensdk.internal.mouse_connectivity.projection_thumbnail.generator`), 203
`SimpleTree` (class in `allensdk.core.simple_tree`), 216
`SimpleTree` (class in `allensdk.core.simpletree`), 267
`simplify_cells_api()` (alias `lensdk.api.queries.cell_types_api.CellTypesApi` method), 51
`simplify_experiment_containers()` (alias `lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 49
`simplify_ophys_experiments()` (alias `lensdk.api.queries.brain_observatory_api.BrainObservatoryApi` method), 49
`simulate()` (`allensdk.model.biophysical.run_simulate.RunSimulate` method), 317
`simulate_neuron()` (in module `lensdk.model.glif.simulate_neuron`), 328
`simulate_sweep()` (in module `al-`

`lensdk.model.glif.simulate_neuron`), 328
`simulate_sweep_from_file()` (in module `al-`
`lensdk.model.glif.simulate_neuron`), 328
`sitk_get_center()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
305
`sitk_get_diagonal_length()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
305
`sitk_get_image_parameters()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
305
`sitk_get_size_parity()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
305
`sitk_np_convert()` (in module `al-`
`lensdk.mouse_connectivity.grid.utilities.image_utilities`),
335
`sitk_paste_into_center()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
305
`sitk_safe_ln()` (in module `al-`
`lensdk.internal.mouse_connectivity.projection_thumbnail.volume_utilities`),
304
`size()` (`allensdk.internal.brain_observatory.mask_set.MaskSet`
method), 260
`slice_arrays()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize`
method), 302
`smooth()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.utilities`),
143
`smooth_STA()` (in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.utilities`),
139
`sobel_grad()` (in module `al-`
`lensdk.internal.brain_observatory.itracker_utils`),
259
`SOMA` (`allensdk.core.swc.Morphology` attribute), 223
`soma` (`allensdk.core.swc.Morphology` attribute), 225
`SOMA` (`allensdk.internal.morphology.morphology.Morphology`
attribute), 291
`soma_root()` (`allensdk.internal.morphology.morphology.Morphology`
method), 294
`sort_data_arrays()` (al-
`lensdk.internal.mouse_connectivity.interval_unionize.interval_unionize`
method), 300
`sort_trials()` (al-
`lensdk.brain_observatory.locally_sparse_noise.LocallySparseNoise`
method), 156
`spacing` (`allensdk.brain_observatory.behavior.image_api.Image`
attribute), 96
`SPACING` (`allensdk.core.swc.Marker` attribute), 223
`SPACING` (`allensdk.internal.core.swc.Marker` attribute),

268
`sparsify()` (`allensdk.core.swc.Morphology` method),
225
`sparsify()` (`allensdk.internal.morphology.morphology.Morphology`
method), 294
`spatial_frequency_to_pix_per_cycle()`
(`allensdk.brain_observatory.stimulus_info.Monitor`
method), 172
`spike_component_of_threshold_exact()`
(in module `al-`
`lensdk.model.glif.glif_neuron_methods`),
327
`spike_component_of_threshold_forward_euler()`
(`lensdk.model.glif.glif_neuron_methods`),
327
`spike_feature()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`spike_visualization_utilities.pages()` (al-
`lensdk.ephys.ephys_extractor.EphysSweepSetFeatureExtractor`
method), 230
`spike_visualization_utilities` (al-
`lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`SPIKE_TIMES` (`allensdk.core.nwb_data_set.NwbDataSet`
attribute), 208
`size_in_bytes()` (`allensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor`
method), 229
`spiral_trials()` (in module `al-`
`lensdk.brain_observatory.circle_plots`), 150
`spiral_trials_polar()` (in module `al-`
`lensdk.brain_observatory.circle_plots`), 150
`square_column()` (in module `al-`
`lensdk.brain_observatory.ecephys.stimulus_table.ephys_pre_spike`
method), 124
`standardize_movie_numbers()` (in module `al-`
`lensdk.brain_observatory.ecephys.stimulus_table.naming_utilities`
method), 125
`start()` (`allensdk.internal.model.biophysical.passive_fitting.output_grab`
method), 275
`START_ROW` (`allensdk.api.queries.rma_api.RmaApi` at-
tribute), 65
`start_specimen()` (al-
`lensdk.internal.model.biophysical.run_optimize.RunOptimize`
method), 275
`StaticGratings` (class in `al-`
`lensdk.brain_observatory.static_gratings`),
160
`stats()` (`allensdk.brain_observatory.sync_dataset.Dataset`
method), 178
`stim_table` (`allensdk.brain_observatory.stimulus_analysis.StimulusAna`
attribute), 171
`stim_table_to_categories()` (in module `al-`

[lensdk.brain_observatory.chisquare_categorical\(\)](#), [attribute](#)), 89
[148](#) [stimulus_times](#) (al-
[stim_timestamps](#) (al- [lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync](#)
[lensdk.internal.brain_observatory.time_sync.OphysTimeAlignment](#)), 313
[attribute](#)), 263 [stimulus_timestamps](#) (al-
[stimuli](#) ([allensdk.brain_observatory.ecephys.file_io.stim_file.CameraStimulusPickleSaveFile](#).[behavior.behavior_data_session.BehaviorDataSession](#)
[attribute](#)), 117 [attribute](#)), 90
[stimuli_in_session\(\)](#) (in module al- [StimulusAnalysis](#) (class in al-
[lensdk.brain_observatory.stimulus_info](#)), [lensdk.brain_observatory.stimulus_analysis](#)),
[174](#) [169](#)
[stimulus_alignment](#) (al- [StimulusSearch](#) (class in al-
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#).[lensdk.brain_observatory.stimulus_info](#)),
[attribute](#)), 313 [172](#)
[stimulus_amplitude\(\)](#) (al- [stitch\(\)](#) ([allensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)
[lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor](#) method), 305
[method](#)), 229 [Stitcher](#) (class in al-
[STIMULUS_CONTENT_TYPE](#) (al- [lensdk.internal.mouse_connectivity.tissuecyte_stitching.stitcher](#)),
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#)
[attribute](#)), 240 [stop\(\)](#) ([allensdk.internal.model.biophysical.passive_fitting.output_grabber](#)
[method](#)), 275
[STIMULUS_CONTENT_TYPE](#) (al- [method](#)), 275
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#) [OPTIMIZE_CONFIG_READERS_PRERELEASE_KEY](#) (al-
[attribute](#)), 243 [lensdk.internal.core.mouse_connectivity_cache_prerelease.MouseConnectivityCachePrerelease](#)
[stimulus_delay](#) (al- [attribute](#)), 266
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#).[lensdk.api.api](#), 75
[attribute](#)), 313
[stimulus_delta](#) (al- [stream_zip_directory_over_http\(\)](#) (in module al-
[lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncOutput](#).[lensdk.api.api](#)), 75
[attribute](#)), 313 [STRING](#) ([allensdk.api.queries.connected_services.ConnectedServices](#)
[stimulus_file_entries\(\)](#) (al- [attribute](#)), 51
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) (al-
[method](#)), 241 [lensdk.core.swc.Morphology](#) method), 226
[stimulus_file_entries\(\)](#) (al- [strip_all_other_types\(\)](#) (al-
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#).[lensdk.morphology.morphology.Morphology](#)
[method](#)), 243 [method](#)), 294
[STIMULUS_MAPPINGS_KEY](#) (al- [strip_type\(\)](#) ([allensdk.core.swc.Morphology](#)
[lensdk.core.brain_observatory_cache.BrainObservatoryCache](#) method), 226
[attribute](#)), 190 [strip_type\(\)](#) ([allensdk.internal.morphology.morphology.Morphology](#)
[stimulus_path\(\)](#) (al- [method](#)), 294
[lensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader](#) (al-
[method](#)), 241 [lensdk.core.ontology.Ontology](#) method),
[stimulus_path\(\)](#) (al- [211](#)
[lensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader](#).[lensdk.core.structure_tree.StructureTree](#)
[method](#)), 243 [method](#)), 222
[stimulus_presentations](#) (al- [method](#)), 222
[lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession](#) (al-
[attribute](#)), 89 [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[stimulus_search](#) (al- [attribute](#)), 214
[lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet](#) (al-
[attribute](#)), 197 [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[STIMULUS_TABLE_TYPES](#) (al- [attribute](#)), 214
[lensdk.core.brain_observatory_nwb_data_set.BrainObservatoryNwbDataSet](#) (al-
[attribute](#)), 194 [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)
[stimulus_templates](#) (al- [attribute](#)), 214
[lensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession](#) (al-
[attribute](#)), 89 [lensdk.core.reference_space_cache.ReferenceSpaceCache](#)

[lensdk.core.mouse_connectivity_cache.MouseConnectivityCache](#) (class in [allensdk.core.mouse_connectivity_cache](#)), 204
[STRUCTURES_KEY](#) (attribute), 204
[lensdk.core.reference_space_cache.ReferenceSpaceCache](#) (class in [lensdk.core.reference_space_cache](#)), 214
[StructureTree](#) (class in [lensdk.core.structure_tree](#)), 220
[stumpify_axon\(\)](#) (method), 226
[stumpify_axon\(\)](#) (method), 229
[SubImage](#) (class in [lensdk.mouse_connectivity.grid.subimage.base_subimage](#)), 331
[SubjectMetadataSchema](#) (class in [lensdk.brain_observatory.behavior.schemas](#)), 101
[subsampling_data\(\)](#) (method), 291
[subsampling_lfp\(\)](#) (method), 118
[subsampling_timestamps\(\)](#) (method), 119
[SUBTHRESH_MAX_AMP](#) (attribute), 227
[sum_pixel_intensity](#) (method), 301
[sum_pixels](#) (method), 301
[sum_projection_pixel_intensity](#) (method), 301
[sum_projection_pixels](#) (method), 301
[summarize\(\)](#) (method), 239
[SUMMARY_STRUCTURE_SET_ID](#) (attribute), 204
[summer_over\(\)](#) (method), 91
[SUPPORTED_PIPELINE_VERSION](#) (attribute), 194
[SupportsStr](#) (class in [allensdk.core.typing](#)), 227
[SvgApi](#) (class in [allensdk.api.queries.svg_api](#)), 70
[SWC_FILE_TYPE](#) (attribute), 90
[sweep_entries\(\)](#) (method), 241
[sweep_entries\(\)](#) (method), 243
[sweep_feature\(\)](#) (method), 229
[sweep_feature_keys\(\)](#) (method), 229
[sweep_features\(\)](#) (method), 230
[sweep_numbers\(\)](#) (method), 241
[sweep_numbers\(\)](#) (method), 243
[sweep_numbers_by_type\(\)](#) (method), 241
[sweep_response](#) (attribute), 157
[sweep_response](#) (attribute), 171
[sweep_response](#) (attribute), 156
[sweep_response](#) (attribute), 157
[sweep_response](#) (attribute), 158
[sweep_response](#) (attribute), 169
[sweep_response](#) (attribute), 230
[SynchronizationApi](#) (class in [allensdk.api.queries.synchronization_api](#)), 70
[synthesize_F\(\)](#) (method), 161
[task_parameters](#) (attribute), 90
[task_parameters](#) (attribute), 65
[task_parameters](#) (attribute), 161
[task_parameters](#) (attribute), 90

tau_m (*allensdk.model.glif.glif_neuron.GlifNeuron* attribute), 322
 temp_dir() (in module *allensdk.test_utilities.temp_dir*), 337
 TEMPLATE_KEY (*allensdk.core.reference_space_cache.ReferenceSpaceCache* attribute), 214
 template_projection() (in module *allensdk.internal.mouse_connectivity.projection_thumbnail_projection_functions*), 304
 template_query() (*allensdk.api.queries.rma_template.RmaTemplate* method), 70
 test() (*allensdk.brain_observatory.ecephys.ecephys_session_api.ecephys_session_api.EcephysSessionApi* method), 115
 test_fit() (in module *allensdk.internal.brain_observatory.fit_ellipse*), 257
 TestNode (class in *allensdk.internal.morphology.validate_swc*), 298
 tfvals (*allensdk.brain_observatory.drifting_gratings.DriftingGratings* attribute), 154
 Tile (class in *allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile*), 306
 timestamps (*allensdk.brain_observatory.running_speed.running_speed.LabelClassifier* attribute), 164
 timestamps (*allensdk.brain_observatory.stimulus_analysis.stimulus_analysis.TrainingMultiLabelClassifier* attribute), 171
 TimeSyncOutputs (class in *allensdk.internal.pipeline_modules.run_ophys_time_sync*), 313
 TimeSyncWriter (class in *allensdk.internal.pipeline_modules.run_ophys_time_sync*), 313
 TissuecyteBaseUnionize (class in *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_base_unionize*), 301
 TissuecyteInjectionUnionize (class in *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_injection_unionize*), 302
 TissuecyteProjectionUnionize (class in *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_projection_unionize*), 302
 TissuecyteUnionizer (class in *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_unionizer*), 302
 to_config_string() (*allensdk.config.app.application_config.ApplicationConfig* method), 181
 to_dict() (*allensdk.internal.model.glif.glif_optimizer.GlifOptimizer* method), 282
 to_dict() (*allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron* method), 284
 to_dict() (*allensdk.internal.morphology.morphology.Morphology* method), 307
 to_dict() (*allensdk.internal.morphology.node.Node* method), 298
 to_dict() (*allensdk.model.glif.glif_neuron.GlifNeuron* method), 322
 to_dict() (*allensdk.model.glif.glif_neuron_methods.GlifNeuronMethod* method), 323
 to_dict() (*allensdk.api.queries.rma_template.RmaTemplate* method), 70
 to_manifest() (*allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader* method), 340
 to_manifest() (*allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader* method), 244
 total_voxel_counts() (*allensdk.core.reference_space.ReferenceSpace* method), 213
 total_voxel_map (*allensdk.core.reference_space.ReferenceSpace* attribute), 213
 TreestitchingTile (class in *allensdk.brain_observatory.circle_plots*), 149
 TrainingMultiLabelClassifier (class in *allensdk.internal.brain_observatory.roi_filter_utils*), 244
 TrainingMultiLabelClassifier (class in *allensdk.internal.brain_observatory.roi_filter_utils*), 262
 transform (*allensdk.mouse_connectivity.grid.image_series_gridded.ImageSeriesGridded* attribute), 337
 translate_image_and_fill() (in module *allensdk.brain_observatory.stimulus_info*), 174
 traverse_h5_file() (in module *allensdk.internal.mouse_connectivity.interval_unionize.tissuecyte_base_unionize*), 202
 tree() (*allensdk.core.swc.Morphology* method), 226
 tree() (*allensdk.internal.morphology.morphology.Morphology* method), 307
 TreeSearchApi (class in *allensdk.api.queries.tree_search_api*), 72
 tree_to_tissuecyte_unionize_record() (in module *allensdk.brain_observatory.behavior.trials_processing*), 108
 tree_to_tissuecyte_unionize_record() (in module *allensdk.brain_observatory.behavior.dprime*), 93
 types() (in module *allensdk.brain_observatory.behavior.trial_masks*), 104
 trials (*allensdk.brain_observatory.behavior.behavior_data_session.BehaviorDataSession* attribute), 90
 trim() (*allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile* method), 307

<code>trim_border_pulses()</code> (in module <code>allensdk.brain_observatory.ecephys.stimulus_sync</code>), 130	<code>update_default_cell_hoc()</code> (allensdk.model.biophysical.utils.Utils method), 319
<code>trim_discontiguous_times()</code> (in module <code>allensdk.brain_observatory.sync_utilities</code>), 145	<code>update_output_sweep_features()</code> (in module <code>allensdk.internal.ephys.core_feature_extract</code>), 169
<code>trim_self()</code> (allensdk.internal.mouse_connectivity.tissuecyte_stitching.tile.Tile method), 307	<code>update_well_known_file()</code> (allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader method), 241
<code>trimmed_stats()</code> (in module <code>allensdk.brain_observatory.ecephys.stimulus_sync</code>), 130	<code>update_well_known_file()</code> (allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader method), 244
<code>TRUE</code> (allensdk.api.queries.rma_api.RmaApi attribute), 65	<code>upsample_image_to_degrees()</code> (in module <code>allensdk.brain_observatory.receptive_field_analysis.utilities</code>), 143
<code>tuple_filters()</code> (allensdk.api.queries.rma_api.RmaApi method), 69	<code>user</code> (allensdk.core.authentication.DbCredentials attribute), 188
<code>two_out_of_three_aint_bad()</code> (in module <code>allensdk.brain_observatory.behavior.criteria</code>), 91	<code>Utils</code> (class in allensdk.model.biophysical.utils), 318
<code>TYPE</code> (allensdk.internal.model.glif.glif_optimizer_neuron.GlifOptimizerNeuron attribute), 282	
<code>TYPE</code> (allensdk.model.glif.glif_neuron.GlifNeuron attribute), 321	

V

U
 validate_epoch_durations() (in module al-
 union() (allensdk.internal.brain_observatory.mask_set.MaskSet
 method), 260
 union_size() (allensdk.internal.brain_observatory.mask_set.MaskSet
 method), 260
 Unionize (class in al-
 allensdk.internal.mouse_connectivity.interval_unionize_mask()), (in module al-
 allensdk.brain_observatory.roi_masks), 164
 unit (allensdk.brain_observatory.behavior.image_api.Image
 attribute), 96
 unknown (allensdk.brain_observatory.argschema_utilities.RaisingSchemaErrors
 attribute), 147
 unknown (allensdk.brain_observatory.behavior.schemas.RaisingSchemaErrors
 attribute), 101
 unpack() (allensdk.config.model.description.Description
 method), 184
 unpack_change_log() (in module al-
 allensdk.brain_observatory.behavior.stimulus_processing), class method), 216
 104
 unpack_manifest() (al-
 allensdk.config.model.description.Description
 method), 184
 unpack_structure_set_ancestors() (al-
 allensdk.api.queries.ontologies_api.OntologiesApi
 method), 63
 unpack_uint32() (in module al-
 allensdk.brain_observatory.sync_dataset),
 178
 update_data() (al-
 allensdk.config.model.description.Description
 method), 184
 validate_epoch_durations() (in module al-
 allensdk.brain_observatory.ecephys.stimulus_table.output_validation), 126
 validate_epoch_order() (in module al-
 allensdk.brain_observatory.ecephys.stimulus_table.output_validation), 126
 validate_max_spontaneous_epoch_duration() (in module al-
 allensdk.brain_observatory.ecephys.stimulus_table.output_validation), 126
 validate_max_trial_duration() (in module al-
 allensdk.brain_observatory.ecephys.stimulus_table.output_validation), 126
 validate_max_trial_duration() (al-
 allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSync
 method), 313
 validate_structure_id() (al-
 allensdk.core.reference_space_cache.ReferenceSpaceCache
 class method), 216
 validate_structure_ids() (al-
 allensdk.core.reference_space_cache.ReferenceSpaceCache
 class method), 216
 validate_structures() (al-
 allensdk.core.reference_space.ReferenceSpace
 method), 214
 validate_swc() (in module al-
 allensdk.internal.morphology.validate_swc),
 298
 validate_swc() (in module al-
 allensdk.morphology.validate_swc), 329
 validate_trial_condition_exclusivity() (in module al-
 allensdk.brain_observatory.behavior.trials_processing),

108
 validate_with_synthetic_F() (in module al-
 lensdk.brain_observatory.r_neuropil), 161
 value_map() (allensdk.core.simple_tree.SimpleTree
 method), 218
 values (allensdk.brain_observatory.running_speed.RunningSpeed
 attribute), 164
 verify_roi_lists_equal() (al-
 lensdk.brain_observatory.session_analysis.SessionAnalysis static method), 222
 method), 166
 VERSION (allensdk.config.manifest.Manifest attribute),
 185
 vin (allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile
 attribute), 117
 visual_degrees_to_pixels() (al-
 lensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor
 method), 171
 visual_degrees_to_pixels() (al-
 lensdk.brain_observatory.stimulus_info.Monitor
 method), 172
 voltage_component_of_threshold_exact() (in
 module al-
 lensdk.model.glif.glif_neuron_methods),
 327
 voltage_component_of_threshold_forward_euler() (in
 module al-
 lensdk.model.glif.glif_neuron_methods),
 328
 voltage_deflection() (al-
 lensdk.ephys.ephys_extractor.EphysSweepFeatureExtractor
 method), 230
 VolumeProjector (class in al-
 lensdk.internal.mouse_connectivity.projection_thumbnails
 module), 305
 VOXEL_RESOLUTION_100_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 63
 VOXEL_RESOLUTION_10_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 63
 VOXEL_RESOLUTION_25_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 63
 VOXEL_RESOLUTION_50_MICRONS (al-
 lensdk.api.queries.reference_space_api.ReferenceSpaceApi
 attribute), 63
 vsig (allensdk.brain_observatory.ecephys.file_io.stim_file.CamStimOnePickleStimFile
 attribute), 117

warp_image() (allensdk.brain_observatory.stimulus_info.BrainObservatoryMonitor
 method), 171
 warp_stimulus_coords() (in module al-
 lensdk.brain_observatory.stimulus_info),
 174
 warp_speed_ring() (in module al-
 lensdk.brain_observatory.circle_plots), 150
 whitelist() (allensdk.core.structure_tree.StructureTree
 static method), 222
 WhitespaceStrippedString (class in al-
 lensdk.test_utilities.custom_comparators),
 337
 write_one_pickle_stim_file() (in module al-
 lensdk.brain_observatory.behavior.criteria),
 91
 write() (allensdk.brain_observatory.stimulus_info.Monitor
 attribute), 172
 window_average() (in module al-
 lensdk.mouse_connectivity.grid.utilities.downsampling_utilities),
 334
 wrap() (allensdk.api.cache.Cache method), 78
 write() (allensdk.config.model.description_parser.DescriptionParser
 method), 184
 write() (allensdk.config.model.formats.hdf5_util.Hdf5Util
 method), 182
 write() (allensdk.config.model.formats.json_description_parser.JsonDescriptionParser
 method), 182
 write() (allensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser
 method), 183
 write() (allensdk.core.swc.Morphology method), 226
 write() (allensdk.internal.brain_observatory.frame_stream.FrameOutputStream
 method), 258
 write() (allensdk.internal.morphology.morphology.Morphology
 method), 295
 write() (allensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncModule
 method), 313
 write() (in module allensdk.core.json_utilities), 203
 write_file() (allensdk.internal.api.queries.biophysical_module_reader.BiophysicalModuleReader
 method), 241
 write_file() (allensdk.internal.api.queries.optimize_config_reader.OptimizeConfigReader
 method), 244
 write_snap_labels() (al-
 lensdk.core.reference_space.ReferenceSpace
 method), 214
 write_json_file() (al-
 lensdk.config.manifest_builder.ManifestBuilder
 method), 187
 write_ndarray_with_sitk() (in module al-
 lensdk.core.sitk_utilities), 219
 write_or_print_outputs() (in module al-
 lensdk.brain_observatory.argschema_utilities),

W

warp_coordinates (al-
 lensdk.brain_observatory.stimulus_info.ExperimentGeometry
 attribute), 171

[147](#)
`write_output()` (in module `al-`
`lensdk.internal.pipeline_modules.run_ophys_eye_calibration`),
[312](#)
`write_output_data()` (`al-`
`lensdk.internal.core.lims_pipeline_module.PipelineModule`
`method`), [265](#)
`write_output_h5()` (`al-`
`lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncWriter`
`method`), [314](#)
`write_output_json()` (`al-`
`lensdk.internal.pipeline_modules.run_ophys_time_sync.TimeSyncWriter`
`method`), [314](#)
`write_receptive_field_to_h5()`
(in module `al-`
`lensdk.brain_observatory.receptive_field_analysis.receptive_field`),
[142](#)
`write_string()` (`al-`
`lensdk.config.model.formats.json_description_parser.JsonDescriptionParser`
`method`), [182](#)
`write_string()` (`al-`
`lensdk.config.model.formats.pycfg_description_parser.PycfgDescriptionParser`
`method`), [183](#)
`write_string()` (in module `al-`
`lensdk.core.json_utilities`), [203](#)
`write_sweep_response()` (in module `al-`
`lensdk.model.glif.simulate_neuron`), [329](#)
`write_volume()` (in module `al-`
`lensdk.mouse_connectivity.grid.utilities.image_utilities`),
[335](#)

Y

`yesterday_was_good()` (in module `al-`
`lensdk.brain_observatory.behavior.criteria`),
[92](#)